



CUSTOMBIOTECH

Custom Biotech Catalog 14th Edition
Expertise You Can Trust



Your Roche Custom Biotech Customer Service _____

**Europe, Middle East, Africa,
and Latin America**

Roche Diagnostics Deutschland GmbH
Sandhofer Straße 116
68305 Mannheim, Germany
Phone +49 621 759 8580
Fax +49 621 759 8610
mannheim.custombiotech@roche.com

United States

Roche Diagnostics Corporation
Roche Applied Science
9115 Hague Road
P.O. Box 50414
Indianapolis, IN 46250-0414, USA
Phone +1 800 428 5433, ext. 14649 (toll-free)
Fax +1 317 521 4065
custombiotech.ussales@roche.com

Canada

Roche Diagnostics
201, Boulevard Armand-Frappier
HV 4A2 Laval, Québec, Canada
Phone +1 450 686 7050
Fax +1 450 686 7012
custombiotech.can@roche.com

Japan

Roche Diagnostics K.K.
Roche Custom Biotech
Applied Science
6-1, Shiba 2-chome
Minato-ku, Tokyo 105-0014, Japan
Phone +81 3 5443 5285
Fax +81 3 5443 7934
japan.custombiotech@roche.com

Asia Pacific

Roche Diagnostics Asia Pacific Pte. Ltd.
Regional Sales and Market Development
298 Tiong Bahru Road
16-01/06 Central Plaza
Singapore, 168730
Phone +65 6371 6638
Fax +65 6371 6601
apac.custombiotech@roche.com

custombiotech.roche.com



Dear Valued Customer,

We are pleased to present you with the new edition of the Roche Custom Biotech Product Catalog. This catalog officially renames our brand, Roche Custom Biotech, replacing Industrial Business. **Custom** underscores Roche's customer focus. **Biotech** is our most significant asset: We are a major global supplier of high quality raw material and reagent chemistries to diagnostics, life sciences and the pharma biotech industry.

Our new name better reflects our mission to provide high-performance products to the diagnostics and pharmaceuticals industry, customized to customer needs. Roche Custom Biotech combines excellence and added value to best serve your business.

Our key account managers are looking forward to supporting your manufacturing success with the finest raw materials and reagents available.

Roche Custom Biotech. Expertise You Can Trust.

This product catalog is divided into four sections:

- Clinical Chemistry: Enzymes, their cofactors and substrates, detergents and biocides
- Immunology: Marker enzymes and substrates, interference eliminating proteins, and biotin/streptavidin reagents
- Molecular Diagnostics: RT-PCR and PCR amplification and detection chemistries
- Pharma Biotech: Enzymes for dissociating cells from primary tissues, downstream proteases, and biocatalysts

Roche Custom Biotech has been manufacturing and supplying industrial chemistries for thirty years, in particular to the diagnostics industry.

With a portfolio of over one thousand products, we also custom manufacture to your needs. Roche is in a unique position to fulfil your manufacturing needs using state-of-the-art innovative and proprietary technologies.

We are committed to the mission of improving the quality of healthcare. Our expanding pharma product portfolio, including Recombinant Trypsin and Liberase Enzyme Blends, are manufactured under GMP regulations, to meet the exacting demands of pharmaceutical manufacturers. We look forward to hearing how we can customize to your individual needs.

Thank you very much for your faith in Roche Custom Biotech providing you with products of quality coupled with Roche Service. This catalog reflects our desire to be a partner in your success. Feel free to use this catalog as the beginning to find the products you need, and as the basis for our planning your custom manufacturing requirements.

Sincerely,

Stefan Schorling
Head of Global Marketing Roche Custom Biotech

Service and Support

Here at Roche Custom Biotech, we make it our constant goal to provide you with close business-to-business collaboration, strict confidentiality, and a strong scientific relationship throughout your products development cycle. This catalog is the first link in a worldwide chain of contact and support we provide for all our customers, whatever your field: Diagnostics, Life Sciences, Pharma Biotech, and beyond.

Need something different?

In addition to the products in this catalog, we can also specifically modify existing products and provide consultancy on the development of completely new items. We offer customization and contract manufacturing in nearly all fields of our portfolio (including OEM of our high-quality IVD reagents), at any scale or production stage. Whether you require raw materials, labeled components, or finished kits, our contract manufacturing service will deliver according to your specifications.

We're here for you.

For quick and easy ordering, your local Roche Custom Biotech consultant is part of a dedicated and responsive order management team.

If you have specific questions or technical support inquiries, rely on our highly trained Key Account Managers to be your dedicated contact at any stage of the process. See this catalog's back cover for contact details.

Discover our scope, service, and support.

This catalog is only the beginning of our commitment to maximizing the efficiency, quality, and profitability of your project.

Thank you for choosing Roche!

Visit us online at custombiotech.roche.com

- Explore Roche products and services.
- Find your local Roche representative.
- Download product literature.
- Obtain comprehensive information on products and applications.
- Access MSDS, pack inserts, certificates of analysis, publications, and more.
- Register to receive product updates and special offers.

The Roche Custom Biotech Team - Your Gateway to the Worldwide Roche Network

The Key Account Manager is the key interface between your company and Roche, accessing products and providing expertise to bring you the best.

The Logistics Specialist maintains and monitors your secure supply.

The Technical Specialist provides product technical support.

The Global Marketing Manager provides project and product management.

Roche Custom Biotech local headquarters


Indianapolis

Laval

Mannheim

Singapore

Tokyo

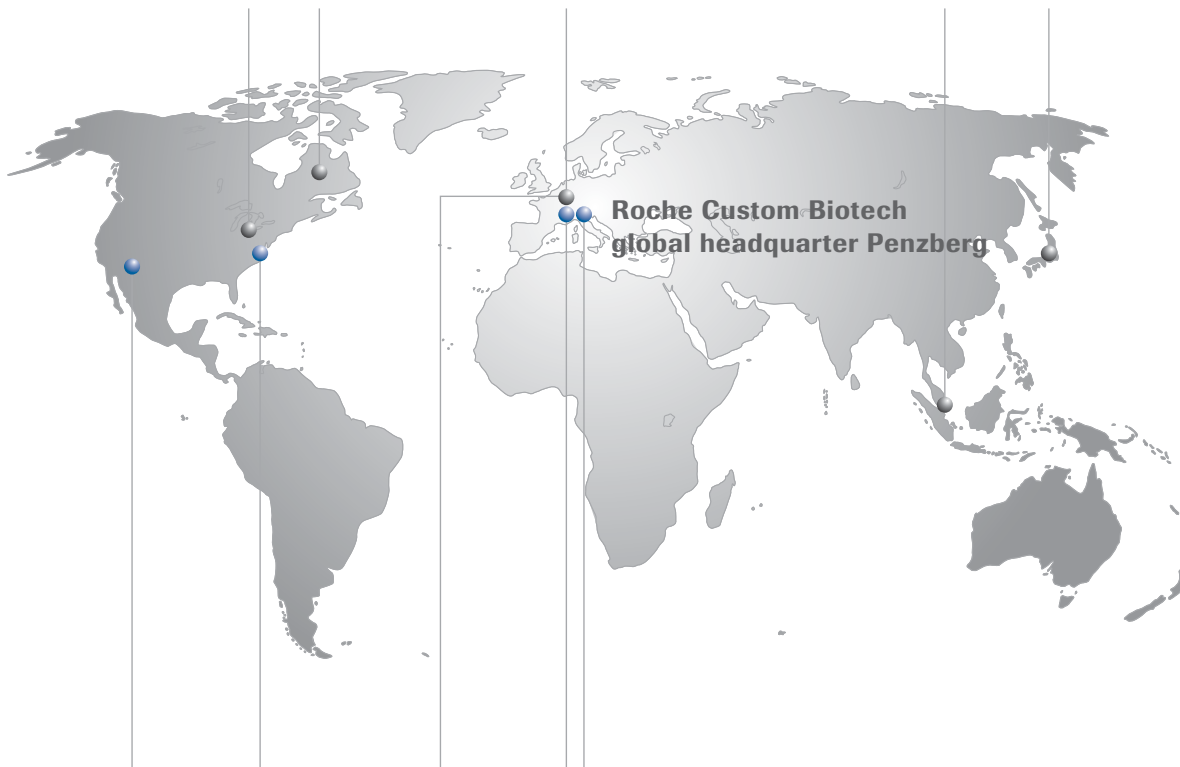

Tuscon Branchburg Mannheim Rotkreuz Penzberg

Roche production facilities

Building on over 30 years of industry experience, Roche Custom Biotech uses the powerful multidisciplinary skills found in Roche facilities across the world. Our research and production complex in

Penzberg, Germany, has become one of the largest biotechnology centers in Europe, and is the centerpiece of a constantly growing international corporate presence that provides solutions to over 500 major customers and partners.

MAB<DD>M-1.2.57 IgG	162
MAB<DD>M-2.1.16 IgG	164
Streptavidin R-Phycoerythrin LumiGrade Ultrasensitive Reagent	182
4-Aminophenyl Phosphate (pAPP), Disodium Salt	214
Proteinase K, recombinant, PCR Grade	223
AptaTaq DNA Polymerase, 5 U/μl	243
AptaTaq DNA Polymerase, 50 U/μl	244
AptaTaq DNA Polymerase LDx, 5 U/μl	245
AptaTaq DNA Polymerase LDx, 50 U/μl	247
AptaTaq Dexo DNA Polymerase, 5 U/μl	243
AptaTaq Dexo DNA Polymerase, 50 U/μl	244
HawkZ05 DNA Polymerase, 40 U/μl	258
HawkTaq DNA Polymerase, 5 U/μl	259
AptaTaq DNA Master	260
AptaTaq DNA Master without Mg ²⁺	261
AptaTaq DNA Master Optimization Kit	261
AptaTaq Genotyping Master	262
AptaTaq Genotyping Master (Rox)	263
EagleTaq Master Mix	264
EagleTaq Master Mix (Rox)	265
AllStart RNA Master	271
AllStart RNA Master (Rox)	272
HawkZ05 Fast One-Step RT-PCR Kit	273
LightCycler® 480 Multiwell Plate 384	296
LightCycler® 480 Multiwell Plate 96	297
T4 Gene 32 Protein, recombinant	300
COT Human DNA, CGH Grade	323
Bovine Serum Albumin, Molecular Biology Grade	334
Liberase MNP-S	340
Liberase MTF C/T, GMP Grade	341
Liberase T-Flex, Research Grade	343
Carboxypeptidase B, recombinant	348
Trypsin, recombinant	358
MycotoOL PCR Mycoplasma Detection Kit	359
CMP-N-Acetylneuraminic Acid	361
UDP-Galactose	363

1 Clinical Chemistry

Biochemicals for Clinical Chemistry	2
Albumin	2
Biocides	5
Buffers	16
Detergents	21
Ionic Detergents.....	21
Non-ionic Detergents.....	25
Zwitterionic Detergents.....	28
Protease Inhibitors	30
Additional Biochemicals	31
Cofactors/Nucleotides for Clinical Chemistry	35
Cofactors	35
Nucleotides	45
Enzymes for Clinical Chemistry	52
Substrates for Clinical Chemistry	139
Colorimetric Substrates.....	139
Non-Colorimetric Substrates.....	147

2 Immunology

Antibodies	160
Monoclonal Antibodies.....	160
Polyclonal Antibodies.....	176
Biotin/Streptavidin System	178
Streptavidin.....	178
Biotin Labels.....	179
Fluorescent Labels.....	181
Solid Phases.....	183
Dyes.....	190
Interference Eliminating Proteins (IEPs)	192
Specific Interference.....	192
Unspecific Interference.....	199
Marker Enzymes and Substrates	204
Enzymes.....	204
Substrates.....	212
Serums.....	216

3 Molecular Diagnostics

Sample Preparation	220
Chaotropic Salts.....	220
Enzymes.....	220
Amplification	260
DNA Polymerases	227
Expand System.....	227
T4 DNA Polymerase.....	232

Taq DNA Polymerase.....	234
Tth DNA Polymerase.....	238
DNA Polymerases, Hot Start	242
ActiQ ^Δ exo DNA Polymerase.....	242
AptaQ ^Δ DNA Polymerase.....	243
EagleQ ^Δ DNA Polymerase.....	252
FastStart DNA Polymerase.....	253
HawkZ05 DNA Polymerase.....	258
HawkQ ^Δ DNA Polymerase.....	259
DNA Master	260
Reverse Transcriptases	267
RNA Master	271
Nucleotides	274
deoxyNTPs.....	274
dideoxyNTPs.....	286
riboNTPs.....	290
Additional Products	296
Labeling and Detection	305
Conjugates	305
Enzymes	308
Labeled Nucleotides	311
Carrier and Competitor Nucleic Acids	323
DNA	323
RNA	326
Glycogen	328
Additional Reagents	329
Enzymes	329
Proteins	334

4 Pharma Biotech

Enzymes	340
Enzymes for cell isolation	340
Glycohydrolase	344
Proteases	348
Industrial Process Control	359
Mycoplasma Testing	359
Biochemicals	361
Activated Sugars	361
Cofactors	365
Additional Reagents	367

Appendix

Disclaimer	372
Trademarks	374
General Information	375
Product Index	376





1 Clinical Chemistry

Biochemicals for Clinical Chemistry	2
Albumin	2
Biocides	5
Buffers	16
Detergents	21
Ionic Detergents	21
Non-ionic Detergents	25
Zwitterionic Detergents	28
Protease Inhibitors	30
Additional Biochemicals	31
Cofactors/Nucleotides for Clinical Chemistry	35
Cofactors	35
Nucleotides	45
Enzymes for Clinical Chemistry	52
Substrates for Clinical Chemistry	139
Colorimetric Substrates	139
Non-Colorimetric Substrates	147

Bovine Serum Albumin (BSA), Fraction V lyophilizate

Serum albumin protein that has numerous biochemical applications.

Application

Use Bovine Serum Albumin (BSA) as a buffering agent, stabilizer, standard and for blending. Bovine Serum Albumin (BSA) is also a versatile tool against non-specific solid phase interference. As blocking reagent Bovine Serum Albumin (BSA) saturates unoccupied binding sites on the solid phase. Use Bovine Serum Albumin (BSA) typically at a concentration of 0.5 to 3% within the reagent buffer.

Benefits

- Keep your diagnostic reagent free from IgG contaminants.
- Take advantage of the strongly reduced concentration of heavy metal ions, sodium and potassium.
- Rely on the proven diagnostic quality of this product.

CAS: 9048-46-8

Properties

Molecular weight: 68 kD

Bovine Serum Albumin (BSA) contains no detectable IgG.

Bovine Serum Albumin (BSA) is controlled for low molecular weight contaminants.

Bovine Serum Albumin (BSA) consists primarily of monomeric albumin.

Specification

Appearance: Slightly yellow lyophilizate

Solubility: Clear, odourless solution in water

A₄₀₅ (against water): ≤0.200

Albumin (gel electrophoresis): ≥98%

Protein (from N according to elementary analysis) : ≥95%

pH value: 6.8-7.2

Water (K. Fischer): ≤5%

Heavy metals (as Pb): ≤0.003%

P_i: ≤0.005%

Chloride (chloride meter): ≤0.15%

Glucose (enzymatically): ≤0.05%

Gluceronol (enzymatically): ≤0.005%

L-Lactate (enzymatically): ≤0.1%

Na (flame photometrically): ≤0.8%

K (flame photometrically): ≤0.015%

Li (flame photometrically): ≤0.0005%

Ca: ≤0.05%

Mg: ≤0.005%

Fe: ≤0.002%

Bioburden: ≤100 CFU/g lyophilizate

Country of origin: New Zealand, USA

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Remarks:

Official veterinary certificate of health of the donor animals is available.

Official certificate of the deactivation of animal material including the method (acid treatment at pH 5 for 3 hours) is available.

Cat. No.

Pack Size

10 738 328 103 custom fill

Will be supplied as "Albumin, Fraction V from Bovine Serum". Unit of Measure is "kg".

For further processing only.

Bovine Serum Albumin (BSA), Fraction V fatty acids ≤ 0.2 mg/g, lyophilizate

Cat. No.	Pack Size
10 774 111 103	custom fill

Highly purified serum albumin protein that has numerous biochemical applications.

Application

Use Bovine Serum Albumin (BSA) as a buffering agent, stabilizer, standard and for blending. Bovine Serum Albumin (BSA) is also a versatile tool against non-specific solid phase interference. As blocking reagent Bovine Serum Albumin (BSA) saturates unoccupied binding sites on the solid phase. Use Bovine Serum Albumin (BSA) typically at a concentration of 0.5 to 3% within the reagent buffer.

Benefits

- Keep your diagnostic reagent free from IgG contaminants.
- Take advantage of the strongly reduced concentration of fatty acids
- Rely on the proven diagnostic quality of this product.

CAS: 9048-46-8

Properties

Molecular weight: 68 kD

Bovine Serum Albumin (BSA) contains no detectable IgG.

Bovine Serum Albumin (BSA) is controlled for low molecular weight contaminants.

Bovine Serum Albumin (BSA) consists primarily of monomeric albumin.

Specification

Appearance: Slightly yellow lyophilizate

Protein (from N, according to elementary analysis): $\geq 97\%$

Water (K. Fischer): $\leq 5\%$

Na (flame photometric): $\leq 0.5\%$

K (flame photometric): $\leq 0.01\%$

Fe (AAS): $\leq 0.001\%$

Cu (AAS): $\leq 0.002\%$

Fatty acids, total (GC): ≤ 0.2 mg/g

Triglycerides (enzymatically): Not detectable

Immunoglobulines (ELISA): Not detectable

Country of origin: USA

Stability: At +2 to +8°C within specification range for 24 months.

Remarks:

Official veterinary certificate of health of the donor animals is available.

Official certificate of the deactivation of animal material including the method (acid treatment at pH 5 for 3 hours) is available.

Bovine Serum Albumin (BSA), reduced sodium and potassium lyophilizate

Serum albumin protein for tests that require a strongly reduced concentration of sodium and potassium.

Application

Use Bovine Serum Albumin (BSA) as a buffering agent, stabilizer, standard and for blending. Bovine Serum Albumin (BSA) is also a versatile tool against non-specific solid phase interference. As blocking reagent Bovine Serum Albumin (BSA) saturates unoccupied binding sites on the solid phase. Use Bovine Serum Albumin (BSA) typically at a concentration of 0.5 to 3% in the reagent buffer.

Benefits

- Keep your diagnostic reagent free from IgG contaminants.
- Take advantage of the strongly reduced concentration of sodium, and potassium.
- Rely on the proven diagnostic quality of this product.

CAS: 9048-46-8

Properties

Molecular weight: 68 kD

Bovine Serum Albumin (BSA) contains no detectable IgG.

Bovine Serum Albumin (BSA) is controlled for low molecular weight contaminants.

Bovine Serum Albumin (BSA) consists primarily of monomeric albumin.

Specification

Appearance: Yellowish lyophilizate

A₄₀₅ (against water): ≤0.200

pH value: 4.3-5.3

Protein (Biuret): ≥80%

Water (K. Fischer): ≤5%

Na (AAS): ≤35 ppm

K (AAS) : ≤4ppm

Li (AS) : ≤50ppm

Ca (AAS): ≤500ppm

Fe (bathophenanthroline): ≤10 ppm

Cu (bathocuproin): ≤15ppm

Heavy metals (as Pb): ≤50ppm

P_i: ≤150ppm

Bioburden: ≤100 CFU/g lyophilizate

NH₄ (enzymatically): ≤10ppm

Glucose (enzymatically): ≤0.02%

Complex creator:

Recovery of Fe: 80-120%

Recovery of Cu: 80-120%

Electrophoresis: Chromatographically homogeneous

Country of origin: Germany

Stability: At +2 to + 8°C within specification range for 36 months.

Remarks:

Official veterinary certificate of health of the donor animals is available.

Official certificate of the deactivation of animal material including the method is available.

Cat. No.

Pack Size

11 297 368 103 custom fill

Unit of Measure is "kg".

For further processing only.

Biocides

Application

Biocides are used for preservation of reagents in diagnostic kits. The concentration recommended for each biocide is based on the highest MIC (minimal inhibition concentration) value to achieve the highest effectiveness. To decrease the biocide concentration and prevent resistance of microorganisms, use different biocides in

combination, according to their efficiency spectrum.

At the recommended concentrations, enzyme performance is not usually influenced. However, this should be verified empirically on a case-by-case basis.

Antimicrobial efficiency – MIC values

This table shows the biocide concentration at which the growth of microorganisms is completely prevented.

	MIC values (mg/ml)				recommended concentration (mg/ml)
	bacteria I	bacteria II	yeast	fungi	
5-Bromo-5-nitro-1.3-dioxane (BND)	0.2	0.08	0.08	0.1	0.2
2-Chloroacetamide (CAA)	>2.6	>2.6	1.3	>2.6	>2.6
2-Hydroxypyridine-N-oxide (Oxy-PYRION)	1.3	1.3	1.3	0.32	1.3
Imidazolinyurea (Germall Grade II)	1.3	>2.6	>2.6	2.6	>2.6
N-Methylisothiazolone (MIT)	0.12	0.2	>0.4	0.2	0.4
In comparison: Na-Azide	2.6	2.6	0.04	<0.02	

MIC = minimal inhibition concentration determined in dilution assays for 72 hours at +28°C

Microorganisms tested

bacteria I	bacteria II	yeast/fungi
<i>Bacillus subtilis</i>	<i>Aeromonas sp.</i>	<i>Candida albicans</i>
<i>Escherichia coli</i>	<i>Alcaligenes sp.</i>	<i>Rhodotorula rubra</i>
<i>Pseudomonas aeruginosa</i>	<i>Flavobacter sp.</i>	<i>Aspergillus oryzae</i>
<i>Staphylococcus aureus</i>	<i>Proteus vulgaris</i>	<i>Mucor racemosus</i>
<i>Streptococcus faecalis</i>	<i>Pseudomonas aeruginosa</i>	<i>Penicillium frequentans</i>
	<i>Pseudomonas fluorescens</i>	
	<i>Pseudomonas putida</i>	

Biocide/enzyme interactions

	AP activity (%)	biocide (mg/ml)	POD activity (%)	biocide (mg/ml)	β -Gal activity (%)	biocide (mg/ml)	Luciferase activity (%)	biocide (mg/ml)
5-Bromo-5-nitro-1.3-dioxane (BND)	100	10	120	10	100	1	66	0.1
2-Chloroacetamide (CAA)	70	3	81	3	90	3	86	3
2-Hydroxypyridine-N-oxide (Oxy-PYRION)	70	3	100	3	90	1.5	30	0.75
Imidazolinyurea (Germall. Grade II)	95	3	89	3	100	3	67	3
N-Methylisothiazolone (MIT)	100	10	100	10	80	1	97	0.5
In comparison:	0.5 mg/ml NaN_3 reduces activity of POD to 70% 1.0 mg/ml NaN_3 reduces activity of POD to 58% 10.0 mg/ml NaN_3 reduces activity of POD to 8%							

Benefits

- Contain no mercury, according to the strict current quality procedures
- Not carcinogenic or mutagenic (Ames-test negative)
- Broad efficiency spectra
- Can be combined for extension of efficiency spectra
- Low MIC values
- No or minor risk of resistance of microorganisms, especially when biocides are combined
- No or minor interferences in assays
- Readily soluble in water (in recommended concentrations)
- Toxically harmless (in recommended concentrations), considerably less toxic than sodium azide
- Highly stable, long lasting efficiency
- No coloring of reagent solution
- No odor
- Applicable in liquid reagents

Remark

In the European Union: For use in medical in vitro diagnostic products only.

2-Chloroacetamide (CAA) crystalline powder

Yeast specialized biocide used as preservative in diagnostic kits.

Application

Use 2-Chloroacetamide (CAA) as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Combine 2-Chloroacetamide (CAA) with non-ionic or ionic detergents.
- Use CAA in your preferred combination in liquid reagents to obtain the most stable reagent.

CAS: 79-07-2

Properties

Formula: C₂H₄ClNO

Molecular weight: 93.51 D

Antimicrobial effect: Highly effective against yeast, effective against bacteria and fungi.

Possible combination: Can be combined broadly; e.g. with N-methylisothiazolone (MIT), 2-hydroxypyridine-N-oxide (Oxy-PYRION).

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentration: ≥2.6 mg/ml

pH optimum: 4.0-8.0

Stability: Stable under most conditions.

Specification

Appearance: White crystalline powder

Solubility: Clear, colorless solution in water (c=10 mg/ml)

Identity (NIR): Corresponds to reference

2-Chloroacetamide (from N; based on anhydrous substance): 98-101%

Water (K. Fischer): ≤0.5%

N (elementary analysis): 14.65-15.05%

Sulfate ash: ≤0.1%

Stability: At +15 to +25°C within specification range for 24 months. Protect from light. Keep tightly sealed.

Cat. No. **Pack Size**

10 623 580 103 custom fill

Will be supplied as "2-Chloroacetamide". Unit of Measure is "g".
For further processing only.

2-Hydroxypyridine-N-oxide (Oxy-PYRION) crystalline powder

Broad range biocide for use in diagnostic kits

Application

Use Oxy-PYRION as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Combine Oxy-PYRION with non-ionic detergents or proteins.
- Use Oxy-PYRION in your preferred combination in liquid reagents to obtain the most stable system.

CAS: 13161-30-3

Properties

Formula: C₅H₅NO₂

Molecular weight: 111.10 D

Antimicrobial effect: Broad effect

Cat. No. **Pack Size**

11 637 126 103 custom fill

Will be supplied as "Hydroxypyridinoxid, pure". Unit of Measure is "kg".
For further processing only.

Possible combination: Can be combined with other biozides.
Interferences: Production of chelates with heavy metal ions.
Toxicity: Harmless in recommended working concentrations; not carcinogenic or mutagenic (Ames-test negative).
Recommended working concentration: 1.3 mg/ml
pH optimum: 5.0-8.0
Stability: Stable

Specification

Appearance: White to beige, crystalline powder
Solubility: Clear, colorless solution in water
A₄₀₅ (aqueous solution, against water): ≤0.010
Identity (NIR): Corresponds to reference
C (elementary analysis): 53.5-54.6%
2-Hydroxypyridine-n-oxide (from C): 99-101%
Water (K. Fischer): ≤0.5%
Heavy metals (as lead): ≤10 ppm
Fe (AAS): ≤10 ppm
Al (AAS): ≤10 ppm
Isopropanol (GC): ≤0.2%
Stability: At +15 to +25°C within specification range for 36 months. Protect from light.

2-Hydroxypyridine-N-oxide (Oxy-PYRION), reduced sodium crystalline powder

Broad range biocide for use in diagnostic kits

Application

Use 2-Hydroxypyridine-N-oxide, reduced sodium as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Combine 2-Hydroxypyridine-N-oxide, reduced sodium with non-ionic detergents or proteins.
- Use 2-Hydroxypyridine-N-oxide, reduced sodium in your preferred combination in liquid reagents to obtain the most stable system.
- Rely on the strongly reduced concentration of sodium.

CAS: 13161-30-3

Properties

Formula: C₅H₅NO₂
Molecular weight: 111.10 D
Antimicrobial effect: Broad effect
Possible combination: Can be combined with other biozides
Interferences: Production of chelates with heavy metal ions
Toxicity: Harmless in recommended working concentrations; not carcinogenic or mutagenic (Ames-test negative).
Recommended working concentration: 1.3 mg/ml
pH optimum: 5.0-8.0
Stability: Stable

Specification

Appearance: White to beige, crystalline powder
Solubility: Clear solution in water (c=1 mg/ml)
A₄₀₅ (c=1 mg/ml water, against water) : ≤0.010
Identity (NIR): Corresponds to reference

8 2-Hydroxypyridine-N-oxide (from C) : 99-101%

Cat. No. **Pack Size**

11 374 559 103 custom fill

Unit of Measure is "kg".
 For further processing only.

C (elementary analysis): 53.5-54.6%

Water (K. Fischer): ≤0.5%

Heavy metals (as Pb): ≤10 ppm

Fe (AAS): ≤10 ppm

Al (AAS): ≤10 ppm

Na (AAS): ≤70 ppm

2-Propanol: ≤0.2%

Stability: At +4°C to +8°C within specification range for 12 months. Protect from light.

5-Bromo-5-Nitro-1,3-Dioxane (BND), Grade I crystalline powder

Broad range biocide for use in diagnostic kits

Application

Use BND, Grade I as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Rely on the tested purity of BND Grade I.
- Combine BND with non-ionic detergents or proteins.
- Use BND in your preferred combination in liquid reagents to obtain the most stable system.
- Use BND in phosphate, HEPES or Tris buffer systems.

CAS: 30007-47-7

Properties

Formula: $C_4H_6BrNO_4$

Molecular weight: 212.0 D

Antimicrobial effect: Very effective against bacteria, yeast and fungi.

Possible combination: Can be combined broadly with other preservatives.

Effect: May oxidize thiol groups in essential enzyme systems.

Interferences: Interferes with cysteine. Reducing agents lower the preservative effectiveness.

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentrations:

0.2 mg/ml, if both the risk of contamination and concentration of protein are low (>0.2%);

0.4 mg/ml, if both the risk of contamination and concentration of protein and substrate are high.

pH optimum: 5.0-7.0

Stability: Not stable at pH < +5 and temperatures > +50°C; corrosive.

Specification

Appearance: White to slightly yellowish crystalline powder

Solubility: Clear, colourless solution in ethanol (c=10 mg/ml)

A₄₀₅ (c=10 mg/ml, ethanol) : ≤0.015

Identity (NIR): Corresponds to reference

5-Bromo-5-nitro-1.3-dioxane (GC) : ≥98.0 area%

Water (K. Fischer): ≤0.5%

Impurities (GC): ≤2.0 area%

Heavy metals (as Pb): ≤10 ppm

Oxidizing substances: ≤1%

Bromine, free: Negative

Stability: At +2 to +8°C within specification range for 18 months.

Cat. No.	Pack Size
11 354 361 103	custom fill

Will be supplied as "5-Bromo-5-nitro-1.3-dioxane, Pure". Unit of Measure is "g".

For further processing only.

5-Bromo-5-Nitro-1,3-Dioxane (BND), Grade II crystalline powder

Broad range biocide for use in diagnostic kits

Application

Use BND, Grade II as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Combine BND with non-ionic detergents or proteins.
- Use BND in your preferred combination in liquid reagents to obtain the most stable system.
- Use BND in phosphate, Hepes or Tris buffer systems.

CAS: 30007-47-7

Properties

Formula: C₄H₆BrNO₄

Molecular weight: 212.0 D

Antimicrobial effect: Very effective against bacteria, yeast and fungi.

Possible combination: Can be combined broadly with other preservatives.

Effect: May oxidize thiol groups in essential enzyme systems.

Interferences: Interferes with cysteine. Reducing agents lower the preservative effectiveness.

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentrations:

0.2 mg/ml, if both the risk of contamination and concentration of protein are low (>0.2%);

0.4 mg/ml, if both the risk of contamination and concentration of protein and substrate are high.

pH optimum: 5.0-7.0

Stability: Not stable at pH <+5 and temperatures >+50°C; corrosive.

Specification

Appearance: White crystalline powder probably slightly yellow

Solubility: Clear, colorless solution in ethanol (c=10 mg/ml)

A₄₀₅ (c=10 mg/ml, ethanol): ≤0.015

Identity (NIR): Corresponds to reference

HPTLC: Corresponds to reference

Water (K. Fischer): ≤0.5%

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No.

11 697 803 103

Pack Size

custom fill

Will be supplied as "Brom-nitro-dioxan,pure". Unit of Measure is "kg".

For further processing only.

Germall 115 crystalline powder

Bacteria specialized biocide used as preservative in diagnostic kits.

Application

Use Germall 115 as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Combine Germall 115 with non-ionic, ionic detergents or proteins.
- Use Germall 115 in your preferred combination in liquid reagents to obtain the most stable system.

CAS: 39236-46-9

Properties

Formula: $C_{11}H_{16}N_8O_8 \times H_2O$

Molecular weight: 406.31 D

Solubility: Easily soluble in water and glycerol; soluble in ethanol 40%; insoluble in ethanol 100%.

Antimicrobial effect: Only moderate effect against fungi and yeast; should be used in combination.

Possible combination: Can be combined with all biozides which do not separate formaldehyde, especially 5-bromo-5-nitro-1,3-dioxane (BND), 2-hydroxypyridine-N-oxide (HPO) or N-methylisothiazolone (MIT).

Effect: Belongs to the "formaldehyde-depot-substances". Inactivating effect due to reactions with $-COOH$, $-NH_2$, $-OH$ and $-SH$ of all proteins.

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentration: 2.6 mg/ml

pH optimum: 4.5-8.5

Stability: At pH 6 and heating $> +60^\circ C$ this product decomposes with separation of formaldehyde.

Specification

Appearance: Fine white crystalline powder

Solubility: Clear, colorless solution in water (c=1%, w/v)

pH value (c=1%, w/v): 6.0-7.5

Identity (NIR-spectrum): Corresponds to reference

Germall 115 (from N): $\geq 94\%$

Water (K. Fischer): $\leq 6\%$

HPLC: Corresponds to reference

Heavy metals (as Pb): ≤ 20 ppm

Reducible components (calculated as H_2O_2): $\leq 0.03\%$

Oxidizable components (calculated as O): $\leq 0.2\%$

Formaldehyde (Nash, photometrically): $\geq 9.0\%$

Stability: At $+15$ to $+40^\circ C$ within specification range for 36 months. Store dry in tightly sealed containers.

Cat. No.	Pack Size
10 235 733 103	custom fill

Will be supplied as "Imidazolidinylurea". Unit of Measure is "kg".
For further processing only.

Germall 115, reduced sodium crystalline powder

Bacteria specialized biocide used as preservative in diagnostic kits

Application

Use Germall 115, reduced sodium as preservative in combination with other biocides to obtain the best possible antimicrobial effect in tests that require a low concentration of sodium.

Cat. No.	Pack Size
11 276 883 103	custom fill

Unit of Measure is "g".
For further processing only.

Benefits

- Take advantage of the strongly reduced concentration of sodium.
- Combine Germall 115, reduced sodium with non-ionic, ionic detergents or proteins.
- Use Germall 115, reduced sodium in your preferred combination in liquid reagents to obtain the most stable system.

CAS: 39236-46-9

Properties

Formula: $C_{11}H_{16}N_8O_8 \cdot x H_2O$

Molecular weight: 406.31 D

Antimicrobial effect: Only moderate effect against fungi and yeast; should be used in combination.

Possible combination: Can be combined with all biozides which do not separate formaldehyde, especially 5-bromo-5-nitro-1,3-dioxane (BND), 2-hydroxypyridine-N-oxide (HPO) or N-methylisothiazolone (MIT).

Effect: Belongs to the "formaldehyde-depot-substances". Inactivating effect due to reactions with $-COOH$, $-NH_2$, $-OH$ and $-SH$ of all proteins.

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentration: 2.6 mg/ml

pH optimum: 4.5-8.5

Stability: At pH 6 and heating up to $>+60^\circ C$ this product decomposes with separation of formaldehyde.

Specification

Appearance: White, crystalline powder

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml, water): 2.8-3.2

Identity (NIR-spectrum): Corresponds to reference

Germall 115 (from N) : $\geq 94\%$

N (elementary analysis): $\geq 26\%$

Water (K. Fischer): $\leq 5\%$

HPTLC: Corresponds to reference

Heavy metals (as Pb): ≤ 20 ppm

Na (AES): ≤ 250 ppm

Reducing substances (as O): $\leq 0.5\%$

Oxidizing substances: $\leq 0.03\%$

Stability: At $+2$ to $+8^\circ C$ within specification range for 36 months.

Micr-O-protect solution

Very effective broad range biocide mix for use in diagnostic kits.

Application

Use Micr-o-protect as effective preservative in reagents, containing buffers, detergents, enzymes and other proteins.

Benefits

- Rely on the highly efficient ready-to-use biozide mix that can be used in most used organic or inorganic buffer systems.
- Combine Micr-O-protect with non-ionic or ionic detergents or proteins.
- Use Micr-O-protect in your preferred combination in liquid reagents to obtain the most stable system.

CAS: 30007-47-7 and 2682-20-4

Cat. No.	Pack Size
11 587 056 103	500 ml

Will be supplied as "Micr-o-protect Biocide Mix". Unit of Measure is "l".

For further processing only.

Properties

Antimicrobial effect: Due to the compatibility of both biocides used, the preparation is extraordinarily efficient at low working concentrations against pro- and eukaryotic microorganisms. Moreover, the combination of both substances provides a higher degree of protection, because development of resistance is prevented.

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentration: 0.1-0.4% (v/v)

pH Optimum: 5.0-7.5

Stability: Stability is lowered at pH >8.5 or <5, at temperatures >+40°C, in the presence of reducing and oxidizing substances and in the presence of strong nucleophilic substances.

Compatibility: Micr-o-protect is well compatible with detergents, all commonly used organic and inorganic buffers and proteins.

Specification

Appearance: Clear, colorless to yellow solution

Bromonitrodioxane (GC): ≥98.0 area%

N-Methylisothiazolon (GC): ≥99.0 area%

Thin layer chromatography (TLC): Corresponds to reference

Stability: At +2 to +8°C within specification range for 24 months.

N-Methylisothiazolone (MIT) crystalline powder

Broad range biocide for use in diagnostic kits

Application

Use MIT as preservative in combination with other biocides to obtain the best possible antimicrobial effect.

Benefits

- Combine MIT with non-ionic or ionic detergents or proteins.
- Use MIT in your favorite combination in liquid reagents to get the most stable system.
- Use MIT in potassium phosphate, Hepes or Tris buffer systems.

CAS: 2682-20-4

Properties

Formula: C₄H₆NOSCl

Nomenclature: 2-methyl-3(2 H)-isothiazolone-hydrochloride

Molecular weight: 151.62 D (MIT: 115.12 D)

Antimicrobial effect: Highly effective against bacteria, yeast and fungi.

Possible combination: Can be well combined with all preservatives. Combination with preservatives that have a good effect against fungi (e.g., 5-bromo-5-nitro-1,3-dioxan (BND), 2-hydroxypyridine-N-oxide (HPO)) results in a very broad spectrum of effectiveness.

Effect: Thought to react with -NH₂ and -SH groups of all proteins.

Interferences: Probably possible with amines and sulfides.

Toxicity: Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).

Recommended working concentrations:

0.4 mg/ml, if both risk of contamination and concentration of protein are low (<0.2%);

1.0 mg/ml, if both risk of contamination and concentration are high.

pH optimum: 4.0-7.0

Stability: Decreases at pH >8.5 and >+50°C; no loss of activity when heated up to +100°C for a short period of time.

Cat. No.	Pack Size
11 085 905 103	custom fill

Will be supplied as "N-Methylisothiazolon-HCl". Unit of Measure is "kg".

For further processing only.

Specification**Appearance:** White to light beige crystalline powder**Solubility:** Clear, colorless solution in water (c=10 mg/ml), free from fuzz**A₄₀₅** (10 mg/ml water, against water): ≤0.025**Melting range** (Büchi): +160 to +180°C**Identity** (NIR): Corresponds to reference**N-Methylisothiazolone** (GC): ≥99 area%**Water** (K. Fischer): ≤0.5%**Impurities** (GC): ≤1 area%**Heavy metals** (as Pb): ≤10 ppm**Oxidizing substances:** ≤2%**Acetone** (GC): ≤0.2%**2-Propanol** (GC): ≤0.2%**Petrolether** (GC): Beyond the limit of detection (approximately 1%); Corresponds to reference**Elecsys® Hbs Ag:** Corresponds**Stability:** At +2 to +8°C within specification range for 24 months.**N-Methylisothiazolone (MIT) for potassium test crystalline powder**

Broad range biocide for use in diagnostic kits

Application

Use MIT for potassium test as preservative in combination with other biocides to obtain the best possible antimicrobial effect in tests that require a low concentration of potassium.

Benefits

- Take advantage of the strongly reduced concentration of potassium.
- Combine MIT with non-ionic or ionic detergents or proteins.
- Use MIT in your favorite combination in liquid reagents to get the most stable system.
- Use MIT in phosphate, HEPES or Tris buffer systems.

CAS: 2682-20-4**Properties****Formula:** C₄H₆NOSCl**Nomenclature:** 2-methyl-3(2H)-isothiazolone-hydrochloride**Molecular weight:** 151.62 D (MIT 115.12 D)**Antimicrobial effect:** Highly effective against bacteria, yeast and fungi.**Possible combination:** Can be well combined with all preservatives. Combination with preservatives that have a good effect against fungi (e.g., 5-bromo-5-nitro-1,3-dioxan (BND), 2-hydroxypyridine-N-oxide (HPO)) results in a very broad spectrum of effectiveness.**Effect:** Thought to react with -NH₂ and -SH groups of all proteins.**Interferences:** Probably possible with amines and sulfides.**Toxicity:** Harmless at recommended concentrations; not carcinogenic or mutagenic (Ames-test negative).**Recommended working concentrations:**

0.4 mg/ml, if both risk of contamination and concentration of protein are low (<0.2%);

1.0 mg/ml, if both risk of contamination and concentration are high.

pH optimum: 4.0-7.0**Stability:** Decreases at pH >8.5 and >+50°C; no loss of activity when heated up to +100°C for a short period of time.**Cat. No.****11 333 917 103****Pack Size**

custom fill

Will be supplied as "N-Methylisothiazolone for Potassium Test". Unit of Measure is "kg".

For further processing only.

Specification

Appearance: White to light beige crystalline powder

Solubility: Clear, colorless solution in water (c=10 mg/ml), free from fuzz

A₄₀₅ (10 mg/ml water, against water): ≤0.025

Melting range (Büchi): +160 to +180°C

Identity (NIR): Corresponds to reference

N-Methylisothiazolone (GC): ≥99 area%

Water (K. Fischer): ≤0.5%

Impurities (GC): ≤1 area%

Na (AAS): ≤200 ppm

K (AAS): ≤10 ppm

Heavy metals (as Pb): ≤10 ppm

Oxidizing substances: ≤2%

Acetone (GC): ≤0.2%

2-Propanol: ≤0.2%

Petrolether (GC): Beyond the limit of detection (approximately 1%); Corresponds to reference

Elecsys® Hbs Ag: Corresponds

Stability: At +2 to +8°C within specification range for 24 months.

Diethanolamine 85% solution

Buffer for enzymatic assays of alkaline phosphatase.

Application

Use Diethanolamine 85% as a buffer in applications that test or use the activity of alkaline phosphatase, such as diagnostic tests for alkaline phosphatase.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 111-42-2

Properties

Formula: C₄H₁₁O₂N

Molecular weight: 105.14 D

Solubility: Miscible with water and ethanol

Suggested pH range: Approximately 10

Specification

Appearance: Colorless, clear liquid

Refractive index: n 20/D: 1.4575-1.4595

Density: D 20/20: 1.092-1.095

Coloration of sample (against water):

A₄₀₅: ≤0.043

A₄₀₅ (10 days at +4°C): ≤0.043

A₄₀₅ (10 days at +35°C): ≤0.051

Coloration of sample in buffer solution (against water):

A₄₀₅: ≤0.01

A₄₀₅ (10 days at +4°C): ≤0.02

A₄₀₅ (10 days at +35°C): ≤0.01

Diethanolamine (HClO₄ titration): 84.0-86.0%

Thin layer chromatography: Corresponds to reference

Mono-ethanolamine (TLC): ≤0.4%

Fe (calculated on Diethanolamine 100%): ≤2 ppm

Heavy metals (as Pb, calculated on Diethanolamine 100%): ≤5 ppm

Water (K. Fischer): 14.0-16.0%

Stability: At +2 to +8°C within specification range 18 months. Store dry. Protect from light. Keep in tightly sealed containers.

Cat. No. **Pack Size**

10 201 294 001 custom fill

Will be supplied as "Diethanolamine 85%". Unit of Measure is "kg". For further processing only.

Glycylglycine crystalline powder

Zwitterionic buffer for diagnostic tests or as a substrate for γ-glutamyltransferase tests.

Application

Use Glycylglycine as a buffer in diagnostic reagents, or as a substrate in γ-glutamyltransferase tests where γ-glutamyltransferase transfers the gamma-glutamyl group of L-γ-glutamyl-3-carboxy-4-nitroanilide (Glupa-C) to glycylglycine.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

Cat. No. **Pack Size**

10 002 887 103 custom fill

Will be supplied as "Glycylglycin". Unit of Measure is "kg". For further processing only.

CAS: 556-50-3

Properties

Formula: $C_4H_8N_2O_3$

Molecular weight: 132.1 D

Solubility: Easily soluble in water

Suggested pH range: 7.5-8.9

Specification

Appearance: White crystalline powder

pH value (c=0.35%, w/v): 5.5-6.5

Glycylglycine (HClO₄ titration, based on anhydrous substance): 99.0-100.5%

TLC: Corresponds to reference

A₄₀₅ (c=10%; w/v, against water): ≤0.01%

Heavy metals (as Pb): ≤5 ppm

Water (K. Fischer): ≤0.5%

Fe: ≤10 ppm

Sulphate ash: ≤0.1%

Glycine (TLC): ≤0.2%

Contaminating amino acids (TLC): Not detectable

Microbiological test: Corresponds to specification

IR Spectrum: Corresponds to reference

Stability: At +15 to +40°C within specification range for 36 months. Store dry in tightly sealed containers.

Hepes crystalline powder

Buffer for diagnostic tests, such as amylase test.

Application

Use Hepes as a buffer in reagents.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 7365-45-9

Properties

Nomenclature: 4-(2-Hydroxyethyl)-1-piperazineethanesulfonic acid

Formula: $C_8H_{18}N_2O_4S$

Molecular weight: 238.3 D

Suggested pH range: 6.8-8.2

Specification

Appearance: White crystalline powder

Solubility: Clear, colorless solution in water (c=1 mol/l)

pK value: 7.21-7.41 (+37°C); 7.45-7.65 (+20°C)

Melting range: +207 to +213°C

Hepes (alcalimetric): ≥97%

Hepes (from N): ≥97%

N (elementary analysis): ≥11.4%

Thin layer chromatography (TLC): Chromatographically homogeneous

A₂₆₀ (against water): ≤0.050

A₄₀₅ (against water): ≤0.030

Cl (chloride meter): ≤0.04%

Hepes mA/min (Purity check, α-amylase contamination): ≤0.1

Exclusion of skin contact and contamination with saliva: Corresponds to specification.

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No.	Pack Size
10 172 944 103	custom fill

Will be supplied as "Hepes". Unit of Measure is "kg".
For further processing only.

Imidazole crystalline powder

Buffer for diagnostic tests, such as creatine kinase test.

Application

Use Imidazole as a buffer for diagnostic tests and other reagents, especially enzymatic reactions such as creatine kinase test.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 288-32-4

Properties

Formula: C₃H₄N₂

Molecular weight: 68.08 D

Suggested pH range: 6.2-7.8

Specification

Appearance: White or slightly yellowish crystallizate

Melting point: +86 to +91°C

Imidazol (titrimetric): ≥99.0%

dA₂₅₀ - dA₃₆₀: ≤0.050

dA₃₃₄: ≤0.050

dA₄₀₅: ≤0.010

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No. **Pack Size**

10 034 428 103 custom fill

Will be supplied as "Imidazole". Unit of Measure is "kg".
For further processing only.

MES crystallizate

Buffer for a variety of diagnostic tests, such as Glucose test.

Application

Use MES as a buffer in reagents that require a pH of approximately 6.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 4432-31-9

Properties

Nomenclature: 4-Morpholineethane Sulphonic Acid

Formula: C₆H₁₃NO₄S

Molecular weight: 195.2 D

Suggested pH range: 5.5-6.7

Specification

Appearance: White crystallized substance

Solubility: Clear, colorless solution in water (c=50 mg/ml)

pK value: 6.2±0.2 (+25°C)

Equivalent point: 8.9±0.3

MES (alkalimetric): ≥98%

MES (from N): ≥98%

N (elementary analysis): ≥7%

Cat. No. **Pack Size**

10 073 571 103 custom fill

Will be supplied as "4-Morpholineethane Sulfonic Acid (Mes)".
Unit of Measure is "kg".
For further processing only.

A₂₆₀ (c=10 mg/ml, neutralized): ≤0.05

Heavy metals (as Pb): ≤10 ppm

Br (chloride meter): ≤0.5%

Thin layer chromatography (TLC): Chromatographically homogeneous

Stability: At +15 to +25°C within specification range for 24 months.

Pipes free acid

Buffer for diagnostic tests, such as cholesterol test.

Application

Use Pipes as a buffer in a variety of diagnostic tests, especially in tests for cholesterol and triglycerides.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 5625-37-6

Properties

Nomenclature: Piperazine-1,4-bis-2-ethane sulfonic acid

Formula: C₈H₁₈N₂O₆S₂

Molecular weight: 302.4 D

Suggested pH range: 6.1-7.5

Specification

Appearance: Colorless crystals

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pK value: 6.8±0.1

Pipes (alkalimetric): ≥98%

Pipes (from N): ≥98%

A₂₆₀ (c=10 mg/ml water) : ≤0.05

Thin layer chromatography (TLC): Chromatographically homogeneous

Br (chloride meter): ≤0.5%

Identity (IR spectrum): Corresponds to reference

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No.	Pack Size
10 239 500 103	custom fill

Will be supplied as "Pipes, Free Acid". Unit of Measure is "kg".
For further processing only.

Pipes disodium salt

Buffer for diagnostic tests, such as cholesterol test.

Application

Use Pipes as a buffer in a variety of diagnostic tests, especially in tests for cholesterol and triglycerides.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 5625-37-6

Properties

Nomenclature: Piperazine-1,4-bis-2-ethane sulfonic acid

Formula: C₈H₁₆N₂O₆S₂Na₂

Molecular weight: 346.3 D

Cat. No.	Pack Size
10 735 361 103	custom fill

Will be supplied as "Pipes, Sodium Salt". Unit of Measure is "kg".
For further processing only.

Suggested pH range: 6.1-7.5

Specification

Appearance: Colorless powder

Solubility: Clear, colorless solution in water (c=10 mg/ml)

Pipes (from N): ≥82%

Na (flame photometric): 11-14%

Water (K. Fischer): ≤5%

N (elementary analysis) : ≥7.6%

A₃₄₀ (c=10 mg/ml water): ≤0.01

HPTLC: Chromatographically homogeneous

Identity (IR spectrum): Corresponds to reference

Stability: At +15 to +25°C within specification range for 24 months.

Tris crystallizate

Buffer for diagnostic tests, such as tests for aminotransferases.

Application

Use Tris as a buffer in diagnostic reagents, especially in tests for aminotransferases or γ-glutamyltransferase.

Benefits

- Rely on the proven quality of this product for the manufacturing of diagnostic tests.

CAS: 77-86-1

Properties

Nomenclature: Tris(hydroxymethyl)-amminomethane

Formula: C₄H₁₁NO₃

Molecular weight: 121.1 D

Suggested pH range: 7.0-9.0

Specification

Appearance: Colorless, odorless crystallizate

Solubility: Clear, colorless solution in water (c=100 mg/ml), free of fuzz

Flow properties: Passes

Melting range: +168 to +171°C

Conductivity (water, 1 μS, +25°C): ≤110 μS

pH value (c=6 mg/ml, in water): 10.0-11.0

Tris (titrimetric, based on dry weight): 99.5-100.5%

Water (K. Fischer): ≤0.2%

Sulfate ash (with concentrated H₂SO₄ at +600°C): ≤0.05%

Fe (AAS): ≤1 ppm

As (AAS): ≤1 ppm

Heavy metals (as Pb): ≤1 ppm

Reducing substances (KMnO₄, 0.002 mol/l): ≤3 ml/100 mg

Acetone (GC): ≤0.05%

Methanol (GC): ≤0.05%

Cl (turbidimetric test with AgNO₃): ≤20 ppm

Bioburden: ≤100 CFU/g

A₃₀₀ (against water, c=100 mg/ml): ≤0.020

A₄₀₅ (against water, c=100 mg/ml): ≤0.004

Stability: At +15 to +25°C within specification range for 24 months. Protect from light.

Cat. No. **Pack Size**

10 153 265 001 custom fill

Will be supplied as "Tris-(hydroxymethyl)-aminomethane". Unit of Measure is "kg".

For further processing only.

Detergents

Application

Detergents are used in diagnostic kits:

- To enhance the solubility of test or sample compounds
- To activate enzymes such as esterases
- To reduce interference from serum lipids
- To reduce carryover effects on analyzers and facilitate dispensing processes
- To reduce non-specific binding to solid phases in immunoassays
- To pretreat samples
- To make use of possible antimicrobial effects

Properties

The choice of a specific detergent depends on the protein to be solubilized, the need of removal, toxicity data, interference with UV-VIS-absorption, interference with subsequent isoelectric focusing or ultracentrifugation.

The following table can therefore only give a first approximation for the optimal choice of a detergent.

	Ability to disperse protein aggregates	Denaturation of protein	Ease of removal
Non-ionic detergents			
n-Octyl- β -D-glycoside	Low	No	Very easy
Polidocanol (Thesit)	Low	No	Difficult
Triton X-100	Low	No	Difficult
n-Dodecyl- β -D-maltoside	Low	No	Difficult
Nonidet P40	Low	No	Difficult
Tween 20	Low	No	Difficult
Ionic detergents			
Cholate	Low	(no)	Very easy
Deoxycholate	High	(no)	easy
Zwitterionic detergents			
CHAPS	High	No	Very easy

Detergents

Ionic Detergents

Dilaurylglycerosulfate powder

Detergent for diagnostic tests.

Application

Use Dilaurylglycerosulfate as a co-emulsifier in the diagnostic test for the determination of lipase.

Benefits

- Enhance the performance of your lipase test.

CAS: 99387-94-7

Properties

Formula: C₂₇H₅₆O₆S

Molecular weight: 508.8 D

Specification

Appearance: White powder

C (elementary analysis): ≥57.3%

H (elementary analysis): ≥10.0%

Water (K. Fischer): ≤5%

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No. **Pack Size**

11 827 294 103 custom fill

Will be supplied as "Dilaurylglycerosulfat". Unit of Measure is "g".
For further processing only.

Cholate ionic detergent, sodium salt

Anionic detergent for diagnostic tests.

Application

Use Cholate in diagnostic reagents, such as for the determination of cholesterol and triglycerides.

Benefits

- Enhance the solubility of your reagent.
- Reduce interferences with serum lipids or carryover effects.
- Rely on the proven diagnostic quality of this product.

CAS: 81-25-4

Properties

Formula: C₂₄H₃₉NaO₅

Molecular weight: 430.6 D

Detergent type: Anionic detergent

Solubility: Limited solubility in the presence of Ca²⁺

Handling advice: Harmful if exposed to skin and if inhaled. Adequate precautions as for handling of irrigating products must be taken.

Specification

Appearance: White crystalline powder

Solubility:

Clear, colorless solution in water (c=10 mg/ml)

Clear, colorless to yellowish solution in water (c=150 mg/ml, +20°C)

Identity (NIR): Corresponds to reference

Cholic acid, Na-salt (HPLC): ≥93 area%

Water (K. Fischer): ≤6%

C (elementary analysis): 62-67%

Na (flame photometric): 5.0-5.5%

Heavy metals (as Pb): ≤10 ppm

Cat. No. **Pack Size**

10 261 084 103 custom fill

Will be supplied as "Cholic Acid Sodium Salt". Unit of Measure is "kg".
For further processing only.

Flame coloration: Positive
A₃₄₀ (against water): ≤0.100
A₅₀₅ (against water): ≤0.005
A₅₄₆ (against water): ≤0.005
A₅₀₅ to A₅₅₀ (against water): ≤0.025
Hydrophilic contaminants (HPLC): ≤15 area%
Lipophilic contaminants (HPLC): ≤4.0 area%
Reducing substances: ≤0.25 ml (KMnO₄, 0.002 mol/l, per 100 mg)
Oxidizing substances: Negative
Bioburden: ≤100 CFU/g, ≤10 moulds/g
Performance: Corresponds to specification
County of origin: USA, Australia
Stability: At +15 to +25°C within specification range for 36 months. Protect from light.

Deoxycholate ionic detergent, sodium salt

Anionic detergent for diagnostic tests.

Application

Use Deoxycholate in diagnostic reagents, such as for the determination of lipase.

Benefits

- Enhance the solubility of your reagent.
- Reduce protein-protein interactions in your reagent.
- Rely on the proven diagnostic quality of this product.

CAS: 83-44-3

Properties

Formula: C₂₄H₃₉O₄Na

Molecular weight: 414.6 D

Detergent type: Anionic detergent

Handling advice: Harmful if exposed to skin and if inhaled. Adequate precautions as for handling of irrigating products must be taken.

Specification

Appearance: White, crystalline powder

Solubility: Clear, colorless solution in water (c=50 mg/ml), free from fuzz

Deoxy cholic acid (sodium deoxycholate, HPLC): ≥83.0 area%

Hydrophilic contaminants: ≤15.0 area%

Lipophilic contaminants: ≤5.0 area%

Water (K. Fischer): ≤10.0%

Acetone (GC): ≤0.5%

Heavy metals (as Pb): ≤5 ppm

Stability: At +15 to +25°C within specification range for 36 months.

Cat. No.	Pack Size
11 434 314 103	custom fill

Will be supplied as "Desoxycholat, Mono-NA, Crystal". Unit of Measure is "kg".

For further processing only.

Taurodesoxycholat sodium salt

Anionic detergent for diagnostic tests.

Application

Use Taurodesoxycholat in diagnostic reagents, such as for the determination of lipase.

Cat. No.	Pack Size
11 332 686 103	custom fill

Will be supplied as "Taurodesoxycholic. acid, Na, pur.". Unit of Measure is "g".

For further processing only.

Benefits

- Enhance the solubility of your reagent.
- Activate the lipase in your reagent.
- Rely on the proven diagnostic quality of this product.

CAS: 1180-95-6

Properties

Formula: C₂₆H₄₄NO₆SNa

Molecular weight: 521.7 D

Detergent type: Anionic detergent

Handling advice: Harmful if exposed to skin and if inhaled. Adequate precautions for handling hazardous products must be used.

Specification

Appearance: White lyophilizate

Taurodesoxy cholate, Na (from C): ≥90%

Taurodesoxy cholate, Na (HPLC): ≥89 area%

C (elementary analysis): 53.9-61.0%

H (elementary analysis): 7.9-8.9%

N (elementary analysis): 2.4-3.0%

Na (flame photometric): 4.4-6.6%

Water (K. Fischer): ≤5%

Stability: At +2 to +8°C within specification range for 24 months.

n-Dodecyl-β-D-maltoside nonionic detergent, powder

Nonionic detergent for diagnostic tests

Application

Use n-Dodecyl-β-D-maltoside as a mild, nondenaturing detergent for the solubilization of proteins, especially antibodies.

Benefits

- Obtain very gentle conditions that stabilize and preserve your reagent activity.
- Rely on the proven diagnostic quality of this product.

CAS: 69227-93-6

Properties

Nomenclature: 1-O-n-Dodecyl-β-D-glucopyranosyl(1-4)α-D-glucopyranoside

Formula: C₂₄H₄₆O₁₁

Molecular weight: 510.62 D

Detergent type: Nonionic alkyl maltoside type

Specification

Appearance: White, crystalline powder

Specific rotation [α]_D²⁵ (in MeOH): +46.0±2.0°

n-Dodecylmaltoside (from C): ≥98%

C (elementary analysis): ≥55.20%

Dodecanol (GC): ≤0.1%

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No. **Pack Size**

10 808 342 103 custom fill

Will be supplied as "n-Dodecyl-β-D-maltoside". Unit of Measure is "g".

For further processing only.

n-Octyl β-D-glucoside nonionic detergent, powder

Nonionic detergent for diagnostic tests

Application

Use n-Octyl β-D-glucoside as a mild, nondenaturing detergent for the solubilization of proteins, especially antibodies. n-Octyl β-D-glucoside can be easily removed by dialysis.

Benefits

- Rely on the proven diagnostic quality of this product.
- Take advantage of the mild conditions that stabilize your valuable reagent.

CAS: 29836-26-8

Properties

Nomenclature: 1-O-Octyl-β-D-glucopyranoside

Formula: C₁₄H₂₈O₆

Molecular weight: 292.4 D

Detergent type: Nonionic alkyl glucoside type

pH stability: Stable in solutions above pH 6.5

Specification

Appearance: White powder

n-Octylglucoside (from C): ≥99%

C (elementary analysis): ≥56.9%

Octanol (GC): ≤0.1%

Stability: At +15 to +25°C within specification range for 36 months. Store dry.

Cat. No. **Pack Size**

10 411 469 103 custom fill

Will be supplied as "n-Octylglucoside". Unit of Measure is "g".

For further processing only.

Detergents

Non-ionic Detergents

Nonidet P40 nonionic detergent, aqueous solution

Nonionic detergent for diagnostic tests.

Application

Use the nondenaturing detergent Nonidet P40 for the solubilization of proteins, especially antibodies.

Benefits

- Rely on the low amount of peroxides and carbonyl groups.
- Stabilize your valuable reagent.

CAS: 9016-45-9

Properties

Formula: $C_{33}H_{60}O_{10}$ (n=9)

Molecular weight: 616,83 D

Detergent type: Nonionic detergent, polyethylene type, nondenaturing

Handling advice: Harmful if exposed to skin and if inhaled. Adequate precautions as for handling of irrigating products must be taken.

Specification

Appearance: Clear, colorless solution

Conductivity: ≤ 100 μ S/cm

Peroxides: ≤ 2 ppm

Carbonyl groups (MBTHS): ≤ 0.02 mg/ml

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No. **Pack Size**

11 333 941 103 custom fill

Will be supplied as "Nonidet P40, cleaned, Solution". Unit of Measure is "l".

For further processing only.

Polidocanol (Thesit)

Nonionic detergent for use in diagnostic reagents.

Application

Use Polidocanol in diagnostic reagents to enhance solubility, minimize interferences and reduce carryover effects.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 9002-92-0

Properties

Nomenclature: Dodecylpolyethyleneglycolether

Formula: $C_{30}H_{62}O_{10}$ (n=approximately 9)

Molecular weight: Approximately 600 D

Detergent type: Nonionic polyoxyethylene type

Handling advice: Polidocanol must be moderately heated (+40 to +50°C) and carefully homogenized by gentle stirring before dispensing.

Specification

Appearance: White, pasty, fatty substance; clear, colorless to slightly yellow liquid at approximately +30°C

Solubility: Clear, colorless solution in water (c=100 mg/ml)

Peroxide (as H₂O₂): ≤ 1 ppm

Stability: At +2 to +8°C within specification range for 24 months. Keep under argon or nitrogen. Protect from light.

Cat. No. **Pack Size**

10 831 620 103 custom fill

Will be supplied as "Polidocanol (PEG Monododecyl Ether)". Unit of Measure is "kg".

For further processing only.

Triton X-100 nonionic detergent, viscous liquid

Nonionic detergent for diagnostic tests.

Application

Use Triton-X 100 for the solubilization of proteins, especially antibodies.

Benefits

- Stabilize your valuable reagent
- Rely on the proven diagnostic quality of this product.

CAS: 9002-93-1

Properties

Nomenclature: Octylphenolpoly(ethyleneglycolether)_n, n=10

Formula: C₃₄H₆₂O₁₁

Molecular weight: 647 D

Detergent type: Nonionic polyethylene type

Handling advice: Triton X-100 must be homogenized carefully before dispensing at +20 to +30°C.

Specification

Appearance: Clear, colorless liquid

Triton X-100 (GC): Corresponds to standard

Peroxide (as H₂O₂): ≤1 ppm

Stability: At +2 to +8°C within specification range for 12 months. Keep under argon or nitrogen. Protect from light. Keep tightly sealed.

Cat. No.	Pack Size
10 743 119 103	custom fill

Will be supplied as "Triton X-100, for Membrane Research". Unit of Measure is "l".

For further processing only.

Tween 20 purified, solution

Non-ionic detergent for diagnostic tests

Application

Use Tween 20 as non-denaturing detergent for the solubilization of proteins, especially antibodies.

Benefits

- Stabilize your valuable reagent
- Rely on the low amount of peroxides and carbonyl groups and salts in this product.

CAS: 9005-64-5

Properties

Formula: C₅₈H₁₁₄O₂₆ (for w+x+y+z=n=20)

Molecular weight: 1228 g/mol

Specification

Appearance: Clear, yellow solution

Conductivity: ≤100 μS/cm

Peroxide (as H₂O₂): ≤2 ppm

Aldehyde: ≤0.02 mg/ml

Stability: At +2 to +8°C within specification range for 24 months. Store under nitrogen. Protect from light.

Cat. No.	Pack Size
11 334 000 103	custom fill

Will be supplied as "Tween 20, gereinigt, Lsg". Unit of Measure is "l".

For further processing only.

Detergents

Zwitterionic Detergents

CHAPS

zwitterionic detergent, crystalline powder

Zwitterionic detergent for diagnostic tests.

Application

Use CHAPS to reduce protein-protein interactions. CHAPS can be easily removed by dialysis.

Benefits

- Minimize protein-protein interactions in your reagent without denaturing the protein or enzyme.
- Rely on the proven diagnostic quality of this product.

CAS: 75621-03-3

Properties

Nomenclature: 3-[[3-(Cholamidopropyl)dimethylammonio]-1-propanesulfonate

Formula: C₃₂H₅₈N₂O₇S

Molecular weight: 614.9 D

Detergent type: Zwitterionic detergent, nondenaturing

Specification

Appearance: White, crystalline powder

CHAPS (from N): ≥99%

N (elementary analysis): ≥4.49%

A₂₆₀ (against water): ≤0.10

A₂₈₀ (against water): ≤0.10

Thin layer chromatography (TLC): Chromatographically homogeneous, corresponds to reference

Stability: At +2 to +8°C within specification range for 24 months. Store dry. Protect from light.

Cat. No.

Pack Size

10 810 681 103 custom fill

Will be supplied as "CHAPS". Unit of Measure is "kg".
For further processing only.

CHAPSO

zwitterionic detergent, crystalline powder

Zwitterionic detergent for diagnostic tests.

Application

Use CHAPSO to reduce protein-protein interactions. CHAPSO can be easily removed by dialysis.

Benefits

- Minimize protein-protein interactions in your reagent without denaturing the protein or enzyme.
- Rely on the proven diagnostic quality of this product.

CAS: 82473-24-3

Properties

Nomenclature: 3-[[3-(Cholamidopropyl)dimethylammonio]-2-hydroxy-1-propansulfonate

Formula: C₃₂H₅₈N₂O₈S

Molecular weight: 630.9 D

Detergent type: Zwitterionic detergent, similar to CHAPS but more soluble.

Handling advice: Harmful if exposed to skin and if inhaled. Adequate precautions as for handling of irrigating products must be taken.

Cat. No.

Pack Size

11 112 392 103 custom fill

Will be supplied as "CHAPSO". Unit of Measure is "g".
For further processing only.

Specification

Appearance: White crystalline powder

CHAPSO (from C): ≥95%

C (elementary analysis): ≥57.6%

Water (K. Fischer): ≤4%

HPTLC: Chromatographically homogeneous

Stability: At +2 to +8°C within specification range for 36 months.

Zwittergent 3-14 zwitterionic detergent, powder

Synthetic zwitterionic detergent for diagnostic tests.

Application

Use Zwittergent 3-14 to reduce protein-protein interactions.

Benefits

- Minimize protein-protein interactions in your reagent without denaturing the protein or enzyme.
- Rely on the proven diagnostic quality of this product.

CAS: 14933-09-6

Properties

Nomenclature: N-Tetradecyl-N,N-dimethyl-3-ammonio-1-propane-sulfonate

Formula: C₁₉H₄₁NO₃S

Molecular weight: 363.65 D

Detergent type: Zwitterionic detergent

Handling advice: Harmful if exposed to skin and if inhaled. Adequate precautions as for handling of irrigating products must be taken.

Specification

Appearance: White powder

UV spectrum (200-400 nm; against water): Corresponds to specification

A₂₂₅ (against water): ≤0.500

A₂₆₀ (against water): ≤0.150

A₂₈₀ (against water): ≤0.150

A₃₂₅ (against water): ≤0.100

N-Tetradecyl-N,N-dimethyl-3-ammonio-1-propane-sulfonate (from C): ≥98.0%

N-Tetradecyl-N,N-dimethyl-3-ammonio-1-propane-sulfonate (HPLC): 95.0 area%

C (elementary analysis): ≥61.40%

H (elementary analysis): ≥11.00%

N (elementary analysis): ≥3.30%

Stability: At +15 to +25°C within specification range for 36 months.

Cat. No. **Pack Size**

11 112 902 103 custom fill

Will be supplied as "N-Tetradec-N,N-dimet-3-am-1-propSulfonat".
Unit of Measure is "g".

For further processing only.

Aprotinin**from bovine lung, lyophilizate**

Protease inhibitor

Application

Use Aprotinin in reagents to inhibit serine proteases, such as kallikrein, plasmin, trypsin and chymotrypsin.

Benefits

- Take advantage of the excellent abilities of Aprotinin to inhibit serine proteases.
- Rely on the proven diagnostic quality of this product.

CAS: 9087-70-1

Specification

Appearance: White lyophilizate

Activity (Chromozym TRY, +25°C): ≥630 inhibitor U/mg lyophilizate

Activity (BAEE, +25°C): ≥200 inhibitor U/mg lyophilizate

Protein (Lowry): 90-100%

Electrophoresis (SDS Page): Corresponds to reference

Country of origin: USA, South Africa, New Zealand, Australia, or Uruguay, respective

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Remarks:

Official veterinary certificate of health of the donor animals is available.

Official certificate of the deactivation of animal material including the method (acid treatment at up to pH 5 for up to 5 h) is available.

Cat. No. **Pack Size**

10 236 632 103 custom fill

Will be supplied as "Aprotinin from Bovine Lung". Unit of Measure is "g".

For further processing only.

Pefabloc SC (AEBSF)**4-(2-Aminoethyl)-benzenesulfonyl fluoride hydrochloride, powder**

Protease inhibitor

Application

Use Pefabloc in reagents to inhibit serine proteases, such as thrombin in serum or plasma.

Benefits

- Take advantage of the excellent abilities of Pefabloc to inhibit thrombin.
- Rely on the low toxicity and enhanced stability of this product compared to PMSF or DFP.

CAS: 34284-75-8

Properties

Molecular weight: 239.5 D

Specification

Appearance: White powder

Pefabloc HCl (HPLC): ≥90 area%

Pefabloc HCl (from C): ≥95%

C (elementary analysis): ≥38.1%

H (elementary analysis): ≥4.2%

N (elementary analysis): ≥5.5%

Thin layer chromatography (TLC): Chromatographically homogeneous

Inhibition chymotrypsin: Corresponds to specification

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No. **Pack Size**

11 427 393 103 custom fill

Will be supplied as "Pefabloc SC AEBSF, Hydrochloride". Unit of Measure is "g".



For further processing only.

1,4-Dithiothreitol (DTT) crystallizate

Reducing agent

Application

Use 1,4-Dithiothreitol primarily to protect free SH-groups from oxidation. Use it routinely in all work with enzymes and proteins during enzyme measurement and the characterization of proteins.

Benefits

- Rely on the higher stability of DTT in comparison to 2-mercaptoethanol in aqueous solution
- Experience more convenience and a less disagreeable odor
- Benefit from the reduced tendency to oxidize in air.

CAS: 3483-12-3

Properties

Nomenclature: Threo-1,4-dimercapto-2,3-butanediol

Formula: C₄H₁₀O₂S₂

Molecular weight: 154.3 D

Specification

Appearance: White to yellowish crystallizate

DTT (with Ellman's reagent): ≥97%

Thin layer chromatography: Chromatographically homogeneous, corresponds to reference

Stability: At +2 to +8°C within specification range for 24 months. Store dry. Protect from light.

Cat. No. **Pack Size**

10 197 785 103 custom fill

Will be supplied as "Dithiothreitol (DTT) Cleland's Reagent". Unit of Measure is "g".

For further processing only.

3-Hydroxy-1,2,3,4-tetrahydrobenzo[h]quinoline crystalline powder

Chemical for dry chemistry diagnostic tests

Application

Use 3-Hydroxy-1,2,3,4-tetrahydrobenzo[h]quinoline in dry chemistry application for the determination of urea.

Benefits

Rely on the proven diagnostic quality of this reagent.

CAS: 5423-67-6

Properties

Nomenclature: 3-Hydroxy-1,2,3,4-tetrahydrobenzo[h]quinoline

Formula: C₁₃H₁₃NO

Molecular weight: 199.25

Toxicity: Harmful

Specification

Appearance: White to greyish, odorless crystalline powder

Solubility: Clear, colorless solution in methanol (c=0.2%; w/v)

Melting range: +150 to +154°C

Loss on drying: ≤0.5%

Sulphate ash: ≤0.2%

3-Hydroxy-1,2,3,4-tetrahydrobenzo[h]quinoline (HClO₄ titration, based on dried substance): 98.0-102.0%

Cat. No. **Pack Size**

10 003 174 001 custom fill

Will be supplied as "BM 32.027". Unit of Measure is "kg".

For further processing only.

UV / VIS spektrum:

Maximum I: 250 to 254 nm (specific absorbance ($A_{1\% / 1\text{cm}}$): 988-1050)

Maximum II: 334 to 338 nm (specific absorbance ($A_{1\% / 1\text{cm}}$): 286-306)

Thin layer chromatography (TLC): Corresponds to reference

IR Spectrum: Corresponds to reference

Stability: At +2 to +8°C within specification range for 24 months. Protect from light.

D-Mannitol reduced sodium

Excipient (inactive substance) for the production of tablets or granulated material

Application

Use D-Mannitol as an excipient (inactive substance) to produce tablets or granulated material that contain reagents or components used in diagnostic applications.

Benefits

- Rely on the proven diagnostic quality of this product.

CAS: 69-65-8

Properties

Formula: $C_6H_{14}O_6$

Molecular weight: 182.2 D

Solubility: Slightly soluble in water and hot ethanol, of low solubility in ethanol.

Specification

Appearance: White, silky crystals or white, crystalline powder

Solubility: Clear, colorless solution in water (c=10%, w/v)

Melting range: +166°C to +168°C

Specific rotation (c=10%, w/v, calculated on dry substance): +23.0° to +25.0°

Purity (HPLC):

D-Mannit: ≥97.5 area%

Sorbit: ≤2.5 area%

Alkaline impurities (calculated as NaOH): ≤80 ppm

Acid impurities (calculated as HCl): ≤45 ppm

As: ≤1 ppm

Cl: ≤50 ppm

Fe: ≤0.1 ppm

Ni: ≤1 ppm

Reducing sugar: Corresponds to reference

Heavy metals (as Pb): ≤2 ppm

Sulphate: ≤100 ppm

Sulphate ash: ≤500 ppm

Loss on drying: ≤0.3%

Na: ≤200 ppm

Microbiological analysis: Corresponds to reference

Stability: At +15 to +25°C within specification range for 36 months. Store dry. Protect from light.

Cat. No. **Pack Size**

11 371 754 103 custom fill

Will be supplied as "D-Mannit, Na-arm". Unit of Measure is "kg".
For further processing only.

Kryptofix 221 solution

Cryptant that binds cations in aqueous solutions

Application

Use Kryptofix 221 in enzymatic tests for potassium to decrease the sodium concentration relative to potassium.

Benefits

- Selectively reduce sodium in your diagnostic reagent.
- Rely on the strongly reduced concentration of potassium and sodium.

CAS: 31364-42-8

Properties

Formula: C₁₆H₃₂N₂O₅

Molecular weight: 332.44 D

Specification

Appearance: Clear, yellow liquid

Solubility: Clear yellowish solution in water (c=4 mg/ml)

A₃₄₀ (aqueous solution): ≤0.100

A₄₀₅ (aqueous solution): ≤0.060

Identity (NIR): Corresponds to reference

¹³C-NMR spectrum: Corresponds to masterlot

Kryptofix 221 agent (HPLC; based on masterlot): ≥75 %

Na (AES): ≤65 ppm

K (AES): ≤5 ppm

NH₃ (evolution in buffer; after 10 days at +55°C): ≤50 μmol/l

Stability: At +2 to +8°C within specification range for 12 months. Keep under nitrogen. Protect from light.

Cat. No.	Pack Size
11 183 958 103	custom fill

Will be supplied as "Kryptofix 221". Unit of Measure is "kg active ingredient".

For further processing only.

Valinomycin crystallizate

Potassium selective ionophoric cyclodepsipeptide

Application

Use Valinomycin in diagnostic tests for potassium where it acts as an ion carrier in potassium selective electrodes.

Benefits

Rely on the proven diagnostic quality of this reagent.

CAS: 2001-95-8

Properties

Formula: C₅₄H₉₀N₆O₁₈

Molecular weight: 1111.4 D

Specification

Appearance: White crystallizate

Solubility: Clear, colorless solution in chloroform (c=10 mg/ml)

Melting point: ≥+183°C

Specific rotation (in chloroform):: +30.0±2.0°

Valinomycin (from N): ≥94%

Valinomycin (HPLC): ≥85.0 area%

C (elementary analysis) : ≥54.86%

H (elementary analysis) : >=7.67%

Cat. No.	Pack Size
10 161 594 103	custom fill

Will be supplied as "Valinomycin". Unit of Measure is "g".

For further processing only.

N (elementary analysis): $\geq 7.10\%$

Thin layer chromatography (TLC):

a) UV: Homogeneous

b) to spray with H_2SO_4 (1%); to vaporize with Iodine: Corresponds to reference

Stability: At +2 to +8°C within specification range for 36 months. Store in safety zone dedicated to poisonous agents.

Acetyl-Coenzyme A trisodium salt

Cofactor for carnitine acetyl transferase

Application

Use Acetyl-Coenzyme A for the determination of L-carnitine.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 72-89-9

Properties

Formula: $C_{23}H_{35}N_7O_{17}P_3SNa_3$

Molecular weight: 875.5 D (Acetyl-CoA: 809.6 D)

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

Acetyl-CoA (enzymatically): 78-97%

Acetyl-CoA (A_{260} , $\epsilon=16.0$ [$l \times mmol^{-1} \times cm^{-1}$]): 80-97%

Na (flame photometric): 6.5-7.5%

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

12 207 273 103 custom fill

Will be supplied as "Acetyl-CoA, Tri-Na". Unit of Measure is "g".



For further processing only.

Acetyl-Coenzyme A trilithium salt

Cofactor for carnitine acetyl transferase

Application

Use Acetyl-Coenzyme A for the determination of L-carnitine.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 72-89-9

Properties

Formula: $C_{23}H_{35}N_7O_{17}P_3SLi_3$

Molecular weight: 827.4 D (Acetyl-CoA: 809.6 D)

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

Acetyl-CoA (enzymatically): $\geq 83\%$

Acetyl-CoA (A_{260} , $\epsilon=16.0$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 85\%$

Li (flame photometric): $2 \pm 0.3\%$

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 150 932 103 custom fill

Will be supplied as "Acetyl Coenzyme A (Acetyl-CoA) Tri-Li". Unit of Measure is "g".



For further processing only.

Coenzyme A, Grade I free acid

Coenzyme A is a cofactor for some enzymes, *e.g.*, citrate lyase.

Application

Use Coenzyme A in diagnostic tests measuring citrate or citrate lyase. It is also used as an enhancer of luciferase light emission or as a precursor for Acetyl-Coenzyme A or in other chemical or enzymatical reactions.

Benefits

- Rely on the enhanced purity of this Grade I product.

CAS: 85-61-0

Properties

Formula: $C_{21}H_{36}N_7O_{16}P_3S$

Molecular weight: 767.6 D (CoA: 767.6 D)

Specification

Appearance: White to slightly yellow lyophilizate

CoA, reduced (enzymatically, 10 U phosphotransacetylase): $\geq 85\%$

CoA (A_{260}^1 , $\epsilon=16.0$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 88\%$

Water (K. Fischer): $\leq 6\%$

Glutathione, reduced (enzymatically): $\leq 1\%$

Stability: At -15 to $-25^\circ C$ within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 151 009 103 custom fill

Will be supplied as "Coenzyme A (CoA), Free Acid, Grade I". Unit of Measure is "g".



For further processing only.

Coenzyme A, Grade I trilithium salt

Coenzyme A is a cofactor for some enzymes, *e.g.*, citrate lyase.

Application

Use Coenzyme A in diagnostic tests measuring citrate or citrate lyase. It is also used as an enhancer of luciferase light emission or as a precursor for Acetyl-Coenzyme A or in other chemical or enzymatical reactions.

Benefits

- Rely on the enhanced purity of this Grade I product.

CAS: 85-61-0

Properties

Formula: $C_{21}H_{33}N_7O_{16}P_3SLi_3$

Molecular weight: 785.4 D (CoA: 767.6 D)

Specification

Appearance: White to slightly yellow lyophilizate

CoA, reduced (enzymatically with 10 U phosphotransacetylase): $\geq 83\%$

CoA (A_{260}^1 , $\epsilon=16.0$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 84\%$

Water (K. Fischer): $\leq 6\%$

Glutathione, reduced (enzymatically): $\leq 1\%$

Stability: At -15 to $-20^\circ C$ within specification range for 12 months.

Cat. No. **Pack Size**

10 121 541 103 custom fill

Will be supplied as "Coenzyme A (CoA), Tri-Li Salt, Grade I". Unit of Measure is "g".



For further processing only.

Coenzyme A, Grade II trilithium salt

Coenzyme A is a cofactor for some enzymes, *e.g.*, citrate lyase.

Application

Use Coenzyme A in diagnostic tests measuring citrate or citrate lyase. It is also used as an enhancer of luciferase light emission or as a precursor for Acetyl-Coenzyme A or in other chemical or enzymatical reactions.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 85-61-0

Properties

Formula: $C_{21}H_{33}N_7O_{16}P_3SLi_3$

Molecular weight: 785.4 D (CoA: 767.6 D)

Specification

Appearance: White to slightly yellow lyophilizate

CoA, reduced (enzymatically with 10 U phosphotransacetylase): $\geq 73\%$

CoA (A_{260} , $\epsilon = 16.0$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 81\%$

Water (K. Fischer): $\leq 8\%$

Stability: At -15 to $-20^\circ C$ within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 155 969 103 custom fill

Will be supplied as "Coenzyme A (CoA), Ttri-Li, Grade II". Unit of Measure is "g".



For further processing only.

FAD disodium salt

Cofactor for dehydrogenases and oxidases.

Application

Use FAD as a cofactor in a variety of enzymatic test or assays to activate enzymes, especially diagnostic tests for triglycerides.

Benefits

- Rely on the proven diagnostic quality of this product.
- Use FAD to design your diagnostic test

CAS: 146-14-5

Properties

Nomenclature: Flavine-adenine dinucleotide

Formula: $C_{27}H_{31}N_9O_{15}P_2Na_2$

Molecular weight: 829.6 D (FAD: 785.7 D)

Specification

Appearance: Yellow powder

FAD (A_{450} , $\epsilon = 11.3$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 86\%$

Na (flame photometric): $5 \pm 1\%$

Water (K. Fischer): $\leq 9\%$

P: $\leq 0.6\%$

Stability: At $+2$ to $+8^\circ C$ within specification range for 24 months. Protect from light.

Cat. No. **Pack Size**

10 154 032 103 custom fill

Will be supplied as "Flavine-adenine Dinucleotide (FAD), Di-Na". Unit of Measure is "g".

For further processing only.

NAD, Grade I free acid

Cofactor for dehydrogenases, *e.g.*, lactate dehydrogenase.

Application

Use NAD, Grade I as a cofactor in a variety of diagnostic tests, such as for the determination of ethanol and lactate dehydrogenase.

Benefits

Rely on the enhanced purity.

CAS: 53-84-9

Properties

Formula: $C_{21}H_{27}N_7O_{14}P_2$

Molecular weight: 663.4 D

Specification

Appearance: Colorless to slightly yellowish lyophilizate

Solubility: Clear, colorless to slightly yellowish solution in water (c=200 mg/ml)

β -NAD (from value found enzymatically, based on dry weight): $\geq 99\%$

β -NAD (enzymatically, A_{340}): $\geq 96.5\%$

β -NAD (A_{260} , $\epsilon=17.6$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 96.5\%$

NAD (HPLC): ≥ 98 area%

Water (K. Fischer): $\leq 3.5\%$

Fe (AA): ≤ 50 ppm

AMP (enzymatically): $\leq 0.1\%$

Ethanol (GC): ≤ 43 ppm

Aceton, isopropanol, methanol (GC): each $\leq 0.05\%$

Reaction rates (LDH) based on NAD II, acid: 95-105%

A_{250}/A_{260} : 0.81-0.85

A_{280}/A_{260} : 0.20-0.24

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 004 618 103 custom fill

Will be supplied as "b-Nicotinamide-adenine Dinucleotide, I". Unit of Measure is "kg".

For further processing only.

NAD, Grade II free acid

Cofactor for dehydrogenases, *e.g.*, lactate dehydrogenase.

Application

Use NAD, Grade II as a cofactor in a variety of diagnostic tests, such as for the determination of ethanol and lactate dehydrogenase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 53-84-9

Properties

Formula: $C_{21}H_{27}N_7O_{14}P_2$

Molecular weight: 663.4 D

Specification

Appearance: Colorless to slightly yellowish lyophilizate

Solubility: Clear, colorless to slightly yellowish solution in water (c=200 mg/ml)

β -NAD (from value found enzymatically, based on dry weight): $\geq 97.5\%$

β -NAD (enzymatically, A_{340}): $\geq 94.5\%$

Cat. No. **Pack Size**

10 004 626 103 custom fill

Will be supplied as "b-Nicotinamide-adenine Dinucleotide, II". Unit of Measure is "kg".

For further processing only.

β -NAD (A_{260} , $\epsilon=17.6$ [$l \times \text{mmol}^{-1} \times \text{cm}^{-1}$]): $\geq 94.5\%$

NAD (HPLC): ≥ 95 area%

Water (K. Fischer): $\leq 3.5\%$

Fe (AA): ≤ 25 ppm

AMP (enzymatically): $\leq 0.1\%$

Ethanol (GC): ≤ 40 ppm

Aceton, isopropanol, methanol (GC): $\leq 0.1\%$, $\leq 0.15\%$, $\leq 0.15\%$

Reaction rates (LDH) based on NAD II, acid: 95-105%

A_{250}/A_{260} : 0.81-0.85

A_{280}/A_{260} : 0.20-0.24

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

NADH, Grade I disodium salt

Cofactor for a variety of dehydrogenases, e.g., malate dehydrogenase and lactate dehydrogenase.

Application

Use NADH, Grade I, as a cofactor in a variety of diagnostic tests, such as for glutamate dehydrogenase, lactate dehydrogenase, α -hydroxybutyrate dehydrogenase, aminotransferases and urea.

Benefits

Rely on the enhanced purity of this Grade I product.

CAS: 58-68-4

Properties

Formula: $C_{21}H_{27}N_7O_{14}P_2Na_2$

Molecular weight: 709.4 D (NADH: 665.4 D)

Specification

Appearance: White to slightly yellowish amorphous powder

Solubility: Clear, colorless to slightly yellowish solution in water (c=50 mg/ml)

NADH- Na_2 (calculated from value found enzymatically, based on dry weight): $\geq 99\%$

NADH (enzymatically, A_{340}): $\geq 85\%$

NADH (A_{340} , $\epsilon=6.3$ [$l \times \text{mmol}^{-1} \times \text{cm}^{-1}$]): $\geq 85\%$

NADH (A_{260} , $\epsilon=14.3$ [$l \times \text{mmol}^{-1} \times \text{cm}^{-1}$]): $\geq 85\%$

Na (flame photometric): $6.5 \pm 0.5\%$

Water (K. Fischer): $\leq 5\%$

NAD (enzymatically): $\leq 0.5\%$

AMP (enzymatically): $\leq 0.2\%$

Ethanol (GC): $\leq 4\%$

Reaction rates (LDH) based on standard: 95-105%

A_{260}/A_{340} : ≤ 2.35

Stability: At +2 to +8°C within specification range for 12 months. Keep under nitrogen. Protect from light.

Cat. No.	Pack Size
10 004 634 103	250g or custom fill

Will be supplied as "b-NADH, Reduced, Disodium Salt, Grade I".
Unit of Measure is "kg".
For further processing only.

NADH, Grade II disodium salt

Cofactor for a variety of dehydrogenases, *e.g.*, malate dehydrogenase and lactate dehydrogenase.

Application

Use NADH, Grade II, as a cofactor in a variety of diagnostic tests, such as for glutamate dehydrogenase, lactate dehydrogenase, α -hydroxybutyrate dehydrogenase, aminotransferases and urea.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 58-68-4

Properties

Formula: $C_{21}H_{27}N_7O_{14}P_2Na_2$

Molecular weight: 709.4 D (NADH: 665.4 D)

Specification

Appearance: White to slightly yellowish amorphous powder

Solubility: Clear, colorless to slightly yellowish solution in water (c=50 mg/ml)

NADH- Na_2 (calculated from value found enzymatically, based on dry weight): $\geq 98\%$

NADH (enzymatically, A_{340}): $\geq 84\%$

NADH (A_{340} , $\epsilon=6.3$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 84\%$

NADH (A_{260} , $\epsilon=14.3$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 85\%$

Na (flame photometric): $6.5 \pm 0.5\%$

Water (K. Fischer): $\leq 6\%$

NAD (enzymatically): $\leq 1\%$

AMP (enzymatically): $\leq 0.2\%$

Ethanol (GC): $\leq 4\%$

Reaction rates (LDH) based on standard: 95-105%

A_{260}/A_{340} : ≤ 2.40

A_{260}/A_{240} : 1.57-2.17

Stability: At +2 to +8°C within specification range for 12 months. Keep under nitrogen. Protect from light.

Cat. No.

Pack Size

10 004 642 103 250g or custom fill

Will be supplied as "b-NADH, Reduced, Disodium Salt, Grade II".
Unit of Measure is "kg".

For further processing only.

NADH, Grade II for potassium test, disodium salt

NADH quality for enzymatic potassium test.

Application

Use this special NADH as a cofactor in a diagnostic test for potassium, together with glutamate dehydrogenase to remove ammonia from the reaction mixture.

Benefits

- Rely on the proven diagnostic quality of this product.
- Rely on the strongly reduced concentration of potassium.

CAS: 58-68-4

Properties

Formula: C₂₁H₂₇N₇O₁₄P₂Na₂

Molecular weight: 709.4 D (NADH: 665.4 D)

Specification

Appearance: White to slightly yellowish amorphous powder

Solubility: Clear, colorless to yellowish solution in water (c=50 mg/ml)

NADH-Na₂ (from content found enzymatically, based on dry weight): ≥98%

NADH (enzymatically): ≥82%

NADH (A₃₄₀, ε=6.3 [l x mmol⁻¹ x cm⁻¹]): ≥84%

NADH (A₂₆₀, ε=14.3 [l x mmol⁻¹ x cm⁻¹]): ≥85%

Na (flame photometric): 6.5±0.5%

K (AAS): ≤250 ppm

Water (K. Fischer): ≤6%

NAD (enzymatically): ≤1%

AMP (enzymatically): ≤0.2%

Ethanol (GC): ≤4%

Reaction rates (Lactate dehydrogenase), based on standard:

Freshly dissolved sample: 95-105%

After 8 days incubation at +45°C: ≥95%

HPLC:

Freshly dissolved sample: ≥96.0 area%

After 8 days incubation at +45 °C: ≥75%

A₂₆₀/A₃₄₀: ≤2.4

A₂₆₀/A₂₄₀: 1.57-2.17

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No. **Pack Size**

11 333 925 103 custom fill

Will be supplied as "NADH Di-Na, Grade II for Potassium Test".
Unit of Measure is "kg".



For further processing only.

NADP disodium salt

Cofactor for dehydrogenases, *e.g.*, glucose-6-phosphate dehydrogenase.

Application

Use NADP as a cofactor in a variety of diagnostic tests, such as for the determination of glucose and creatine kinase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 53-59-8

Properties

Formula: $C_{21}H_{26}N_7O_{17}P_3Na_2$

Molecular weight: 787.4 D (NADP: 743.4 D)

Specification

Appearance: White to yellowish amorphous powder

Solubility: Clear, colorless solution in water (c=50 mg/ml)

NADP, Na-salt (from value found enzymatically, based on dry weight): $\geq 97\%$

NADP (enzymatically, A_{340} , $\epsilon=6.3$ [l x mmol⁻¹ x cm⁻¹]): $\geq 85\%$

NADP (A_{260} , $\epsilon=18$ [l x mmol⁻¹ x cm⁻¹]): $\geq 85\%$

Na (flame photometric): $4.5 \pm 0.5\%$

Water (K. Fischer): $\leq 6\%$

NAD (enzymatically): $\leq 0.5\%$

Methanol (GC): $\leq 3\%$

A₃₆₀ (c=10 mg/ml water, against water): ≤ 0.600

A₃₄₀ (c=0.01 mg/ml phosphate buffer, pH 7.0): ≤ 0.005

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No. **Pack Size**

10 004 669 103 custom fill

Will be supplied as "b-NADP, Disodium Salt". Unit of Measure is "kg".

For further processing only.

NADP monopotassium salt

Cofactor for dehydrogenases, *e.g.*, glucose-6-phosphate dehydrogenase.

Application

Use NADP as a cofactor in a variety of diagnostic tests, such as for the determination of glucose and creatine kinase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 53-59-8

Properties

Formula: $C_{21}H_{27}N_7O_{17}P_3K \times 2 H_2O$

Molecular weight: 817.4 D (NADP: 743.4 D)

Specification

Appearance: White crystalline powder

NADP, K-salt (calculated from value determined enzymatically, based on dry weight): $\geq 97\%$

NADP (enzymatically, A_{340} , $\epsilon=6.3$ [l x mmol⁻¹ x cm⁻¹]): $\geq 88\%$

NADP (A_{260} , $\epsilon=18$ [l x mmol⁻¹ x cm⁻¹]): $\geq 88\%$

K (flame photometric): 4.0-5.0%

Water (K. Fischer): $\leq 4.5 \pm 1.0\%$

NAD (enzymatically): $\leq 0.2\%$

Methanol (GC): $\leq 2\%$

Mg (AAS): ≤ 40 ppm

A₃₆₀ (c=10 mg/ml water, against water): ≤ 0.600

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. **Pack Size**

10 233 536 103 custom fill

Will be supplied as "b-NADP, Monopotassium Salt". Unit of Measure is "g".



For further processing only.

NADPH

tetrasodium salt

Cofactor for glutamate dehydrogenase.

Application

Use NADPH as a cofactor in diagnostic tests for ammonia, urea and creatinine.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 53-57-6

Properties

Formula: $C_{21}H_{26}N_7O_{17}P_3Na_4$

Molecular weight: 833.4 D (NADPH: 745.4 D)

Specification

Appearance: White to slightly yellowish amorphous powder

Particle size (screen analysis according to Ph.Eur., US mesh <30): $\geq 95\%$

Solubility: Clear, colorless to slightly yellowish solution in water (c=50 mg/ml)

NADPH- Na_4 (calculated from content found enzymatically, based on dry weight): $\geq 97\%$

NADPH (enzymatically, A_{340}): $\geq 79\%$

NADPH (A_{340} , $\epsilon=6.3$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 0.79\%$

NADPH (A_{260} , $\epsilon=15$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 80\%$

NADPH (HPLC): ≥ 95 area%

Na (flame photometric): $11 \pm 1\%$

Water (K. Fischer): $\leq 6\%$

NADH (HPLC): ≤ 0.5 area%

NADP (HPLC): ≤ 0.5 area%

NADP (enzymatically): $\leq 0.5\%$

Nicotinic acid amide: ≤ 2 area%

Ethanol (GC): $\leq 3\%$

A_{260}/A_{340} : 2.32-2.65

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 041 939 103 custom fill

Will be supplied as "b-NADPH, Reduced, Tetrasodium Salt". Unit of Measure is "g".

For further processing only.

Thio-NAD

free acid

NAD analog with an absorbance maximum of thio-NADH at 405 nm.

Application

Use Thio-NAD instead of NAD in enzymatic reactions to measure the reaction kinetics at 405 nm, e.g., kinetic enzyme cycling methods using Thio-NAD and NADH.

Benefits

- Use Thio-NAD in your sensitive enzymatic diagnostic test.
- Take advantage of the increased absorption maximum.

CAS: 4090-29-3

Properties

Formula: $C_{21}H_{27}N_7O_{13}SP_2$

Molecular weight: 679.5 D

Specification

Appearance: Yellowish powder

Solubility: Clear, colorless to slightly yellowish solution in water (c=37 mg/ml)

Cat. No. **Pack Size**

04 635 396 103 custom fill

Will be supplied as "Thio-NAD free acid". Unit of Measure is "g".



For further processing only.

Thio-NAD ($A_{259} \epsilon = 19.7 \text{ l x mmol}^{-1} \text{ x cm}^{-1}$): $\geq 95\%$

Water (K. Fischer): $\leq 4\%$

A_{398}/A_{340} (against water) : ≤ 0.124

A_{236}/A_{259} (against water) : 0.640-0.670

A_{296}/A_{259} (against water) : 0.244-0.264

Revision of absorption in presence of 0.4 mmol/l NADH: $\leq 0.001 \text{ abs/min}$

Ca (AAS): No limit

Mg (AAS): No limit

Stability: At -15 to -25°C within specification range for 12 months.

Adenosine-5'-O-(2-thiodiphosphate) trilithium salt

Non-hydrolyzable ADP analog.

Application

Use Adenosine-5'-O-(2-thiodiphosphate) to inhibit ADP-binding enzymes.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 73536-95-5

Properties

Formula: C₁₀H₁₂N₅O₉P₂SLi₃

Molecular weight: 461.0 D (ATP-β-S: 443.2 D)

Specification

Appearance: White powder

ATP-β-S, Li (A₂₆₀): ≥81%

ATP-β-S (A₂₆₀, ε= 15.0 [l x mmol⁻¹ x cm⁻¹]): ≥78%

ATP-β-S (HPLC): ≥90 area%

Li (flame photometric): 3-5%

Water (K. Fischer): ≤12%

ATP (HPLC): ≤5 area%

ADP (HPLC): ≤1 area%

AMP (HPLC): ≤4 area%

A₂₅₀/A₂₆₀: 0.75-0.83

A₂₈₀/A₂₆₀: 0.14-0.18

A₂₉₀/A₂₆₀: 0.00-0.01

Stability: At -15 to -25°C within specification range for 12 months. Protect from light.

Cat. No. **Pack Size**

10 200 166 103 custom fill

Will be supplied as "Adenosine-5'-O-(2-thiodiphosphate) Tri-Li".
Unit of Measure is "g".



For further processing only.

Adenosine-5'-O-(3-thiotriphosphate) tetralithium salt

Non-hydrolyzable ATP analog.

Application

Use Adenosine-5'-O-(3-thiotriphosphate) to inhibit ATP-binding enzymes.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 35094-45-2

Properties

Formula: C₁₀H₁₂N₅O₁₂P₃SLi₄

Molecular weight: 547.0 D (ATP-γ-S: 523.2 g/mol)

Specification

Appearance: White powder

ATP-γ-S, Li (A₂₆₀): ≥78%

ATP-γ-S (A₂₆₀, ε=15.0 [l x mmol⁻¹ x cm⁻¹]): ≥74%

ATP-γ-S (HPLC): ≥85 area%

Li (flame photometric): 3-5%

Water (K. Fischer): ≤12%

ADP (HPLC): ≤12 area%

AMP (HPLC): ≤3 area%

A₂₅₀/A₂₆₀: 0.79±0.02

Cat. No. **Pack Size**

10 122 734 103 custom fill

Will be supplied as "Adenosine-5'-O-(3-thiotriphosphate), Li₄".
Unit of Measure is "g".



For further processing only.

A_{280}/A_{260} : 0.16±0.01

A_{290}/A_{260} : ≤0.05

Stability: At -15 to -25°C within specification range for 6 months. Protect from light.

ADP potassium salt

Cofactor for diagnostic tests.

Application

Use ADP in variety of diagnostic tests, such as for the determination of creatine kinase and pyruvate kinase. Use it also for the activation of glutamate dehydrogenase in the determination of for example urea or ammonia.

Benefits

- Enhance the activation of your glutamate dehydrogenase containing reagents and rely on the proven diagnostic quality of this product.

CAS: 58-64-0

Properties

Formula: $C_{10}H_{14}N_5O_{10}P_2K \times 2 H_2O$

Molecular weight: 501.3 D (ADP: 427.2 D)

Remark: Crystalline ADP-K x 2 H₂O is the purest and most stable form of ADP available.

Specification

Appearance: Colorless crystals

Solubility: Clear, colorless solution in water (c=50 mg/ml)

ADP-K x 2 H₂O (based on value found enzymatically): ≥98%

ADP (enzymatically): ≥84%

ADP (A_{260} , $\epsilon = 15 [l \times mmol^{-1} \times cm^{-1}]$): ≥84%

K (flame photometric): 7.8±0.5%

Water (K. Fischer): 7.2±1%

P_i (Fiske and Subbarow): ≤0.3%

AMP (enzymatically): ≤1%

ATP (enzymatically): ≤0.2%

NH₄ (enzymatically): ≤0.005%

A_{250}/A_{260} : 0.78±0.02

A_{280}/A_{260} : 0.16±0.01

A_{290}/A_{260} : ≤0.01

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No.	Pack Size
10 233 528 103	custom fill

Will be supplied as "Adenosine-5'-diphosphate (ADP), K-Salt".
Unit of Measure is "kg".
For further processing only.

ADP disodium salt

Cofactor for diagnostic tests.

Application

Use ADP in variety of diagnostic tests, such as for the determination of creatine kinase and pyruvate kinase. Use it also for the activation of glutamate dehydrogenase in the determination of for example urea or ammonia.

Benefits

- Enhance the activation of your glutamate dehydrogenase containing reagents
- Rely on the proven diagnostic quality of this product.

Cat. No.	Pack Size
10 129 062 103	custom fill

Will be supplied as "ADP, Di-Na". Unit of Measure is "kg".
For further processing only.

CAS: 58-64-0

Properties

Formula: $C_{10}H_{13}N_5O_{10}P_2Na_2$

Molecular weight: 471.2 D (ADP: 427.2 D)

Remark: ATP and AMP may form during storage.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=50 mg/ml)

ADP-Na-salt (calculated on value found enzymatically): $\geq 90\%$

ADP (enzymatically): $\geq 82\%$

ADP (A_{260} , $\epsilon = 15$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 82\%$

Na (flame photometric): $9 \pm 1\%$

Water (K. Fischer): $\leq 7\%$

P_i (Fiske and Subbarow): $\leq 0.6\%$

AMP (enzymatically): $\leq 3\%$

ATP (enzymatically): $\leq 1\%$

NH₄ (enzymatically): $\leq 0.01\%$

A₂₅₀/A₂₆₀: 0.78 ± 0.02

A₂₈₀/A₂₆₀: 0.16 ± 0.01

A₂₉₀/A₂₆₀: ≤ 0.01

Stability: At -15 to -25°C within specification range for 6 months. Store dry.

ADP for potassium test free acid

Application

Use ADP for potassium testing as an activator for glutamate dehydrogenase in enzymatic potassium tests.

Benefits

- Take advantage of the strongly reduced potassium concentration.
- Enhance your enzymatic potassium test
- Rely on the proven diagnostic quality of this product.

CAS: 58-64-0

Properties

Formula: $C_{10}H_{15}N_5O_{10}P_2$

Molecular weight: 427.2 D

Remark: ATP and AMP may form during storage.

Specification

Appearance: White crystallizate

Solubility: Clear, colorless solution in NaOH, 0.1 mol/l (c=50 mg/ml)

IR-spectrum (KBr-pellet): Corresponds to reference

ADP (enzymatically): $\geq 97\%$

ADP (A_{260} , $\epsilon = 15$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 97\%$

Water (K. Fischer): $\leq 2\%$

P_i (Fiske and Subbarow): $\leq 0.6\%$

AMP (enzymatically): $\leq 3\%$

ATP (enzymatically, HK/G6P-DH): $\leq 0.3\%$

Na (AAS): ≤ 750 ppm

K (AAS): ≤ 20 ppm

NH₄ (enzymatically): $\leq 0.01\%$

A₂₅₀/A₂₆₀: 0.78 ± 0.02

A₂₈₀/A₂₆₀: 0.16 ± 0.01

A₂₉₀/A₂₆₀: ≤ 0.01

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

11 333 879 103 custom fill

Will be supplied as "ADP Free Acid for Potassium Test". Unit of Measure is "kg".



For further processing only.

AMP free acid

Nucleotide for diagnostic tests.

Application

Use AMP for the determination of 5'-nucleotidase and in diagnostic tests for the determination of creatine kinase.

Benefits

- Enhance the performance of your creatine kinase test.
- Rely on the proven diagnostic quality of this product.

CAS: 61-19-8

Properties

Formula: $C_{10}H_{14}N_5O_7P \times H_2O$

Molecular weight: 365.2 D (AMP: 347.2 D)

Specification

Appearance: White, crystalline powder

Solubility: Clear, colorless solution in NaOH, 1 mol/l (c=50 mg/ml)

AMP x H₂O (based on value found enzymatically): ≥98%

AMP (enzymatically): ≥93%

AMP (A_{260} , $\epsilon=15.0$ [$l \times mmol^{-1} \times cm^{-1}$]): ≥93%

Water (K. Fischer): ≤5±2%

P_i: ≤0.3%

Fe (bathophenanthroline): ≤10 ppm

Heavy metals (as Pb): ≤10 ppm

A_{250}/A_{260} : 0.78±0.02

A_{280}/A_{260} : 0.15±0.01

A_{290}/A_{260} : ≤0.01

Stability: At +15 to +25°C within specification range for 36 months. Store dry.

Cat. No. **Pack Size**

10 000 086 103 custom fill

Will be supplied as "Adenosine-5'-monophosphoric Acid (AMP)".
Unit of Measure is "kg".

For further processing only.

AMP disodium salt

Nucleotide for diagnostic tests.

Application

Use AMP for the determination of 5'-nucleotidase and in diagnostic tests for the determination of creatine kinase.

Benefits

- Enhance the performance of your creatine kinase test
- Rely on the proven diagnostic quality of this product.

CAS: 61-19-8

Properties

Formula: $C_{10}H_{12}N_5O_7PNa_2$

Molecular weight: 391.2 D (AMP: 347.2 D)

Specification

Appearance: White crystals

Solubility: Clear, colorless solution in water (c=50 mg/ml)

AMP (enzymatically): $\geq 76\%$

AMP (A_{260}): $\geq 76\%$

Na (flame photometric): 10-12%

Water (K. Fischer): $\leq 12\%$

P_i: $\leq 0.3\%$

A₂₅₀/A₂₆₀: 0.78 ± 0.02

A₂₈₀/A₂₆₀: 0.16 ± 0.01

A₂₉₀/A₂₆₀: ≤ 0.01

Stability: At +15 to +25°C within specification range for 36 months. Store dry.

Cat. No. **Pack Size**

10 000 094 103 custom fill

Will be supplied as "Adenosine-5'-monophosphate (AMP), Di-Na". Unit of Measure is "kg".

For further processing only.

ATP, Grade I disodium salt

Cofactor for kinases, e.g. glycerokinase and hexokinase.

Application

Use ATP, Grade I, in a variety of diagnostic tests, such as for the determination of triglycerides, creatine kinase and glucose.

Benefits

Rely on the enhanced purity of this grade I product.

CAS: 56-65-5

Properties

Formula: $C_{10}H_{14}N_5O_{13}P_3Na_2 \times 3 H_2O$

Molecular weight: 605.2 D (ATP: 507.2 D)

Specification

Appearance: White crystals

Solubility: Clear, colorless solution in water (c=50 mg/ml)

ATP-Na₂ x 3 H₂O (based on value found enzymatically): $\geq 99\%$

ATP (enzymatically): $\geq 84\%$

ATP (A_{260} , $\epsilon=15.0$ [l x mmol⁻¹ x cm⁻¹]): $\geq 84\%$

Na (flame photometric): $7.5 \pm 0.5\%$

Water (K. Fischer): $\leq 8 \pm 1\%$

P_i: $\leq 0.15\%$

ADP, AMP (enzymatically): $\leq 0.5\%$

Cat. No. **Pack Size**

10 422 495 103 custom fill

Will be supplied as "ATP, Di-Na, Special Quality". Unit of Measure is "kg".

For further processing only.

GTP (HPLC): ≤ 0.01 area%

Fe (AAS): ≤ 10 ppm

Mg (AAS): ≤ 10 ppm

Ca (AAS): ≤ 20 ppm

Zn (AAS): ≤ 5 ppm

V (AAS): ≤ 1 ppm

A_{250}/A_{260} : 0.79 ± 0.02

A_{280}/A_{260} : 0.15 ± 0.01

A_{290}/A_{260} : ≤ 0.01

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

ATP, Grade II disodium salt

Cofactor for kinases, *e.g.* glycerokinase and hexokinase.

Application

Use ATP, Grade II in a variety of diagnostic tests, such as for the determination of triglycerides and glucose.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 56-65-5

Properties

Formula: $C_{10}H_{14}N_5O_{13}P_3Na_2 \times 3 H_2O$

Molecular weight: 605.2 D (ATP: 507.2 D)

Specification

Appearance: White crystals

Solubility: Clear, colorless solution in water (c=50 mg/ml)

ATP-Na₂ x 3 H₂O (based on value found enzymatically): $\geq 98\%$

ATP (enzymatically): $\geq 82\%$

ATP (A_{260} , $\epsilon = 15.0$ [$l \times mmol^{-1} \times cm^{-1}$]): $\geq 82\%$

Na (flame photometric): $7.5 \pm 0.5\%$

Water (K. Fischer): $\leq 10\%$

P_i: $\leq 0.3\%$

ADP, AMP (enzymatically): $\leq 0.5\%$

GTP (HPLC): ≤ 0.01 area%

Fe (AAS): ≤ 15 ppm

Heavy metals (as Pb): ≤ 30 ppm

A_{250}/A_{260} : 0.79 ± 0.02

A_{280}/A_{260} : 0.15 ± 0.01

A_{290}/A_{260} : ≤ 0.01

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No. **Pack Size**

10 000 116 103 custom fill

Will be supplied as "Adenosine-5'-triphosphate (ATP), Di-Na". Unit of Measure is "kg".

For further processing only.

Guanosine-5'-O-(3-thiodiphosphate) tetralithium salt

Non-hydrolyzable GTP analog.

Application

Use Guanosine-5'-O-(3-thiodiphosphate) to inhibit GTP-binding enzymes and activate G-proteins.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 37589-80-3

Properties

Formula: C₁₀H₁₂N₅O₁₃P₃SLi₄

Molecular weight: 563.0 D (GDP-γ-S: 539.2 D)

Specification

Appearance: White powder

GTP-γ-S-Li (A₂₅₄): ≥77%

GTP-γ-S (A₂₅₄, ε=13.5 [l x mmol⁻¹ x cm⁻¹]): ≥73%

GDP-γ-S (HPLC): ≥85 area%

Li (flame photometric): 5±1%

Water (K. Fischer): ≤12%

GDP (HPLC): ≤12 area%

GMP (HPLC): ≤1 area%

GTP (HPLC): ≤2 area%

A₂₅₀/A₂₆₀: 1.14±0.05

A₂₈₀/A₂₆₀: 0.65±0.04

A₂₉₀/A₂₆₀: 0.27±0.03

Stability: At -15 to -25°C within specification range for 6 months. Store dry.

Cat. No. **Pack Size**

10 220 655 103 custom fill

Will be supplied as "Guanosine-5'-O-(3-thiotriphosphate), Li⁴".
Unit of Measure is "g".



For further processing only.

Guanosine-5'-O-(2-thiodiphosphate) trilithium salt

Non-hydrolyzable GDP analog.

Application

Use Guanosine-5'-O-(2-thiodiphosphate) to completely inhibit G protein activation by GTP and GTP analogs.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 71376-97-1

Properties

Formula: C₁₀H₁₂N₅O₁₀P₂SLi₃

Molecular weight: 477.0 D (GDP-β-S: 459.3 D)

Specification

Appearance: White powder

GDP-β-S (HPLC): ≥85 area%

Li (flame photometric): 4±1%

Water (K. Fischer): ≤12%

GMP (HPLC): ≤10 area%

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. **Pack Size**

10 526 134 103 custom fill

Will be supplied as "Guanosine-5'-O-(2-thiodiphosphate) Tri-Li".
Unit of Measure is "g".




For further processing only.

Enzymes for Clinical Chemistry

Colorimetric Tests	53	Glycerol Kinase (GK), concentrated	103
UV Tests	55	Glycerol-3-phosphate Dehydrogenase	104
Acetate-CoA Ligase (Acetyl-CoA Synthetase)	58	Glycerol-3-phosphate Oxidase	105
Acid Phosphatase	59	Glycerol-3-phosphate Oxidase, chemically modified	106
Acyl-CoA Oxidase	59	Hexokinase (HK)	107
Adenosine Deaminase	59	Hexokinase (HK), chemically modified	108
Alanine Aminotransferase (ALT) (GPT)	60	D-Lactate Dehydrogenase (D-LDH)	110
Alcohol Dehydrogenase	61	D-Lactate Dehydrogenase (D-LDH), Grade I	111
Alcohol Dehydrogenase, chemically modified	62	D-Lactate Dehydrogenase (D-LDH), Grade II	112
Aldehyde Dehydrogenase	64	L-Lactate Dehydrogenase (L-LDH)	113
Aldose 1-Epimerase (Mutarotase)	65	L-Lactate Dehydrogenase (L-LDH), chemically modified	116
Ascorbate Oxidase	66	Lactate 2-Monooxygenase (Lactate oxidase), Grade I	118
Ascorbate Oxidase, chemically modified	68	Lactate 2-Monooxygenase (Lactate oxidase), Grade II	119
Aspartate Aminotransferase (AST) (GOT)	69	Lactate 2-Monooxygenase (Lactate oxidase)	119
N-Carbamoylsarcosine Amidase	70	Lipase	120
Cholesterol Esterase	71	Lipoprotein Lipase	120
Cholesterol Esterase, chemically modified	71	Lipoprotein Lipase, chemically modified	121
Cholesterol Oxidase	75	Lysozyme	122
Citrate Lyase	80	Malate Dehydrogenase	123
Citrate Synthase	80	Malate Dehydrogenase, chemically modified	123
Colipase	81	Malate Dehydrogenase, IFCC Quality	124
Creatinase	82	N-Methylhydantoinase (ATP-hydrolyzing)	126
Creatininase	82	NAD(P)H Dehydrogenase (quinone) (Diaphorase)	126
Creatinine Deaminase	83	Nitrate Reductase	127
Formate Dehydrogenase	84	Oxalate Oxidase	128
Galactose 1-Dehydrogenase	84	Peroxidase (POD), Grade II	128
Glucose Oxidase (GOD), Grade I	86	Phosphogluconate Dehydrogenase (decarboxylating)	129
Glucose Oxidase (GOD), Grade II	87	6-Phosphogluconolactonase	130
Glucose Oxidase (GOD), chemically modified	88	Pyruvate Kinase	133
Glucose-6-phosphate Dehydrogenase (G6P-DH)	89	Pyruvate Oxidase	134
Glucose-6-phosphate Isomerase	95	Sarcosine Oxidase	135
α -Glucosidase	96	Triose-phosphate Isomerase	136
β -Glucuronidase	97	Thrombin	136
Glutamate Dehydrogenase (NAD(P))	97	Urease	137
γ -Glutamyltransferase	101	Uricase	138
Glycerol Kinase (GK)	102		

Colorimetric Tests

*Indicator reaction

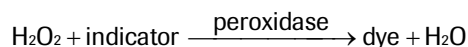


Peroxidase (POD), Grade II (128)

4-Aminoantipyrine (4-APP) (139)

TOOS (145)

Alkaline phosphatase (ALP)

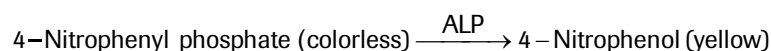

4-Nitrophenyl phosphate (4-NPP) (141)




Cholesterol



Cholesterol Esterase (71)

Cholesterol Oxidase (75)





α -Amylase

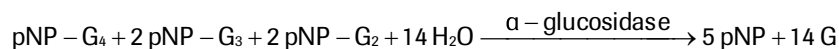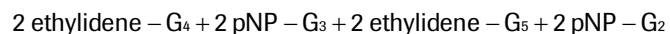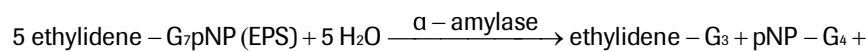

Ethylidene-4-NP-G7 (144)

α -Glucosidase (96)

Creatinine

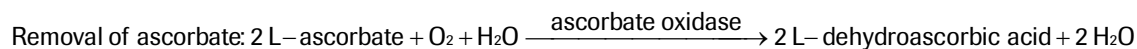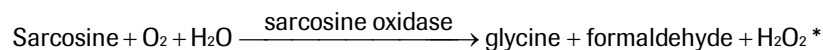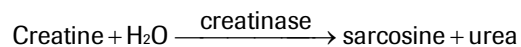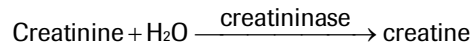



Creatininase (82)

Creatinase (82)

Sarcosine oxidase (135)

Ascorbate oxidase (66)

Glucose

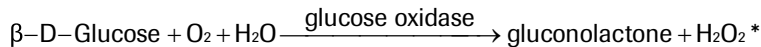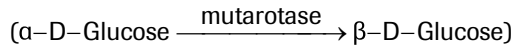


Glucose oxidase (GOD), Grade II (87)

Glucose oxidase (GOD), chemically modified (88)

Aldose 1-epimerase (Mutarotase) (65)

γ -Glutamyltransferase


Gluta - carboxylate (donor) + glycylglycine (acceptor)

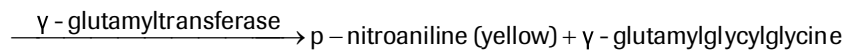

Gluta-carboxylate (144)


Glycylglycine (16)

Lactate

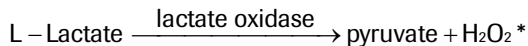

Lactate 2-monooxygenase (Lactate oxidase), Grade II (119)




Lipase

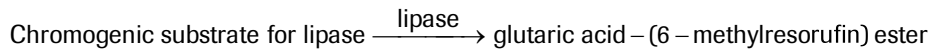

Chromogenic Substrate for Lipase (142)

Colipase (81)

Triglycerides

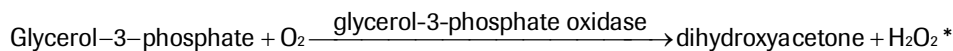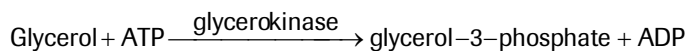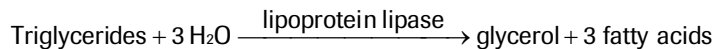


Lipoprotein lipase (120)

Glycerokinase (GK) (102)

Glycerol-3-phosphate oxidase (104)

ATP, Grade I (49)

Uric acid


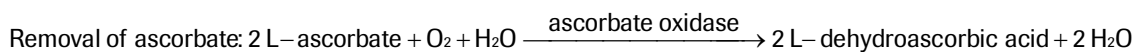


Uricase (138)

Ascorbate oxidase (66)

UV Tests

Alanine aminotransferase (ALT)

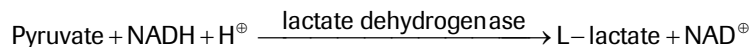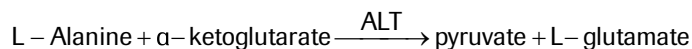



L-Alanine (147)

α -Ketoglutarate (2-Oxoglutarate) (147)

D-Lactate dehydrogenase (D-LDH) (110)

NADH, Grade I (39)

Aspartate aminotransferase (AST)

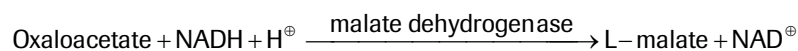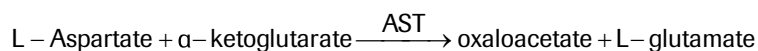




α -Ketoglutarate (2-Oxoglutarate) (147)

Malate dehydrogenase (123)

D-Lactate dehydrogenase (D-LDH), Grade II (112)

NADH, Grade I (39)

Creatine kinase

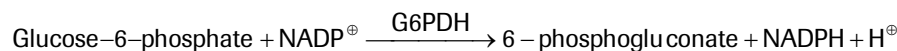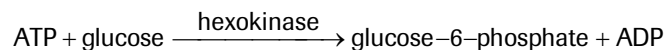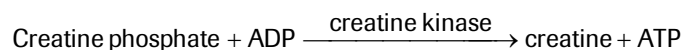


Creatine phosphate (150)

Hexokinase (107)

Glucose-6-phosphate dehydrogenase (G6P-DH) (89)

NADP (42)



Ethanol

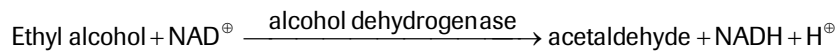

Alcohol dehydrogenase (63)

NAD, Grade I (38)

Glucose

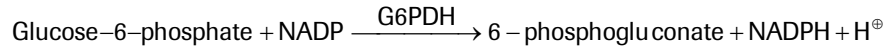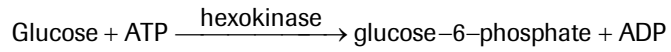


Hexokinase (107)

Glucose-6-phosphate dehydrogenase (G6P-DH) (89)

NADP (42)



Lactate dehydrogenase

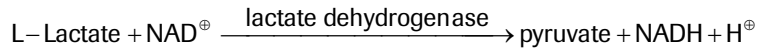

D(-)-Lactate (151)

NAD, Grade I (38)

Urea

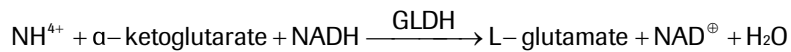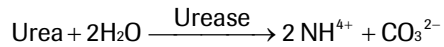

Urease (137)

α -Ketoglutarate (2-Oxoglutarate) (147)

NADH, Grade I (39)

Acetate-CoA Ligase (Acetyl-CoA Synthetase) from yeast, lyophilizate

Ligase that catalyzes the synthesis of acetyl-CoA from acetate and coenzyme A.

Application

Use Acetate-CoA Ligase (Acetyl-CoA Synthetase) in diagnostic tests for the determination of free fatty acids in combination with Acyl-CoA Oxidase, Catalog No. 10 885 550 103 or for the determination of acetic acid in combination with Citrate Synthase, Catalog No.10 153 605 103 and Malate Dehydrogenase, Catalog No. 10 200 387 103

Benefits

Rely on the proven diagnostic quality of this product.

EC 6.2.1.1

Specification

Appearance: White lyophilizate

pH value (c=20 mg/ml in water): 6.8-7.8

Specific activity (+37°C, acetate): ≥4 U/mg protein

Protein (Biuret): ≤0.25 mg/mg lyophilizate

Acetate (enzymatically): ≤0.1%

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 128 180 103 custom fill

Will be supplied as "Acetyl-CoA Synthetase from Yeast". Unit of Measure is "kU".

For further processing only.

Acetate-CoA Ligase (Acetyl-CoA Synthetase) from microorganism, lyophilizate

Ligase that catalyzes the synthesis of acetyl-CoA from acetate and coenzyme A.

Application

Use Acetate-CoA Ligase (Acetyl-CoA Synthetase) in diagnostic tests for the determination of free fatty acids in combination with Acyl-CoA Oxidase, Catalog No. 10 885 550 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 6.2.1.1

Specification

Appearance: White lyophilizate

Activity (+37°C, enzymatically): ≥1.5 U/mg lyophilizate

Specific Activity (enzymatically): ≥1.5 U/mg protein

Absorbance of the solution (2 U/ml):

A_{400} : ≤0.05

A_{500} : ≤0.025

A_{650} : ≤0.012

Contaminants:

Lipase (indirect): -5% bis +15%

Stability: At -15 to -25 °C within specification range for 12 months. Store dry in tightly sealed containers.

Cat. No. **Pack Size**

10 885 568 103 custom fill

Will be supplied as "Acyl-CoA-Synthetase, Lyo.". Unit of Measure is "kU".



For further processing only.

Acid Phosphatase from potato, lyophilizate

Application

Use Acid Phosphatase in your controls or calibrators.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.1.3.2

Specification

Appearance: Pale brown lyophilizate

Activity (+25°C, 4-nitrophenyl phosphate): ≥2 U/mg lyophilizate

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.	Pack Size
10 154 393 103	custom fill

Will be supplied as "Phosphatase, Acid, Grade II from Potato". Unit of Measure is "kU".

For further processing only.

Acyl-CoA Oxidase from microorganisms, lyophilizate

Oxidoreductase that catalyzes the interconversion of acyl-CoA to trans-2,3-dehydroacyl-CoA.

Application

Use Acyl-CoA Oxidase in diagnostic tests for the determination of free fatty acids in combination with Acetate-CoA Ligase (Acetyl-CoA Synthetase), Catalog Nos. 10 885 568 103 or 10 128 180 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.3.3.6

Specification

Appearance: Yellow lyophilizate

Activity (+37°C, enzymatically): ≥20 U/mg lyophilizate

Specific activity (enzymatically): ≥20 U/mg protein

Absorbance of the solution (20 U/ml):

A_{400} : ≤0.08

A_{500} : ≤0.04

A_{650} : ≤0.02

Contaminants:

Catalase: ≤12 U/U Acyl-CoA oxidase

Stability: At -15 to -25 °C within specification range for 12 months. Store dry in tightly sealed containers.

Cat. No.	Pack Size
10 885 550 103	custom fill

Will be supplied as "Acyl-CoA-Oxidase, Lyo.". Unit of Measure is "kU".



For further processing only.

Adenosine Deaminase from calf intestine, suspension

Application

Use Adenosine Deaminase for the deamination of adenosine analogs to the corresponding inosine analogs.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.4.4

Cat. No.	Pack Size
10 153 460 103	custom fill

Will be supplied as "Adenosine Deaminase (ADA), Calf Intestine". Unit of Measure is "g".

For further processing only.

Specification**Appearance:** White suspension in ammonium sulfate**Specific activity:** ≥ 200 U/mg protein**Protein** (Biuret): ≥ 10 mg/ml**Ammonium sulphate:** 3.2 ± 0.2 mol/l**Contaminants** (expressed as percentage of Adenosine Deaminase specific activity):AMP deaminase: ≤ 0.01 Guanase: ≤ 0.01 Nucleoside phosphorylase: ≤ 0.01 Phosphatase, alkaline: ≤ 0.01 **pH 5.5 treatment** (30 minutes): Corresponds to specification**Stability:** At +2 to +8°C within specification range for 24 months.**Alanine Aminotransferase (ALT) (GPT)
from pig heart, lyophilizate**

Cat. No.	Pack Size
10 170 674 103	custom fill

ApplicationUse Alanine Aminotransferase for designing your calibrator/control reagent and for the synthesis of unnatural L-amino acids from α -keto acids.Will be supplied as "GPT from Pig Heart". Unit of Measure is "kU".
For further processing only.**Benefits**

Rely on the proven diagnostic quality of this product.

EC 2.6.1.2

Specification**Appearance:** Slightly yellow lyophilizate**pH value** (c=10 mg/ml in water): 7.0-8.0**Activity** (+25°C, L-alanine, α -oxoglutarate): ≥ 3 U/mg lyophilizate**Activity** (+37°C, ALT (ALAT/GPT)-kit): ≥ 4.8 U/mg lyophilizate**Contaminants** (expressed as percentage of Alanine Aminotransferase activity):Contaminating oxidases (FOX): ≤ 0.7 Glutamate dehydrogenase: ≤ 0.01 Aspartate Aminotransferase (AST/GOT): ≤ 0.135 Lactate dehydrogenase: ≤ 0.01 Malate dehydrogenase: ≤ 0.01 **SVD free:** Corresponds to specification**pH 5.5 treatment** (30 minutes): Corresponds to specification**Stability:** At +2 to +8°C within specification range for 12 months. Store dry.**Alanine Amionotransferase (ALT) (GPT)
from pig heart, suspension**

Cat. No.	Pack Size
10 153 443 103	custom fill

ApplicationUse Alanine Aminotransferase for designing your calibrator/control reagent and for the synthesis of unnatural L-amino acids from α -keto acids.Will be supplied as "GPT from Pig Heart". Unit of Measure is "kU".
For further processing only.**Benefits**

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 2.6.1.2

Specification**Appearance:** Slightly yellow suspension in ammonium sulfate, 3.2 mol/l

pH value: 5.5-6.5

Specific activity (+25 °C; L-alanine, α-ketoglutarate): ≥80 U/mg protein

Protein (Biuret): ≥10 mg/ml (standardized to 10±1 mg/ml)

Ammonium sulphate: 3.2±0.2 mol/l

Contaminants (expressed as percentage of Alanine Aminotransferase activity):

Glutamate dehydrogenase: ≤0.01

Aspartate Aminotransferase (AST/GOT) : ≤0.03

Lactate dehydrogenase: ≤0.01

Malate dehydrogenase: ≤0.01

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Alcohol Dehydrogenase from yeast, lyophilizate

Dehydrogenase that catalyzes the interconversion of alcohols to the corresponding aldehydes.

Application

Use Alcohol Dehydrogenase in diagnostic tests for the determination of alcohol or aldehyde (formate).

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.1

Properties

Nomenclature: Alcohol:NAD⁺ oxidoreductase

Molecular weight: 141 kD (pH 7.0)

Isoelectric point: 5.4-5.8

Michaelis constants (Phosphate buffer, pH 7.15, +25°C):

Ethanol: 1.3×10^{-2} mol/l

NAD: 7.4×10^{-5} mol/l

Acetaldehyde: 7.8×10^{-4} mol/l

NADH: 1.1×10^{-5} mol/l

Inhibitor constants (Phosphate buffer, pH 7.15, +25°C):

Ethanol: 4.3×10^{-2} mol/l

NAD: 6.1×10^{-4} mol/l

Acetaldehyde: 6.7×10^{-4} mol/l

NADH: 1.8×10^{-5} mol/l

Inhibitors:

-SH-reagents and heavy metals, such as derivatives, 4-chloromercuribenzoate, iodoacetic acid, N-substituted maleinimides, Hg²⁺, Ag⁺ and Cu²⁺.

-Complexing agents, e.g., o-phenanthroline, EDTA, oxalate.

-NAD analogs and NAD partial structures, e.g., NADP, NADH, ADP, ADP-ribose.

-Substances, which react with enzyme bound NAD, e.g., sulfite, hydroxylamine, cyanide.

-Substrate analogs, e.g., fluoroethanol.

-Oxidizers, e.g., H₂O₂ and aerial oxygen inactivate by oxidation of essential groups.

pH optimum: 9.0 (see figure)

Temperature dependence: See figure


pH stability: 6.0-8.0 (see figure)

Thermal stability: Up to +50°C (see figure)


Specificity: Alcohol dehydrogenase oxidizes primary alcohols. Isopropanol and secondary butanol are slowly oxidized, while higher secondary and

Cat. No. 11 452 541 103 **Pack Size** custom fill


Will be supplied as "Alcohol Dehydrogenase, Yeast". Unit of Measure is "MU".

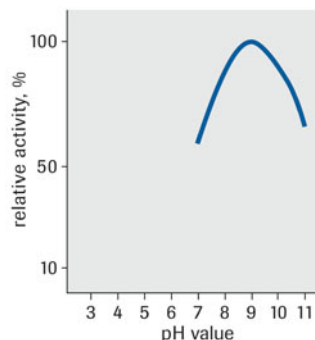

pH optimum

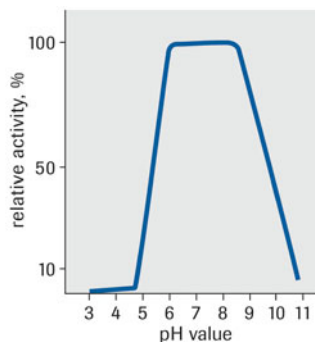

pH stability

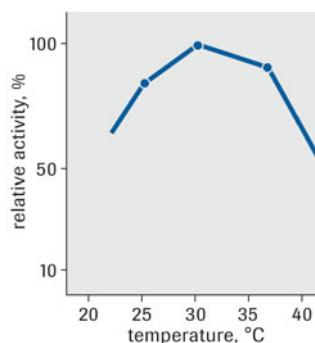

Temperature dependence

Incubation:
25°C, 120 min
pH 3.0 – 5.0:
citrate buffer, 0.2 mol/l
pH 6.0 – 8.0:
phosphate buffer,
0.2 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.2 mol/l
180 U ADH/ml

tertiary alcohols do not react. Numerous aldehydes are reduced in the reverse reaction. The enzyme does not react with NADP.

Remarks: Alcohol dehydrogenase tends to show turbidity in solution at +37°C storage. Modified Alcohol dehydrogenase shows no turbidity for at least 4 weeks in solution at +37°C.

Specification

Appearance: White lyophilizate (50 mg lyophilizate contain approximately 30 mg enzyme protein, 15 mg sucrose, 5 mg phosphate)

pH value (c=50 mg/ml in water): 7.0-8.0

Activity (+25°C, ethanol): No limit

Specific activity: ≥400 U/mg (protein)


Protein (Biuret): ≥0.5 mg/mg lyophilizate

Contaminants (expressed as percentage of Alcohol Dehydrogenase activity):

Lactate dehydrogenase: ≤0.01

Malate dehydrogenase ≤0.01

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

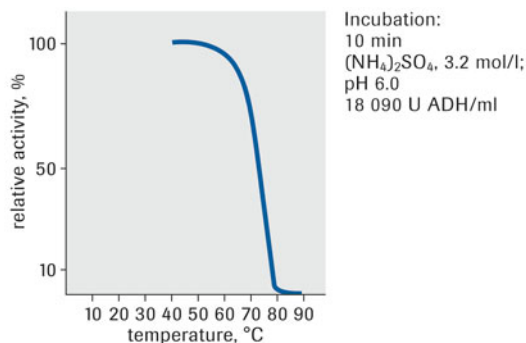

Thermal stability

For further processing only.

Alcohol Dehydrogenase, chemically modified from yeast, lyophilizate

Dehydrogenase that catalyzes the interconversion of alcohols to the corresponding aldehydes.

Application

Use Alcohol Dehydrogenase in diagnostic tests for the determination of alcohol or aldehyde (formate).

Benefits

- Take advantage of the enhanced liquid stability of this enzyme.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.1

Properties

Nomenclature: Alcohol:NAD⁺ oxidoreductase

Molecular weight: 141 kD (pH 7.0)

Isoelectric point: 5.4-5.8

Michaelis constants (Phosphate buffer, pH 7.15, +25°C):

Ethanol: 1.3×10^{-2} mol/l

NAD: 7.4×10^{-5} mol/l

Acetaldehyde: 7.8×10^{-4} mol/l

NADH: 1.1×10^{-5} mol/l

Inhibitor constants (Phosphate buffer, pH 7.15, +25°C):

Ethanol: 4.3×10^{-2} mol/l

NAD: 6.1×10^{-4} mol/l

Acetaldehyde: 6.7×10^{-4} mol/l

NADH: 1.8×10^{-5} mol/l

Inhibitors:

-SH-reagents and heavy metals, such as derivatives, 4-chloromercuribenzoate, iodoacetic acid, N-substituted maleinimides, Hg²⁺, Ag⁺ and Cu²⁺.

-Complexing agents, e.g., o-phenanthroline, EDTA, oxalate.

-NAD analogs and NAD partial structures, e.g., NADP, NADH, ADP, ADP-ribose.

-Substances, which react with enzyme bound NAD, e.g., sulfite, hydroxylamine, cyanide.

-Substrate analogs, e.g., fluoroethanol.

-Oxidizers, e.g., H₂O₂ and aerial oxygen inactivate by oxidation of essential groups.

Cat. No.

Pack Size

11 644 980 103 custom fill

Will be supplied as "Alcohol Dehydrogenase, Yeast, Modified".

Unit of Measure is "MU".

For further processing only.

pH optimum: 9.0 (see figure)
Temperature dependence: See figure
pH stability: 6.0-8.0 (see figure)
Thermal stability: Up to +50°C (see figure)
Specificity: Alcohol dehydrogenase oxidizes primary alcohols. Isopropanol and secondary butanol are slowly oxidized, while higher secondary and tertiary alcohols do not react. Numerous aldehydes are reduced in the reverse reaction. The enzyme does not react with NADP.

Specification

Appearance: White lyophilizate
Solubility: Clear, colorless solution in water (c = 50 mg/ml)
pH value: 6.5-8.0
Activity (+25°C, ethanol): ≥25 U/mg lyophilizate
Contaminants (expressed as percentage of Alcohol Dehydrogenase activity):
 Lactate dehydrogenase: ≤0.01
 Malate dehydrogenase: ≤0.01
Stability: At +2 to +8°C within specification range for 12 months.

Alcohol Dehydrogenase from yeast, suspension

Cat. No.	Pack Size
11 531 034 103	custom fill

Dehydrogenase that catalyzes the interconversion of alcohols to the corresponding aldehydes.

Will be supplied as "ADH-Y, As, new". Unit of Measure is "g".
For further processing only.

Application

Use Alcohol Dehydrogenase in diagnostic tests for the determination of alcohol or aldehyde (formate).

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.1

Properties

Nomenclature: Alcohol:NAD⁺ oxidoreductase
Molecular weight: 141 kD (pH 7.0)
Isoelectric point: 5.4-5.8
Michaelis constants (Phosphate buffer, pH 7.15, +25°C):
 Ethanol: 1.3×10^{-2} mol/l
 NAD: 7.4×10^{-5} mol/l
 Acetaldehyde: 7.8×10^{-4} mol/l
 NADH: 1.1×10^{-5} mol/l
Inhibitor constants (Phosphate buffer, pH 7.15, +25°C):
 Ethanol: 4.3×10^{-2} mol/l
 NAD: 6.1×10^{-4} mol/l
 Acetaldehyde: 6.7×10^{-4} mol/l
 NADH: 1.8×10^{-5} mol/l

Inhibitors:

- SH-reagents and heavy metals, such as derivatives, 4-chloromercuribenzoate, iodoacetic acid, N-substituted maleinimides, Hg²⁺, Ag⁺ and Cu²⁺.
- Complexing agents, e.g., o-phenanthroline, EDTA, oxalate.
- NAD analogs and NAD partial structures, e.g., NADP, NADH, ADP, ADP-ribose.
- Substances, which react with enzyme bound NAD, e.g., sulfite, hydroxylamine, cyanide.
- Substrate analogs, e.g., fluoroethanol.
- Oxidizers, e.g., H₂O₂ and aerial oxygen inactivate by oxidation of essential groups.

pH optimum: 9.0 (see figure)

Temperature dependence: See figure

pH stability: 6.0-8.0 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Alcohol dehydrogenase oxidizes primary alcohols. Isopropanol and secondary butanol are slowly oxidized, while higher secondary and tertiary alcohols do not react. Numerous aldehydes are reduced in the reverse reaction. The enzyme does not react with NADP.

Remark: Alcohol dehydrogenase tends to show turbidity in solution at +37°C storage.

Modified Alcohol dehydrogenase shows no turbidity for at least 4 weeks in solution at +37°C.

Specification

Appearance: White to yellow-brown crystalline suspension in ammonium sulfate solution, 3.2 mol/l, pH approximately 6.5

pH value: 6.0-7.0

Activity: 9,000 to 18,000 U/ml

Specific activity: ≥ 300 U/mg (protein)

Protein (Biuret): 30 ± 3 mg/ml lyophilizate

Contaminants (expressed as percentage of Alcohol Dehydrogenase activity):

Lactate dehydrogenase: ≤ 0.01

Malate dehydrogenase: ≤ 0.01

Stability: At +2 to +8°C within specification range for 9 months. Store under nitrogen.

Aldehyde Dehydrogenase from yeast, lyophilizate

Dehydrogenase that catalyzes the oxidation of aldehydes using NAD(P)⁺ as acceptor.

Application

Use Aldehyde Dehydrogenase in diagnostic tests that use an NADH/NADPH recycling system.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.2.1.5

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=20 mg/ml)

Activity (+25°C, acetaldehyde): ≥ 2.0 U/mg lyophilizate

Specific activity: ≥ 20 U/mg protein

Protein (Biuret): No limit (approximately 10%)

Contaminants (expressed as percentage of Aldehyde Dehydrogenase activity):

Alcohol dehydrogenase: ≤ 0.01

Lactate dehydrogenase: ≤ 0.01

"NADH oxidase": ≤ 0.01

"NADPH oxidase": ≤ 0.01

Stability: At +2 to +8°C within specification range for 12 months. Store dry. Store under nitrogen.

Cat. No.	Pack Size
10 145 947 103	custom fill

Will be supplied as "Aldehyde Dehydrogenase from Yeast". Unit of Measure is "kU".

For further processing only.

Aldose 1-Epimerase (Mutarotase) from pig kidney, suspension

Enzyme for mutarotation of sugars.

Application

Use Aldose 1-Epimerase (Mutarotase) in diagnostic tests for the determination of glucose anomers.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 5.1.3.3

Specification

Appearance: White suspension in ammonium sulfate, 3.2 mol/l, pH approximately 6

pH value: 5.5-6.5

Specific activity (+25°C, α-D-glucose): ≥5,000 U/mg protein

Protein (Biuret): 5±0.5 mg/ml

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
10 152 331 103	custom fill

Will be supplied as "Mutarotase from Hog Kidney". Unit of Measure is "MU".

For further processing only.

Ascorbate Oxidase from *Cucurbita* species, lyophilizate

Oxidoreductase that oxidizes ascorbic acid to dehydroascorbate.

Application

Use Ascorbate Oxidase in a variety of diagnostic tests to eliminate the interference of ascorbic acid, since ascorbic acid interferes with the Trinder reaction that is widely used for the colorimetric determination of analytes. It is useful in liquid as well as dry chemistry test, *e.g.*, for the determination of uric acid, lactate, creatine kinase or transaminases.

Benefits

- Rely on the proven diagnostic quality of this product.

EC 1.10.3.3

Properties

Nomenclature: L-ascorbate:oxygen oxidoreductase

Molecular weight: Approximately 140 kD

Isoelectric point: 5.0-6.0

Michaelis constant (Phosphate buffer, pH 5.6, +25°C):

L-ascorbate: 3×10^{-4} mol/l

Inhibitors: 4-chloromercuribenzoate, CN⁻, Na₂S, diethyl-dithiocarbamate, 8-hydroxyquinoline, K-ethylxanthate

pH optimum: 5.6-7.0 (see figure)

Temperature dependence: See figure

pH stability: 6.5-9.0 (see figure)

Thermal stability: Up to +70°C (see figure)

Stability of the lyophilizate: Stable at +35°C for at least 3 weeks (see figure).

Stability in solution: See figure

Specificity: Several analogs of ascorbate react.

Remark: A decrease in activity of approximately 10% may occur.

Specification

Appearance: Turquoise lyophilizate

Solubility: Clear, slightly turquoise solution in water (c=50 mg/ml)

pH value (c=50 mg/ml in water): 7.0-8.0

Activity (+25°C, L-ascorbate): ≥170 U/mg lyophilizate

Specific activity (+25°C): ≥1,700 U/mg protein

Activity (+37°C, L-ascorbate): ≥180 U/mg lyophilizate

Specific activity (+37°C): ≥1,800 U/mg protein

Protein (BCA): 0.07-0.14 mg/mg lyophilizate

Contaminants (+25°C; expressed as percentage of Ascorbate Oxidase activity):

Catalase: ≤0.2

Aspartate aminotransferase (AST/GOT): No limit

Alanine aminotransferase (ALT/GPT): No limit

Peroxidase: ≤0.005


Stability: At -15 to -25°C within specification range for 12 months. Store dry. Keep tightly sealed.

Cat. No.


Pack Size

10 199 605 103 custom fill


Will be supplied as "Ascorbate Oxidase from *Cucurbita* species". Unit of Measure is "kU".

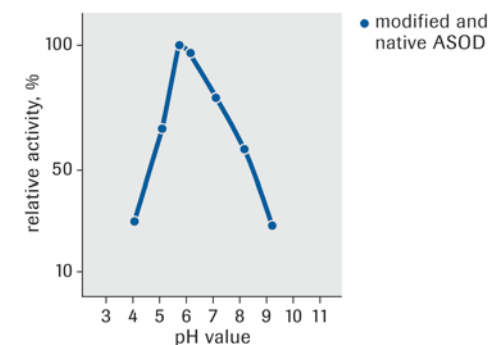


pH optimum

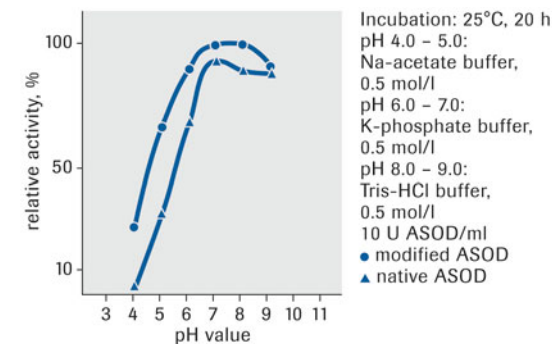


pH stability

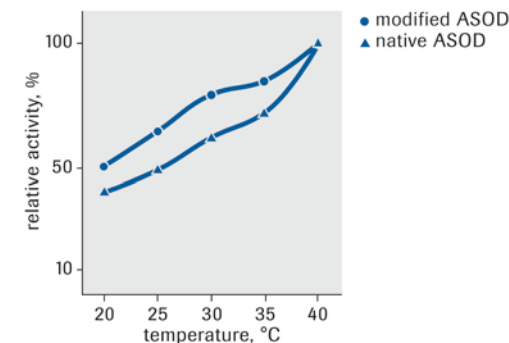


Thermal stability

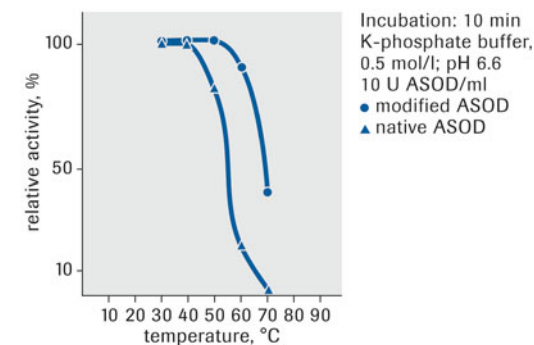

pH optimum

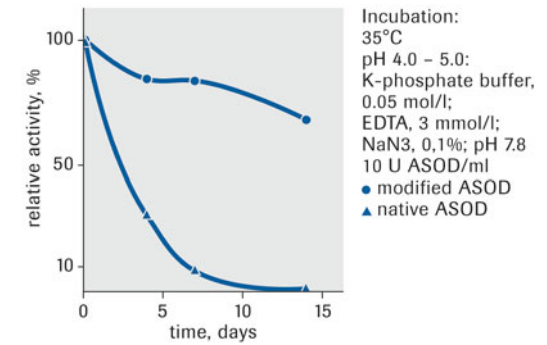

pH stability

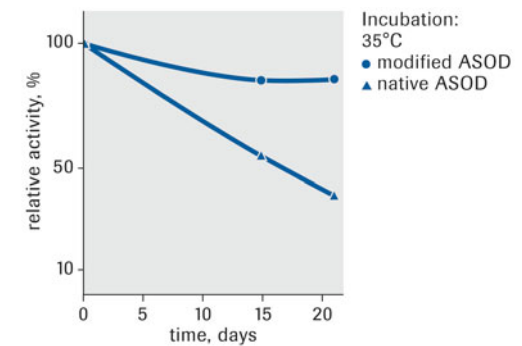

Thermal stability

For further processing only.

Ascorbate Oxidase

from *Cucurbita* species, poor of Aspartate aminotransferase (AST/GOT), lyophilizate

Oxidoreductase that oxidizes ascorbic acid to dehydroascorbate.

Application

Use Ascorbate Oxidase in a variety of diagnostic tests to eliminate the interference of ascorbic acid, since ascorbic acid interferes with the Trinder reaction that is widely used for the colorimetric determination of analytes. It is useful in liquid as well as dry chemistry test, *e.g.*, for the determination of uric acid, lactate, creatine kinase or transaminases.

Benefits

- Rely on the tested deficiency of aspartate transaminase.
- Rely on the proven diagnostic quality of this product.

EC 1.10.3.3

Properties

Nomenclature: L-ascorbate:oxygen oxidoreductase

Molecular weight: Approximately 140 kD

Isoelectric point: 5.0-6.0

Michaelis constant (Phosphate buffer, pH 5.6, +25°C):

L-ascorbate: 3×10^{-4} mol/l

Inhibitors: 4-chloromercuribenzoate, CN⁻, Na₂S, diethyl-dithiocarbamate, 8-hydroxyquinoline, K-ethylxanthate

pH optimum: 5.6-7.0 (see figure)

Temperature dependence: See figure

pH stability: 6.5-9.0 (see figure)

Cat. No. **Pack Size**

11 136 364 103 custom fill

Will be supplied as "Ascorbate Oxidase GOT-deficient". Unit of Measure is "MU".



For further processing only.

Thermal stability: Up to +70°C (see figure)
Stability of the lyophilizate: Stable at +35°C for at least 3 weeks (see figure).
Stability in solution: See figure
Specificity: Several analogs of ascorbate react.
Remark: A decrease in activity of approximately 10% may occur.

Specification

Appearance: Turquoise lyophilizate
Solubility: Clear, slightly turquoise solution in water (c=50 mg/ml)
pH value (c=50 mg/ml in water): 7.0-8.0
Activity (+25°C, L-ascorbate): ≥170 U/mg lyophilizate
Specific activity (+25°C): ≥1,700 U/mg protein
Activity (+37°C, L-ascorbate): ≥180 U/mg lyophilizate
Specific activity (+37°C): ≥1,800 U/mg protein
Protein (BCA): 0.08-0.14 mg/mg lyophilizate
Contaminants (+25°C, expressed as percentage of Ascorbate Oxidase activity):
 Catalase: ≤0.2
 Aspartate aminotransferase (AST/GOT): ≤0.0003
 Alanine aminotransferase (ALT/GPT): ≤0.0005
 Contaminating oxidases (FOX): ≤0.0002
Stability: At -15 to -25°C within specification range for 12 months. Store dry. Keep tightly sealed.

Ascorbate Oxidase, chemically modified from *Cucurbita* species, lyophilizate

Oxidoreductase that oxidizes ascorbic acid to dehydroascorbate.

Application

Use Ascorbate Oxidase, chemically modified, in a variety of diagnostic tests to eliminate the interference of ascorbic acid, since ascorbic acid interferes with the Trinder reaction that is widely used for the colorimetric determination of analytes. It is useful in liquid as well as dry chemistry test, *e.g.*, for the determination of uric acid, lactate, creatine kinase or transaminases.

Benefits

- Take advantage of the improved stability in liquid reagents.
- Rely on the proven diagnostic quality of this product.

EC 1.10.3.3

Properties

Nomenclature: L-ascorbate:oxygen oxidoreductase
Molecular weight: Approximately 140 kD
Isoelectric point: 5.0-6.0
Michaelis constant (Phosphate buffer, pH 5.6, +25°C):
 L-ascorbate: 3×10^{-4} mol/l
Inhibitors: 4-chloromercuribenzoate, CN⁻, Na₂S, diethyl-dithiocarbamate, 8-hydroxyquinoline, K-ethylxanthate
pH optimum: 5.6-7.0 (see figure)
Temperature dependence: See figure
pH stability: 6.5-9.0 (see figure)
Thermal stability: Up to +70°C (see figure)
Stability of the lyophilizate: Stable at +35°C for at least 3 weeks (see figure).
Stability in solution: See figure
Specificity: Several analogs of ascorbate react.
Remark: The modified enzyme is especially suited for liquid stable applications with extended shelf life requirements.

Cat. No.	Pack Size
11 558 668 103	custom fill

Will be supplied as "AOD, modified". Unit of Measure is "MU".
 For further processing only.

Specification

Appearance: Turquoise lyophilizate
Solubility: Clear, slightly turquoise solution in water (c=50 mg/ml)
pH value (c=50 mg/ml in water): 7.2-8.2
Activity (+25°C, L-ascorbate): ≥ 120 U/mg lyophilizate
Specific activity (+25°C): ≥ 1,200 U/mg protein
Activity (+37°C, L-ascorbate): ≥ 180 U/mg lyophilizate
Specific activity (+37°C): ≥ 1,800 U/mg protein
Protein (BCA): 0.04-0.10 mg/mg lyophilizate
Contaminants (+25°C, expressed as percentage of Ascorbate Oxidase activity):
 Catalase: ≤ 0.2
 Glutamate oxalacetate transaminase (AST): ≤ 0.0003
 Glutamate pyruvate transaminase (ALT): ≤ 0.0005
 Contaminating oxidases (FOX): ≤ 0.0002
Stability: At -15 to -25°C within specification range for 12 months. Store dry. Keep tightly sealed.

Aspartate Aminotransferase (AST) (GOT) from pig heart, lyophilizate

Cat. No.	Pack Size
10 170 666 103	custom fill

Will be supplied as "GOT from Pig Heart". Unit of Measure is "kU".
For further processing only.

Application

Use Aspartate Aminotransferase for designing your calibrator/control reagent and for the synthesis of unnatural L-amino acids from α-keto acids.

Benefits

- Rely on the proven diagnostic quality of this product.

EC 2.6.1.1

Specification

Appearance: Slightly yellow lyophilizate
pH value (c=10 mg/ml in water): 7.0-8.0
Activity (+37°C, AST (ASAT/GOT)-kit): ≥ 45 U/mg lyophilizate
Contaminants (expressed as percentage of Aspartate Aminotransferase activity):
 Contaminating oxidases (FOX): ≤ 0.7
 Glutamate dehydrogenase: ≤ 0.01
 Alanine Aminotransferase (ALT/GPT) : ≤ 0.01
 Lactate dehydrogenase: ≤ 0.01
 Malate dehydrogenase: ≤ 0.01
 Oxaloacetate decarboxylase: ≤ 0.01
SDV free: Corresponds to specification
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Aspartate Aminotransferase (AST) (GOT) from pig heart, suspension

Cat. No.	Pack Size
10 153 354 103	custom fill

Will be supplied as "GOT from Pig Heart". Unit of Measure is "MU".
For further processing only.

Application

Use Aspartate Aminotransferase for designing your calibrator/control reagent and for the synthesis of unnatural L-amino acids from α-keto acids.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 2.6.1.1

Specification**Appearance:** Yellow suspension in ammonium sulfate**pH value:** 5.5-6.5**Specific activity** (+25°C; L-aspartate, α-ketoglutarate): ≥200 U/mg protein**Protein:** ≥10 mg/ml (standardized to 10±1 mg/ml)**Ammonium sulphate:** 3.2 0.2 mol/l**Contaminants** (expressed as percentage of Aspartate Aminotransferase activity):

Glutamate dehydrogenase: ≤0.01

Alanine Aminotransferase (ALT/GPT) : ≤0.01

Lactate dehydrogenase: ≤0.01

Malate dehydrogenase: ≤0.01

Oxaloacetate decarboxylase: ≤0.01

SDV free: Corresponds**pH 5.5 treatment** (30 minutes): Corresponds to specification**Stability:** At +2 to +8°C within specification range for 24 months.**N-Carbamoylsarcosine Amidase
from *E.coli* overproducer, lyophilizate**

Hydrolase that catalyzes the interconversion of N-carbamoylsarcosine to sarcosine.

Application

Use N-Carbamoylsarcosine Amidase in diagnostic tests for the determination of creatinine in combination with Creatinine Deaminase Catalog No. 11 330 764 103, N-Methylhydantoinase (ATP-hydrolysing), Catalog No. 11 288 555 103, and Sarcosine Oxidase, Catalog No. 11 378 856 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.1.59

Specification**Appearance:** White lyophilizate**Solubility:** Clear, colorless solution in water (c=10 mg/ml)**pH value** (c=10 mg/ml in water): 7.3-8.3**Activity** (+25°C, carbamoylsarcosine): 0.80-1.30 U/mg lyophilizate**Protein** (Biuret): 30-50 mg/100 mg lyophilizate**Contaminants** (expressed as percentage of Carbamoylsarcosine Amidase activity):

Creatinase: ≤0.013

Creatininase: ≤0.01

Catalase: ≤30

Uricase: ≤0.01

Stability: At -15 to -25°C within specification range for 12 months. Store dry. Protect from light.**Cat. No.** **Pack Size****11 248 847 103** custom fill

Will be supplied as "N-Carbamoyl-sarcosine Hydrolase". Unit of Measure is "kU".



For further processing only.

Cholesterol Esterase, Grade I from *Pseudomonas species*, lyophilizate

Hydrolase that splits fatty acids from sterols.

Application

Use Cholesterol Esterase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Oxidase, Catalog Nos. 10 634 522 103, 10 129 054103 or 11 479 709 103.

Benefits

- Use this enzyme for you liquid applications.
- Rely on the proven diagnostic quality of this product.

EC 3.1.1.13

Properties

Nomenclature: Sterol-ester acylhydrolase

Molecular weight: ~129 kD

Isoelectric point: 4.5

Michaelis constant (Phosphate buffer, pH 7.5):

Cholesterol oleate: 7×10^{-5} mol/l

Inhibitors: Heavy metals such as Cu^{2+} , Ag^{+} , Zn^{2+}

Activators: Detergents

pH optimum: 6.0-8.0; (maximum at pH 7.6) (see figure)

Temperature dependence: Not possible to determine under assay conditions due to turbidity of Thesit at temperatures above +27°C.

pH stability: 6.0-6.5 (see figure)

Thermal stability: Below +20°C (see figure)

Specificity: Cholesterol esterase is an enzyme of lipid metabolism and gives complete cleavage of all serum cholesterol esters.

Remark: This Cholesterol esterase is especially suited for liquid stable applications with extended shelf life requirements.

Specification

Appearance: Brownish lyophilizate

Solubility: Clear, brown solution in water (c=50 mg/ml)

pH value (c=50 mg/ml in water): 7.0-8.0

Activity (+25°C, cholesterol oleate): ≥ 100 U/mg lyophilizate

Specific Activity: ≥ 100 U/mg protein

Protein (Biuret): No limit

Contaminants (expressed as percentage of Cholesterol Esterase activity):

ATPase: ≤ 0.005

Catalase: ≤ 1.00

Glycerokinase: ≤ 0.001

Glucose oxidase: ≤ 0.001

Hexokinase: ≤ 0.005

"NADH oxidase": ≤ 0.001


Uricase: ≤ 0.005

Stability: At +2 to +8°C within specification range for 12 months. Store dry.


Cat. No. **Pack Size**

11 520 857 103 custom fill

Will be supplied as "CE, Ps.species, Lyo., SQ". Unit of Measure is "MU".

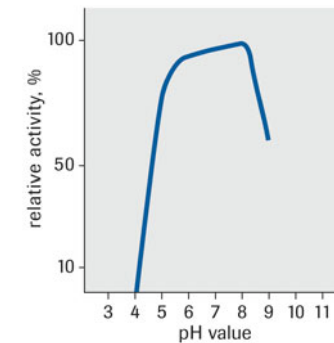


pH optimum

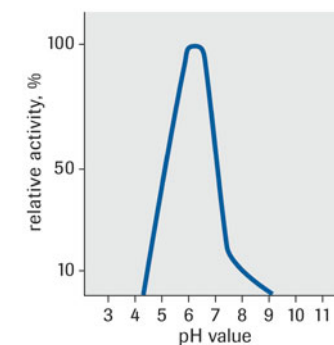

pH stability

Incubation:
25°C, 25 h
K-phosphate buffer,
0.7 mol/l
42.9 U CE/ml

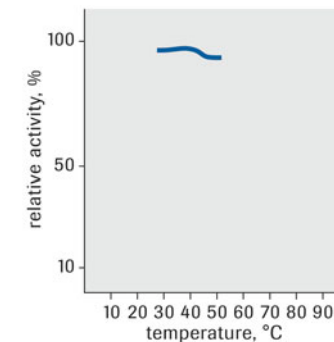

Thermal stability

Incubation:
10 min
K-phosphate buffer,
0.05 mol/l;
pH 6.5
18 U CE/ml

For further processing only.

Cholesterol Esterase, chemically modified from *Pseudomonas species*, lyophilizate

Hydrolase that splits fatty acids from sterols.

Cat. No. **Pack Size**

11 641 735 103 custom fill

Will be supplied as "Cholesterol Esterase Modified". Unit of Measure is "MU".

For further processing only.

Application

Use Cholesterol Esterase, chemically modified in diagnostic tests for the determination of cholesterol in combination with Cholesterol Oxidase, Catalog Nos. 11 479 709 103, 10 634 522 103 or 10 129 054.

Benefits

- Take advantage of the enhanced stability of this enzyme in liquid reagents.
- Rely on the proven diagnostic quality of this product.

EC 3.1.1.13

Properties

Nomenclature: Sterol-ester acylhydrolase

Molecular weight: ~129 kD

Isoelectric point: 4.5

Michaelis constant (Phosphate buffer, pH 7.5):

Cholesterol oleate: 7×10^{-5} mol/l

Inhibitors: Heavy metals such as Cu^{2+} , Ag^{+} , Zn^{2+}

Activators: Detergents

pH optimum: 6.0-8.0; (maximum at pH 7.6) (see figure)

Temperature dependence: Not possible to determine under assay conditions due to turbidity of Thesit at temperatures above +27°C.

pH stability: 6.0-6.5 (see figure)

Thermal stability: Below +20°C (see figure)

Specificity: Cholesterol esterase is an enzyme of lipid metabolism and gives complete cleavage of all serum cholesterol esters.

Specification

Appearance: Brownish lyophilizate

Solubility: Clear, brown solution in water (c=50 mg/ml)

pH value (c=50 mg/ml, in water): 7.0-8.0

Activity (+25°C, cholesterol oleate): ≥ 10 U/mg lyophilizate

Specific Activity: ≥ 100 U/mg protein

Volume Activity: For information only

Contaminants (expressed as percentage of Cholesterol Esterase activity):

ATPase: ≤ 0.005

Catalase: ≤ 1.00

Glycerokinase: ≤ 0.001

Hexokinase: ≤ 0.005

"NADH oxidase": ≤ 0.001

Uricase: ≤ 0.005

NaCl: 3 ± 0.2 mol/l

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cholesterol Esterase, Grade II from *Pseudomonas* species, lyophilizate

Hydrolase that splits fatty acids from sterols.

Application

Use Cholesterol Esterase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Oxidase, Catalog Nos. 10 634 522 103, 10 129 054 103 or 11 479 709 103.

Benefits

- Use this enzyme for your dry chemistry applications.
- Rely on the proven diagnostic quality of this product.

Cat. No.	Pack Size
11 015 923 103	custom fill

Will be supplied as "Cholesterol Esterase from Pseudom.species".

Unit of Measure is "MU".

For further processing only.

EC 3.1.1.13

Properties

Nomenclature: Sterol-ester acylhydrolase

Molecular weight: ~129 kD

Isoelectric point: 4.5

Michaelis constant (Phosphate buffer, pH 7.5):

Cholesterol oleate: 7×10^{-5} mol/l

Inhibitors: Heavy metals such as Cu^{2+} , Ag^+ , Zn^{2+}

Activators: Detergents

pH optimum: 6.0-8.0; (maximum at pH 7.6) (see figure)

Temperature dependence: Not possible to determine under assay conditions due to turbidity of Thesit at temperatures above +27°C.

pH stability: 6.0-6.5 (see figure)

Thermal stability: Below +20°C (see figure)

Specificity: Cholesterol esterase is an enzyme of lipid metabolism and gives complete cleavage of all serum cholesterol esters.

Specification

Appearance: Slightly yellowish lyophilizate

Solubility: Clear, colorless solution in water (c=50 mg/ml)

pH value (c=50 mg/ml in water): 7.0-8.0

Protein (Biuret): 0.14-0.20 mg/mg lyophilizate

Activity (+25°C, cholesterol oleate): ≥ 25 U/mg lyophilizate

Specific Activity: ≥ 100 U/mg protein

Contaminants (expressed as percentage of Cholesterol Esterase activity):

ATPase: ≤ 0.005

Catalase: ≤ 200 U/mg lyophilizate

Glycerokinase: ≤ 0.001

GOD: ≤ 0.001

Hexokinase: ≤ 0.005

"NADH oxidase": ≤ 0.005

Uricase: ≤ 0.005

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cholesterol Esterase from *Candida cylindracea*, lyophilizate

Hydrolase that splits fatty acids from sterols.

Application

Use Cholesterol Esterase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Oxidase, Catalog Nos. 10 634 522 103, 10 129 054 103 or 11 479 709 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.1.1.13

Specification

Appearance: Almost white lyophilizate

Solubility: Clear, colorless solution in phosphate buffer, 0.05 mol/l, pH 6.0 (c=10 mg/ml)

pH value: 5.5-6.5

Activity (+25°C; cholesterol oleate): ≤ 10.5 U/mg lyophilizate


Protein (Lowry): 0.20-0.30 mg/mg lyophilizate

Contaminants (expressed as percentage of Cholesterol Esterase activity):

Cat. No. **Pack Size**


10 129 046 103 custom fill

Will be supplied as "Cholesterol Esterase, *Candida cylindracea*".
Unit of Measure is "MU".

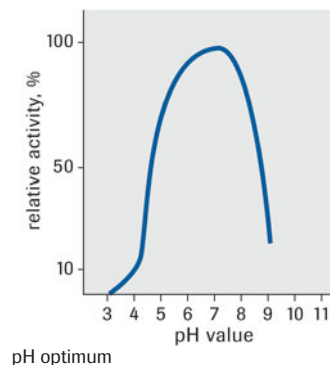


ATPase: ≤ 0.005
 Glucose oxidase: ≤ 0.001
 Glycerokinase: ≤ 0.001
 Hexokinase: ≤ 0.005
 Catalase: ≤ 1 U/mg lyophilizate
 "NADH oxidase": ≤ 0.005
 Proteases: No limit
 Uricase: ≤ 0.005

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

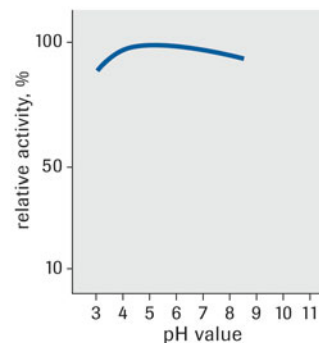

Incubation:
 25°C, 20 h
 pH 3.0 – 5.0:
 Na-acetate buffer,
 0.7 mol/l
 pH 4.0 – 8.2:
 K-phosphate buffer,
 0.7 mol/l
 48.7 U CE/ml

pH stability

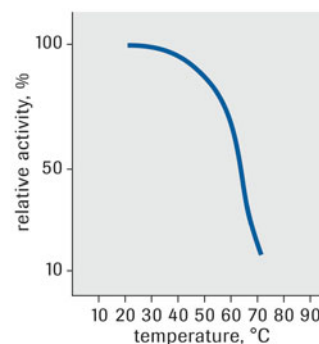

Incubation:
 10 min
 K-phosphate buffer,
 0.05 mol/l;
 pH 6.5
 48.7 U CE/ml

Thermal stability

For further processing only.

Cholesterol Esterase from *Candida cylindracea*, solution

Hydrolase that splits fatty acids from sterols.

Application

Use Cholesterol Esterase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Oxidase, Catalog Nos. 10 634 522 103, 10 129 054 103 or 11 479 709 103. Apply this ready-to-use enzyme directly in your diagnostic test.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.1.1.13

Specification

Appearance: Clear to turbid, brownish-yellow solution in NaCl

pH value: 5.7-6.3

Specific activity (+25°C; cholesterol oleate): ≥ 26 U/mg

Protein (Lowry): ≥ 2 mg/ml

NaCl (chloride meter): 3 ± 0.2 mol/l


Contaminants (expressed as percentage of Cholesterol Esterase activity):

ATPase: ≤ 0.005
 Glucose oxidase: ≤ 0.001
 Glycerokinase: ≤ 0.001
 Hexokinase: ≤ 0.005
 Catalase: ≤ 200 U/mg
 "NADH oxidase": ≤ 0.005
 Uricase: ≤ 0.005


Cat. No. **Pack Size**

10 262 609 103 custom fill

Will be supplied as "Cholesterol Esterase, *Candida cylindracea*".
 Unit of Measure is "MU".

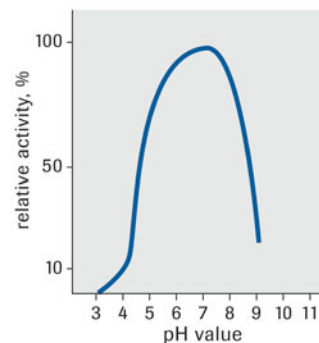


pH optimum

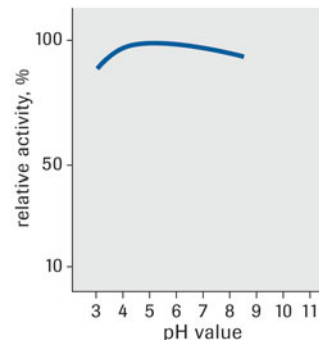

Incubation:
 25°C, 20 h
 pH 3.0 – 5.0:
 Na-acetate buffer,
 0.7 mol/l
 pH 4.0 – 8.2:
 K-phosphate buffer,
 0.7 mol/l
 48.7 U CE/ml

Stability: At +2 to +8°C within specification range for 12 months.

pH stability

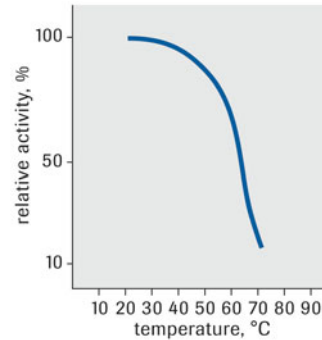

Incubation:
10 min
K-phosphate buffer,
0.05 mol/l;
pH 6.5
48.7 U CE/ml

Thermal stability

For further processing only.

Cholesterol Oxidase

from *Brevibacterium* species, expressed in *E.coli*, lyophilizate

Oxidoreductase that catalyzes the interconversion of cholesterol to cholest-4-en-3-one.

Application

Use Cholesterol Oxidase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Esterase, Catalog Nos. 10 129 046 103, 10 262 609 103, 11 015 923 103, 11 520 857 103 or 11 641 735 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.3.6

Properties

Nomenclature: Cholesterol:oxygen oxidoreductase

Molecular weight: 60 kD (native and SDS)

Isoelectric point: ~5.0

Michaelis constant (Phosphate buffer, 0.5 mol/l, pH 7.5; +25°C):

Cholesterol: 1×10^{-4} mol/l

Inhibitors: Hg²⁺, ZnCl₂, SDS

Activators: Non ionic detergents

pH optimum: 5.5-8.0 (see figure)

Temperature dependence: See figure

pH stability: 5.0-10.0 (see figure)

Thermal stability: Up to +55°C (see figure)

Storage and Stability: No decrease in activity over 6 weeks at +35°C (see figure)

Specificity:

cholesterol 100%

pregnenolon 52%

stigmasterol 17%

dehydroandrosterone 0.5%

androsterone 0%

estradiol 0%

cholate 0%

Specification


Appearance: Yellow lyophilizate

Solubility: Clear, yellowish solution in water (c=10 mg/ml)

Cat. No. **Pack Size**


11 479 709 103 custom fill

Will be supplied as "ChOD, Brevibacterium rec.". Unit of Measure is "MU".

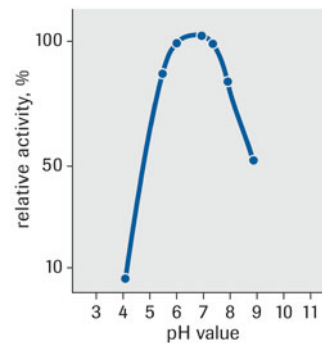

pH 4.0:
citrate buffer, 0.5 mol/l
pH 5.5 – 8.0:
phosphate buffer,
0.5 mol/l;
pH 9.0:
glycine buffer,
0.5 mol/l
0.22 U CO/ml

pH optimum

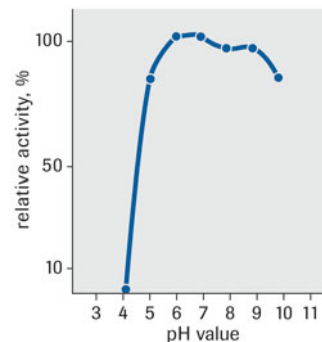

Incubation:
25°C, 24 h
pH 4.0 – 5.0:
citrate buffer, 0.5 mol/l
pH 5.5 – 8.0:
phosphate buffer,
0.5 mol/l
pH 9.0 – 10.0:
glycine buffer,
0.5 mol/l
3 U CO/ml

pH stability

pH value: 6.0-7.0

Protein (Biuret): 0.1-0.3 mg/mg lyophilizate

Activity (+25°C, cholesterol): 10-20 U/mg lyophilizate

Contaminants (expressed as percentage of Cholesterol Oxidase activity):

Catalase: ≤6.0

Glucose oxidase: ≤0.01


"NADH oxidase": ≤0.01

Uricase: ≤0.01


Stability: At -15 to -25°C within specification range for 12 months. Store dry. Protect from light.

Literature


- 1) T. Ohta, K. Fujishiro, K. Yamaguchi, Y. Tamura, K. Aisaka, T. Unajima, M. Hasegawa, *Gene* **103**, 93 (1991)
- 2) A. Vrielink, L.F. Lloyd, D.M. Blow, *J. Mol. Biol.* **219**, 533 (1991)
- 3) K. Fujishiro, T. Ohta, M. Hasegawa, K. Yamaguchi, T. Mizukami, T. Uwajima, *Biochem. Biophys. Res. Comm.* **172**, 721 (1990)

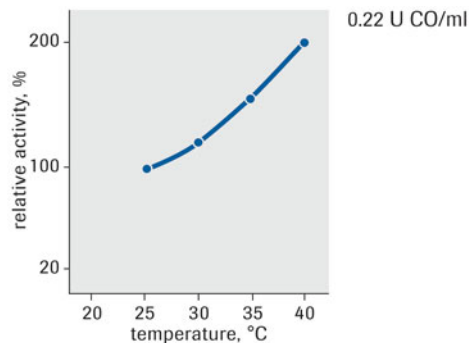

Temperature dependence

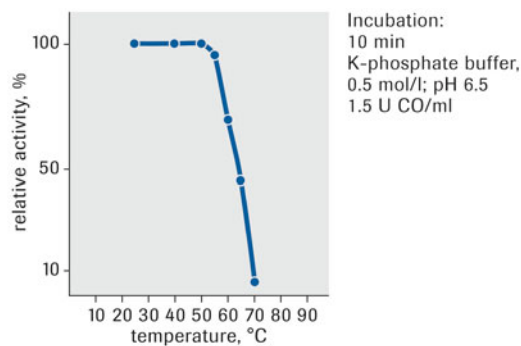

Thermal stability

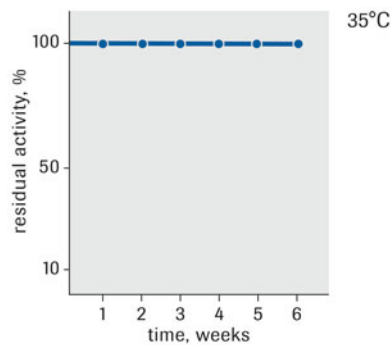

Stability of the lyophilizate

For further processing only.

Cholesterol Oxidase from *Nocardia erythropolis*, lyophilizate

Oxidoreductase that catalyzes the interconversion of cholesterol to cholest-4-en-3-one.

Application

Use Cholesterol Oxidase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Esterase, Catalog Nos. 10 129 046 103, 10 262 609 103, 11 015 923 103, 11 520 857 103 or 11 641 735 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.3.6

Cat. No.	Pack Size
10 129 054 103	custom fill

Will be supplied as "Cholesterol Oxidase, *Nocardia erythropolis*".
Unit of Measure is "MU".



Properties

Nomenclature: Cholesterol:oxygen oxidoreductase

Molecular weight: ~59 kD

Isoelectric point: 4.85

Michaelis constants (Cholesterol):

Phosphate buffer, pH 7.0: 1×10^{-6} mol/l

Triton X-100/isopropanol, pH 7.0: 7×10^{-6} mol/l

Inhibitors: $ZnCl_2$, Brij 35, Tween 40, Tween 60, Hg^{2+} , GSH, sodium dodecyl sulfate, laurylbenzenesulfonate

Activators: Hydroxypolyethoxydodecane, Triton X-100, DOC

pH optimum: 7.0-8.0 (see figure)

Temperature dependence: See figure

pH stability: 5.5-6.5 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Cholesterol (100%), b-cholestanol, pregnenolone, b-sitosterol and stigmasterol react as substrates (65-70%). Dehydroisoandrosterone and ergosterol show 5% relative activity. Androsterone, testosterone, β -estradiol, vitamin D3, cholic acid and cholesterol acetate do not react.

Specification

Appearance: Yellow-brown lyophilizate

Solubility: Clear, yellowish solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 5.5-6.5

Activity (+25°C, cholesterol): ≥ 2.5 U/mg lyophilizate

Protein (Lowry): 0.05-0.12 mg/mg lyophilizate

Contaminants (expressed as percentage of Cholesterol Oxidase activity):

Catalase: ≤ 6.0


Cholesterol esterase: No limit

Glucose oxidase: ≤ 0.01


"NADH oxidase": ≤ 0.01

Uricase: ≤ 0.01


Stability: At -15 to -25°C within specification range for 15 months. Store dry.

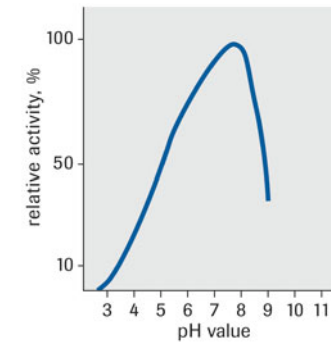


pH optimum

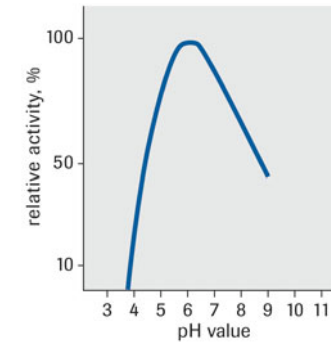

pH stability

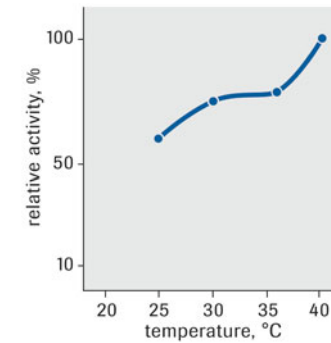

Temperature dependence

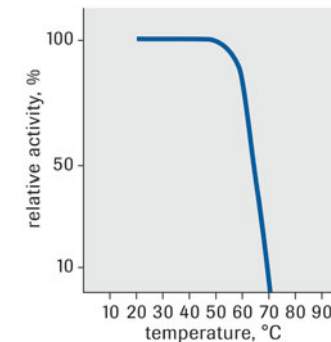

Thermal stability

Incubation:
25°C, 25 h
pH 4.2- 8.8:
K-phosphate buffer,
0.5 mol/l
10.9 U CO/ml

Incubation:
10 min
K-phosphate buffer,
50 mmol/l; pH 6.5
2.24 U CO/ml

For further processing only.

Cholesterol Oxidase

from *Nocardia erythropolis*, CE ≤0.005%, solution

Cat. No.	Pack Size
10 262 595 103	custom fill

Oxidoreductase that catalyzes the interconversion of cholesterol to cholest-4-en-3-one.

Application

Use Cholesterol Oxidase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Esterase, Catalog Nos. 10 129 046 103, 10 262 609 103, 11 015 923 103, 11 520 857 103 or 11 641 735 103.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.3.6

Properties

Nomenclature: Cholesterol:oxygen oxidoreductase

Molecular weight: ~59 kD

Isoelectric point: 4.85

Michaelis constants (Cholesterol):

Phosphate buffer, pH 7.0: 1×10^{-6} mol/l

Triton X-100/isopropanol, pH 7.0: 7×10^{-6} mol/l

Inhibitors: ZnCl₂, Brij 35, Tween 40, Tween 60, Hg²⁺, GSH, sodium dodecyl sulfate, laurylbenzenesulfonate

Activators: Hydroxypolyethoxydodecane, Triton X-100, DOC

pH optimum: 7.0-8.0 (see figure)

Temperature dependence: See figure

pH stability: 5.5-6.5 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Cholesterol (100%), b-cholestanol, pregnenolone, b-sitosterol and stigmasterol react as substrates (65-70%). Dehydroisoandrosterone and ergosterol show 5% relative activity. Androsterone, testosterone, β-estradiol, vitamin D₃, cholic acid and cholesterol acetate do not react.

Specification

Appearance: Brownish-yellow, slightly turbid solution in NaCl solution, 3 mol/l

pH value: 5.7-6.3

Specific Activity: 25 U/mg protein

Protein: 1±0.1 mg/ml

NaCl (chloride meter): 3±0.2 mol/l

Contaminants (expressed as percentage of Cholesterol Oxidase activity):

Cholesterol esterase: ≤0.005

Glucose oxidase: ≤0.01

"NADH oxidase": ≤0.01

Uricase: ≤0.01

Stability: At +2 to +8°C within specification range for 12 months.

Will be supplied as "Chol. Oxidase, Nocardia erythropolis". Unit of Measure is "kU".

For further processing only.

Cholesterol Oxidase from *Streptomyces* species, lyophilizate

Oxidoreductase that catalyzes the interconversion of cholesterol to cholest-4-en-3-one.

Application

Use Cholesterol Oxidase in diagnostic tests for the determination of cholesterol in combination with Cholesterol Esterase, Catalog Nos. 10 129 046 103, 10 262 609 103, 11 015 923 103, 11 520 857 103 or 11 641 735 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.3.6

Specification

Appearance: Yellow lyophilizate

Solubility: Clear, yellow solution in water (c=20 mg/ml)

pH value (c=20 mg/ml): 7.0-8.0

Activity (+25°C, cholesterol): ≥3.0 to 4.6 U/mg lyophilizate

Specific activity: ≥40.0 U/mg protein

Protein (Biuret): No limit

Contaminants (expressed as percentage of Cholesterol Oxidase activity):

Glucose oxidase: ≤0.01

Catalase: ≤1.00

Uricase: ≤0.01

Stability: At -15 to -20°C within specification range for 12 months. Store dry.

Cat. No.	Pack Size
10 634 522 103	custom fill

Will be supplied as "Cholesterol Oxidase, Streptomyces species".
Unit of Measure is "kU".



For further processing only.

Citrate Lyase**from *Klebsiella pneumoniae*, lyophilizate**

Enzyme that catalyzes the interconversion of oxalacetate and acetate to citrate.

Application

Use Citrate lyase in tests for citric acid in combination with Malate Dehydrogenase, Catalog Nos. 11 866 109 103 and 10 200 387 103, and Lactate Dehydrogenase, Catalog Nos. 11 291 416 103, 12 235 650 103 or 10 003 557 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 4.1.3.6

Specification

Appearance: Slightly beige lyophilizate

pH value (hydrous solution): 6.5-7.5

Activity (+25 °C, citrate): ≥0.25 U/mg lyophilizate

Contaminants (expressed as percentage of Citrate Lyase activity):

Isocitrate dehydrogenase (NAD specific): ≤0.05

"NADH-oxidase": ≤0.05

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.

10 354 066 103 custom fill

Pack Size

Will be supplied as "Citrate Lyase (CL), *Aerobacter aerogenes*". Unit of Measure is "kU".

For further processing only.

Citrate Synthase**from pig heart, suspension**

Enzyme that catalyzes the formation of citrate from acetyl-CoA and oxalacetate.

Application

Use Citrate Synthase in reagents for acetic acid testing in combination with Acetate-CoA Ligase (Acetyl-CoA Synthetase), Catalog Nos. 10 128 180 103 and 10 885 568 103, and L-Malate Dehydrogenase, Catalog Nos. 11 866 109 103 and 10 200 387 103 and 10 200 387 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 2.3.3.1

Specification

Appearance: Slightly grey-brown suspension in ammonium sulfate, 3.2 mol/l; potassium phosphate, 0.02 mol/l; pH approximately 7

Specific activity (+25°C, oxaloacetic acid): ≥110 U/mg protein

Protein (Biuret): ≥10 mg/ml

Contaminants (expressed as percentage of Citrate Synthase activity):

Oxaloacetate decarboxylase: ≤0.1

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No.

10 153 605 103 custom fill

Pack Size

Will be supplied as "Citrate Synthase (CS) from Pig Heart". Unit of Measure is "g".

For further processing only.

Colipase from porcine pancreas, lyophilizate

Application

Use Colipase as a co-emulsifier in diagnostic tests for the determination lipase activity in combination with chromogenic Lipase Substrate, Catalog No. 11 034 618 103 or for the determination of triglycerides in combination with Lipase, Catalog no. 10 410 551 103.

Benefits

- Activate lipase activity in your reagent mix.
- Rely on the proven diagnostic quality of this product.

Properties

Molecular weight: Approximately 10 kD

Isoelectric point: 5.0

pH optimum: 8.8 (see figure)

Temperature dependence: See figure

pH stability: 3.5-11.5 (see figure)

Thermal stability: +25 to +80°C (see figure)

Specificity: Pancreatic colipase consists of 3 forms, colipase101 (procolipase), colipase96 and colipase85 (numbers stand for amino acid residues). colipase96 and colipase85 are trypsin digestion products of colipase101.

Specification

Appearance: White lyophilizate, TEA buffer, pH 6.0

Solubility: Clear, colorless solution in water (c=1 mg/ml)

pH value (c=1 mg/ml in water): 5.0-7.0

Protein (Lowry): 0.8±0.2 mg/mg lyophilizate

Colipase: ≥0.6 mg/mg lyophilizate

Activity (+25°C, tributyrin): ≥70 000 U/mg lyophilizate

Contaminants:

Lipase: ≤0.0005 U/mg lyophilizate

Proteases: ≤180 U/mg lyophilizate

SVD free: Corresponds to specification


pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 24 months. Store dry.


Cat. No. **Pack Size**

10 204 307 103 custom fill

Will be supplied as "Colipase from Porcine Pancreas". Unit of Measure is "g active ingredient".

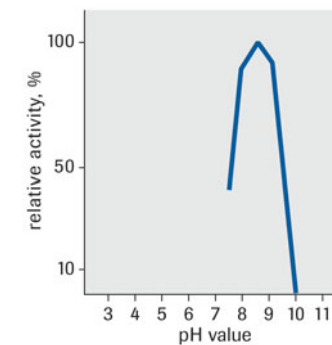


pH optimum

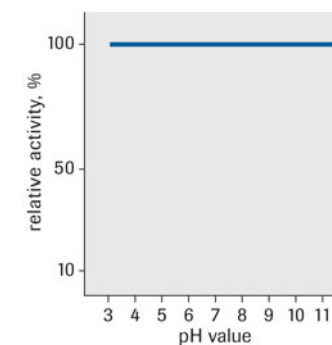


pH stability

Incubation:
25°C, 60 min
pH 3.0 – 5.0:
citrate buffer, 0.1 mol/l
pH 6.0 – 8.0
phosphate buffer,
0.1 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.1 mol/l
0.7 mg colipase/ml

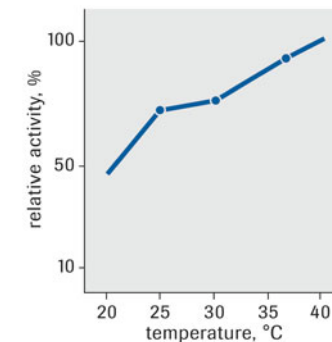

Temperature dependence

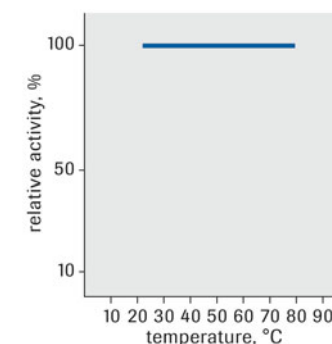

Thermal stability

Incubation:
10 min
water;
pH 4.5
0.7 mg colipase/ml

For further processing only.

Creatinase

from microorganism, lyophilizate

Hydrolase for creatinine determination that catalyzes the conversion of creatine to sarcosine and urea.

Application

Use Creatinase in diagnostic reagent for the determination of creatinine in combination with Creatininase, Catalog No. 11 865 471 103 and Sarcosine Oxidase, Catalog No. 11 378 856 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.3.3

Specification

Appearance: White to slightly yellowish lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 5.5-6.5

Activity (+25°C, creatine, POD/PAP method): ≥4 U/mg lyophilizate

Specific activity: ≥9 U/mg protein

Protein (Biuret): 0.3-0.5 mg/mg lyophilizate

Contaminants (expressed as percentage of Creatinase activity):

Creatininase: ≤0.01

Catalase: ≤2

Creatinine deaminase: ≤0.01

Proteases (casein/resorufin, 2 hours stress duration): ≤0.001

Contaminating oxidases (FOX): ≤0.001

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. **11 799 142 103** Pack Size **custom fill**

Will be supplied as "Creatinase, Microbial Lyophil. Substance". Unit of Measure is "MU".



For further processing only.

Creatininase

from *Pseudomonas* species, expressed in *E.coli*, lyophilizate

Hydrolase for creatinine determination that catalyzes the conversion of creatinine to creatine.

Application

Use Creatininase as a diagnostic reagent for the determination of creatinine in combination with Creatinase, Catalog No. 11 799 142 103 and Sarcosine Oxidase, Catalog No. 11 378 856 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.2.10

Properties

Nomenclature: Creatinine amidohydrolase

Molecular weight (gel filtration): 175 kD

Structure (SDS PAGE): 8 equal subunits (23 kD + zinc)

Isoelectric point (IEF): 4.7

Michaelis constants (Glycylglycine buffer, pH 8.0, +25°C):


Creatinine: 3×10^{-2} mol/l

Creatine: 6×10^{-2} mol/l

Inhibitors: Hg²⁺, Fe³⁺, Cu²⁺ (1 mmol/l), N-bromosuccinimide, o-phenanthroline, 4-chloromercuribenzoate. The enzyme is sensitive against photooxidation. A stable, inactive apoenzyme free of zinc can be obtained after EDTA

Cat. No. **11 865 471 103** Pack Size **custom fill**

Will be supplied as "Creatininase, Recombinant Lyo". Unit of Measure is "MU".

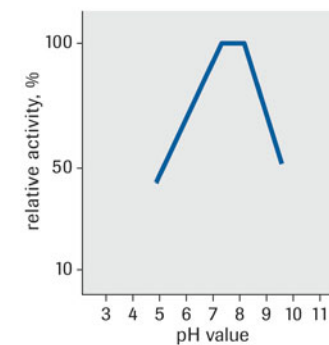

pH optimum

incubation which can be reactivated completely with Zn²⁺, Mn²⁺, Mg²⁺, Co²⁺, Fe²⁺ or Ni²⁺ (1 mmol/l).

Activators: Mn²⁺, Mg²⁺. The enzyme requires metal ions. Phenylmethylsulfonyl fluoride and iodoacetamide do not react.

pH optimum: 7.8 (see figure)

Temperature dependence: See figure

pH stability: 7.5-9.0 (see figure)

Thermal stability: Up to +65°C (see figure)

Specificity: Creatininase is specific for creatinine. It also reacts with glycohydantoin and glycohydantoin. It does not react with hydantoin and its derivatives.

Specification

Appearance: White to slightly yellowish lyophilizate

Solubility: Clear, colorless to slightly yellowish solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 7.0-8.0

Activity (+25°C, creatinine): ≥55 U/mg lyophilizate

Specific activity: ≥250 U/mg protein

Protein (Biuret): 0.10-0.35 mg/mg lyophilizate

Contaminants (expressed as percentage of Creatininase activity):

ATPase: ≤0.01

Catalase: ≤2.0

Contaminating oxidases (FOX:) ≤0.001

Creatinine deaminase: ≤0.0015

Kinase test: ≤0.01


Proteases (casein/resorufin, 2 hours stress duration): ≤0.005

Stability: At +2 to +8°C within specification range for 12 months. Protect from light.

Literature


1) A.W. Wahlefeld, J. Siedel, in H.U. Bergmeyer, Methods of Enzymatic Analysis **8**, 488 – 500 (1985)

2) J. Siedel, H. Möllering, J. Ziegenhorn, Clin. Chem. **30**, 968 (1984)

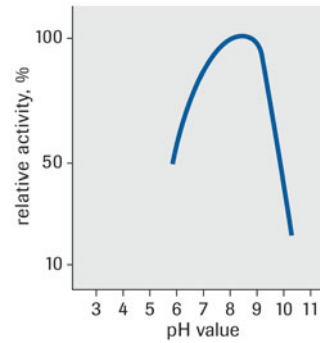


Incubation:
25°C, 120 min
pH 6.0 – 10.0:
K-phosphate
buffer, 20 mmol/l
40 U creatininase/ml

pH stability

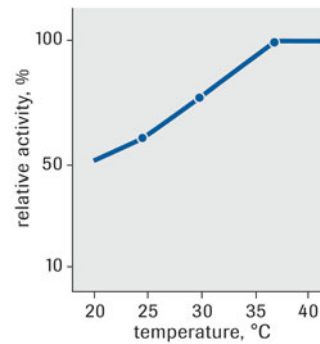

Temperature dependence

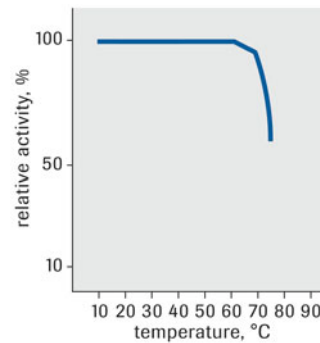

Thermal stability

Incubation:
60 min
K-phosphate
buffer, 20 mmol/l;
pH 8.0
40 U creatininase/ml

For further processing only.

Creatinine Deaminase
from *Corynebacterium lilium*, lyophilizate

Hydrolase for creatinine determination that catalyzes the conversion of creatinine to N-methylhydantoin and ammonia.

Application

Use Creatinine Deaminase in diagnostic tests for the determination of creatinine in combination with N-Carbamoylsarcosine Amidase, Catalog No. 11248 847 103, N-Methylhydantoinase (ATP-hydrolysing), Catalog No. 11 288 555 103 and Sarcosine Oxidase, Catalog No. 11 378 856 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.4.21

Cat. No.	Pack Size
11 330 764 103	custom fill

Will be supplied as "Creatinine Deiminase". Unit of Measure is "MU".

For further processing only.

Specification**Appearance:** Beige lyophilizate**Solubility:** Clear, yellowish solution in water (c=10 mg/ml)**pH value** (c=10 mg/ml in water): 8.0-9.0**Activity** (+25°C, creatinine, via N-methylhydantoin, UV): 45.0-90.0 U/mg lyophilizate**Activity** (+25°C, creatinine, via NH₃, UV): 35.0-70.0 U/mg lyophilizate**Protein** (BCA): 10-30 mg/100 mg lyophilizate**Contaminants** (expressed as percentage of Creatinine Desaminase activity)

ATPase: ≤0.1

Creatinase: ≤0.013

Creatininase: ≤0.01

Catalase: ≤10.0

Urease: ≤0.007

Uricase: ≤0.01

NH₃: ≤0.01µg/U**Stability:** At +2 to +8°C within specification range for 12 months. Store dry.

Protect from light.

**Formate Dehydrogenase
from yeast, lyophilizate**

Dehydrogenase that catalyzes the interconversion of formate to carbon dioxide.

Application

Use Formate Dehydrogenase in diagnostic tests for the determination of oxalate in combination with Oxalate Oxidase, Catalog No. 10 570 524 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.2.1.2

Specification**Appearance:** White lyophilizate**pH value** (c=10 mg/ml in water): Approximately 7.5**Activity** (+25°C, formiate): ≥0.40 U/mg lyophilizate**Specific activity:** ≥3.0 U/mg protein**Contaminants** (expressed as percentage of Formate Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.05

Lactate dehydrogenase: ≤0.05

Malate dehydrogenase: ≤0.1

Stability: At +2 to +8°C within specification range for 12 months. Store dry.**Galactose 1-Dehydrogenase
from *E.coli* overproducer, lyophilizate**

Dehydrogenase that catalyzes the oxidation of galactose to D-galactono-1,4-lactone.

Application

Use Galactose 1-Dehydrogenase in diagnostic tests for the determination of total galactose.

Benefits

- Rely on the proven diagnostic quality of this recombinant enzyme.

84 EC 1.1.1.48

Cat. No. Pack Size**10 204 226 103** custom fill

Will be supplied as "Formate Dehydrogenase from Yeast". Unit of Measure is "kU".

For further processing only.

Cat. No. Pack Size**11 290 983 103** custom fill

Will be supplied as "b-Galactose Dehydrogenase S". Unit of Measure is "kU".

For further processing only.

Specification

Appearance: White lyophilizate

Specific activity (+25°C, galactose): ≥50 U/mg protein

Protein (Biuret): ≥0.3-0.7 mg/mg lyophilizate

Contaminants (expressed as percentage of Galactose 1-Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.01

β-Galactosidase: ≤0.01

Glutamate dehydrogenase (standard): ≤0.5

Lactate dehydrogenase: ≤0.1

Malate dehydrogenase: ≤1.0

"NADH-oxidase": ≤0.05

Stability: At +2 to +8°C within specification range for 12 months.

Galactose 1-Dehydrogenase from *E.coli* overproducer, suspension

Dehydrogenase that catalyzes the oxidation of galactose to D-galactono-1,4-lactone.

Application

Use Galactose 1-Dehydrogenase in diagnostic tests for the determination of total galactose.

Benefits

- Apply this ready-to-use recombinant enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.48

Specification

Appearance: White suspension in ammonium sulfate solution, 3.2 mol/l, pH approximately 6

Specific activity (+25°C, D-galactose): ≥100 U/mg protein

Protein (Biuret): ≥1 mg/ml

Contaminants (expressed as percentage of Galactose 1-Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.01

β-Galactosidase: ≤0.01

Lactate dehydrogenase: ≤0.1

Malate dehydrogenase: ≤1.0

"NADH-oxidase": ≤0.05

Stability: At +2 to +8°C within specification range for 12 months.

Galactose 1-Dehydrogenase from *Pseudomonas fluorescens*, suspension

Dehydrogenase that catalyzes the oxidation of galactose to D-galactono-1,4-lactone.

Application

Use Galactose 1-Dehydrogenase in diagnostic tests for the determination of total galactose.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

Cat. No.	Pack Size
10 633 313 103	custom fill

Will be supplied as "b-Galactose Dehydrogenase S. E. coli". Unit of Measure is "kU".

For further processing only.

Cat. No.	Pack Size
10 150 959 103	custom fill

Will be supplied as "b-Gal-DH from Pseudomonas fluorescens". Unit of Measure is "kU".

For further processing only.

EC 1.1.1.48

Specification

Appearance: White suspension in ammonium sulfate solution, 3.2 mol/l; EDTA, 1 mmol/l; pH approximately 6

Specific activity (+25°C, D-galactose): ≥5 U/mg protein

Protein (Biuret): 5±0.5 mg/ml

Contaminants (expressed as percentage of Galactose 1-Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.01

β-Galactosidase: ≤0.01

Lactate dehydrogenase: ≤0.5

"NADH-oxidase": ≤0.1

Stability: At +2 to +8°C within specification range for 12 months.

Glucose Oxidase (GOD), Grade I from *Aspergillus niger* overproducer, lyophilizate

Oxidoreductase that catalyzes the conversion of D-glucose to D-glucono-1,5-lactone which hydrolyzes spontaneously to gluconate.

Application

Use Glucose Oxidase (GOD), Grade I for the determination of α-amylase and D-glucose or O₂.

Benefits

- Take advantage of the tested low conductivity.
- Rely on the proven diagnostic quality of this product.

EC 1.1.3.4

Properties

Nomenclature: β-D-glucose: oxygen 1-oxidoreductase

Molecular weight: 79 kD

Isoelectric point: 4.3

Michaelis constants (Glucose):

Acetate buffer, pH 5.0, +25°C: 3.6×10^{-2} mol/l

Potassium phosphate buffer, 0.2 mol/l, pH 7.5, +25°C: 4.8×10^{-2} mol/l

Inhibitors: Ag⁺, Hg²⁺, Cu²⁺, 4-chloromercuribenzoate, D-arabinose (50%).

FAD binding is inhibited by several nucleotides.

pH optimum: 7.0 (see figure)

Temperature dependence: See figure

pH stability: See figure

Thermal stability: See figure

Specificity: Glucose oxidase is specific for β-D-glucose. O₂ can be replaced by hydrogen acceptors such as 2,6-dichlorophenol indophenol.

Specification

Appearance: Yellowish lyophilizate

Conductivity (1%, w/v): ≤250 μS/cm

Activity (+25°C, glucose): ≥300 U/mg lyophilizate

Contaminants (expressed as percentage of Glucose Oxidase activity):

Amylase: ≤0.01


Catalase: ≤0.5


Saccharase: ≤0.01

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No.	Pack Size
12 158 566 103	custom fill

Will be supplied as "GOD, RG I, rec., Lyo.". Unit of Measure is "MU".

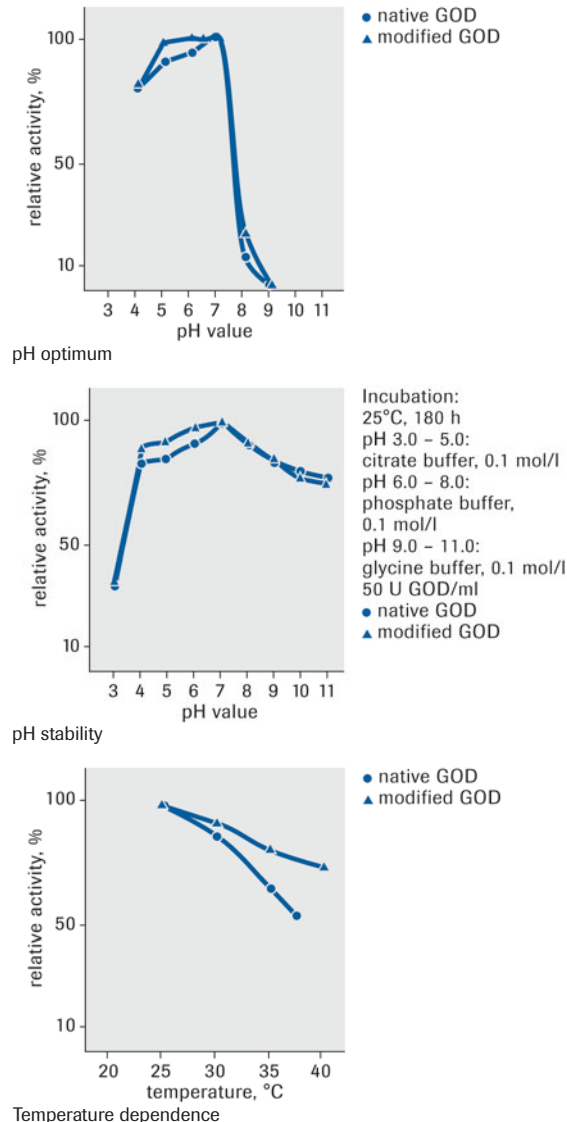

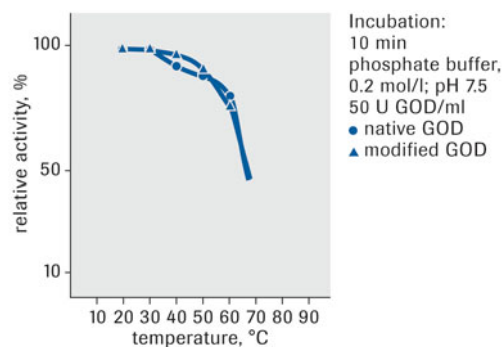

Thermal stability

For further processing only.

Glucose Oxidase (GOD), Grade II from *Aspergillus niger* overproducer, lyophilizate

Oxidoreductase that catalyzes the conversion of D-glucose to D-glucono-1,5-lactone which hydrolyzes spontaneously to gluconate.

Application

Use Glucose Oxidase (GOD), Grade II for the determination of α -amylase and D-glucose or O_2 .

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.3.4

Properties

Nomenclature: β -D-glucose: oxygen 1-oxidoreductase

Molecular weight: 79 kD

Isoelectric point: 4.3

Michaelis constants (Glucose):

Acetate buffer, pH 5.0, +25°C: 3.6×10^{-2} mol/l

Potassium phosphate buffer, 0.2 mol/l, pH 7.5, +25°C: 4.8×10^{-2} mol/l

Inhibitors: Ag^+ , Hg^{2+} , Cu^{2+} , 4-chloromercuribenzoate, D-arabinose (50%).

FAD binding is inhibited by several nucleotides.

pH optimum: 7.0 (see figure)

Temperature dependence: See figure

pH stability: See figure

Thermal stability: See figure

Specificity: Glucose oxidase is specific for β -D-glucose. O_2 can be replaced by hydrogen acceptors such as 2,6-dichlorophenol indophenol.

Specification

Appearance: Yellow brown lyophilizate

Solubility: Clear, yellow solution in phosphate buffer, 0.1 mol/l, pH 7.0 (c=5 mg/ml)

pH value (c=10 mg/ml in water): 6.8-7.8

Protein (Pierce): No limit

Activity (+25°C, glucose): ≥ 250 U/mg lyophilizate

Contaminants (expressed as percentage of Glucose Oxidase activity):

Amylase: ≤ 0.1

Catalase: ≤ 5 U/mg lyophilizate

Saccharase: ≤ 0.1

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No.	Pack Size
11 939 998 103	custom fill

Will be supplied as "GOD, rec., Lyo.". Unit of Measure is "MU".
For further processing only.

Glucose Oxidase (GOD), chemically modified

from *Aspergillus niger* overproducer, lyophilizate

Oxidoreductase that catalyzes the conversion of D-glucose to D-glucono-1,5-lactone which hydrolyzes spontaneously to gluconate.

Application

Use Glucose Oxidase (GOD), chemically modified for the determination of α -amylase and D-glucose or O₂.

Benefits

- Take advantage of the enhanced liquid stability.
- Rely on the proven diagnostic quality of this product.

EC 1.1.3.4

Properties

Nomenclature: β -D-glucose:oxygen 1-oxidoreductase

Molecular weight: 79 kD

Isoelectric point: 4.3

Michaelis constants (Glucose):

Acetate buffer, pH 5.0, +25°C: 3.6×10^{-2} mol/l

Potassium phosphate buffer, 0.2 mol/l, pH 7.5, +25°C: 4.8×10^{-2} mol/l

Inhibitors: Ag⁺, Hg²⁺, Cu²⁺, 4-chloromercuribenzoate, D-arabinose (50%).

FAD binding is inhibited by several nucleotides.

pH optimum: 7.0 (see figure)

Temperature dependence: See figure

pH stability: See figure

Thermal stability: See figure

Specificity: Glucose oxidase is specific for β -D-glucose. O₂ can be replaced by hydrogen acceptors such as 2,6-dichlorophenol indophenol.

Remark: The modified enzyme is especially suited for liquid stable applications with extended shelf life requirements.

Specification

Appearance: Yellowish white lyophilizate

pH value (c=40 mg/ml in water): 6.5-7.5

Activity (+25°C, glucose): ≥ 20 U/mg lyophilizate

Contaminants (expressed as percentage of Glucose Oxidase activity):

Catalase: ≤ 20 U/mg lyophilizate

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

11 485 938 103 custom fill

Will be supplied as "GOD, Asp.niger, Bound to Dextran". Unit of Measure is "MU".

For further processing only.

Glucose-6-phosphate Dehydrogenase (G6P-DH)

from *Leuconostoc mesenteroides*, expressed in *E. coli*, lyophilizate

Dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase for the determination of blood glucose or creatine kinase.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Properties

Nomenclature: D-glucose-6-phosphate:NAD(P)⁺ 1-oxidoreductase

Molecular weight: 110 kD (1) (2 identical subunits 55,000 D)

Isoelectric point: pH 4.6

Michaelis constants (Tris: 0.1 mol/l; pH 7.8, +25°C):

NAD: 1.4×10^{-4} mmol/l

NADP: 3.7×10^{-5} mmol/l

Glucose-6-P: 3.7×10^{-4} mmol/l (NAD as coenzyme)

Glucose-6-P: 2.0×10^{-4} mmol/l (NADP as coenzyme)

Activators/inhibitors:

Phosphate, 5 mmol/l: 100% (NAD), 80% (NADP)

Phosphate, 50 mmol/l: 100% (NAD), 80% (NADP)

Without Mg²⁺: 90% (NAD), 80% (NADP)

Mg²⁺, 3 mmol/l: 100% (NAD), 100% (NADP)

Mg²⁺, 30 mmol/l: 100% (NAD), 100% (NADP)

HCO³⁻, 3 mmol/l: 100% (NAD), 100% (NADP)

Inhibitors: NADPH is a competitive inhibitor in the NAD-dependent reaction. Unlike the yeast enzyme, myristic acid, dehydroepiandrosterone and palmitoyl CoA do not inhibit.

pH optimum: 7.8 (see figure)

Temperature dependence: See figure

pH stability: 5.0-10.0 (see figure)

Thermal stability: Up to +40°C for native G6P-DH, up to +50°C for modified G6P-DH (see figure)

Buffer stability: Temperature stability can be significantly improved by the increase of ionic strength (see figure).

Stability of the lyophilizate: 100% residual activity after 3 weeks at +35°C

Specificity: G6P-DH is highly specific for glucose-6-phosphate and does not react with fructose-6-P, fructose-1,6-P₂ or glucose-1-P. 2-Deoxyglucose-6-P is slowly oxidized with NAD (5%) and with NADP (4%).

Specification

Appearance: White or slightly yellowish lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml, water): 6.5-7.5

Activity (+25°C, glucose-6-P, NAD): ≥600 U/mg lyophilizate

Specific activity (+25°C): ≥800 U/mg protein

Activity (+30°C): ≥750 U/mg lyophilizate

Activity (+37°C): ≥1,000 U/mg lyophilizate

Activity (+25°C, glucose-6-P, NADP): No limit

Protein (Biuret): 0.7-0.9 mg/mg lyophilizate

Contaminants (expressed as percentage of Glucose-6-phosphate dehydrogenase activity):


ATPase: ≤0.02

Creatine kinase: ≤0.001


Cat. No. Pack Size

11 293 206 103 custom fill


Will be supplied as "G6P-DH, rec., Lyo". Unit of Measure is "MU".

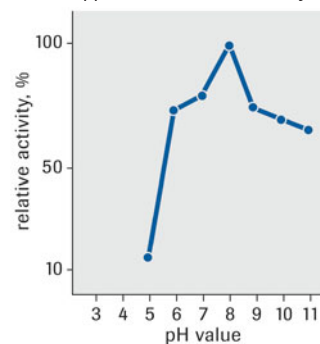


pH optimum

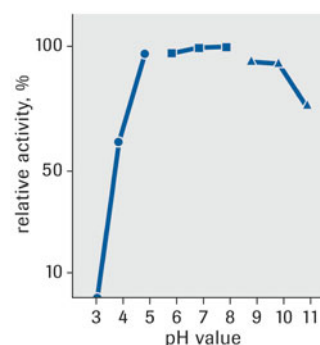

pH stability

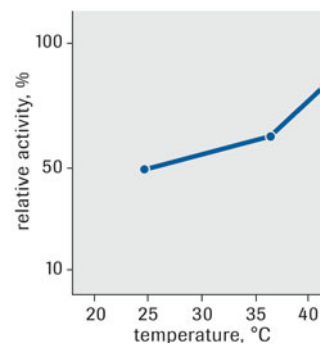

Temperature dependence

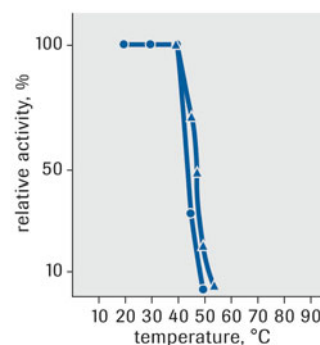


Thermal stability

Incubation:
25°C, 180 min
● pH 3.0 – 5.0:
citrate buffer, 0.1 mol/l
■ pH 6.0 – 8.0:
phosphate buffer,
0.1 mol/l
▲ pH 9.0 – 11.0:
glycine buffer, 0.1 mol/l
500 U G6P-DH/ml

Incubation:
20 min
phosphate buffer,
0.02 mol/l; pH 7.5
50 U G6P-DH/ml
● native G6P-DH
▲ modified G6P-DH

Glutamate dehydrogenase: ≤ 0.01
 Glutathione reductase: ≤ 0.001
 Hexokinase and Glucose dehydrogenase: ≤ 0.05
 Myokinase: ≤ 0.001
 "NADH oxidase": ≤ 0.02
 "NADPH oxidase": ≤ 0.0005
 6-Phosphogluconate dehydrogenase: ≤ 0.001
 Phosphoglucose isomerase : ≤ 0.01
 Phosphoglucomutase: ≤ 0.001

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

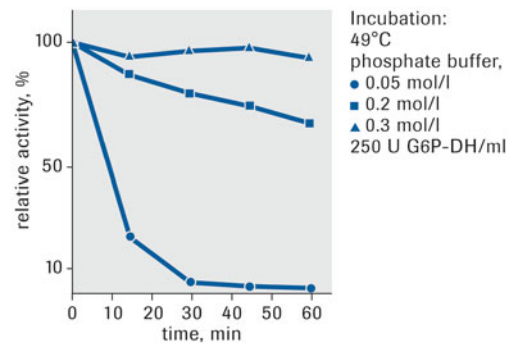

Buffer stability

For further processing only.

Glucose-6-phosphate Dehydrogenase (G6P-DH), chemically modified

from *Leuconostoc mesenteroides*, expressed in *E. coli*, lyophilizate

Recombinant dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase for the determination of blood glucose or creatine kinase.

Benefits

- Take advantage of the improved stability in liquid reagents.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Properties

Nomenclature: D-glucose-6-phosphate:NAD(P)⁺ 1-oxidoreductase

Molecular weight: 110 kD (2 identical subunits 55,000 D)

Isoelectric point: pH 4.6

Michaelis constants (Tris: 0.1 mol/l; pH 7.8, +25°C):

NAD: 1.4×10^{-4} mmol/l

NADP: 3.7×10^{-5} mmol/l

Glucose-6-P: 3.7×10^{-4} mmol/l (NAD as coenzyme)

Glucose-6-P: 2.0×10^{-4} mmol/l (NAD as coenzyme)

Activators/inhibitors:

Phosphate, 5 mmol/l: 100% (NAD), 80% (NADP)

Phosphate, 50 mmol/l: 100% (NAD), 80% (NADP)

Without Mg²⁺: 90% (NAD), 80% (NADP)

Mg²⁺, 3 mmol/l: 100% (NAD), 100% (NADP)

Mg²⁺, 30 mmol/l: 100% (NAD), 100% (NADP)

HCO³⁻, 3 mmol/l: 100% (NAD), 100% (NADP)

Inhibitors: NADPH is a competitive inhibitor in the NAD-dependent reaction.

Unlike the yeast enzyme, myristic acid, dehydroepiandrosterone and palmitoyl CoA do not inhibit.

pH optimum: 7.8 (see figure)

Temperature dependence: See figure

pH stability: 5.0-10.0 (see figure)

Thermal stability: Up to +40°C for native G6P-DH, up to +50°C for modified G6P-DH (see figure)

Buffer stability: Temperature stability can be significantly improved by the

Cat. No.

Pack Size

11 389 343 103 custom fill

Will be supplied as "G6P-DH, rec., Lyo., mod.". Unit of Measure is "MU".

For further processing only.

increase of ionic strength (see figure).

Stability of the lyophilizate: 100% residual activity after 3 weeks at +35°C

Specificity: G6P-DH is highly specific for glucose-6-phosphate and does not react with fructose-6-P, fructose-1,6-P2 or glucose-1P. 2-Deoxyglucose-6-P is slowly oxidized with NAD (5%) and with NADP (4%).

Remark: The modified enzyme is especially suited for liquid stable applications with extended shelf life requirements.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=40 mg/ml)

pH value (c=40 mg/ml in water): 6.5-7.5

Activity (+25°C, glucose-6-P, NAD): ≥30 U/mg lyophilizate

Activity (+30°C): ≥39 U/mg lyophilizate

Activity (+37°C): ≥54 U/mg lyophilizate

Contaminants (expressed as percentage of Glucose-6-phosphate Dehydrogenase activity):

ATPase: ≤0.02

Creatine kinase: ≤0.001

Glutamate dehydrogenase: ≤0.01

Glutathione reductase: ≤0.001

Hexokinase and Glucose dehydrogenase: ≤0.05

Myokinase: ≤0.05

"NADH oxidase": ≤0.02

"NADPH oxidase": ≤0.0005

6-Phosphogluconate dehydrogenase: ≤0.001

Phosphoglucose isomerase : ≤0.01

Phosphoglucomutase: ≤0.001

Glucose: ≤0.3 µg/mg lyophilizate

Stability: At +2 to +8°C within specification range for 18 months. Store dry.

Glucose-6-phosphate Dehydrogenase (G6P-DH)

from *Leuconostoc mesenteroides*, expressed in *E. coli*, solution

Recombinant dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase for the determination of blood glucose or creatine kinase.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Properties

Nomenclature: D-glucose-6-phosphate:NAD(P)⁺ 1-oxidoreductase

Molecular weight: 110 kD (1) (2 identical subunits 55,000 D)

Isoelectric point: pH 4.6

Michaelis constants (Tris: 0.1 mol/l; pH 7.8, +25°C):

NAD: 1.4 x 10⁻⁴ mmol/l

NADP: 3.7 x 10⁻⁵ mmol/l

Glucose-6-P: 3.7 x 10⁻⁴ mmol/l (NAD as coenzyme)

Glucose-6-P: 2.0 x 10⁻⁴ mmol/l (NAD as coenzyme)

Activators/inhibitors:

Cat. No.	Pack Size
11 650 742 103	custom fill

Will be supplied as "G6P-DH, Recombinant (E. coli)". Unit of Measure is "MU".

For further processing only.

Phosphate, 5 mmol/l: 100% (NAD), 80% (NADP)
 Phosphate, 50 mmol/l: 100% (NAD), 80% (NADP)
 Without Mg²⁺: 90% (NAD), 80% (NADP)
 Mg²⁺, 3 mmol/l: 100% (NAD), 100% (NADP)
 Mg²⁺, 30 mmol/l: 100% (NAD), 100% (NADP)
 HCO³⁻, 3 mmol/l: 100% (NAD), 100% (NADP)

Inhibitors: NADPH is a competitive inhibitor in the NAD-dependent reaction. Unlike the yeast enzyme, myristic acid, dehydroepiandrosterone and palmitoyl CoA do not inhibit.

pH optimum: 7.8 (see figure)

Temperature dependence: See figure

pH stability: 5.0-10.0 (see figure)

Thermal stability: Up to +40°C for native G6P-DH, up to +50°C for modified G6P-DH (see figure)

Buffer stability: Temperature stability can be significantly improved by the increase of ionic strength (see figure).

Stability of the lyophilizate: 100% residual activity after 3 weeks at +35°C

Specificity: G6P-DH is highly specific for glucose-6-phosphate and does not react with fructose-6-P, fructose-1,6-P₂ or glucose-1P. 2-Deoxyglucose-6-P is slowly oxidized with NAD (5%) and with NADP (4%).

Specification

Appearance: Clear, yellowish solution in glycerol

pH value: 6.0-7.0

Activity (+25°C, glucose-6-P): ≥2,500 U/ml

Activity (+30°C): ≥3,000 U/ml

Contaminants (expressed as percentage of Glucose-6-phosphate Dehydrogenase activity):

ATPase: ≤0.0200

Creatine kinase: ≤0.001

Glutamate dehydrogenase: ≤0.01

Glutathione reductase: ≤0.001

Hexokinase and Glucose dehydrogenase: ≤0.05

Myokinase: ≤0.01

"NADH oxidase": ≤0.02

"NADPH oxidase": ≤0.0005

Phosphoglucomutase: ≤0.001

6-Phosphogluconate dehydrogenase (NAD): ≤0.001

6-Phosphogluconate dehydrogenase (NADP): ≤0.001

Phosphoglucose isomerase : ≤0.01

Glycerol (enzymatically): 45-55% (v/v)

Stability: At +2 to +8°C within specification range for 12 months.

Glucose-6-phosphate Dehydrogenase (G6P-DH)

from *Leuconostoc mesenteroides*, lyophilizate

Dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase for the determination of blood glucose or creatine kinase.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Cat. No.	Pack Size
10 186 783 103	custom fill

Will be supplied as "G6P-DH from *Leuconostoc mesenteroides*".
 Unit of Measure is "kU".
 For further processing only.

Specification

Appearance: White or slightly yellowish lyophilizate
Solubility: Clear, colorless solution in water (c=10mg/ml)
pH value (c=10mg/ml in water): 6.5-7.5
Activity (+25°C, glucose-6-P): ≥400 U/mg lyophilizate
Activity (+30°C): ≥520 U/mg lyophilizate
Activity (+37°C): ≥720 U/mg lyophilizate
Contaminants (expressed as percentage of Glucose-6-phosphate Dehydrogenase activity):
 ATPase: ≤0.02
 Creatine kinase: ≤0.001
 Glutamate dehydrogenase: ≤0.01
 Glutathione reductase: ≤0.001
 Hexokinase and Glucose dehydrogenase: ≤0.05
 Myokinase: ≤0.05
 "NADH oxidase": ≤0.02
 "NADPH oxidase": ≤0.0005
 Phosphoglucomutase: ≤0.001
 6-Phosphogluconate dehydrogenase (NAD): ≤0.001
 6-Phosphogluconate dehydrogenase (NADP): ≤0.001
 Phosphoglucose isomerase: ≤0.01
Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Glucose-6-phosphate Dehydrogenase (G6P-DH)

from *Leuconostoc mesenteroides*, expressed in *E. coli*, reduced phosphate, lyophilizate

Recombinant dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase, reduced phosphate for the determination of inorganic phosphate in an colorimetric enzymatic reaction.

Benefits

- Take advantage of the strongly reduced concentration of phosphate.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Specification

Appearance: White to slightly yellowish lyophilizate
Solubility: Clear, colorless solution in water (c=10 mg/ml)
pH value (c=10mg/ml in water): 6.5-7.5
Activity (+25°C, glucose-6-P): ≥600 U/mg lyophilizate
Activity (+30°C): ≥750 U/mg lyophilizate
Activity (+37°C): ≥1,000 U/mg lyophilizate
Specific activity (+25°C, glucose-6-P): ≥800 U/mg
Protein (Biuret): 0.8-1.0 mg/mg lyophilizate
Contaminants (expressed as percentage of Glucose-6-phosphate Dehydrogenase activity):
 ATPase: ≤0.00005
 Creatine kinase: ≤0.001
 Glutamate dehydrogenase: ≤0.01
 Glutathione reductase: ≤0.001
 α-Glucosidase: ≤0.00010
 Hexokinase and glucose dehydrogenase: ≤0.05
 Myokinase: ≤0.001
 "NADH oxidase": ≤0.02

Cat. No.	Pack Size
11 650 734 103	custom fill

Unit of Measure is "MU".
For further processing only.

"NADPH oxidase": ≤ 0.0005
 Phosphoglucomutase: ≤ 0.001
 6-Phosphogluconate dehydrogenase (NAD): ≤ 0.001
 6-Phosphogluconate dehydrogenase (NADP): ≤ 0.0001
 Phosphoglucose isomerase: ≤ 0.01
Phosphate (as P_i): $\leq 10 \mu\text{g}/\text{mg}$ lyophilizate
Stability: At +2 to +8°C within specification range for 18 months. Store dry.

Glucose-6-phosphate Dehydrogenase (G6P-DH)

from *Leuconostoc mesenteroides*, suspension

Dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase for the determination of blood glucose or creatine kinase.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Specification

Appearance: Yellowish suspension in ammonium sulfate

pH value: 5.5-6.5

Specific activity (+25°C, glucose-6-P): $\geq 550 \text{ U}/\text{mg}$

Specific activity (+30°C): $\geq 650 \text{ U}/\text{mg}$

Specific activity (+37°C): $\geq 800 \text{ U}/\text{mg}$

Protein (Biuret): $\geq 5 \text{ mg}/\text{ml}$

Contaminants (expressed as percentage of Glucose-6-phosphate Dehydrogenase activity):

ATPase: ≤ 0.02

Creatine kinase: ≤ 0.001

Glutamate dehydrogenase: ≤ 0.01

Glutathione reductase (NADH): ≤ 0.001

Hexokinase and Glucose dehydrogenase: ≤ 0.05

Myokinase: ≤ 0.05

"NADH oxidase": ≤ 0.02

"NADPH oxidase": ≤ 0.0005

Phosphoglucomutase: ≤ 0.001

6-Phosphogluconate dehydrogenase (NAD): ≤ 0.001

6-Phosphogluconate dehydrogenase (NADP): ≤ 0.001

Phosphoglucose isomerase: ≤ 0.01

Stability: At +2 to +8°C within specification range for 18 months.

Cat. No. **Pack Size**

10 128 171 103 custom fill

Will be supplied as "G6P-DH from *Leuconostoc mesenteroides*".
 Unit of Measure is "MU".

For further processing only.

Glucose-6-phosphate Dehydrogenase (G6P-DH)

from yeast, lyophilizate

Dehydrogenase, that catalyzes the oxidation of Glucose-6-phosphate.

Application

Use Glucose-6-phosphate Dehydrogenase for the determination of blood glucose or creatine kinase.

Cat. No. **Pack Size**

10 190 454 103 custom fill

Will be supplied as "Glucose-6-phosphate Dehydrogenase, Yeast".
 Unit of Measure is "MU".

For further processing only.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.49

Specification

Appearance: White to slightly yellowish lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml): 6.0-7.0

Activity (+25°C, glucose-6-P): ≥15.0 U/mg lyophilizate

Contaminants (expressed as percentage of Glucose-6-phosphate Dehydrogenase activity):

Creatine kinase: ≤0.001

Glutathione reductase: ≤0.05

Hexokinase: ≤0.02

Phosphoglucomutase: ≤0.01

6-Phosphogluconate dehydrogenase: ≤0.01

Phosphoglucose isomerase: 0.002

Bioburden: ≤10,000 CFU/g

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

**Glucose-6-phosphate Isomerase
from yeast, suspension**

Isomerase, that catalyzes the conversion of glucose-6-phosphate into fructose 6-phosphate.

Application

Use Glucose-6-phosphate Isomerase for the isomerization of ketoses to al-doses and can be used for the determination of fructose-6-phosphate.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 5.3.1.9

Specification

Appearance: White suspension in ammonium sulfate, 3.2 mol/l

Specific activity (+25°C, fructose-6-P): ≥350 U/mg protein

Protein (Biuret): 10±1 mg/ml

Contaminants (expressed as percentage of Glucose-6-phosphate Isomerase activity):

Fructose-6-phosphate kinase: ≤0.01

β-Fructosidase: ≤0.2

Glutathione reductase: ≤0.01

Phosphoglucomutase: ≤0.01

6-Phosphogluconate dehydrogenase: ≤0.01

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No.	Pack Size
10 154 334 103	custom fill

Will be supplied as "Phosphoglucose Isomerase (PGI) from Yeast".
Unit of Measure is "MU".
For further processing only.

α -Glucosidase from yeast overproducer, multifunctional, lyophilizate

Recombinant glucosidase, that hydrolyzes 1,4-linked α -D-glucose residues with release of α -D-glucose.

Application

Use α -Glucosidase in diagnostic tests for the determination of α -amylase and pancreatic α -amylase activity according to the recommendations of the International Federation of Clinical Chemistry and Laboratory Medicine (IFCC). It catalyzes the 100% liberation of p-nitrophenol residues from the amylase substrate EPS (Catalog No. 10 880 078 103) once it has been cleaved by α -amylase.

Benefits

Use this recombinant enzyme in your amylase reagent mix and rely on the 100% chromophore liberation and the proven diagnostic quality of this product.

EC 3.2.1.20

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml): 6.8-7.4

Activity (+37°C, 4-NP- α -D-glucoside): ≥ 60 U/mg lyophilizate


Specific activity: ≥ 130 U/mg protein

Protein (Biuret): 25-45 mg/100 mg lyophilizate


Contaminants (expressed as percentage of α -Glucosidase activity):

α -Amylase: ≤ 0.00001

Stability: At -15 to -20°C within specification range for 12 months. Store dry.

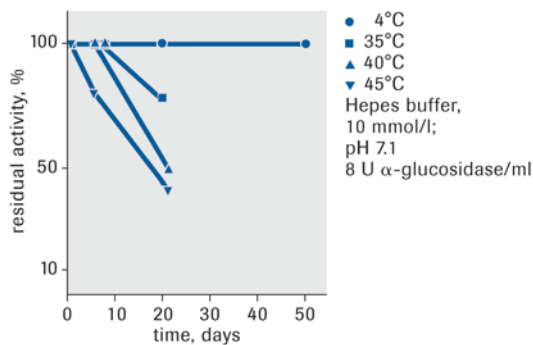


Stability in solution

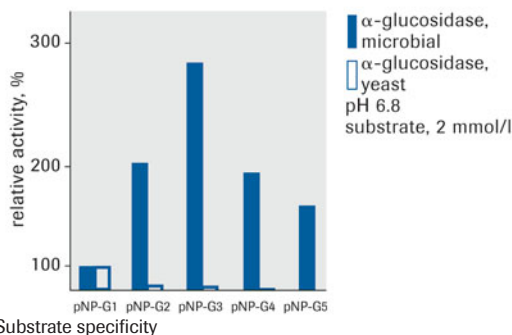


Substrate specificity

Cat. No. 11 626 329 103
Pack Size custom fill


Will be supplied as "a-Glucosidase Multifunctional". Unit of Measure is "MU".

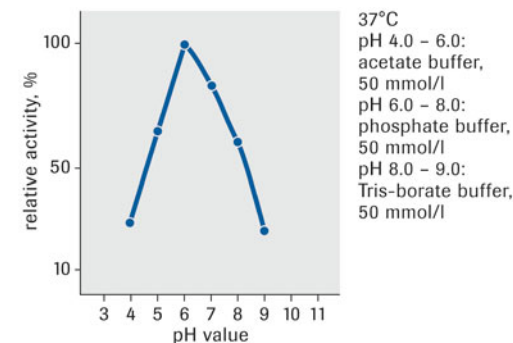


pH optimum

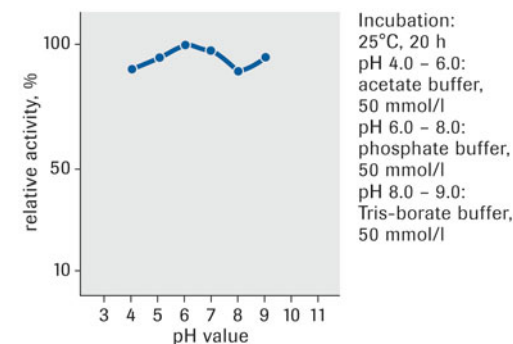

pH stability

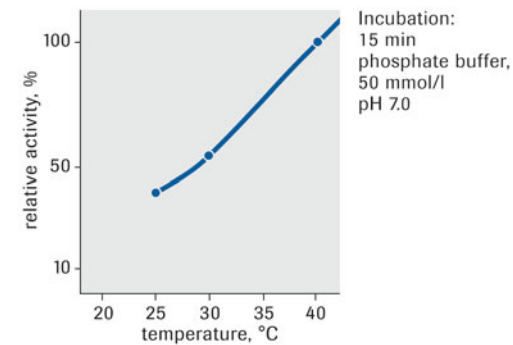

Temperature dependence

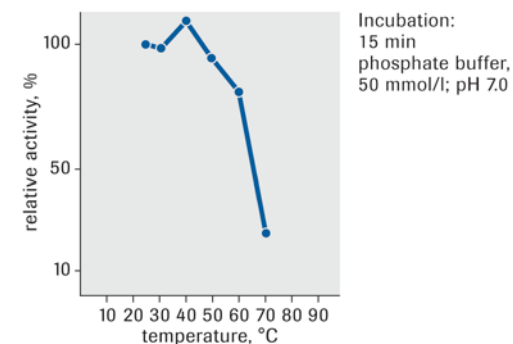

Thermal stability

For further processing only.

β -Glucuronidase from *E.coli*, solution

Hydrolase that cleaves β -linked terminal glucuronic acid.

Application

Use β -Glucuronidase in reagents for drug monitoring and doping analysis where it catalyzes the hydrolysis of steroid conjugates to detect various steroids in urine.

Benefits

- Take advantage of the high specific activity of the *E.coli* enzyme and its great affinity for the various β -glucuronides.
- Save time. Develop your procedure without the need of cleaning up the reaction and buffering the urine.

EC 3.2.1.31

Specification

Appearance: Clear, colorless solution

Specific activity (+25°C, 4-NP-glucuronide): ≥ 80 U/mg protein

Specific activity (+37°C, 4-NP-glucuronide): ≥ 140 U/mg protein

Protein (Biuret): ≥ 0.5 mg/ml solution

Stability: At +2 to +8°C within specification range for 18 months.

Cat. No. **Pack Size**

03 708 446 103 custom fill

Will be supplied as "beta-Glucuronidase E.coli K12 Glycerol". Unit of Measure is "MU".

For further processing only.

Glutamate Dehydrogenase (NAD(P)) from *E.coli* overproducer, lyophilizate

Recombinant glutamate dehydrogenase.

Application

Use recombinant Glutamate Dehydrogenase in diagnostic tests for the determination of ammonia, urea, L-glutamate, glutamate pyruvate transaminase and leucine aminopeptidase.

Benefits

- Rely on the proven diagnostic quality of this product.

EC 1.4.1.3

Properties

Nomenclature: L-glutamate:NAD(P)⁺ oxidoreductase (deaminating)

Molecular weight: ~ 2 200 kD for the associated enzyme with 8 subunits; 280 kD for one subunit.

Michaelis constants (Tris buffer, pH 8.0, +23°C):

L-glutamate: 1.8×10^{-3} mol/l

NADP: 4.7×10^{-5} mol/l

α -ketoglutarate: 7.0×10^{-4} mol/l

NH₄⁺: 3.2×10^{-3} mol/l

NADPH: 2.6×10^{-5} mol/l

Km values for NAD or NADH are difficult to obtain due to their inhibitory action.

Inhibitors: 4-chloromercuribenzoate, Na₂S, diethyldithiocarbamate, 1,10-phenanthroline,

8-hydroxyquinoline, NaN₃, thyroxine, heparin, sulfonylcarbamides, Cu²⁺, Hg²⁺, Ag²⁺, Fe³⁺, Zn²⁺, K⁺, PO₄²⁻, NO₃⁻

Activators: Thioglycolic acid, b-mercaptoethylamine, EDTA, α , α' -dipyridyl

pH optimum: 8.0 (see figure)

Temperature dependence: See figure

Cat. No. **Pack Size**

11 745 727 103 custom fill

Will be supplied as "GIDH, Lyo., rec.". Unit of Measure is "MU".

For further processing only.

pH stability: 5.5-6.5 (see figure)

Thermal stability: Up to +60°C (see figure)

Specificity: The oxidation of L-glutamate is stimulated by ADP and inhibited by GTP. In contrast, the oxidation of alanine, leucine, isoleucine, methionine, valine, norleucine, norvaline and 2-aminobutyrate is stimulated by GTP and inhibited by ADP.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=20 mg/ml)

pH value (c=20 mg/ml in water): 6.5-7.5

Activity (+25°C, α-ketoglutarat): ≥80 U/mg lyophilizate

Contaminants (expressed as percentage of Glutamate Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.005

Lactate dehydrogenase: ≤0.005

Malate dehydrogenase: ≤0.005

"NADH-Oxidase": ≤0.005

NH₄⁺: ≤0.05 µg/mg lyophilizate

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Glutamate Dehydrogenase (NAD(P)) from beef liver, lyophilizate

Dehydrogenase that catalyzes the conversion of glutamate to α-ketoglutarate.

Application

Use Glutamate Dehydrogenase in diagnostic tests for the determination of ammonia, urea, L-glutamate, glutamate pyruvate transaminase and leucine aminopeptidase.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.4.1.3

Properties

Nomenclature: L-glutamate:NAD(P)⁺ oxidoreductase (deaminating)

Molecular weight: ~2 200 kD for the associated enzyme with 8 subunits; 280 kD for one subunit.

Michaelis constants (Tris buffer, pH 8.0, +23°C):

L-glutamate: 1.8×10^{-3} mol/l

NADP: 4.7×10^{-5} mol/l

α-ketoglutarate: 7.0×10^{-4} mol/l

NH₄⁺: 3.2×10^{-3} mol/l

NADPH: 2.6×10^{-5} mol/l

K_m values for NAD or NADH are difficult to obtain due to their inhibitory action.

Inhibitors: 4-chloromercuribenzoate, Na₂S, diethyldithiocarbamate, 1,10-phenanthroline, 8-hydroxyquinoline, NaN₃, thyroxine, heparin, sulfonylcarbamides, Cu²⁺, Hg²⁺, Ag²⁺, Fe³⁺, Zn²⁺, K⁺, PO₄²⁻, NO₃⁻

Activators: Thioglycolic acid, b-mercaptoethylamine, EDTA, α, α'-dipyridyl

pH optimum: 8.0 (see figure)

Temperature dependence: See figure

pH stability: 5.5-6.5 (see figure)


Thermal stability: Up to +60°C (see figure)

Specificity: The oxidation of L-glutamate is stimulated by ADP and inhibited by GTP. In contrast, the oxidation of alanine, leucine, isoleucine, methionine,


Cat. No. **Pack Size**

10 190 462 103 custom fill

Will be supplied as "Glutamate Dehydrogenase from Beef Liver".
Unit of Measure is "MU".

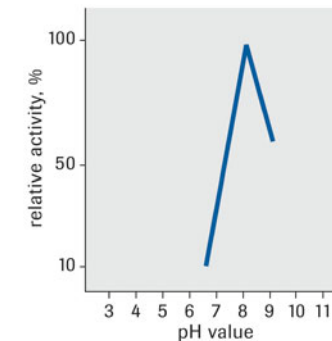

pH optimum

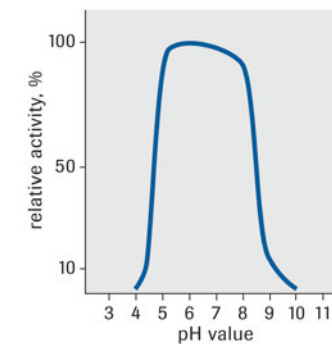

pH stability

Incubation:
25°C, 180 min
pH 3.0 – 5.0:
citrate buffer, 0.1 mol/l
pH 6.0 – 8.0:
phosphate buffer,
0.1 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.1
mol/l
120 U GIDH/ml

valine, norleucine, norvaline and 2-aminobutyrate is stimulated by GTP and inhibited by ADP.

Remarks:

- Glutamate dehydrogenase suspension or solution can be dialyzed against phosphate buffer, 10 mmol/l. Glutamate dehydrogenase molecules have the tendency to associate in some test formulations, modified Glutamate dehydrogenase minimizes this effect.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless to slightly opalescent solution in water (c=20 mg/ml)

pH value (c=20 mg/ml in water): 6.5-7.5

Activity (+25°C, α-oxoglutarat): ≥10 U/mg lyophilizate

Contaminants (expressed as percentage of Glutamate Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.005


Lactate dehydrogenase: ≤0.005

Malate dehydrogenase: ≤0.005


NH₄⁺: ≤0.1 µg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 18 months. Store dry.

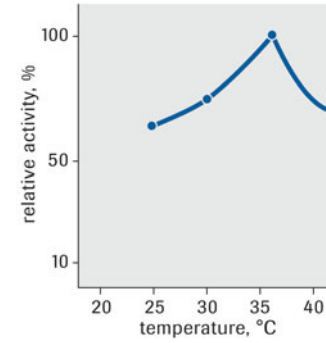

Temperature dependence

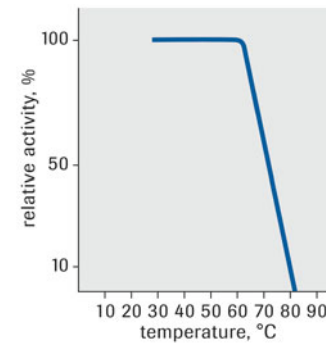

Thermal stability

Incubation:
10 min
(NH₄)₂ SO₄, 2.0 mol/l,
pH 7.0
2400 U GIDH/ml

For further processing only.

**Glutamate Dehydrogenase (NAD(P))
from beef liver, chemically modified, lyophilizate**

Dehydrogenase that catalyzes the conversion of glutamate to α-ketoglutarate.

Application

Use modified Glutamate Dehydrogenase in diagnostic tests for the determination of ammonia, urea, L-glutamate, glutamate pyruvate transaminase and leucine aminopeptidase.

Benefits

- Take advantage of the enhanced liquid stability.
- Rely on the proven diagnostic quality of this product.

EC 1.4.1.3

Properties

Nomenclature: L-glutamate:NAD(P)⁺ oxidoreductase (deaminating)

Molecular weight: ~2 200 kD for the associated enzyme with 8 subunits; 280 kD for one subunit.

Michaelis constants (Tris buffer, pH 8.0, +23°C):

L-glutamate: 1.8 x 10⁻³ mol/l

NADP: 4.7 x 10⁻⁵ mol/l

α-ketoglutarate: 7.0 x 10⁻⁴ mol/l

NH₄⁺: 3.2 x 10⁻³ mol/l

NADPH: 2.6 x 10⁻⁵ mol/l

Km values for NAD or NADH are difficult to obtain due to their inhibitory action.

Inhibitors: 4-chloromercuribenzoate, Na₂S, diethyldithiocarbamate, 1,10-phenanthroline,

Cat. No. **Pack Size**

11 434 993 103 custom fill

Will be supplied as "GIDH, Modified from Beef Liver". Unit of Measure is "MU".

For further processing only.

8-hydroxyquinoline, NaN_3 , thyroxine, heparin, sulfonylcarbamides, Cu^{2+} , Hg^{2+} , Ag^{2+} , Fe^{3+} , Zn^{2+} , K^+ , PO_4^{2-} , NO_3^-

Activators: Thioglycolic acid, b-mercaptoethylamine, EDTA, α , α' -dipyridyl

pH optimum: 8.0 (see figure)

Temperature dependence: See figure

pH stability: 5.5-6.5 (see figure)

Thermal stability: Up to $+60^\circ\text{C}$ (see figure)

Specificity: The oxidation of L-glutamate is stimulated by ADP and inhibited by GTP. In contrast, the oxidation of alanine, leucine, isoleucine, methionine, valine, norleucine, norvaline and 2-aminobutyrate is stimulated by GTP and inhibited by ADP.

Remarks:

- Glutamate dehydrogenase suspension or solution can be dialyzed against phosphate buffer, 10 mmol/l. Glutamate dehydrogenase molecules have the tendency to associate in some test formulations, modified Glutamate dehydrogenase minimizes this effect.

-The modified enzyme is especially suited for liquid stable applications with extended shelf life requirements.

Specification

Appearance: White lyophilizate, stabilized with RPLA 4 and ADP

Solubility: Clear, slightly opalescent solution in water ($c=40$ mg/ml)

Activity ($+25^\circ\text{C}$, α -oxoglutarat): ≥ 7 U/mg lyophilizate

Contaminants (expressed as percentage of Glutamate Dehydrogenase activity):

Alcohol dehydrogenase: ≤ 0.005

Lactate dehydrogenase: ≤ 0.005

Malate dehydrogenase: ≤ 0.005

NH_4^+ : ≤ 0.16 $\mu\text{mol/KU}$ Glutamate Dehydrogenase

K (flame photometric): ≤ 0.1 $\mu\text{mol/KU}$ Glutamate Dehydrogenase

Na (flame photometric): ≤ 2.0 $\mu\text{mol/KU}$ Glutamate Dehydrogenase

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At $+2$ to $+8^\circ\text{C}$ within specification range for 12 months. Store dry.

Glutamate Dehydrogenase (NAD(P)) from beef liver, solution

Dehydrogenase that catalyzes the conversion of glutamate to α -ketoglutarate.

Application

Use Glutamate Dehydrogenase in diagnostic tests for the determination of ammonia, urea, L-glutamate, glutamate pyruvate transaminase and leucine aminopeptidase.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.4.1.3

Properties

Nomenclature: L-glutamate:NAD(P)⁺ oxidoreductase (deaminating)

Molecular weight: ~ 2200 kD for the associated enzyme with 8 subunits; 280 kD for one subunit.

Michaelis constants (Tris buffer, pH 8.0, $+23^\circ\text{C}$):

L-glutamate: 1.8×10^{-3} mol/l

NADP: 4.7×10^{-5} mol/l

α -ketoglutarate: 7.0×10^{-4} mol/l

NH_4^+ : 3.2×10^{-3} mol/l

NADPH: 2.6×10^{-5} mol/l

Cat. No. **Pack Size**

10 120 847 103 custom fill

Will be supplied as "Glutamate Dehydrogenase, Beef Liver". Unit of Measure is "IU".

For further processing only.

K_m values for NAD or NADH are difficult to obtain due to their inhibitory action.

Inhibitors: 4-chloromercuribenzoate, Na₂S, diethyldithiocarbamate, 1,10-phenanthroline, 8-hydroxyquinoline, NaN₃, thyroxine, heparin, sulfonylcarbamides, Cu²⁺, Hg²⁺, Ag²⁺, Fe³⁺, Zn²⁺, K⁺, PO₄²⁻, NO₃⁻

Activators: Thioglycolic acid, b-mercaptoethylamine, EDTA, α, α'-dipyridyl

pH optimum: 8.0 (see figure)

Temperature dependence: See figure

pH stability: 5.5-6.5 (see figure)

Thermal stability: Up to +60°C (see figure)

Specificity: The oxidation of L-glutamate is stimulated by ADP and inhibited by GTP. In contrast, the oxidation of alanine, leucine, isoleucine, methionine, valine, norleucine, norvaline and 2-aminobutyrate is stimulated by GTP and inhibited by ADP.

Remarks: GIDH suspension or solution can be dialyzed against phosphate buffer, 10 mmol/l. GIDH molecules have the tendency to associate in some test formulations, modified GIDH minimizes this effect.

Specification

Appearance: Clear, colourless solution in glycerol

pH value: 7.0-7.8

Specific Activity: ≥ 120 U/mg

Protein (Biuret): 30±3 mg/ml

Contaminants (expressed as percentage of Glutamate Dehydrogenase activity):

Alcohol dehydrogenase: ≤0.01

Lactate dehydrogenase: ≤0.01

Malate dehydrogenase: ≤0.01

NH₄⁺: ≤0.16 µg/mg protein

Glycerol: 560-680 mg/ml (45-55% (v/v))

EDTA (complexometric): 12.2-13.4 mmol/l

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 18 months.

γ-Glutamyltransferase from hog kidney, lyophilizate

Application

Use γ-Glutamyltransferase for designing your calibrator or control reagent.

Benefits

Rely on the proven diagnostic quality of this product.

EC 2.3.2.2

Specification

Appearance: White/off white to buff lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

Activity (+37°C, with γ-GT kit): >23 U/mg lyophilizate

Contaminants (expressed as percentage of γ-Glutamyltransferase activity):

Leucine aminopeptidase: <0.10

Phosphatase, alkaline: <2

"NADH oxidase": ≤0.01

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.	Pack Size
10 445 363 103	custom fill

Will be supplied as "γ-Glutamyltransferase from Hog Kidney". Unit of Measure is "kU".

For further processing only.

Glycerol Kinase (GK)

from *Bacillus stearothermophilus*, lyophilizate

Enzyme that catalyzes the phosphorylation of glycerol to glycerol-3-phosphate.

Application

Use Glycerol Kinase for diagnostic tests for the determination of triglycerides together with Glycerol-3-phosphate Oxidase, Catalog No.11 582 003 103 and 11 654 730 103, and Lipoprotein Lipase, Catalog No. 11 145 991 103.

Benefits

- Rely on the proven diagnostic quality of this product.

EC 2.7.1.30

Properties

Nomenclature: ATP:glycerol 3-phosphotransferase**Molecular weight:** 230 000 D (Sephadex G 200), 4 x 58 000 D (SDS-gel electrophoresis)**Michaelis constants** (Glycine buffer, pH 9.8; +30°C):Glycerol: 4.4×10^{-5} mol/l**Inhibitors:** Unknown; Inhibitors of glycerokinase from *Candida mycoderma* do not inhibit the glycerokinase from *Bacillus stearothermophilus*.**pH optimum:** 10.0-10.5 (see figure)**Temperature dependence:** See figure**pH stability:** 5.0-11.0 (see figure)**Thermal stability:** Up to +60° C (see figure)

Specification


Appearance: White to slightly yellowish lyophilizate**Solubility:** Clear, colorless solution in water (c=10 mg/ml)**pH value** (c=10 mg/ml in water): 5.0-7.0**Activity** (+25°C, glycerol): 18-25 U/mg lyophilizate**Specific activity:** ≥ 80 U/mg protein**Protein** (Biuret): 0.18-0.26 mg/mg lyophilizate**Contaminants** (expressed as percentage of Glycerol Kinase activity):ATPase: ≤ 0.005 Hexokinase: ≤ 0.01 "NADH oxidase": ≤ 0.005 Glycerol: ≤ 4 $\mu\text{g}/100$ U (40 Mg/kU)**Stability:** At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.


Pack Size

11 499 530 103 custom fill


Will be supplied as "GK, B.stearot., Lyo., w. Lactose". Unit of Measure is "MU".

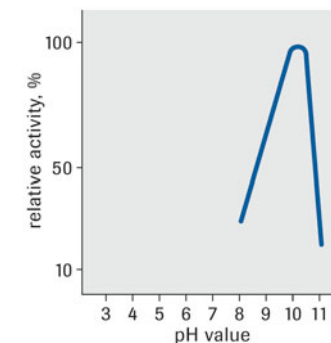


pH optimum

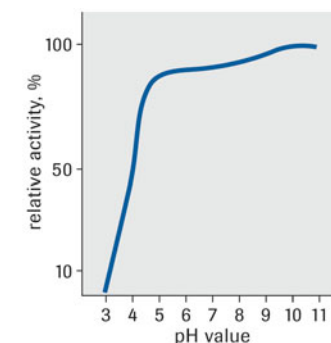

pH stability

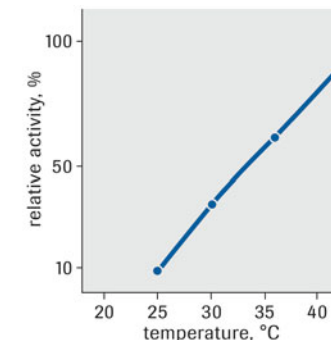

Temperature dependence

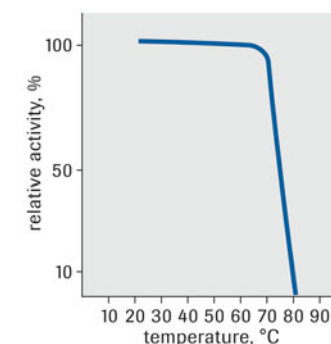

Thermal stability

Incubation:
25°C, 180 min
pH 3.0 – 5.0:
citrate buffer, 0.1 mol/l
pH 6.0 – 8.0:
phosphate buffer,
0.1 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.1
mol/l
85 U GK/ml

Incubation:
10 min
Tris buffer, 0.1 mol/l;
pH 9.0
510 U GK/ml

For further processing only.

Glycerol Kinase (GK) from *Bacillus stearothermophilus*, solution

Enzyme that catalyzes the phosphorylation of glycerol to glycerol-3-phosphate.

Application

Use Glycerol Kinase for diagnostic tests for the determination of triglycerides together with Glycerol-3-phosphate Oxidase, Catalog No. 11 582 003 103 and 11 654 730 103, and Lipoprotein Lipase, Catalog No. 11 145 991 103.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 2.7.1.30

Properties

Nomenclature: ATP:glycerol 3-phosphotransferase

Molecular weight: 230 000 D (Sephadex G 200), 4 x 58 000 D (SDS-gel electrophoresis)

Michaelis constants (Glycine buffer, pH 9.8; +30°C):

Glycerol: 4.4×10^{-5} mol/l

Inhibitors: Unknown; Inhibitors of glycerokinase from *Candida mycoderma* do not inhibit the glycerokinase from *Bacillus stearothermophilus*.

pH optimum: 10.0-10.5 (see figure)

Temperature dependence: See figure

pH stability: 5.0-11.0 (see figure)

Thermal stability: Up to +60° C (see figure)

Specification

Appearance: Clear, colorless to slightly yellowish solution in Tris buffer; 1 mol/l

Activity (+25°C, glycerol): ≥ 500 U/ml

Specific activity: ≥ 85 U/mg protein

Protein (Biuret): No limit

Contaminants (expressed as percentage of Glycerol Kinase activity):

Hexokinase: ≤ 0.01

"NADH oxidase": ≤ 0.005

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
10 691 666 103	custom fill

Will be supplied as "Glycerokinase from Bac.stearothermophil.". Unit of Measure is "kU".

For further processing only.

Glycerol Kinase (GK), concentrated from *Bacillus stearothermophilus*, solution

Enzyme that catalyzes the phosphorylation of glycerol to glycerol-3-phosphate.

Application

Use Glycerol Kinase for diagnostic tests for the determination of triglycerides together with Glycerol-3-phosphate Oxidase, Catalog No. 11 582 003 103 and 11 654 730 103, and Lipoprotein Lipase, Catalog No. 11 145 991 103.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 2.7.1.30

Properties

Nomenclature: ATP:glycerol 3-phosphotransferase

Molecular weight: 230 000 D (Sephadex G 200), 4 x 58 000 D (SDS-gel electrophoresis)

Cat. No.	Pack Size
10 539 937 103	custom fill

Will be supplied as "Glycerokinase, Bac. stearothermophilus". Unit of Measure is "MU".

For further processing only.

Michaelis constants (Glycine buffer, pH 9.8; +30°C):

Glycerol: 4.4×10^{-5} mol/l

Inhibitors: Unknown; Inhibitors of glycerokinase from *Candida mycoderma* do not inhibit the glycerokinase from *Bacillus stearothermophilus*.

pH optimum: 10.0-10.5 (see figure)

Temperature dependence: See figure

pH stability: 5.0-11.0 (see figure)

Thermal stability: Up to +60° C (see figure)

Specification

Appearance: Clear, colorless to slightly yellowish solution in Tris buffer; pH approx. 7.3; stabilized. Particles as result of recrystallized salts.

pH value (c=10 mg/ml, in water): 7.1-7.5

Activity (+25°C, glycerol): ≥ 2200 U/ml

Specific activity: ≥ 85 U/mg protein

Protein (Biuret): No limit

Contaminants (expressed as percentage of Glycerol Kinase activity):

Hexokinase: ≤ 0.01

"NADH oxidase": ≤ 0.005

Stability: At +2 to +8°C within specification range for 12 months.

Glycerol-3-phosphate Dehydrogenase from rabbit muscle, suspension

Dehydrogenase that catalyzes the interconversion of dihydroxyacetone phosphate to glycerol 3-phosphate.

Application

Use Glycerol-3-phosphate Dehydrogenase in diagnostic reagents for the determination of aldolase in combination with Triose-phosphate Isomerase, Catalog No. 10 153 338 103, and Fructose-1,6-diphosphate, Catalog No. 10 041 793 103.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.8

Specification

Appearance: White suspension in ammonium sulfate, 3.2 mol/l, pH approximately 6

pH value: 5.5-6.5

Activity: ≥ 2000 U/ml solution

Specific activity (+25°C, glyceraldehyde-3-phosphate): ≥ 170 U/mg protein

Protein (Biuret): 10 mg/ml

Ammonium sulphate: 3.2 ± 0.2 mol/l

Contaminants (expressed as percentage of Glycerol-3-phosphate Dehydrogenase activity):

Aldolase: ≤ 0.001

Glyceraldehyde-3-phosphate dehydrogenase: ≤ 0.001

Lactate dehydrogenase: ≤ 0.01

Triose-phosphate isomerase: ≤ 0.01

Stability: At +2 to +8°C within specification range for 18 months.

Cat. No.

Pack Size

10 151 351 103 custom fill

Will be supplied as "GDH from Rabbit Muscle". Unit of Measure is "MU".

For further processing only.

Glycerol-3-phosphate Oxidase from *E.coli* overproducer, lyophilizate

Recombinant oxidoreductase that catalyzes the interconversion of glycerol 3-phosphate to dihydroxyacetone phosphate.

Application

Use Glycerol-3-phosphate Oxidase in diagnostic tests for the determination of triglycerides together with Glycerol Kinase, Catalog No. 10 539 937 103 or 11 499 530 103 and Lipoprotein Lipase, Catalog No. 11 145 991 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.3.21

Properties

Nomenclature: Glycerol-3-phosphate : oxygen oxidoreductase

Molecular weight: 75 kD (SDS-PAGE); 74 kD (gel filtration, Sephadex G 150)

Isoelectric point: ~4.2 (230 000 D (Sephadex G 200), 4 x 58 000 D (SDS-gel electrophoresis))

Michaelis constants (L-glycerol phosphate):

K-phosphate buffer, 0.1 mol/l; pH 7.5: 1.36×10^{-2} mol/l (o-dianisidine assay)

Tris buffer, 0.1 mol/l; pH 7.6: 2.90×10^{-3} mol/l (o-dianisidine assay)

Tris buffer, 0.1 mol/l; pH 8.1: 1.40×10^{-3} mol/l (PAP assay)

Structure: Monomeric protein with FAD as cofactor

Inhibitors: Ag, Hg-salts and SDS

pH optimum: 8.0-8.5 (see figure)

Temperature dependence: See figure

pH stability: 6.5-8.5 (see figure)

Thermal stability: See figure

Specificity: Glycerol phosphate oxidase reacts highly specific with L- α -glycerol phosphate.

Specification

Appearance: Greenish yellow lyophilizate

Solubility: Clear yellow solution in Tris/HCl, 150 mmol/l, pH 7.6 (c=10 mg/ml)

Activity (+25°C, L- α -glycerol-3-phosphate): ≥ 50 U/mg lyophilizate

Activity (+37°C): ≥ 90 U/mg lyophilizate

Contaminants (expressed as percentage of Glycerol-3-phosphate Oxidase activity):

Cholesterol oxidase: ≤ 0.001

Lactate oxidase: ≤ 0.002

Uricase: ≤ 0.001

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Literature





1) J. Siedel, M. Town, W. Hölzel, Long-term stable, liquid, ready-to-use monoreagent for the enzymatic assay of serum or plasma triglycerides (GPO-PAP method), Poster presented at the XV Int. Congress of Clin. Chem., Melbourne/Australia, Nov. 1993

Cat. No. Pack Size

11 654 730 103 custom fill

Will be supplied as "L- α -Glycerol-phosphate Oxidase, rec., Lyo". Unit of Measure is "MU".

Additional formulation: Lyophilizate, Catalog No. 11 582 003 103

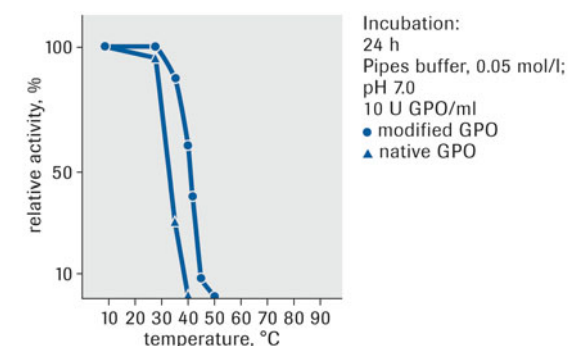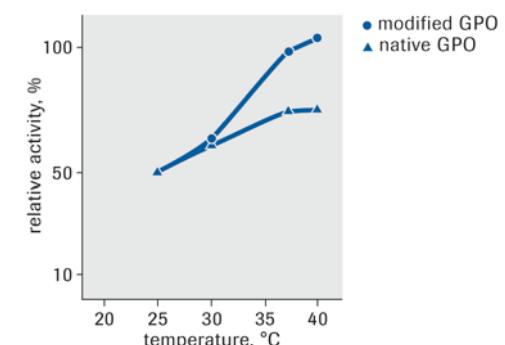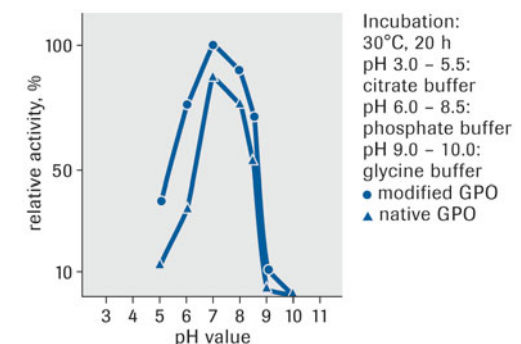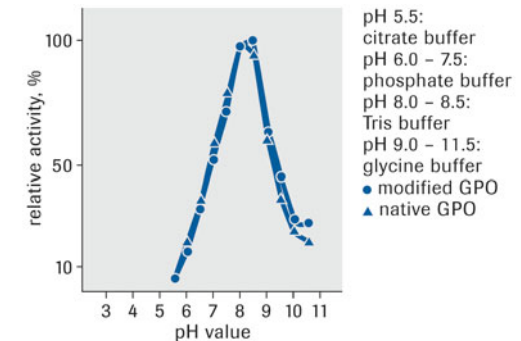

For further processing only.

Glycerol-3-phosphate Oxidase, chemically modified

from *E.coli* overproducer, lyophilizate

Recombinant oxidoreductase that catalyzes the interconversion of glycerol 3-phosphate to dihydroxyacetone phosphate.

Application

Use Glycerol-3-phosphate Oxidase in diagnostic tests for the determination of triglycerides together with Glycerol Kinase, Catalog No. 10 539 937 103 or 11 499 530 103 and Lipoprotein Lipase, Catalog No. 11 145 991 103.

Benefits

- Take advantage of the enhanced liquid stability of this enzyme.
- Rely on the proven diagnostic quality of this product.

EC 1.1.3.21

Properties

Nomenclature: Glycerol-3-phosphate : oxygen oxidoreductase

Molecular weight: 75 kD (SDS-PAGE); 74 kD (gel filtration, Sephadex G 150)

Isoelectric point: ~4.2 (230 000 D (Sephadex G 200), 4 x 58 000 D (SDS-gel electrophoresis))

Michaelis constants (L-glycerol phosphate):

K-phosphate buffer, 0.1 mol/l; pH 7.5: 1.36×10^{-2} mol/l (o-dianisidine assay)

Tris buffer, 0.1 mol/l; pH 7.6: 2.90×10^{-3} mol/l (o-dianisidine assay)

Tris buffer, 0.1 mol/l; pH 8.1: 1.40×10^{-3} mol/l (PAP assay)

Structure: Monomeric protein with FAD as cofactor

Inhibitors: Ag, Hg-salts and SDS

pH optimum: 8.0-8.5 (see figure)

Temperature dependence: See figure

pH stability: 6.5-8.5 (see figure)

Thermal stability: See figure

Specificity: Glycerol phosphate oxidase reacts highly specific with L- α -glycerol phosphate.

Specification

Appearance: Green yellow amorphous lyophilizate

Solubility: Clear yellow solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 6.8-7.8

Activity (+25°C, L- α -glycerol phosphate): ≥ 5 U/mg lyophilizate

Specific activity (+25°C): ≥ 40 U/mg protein

Activity (+37°C, L- α -glycerol phosphate): ≥ 10 U/mg lyophilizate

Protein (BCA): ≥ 0.1 mg/mg lyophilizate

Contaminants (expressed as percentage of Glycerol-3-phosphate Oxidase activity):

Cholesterol oxidase: ≤ 0.001

Lactate oxidase: ≤ 0.002

Uricase: ≤ 0.001

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Literature

1) J. Siedel, M. Town, W. Hölzel, Long-term stable, liquid, ready-to-use monoreagent for the enzymatic assay of serum or plasma triglycerides (GPO-PAP method), Poster presented at the XV Int. Congress of Clin. Chem., Melbourne/Australia, Nov. 1993

Cat. No.

Pack Size

11 582 003 103 custom fill

Will be supplied as "L-Glycerol-3-phosphate Oxidase rec. mod.". Unit of Measure is "MU".

Additional formulation: Lyophilizate, Catalog No. 11 654 730 For further processing only.

Hexokinase (HK) from yeast overproducer, lyophilizate

Recombinant enzyme that converts hexose to hexose-6-phosphate.

Application

Use Hexokinase in diagnostic tests for blood glucose using the hexokinase method.

Benefits

Rely on the proven diagnostic quality of this product.

EC 2.7.1.1

Properties

Nomenclature: ATP:D-hexose 6-phosphotransferase

Molecular weight: 57 kD (SDS-PAGE)

Isoelectric point: 4.5-5.0

Michaelis constants (D-glucose):

Phosphate buffer, 0.1 mol/l, pH 7.0; +25°C: 3.05×10^{-4} mol/l

Phosphate buffer, 0.1 mol/l, pH 7.4; +30°C: 1.90×10^{-4} mol/l

Tea buffer, 0.1 mol/l, pH 7.6; +25°C: 2.30×10^{-4} mol/l

Michaelis constants (ATP):

Tris buffer, 0.1 mol/l, pH 7.6; +28°C: 1.60×10^{-4} mol/l

Tea buffer, 0.1 mol/l, pH 7.6; +25°C: 1.90×10^{-4} mol/l

Inhibitors: EDTA, SH-blocking compounds, sorbose-1-phosphate, polyphosphates, 6-deoxy-6-fluoroglucose, 2-C-hydroxymethylglucose, lyxose.

Activators: Mg^{2+} , catecholamines

pH optimum: 7.0-10.0 (see figure)

Temperature dependence: See figure

pH stability: 5.0-9.0 (see figure)

Thermal stability: 100% for 20 minutes at +37°C, 50% for 20 minutes at +45°C, 75% for 5 hours at +37°C (see figures)

Specificity: Hexokinase phosphorylates D-glucose, D-fructose, D-mannose, D-glucosamin, 2-deoxyglucose. L-Arabinose, D-xylose, L-rhamnose, D-galactose, D-lactose, sucrose, maltose, trehalose, raffinose, N-acetyl glucosamine do not react. ATP can be partially replaced by other nucleotides.

Specification

Appearance: Yellowish lyophilizate

Solubility: Clear, colorless solution in water ($c=10$ mg/ml)

pH value ($c=10$ mg/ml in water): 6.5-7.5

Activity (+25°C, glucose): ≥ 70 U/mg lyophilizate

Activity (+30°C): ≥ 98 U/mg lyophilizate

Activity (+37°C): ≥ 115 U/mg

Protein (Biuret): 0.15 ± 0.05 mg/mg lyophilizate

Contaminants (expressed as percentage of Hexokinase activity):

Alcohol dehydrogenase: ≤ 0.001

ATPase: ≤ 0.05

Creatine kinase: ≤ 0.001

G6P-DH: ≤ 0.005

Glutamate dehydrogenase: ≤ 0.05

Glutathione reductase: ≤ 0.005

Myokinase: ≤ 0.001

"NADH oxidase": ≤ 0.001

"NADPH oxidase": ≤ 0.001

6-Phosphogluconate dehydrogenase: ≤ 0.001

Phosphoglucose isomerase: ≤ 0.002


Phosphoglucomutase: ≤ 0.02

Stability: At +2 to +8°C within specification range for 18 months. Store dry.


Cat. No. **Pack Size**

11 119 796 103 custom fill


Will be supplied as "Hexokinase (HK) from Yeast Overproducer".
Unit of Measure is "MU".

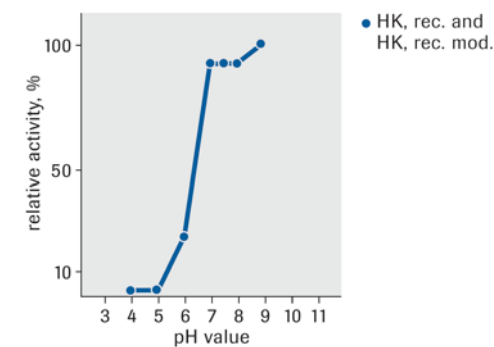


pH optimum

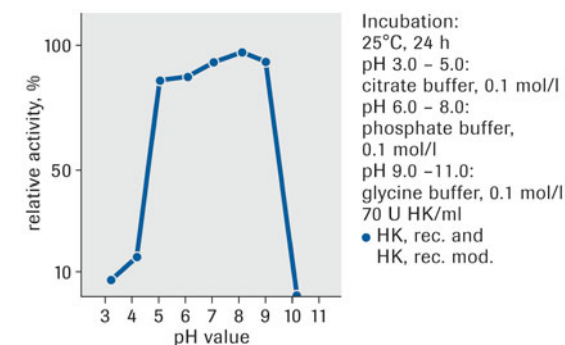


pH stability

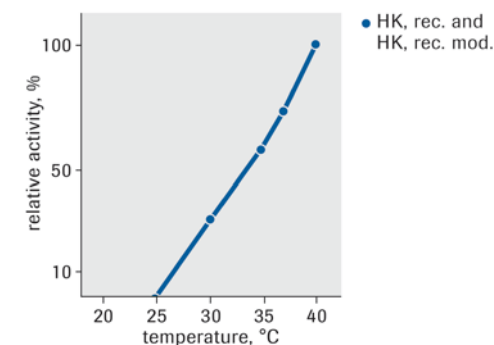

Temperature dependence

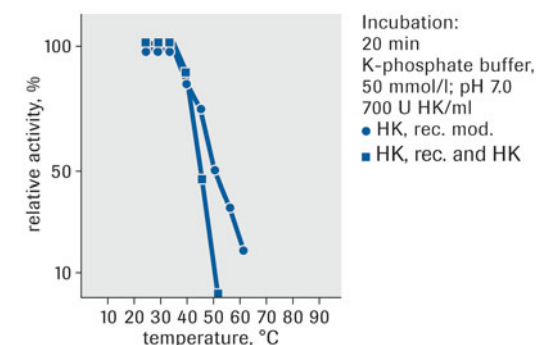

Thermal stability

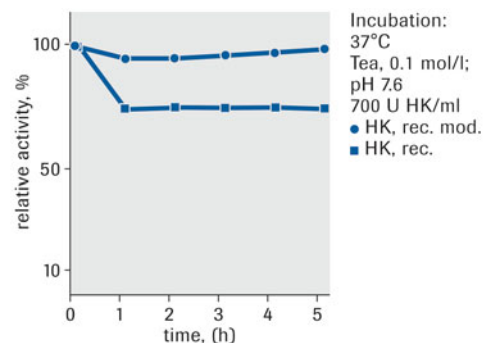

Thermal stability

For further processing only.

Hexokinase (HK), chemically modified from yeast overproducer, lyophilizate

Recombinant enzyme that converts hexose to hexose-6-phosphate.

Application

Use Hexokinase in diagnostic tests for blood glucose using the hexokinase method.

Benefits

- Take advantage of the improved stability in liquid reagents.
- Rely on the proven diagnostic quality of this product.

EC 2.7.1.1

Properties

Nomenclature: ATP:D-hexose 6-phosphotransferase

Molecular weight: 57 kD (SDS-PAGE)

Isoelectric point: 4.5-5.0

Michaelis constants (D-glucose):

Phosphate buffer, 0.1 mol/l, pH 7.0; +25°C: 3.05×10^{-4} mol/l

Phosphate buffer, 0.1 mol/l, pH 7.4; +30°C: 1.90×10^{-4} mol/l

Tea buffer, 0.1 mol/l, pH 7.6; +25°C: 2.30×10^{-4} mol/l

Michaelis constants (ATP):

Tris buffer, 0.1 mol/l, pH 7.6; +28°C: 1.60×10^{-4} mol/l

Tea buffer, 0.1 mol/l, pH 7.6; +25°C: 1.90×10^{-4} mol/l

Inhibitors: EDTA, SH-blocking compounds, sorbose-1-phosphate, polyphosphates, 6-deoxy-6-fluoroglucose, 2-C-hydroxymethylglucose, lyxose.

Activators: Mg^{2+} , catecholamines

pH optimum: 7.0-10.0 (see figure)

Temperature dependence: See figure

pH stability: 5.0-9.0 (see figure)

Thermal stability: see figure p. 108.

Specificity: Hexokinase phosphorylates D-glucose, D-fructose, D-mannose, D-glucosamin, 2-deoxyglucose. L-Arabinose, D-xylose, L-rhamnose, D-galactose, D-lactose, sucrose, maltose, trehalose, raffinose, N-acetyl glucosamine do not react. ATP can be partially replaced by other nucleotides.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=40 mg/ml)

pH value (c=40 mg/ml in water): 6.5-7.5

Activity (+25°C, glucose): ≥ 40 U/mg lyophilizate

108 Contaminants (expressed as percentage of hexokinase activity):

Cat. No. **Pack Size**

11 370 600 103 custom fill

Will be supplied as "Hexokinase (HK) from Rec.Yeast, Modif.". Unit of Measure is "MU".

For further processing only.

Alcohol dehydrogenase: ≤ 0.001
 ATPase: ≤ 0.05
 Creatine kinase: ≤ 0.001
 G6P-DH: ≤ 0.005
 Glutamate dehydrogenase: ≤ 0.05
 Glutathione reductase: ≤ 0.005
 Myokinase: ≤ 0.001
 "NADH oxidase": ≤ 0.001
 6-Phosphogluconate dehydrogenase: ≤ 0.001
 Phosphoglucose isomerase: ≤ 0.002
 Phosphoglucomutase: ≤ 0.02
 Glucose: $\leq 0.3 \mu\text{g}/\text{mg}$ lyophilizate
Stability: At +2 to +8°C within specification range for 18 months. Store dry.
Remark: This enzyme is especially suited for liquid stable applications with extended shelf life requirements.

Hexokinase (HK) from yeast overproducer, solution

Recombinant enzyme that converts hexose to hexose-6-phosphate.

Application

Use Hexokinase in diagnostic tests for blood glucose using the hexokinase method.

Benefits

- Apply this ready-to-use recombinant enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 2.7.1.1

Properties

Nomenclature: ATP:D-hexose 6-phosphotransferase

Molecular weight: 57 kD (SDS-PAGE)

Isoelectric point: 4.5-5.0

Michaelis constants (D-glucose):

Phosphate buffer, 0.1 mol/l, pH 7.0; +25°C: 3.05×10^{-4} mol/l

Phosphate buffer, 0.1 mol/l, pH 7.4; +30°C: 1.90×10^{-4} mol/l

Tea buffer, 0.1 mol/l, pH 7.6; +25°C: 2.30×10^{-4} mol/l

Michaelis constants (ATP):

Tris buffer, 0.1 mol/l, pH 7.6; +28°C: 1.60×10^{-4} mol/l

Tea buffer, 0.1 mol/l, pH 7.6; +25°C: 1.90×10^{-4} mol/l

Inhibitors: EDTA, SH-blocking compounds, sorbose-1-phosphate, polyphosphates, 6-deoxy-6-fluoroglucose, 2-C-hydroxymethylglucose, lyxose.

Activators: Mg^{2+} , catecholamines

pH optimum: 7.0-10.0 (see figure)

Temperature dependence: See figure

pH stability: 5.0-9.0 (see figure)

Thermal stability: 100% for 20 minutes at +37°C, 50% for 20 minutes at +45°C, 75% for 5 hours at +37°C (see figures)

Specificity: Hexokinase phosphorylates D-glucose, D-fructose, D-mannose, D-glucosamin, 2-deoxyglucose. L-Arabinose, D-xylose, L-rhamnose, D-galactose, D-lactose, sucrose, maltose, trehalose, raffinose, N-acetyl glucosamine do not react. ATP can be partially replaced by other nucleotides.

Specification

Appearance: Clear, yellowish solution, in 50% glycerol (v/v)

pH value: 6.0-7.0

Cat. No.	Pack Size
11 149 130 103	custom fill

Will be supplied as "Hexokinase (HK) from Recombinant Yeast".
 Unit of Measure is "MU".
 For further processing only.

Activity (+25°C, glucose): $\geq 1,200$ U/ml
Activity (+30°C): $\geq 1,680$ U/ml
Protein (Biuret): ≥ 75 mg/ml
Contaminants (expressed as percentage of Hexokinase activity):
 Alcohol dehydrogenase: ≤ 0.001
 ATPase: ≤ 0.05
 Creatine kinase: ≤ 0.001
 G6P-DH: ≤ 0.005
 Glutamate dehydrogenase: ≤ 0.01
 Glutathione reductase: ≤ 0.005
 Myokinase: ≤ 0.001
 6-Phosphogluconate dehydrogenase: ≤ 0.001
 Phosphoglucose isomerase: ≤ 0.002
 Phosphoglucomutase: ≤ 0.02
 Glucose: ≤ 0.125 $\mu\text{g}/\text{mg}$ lyophilizate
 Glycerol (enzymatically): 45-55%
Stability: At +2 to +8°C within specification range for 18 months. Store dry.

Hexokinase (HK) from yeast, lyophilizate

Enzyme that converts hexose to hexose-6-phosphate.

Application

Use Hexokinase in diagnostic tests for blood glucose using the hexokinase method.

Benefits

Rely on the proven diagnostic quality of this product.

EC 2.7.1.1

Specification

Appearance: Yellowish lyophilizate
Solubility: Clear, colorless solution in water (c=10 mg/ml)
pH value (c=10 mg/ml): 6.5-7.5
Activity (+25°C; Glucose): ≥ 40 U/mg lyophilizate
Protein (Biuret): 0.4 ± 0.1 mg/mg lyophilizate
Contaminants (expressed as percentage of Hexokinase activity):
 ATPase: ≤ 0.05
 Alcohol dehydrogenase: ≤ 0.001
 Creatine kinase: ≤ 0.001
 G6P-DH: ≤ 0.005
 Glutamate dehydrogenase: ≤ 0.05
 Glutathione reductase: ≤ 0.005
 Myokinase: ≤ 0.001
 Phosphoglucomutase: ≤ 0.02
 6-Phosphogluconate dehydrogenase: ≤ 0.001
 Phosphoglucose isomerase: ≤ 0.002
Stability: At +2 to +8°C within specification range for 18 months. Store dry.

Cat. No.	Pack Size
10 152 676 103	custom fill

Will be supplied as "Hexokinase (HK) from Yeast". Unit of Measure is "MU".

For further processing only.

D-Lactate Dehydrogenase (D-LDH) from microorganism, lyophilizate

Recombinant dehydrogenase that catalyzes the interconversion of D(-)-lactate to pyruvate.

Cat. No.	Pack Size
12 235 650 103	custom fill

Will be supplied as "LDH". Unit of Measure is "MU".



For further processing only.

Application

Use D-Lactate Dehydrogenase in a variety of diagnostic tests, *e.g.*, in the determination of alanine aminotransferases, lactate or pyruvate. Used for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aspartate aminotransferases, glutamate dehydrogenase).

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.28

Specification

Appearance: White to slightly yellow powder or lyophilizate

Solubility: Soluble in water

Activity (+25°C, lyophilizate): ≥340 U/mg

Stability: At -15 to -25°C within specification range for 18 months. Store dry in tightly sealed containers.

D-Lactate Dehydrogenase (D-LDH), Grade I

from *Lactobacillus delbrückii*, lyophilizate

Dehydrogenase that catalyzes the interconversion of D(-)-lactate to pyruvate.

Application

Use D-Lactate Dehydrogenase (D-LDH), Grade I, in a variety of diagnostic tests, *e.g.*, in the determination of alanine aminotransferases, lactate or pyruvate. Used for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aspartate aminotransferases, glutamate dehydrogenase).

Benefits

- Rely on the proven diagnostic quality of this product.
- Benefit from the extended shelf life of this enzyme.

EC 1.1.1.28

Properties

Nomenclature: D-lactate:NAD⁺ oxidoreductase

Michaelis constants (Tris maleate buffer, pH 8.0, +25°C):

D-lactate: 0.7×10^{-1} mol/l (NAD, 2 mmol/l)

Pyruvate: 1.2×10^{-3} mol/l (NADH, 0.1 mmol/l)

NADH: 7.1×10^{-5} mol/l (pyruvate, 20 mmol/l)

pH optimum: 7.0 (see figure)

Temperature dependence: See figure

pH stability: 4.0-10.0 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Lactate dehydrogenase is specific for D(-)-lactate, L(+)-lactate does not react.

Remark: Lactate dehydrogenase, Grade I is especially suited for liquid stable applications with extended shelf life requirements.

Specification

Appearance: White to yellowish lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 6.0-7.0

Activity (+25°C, pyruvate): ≥180 U/mg lyophilizate

Specific activity: ≥450 U/mg protein


Protein (Biuret) : No limit, 0.3-0.8 mg/mg lyophilizate

Contaminants (expressed as percentage of D-Lactate Dehydrogenase


Cat. No. **Pack Size**

11 291 416 103 custom fill

Will be supplied as "D(-)-LDH, Special Quality". Unit of Measure is "MU".

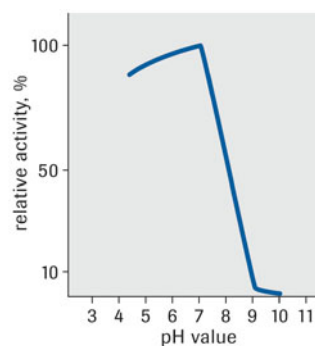

pH optimum

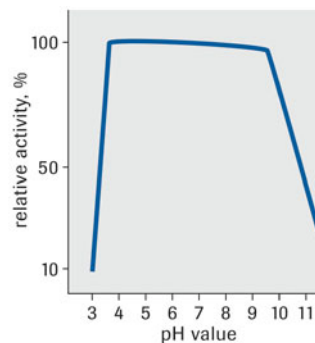

pH stability

Incubation:
25°C, 60 min
pH 3.0 – 5.0:
citrate buffer, 0.2 mol/l
pH 6.0 – 8.0:
phosphate buffer,
0.2 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.2 mol/l
2400 U D-LDH/ml

activity):

Alcohol dehydrogenase: ≤ 0.01

Malate dehydrogenase: ≤ 0.1

"NADH oxidase": ≤ 0.0005


Succinate dehydrogenase: ≤ 0.01

NH₄⁺: ≤ 0.01 $\mu\text{mol/KU}$


Na (flame photometric): ≤ 0.5 $\mu\text{mol/KU}$

K (flame photometric): ≤ 0.007 $\mu\text{mol/KU}$

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

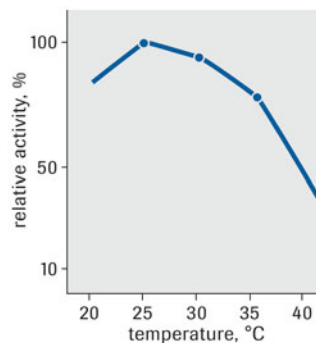

Temperature dependence

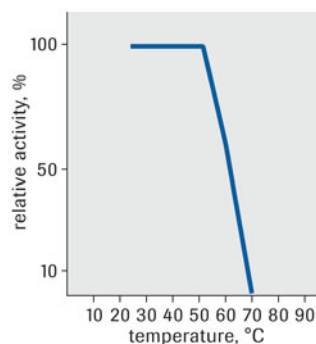

Thermal stability

Incubation:
10 min
water; pH 5.2
2400 U D-LDH/ml

For further processing only.

D-Lactate Dehydrogenase (D-LDH), Grade II

from *Lactobacillus delbrückii*, lyophilizate

Dehydrogenase that catalyzes the interconversion of D(-)-lactate to pyruvate.

Application

Use D-Lactate Dehydrogenase (D-LDH), Grade II, in a variety of diagnostic tests, *e.g.*, in the determination of alanine aminotransferases, lactate or pyruvate. Used for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aspartate aminotransferases, glutamate dehydrogenase).

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.28

Properties

Nomenclature: D-lactate:NAD⁺ oxidoreductase

Michaelis constants (Tris maleate buffer, pH 8.0, +25°C):

D-lactate: 0.7×10^{-1} mol/l (NAD, 2 mmol/l)

Pyruvate: 1.2×10^{-3} mol/l (NADH, 0.1 mmol/l)

NADH: 7.1×10^{-5} mol/l (pyruvate, 20 mmol/l)

pH optimum: 7.0 (see figure)

Temperature dependence: See figure

pH stability: 4.0-10.0 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Lactate dehydrogenase is specific for D(-)-lactate, L(+)-lactate does not react.

Cat. No.

Pack Size

10 679 666 103 custom fill

Will be supplied as "D(-)-Lactate Dehydrogenase (D-LDH)". Unit of Measure is "MU".

For further processing only.

Specification

Appearance: White to yellowish lyophilizate
Solubility: Clear, colorless solution in water (c=10 mg/ml)
pH value (c=10 mg/ml in water): 5.7-6.7
Activity (+25°C, pyruvate): ≥150 U/mg lyophilizate
Specific activity: ≥300 U/mg protein
Protein (Biuret): No limit, approximately 0.4-0.7 mg/mg lyophilizate
Contaminants (expressed as percentage of D-Lactate Dehydrogenase activity):
 Alcohol dehydrogenase: ≤0.01
 Glucose dehydrogenase: ≤0.01
 Malate dehydrogenase: ≤0.1
 Succinate dehydrogenase: ≤0.01
Stability: At +2 to +8°C within specification range for 12 months. Store dry.

L-Lactate Dehydrogenase (L-LDH) from pig muscle, for use of AST/GOT-Determination according to IFCC recommendations, lyophilizate

Dehydrogenase that catalyzes the interconversion of specific for L(+)-lactate to pyruvate.

Application

Use L-Lactate Dehydrogenase in a variety of diagnostic tests for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aminotransferases, glutamate dehydrogenase).

Benefits

- Rely on the proven diagnostic quality of this product.
- Tested according to the recommendations of the International Federation of Clinical Chemistry (IFCC).

EC 1.1.1.27

Properties


Nomenclature: L-lactate:NAD⁺ oxidoreductase
Molecular weight: 140 kD
Isoelectric point: 4.6
Michaelis constants (Phosphate buffer, pH 7.5; +25°C):
 Pyruvate: 1.5×10^{-4} mol/l (NADH: 0.18 mmol/l)
 L-lactate: 3.3×10^{-3} mol/l (NAD: 0.5 mmol/l)
 NADH: 1.1×10^{-5} mol/l (Pyruvate: 0.6 mmol/l)
 NAD: 6.7×10^{-5} mol/l (L-lactate: 34 mmol/l)
Inhibitors: Oxamate, pyruvate (excess), oxalate, Ag⁺, Hg²⁺, Cu²⁺
pH optimum: 3.0-7.0 (see figure)
Temperature dependence: See figure
pH stability: 5.5-8.5 (see figure)
Thermal stability: Up to +40°C (see figure)
Specificity: Lactate dehydrogenase is specific for L(+)-lactate, D(-)-lactate does not react. Glyoxylate is also a lactate dehydrogenase substrate. Apart from pyruvate some 2-oxoacids are reduced. NAD analogs (*e.g.*, APAD) react at similar rates.

Specification


Appearance: White lyophilizate
Solubility: Clear colorless solution in water (c=10 mg/ml)
pH value (c=10 mg/ml in water): 6.0-7.0
Activity (+30°C, pyruvate; according to IFCC recommendations):
 for Aspartate aminotransferase (AST/GOT) determination: ≥50 U/mg

Cat. No.	Pack Size
10 254 754 103	custom fill

Will be supplied as "LDH IFCC-quality from Hog Muscle". Unit of Measure is "MU".
 Additional formulation: Suspension in glycerol solution, Catalog No. 10 417 718

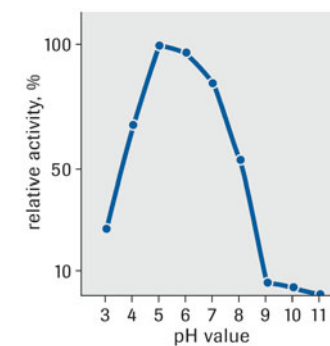

pH optimum

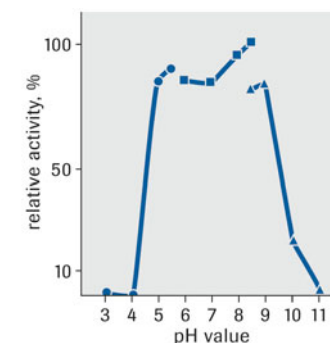

pH stability

Incubation:
 25°C, 6 h resp. 24 h
 ● pH 3.0 – 5.5:
 citrate buffer, 0.1 mol/l
 ■ pH 5.5 – 8.5:
 phosphate buffer,
 0.1 mol/l
 ▲ pH 8.5 – 11.0:
 glycine buffer, 0.1 mol/l
 10 U LDH/ml

lyophilizate

for Alanine aminotransferase (ALT/GPT) determination: ≥ 50 U/mg lyophilizate

Activity (mean value of both determinations): ≥ 50 U/mg lyophilizate

Contaminants (expressed as percentage of Lactate Dehydrogenase activity, assayed according to the IFCC recommendations):

Aspartate aminotransferase (AST/GOT): ≤ 0.001

unspecificity of Lactate dehydrogenase: ≤ 0.005


Alanine aminotransferase (ALT/GPT): ≤ 0.001

unspecificity of Lactate dehydrogenase: ≤ 0.005


SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

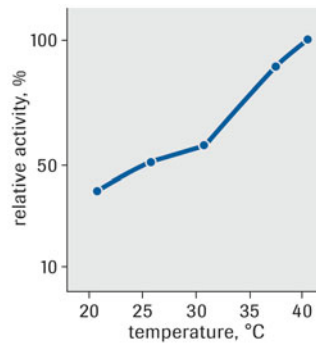

Temperature dependence

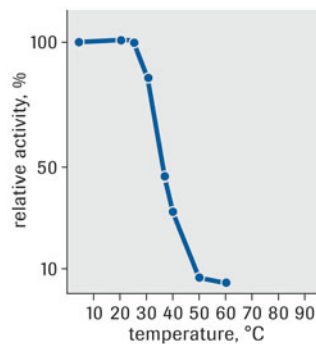

Thermal stability

For further processing only.

Incubation:
10 min
Tris buffer, 0.1 mol/l;
pH approx. 7.0
834 U LDH/ml

L-Lactate Dehydrogenase (L-LDH) from pig muscle, for use of AST/GOT-Determination according to IFCC recommendations, solution

Dehydrogenase that catalyzes the interconversion of specific for L(+)-lactate to pyruvate.

Application

Use L-Lactate Dehydrogenase in a variety of diagnostic tests for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aminotransferases, glutamate dehydrogenase).

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.
- Tested according to the recommendations of the International Federation of Clinical Chemistry (IFCC).

EC 1.1.1.27

Properties

Nomenclature: L-lactate:NAD⁺ oxidoreductase

Molecular weight: 140 kD

Isoelectric point: 4.6

Michaelis constants (Phosphate buffer, pH 7.5; +25°C):

Pyruvate: 1.5×10^{-4} mol/l (NADH: 0.18 mmol/l)

L-lactate: 3.3×10^{-3} mol/l (NAD: 0.5 mmol/l)

NADH: 1.1×10^{-5} mol/l (Pyruvate: 0.6 mmol/l)

NAD: 6.7×10^{-5} mol/l (L-lactate: 34 mmol/l)

Inhibitors: Oxamate, pyruvate (excess), oxalate, Ag⁺, Hg²⁺, Cu²⁺

Cat. No. **Pack Size**

10 417 718 103 custom fill

Will be supplied as "LDH, IFCC-quality from Hog Muscle". Unit of Measure is "MU".

For further processing only.

pH optimum: 3.0-7.0 (see figure)
Temperature dependence: See figure
pH stability: 5.5-8.5 (see figure)
Thermal stability: Up to +40°C (see figure)
Specificity: Lactate dehydrogenase is specific for L(+)-lactate, D(-)-lactate does not react. Glyoxylate is also a lactate dehydrogenase substrate. Apart from pyruvate some 2-oxoacids are reduced. NAD analogs (*e.g.*, APAD) react at similar rates.

Specification

Appearance: Clear, colorless solution in glycerol, 50% (v/v), pH approximately 7
pH value: 6.5-7.5
Activity (+30°C, pyruvate, according to the IFCC recommendations): ≥9,600 U/ml solution
Specific activity: ≥480 U/mg protein
Protein (Biuret): ≥20 mg/ml solution
Glycerol: 45-55 % (v/v)
Contaminants (expressed as percentage of Lactate Dehydrogenase specific activity):
 Glutamate dehydrogenase: ≤0.003
 Aspartate aminotransferase (AST/GOT): ≤0.005
 Alanine aminotransferase (ALT/GPT): ≤0.005
Reagent blank for determination of aspartate aminotransferase (AST/GOT): ≤0.9 mA/min
SVD free: Corresponds to specification
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At +2 to +8°C within specification range for 12 months.

L-Lactate Dehydrogenase (L-LDH) from pig muscle, suspension

Dehydrogenase that catalyzes the interconversion of specific for L(+)-lactate to pyruvate.

Application

Use L-Lactate Dehydrogenase in a variety of diagnostic tests for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aminotransferases, glutamate dehydrogenase).

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.27

Specification

Appearance: White suspension in ammonium sulfate, 3.2 mol/l; Tris, 10 mmol/l, pH approximately 6.5
pH value: 6.0-7.0
Specific activity (+25°C, pyruvate): ≥550 U/mg protein
Protein (Biuret): ≥10 mg/ml
Ammonium sulphate: 3.2±0.2 mol/l
Contaminants (expressed as percentage of Lactate Dehydrogenase activity):
 Aldolase: ≤0.001
 Glutamate dehydrogenase: ≤0.01
 Aspartate aminotransferase (AST/GOT): ≤0.005
 Alanine aminotransferase (ALT/GPT): ≤0.005
 Malate dehydrogenase: ≤0.01
 Myokinase: ≤0.01
 Pyruvate kinase: ≤0.001
SVD free: Corresponds to specification
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
10 021 415 103	custom fill

Will be supplied as "Lactate Dehydrogenase (LDH), Hog Muscle".
 Unit of Measure is "MU".
 For further processing only.

L-Lactate Dehydrogenase (L-LDH), chemically modified from pig heart, lyophilizate

Dehydrogenase that catalyzes the interconversion of L(+)-lactate to pyruvate.

Application

Use L-Lactate Dehydrogenase (L-LDH), chemically modified, in a variety of diagnostic tests for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aminotransferases, glutamate dehydrogenase).

Benefits

- Take advantage of the enhanced liquid stability of this enzyme.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.27

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value: 7.1-8.1

Activity (+25°C, pyruvate): ≥25 U/mg lyophilizate

Specific activity: ≥150 U/mg protein

Protein (BCA): 0.15-0.25 mg/mg lyophilizate

Contaminants (expressed as percentage of Lactate Dehydrogenase activity):

Aspartate aminotransferase (AST/GOT): ≤0.005

Unspecificity of Lactate dehydrogenase: ≤0.05

Alanine aminotransferase (ALT/GPT): ≤0.01

Unspecificity of Lactate dehydrogenase: ≤0.05

"NADH-Oxidase": ≤0.001

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability in CAPSO (pH 9.4, at +60°C for 1 hour): ≥50%

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.

Pack Size

11 866 117 103 custom fill

Will be supplied as "Lactate Dehydrogenase (LDH)". Unit of Measure is "MU".

For further processing only.

L-Lactate Dehydrogenase (L-LDH) from pig heart, suspension

Dehydrogenase that catalyzes the interconversion of L(+)-lactate to pyruvate.

Application

Use L-Lactate Dehydrogenase (L-LDH), in a variety of diagnostic tests for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aminotransferases, glutamate dehydrogenase).

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality.

EC 1.1.1.27

Specification

Appearance: White suspension in ammonium sulfate solution, 3.2 mol/l, pH approximately 6

pH value: 5.5-6.5

Specific activity (+25°C, pyruvate): ≥300 U / mg protein

Protein (Biuret): ≥10 mg/ml

Ammonium sulphate: 3.2±0.2 mol/l

Cat. No.

Pack Size

10 153 729 103 custom fill

Will be supplied as "Lactate Dehydrogenase (LDH), Pig Heart". Unit of Measure is "MU".

For further processing only.

Contaminants (expressed as percentage of Lactate Dehydrogenase activity):

ATPase: ≤ 0.001

Glutamate dehydrogenase: ≤ 0.02

Aspartate aminotransferase (AST/GOT): ≤ 0.005

Malate dehydrogenase: ≤ 0.01

Myokinase: ≤ 0.01

"NADH-Oxidase": ≤ 0.001

Pyruvate kinase: ≤ 0.01

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months.

L-Lactate Dehydrogenase (L-LDH) from rabbit muscle, suspension

Dehydrogenase that catalyzes the interconversion of specific for L(+)-lactate to pyruvate.

Application

Use L-Lactate Dehydrogenase in a variety of diagnostic tests for the removal of pyruvate in determinations working with NADH (*i.e.*, triglycerides, lipase, aldolase, aminotransferases, glutamate dehydrogenase).

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.27

Specification

Appearance: White suspension in ammonium sulfate, 3.2 mol/l; Tris, 10 mmol/l; pH approximately 6.5

pH value: 6.0-7.0

Specific activity (+25°C, pyruvate): ≥ 550 U/mg protein

Protein (Biuret): ≥ 10 mg/ml

Ammonium sulphate: 3.2 ± 0.2 mol/l

Contaminants (expressed as percentage of Lactate Dehydrogenase activity):

Aldolase: ≤ 0.001

Aspartate aminotransferase (AST/GOT): ≤ 0.01

Alanine aminotransferase (ALT/GPT): ≤ 0.01

Malate dehydrogenase: ≤ 0.01

Myokinase: ≤ 0.01

Pyruvate kinase: ≤ 0.00

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
10 003 557 103	custom fill

Will be supplied as "Lactate Dehydrogenase, Rabbit Muscle". Unit of Measure is "MU".

For further processing only.

Lactate 2-Monooxygenase (Lactate oxidase), Grade I

from *Aerococcus viridans*, expressed in *E. coli*, lyophilizate

Recombinant oxidoreductase that catalyzes the conversion of lactate to pyruvate.

Application

Use Lactate 2-Monooxygenase, Grade I in diagnostic tests for the determination of L-lactate.

Benefits

- Benefit from the high activity in this Grade I enzyme.
- Rely on the proven diagnostic quality of this product.

EC 1.13.12.4

Properties

Nomenclature: L-lactate:oxigen oxidoreductase

Michaelis constant: L-lactate: 5×10^{-4} mol/l

V_{maximum} : L-lactate: 0.2 mol/l

pH optimum: 6.5-7.5 (see figure)

Temperature dependence: See figure

pH stability: 6.0-9.0 (see figure)

Thermal stability: Up to +65°C (see figure)

Specification

Appearance: Yellow lyophilizate

Activity: (+25°C, L-lactate): ≥ 40 U/mg lyophilizate

Specific activity: ≥ 55 U/mg protein

Contaminants (expressed as percentage of Lactate 2-Monooxygenase activity):

Catalase: ≤ 0.2

Glucose oxidase: ≤ 0.001

Pyruvate oxidase: ≤ 0.001

Uricase: ≤ 0.001





Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No.

Pack Size

04 822 277 103 custom fill

Will be supplied as "Lactat-OD, SQ, rec., lyo". Unit of Measure is "MU".

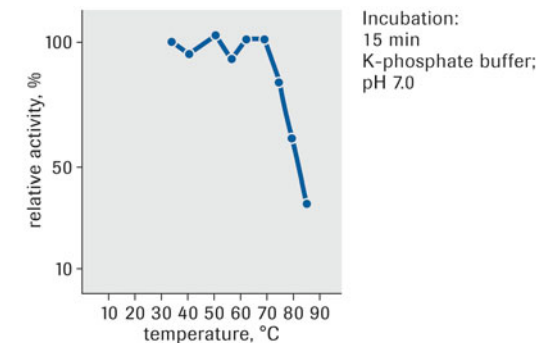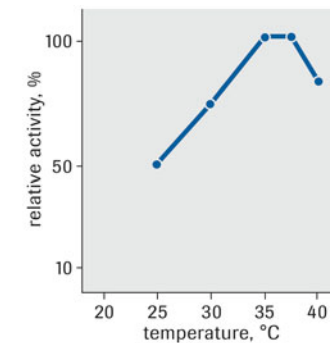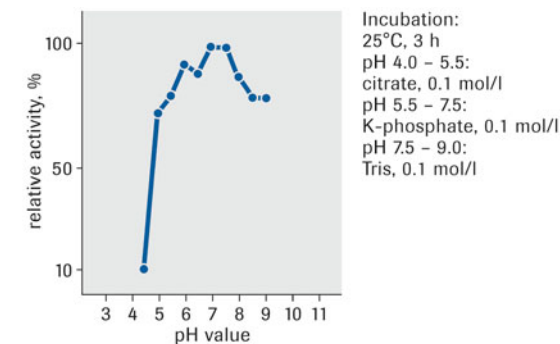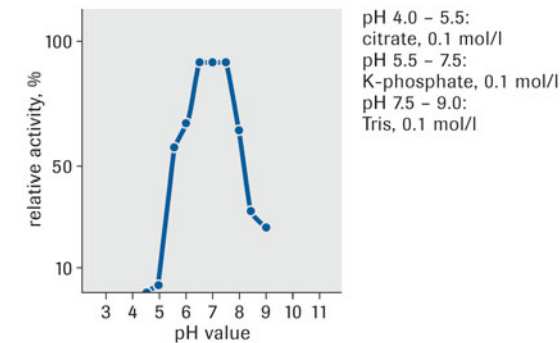

For further processing only.

Lactate 2-Monooxygenase (Lactate oxidase), Grade II from *Aerococcus viridans*, expressed in *E. coli*, lyophilizate

Recombinant oxidoreductase that catalyzes the conversion of lactate to pyruvate.

Application

Use Lactate 2-Monooxygenase, Grade II in diagnostic tests for the determination of L-lactate.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.13.12.4

Properties

Nomenclature: L-lactate:oxigen oxidoreductase

Michaelis constant: L-lactate: 5×10^{-4} mol/l

V_{maximum}: L-lactate: 0.2 mol/l

pH optimum: 6.5-7.5 (see figure)

Temperature dependence: See figure

pH stability: 6.0-9.0 (see figure)

Thermal stability: Up to +65°C (see figure)

Specification

Appearance: Yellow lyophilizate

Activity (+25°C, L-lactate): ≥ 20 U/mg lyophilizate

Specific activity: ≥ 55 U/mg protein

Protein (BCA): 0.3-0.7 mg/mg lyophilizate

Contaminants (expressed as percentage of Lactate 2-Monooxygenase activity):

Catalase: ≤ 0.2

Glucose oxidase: ≤ 0.001

Pyruvate oxidase: ≤ 0.001

Uricase: ≤ 0.001

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

11 798 197 103 custom fill

Will be supplied as "Lactat-OD, rec., Lyo.". Unit of Measure is "MU".

For further processing only.

Lactate 2-Monooxygenase (Lactate oxidase) from *Pediococcus* species, lyophilizate

Oxidoreductase that catalyzes the conversion of lactate to pyruvate.

Application

Use Lactate 2-Monooxygenase in diagnostic tests for the determination of lactate.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.13.12.4

Specification

Appearance: Yellow lyophilizate

Activity (25°C, L-lactate): ≥ 20 U/mg lyophilizate

Specific activity: ≥ 55 U/mg protein

Protein (Lowry): 0.2-0.4 mg/mg lyophilizate

Stability: At -15 to -25°C within specification range for 18 months. Store dry.

Cat. No. **Pack Size**

10 980 927 103 custom fill

Will be supplied as "Lactate Qxidase from *Pediococcus* species". Unit of Measure is "kU".



For further processing only.

Lipase from porcine pancreas, lyophilizate

Lipolytic enzyme that hydrolyzes triglycerides.

Application

Use Lipase in diagnostic tests for the determination of triglycerides in combination with Colipase, Catalog No. 10 204 307 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.1.1.3

Specification

Appearance: White lyophilizate, stabilized

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value: 7.0-8.0

Activity: ≥300 U/mg lyophilizate

Contaminants (expressed as percentage of Lipase activity):

Proteases: ≥85

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 410 551 103 custom fill

Will be supplied as "Lipase from Porcine Pancreas". Unit of Measure is "kU".

For further processing only.

Lipoprotein Lipase from *Pseudomonas* species, lyophilizate

Enzyme that hydrolyzes triglycerides into three free fatty acids and glycerol.

Application

Use Lipoprotein Lipase in diagnostic tests for the determination of triglycerides together with Glycerol Kinase, Catalog No. 10 539 937 103 or 11 499 530 103 and Glycerol-3-phosphate Dehydrogenase, Catalog No. 11 654 730 103 or 11 582 003 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.1.1.34

Properties

Nomenclature: Triacylglycerol-protein acylhydrolase

Molecular weight: 47 kD

Effectors: Hg²⁺, Ag⁺, Cr²⁺, Sn²⁺, Cu²⁺ and ionic detergents inhibit. Mg²⁺, sodium cholate and BSA stabilize the enzyme. 4-Chloromercuribenzoate (2 mmol/l), monoiodoacetate (2 mmol/l), NaF (20 mmol/l), NaN₃ (20 mmol/l), EDTA (5 mmol/l) and 2-phenanthroline (2 mmol/l) do not affect the enzyme activity while SDS (0.1% (w/v)) is inactivating.

pH optimum: 7.5 (see figure)

pH stability: 6.0-10.0 (see figure)


Thermal stability: Up to +50°C (see figure)

Specificity: Lipoprotein Lipase has both lipolytic and sterol ester hydrolytic activities. It hydrolyzes triacylglycerols in chylomicrons, lipoproteins and diacylglycerols. With human plasma as substrate triglycerides are hydrolyzed more rapidly than cholesterol esters. The effects of pH and ionic strength on the enzymatic activity are somewhat different between the hydrolysis of triglyceride

Cat. No. **Pack Size**

10 734 284 103 custom fill

Will be supplied as "Lipoprotein Lipase from Pseudomon.spec.". Unit of Measure is "MU".

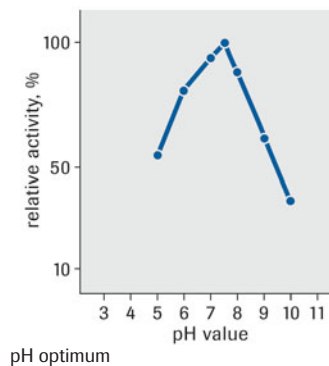

and of cholesterol ester depending on the different states of these substrates in the plasma or the transfer of the reaction products at the interface of substrates.

Lipolytic activity (Substrate, Number of C-atoms to number of double bonds, Relative rate):

- olive oil: 94%
- triolein (18:1): 100%
- tripalmitin (16:0): 2%
- trimyristin (14:0): 7%
- trilaurin (12:0): 4%
- tricaprin (10:0): 17%
- tricaprylin (8:0): 64%
- tricaproin (6:0): 2%
- tributylin (4:0): 2%
- tripropionin (3:0): 2%
- triacetin (2:0): 1%

Specification

Appearance: Brownish lyophilizate


Solubility: Clear, brown solution in water (c=50 mg/ml)

Activity (+25°C, cholesterol oleate): ≥100 U/mg lyophilizate

Contaminants (expressed as percentage of Lipoprotein Lipase activity):


- ATPase: ≤0.005
- Catalase: ≤1.0
- Glycerokinase: ≤0.001
- Glucose oxidase: ≤0.001
- Hexokinase: ≤0.005
- "NADH oxidase": ≤0.001
- Uricase: ≤0.005

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

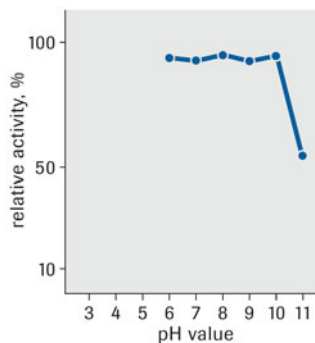

Incubation:
25°C, 60 min
pH 6.0 – 8.0:
phosphate buffer,
0.7 mol/l
pH 9.0 – 11.0:
glycine buffer, 2.1 mol/l
50 U LPL/ml

pH stability

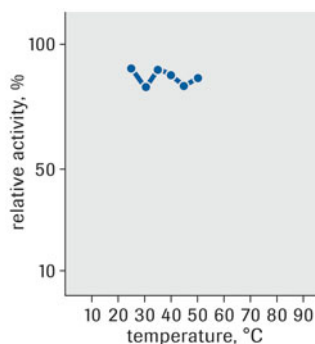

Incubation:
30 min
Tris buffer, 0.1 mol/l;
pH 7.7
50 U LPL/ml

Thermal stability

For further processing only.

Lipoprotein Lipase, chemically modified from *Pseudomonas* species, lyophilizate

Enzyme that hydrolyzes triglycerides into three free fatty acids and glycerol.

Application

Use Lipoprotein lipase in diagnostic tests for the determination of triglycerides together with Glycerol Kinase, Catalog Nos. 10 539 937 103 or 11 499 530 103 and Glycerol-3-phosphate Dehydrogenase, Catalog Nos. 11 654 730 103 or 11 582 003 103.

Benefits

- Take advantage of the enhanced liquid stability of this enzyme.
- Rely on the proven diagnostic quality of this product.

EC 3.1.1.34

Properties

Nomenclature: Triacylglycerol-protein acylhydrolase


Molecular weight: 47 kD

Effectors: Hg²⁺, Ag⁺, Cr²⁺, Sn²⁺, Cu²⁺ and ionic detergents inhibit. Mg²⁺, sodium cholate and BSA stabilize the enzyme. 4-Chloromercuribenzoate (2 mmol/l), monoiodoacetate (2 mmol/l), NaF (20 mmol/l), NaN₃ (20 mmol/l), EDTA (5 mmol/l) and 2-phenanthroline (2 mmol/l) do not affect the enzyme activity while SDS (0.1% (w/v)) is inactivating.

Cat. No. **Pack Size**

11 145 991 103 custom fill

Will be supplied as "Lipoprotein Lipase Modified". Unit of Measure is "MU".

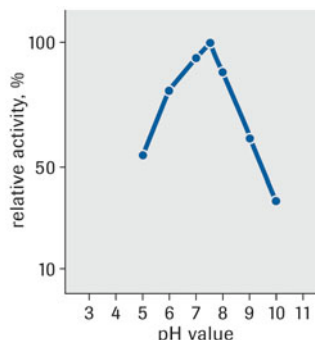

pH optimum

pH optimum: 7.5 (see figure)

pH stability: 6.0-10.0 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Lipoprotein Lipase has both lipolytic and sterol ester hydrolytic activities. It hydrolyzes triacylglycerols in chylomicrons, lipoproteins and diacylglycerols. With human plasma as substrate triglycerides are hydrolyzed more rapidly than cholesterol esters. The effects of pH and ionic strength on the enzymatic activity are somewhat different between the hydrolysis of triglyceride and of cholesterol ester depending on the different states of these substrates in the plasma or the transfer of the reaction products at the interface of substrates.

Lipolytic activity (Substrate, Number of C-atoms to number of double bonds, Relative rate):

olive oil: 94%

triolein (18:1): 100%

tripalmitin (16:0): 2%

trimyristin (14:0): 7%

trilaurin (12:0): 4%

tricaprin (10:0): 17%

tricaprylin (8:0): 64%

tricaproin (6:0): 2%

tributylin (4:0): 2%

tripropionin (3:0): 2%

triacetin (2:0): 1%

Remark: Chemically modified Lipoprotein Lipase (LPL) is more hydrophilic than native LPL. Carryover effect is therefore reduced.

Specification

Appearance: Brownish lyophilizate

Solubility: Clear, brown solution in water (c=50 mg/ml)

Activity (+25°C, cholesterol oleate): ≥10 U/mg lyophilizate

Contaminants (expressed as percentage of Lipoprotein Lipase activity):

ATPase: ≤0.005

Catalase: ≤1.0


Glycerokinase: ≤0.001

Hexokinase: ≤0.005


"NADH oxidase": ≤0.001

Uricase: ≤0.005

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

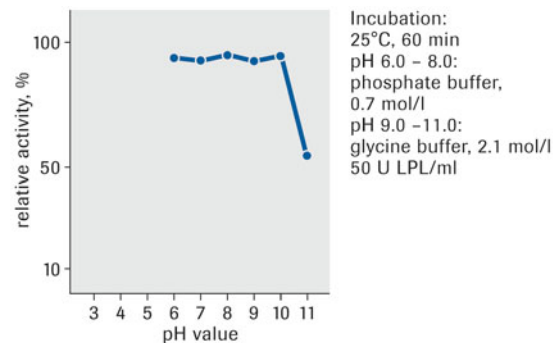

pH stability

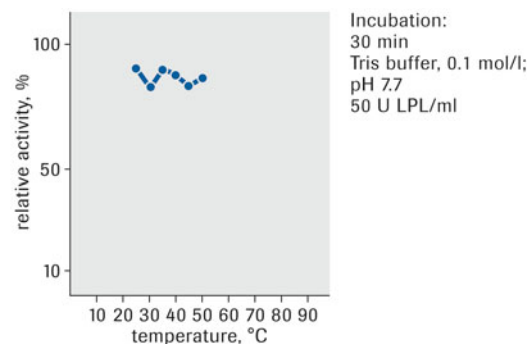

Thermal stability

For further processing only.

Lysozyme from hen egg white, crystalline powder

Glucosidic bond hydrolyzing enzyme

Application

Use Lysozyme for bacteriolysis, preparation of protoplasts and sample preparation prior to isolation of nucleic acids. It can also be used in pharmaceutical applications and food industry.

Benefits

Rely on the proven diagnostic quality of this product.

Product Description

Lysozyme from chicken egg. During purification and processing of this enzyme, steps included pH treatment at pH 3.5 or less for at minimum 30 minutes.

EC 3.2.1.17

Cat. No. **Pack Size**

10 153 516 103 custom fill

Will be supplied as "Lysozyme (Muramidase) from Hen Egg White". Unit of Measure is "g".

For further processing only.

Specification

Appearance: White, crystalline powder
Activity (+25°C, with *Micrococcus luteus*): ≥ 12,200 U/mg substance
Activity (+25°C, with *M. luteus*, previous Roche-substrate; calculated): ≥ 50,000 U/mg substance
Proteases: ≤ 0.5 U/mg substance
Stability: At +2 to +8°C within specification range for 36 months. Store dry.

Malate Dehydrogenase from pig heart, lyophilizate

Dehydrogenase that catalyzes the interconversion of malate to oxaloacetate.

Application

Use Malate Dehydrogenase in diagnostic tests for the determination of aspartate aminotransferase or in applications for citric and acetic acid testing.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.37

Specification

Appearance: White lyophilizate, stabilized
Solubility: Clear, colorless solution in water (c=10 mg/ml)
pH value: 7.0-8.0
Activity (+25°C, oxaloacetate): ≥ 70 U/mg lyophilizate
Contaminants (expressed as percentage of Malate Dehydrogenase activity):
 Fumarase: ≤ 0.01
 Aspartate aminotransferase (AST/GOT): ≤ 0.002
 Lactate dehydrogenase: ≤ 0.01
SVD free: Corresponds to specification
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 200 387 103 custom fill

Will be supplied as "Malate Dehydrogenase, Pig Heart (Mitochond.)". Unit of Measure is "MU".
 For further processing only.

Malate Dehydrogenase, chemically modified from pig heart, lyophilizate

Dehydrogenase that catalyzes the interconversion of malate to oxaloacetate.

Application

Use Malate Dehydrogenase in diagnostic tests for the determination of aspartate aminotransferase or in applications for citric and acetic acid testing.

Benefits

- Take advantage of the enhanced liquid stability
- Rely on the proven diagnostic quality of this product.

EC 1.1.1.37

Specification

Appearance: White lyophilizate
pH value: 7.5-8.5
Activity (+25°C, oxaloacetate): ≥ 20 U/mg lyophilizate
Specific activity: ≥ 400 U/mg protein
Protein (BCA): ≥ 0.02 mg/mg lyophilizate
Contaminants (expressed as percentage of Malate Dehydrogenase activity):

Cat. No. **Pack Size**

11 866 109 103 custom fill

Will be supplied as "MDH, Lyo., mod.". Unit of Measure is "MU".



For further processing only.

Aspartate aminotransferase (AST/GOT): ≤ 0.01

"NADH-Oxidase": ≤ 0.005

SVD free: Corresponds to specification

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 18 months.

Malate Dehydrogenase, IFCC Quality from pig heart, lyophilizate

Dehydrogenase that catalyzes the interconversion of malate to oxaloacetate.

Application

Use Malate Dehydrogenase in diagnostic tests for the determination of aspartate aminotransferase or in applications for citric and acetic acid testing.

Benefits

- Rely on the proven diagnostic quality of this product.
- Tested according to the recommendations of the International Federation of Clinical Chemistry (IFCC).

EC 1.1.1.37

Properties

Nomenclature: L-malate:NAD⁺ oxidoreductase

Molecular weight: 70 kD

Isoelectric point: 6.1-6.4

Michaelis constants (Phosphate buffer, 95 mmol/l, pH 8.3, $+25^{\circ}\text{C}$):

L-malate: 4.0×10^{-4} mol/l

L-tartrate: 9.0×10^{-3} mol/l

meso-tartrate: 1.2×10^{-3} mol/l

oxaloacetate: 3.3×10^{-5} mol/l

Inhibitors: Iodinated compounds such as iodine cyanide, thyroxine and molecular iodine, phenols, 1,10-phenanthroline, 8-hydroxyquinoline, sulfide, nicotinic acidamide, adenine, AMP, ATP; oxaloacetate (excess).

Activators: Phosphate, arsenate, Zn^{2+}

pH optimum: 7.5 (see figure)

Temperature dependence: See figure

pH stability: 7.0-9.0 (see figure)

Thermal stability: Up to $+40^{\circ}\text{C}$ (see figure)

Specificity: L-configuration of malate and tartrate. NAD can be replaced by its analogs, but not by NADP.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water ($c=10$ mg/ml)

Activity ($+30^{\circ}\text{C}$, oxaloacetate; according to the IFCC recommendations): ≥ 70 U/mg lyophilizate

Contaminants (expressed as percentage of Malat Dehydrogenase activity; assayed according to the IFCC recommendations):

Aspartate aminotransferase (AST/GOT): ≤ 0.001

Alanine aminotransferase (ALT/GPT): ≤ 0.001

Glutamate dehydrogenase: ≤ 0.005

SVD free: Corresponds to specification


pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At $+2$ to $+8^{\circ}\text{C}$ within specification range for 12 months. Store dry.


Cat. No. **Pack Size**

10 267 155 103 custom fill


Will be supplied as "MDH IFCC-quality, Pig Heart (Mitochond.)".
Unit of Measure is "MU".

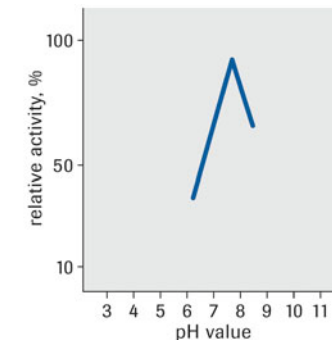

pH optimum

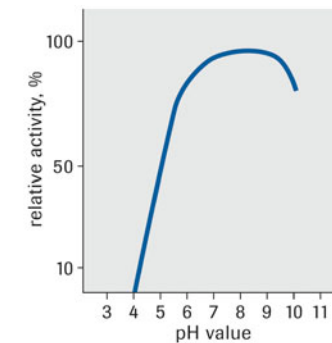


pH stability

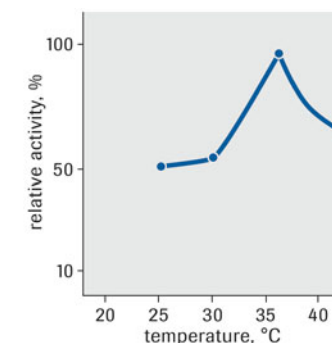

Temperature dependence

Incubation:
 25°C , 180 min
pH 3.0 – 5.0:
citrate buffer, 0.1 mol/l
pH 6.0 – 8.0:
phosphate buffer,
0.1 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.1 mol/l
1200 U MDH/ml

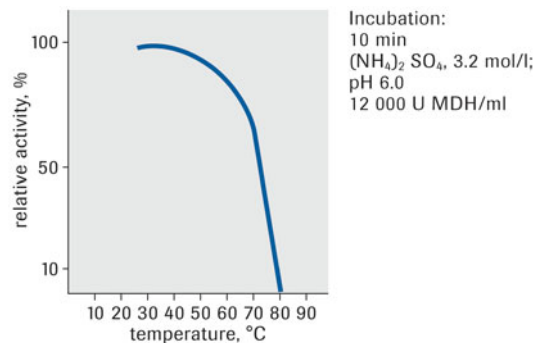

Thermal stability

For further processing only.

Malate Dehydrogenase, IFCC Quality from pig heart, solution

Dehydrogenase that catalyzes the interconversion of malate to oxaloacetate.

Application

Use Malate Dehydrogenase in diagnostic tests for the determination of aspartate aminotransferase or in applications for citric and acetic acid testing.

Benefits

- Apply this ready-to-use enzyme directly in your diagnostic test.
- Rely on the proven diagnostic quality of this product.
- Tested according to the recommendations of the International Federation of Clinical Chemistry (IFCC)

EC 1.1.1.37

Properties

Nomenclature: L-malate:NAD⁺ oxidoreductase

Molecular weight: 70 kD

Isoelectric point: 6.1-6.4

Michaelis constants (Phosphate buffer, 95 mmol/l, pH 8.3, +25°C):

L-malate: 4.0 x 10⁻⁴ mol/l

L-tartrate: 9.0 x 10⁻³ mol/l

meso-tartrate: 1.2 x 10⁻³ mol/l

oxaloacetate: 3.3 x 10⁻⁵ mol/l

Inhibitors: Iodinated compounds such as iodine cyanide, thyroxine and molecular iodine, phenols, 1,10-phenanthroline, 8-hydroxyquinoline, sulfide, nicotinic acidamide, adenine, AMP, ATP; oxaloacetate (excess).

Activators: Phosphate, arsenate, Zn²⁺

pH optimum: 7.5 (see figure)

Temperature dependence: See figure

pH stability: 7.0-9.0 (see figure)

Thermal stability: Up to +40°C (see figure)

Specificity: L-configuration of malate and tartrate. NAD can be replaced by its analogs, but not by NADP.

Specification

Appearance: Clear, colorless solution in glycerol (50% (v/v))

pH value: 6.5-7.5

Specific activity (+30°C, oxaloacetate): ≥900 U/mg protein

Protein (Biuret): ≥10 mg/ml

Contaminants (expressed as percentage of Malate Dehydrogenase activity):

Aspartate aminotransferase (AST/GOT): ≤0.005

Cat. No. **Pack Size**

10 417 726 103 custom fill

Will be supplied as "MDH, IFCC-quality, Pig Heart (Mitochon)".
Unit of Measure is "MU".

For further processing only.

Alanine aminotransferase (ALT/GPT): ≤ 0.005

Glutamate dehydrogenase: ≤ 0.003

Reagent blank for determination of aspartate aminotransferase (AST/GOT): ≥ 0.009 ($\Delta A_{334}/10$ minutes)

SVD free: Corresponds to specification

pH 5.5 treatment (for at minimum 30 minutes): Corresponds to specification

Stability: At +2 to +8°C within specification range for 12 months.

N-Methylhydantoinase (ATP-hydrolyzing) from *Arthrobacter* species, expressed in *E. coli*, lyophilizate

Hydrolase for creatinine determination that uses ATP to catalyze the conversion of N-methylhydantoin to N-carbamoylsarcosine.

Application

Use N-Methylhydantoinase (ATP-hydrolyzing) in diagnostic tests for the determination of creatinine in combination with Creatinine Deaminase Catalog No. 11 330 764 103, N-Carbamoylsarcosine Amidase, Catalog No. 11 248 847 103 and Sarcosine Oxidase, Catalog No. 11 378 856 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.2.14

Specification

Appearance: White lyophilizate

Solubility. Clear, colorless solution in water (c=10 mg/ml, +25 °C)

pH value (c=100 mg/ml in water): 7.8-8.8

Activity (+25°C, N-methylhydantoin): 0.6-1.0 U/ mg lyophilizate

Protein (Biuret): 20-43 mg/100 mg lyophilizate

Contaminants (expressed as percentage of N-Methylhydantoinase activity):

Creatinase: ≤ 0.013

Creatininase: ≤ 0.01

Catalase: ≤ 100

Uricase: ≤ 0.01

Stability: At -15 to -25°C within specification range for 12 months. Store dry. Protect from light.

Cat. No. **Pack Size**

11 288 555 103 custom fill

Will be supplied as "N-Methylhydantoin Hydrolase". Unit of Measure is "kU".



For further processing only.

NAD(P)H Dehydrogenase (quinone) (Diaphorase) from pig heart, suspension

Dehydrogenase that catalyzes the oxidation of dihydrolipoyl groups and has diaphorase activity.

Application

Use the diaphorase activity of NAD(P)H Dehydrogenase (quinone) for the determination of NAD(P)H and many dehydrogenases when coupled with various dyes which act as hydrogen acceptors from NAD(P)H, e.g. tetrazolium salts.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.6.5.2

Cat. No. **Pack Size**

10 153 427 103 custom fill

Will be supplied as "Diaphorase, Grade I from Pig Heart". Unit of Measure is "g".

For further processing only.

Specification

Appearance: Yellow suspension in ammonium sulfate, 3.2 mol/l
pH value: 5.5-6.5
Specific activity (+25°C, lipoate): ≥25 U/mg protein
Protein (Biuret): 10±1 mg/ml
Ammonium sulphate: 3.2±0.2 mol/l
SVD free: Corresponds to specification
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At +2 to +8°C within specification range for 12 months.

Nitrate Reductase
from *Aspergillus* species, lyophilizate

Oxidoreductase that catalyzes the reduction of nitrate to nitrite.

Application

Use Nitrate Reductase in diagnostic tests for the determination of nitrate.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.6.6.2

Specification

Appearance: Yellow lyophilizate
Activity (+25°C, nitrate): ≥0.4 U/mg lyophilizate
Specific activity: ≥10 U/mg protein
Protein (Biuret): No limit
Contaminants (expressed as percentage of Nitrate Reductase activity):
 Alcohol dehydrogenase (NADPH dependent): ≤0.8
 "NADPH oxidase": ≤0.5
 Nitrite reductase: ≤0.15
Stability: At -15 to -25°C within specification range for 12 months. Store dry.
 Protect from light.

Cat. No.	Pack Size
10 918 202 103	custom fill

Will be supplied as "Nitrate Reductase (*Aspergillus* species)". Unit of Measure is "kU".



For further processing only.

Oxalate Oxidase from barley seedings, lyophilizate

Oxidoreductase that catalyzes the interconversion of oxalate to carbon dioxide and hydrogen peroxide.

Application

Use Oxalate Oxidase in diagnostic tests for the determination of oxalate in combination with Formate Dehydrogenase, Catalog No. 10 204 226 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.2.3.4

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=40 mg/ml)

Activity (+37°C, oxalate): ≥0.25 U/mg lyophilizate

Specific activity: ≥5 U/mg protein

Protein (Biuret): No limit

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

10 570 524 103 custom fill

Will be supplied as "Oxalate Oxidase from Barley Seedlings". Unit of Measure is "U".



For further processing only.

Peroxidase (POD), Grade II from horse radish, lyophilizate

Application

Use Peroxidase (POD), Grade II, for the oxidation of reduced dyes in the indicator reaction of many diagnostic tests, e.g., for the determination of blood glucose, triglycerides or lactate. It may also be used as a marker enzyme for enzyme immunoassays (EIA).

Benefits

Rely on the proven diagnostic quality of this enzyme.

EC 1.11.1.7

Properties

Nomenclature: Donor:hydrogen-peroxide oxidoreductase

Molecular weight: 40 kD

Structure: Glycoprotein with one mole of protoheme IX

Isoelectric point: 7.2

Rate constants:

a) Hydrogen acceptors:

H₂O₂ 9 x 10⁸ [l x mol⁻¹ x s⁻¹]

methyl peroxide 1.5 x 10⁶ [l x mol⁻¹ x s⁻¹]

ethyl peroxide 3.6 x 10⁶ [l x mol⁻¹ x s⁻¹]

b) Hydrogen donors: Many

Inhibitors: Cyanide, sulfide, fluoride, azide, hydroxylamine, hydroxyl ions

Activators: Peroxidation of o-dianisidine is accelerated by ammonia, pyridine, imidazole at pH values >7.0

pH optimum: 6.0-6.5 (see figure)

Temperature dependence: See figure

pH stability: 4.0-10.0 (see figure)


Thermal stability: Up to +40°C (see figure)

Specificity: Peroxidase is specific for the hydrogen acceptor; only H₂O₂, methyl- and ethylperoxides are active. In contrast the enzyme is not specific for the hydrogen donor. A large number of phenols, aminophenols, diamines, indophenols, leucocyte dyes, ascorbate and several amino acids react.


Cat. No. **Pack Size**

11 378 783 103 custom fill

Will be supplied as "Peroxidase (POD), Grade II, Horse-radish". Unit of Measure is "MU".

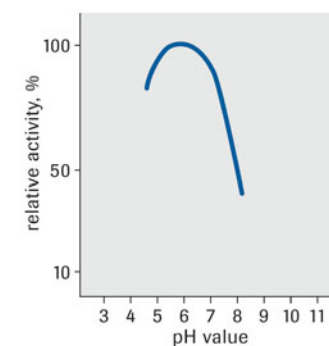

pH optimum

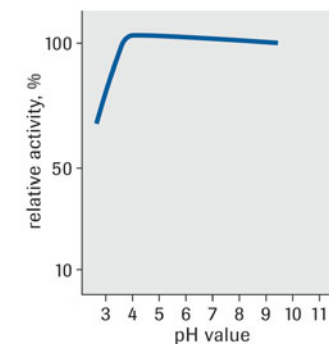


pH stability


Incubation:
25°C, 180 min
pH 3.0 – 5.0:
citrate buffer, 0.1 mol/l
pH 6.0 – 8.0:
phosphate buffer,
0.1 mol/l
pH 9.0 – 11.0:
glycine buffer, 0.1 mol/l
200 U POD/ml

Specification

Appearance: Red-brown lyophilizate
Solubility: Clear, red-brown solution in water (c=10 mg/ml)
pH value (c=10 mg/ml in water): 6.6-7.6
Activity (+25°C, guaiacol, H₂O₂): ≥200 U/mg lyophilizate
Purity number (A₄₀₃/A₂₇₅): 2.0-3.5
A₅₀₀ (100 U/ml): ≤0.120
Contaminants (expressed as percentage of Peroxidase activity):
 ATPase: ≤0.001
 Catalase: ≤0.7
 Contaminating oxidases: ≤0.00005
 Phosphatase, acidic: ≤0.001
Glucose: ≤0.25 µg/mg lyophilizate
Stability: At +2 to +8°C within specification range for 24 months. Store dry. Keep tightly sealed.

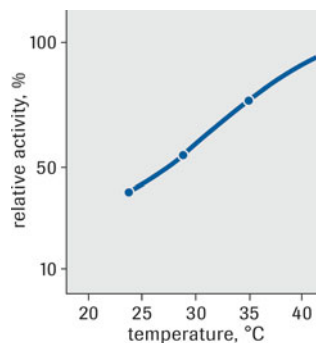

Temperature dependence

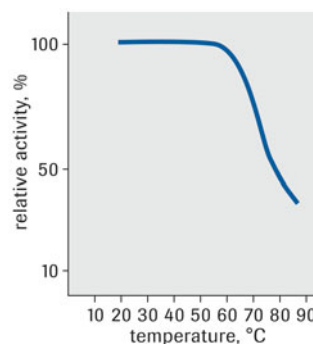

Thermal stability

Incubation:
 10 min
 phosphate buffer,
 0.1 mol/l; pH 7.0
 2 000 U POD/ml

For further processing only.

Phosphogluconate Dehydrogenase (decarboxylating) from yeast, lyophilizate

Dehydrogenase that catalyzes the formation of ribulose 5-phosphate from 6-phosphogluconate.

Application

Use Phosphogluconate Dehydrogenase in diagnostic tests for the determination of creatine kinase or glucose in the combination with Hexokinase, Catalog Nos. 11 119 796 103 or 10152 676 103, Glucose-6-phosphate Dehydrogenase, Catalog Nos. 10 186 783 103, 11 389 343 103, 11 293 206 103 or 10 190 454, and 6-Phosphogluconolactonase, Catalog No. 11 373 129 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 1.1.1.44

Properties

Nomenclature: 6-phospho D-gluconate:NADP⁺ 2-oxidoreductase (decarboxylating)

Molecular weight: 150 kD (native), 47 kD (SDS-PAGE)

Michaelis constants (TEA, 0.1 mol/l, pH 7.6, +25°C):

6-Phosphogluconate: 7.1 x 10⁻⁵ mol/l

NADP: 1.3 x 10⁻⁴ mol/l

Inhibitor constant (Phosphate buffer, pH 7.5):

Pyridoxal-5-P: 4.3 x 10⁻⁵ mol/l competitive


Inhibitors: Pyridoxal-5-P, iodoacetate and 4-hydroxymercuribenzoate

Activators: Chelators (EDTA, cysteine) plus metal ions (Mg²⁺); NaCl (0.2

Cat. No. **Pack Size**

11 126 482 103 custom fill

Will be supplied as "6-PGDH from Yeast, Lyophilizate". Unit of Measure is "kU".

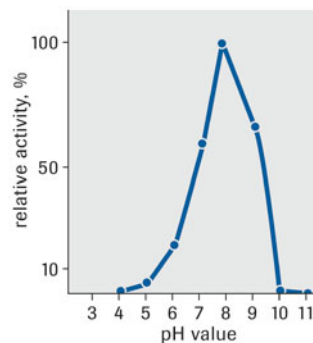

pH optimum

mol/l), KCl (0.2 mol/l).

pH optimum: 7.8 (see figure)

Temperature dependence: See figure

pH stability: 5.0-8.0 (see figure)

Thermal stability: Up to +45°C (see figure)

Specificity: Phosphogluconate dehydrogenase is specific for NADP; NAD does not react.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 7.0-8.0

Protein (Biuret): 0.08-0.16 mg/mg lyophilizate

Activity (+25°C, gluconate-6-P): ≥2 U/mg lyophilizate

Specific activity: ≥12 U/mg protein

Contaminants (expressed as percentage of 6-Phosphogluconate Dehydrogenase activity):

Creatine kinase: ≤0.006


G6P-DH: ≤0.01

Glutathione reductase: ≤0.01


Hexokinase: ≤0.01

Phosphoglucose isomerase: ≤0.03


Stability: At +2 to +8°C within specification range for 12 months. Store dry.

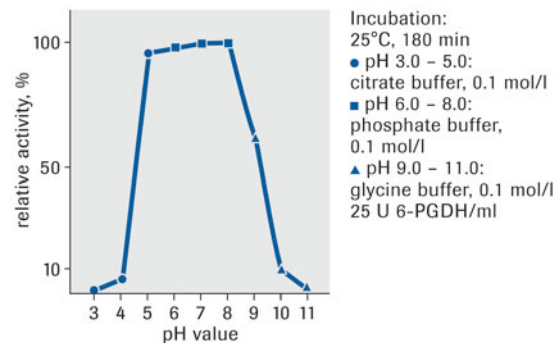

pH stability

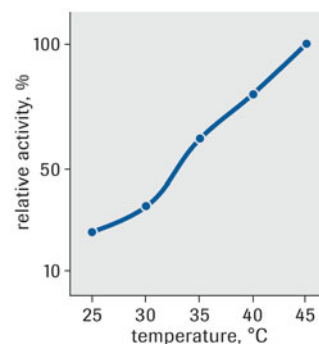

Temperature dependence

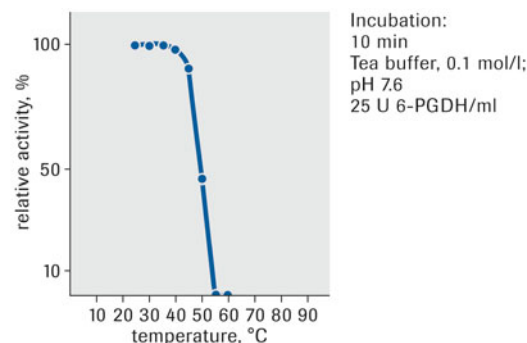

Thermal stability

For further processing only.

6-Phosphogluconolactonase from *Leuconostoc mesenteroides*, lyophilizate

Hydrolase that catalyzes the conversion of 6-phosphogluconolactone to 6-phosphogluconate.

Application

Use 6-Phosphogluconolactonase in diagnostic tests for the determination of creatine kinase or glucose in the combination with Hexokinase, Catalog Nos. 11 119 796 103 or 10152 676 103, Glucose-6-phosphate Dehydrogenase, Catalog Nos. 10 186 783 103, 11 389 343 103, 11 293 206 103 or 10 190 454 103, and Phosphogluconate Dehydrogenase, Catalog No. 11 126 482 103.

Cat. No. **Pack Size**

11 373 129 103 custom fill

Will be supplied as "6-Phosphogluconolactonase". Unit of Measure is "kU".



Benefits

Rely on the proven diagnostic quality of this product.

EC 3.1.1.31

Properties

Nomenclature: 6-phosphogluconolactonase

Molecular weight: 38 kD (SDS)

Isoelectric point: 6.0

Michaelis constants (Mes buffer, pH 6.5; +25°C):

6-Phosphogluconolactone: $<1 \times 10^{-7}$ mol/l

Inhibitors: $(\text{NH}_4)_2\text{SO}_4$ (> 20 mmol/l), Mg^{2+} (> 10 mmol/l), NaCl (> 10 mmol/l). The enzyme is not inhibited by Cu^{2+} , Zn^{2+} , EDTA, 5.5'-dithiobis-2-nitrobenzoic acid, octanol (0.01%), Triton X-100 (1%) and Thesit (1%).

pH optimum: 6.0-7.5 (see figure)

Temperature dependence: See figure

pH stability: 7.0-9.0 (at +4°C, see figure)

Thermal stability: Up to +50°C (see figure)

Stability at different ionic strength: See figure

Specificity: 6-Phosphogluconolactone 100%, gluconolactone 0.5%

Specification

Appearance: White lyophilizate

Activity (+25°C, 6-phosphogluconolactone): ≥ 50 U/mg lyophilizate

Contaminants (expressed as percentage of 6-Phosphogluconolactonase activity):

Creatine kinase: ≤ 0.001

G6P-DH: ≤ 0.02

Myokinase: ≤ 0.001

"NADPH oxidase": ≤ 0.001


6- Phosphogluconate dehydrogenase: ≤ 0.01

Function testing (with G6P-DH, reaction time up to 5 minutes): $\geq 98\%$


Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Literature

- 1) H.P. Bauer, T. Shihari, J.C. Jochim, H.W. Hofer, *Europ. J. Biochem.* **133**, 163-168 (1983)
- 2) R.K. Scopes, *FEBS Lett.* **193/2**, 185-188 (1985)
- 3) R. Khanna, A.R. Data, J.L. Rosner, *Plasmid* **17**, 76-82 (1987)
- 4) R.D. Moir, G.B. Stokes, *Biochem. J.* **256**, 69-73 (1988)
- 5) R. Vormbrock, R. Helger, *Enzyme* **38/1**, 20-21 (1987)

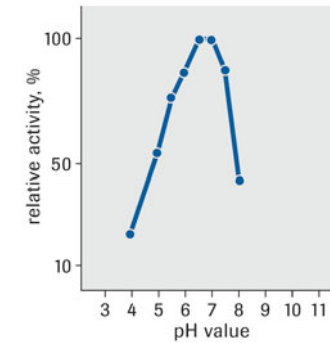


pH optimum

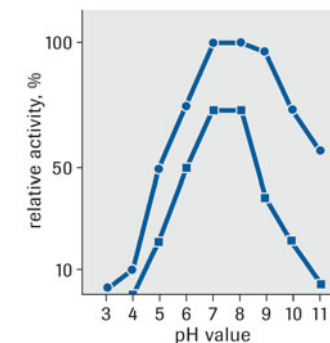


pH stability


Incubation:
 ● 4°C, 24 h
 ■ 25°C, 24 h
 pH 3.0 - 5.0: citrate buffer, 0.05 mol/l
 pH 6.0 - 8.0: phosphate buffer, 0.05 mol/l
 pH 9.0 - 11.0: glycine buffer, 0.05 mol/l
 1 U 6-PGL/ml

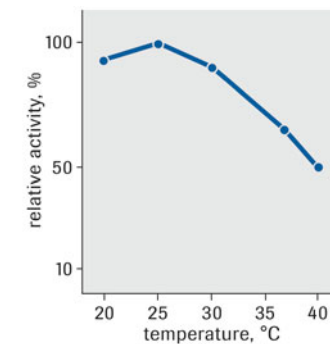

Temperature dependence

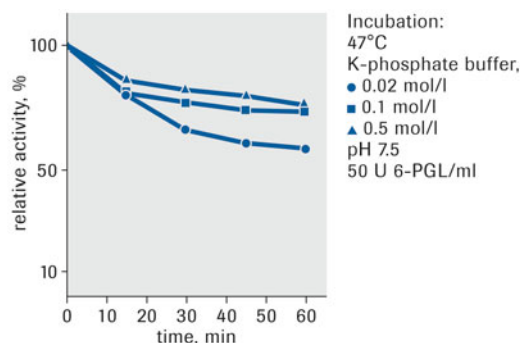

Stability at different ionic strength

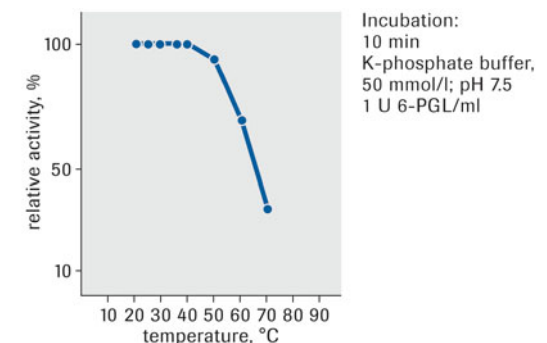

Thermal stability

For further processing only.

Pyruvate Kinase

from *Bacillus stearothermophilus*, lyophilizate

Application

Use Pyruvate Kinase to catalyze the transfer of a phosphate group from phosphoenolpyruvate (PEP) to ADP, *e.g.*, for the enzymatic determination of potassium or triglycerides.

Benefits

Rely on the proven diagnostic quality of this product.

EC 2.7.1.40

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=10 mg/ml)

pH value (+25°C, c=10 mg/ml): 7.6-8.7

Activity (+37°C, PEP): ≥120 U/mg lyophilized material

Specific activity: ≥150 U/mg protein

Protein (Biuret): 0.55-0.95 mg/mg lyophilizate

Contaminants (expressed as percentage of Pyruvate Kinase activity):

"NADH oxidase" (dA_{365} , 48 hours): ≤0.060

Na (flame photometric): ≤0.60 mmol/KU

K (flame photometric): ≤30 mmol/KU

NH₄ (enzymatically): ≤0.13 mg/KU

Mg (AAS): ≤50 mmol/KU

Mn (AAS): ≤2.4 mmol/KU

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.

11 462 652 103

Pack Size

custom fill

Will be supplied as "PK, Bacillus stearothermoph., Lyo.". Unit of Measure is "MU".

Contents

100 mg lyophilized material contains 60 mg protein and 40 mg Tris

For further processing only.

Pyruvate Kinase from rabbit muscle, suspension

Cat. No.	Pack Size
10 005 533 103	custom fill

Application

Use Pyruvate Kinase to catalyze the transfer of a phosphate group from phosphoenolpyruvate (PEP) to ADP, *e.g.*, for the enzymatic determination of potassium or triglycerides.

Benefits

Rely on the proven diagnostic quality of this product.

EC 2.7.1.40

Specification

Appearance: White suspension in ammonium sulfate, 3.2 mol/l, pH approximately 6

pH value: 5.5 to 6.5

Activity (+25°C, PEP): ≥2000 U/ml solution

Specific activity: ≥200 U/mg protein

Protein (Biuret): ≥10 mg/ml

Ammonium sulphate: 3.2±0.2 mol/l

Contaminants (expressed as percentage of Pyruvate Kinase activity):

ATPase: ≤0.002

Enolase: ≤0.01

Glycerokinase: ≤0.001

Hexokinase: ≤0.002

Lactate dehydrogenase: ≤0.01

Myokinase: ≤0.01

"NADH oxidase": ≤0.002

Glycerol (enzymatically): ≤10 µg/10 mg

Stability: At +2 to +8°C within specification range for 24 months.

Will be supplied as "Pyruvate Kinase (PK) from Rabbit Muscle".
Unit of Measure is "MU".
For further processing only.

Pyruvate Oxidase from *E.coli* overproducer, lyophilizate

Recombinant oxidoreductase that catalyzes the interconversion of pyruvate to acetyl phosphate.

Application

Use Pyruvate Oxidase in a variety of diagnostic tests, such as for the determination of pyruvate, lactate or aminotransferases.

Benefits

Rely in the proven diagnostic quality of this product.

EC 1.2.3.3

Specification

Appearance: Yellow lyophilizate

Solubility: Clear, yellowish solution in potassium phosphate buffer, 0.1 M, pH 6.5 (c=10 mg/ml)

Activity (+25°C, pyruvate, O₂, P_i): ≥1.5 U/mg lyophilizate

Specific activity: ≥3 U/mg protein

Protein (Biuret): ≥0.4 mg/mg lyophilizate

Contaminants (expressed as percentage of Pyruvate Oxidase activity):

ATPase: No limit

Glucose oxidase: ≤0.001

Aspartate aminotransferase (AST/GOT): ≤0.01

Alanine aminotransferase (ALT/GPT): ≤0.01

apo-Alanine aminotransferase (apo-ALT/apo-GPT) : ≤0.005

Impurities, total: ≤0.02

Lactate oxidase: ≤0.002

α-Ketoglutarate oxidase: ≤0.02

"NADH oxidase": ≤0.02

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

11 418 912 103 custom fill

Will be supplied as "Pyruvate Oxidase Recombinant (E. coli)". Unit of Measure is "kU".



For further processing only.

Sarcosine Oxidase from *E.coli* overproducer, lyophilizate

Oxidoreductase that catalyzes the demethylation of sarcosine to glycine.

Application

Use Sarcosine Oxidase in diagnostic tests for the determination of creatinine. This can be done using one of two methods:

- (1) In combination with Creatinase, Catalog No. 11 799 142 103 and Creatininase, Catalog No. 11 865 471 103.
- (2) In combination with Creatinine Deaminase, Catalog No. 11 330 764 103, N-Carbamoylsarcosine Amidase, Catalog No. 11 248 847 103 and N-Methylhydantoinase (ATP-hydrolysing), Catalog No. 11 288 555 103.

Benefits

Use Sarcosine Oxidase in your preferred creatinine reagent mix and rely on the proven diagnostic quality of this product.

EC 1.5.3.1

Properties

Nomenclature: Sarcosine:oxygen oxidoreductase (demethylating)

Molecular weight: 40 kD (PAGE, native Phast®-System)

Isoelectric point: 5.3 (Phast®-System)

Michaelis constants (Tris buffer, 0.1 mol/l, pH 8.0; Sarcosine):
at +25°C: 3.7×10^{-3} mol/l

at +37°C: 6.3×10^{-3} mol/l

Inhibitors: Completely inhibited by $ZnCl_2$ (7 mmol/l), $CdCl_2$ (7 mol/l), heavy metals and NaN_3 . Chloroacetic amine (0.2%) does not inhibit.

pH optimum: 8.0 (see figure)

Temperature dependence: See figure

pH stability: 7.0-10.0 (see figure)

Thermal stability: Up to +50°C (see figure)

Specificity: Sarcosine Oxidase reacts with sarcosine (100%), N-ethylglycine, 2 mmol/l (4%), L(-)-proline (0.28%), carbamoylsarcosine (0%), and glycine (0%).

Specification

Appearance: Yellow lyophilizate

Solubility: Clear, yellow solution in water (c=10 mg/ml)

pH value (c=10 mg/ml in water): 7.5-8.5

Protein (Biuret): 0.4-0.6 mg/mg lyophilizate

Activity (+25°C, sarcosine): 22-40 U/mg lyophilizate

Specific activity: ≥ 45 U/mg protein

Contaminants (expressed as percentage of Sarcosine Oxidase activity):

ATPase: ≤ 0.01

Catalase: ≤ 10.0

Contaminating oxidases (FOX): ≤ 0.005

Creatinase: ≤ 0.001

Creatininase: ≤ 0.01

Creatinine deaminase: ≤ 0.001

N-Carbamoylsarcosine amidohydrolase: ≤ 0.001

N-Methylhydantoinase: ≤ 0.001

Stability: At -15 to -25°C within specification range for 12 months. Store dry. Protect from light.


Literature

- 1) N. Mori, M. Sano, Y. Tani, H. Yamada, *Agr. Biol. Chem.* **44** (6), 1391 (1980)
- 2) M. Suzuki, *J. Biochem.* **89**, 599 (1981)
- 3) S. Hayashi, M. Suzuki, S. Nakamura, *Biochem. Int.* **4**, 617 (1982)
- 4) S. Hayashi, S. Nakamura, M. Suzuki, *Biochem. Biophys. Res. Com.* **96**, 924 (1980)


Cat. No. **Pack Size**

11 378 856 103 custom fill


Will be supplied as "Sarcosine Oxidase, Recombinant (E. coli)".
Unit of Measure is "kU".

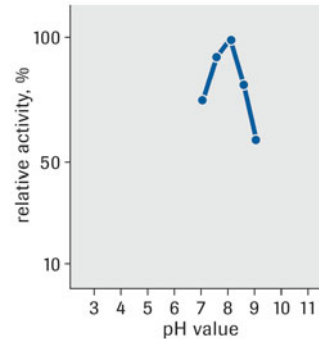


pH optimum

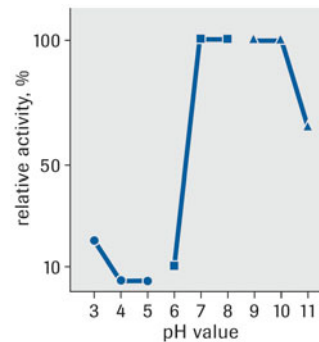

pH stability

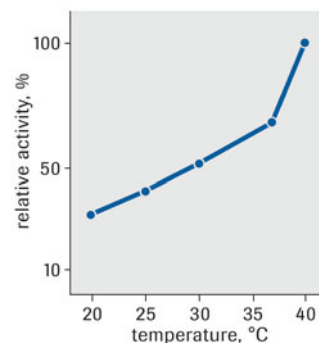

Temperature dependence

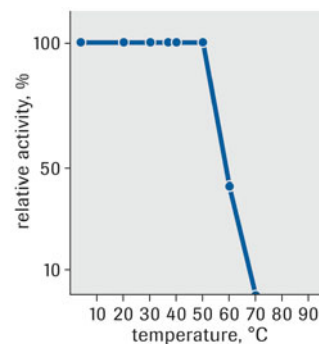

Thermal stability

For further processing only.

Incubation:
25°C, 360 min
● pH 3.0 – 5.0: citrate buffer, 50 mmol/l
■ pH 6.0 – 8.0: phosphate buffer, 50 mmol/l
▲ pH 9.0 – 11.0: glycine buffer, 50 mmol/l
10 U sarcosine OD/ml

Incubation:
10 min
phosphate buffer, 0.1 mol/l; pH 8.0
10 U sarcosine OD/ml

Triose-phosphate Isomerase from rabbit muscle, suspension

Isomerase that interconverts dihydroxyacetone phosphate and D-glyceraldehyde 3-phosphate.

Application

Use Triose-phosphate Isomerase in diagnostic reagents for the determination of aldolase in combination with Glycerol-3-phosphate Dehydrogenase, Catalog No. 10 151 351 103, and Fructose-1,6-diphosphate, Catalog. No. 10 041 793 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 5.3.1.1

Specification

Appearance: White suspension in ammonium sulfate

pH value: 5.5-6.5

Specific activity (+25°C, glyceraldehyde-3-phosphate): ≥5,000 U/mg protein

Protein (Biuret): 10±1 mg/ml

Contaminants (expressed as percentage of Triose-phosphate Isomerase activity):

Aldolase: ≤0.01

Glyceraldehyde-3-phosphate dehydrogenase: ≤0.001

Glycerol-phosphate dehydrogenase: ≤0.01

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No. **Pack Size**

10 153 338 103 custom fill

Will be supplied as "Triosephosphate Isomerase, Rabbit Muscle". Unit of Measure is "MU".

For further processing only.

Thrombin from human plasma, lyophilizate

Plasma derived coagulation factor II a that selectively cleaves the Arg--Gly bonds of fibrinogen to form fibrin.

Application

Use Thrombin to generate reference antigens for anti D-dimer antibodies.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.4.21.5

Specification

Appearance: White lyophilizate

Specific activity (+25°C, Chromozym TH): ≥120 U/mg protein

Protein (Lowry): Approximately 0.004 mg/mg lyophilizate

Factor Xa: ≤3 %

Anti HIV: Negative

HBsAg: Negative

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No. **Pack Size**

10 582 514 103 custom fill

Will be supplied as "Thrombin (Coagulation Factor II a)". Unit of Measure is "U".

For further processing only.

Urease from jack bean, lyophilizate

Hydrolase that catalyzes the breakdown of urea in carbon dioxide and ammonia.

Application

Use Urease in diagnostic tests for the determination of urea in combination with Glutamate Dehydrogenase, Catalog Nos. 10 190 462 103 or 11 434 993 103.

Benefits

Rely on the proven diagnostic quality of this product.

EC 3.5.1.5

Properties

Nomenclature: Urea amidohydrolase

Molecular weight: 480 kD

Isoelectric point: 5.0-5.1

Michaelis constant (Phosphate buffer, pH 7.0; +25°C):

Urea: 1.05×10^{-2} mol/l

Inhibitors: Na⁺, K⁺, NH₄⁺; suramin and thiourea are competitive inhibitors.

Activators: P_i

pH optimum: 7.5 (see figure)

Temperature dependence: See figure

pH stability: 6.0-9.5 (see figure)

Thermal stability: Up to +70° C (see figure)

Specificity: Urease is specific for urea.

Specification

Appearance: Almost white lyophilizate

Solubility: Clear, colorless solution in water (c=20 mg/ml)

pH value (c=20 mg/ml in water): 6.0-7.0

Activity (+25° C, urea): ≥45 U/mg lyophilizate

Specific activity: ≥600 U/mg protein

Protein (Biuret): ≤0.15 mg/mg lyophilizate

Contaminants (expressed as percentage of Urease activity):

L-Arginase: ≤0.002


NH₄⁺: ≤1.5 µg/KU

Stability: At +2 to +8°C within specification range for 12 months. Store dry.


Cat. No. **Pack Size**

11 759 132 103 custom fill


Will be supplied as "Urease, Lyo., SQ". Unit of Measure is "MU".

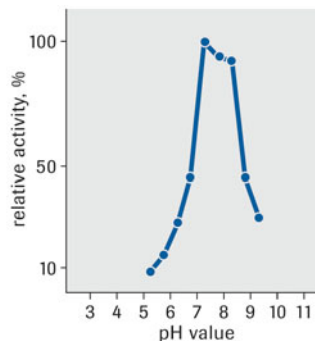


pH optimum

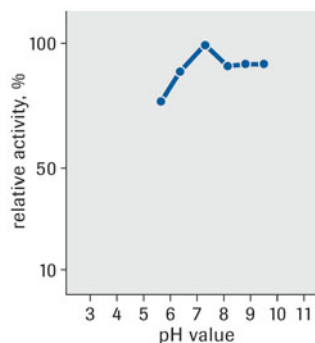

pH stability

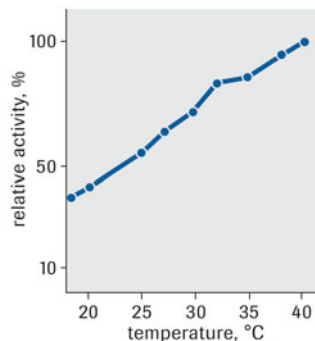

Temperature dependence

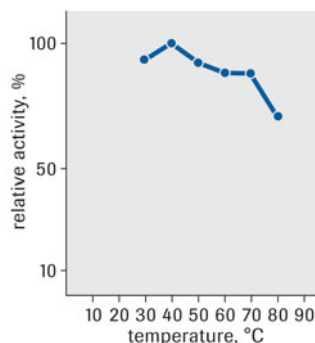

Thermal stability

For further processing only.

Incubation:
25°C, 24 h
K-phosphate buffer,
0.1 mol/l
urease,
approx. 45 U/ml

Incubation:
15 min
K-phosphate buffer,
0.1 mol/l; pH 7.0
urease,
approx. 45 U/ml

Uricase

from *Arthrobacter protophormiae*, lyophilizate

Oxidase that catalyzes the oxidation of uric acid to 5-hydroxyisourate which decomposes to allantoin under *in vitro* conditions.

Application

Use Uricase in diagnostic tests for the determination of uric acid and for the elimination of uric acid interferences.

Benefits

- Eliminate uric acid interferences.
- Rely on the proven diagnostic quality.

EC 1.7.3.3

Properties

Nomenclature: Urate:oxygen oxidoreductase

Molecular weight: ~170 kD, with four subunits of ~40 kD

Michaelis constant (Phosphate buffer, 0.1 mol/l, pH 8.0; +25°C):

Urate: 6.6×10^{-5} mol/l

Stabilizer/activators: EDTA is good for stabilization. DTT or DTE may show a stabilizing effect depending on reagent composition. Triton X-100 (1-2 ml/l) may show an activating effect.

Inhibitors: Zn^{2+} , Cl^- (Tris-HCl buffer is not suitable) and borate inhibit strongly. NaN_3 , 0.1% does not inhibit.

pH optimum: 9.0 (see figure). Roche uric acid reagent contains phosphate buffer, pH 7.8. For these conditions, high activity and higher stability of the Uricase are achieved.

Temperature dependence: See figure, above +50°C there is a decrease in activity due to lower O_2 concentration.

pH stability: 6.5-10.0 (see figure)

Thermal stability: +20°C to +60°C (see figure)

Specificity: Uricase is specific for urea.

Specification

Appearance: White lyophilizate

pH value (c=10 mg/ml in water): 6.7-7.5

Protein (Biuret): No limit

Activity (+25°C, urate, saturated with O_2 , pH 8.5): ≥ 20 U/mg lyophilizate

Specific activity: ≥ 50 U/mg protein

Contaminant (expressed as percentage of Uricase activity):

Catalase: ≤ 0.5


Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.


Pack Size

10 828 475 103 custom fill

Will be supplied as "Uricase from *Arthrobact. protophormiae*". Unit of Measure is "MU".

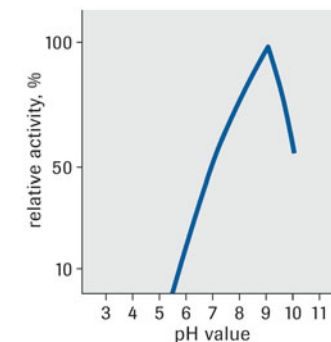


pH optimum

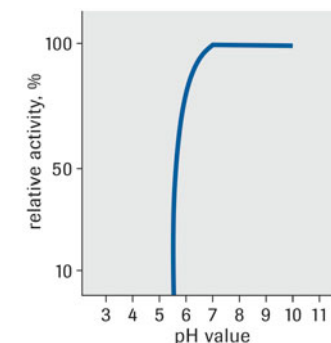


pH stability

Incubation:
25°C, 24 h
pH 5.5 - 10.0:
potassium phosphate
buffer, 0.1 mol/l
40 U uricase/ml

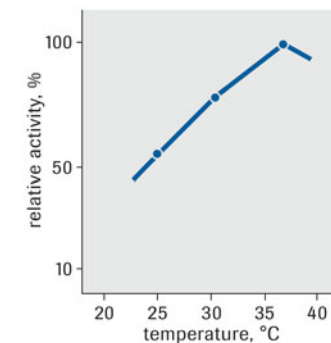

Temperature dependence

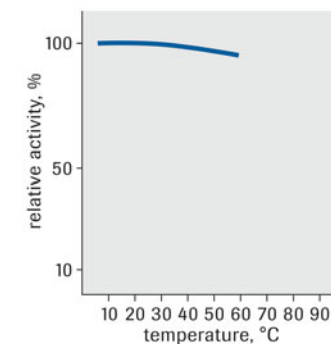

Thermal stability

Incubation:
15 min
phosphate buffer,
0.1 mol/l;
pH 8.0
40 U uricase/ml

For further processing only.

3,5-Dichlorophenolsulfonic Acid disodium salt

Color reagent for diagnostic tests

Application

Use 3,5-Dichlorophenolsulfonic Acid instead of phenol as a component in the trinder reaction.

Benefits

- Take advantage of the higher molar absorptivity compared to phenol.

CAS: 95041-38-6

Properties

Nomenclature: 3,5-Dichloro-2-hydroxy-benzosulphonic acid disodium salt

Formula: $C_6H_2O_4C_{12}SNa_2$

Molecular weight: 287.0 D

Specification

Appearance: White powder

Solubility: Clear, colorless solution in water (c=20 mg/ml)

Dichlorophenolsulphonic acid-Na (from C): $\geq 98.0\%$

C (elementary analysis): 24.6-25.6%

H (elementary analysis): 0.69-0.80%

Thin layer chromatography (TLC): Chromatographically homogeneous; corresponds to reference

Stability: At +15 to +25°C within specification range for 36 months.

Cat. No.	Pack Size
10 667 536 103	custom fill

Will be supplied as "3,5-Dichlorophenol Sulfonic Acid, Di-Na". Unit of Measure is "g".

For further processing only.

4-Aminoantipyrine (4-APP) crystalline powder

Substrate for peroxidase

Application

Use 4-Aminoantipyrine in a variety of diagnostic tests that use the Trinder reaction for the colorimetric determination of analytes, such as for the determination of cholesterol, glucose, creatinine or uric acid.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 83-07-8

Properties

Formula: $C_{11}H_{13}N_3O$

Molecular weight: 203.25 D

Specification

Appearance: Yellow to redish brown, crystalline powder

Solubility: Clear, colorless solution in water (c=0.1%, w/v)

Melting range: +106 to +109°C

A 370 (c=0.1%, w/v): ≤ 0.04

UV-spectrum (c=0.002%, w/v): Corresponds to reference

Maximum: 242-246 nm

Minimum: 217-219 nm

Shoulder: 274 nm

Heavy metals (as Pb): ≤ 5 ppm \pm 0.0005%

IR-spectrum: Corresponds to reference

4-Aminoantipyrine (HClO₄-titration, based on undried substance): $\geq 98.0\%$

Cat. No.	Pack Size
10 073 474 001	custom fill

Will be supplied as "4-Aminoantipyrine". Unit of Measure is "kg".

For further processing only.

Purity (HPLC): ≥99.0 area%

Stability: At +15 to +40°C within specification range for 36 months. Store dry in tightly closed containers.

4-Nitrophenyl-α-D-maltohexaoside powder

Nitrophenyl substrate

Application

Use 4-Nitrophenyl-α-D-maltohexaoside in diagnostic tests for the determination of α-amylase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 74173-30-1

Properties

Formula: C₄₂H₆₅NO₃₃

Molecular weight: 1112.1 D

Specification

Appearance: White to slightly yellowish, amorphous powder

Solubility: Clear, slightly yellowish solution in water (c=70 mg/ml)

4-Nitrophenyl-maltohexaoside (enzymatically): ≥90%

4-Nitrophenyl-maltohexaoside (HPLC): ≥96.0 area%

Water (K. Fischer): ≤3.0%

4-Nitrophenyl-maltopentaoside (HPLC): ≤1.0 area%

4-Nitrophenyl-maltoheptaoside (HPLC): ≤2.0 area%

4-Nitrophenol, free: ≤0.05%

2-Propanol (GC): ≤6%

Stability: At +2 to +8°C within specification range for 18 months. Store dry.

Cat. No. **Pack Size**

10 691 682 103 custom fill

Will be supplied as "4-Nitrophenyl-α-D-maltohexaoside". Unit of Measure is "g".

For further processing only.

4-Nitrophenyl-α-D-maltopentaoside powder

Nitrophenyl substrate

Application

Use 4-Nitrophenyl-α-D-maltopentaoside in diagnostic tests for the determination of α-amylase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 66068-38-0

Properties

Formula: C₃₆H₅₅NO₂₈

Molecular weight: 949.9 D

Specification

Appearance: White or slightly yellowish, amorphous powder

Solubility: Clear, slightly yellowish solution in water (c=70 mg/ml)

4-Nitrophenyl-maltopentaoside (enzymatically): ≥90%

4-Nitrophenyl-maltopentaoside (HPLC): ≥98.0 area%

Water (K. Fischer): ≤2.0%

4-Nitrophenyl-maltotetraoside (HPLC): ≤0.4 area%

Cat. No. **Pack Size**

10 691 747 103 custom fill

Will be supplied as "4-Nitrophenyl-α-D-maltopentaoside". Unit of Measure is "g".

For further processing only.

4-Nitrophenyl-maltohexaoside (HPLC): ≤0.6 area%
4-Nitrophenol, free: ≤0.05%
2-Propanol (GC): ≤6%
Stability: At +2 to +8°C within specification range for 18 months. Store dry.

4-Nitrophenyl Phosphate (4-NPP) di-Tris salt

Substrate for alkaline phosphatase

Application

Use 4-Nitrophenyl Phosphate in diagnostic test for the determination of alkaline phosphatase according to the recommendations of the International Federation of Clinical Chemistry and Laboratory Medicine (IFCC).

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 68189-42-4

Properties

Formula: $C_6H_4NO_6P(C_4H_{12}NO_3)_2$

Molecular weight: 461.3 D (4-NPP: 219.1 D)

Specification

Appearance: White to slightly yellow, crystalline powder

Solubility: Clear, colorless to slightly yellow solution in water

4-NPP, di-Tris (from content found enzymatically): ≥88%

4-NPP (enzymatically): ≥42%

Tris (titrimetric): ≥46%

Water (K. Fischer): ≤6%

4-NP free: ≤0.07%

Reaction rates (alkaline phosphatase): 100±5%

Stability: At +2 to +8°C within specification range for 24 months. Store dry. Protect from light.

Cat. No.	Pack Size
10 270 857 103	custom fill

Will be supplied as "4-Nitrophenyl Phosphate, Di-Tris Salt". Unit of Measure is "kg".

For further processing only.

ABTS

2,2'-Azino-di-[3-ethylbenzthiazoline sulfonate (6)] diammonium salt

Substrate for peroxidase

Application

Use ABTS for activity measurements of peroxidase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 30931-67-0

Properties

Formula: $C_{18}H_{16}N_4O_6S_4(NH_4)_2$

Molecular weight: 548.5 D (ABTS: 514.6 D)

POD changes absorbance with ABTS from colorless to green at 405 nm.

Specification

Appearance: Green crystals

Solubility: Clear, slightly green solution in water (c=20 mg/ml)

ABTS-(NH₄)₂ (A_{340} , $\epsilon=40.28$ [l x mmol⁻¹ x cm⁻¹]): ≥97.5%

Cat. No.	Pack Size
10 122 661 103	custom fill

Will be supplied as "ABTS, Diammonium Salt". Unit of Measure is "kg".

For further processing only.

Thin layer chromatography: Chromatographically homogeneous
Stability: At +2 to +8°C within specification range for 24 months. Keep under nitrogen or argon. Protect from light.

Benzylidene-4-NP-G7

4,6-Benzylidene-4-nitrophenyl- α -D-maltoheptaoside, lyophilizate

Nitrophenyl substrate

Application

Use Benzylidene-4-NP-G7 in diagnostic tests for the determination of α -amylase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 109055-07-4

Properties

Formula: C₅₅H₇₉NO₃₈

Molecular weight: 1362.1 D

Specification

Appearance: White to slightly yellowish lyophilizate

Benzylidene-4-NP-G7 (enzymatically): $\geq 90\%$

Water (K. Fischer): $\leq 3\%$

4-NP-maltoheptaoside (HPLC): ≤ 1.0 area%

4-Nitrophenol, free: $\leq 0.01\%$

Reaction rates (α -amylase):

In Precinorm (R) U: $100 \pm 5\%$

In Precipath (R) U: $100 \pm 5\%$

Stability: At -15 to -25°C within specification range for 18 months. Store dry.

Cat. No. **Pack Size**

11 378 872 103 custom fill

Will be supplied as "Benzylidene-4-NP-G7". Unit of Measure is "kg".



For further processing only.

Chromogenic Substrate for Lipase

Substrate for lipase

Application

Use Chromogenic Substrate for Lipase in diagnostic tests for the determination of lipase activity.

Benefits

- Design the best possible diagnostic test for lipase activity.
- Rely on the proven diagnostic quality of this product.

CAS: 195833-46-6

Properties

Nomenclature: 1,2-O-Dilauryl-rac-glycero-3-glutaric acid-(6'-methylresorufin)ester

Formula: C₄₅H₆₉NO₈

Molecular weight: 752.05 D

$\lambda_{\text{max,substrate}}$: 470 nm (Tris-HCl, pH 8.4)

ϵ_{470} : 57.94 [l x mmol⁻¹ x cm⁻¹]

$\lambda_{\text{max,methylresorufin}}$: 581 nm (Tris-HCl, pH 8.4)

Melting range: +29 to +31°C

pH optimum: 7.0-9.5

Cat. No. **Pack Size**

11 034 618 103 custom fill

Will be supplied as "Chromogenic Substrate for Lipase". Unit of Measure is "g".

For further processing only.

Solubility: Soluble in polar organic solvents, *e.g.*, n-propanol, ethyl acetate, dioxane, methanol, dimethyl sulfoxide. The limit of solubility in n-propanol is 42.9 mg/ml.

Specification

Appearance: Red, smear substance

Chromogenic Lipase Substrate (from C): $\geq 95\%$

Chromogenic Lipase Substrate (HPLC): ≥ 95 area%

C (elementary analysis): 68.2-72.5%

H (elementary analysis): 8.7-9.7%

N (elementary analysis): 1.4-2.4%

Methylresorufin, free (HPLC): ≤ 0.5 area%

Isomer (HPLC): ≤ 2 area%

Stability: At +2 to +8°C within specification range for 36 months.

Background Information

Reagent proposal for lipase test

The sensitivity of this lipase test is especially influenced by the extinction of the substrate solution, the concentration of taurodesoxycholate, the pH value, and the molarity of the Tris buffer.

Final test concentration of the substrate solution:

chrom. Lipase Substrate: 0.24 mmol/l (=180 mg/l)

colipase: 0.98 mg/l

taurodesoxycholate *(see buffer solution)

CaCl₂: 0.1 mol/l

tartrate buffer, pH 4.0: 1.6 mmol/l

stabilizers

The lipase substrate has to be dissolved in a small quantity of an organic solvent (*e.g.* n-propanol) first. Under vigorous stirring this organic solution should be injected into the tartare buffered aqueous solution with a thin beam. (The lipase substrate starts to hydrolyze at alkaline pH values.) The lipase substrate containing solution should be a micro-emulsion with an extinction of about 0.5 E. (Lower extinction of the reagent results in measurement of nonspecific serum esterases.) Stabilizers like mannitol, polywax 4000 and co-emulsifiers like lecithin, phosphoryl choline or dilauryl-glycerol-sulfate improve the stability of the micro-emulsion.

Final test concentration of the buffer solution:

Tris-HCl, pH 8.4: 41 mmol/l

taurodesoxycholate: *7.2 mmol/l (total)

desoxycholate: 1.77 mmol/l

Wavelength: 578 nm or 580 nm

Temperature: +25°C, +30°C or +37°C, respective

Buffer solution: 1 ml

Substrate solution: 0.2 ml

Sample volume: 0.02 ml

Assay time: 2 to 10 minutes

Start of reaction: with substrate

Literature

- 1) Neumann U, Junius M, Maier B. A sensitive colorimetric assay for the kinetic Lipase determination in serum (Boehringer Mannheim/now Roche Diagnostics). Abstract 13th Int. Congress for Clin Chem (ICCC), Den Haag, Netherlands, 28.6-3.7; 1987
- 2) Panteghini M, Bonora R, Pagani F. Measurement of pancreatic lipase activity in serum by a kinetic colorimetric assay using a new chromogenic substrate, Ann. Clin. Biochem. (2001) Jul ;**38** (Pt 4), 365-70.

Ethylidene-4-NP-G7**Ethyliden-4-nitrophenyl- α -D-maltoheptaosid (EPS), powder**

Nitrophenyl substrate

Application

Use Ethylidene-4-NP-G7 in diagnostics tests for the determination of α -amylase and pancreatic α -amylase according to the recommendations of the International Federation of Clinical Chemistry and Laboratory Medicine (IFCC).

Benefits

- Minimize the value of your blank for the liquid α -amylase reagent.
- Rely on the proven purity and the very low content of free maltoheptaoside.

CAS: 96597-16-9

Properties

Formula: C₅₀H₇₇NO₃₈

Molecular weight: 1300.1 D

Specification

Appearance: White to slightly yellowish, amorphous powder

Solubility: Clear, slightly yellowish solution in water (c=70 mg/ml)

EPS (enzymatically): $\geq 90\%$

Water (K. Fischer): $\leq 3\%$

pNP-G7 (enzymatically): $\leq 0.1\%$

pNP, free: $\leq 0.01\%$

Reaction rates (α -amylase):

In Precinorm U: 95-105%

In Precipath U: 95-105%

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Literature

E. Rauscher *et al.*, Fresenius Z. Analyt. Chem. **324**, 304 (1986)

Cat. No.

10 880 078 103 custom fill

Pack Size

Will be supplied as "Ethylidene-4-NP-G7". Unit of Measure is "kg".
Additional products: OEM reagents for the determination of α -amylase and pancreatic amylase, as well as specific inhibitory antibodies. Catalog Nos. 11 543 598 103 and 11 543 601 103 for the detection of pancreatic amylase.

For further processing only.

Glupa-carboxylate monoammonium salt

Substrate for γ -glutamyltransferase

Application

Use Glupa-carboxylate in diagnostic tests for the determination of γ -glutamyltransferase, according to the recommendations of the International Federation of Clinical Chemistry and Laboratory Medicine (IFCC).

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 63699-78-5

Properties

Formula: C₁₂H₁₂N₃O₇NH₄

Molecular weight: 328.3 D

Specification

Appearance: White to yellowish crystalline powder

Solubility: Clear, yellow solution in water (c=100 mg/ml), tested for insolubles

pH value: 4.0-6.0

Molar rotation: $[\alpha]_{25/D} +32.0 \pm 2.0^\circ$

Cat. No.

10 413 151 103 custom fill

Pack Size

Will be supplied as "Glupa-carboxylate, Monoammonium Salt".
Unit of Measure is "kg".

For further processing only.

Melting range (Kofler): Approximately +170 to +180°C

Glupa-carboxylate, free acid (enzymatically): ≥87%

Glupa-carboxylate (HPLC): ≥99 area%

Water (K. Fischer): ≤6.2%

NH₄ (Nessler's reagent): 5.2±1%

5-Amino-2-nitrobenzoate (HPLC): ≤0.1 area%

α-Glupa-carboxylate (HPLC): ≤0.4 area%

Thin layer chromatography (silica gel F; n-butanol/glacial acetic acid/H₂O = 50/15/25; UV, with Nihydrin): Chromatographically homogeneous

A₄₀₅ (Glupa-carboxylate, 6 mmol/l): 0.65-0.80

Stability: At +2 to +8°C within specification range for 24 months. Protect from light.

Additional formulations (tablets) are available on request.

Literature

Schumann G, Bonora R, Ceriotti F, Féraud G, Ferrero CA, et al. IFCC primary reference procedures for the measurement of catalytic activity concentrations of enzymes at 37 degrees C. International Federation of Clinical Chemistry and Laboratory Medicine. Part 6. Reference procedure for the measurement of catalytic concentration of gamma-glutamyltransferase. Clin Chem Lab Med. 2002 Jul, **40** (7): 734-8.

TOOS

(N-Ethyl-N-(2-hydroxy-3-sulfopropyl)-m-toluidine, monosodium salt, dihydrate

Substrate for peroxidase

Application

Use TOOS together with 4-Aminoantipyrine in an indicator reaction using peroxidase to form a quinoneimine dye.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 679787-10-1

Properties

Nomenclature: Dihydrate (N-Ethyl-N-(2-hydroxy-3-sulfopropyl)-m-toluidine

Formula: C₁₂H₁₈NO₄SNa x 2 H₂O

Molecular weight: 331.37 D

Specification

Appearance: White to slightly bluish crystallize

TOOS, mono-Na x 2 H₂O (from C): ≥98.0%

C (elementary analysis): 42.6-46.5%

H (elementary analysis): 6.3-6.8%

N (elementary analysis): 3.8-4.7%

Water (K. Fischer): 7.0-12.0%

Heavy metals (as Pb): ≤20 ppm

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No.	Pack Size
11 650 670 103	custom fill

Will be supplied as "TOOS". Unit of Measure is "kg".

For further processing only.

Tribromo-hydroxybenzoic acid crystallizate

Color reagent for diagnostic tests

Application

Use Tribromo-hydroxybenzoic acid instead of phenol as a component in the trinder reaction.

Benefits

- Take advantage of the higher molar absorptivity compared to phenol.

CAS: 14348-40-4

Properties

Formula: C₇H₃O₃Br₃

Molecular weight: 374.8 D

Specification

Appearance: White crystallizate

Solubility: Clear, colorless solution in NaOH, 0.1 mol/l (c=2.09%, w/v)

Dissolving time (c=2.09%, w/v): 10-20 minutes in NaOH, 0.1 mol/l

A₄₀₅ (c=2.09%, w/v, in NaOH, 0.1 mol/l): ≤0.020

Melting range: +143 to +148°C

Thin-layer chromatography (TLC): Corresponds to reference

2,4,6-Tribromo-hydroxybenzoic acid (alkalimetrically): ≥98.0%

Stability: At +15 to +40°C within specification range for 24 months. Store dry in tightly sealed containers. Protect from light.

Cat. No.

Pack Size

10 755 745 103 custom fill

Will be supplied as "Tribrom-Hydroxybenzoic acid". Unit of Measure is "g".

For further processing only.

L(+)-Alanine crystalline powder

Substrate for alanine aminotransferase

Application

Use L(+)-Alanine in diagnostic tests for the determination of alanine aminotransferase (ALT).

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 56-41-7

Properties

Formula: C₃H₇NO₂

Molecular weight: 89.09 D

Solubility: Easily soluble in water and mineral acids, insoluble in organic solvents.

Specification

Appearance: White, crystalline powder or crystals

Solubility: Clear, colorless solution in phosphate buffer (c=0.9%, w/v, pH 7.4)

Microbiological test: Corresponds

Heavy metals (as Pb): ≤20 ppm ± 0.002%

Sulfate ash: <00.1%

Thin layer chromatography: Corresponds to reference

Water (K.Fischer): ≤1.0%

L-alanine (HClO₄ titration, based on anhydrous substance): 98.5-100.5%

L-alanine (enzymatically, based on anhydrous substance): 97.0-105.0%

Stability: At +15 to +40°C within specification range for 36 months. Store dry in tightly closed containers.

Cat. No.	Pack Size
10 136 921 103	custom fill

Will be supplied as "L(+)-Alanin". Unit of Measure is "kg".
For further processing only.

α-Ketoglutarate (2-Oxoglutarate) free acid

Substrate for transaminases and glutamate dehydrogenase

Application

Use α-Ketoglutarate in a variety of diagnostic tests, such as for the determination of alanine aminotransferase, aspartate aminotransferase, ammonia, urea and glutamate dehydrogenase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 328-50-7

Properties

Formula: C₅H₆O₅

Molecular weight: 146.1 D

Specification

Appearance: White crystallizate

Solubility: Clear, colorless solution in water (c=50 mg/ml)

A₄₀₅ (c=50 mg/ml in water, against water): ≤0.020

Melting range: +113 to +117°C

α-Ketoglutaric acid (enzymatically): ≥98%

Water (K. Fischer): ≤1%

NH₄ (enzymatically): ≤0.1%

Cat. No.	Pack Size
10 156 736 103	custom fill

Will be supplied as "α-Ketoglutaric Acid, Free Acid". Unit of Measure is "kg".

For further processing only.

Heavy metals (as Pb): ≤10 ppm

Bioburden: ≤100 CFU/g

Reaction rates (Glutamate pyruvate transaminase (ALT)): ≥95%

Reaction rates (Glutamate oxalacetate transaminase (AST)): ≥95%

Stability: At +15 to +25°C within specification range for 36 months.

α-Ketoglutarate (2-Oxoglutarate) disodium salt, dihydrate

Substrate in enzymatic reactions with glutamate dehydrogenase or transaminases

Application

Use α-Ketoglutarate in a variety of diagnostic tests, such as for the determination of glutamate dehydrogenase, ammonia, alanine- and aspartate aminotransferases and urea. The dihydrate formulation is well suited for dry chemistry tests.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 305-72-6

Properties

Formula: C₅H₈O₇Na₂

Molecular weight: 226.1 D (α-KG: 146.1 D)

Specification

Appearance: White, crystalline powder

Solubility: Clear, colorless solution in water, pH 7.3 (c=200 mg/ml)

A₄₀₅ (c=10 mg/ml, water; against water): ≤0.020

α-Ketoglutarate, salt (based on value found enzymatically): ≥97%

α-Ketoglutarate, free acid (enzymatically): ≥63%

Na (flame photometric): 20.5±1%

Water (K. Fischer): 15±2%

Stability: At +15 to +25°C within specification range for 24 months. Store dry.

Cat. No. **Pack Size**

10 040 584 103 custom fill

Will be supplied as "α-Ketoglutarate (α-Oxoglutarate), Di-Na". Unit of Measure is "kg".

Additional formulation: Crystallized free acid, Catalog No. 10 156 736

For further processing only.

α-Ketoglutarate (2-Oxoglutarate) disodium salt

Substrate in enzymatic reactions with glutamate dehydrogenase or transaminases

Application

Use α-Ketoglutarate in a variety of diagnostic tests, such as for the determination of glutamate dehydrogenase, ammonia, alanine- and aspartate aminotransferases and urea. The disodium formulation is well suited for liquid tests.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 305-72-6

Properties

Formula: C₅H₄O₅Na₂

Molecular weight: 190.1 D (α-KG: 146.1 D)

Cat. No. **Pack Size**

10 266 400 103 custom fill

Will be supplied as "α-Ketoglutarate, Di-Na, (M 190.1 g/mol)". Unit of Measure is "kg".

For further processing only.

Specification**Appearance:** White, crystalline powder**Solubility:** Clear, colorless solution in water, pH 7.3 (c=200 mg/ml)**A₄₀₅** (c=10 mg/ml in water, against water): ≤0.020**α-Ketoglutarate, salt** (based on value found enzymatically): ≥97.5%**α-Ketoglutarate, free acid** (enzymatically): ≥74%**Na** (flame photometric): 24±2%**Water** (K. Fischer): ≤2%**Heavy metals** (as Pb): ≤0.002%**Stability:** At +15 to +25°C within specification range for 24 months. Store dry.

Additional formulation crystallized free acid, Catalog No. 10 156 736

α-Ketoglutarate (2-Oxoglutarate) for potassium test free acid

Substrate for transaminases and glutamate dehydrogenase

Application

Use α-Ketoglutarate for enzymatic potassium tests especially to remove ammonia from the reaction.

Benefits

- Rely on the proven diagnostic quality of this product.
- Rely on the strongly reduced concentration of potassium.

CAS: 328-50-7**Properties****Formula:** C₅H₆O₅**Molecular weight:** 146.1 D**Specification****Appearance:** White crystallize**Solubility:** Clear, colorless solution in water (c=50 mg/ml)**A₄₀₅** (c=50 mg/ml in water, against water): ≤0.020**Melting range:** +113 to +117°C**α-Ketoglutaric acid** (enzymatically): ≥98%**Water** (K. Fischer): ≤1%**NH₄** (enzymatically): ≤0.1%**Na** (AES): ≤500 ppm**K** (AES): ≤10 ppm**Heavy metals** (as Pb): ≤10 ppm**Bioburden:** ≤100 CFU/g**Reaction rates** (Glutamate pyruvate transaminase (ALT)): ≥95%**Reaction rates** (Glutamate oxalacetate transaminase (AST)): ≥95%**Stability:** At +15 to +25°C within specification range for 36 months.**Cat. No.** **Pack Size****11 332 775 103** custom fill

Will be supplied as "α-Ketoglutaric Acid for Potassium Test". Unit of Measure is "kg".

For further processing only.

Creatine Phosphate disodium salt

Substrate for creatine kinase (reverse reaction)

Application

Use Creatine Phosphate in diagnostic tests for the determination of creatine kinase, according to the recommendations of the International Federation of Clinical Chemistry and Laboratory Medicine (IFCC).

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 922-32-7

Properties

Molecular weight: 327.2 D (CP: 211.1 D)

Specification

Appearance: White crystals

Solubility: Clear, colorless solution in water (c=150 mg/ml), free from fuzz

pH value (c=10 mg/ml in water): 7.7-8.7

Creatine-P-Na₂ x 4 H₂O (based on value found enzymatically): ≥97%

Creatine-P (enzymatically): ≥63%

Creatine-P (from P_{organic}): ≥63%

Na (flame photometric): 14±1%

Water (K. Fischer): 22±2%

P_{organic} (P_{total} - P_i): ≥9.25%

P_{total}: ≥9.25%

P_i (acid labile): ≤0.5%

P_i (Fiske and Subbarow): ≤1.5%

PP_i (enzymatically): ≤0.02%

ATP (enzymatically with hexokinase/G6P-DH): ≤0.002%

Sulfate (qualitative): Negative

Creatine, free: ≤0.5%

Glucose-6-P (enzymatically): ≤0.006%

PEP (enzymatically): ≤0.05%

Pyruvate (enzymatically): ≤0.02%

Kinetic of creatine kinase reaction: Corresponds to standard

Reaction rates (creatine kinase): 95-105%

A₃₃₄ (c=9 ml/ml water): ≤0.005

A₃₃₄ (against reaction mixture CK NAC active): ≤0.040

A₃₄₀ (hydrous solution): ≤0.120

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.

Pack Size

10 003 506 103 custom fill

Will be supplied as "Creatine Phosphate, Disodium Salt". Unit of Measure is "kg".

For further processing only.

D(-)-Lactate monolithium salt

Substrate for D-lactate dehydrogenase

Application

Use D(-)-Lactate as a standard in tests for lactic acid.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 27848-80-2

Properties

Formula: C₃H₅O₃Li

Molecular weight: 96.0 D (Lactate: 89.1 D, Lactic acid: 90.1 D)

Specification

Appearance: White, crystalline powder

D(-)-Lactate (enzymatically, as anion): ≥91%

Li (flame photometric): 7.0±1.0%

L(+)-Lactate (enzymatically, as anion): ≤0.2%

Stability: At +15 to +25°C within specification range for 36 months.

Cat. No. **Pack Size**

10 151 874 103 custom fill

Will be supplied as "D(-)-Lactate, Monolithium Salt". Unit of Measure is "g".

For further processing only.

Di(adenosine-5'-)penta-phosphate trilithium salt

Inhibitor of adenylate kinase

Application

Use Di(adenosine-5'-)penta-phosphate in diagnostic reagents for the determination of creatine kinase to inhibit adenylate kinase in the reaction.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 75522-97-3

Properties

Formula: C₂₀H₂₆N₁₀O₂₂P₅Li₃

Molecular weight: 934.2 D (Ap5A: 916.4 D)

Specification

Appearance: White to slightly yellowish, amorphous powder

Solubility: Clear, colorless to slightly yellowish solution in water (c=10 mg/ml)

Ap5A-Li₃ (from P_{organic}): ≥91%

Ap5A (A₂₆₀, ε=26.4 [l x mmol⁻¹ x cm⁻¹]): ≥90%

Ap5A (from P_{organic}): ≥90%

Ap5A (HPLC): ≥95 area%

Li (flame photometric): 2.1±0.3%

Water (K. Fischer): ≤5%

P_{total} (ammonium vanadate): ≥15.2%

P_i (ammonium vanadate): ≤1.5%

Thin layer chromatography (PEI-cellulose, KH₂PO₄, 0.75 mol/l): Chromatographically homogeneous

A₂₅₀/A₂₆₀: 0.79±0.04

A₂₈₀/A₂₆₀: 0.21±0.03

A₂₉₀/A₂₆₀: 0.02±0.02

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No. **Pack Size**

10 161 624 103 custom fill

Will be supplied as "Di(adenosine-5'-)penta-phosphate, Tri-Li". Unit of Measure is "g".

For further processing only.

Fructose-1,6-diphosphate trisodium salt

Substrate for aldolase and phosphatases

Application

Use Fructose-1,6-diphosphate for the determination of aldolase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 23558-08-9

Properties

Formula: $C_6H_{11}O_{12}P_2Na_3 \times 8 H_2O$

Molecular weight: 550.2 D (Fructose-1,6- P_2 : 340.1 D)

Specification

Appearance: White to slightly yellowish crystallizate

Solubility: Clear, colorless to slightly yellow solution in water (c=50 mg/ml)

Fructose-1,6- P_2 - $Na_3 \times 8 H_2O$ (from content enzymatically): $\geq 97\%$

Fructose-1,6- P_2 (enzymatically): $\geq 60\%$

Na (flame photometric): 11-15%

Water (K. Fischer): 23-29%

P_i : $\leq 0.6\%$

Heavy metals (as Pb): ≤ 10 ppm

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No. **Pack Size**

10 041 793 103 custom fill

Will be supplied as "Fructose-1,6-diphosphate, Trisodium Salt".
Unit of Measure is "kg".

For further processing only.

Glucose-1,6-diphosphate tetra(cyclohexylammonium) salt

Substrate in diagnostic tests

Application

Use Glucose-1,6-diphosphate in diagnostic tests for the determination of inorganic phosphate and sucrose.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 10139-18-1

Properties

Formula: $C_6H_{14}O_{12}P_2 \times (C_6H_{14}N)_4 \times 4 H_2O$

Molecular weight: 808.9 D (Glucose-1,6- P_2 : 340.1 D)

Specification

Appearance: Yellowish crystallizate

Glucose-1,6- P_2 (CHA) $_4 \times 4 H_2O$: 93.0-105.0%

Glucose-1,6- P_2 (from P_{organic}): 39.0-44.0%

CHA (titrimetric): 46.-50.0%

Water (K. Fischer): 6.0-10.0%

P_{organic} ($P_{\text{total}} - P_i$): 7.10-8.00%

P_i : $\leq 0.30\%$

Stability: At +15 to +25°C within specification range for 36 months

Cat. No. **Pack Size**

10 150 827 103 custom fill

Will be supplied as "Glucose-1,6-diphosphate, Tetra-CHA Salt".
Unit of Measure is "g".

For further processing only.

Glucose-6-phosphate disodium salt

Substrate for glucose-6-phosphate dehydrogenase

Application

Use Glucose-6-phosphate in diagnostic tests for the determination of glucose-6-phosphate dehydrogenase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 3671-99-6

Properties

Formula: $C_6H_{11}O_9PNa_2$

Molecular weight: 304.2 D (Glucose-6-P: 260.2 D)

Specification

Appearance: White lyophilizate

Solubility: Clear solution in water (c=50 mg/ml)

Glucose-6-P (enzymatically): $\geq 77\%$

Glucose-6-P (from P_{organic}): $\geq 77\%$

Na (flame photometric): $12.5 \pm 1\%$

Water (K. Fischer): $8 \pm 2\%$

P_{organic} ($P_{\text{total}} - P_i - P_{\text{fructose-6-P}}$): $\geq 8.9\%$

P_i : $\leq 0.6\%$

Fructose-6-P (enzymatically): $\leq 2\%$

Glucose (enzymatically): $\leq 0.2\%$

Stability: At +15 to +25°C within specification range for 24 months. Store dry.

Cat. No. **Pack Size**

10 153 079 103 custom fill

Will be supplied as "Glucose-6-phosphate, Disodium Salt". Unit of Measure is "g".

For further processing only.

N-Acetyl-L-Cysteine crystallizate

Activator of creatine kinase

Application

Use N-Acetyl-L-Cysteine in diagnostic tests for the determination of creatine kinase, where it is used to reactivate creatine kinase as recommended by the International Federation of Clinical Chemistry and Laboratory Medicine (IFCC).

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 616-91-1

Properties

Formula: $C_5H_9NO_3S$

Molecular weight: 163.19 D

Specification

Appearance: White crystals

Solubility: Clear, colorless solution in water (c=5%, w/v)

Melting range: +107 to +113°C

Specific rotation (c=2%, w/v, based on undried substance): $[\alpha]_{25/D}$: +3.0° to +5.0°

Heavy metals (as Pb): $\leq 5 \text{ ppm} \pm 0.0005\%$

Screening analysis:

Particle size $\geq 250 \mu\text{m}$: $\leq 15\%$

Particle size $\geq 100 \mu\text{m}$: $\leq 50\%$

Cat. No. **Pack Size**

10 068 365 103 custom fill

Will be supplied as "N-Acetyl-L-Cystein". Unit of Measure is "kg".

For further processing only.

Fe (AAS): ≤ 2.0 ppm \pm 0.0002%

Cu (AAS): ≤ 1.0 ppm \pm 0.0001%

Mn (AAS): ≤ 1.0 ppm \pm 0.0001%

Microbiological test: Corresponds to specification

IR-spectrum: Corresponds to reference

Purity (HPLC): ≥ 99.0 area%

N-Acetyl-L-Cysteine (Ellmann's reagent, based on undried substance): $\geq 99.0\%$

N-Acetyl-L-Cysteine (alkalimetric, based on undried substance): $\geq 99.0\%$

Content of nitrogen (elementary analysis, based on undried substance): $\geq 8.5\%$

Content from nitrogen (elementary analysis, based on undried substance): $\geq 99.0\%$

Stability: At +15 to +40°C within specification range for 24 months. Store dry in tightly closed containers.

Phosphoenolpyruvate (PEP), for potassium test tri(cyclohexylammonium) salt

Substrate for phosphoenolpyruvate carboxylase

Application

Use Phosphoenolpyruvate as a substrate for pyruvate kinase, stimulated by potassium, for the enzymatic determination of potassium.

Benefits

- Rely on the proven diagnostic quality of this product.
- Take advantage of the strongly reduced concentration of potassium.

CAS: 138-08-9

Properties

Formula: $C_3H_2O_6P(C_6H_{14}N)_3 \times H_2O$

Molecular weight: 483.3 D (PEP: 168.0 D)

Specification

Appearance: Colorless crystallizate

PEP-(CHA)₃ (from content found enzymatically): $\geq 96\%$

PEP (enzymatically): $\geq 34.5\%$

CHA (titrimetric): 57-67%

Water (K. Fischer): $\leq 4.5\%$

P_i: $\leq 0.6\%$

Pyruvate (enzymatically): $\leq 0.1\%$

Na (AES) : ≤ 100 ppm

K: ≤ 10 ppm

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No.	Pack Size
11 333 968 103	custom fill

Will be supplied as "PEP, tri-CHA for Potassium Test". Unit of Measure is "kg".

Additional formulation: Crystallized monosodium salt, Catalog No. 10 152 960 103

For further processing only.

Phosphoenolpyruvate (PEP) tri(cyclohexylammonium) salt

Substrate for phosphoenolpyruvate carboxylase

Application

Use Phosphoenolpyruvate in diagnostic tests for the determination of carbon dioxide, creatine or pyruvate kinase.

Cat. No.	Pack Size
10 005 185 103	custom fill

Will be supplied as "Phosphoenolpyruvate (PEP), CHA-Salt". Unit of Measure is "kg".

For further processing only.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 138-08-9

Properties

Formula: $C_3H_2O_6P(C_6H_{14}N)_3 \times H_2O$

Molecular weight: 483.3 D (PEP: 168.0 D)

Specification

Appearance: Colorless, crystalline powder

PEP salt (based on value found enzymatically): $\geq 96\%$

PEP (enzymatically): $\geq 34.5\%$

CHA (titrimetric with perchloric acid): 57-67%

Water (K. Fischer): $\leq 4.5\%$

P_i (Fiske and Subbarow): $\leq 0.6\%$

Pyruvate (enzymatically): $\leq 0.1\%$

Stability: At +2 to +8°C within specification range for 24 months.

Additional formulation crystallized monosodium salt, Catalog No. 10 152 960

Phosphoenolpyruvate (PEP) monosodium salt

Substrate for phosphoenolpyruvate carboxylase

Application

Use Phosphoenolpyruvate in diagnostic tests for the determination of carbon dioxide, creatine or pyruvate kinase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 138-08-9

Properties

Formula: $C_3H_4O_6PNa \times H_2O$

Molecular weight: 208.0 D (PEP: 168.0 D)

Specification

Appearance: White, crystalline powder

PEP-Na x H₂O (based on value found enzymatically): $\geq 94\%$

PEP (enzymatically): $\geq 76.0\%$

Na (flame photometric): 9-13%

Water (K. Fischer): 8-10%

P_i: $\leq 0.6\%$

Pyruvate (enzymatically): $\leq 0.1\%$

Stability: At +2 to +8°C within specification range for 24 months.

Pyruvate monosodium salt

Substrate for many enzymes, such as lactate dehydrogenase and pyruvate kinase.

Application

Use Pyruvate in diagnostic tests for the determination of lactate dehydrogenase.

Cat. No.	Pack Size
10 152 960 103	custom fill

Will be supplied as "Phosphoenolpyruvate (PEP), Mono-Na Salt".

Unit of Measure is "g".

For further processing only.

Cat. No.	Pack Size
10 005 525 103	custom fill

Will be supplied as "Pyruvate Monosodium Salt". Unit of Measure is "g".

For further processing only.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 57-60-3

Properties

Formula: $C_3H_3O_3Na$

Molecular weight: 110.0 D

Specification

Appearance: White, crystalline powder

Pyruvate-Na (from content found enzymatically): 96-103%

Pyruvate (enzymatically, based on anion): 77-81%

Na (flame photometric): 20.5-21.5%

Bioburden: ≤ 100 CFU/g

Heavy metals (as Pb): ≤ 10 ppm

Stability: At +15 to +25°C within specification range for 18 months.

S-Butyrylthiocholine Iodide crystallizate

Substrate for cholinesterase

Application

Use S-Butyrylthiocholine Iodide in diagnostic tests for the determination of cholinesterase.

Benefits

Rely on the proven diagnostic quality of this product.

CAS: 1866-16-6

Properties

Formula: $C_9H_{20}NOSJ$

Molecular weight: 317.2 D

Specification

Appearance: Colorless crystallizate

Solubility: Clear, colorless solution in water (c=70 mg/ml)

Melting range: +172 to +174°C

Butyrylthiocholine iodide (titrimetric) : $\geq 98.0\%$

Thiocholine iodide, free: $\leq 0.15\%$

Test for inhibitors of choline esterase; reaction rates: $100 \pm 5\%$

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No.

10 034 614 103

Pack Size

custom fill

Will be supplied as "S-Butyrylthiocholine Iodide". Unit of Measure is "kg".

For further processing only.





2 Immunology

Antibodies	160
Monoclonal Antibodies	160
Polyclonal Antibodies	176
Biotin/Streptavidin System	178
Streptavidin	178
Biotin Labels	179
Fluorescent Labels	181
Solid Phases	183
Dyes	190
Interference Eliminating Proteins (IEPs)	192
Specific Interference	192
Unspecific Interference	199
Marker Enzymes and Substrates	204
Enzymes	204
Substrates	212
Serums	216

MAB<CK-MB>M-7.4.5 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<CK-MB>M-7.4.5 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<CK-MB>M-6.12.47 IgG.

Product Description

Antibody class: IgG 1, kappa

Properties

MAB<CK-MB>M-7.4.5 IgG is a monoclonal antibody directed to creatine kinase MB. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

Cat. No.

11 719 815 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<CK-MB>M-7.4.5 IGG". Unit of Measure is "mg".

For further processing only.

MAB<CK-MM>Mix frozen solution

For measurement of human creatine kinase isoenzyme (CK-MB), which is a well established tool in confirming the diagnosis of acute myocardial infarction.

Application

MAB<CK-MM>Mix is a main ingredient of the CK-MB assay.

Benefits

- Inhibit the CK-M subunit with more than 99% efficiency.

Product Description

Immunogen: Human creatinine kinase isoenzyme MM (h-CK-MM)

Inhibitor capacity (for information only): Determined by Roche test CK liquid (+37°C):

25 µg (=2.1 µl) antibody solution/ml Hitachi reagent inhibits

4500 U/l CK-MM: ≥99.6%

300 U/l CK-BB: ± 5%

Properties

The MAB<CK-MM>Mix consists of four highly specific monoclonal mouse antibodies directed to Human creatinin kinase isoenzyme MM (h-CK-MM). The frozen solution contains protein, potassium phosphate buffer and sodium chloride. No preservative are added.

Specification

Appearance: Slightly opalescent colorless solution

Protein (BCA): 10.8-13.2

Cat. No.

04 688 457 103

Pack Size

1, 10, 50, 100, 1000 ml

Will be supplied as "Mab<CK-MM>Mix". Unit of Measure is "l".



For further processing only.

pH value (+25°C): 7.4-7.6
Purity (TSK 3000): ≥ 90 area%
HPLC (Mono Q basic material): ≥80%
Aggregates (HPLC / TSK3000): ≤10%
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 36 months. Avoid repeated freezing and thawing.

Background Information

Human creatine kinase isoenzyme CK-MB consist of two subunits: CK-M and CK-B. By inhibiting the CK-M subunit, the creatine kinase reaction is triggered exclusively by the β subunit of creatine kinase, which accounts for one-half of the activity of CK-MB.

Literature

Roche Applied Science, MABMix, May 2006

MAB<AFP>M-LJ738 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<AFP>M-LJ738 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<AFP>M-TU11 IgG.

Product Description

Antibody class: IgG 1, kappa

Properties

MAB<AFP>M-LJ738 IgG is a monoclonal antibody directed to alpha fetoprotein. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate
Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
Purity (HPLC / Mono Q): ≥90 area%
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

Cat. No.	Pack Size
11 492 101 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<AFP>M-LJ738-IGG(DE)". Unit of Measure is "g active ingredient".
 For further processing only.

MAB<AFP>M-TU11 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<AFP>M-TU11 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<AFP>M-LJ738 IgG.

Product Description

Antibody class: IgG 2a, kappa

Cat. No.	Pack Size
11 492 080 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<AFP>M-TU11-IGG(DE)". Unit of Measure is "g active ingredient".
 For further processing only.

Properties

MAB<AFP>M-TU11 IgG is a monoclonal antibody directed to alpha fetoprotein. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

MAB<CEA>M-TU2 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<CEA>M-TU2 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<CEA>M-TU3 IgG.

Product Description

Immunogen: Carcinoembryonic antigen (CEA)

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, kappa

Properties

MAB<CEA>M-TU2 IgG is a monoclonal antibody directed to carcinoembryonic antigen. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride; pH 7.5. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear, to slightly opalescent solution in NaCl, 0.9% (c=5mg/ml)

Purity (HPLC / Mono Q): ≥90 area% of total protein

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<CEA>M-TU3 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<CEA>M-TU3 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<CEA>M-TU2 IgG.

Cat. No.	Pack Size
11 353 713 103	5, 10, 100 mg

Will be supplied as "MAK<CEA>M-TU2-IGG *SQ". Unit of Measure is "mg active ingredient".
For further processing only.

Cat. No.	Pack Size
10 777 498 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<CEA>M-TU3-IGG(DE)". Unit of Measure is "mg active ingredient".
For further processing only.

Product Description

Antibody class: IgG 1, kappa

Properties

MAB<CEA>M-TU3 IgG is a monoclonal antibody directed to carcinoembryonic antigen. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate
Solubility: Clear, to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
Purity (HPLC / Mono Q): ≥90 area% of total protein
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<CK-MB>M-6.12.47 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<CK-MB>M-6.12.47 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<CK-MB>M-7.4.5 IgG.

Product Description

Antibody class: IgG 1, kappa

Properties

MAB<CK-MB>M-6.12.47 IgG is a monoclonal antibody directed to creatine kinase MB. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate
Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
Purity (HPLC / Mono Q): ≥90 area%
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

MAB<DD>M-1.2.57 IgG lyophilizate

Qualified for the Cobas® Core Modular Platform.

Application

MAB<DD>M-1.2.57 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<DD>M-2.1.16 IgG.

Cat. No.	Pack Size
11 719 823 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<CK-MB>M-6.12.47 IGG". Unit of Measure is "mg".
 For further processing only.

Cat. No.	Pack Size
12 156 903 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<DD>M-1.2.57-IgG(SP/Q)". Unit of Measure is "g active ingredient".
 For further processing only.

Product Description**Immunogen:** Human fibrinogen cleavage product D-Dimer**Spleen donor:** Mouse Balb/c**Antibody class:** IgG 1, kappa**Properties**

MAB<DD>M-1.2.57 IgG is a monoclonal antibody directed to human fibrinogen cleavage product D-Dimer. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification**Appearance:** White lyophilizate**Solubility:** Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)**Protein (Biuret) :** ≥0.6 mg/mg lyophilizate**Aggregated IgG (HPLC / TSK3000):** ≤10 area%**pH 5.5 treatment (30 minutes):** Corresponds to specification**Stability:** At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<DD>M-2.1.16 IgG lyophilizate

Qualified for the Cobas® Core Modular Platform.

Application

MAB<DD>M-2.1.16 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<DD>M-1.2.57 IgG.

Product Description**Immunogen:** Human fibrinogen cleavage product D-Dimer**Spleen donor:** Mouse Balb/c**Antibody class:** IgG 1, kappa**Properties**

MAB<DD>M-2.1.16 IgG is a monoclonal antibody directed to human fibrinogen cleavage product D-Dimer. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification**Appearance:** White lyophilizate**Solubility:** Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)**Protein (Biuret):** ≥0.6 mg/mg lyophilizate**Aggregated IgG (HPLC / TSK3000):** ≤10 area%**pH 5.5 treatment (30 minutes):** Corresponds to specification**Stability:** At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.**Cat. No.****12 045 206 103****Pack Size**

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<DD>M-2.1.16-IgG(SP/Q)". Unit of Measure is "g active ingredient".
For further processing only.

MAB<Ferr>M-3.170 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<Ferr>M-3.170 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<Ferr>M-4.184 IgG.

Product Description

Immunogen: Human liver ferritin

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, kappa

Cross reactivity to: Spleen ferritin 74%; heart ferritin 11%

Properties

MAB<Ferr>M-3.170 IgG is a monoclonal antibody directed to human liver ferritin. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Cat. No.

11 547 089 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<Ferr>M-3.170-IgG". Unit of Measure is "g active ingredient".

For further processing only.

MAB<Ferr>M-4.184 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<Ferr>M-3.170 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<Ferr>M-3.170 IgG.

Product Description

Immunogen: Human liver ferritin

Spleen donor: Mouse balb/c

Antibody class: IgG 2a, kappa

Cross reactivity to: Spleen ferritin 60%; heart ferritin 11%

Properties

MAB<Ferr>M-4.184 IgG is a monoclonal antibody directed to human liver ferritin. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Cat. No.

11 547 119 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<Ferr>M-4.184-IgG". Unit of Measure is "g active ingredient".

For further processing only.

MAB<FSH>M-1.303 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<FSH>M-1.303 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<FSH>M-W3 IgG.

Product Description

Immunogen: Human follicle stimulating hormone (FSH)

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, kappa

Cross reactivity to: Human chorionic gonadotropin (HCG) <0.1%; Luteinizing hormone (LH) <1.5%; Thyroid stimulating hormone (TSH) <2.5%

Properties

MAB<FSH>M-1.303-IgG is a monoclonal antibody directed to human follicle stimulating hormone. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Cat. No.

11 493 540 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<FSH>M-1.303-IGG". Unit of Measure is "mg active ingredient".

For further processing only.

MAB<FSH>M-W3g IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<FSH>M-W3 IgG highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<FSH>M-1.303 IgG.

Product Description

Immunogen: Human follicle stimulating hormone (FSH)

Spleen donor: Mouse Balb/c

Antibody class: IgG 2a, kappa

Cross reactivity to: Human chorionic gonadotropin (HCG) not detectable; Luteinizing hormone (LH) <1.0%; Thyroid stimulating hormone (TSH) <3.5%

Properties

MAB<FSH>M-W3 IgG is a monoclonal antibody directed to human follicle stimulating hormone. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Cat. No.

11 493 531 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<FSH>M-W3-IGG". Unit of Measure is "mg active ingredient".

For further processing only.

Specification

Appearance: White lyophilizate
Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<HCG>M-INN2 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<HCG>M-INN2 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<HCG>M-INN2.

Product Description

Immunogen: Human chorionic gonadotropin b-chain (HCG)
Spleen donor: Mouse Balb/c
Antibody class: IgG 1, kappa
Cross reactivity to: Luteinizing hormone (LH) <0,3%; Follicle stimulating hormone (FSH) <0.1%; Thyroid stimulating hormone (TSH) <0.1%

Properties

MAB<HCG>M-INN2 IgG is a monoclonal antibody directed to human chorionic gonadotropin β-chain. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate
Solubility: Clear, to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<HCG>M-INN22 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<HCG>M-INN22 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAK<HCG>M-INN2.

Product Description

Immunogen: Human chorionic gonadotropin β-chain (HCG)
Spleen donor: Mouse Balb/c
Antibody class: IgG 1, kappa
Cross reactivity to: Luteinizing hormone (LH) <5.0%; Follicle stimulating hormone (FSH) <0.2%; Thyroid stimulating hormone (TSH) <0.5%

Cat. No.	Pack Size
03 116 263 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<HCG>M-INN2-IgG". Unit of Measure is "mg active ingredient".
 For further processing only.

Cat. No.	Pack Size
11 812 564 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<HCG>M-INN22-IgG". Unit of Measure is "mg active ingredient".



For further processing only.

Properties

MAB<HCG>M-INN22 IgG is a monoclonal antibody directed to human chorionic gonadotropin β -chain. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear, to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥ 0.6 mg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<H-S-Amy>M-Tu66C7 IgG lyophilizate

For measurement of acute pancreatitis (pancreatic α -amylase) in human serum and urine selective blocking of salivary α -amylase isoenzyme is achieved in the presence of the pancreatic h- α -amylase.

Application

The combination of MAB<H-S-Amy>Tu88E8 and MAB<H-S-Amy>Tu66C7 inhibits the human salivary α -amylase $\geq 97\%$ while maintaining the activity of the pancreatic h- α -amylase.

Benefits

Profit from the highly selective blocking of salivary α -amylase isoenzyme.

Product Description

Immunogen: Human salivary amylase

Spleen donor: Mouse balb/c

Antibody class: IgG 1, kappa

Clone: Tu66C7

Properties

MAB<H-S-Amy>M-Tu66C7 IgG fraction is purified by chromatography and lyophilized from a solution containing protein (≥ 20 mg/ml), potassium-phosphate buffer and NaCl. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Reconstitute with 0.9% saline solution (c=10 mg/ml)

Protein (Biuret): ≥ 0.7 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥ 90 area% IgG

Cross reactivity to h-pancreas α -amylase: ≤ 1 U/gW

Function testing (synergetic effects at +37°C):

h salivary amylase + MAB <S-AMY>: $\leq 3\%$ amylase activity

h pancreas amylase + MAB <S-AMY>: $\geq 98\%$ amylase activity

Stability: At -15 to -25°C within specification range for 36 months. Avoid repeated freezing and thawing.

Cat. No.

11 543 601 103

Pack Size

custom fill

Will be supplied as "MAK<H-S-Amy>M-Tu66C7-IgG(BR)SQ". Unit of Measure is "g active ingredient".

For further processing only.

MAB<H-S-Amy>M-Tu88E8 IgG lyophilizate

For measurement of acute pancreatitis (pancreatic α -amylase) in human serum and urine selective blocking of salivary α -amylase isoenzyme is achieved in the presence of the pancreatic h- α -amylase.

Application

The combination of MAB<H-S-Amy>Tu88E8 and MAB<H-S-Amy>Tu66C7 inhibits the human salivary α -amylase $\geq 97\%$ while maintaining the activity of the pancreatic h- α -amylase.

Benefits

Profit from the highly selective blocking of salivary α -amylase isoenzyme.

Product Description

Immunogen: Human salivary amylase

Spleen donor: Mouse balb/c

Antibody class: IgG 2a, kappa

Clone: Tu88E8

Properties

MAB<H-S-Amy>M-Tu88E8 IgG fraction is purified by chromatography and lyophilized from a solution containing protein (≥ 20 mg/ml), potassium-phosphate buffer and NaCl. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear solution in NaCl, 0.9% (c=10 mg/ml)

Protein (Biuret): ≥ 0.7 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥ 90 area% IgG

Cross reactivity to h-pancreas α -amylase: ≤ 1 U/gW

Function testing (synergetic effects at $+37^\circ\text{C}$):

h salivary amylase + MAB : $\leq 3\%$ amylase activity

h pancreas amylase + MAB : $\geq 98\%$ amylase activity

Stability: At -15 to -25°C within specification range for 36 months. Avoid repeated freezing and thawing.

Literature

1) M. Gerber, K. Naujoks, H. Lenz, W. Gerhard, K. Wulff, Clin. Chem. **31**, 1331 (1985)

2) M. Gerber, K. Naujoks, H. Lenz, K. Wulff, Clin. Chem. **33**, 1158 (1987)

MAB<IGE>M-323 IgG lyophilizate

Qualified for the Cobas[®] Core / Elecsys[®] Modular Platform.

Application

MAB<IGE>M-323 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<IGE>M-7H8.

Product Description

Immunogen: Human IgE

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, light chain kappa

Cross reactivity to: Human IgM, IgG and IgA: Not detectable

Cat. No.	Pack Size
11 543 598 103	custom fill

Will be supplied as "MAK<H-S-Amy>M-Tu88E8-IgG(BR)SQ". Unit of Measure is "g active ingredient".

For further processing only.

Cat. No.	Pack Size
11 543 393 103	5 mg (samples), ≥ 50 mg (custom fill)

Will be supplied as "MAK<IGE>M-323-IgG". Unit of Measure is "g active ingredient".

For further processing only.

Properties

MAB<IGE>M-323 IgG is a monoclonal antibody directed to human IgE. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<IGE>M-7H8 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<IGE>M-7H8 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<IGE>M-323.

Product Description

Immunogen: Human IgE

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, kappa

Cross reactivity to: Human IgM, IgG and IgA not detectable

Properties

MAB<IGE>M-7H8 IgG is a monoclonal antibody directed to human IgE. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

MAB<INSULIN>M-BM1 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<INSULIN>M-BM1 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<INSULIN>M-ST3.

Product Description

Antibody class: IgG 1, kappa

Cat. No.

11 988 204 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<IGE>M-7H8-IGG". Unit of Measure is "g active ingredient".

For further processing only.

Cat. No.

12 208 725 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<INSULIN>M-BM1-IgG". Unit of Measure is "mg active ingredient".

For further processing only.

Properties

MAB<INSULIN>M-BM1 IgG is a monoclonal antibody to insulin. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate
Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
Purity (HPLC / Mono Q): ≥90 area%
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

MAB<INSULIN>M-ST3 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<INSULIN>M-ST3 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<INSULIN>M-BM1.

Product Description

Antibody class: IgG 1, kappa
Cross reactivity to: Human pro insulin 1.2%; porcine insulin (strong recognition, determined by radio immuno assay(RIA)); bovine insulin (weak recognition, determined by RIA)

Properties

MAB<INSULIN>M-ST3 IgG is a monoclonal antibody directed to insulin. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate
Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)
Protein (Biuret): ≥0.6 mg/mg lyophilizate
Purity (HPLC / Mono Q): ≥90 area%
pH 5.5 treatment (30 minutes): Corresponds to specification
Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

MAB<LH>M-11412 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<LH>M-11412 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<LH>M-2.406-IgG.

Cat. No.	Pack Size
12 208 750 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<INSULIN>M-ST3-IgG". Unit of Measure is "mg active ingredient".
 For further processing only.

Cat. No.	Pack Size
11 547 925 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<LH>M-11412-IgG". Unit of Measure is "g active ingredient".
 For further processing only.

Product Description**Immunogen:** Luteinizing hormone (LH)**Antibody class:** IgG 1, kappa**Cross reactivity to:** Human chorionic gonadotropin <0.1%**Properties**

MAB<LH>M-11412 IgG is a monoclonal antibody directed to luteinizing hormone. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification**Appearance:** White lyophilizate**Solubility:** Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)**Protein** (Biuret): ≥0.6 mg/mg lyophilizate**Purity** (HPLC / Mono Q): ≥90 area%**pH 5.5 treatment** (30 minutes): Corresponds to specification**Stability:** At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.**MAB<LH>M-2.406 IgG lyophilizate**

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<LH>M-2.406 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<LH>M-11412 IgG.

Product Description**Immunogen:** Luteinizing hormone (LH)**Antibody class:** IgG 1, kappa**Cross reactivity to:** Human chorionic gonadotropin (HCG) not detectable; Follicle stimulating hormone (FSH) <0.3%; Thyroid stimulating hormone (TSH) <0.3%**Properties**

MAB<LH>M-2.406-IgG is a monoclonal antibody directed to luteinizing hormone. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification**Appearance:** White lyophilizate**Solubility:** Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)**Protein** (Biuret): ≥0.6 mg/mg lyophilizate**Purity** (HPLC / Mono Q): ≥90 area%**pH 5.5 treatment** (30 minutes): Corresponds to specification**Stability:** At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.**MAB<PRL>M-C4E4 IgG lyophilizate**

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<PRL>M-C4E4 IgG is highly qualified for heterogeneous immunoassays.

Cat. No.**11 547 038 103****Pack Size**5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<LH>M-2.406-IgG". Unit of Measure is "g active ingredient".

For further processing only.

Cat. No.**11 458 701 103****Pack Size**5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<PRL>M-C4E4-IGG". Unit of Measure is "g active ingredient".

For further processing only.

Benefits

Rely on optimal results obtained with sandwich partner MAB<PRL>M-H12G10.

Product Description

Antibody class: IgG 1, kappa

Cross reactivity to: Human chorionic gonadotropin (HCG) not detectable; Luteinizing hormone (LH) not detectable; Thyroid stimulating hormone (TSH) not detectable.

Properties

MAB<PRL>M-C4E4 IgG is a monoclonal antibody directed to prolactin. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

MAB<PRL>M-H12G10 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<PRL>M-H12G10 IgG is highly qualified for heterogeneous Immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<PRL>M-C4E4 IgG.

Product Description

Antibody class: IgG 2a, kappa

Cross reactivity to: Human chorionic gonadotropin (HCG) not detectable; Luteinizing hormone (LH) not detectable; Thyroid stimulating hormone (TSH) not detectable.

Properties

MAB<PRL>M-H12G10 IgG is a monoclonal antibody directed to prolactin. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 60 months. Avoid repeated freezing and thawing.

Cat. No.	Pack Size
11 027 689 103	5 mg (samples), ≥50 mg (custom fill)

Will be supplied as "MAK<PRL>M-H12G10-IGG". Unit of Measure is "mg active ingredient".
For further processing only.

MAB<TSH>M-A8 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<TSH>M-A8 IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<TSH>M-Tu1.20.

Product Description

Immunogen: Human thyroid stimulating hormone

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, kappa

Cross reactivity to: Luteinizing hormone (LH) <0.1%; Follicle stimulating hormone (FSH) <0.1%

Properties

MAB<TSH>M-A8 IgG is a monoclonal antibody directed to thyroid stimulating hormone. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

Purity (HPLC / Mono Q): ≥90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Cat. No.

11 367 978 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<TSH>M-A8-IGG(BR)". Unit of Measure is "mg active ingredient".

For further processing only.

MAB<TSH>M-TU1.20 IgG lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

MAB<TSH>M-Tu1.20 IgG is highly qualified for heterogenous immunoassays.

Benefits

Rely on optimal results obtained with sandwich partner MAB<TSH>M-A8.

Product Description

Immunogen: Human thyroid stimulating hormone

Spleen donor: Mouse Balb/c

Antibody class: IgG 1, kappa

Cross reactivity to: Luteinizing hormone (LH) <0.2%; Follicle stimulating hormone (FSH) <3.0%.

Properties

MAB<TSH>M-Tu1.20 IgG is a monoclonal antibody directed to thyroid stimulating hormone. It is lyophilized from a solution containing protein, potassium phosphate buffer and sodium chloride. No preservatives are added.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c=10mg/ml)

Protein (Biuret): ≥0.6 mg/mg lyophilizate

Cat. No.

10 767 778 103

Pack Size

5 mg (samples),
≥50 mg (custom fill)

Will be supplied as "MAK<TSH>M-TU1.20-IGG(DE),LYO.". Unit of Measure is "mg active ingredient".

For further processing only.

Purity (HPLC / Mono Q): ≥ 90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

PAB<CRP>S IgG

frozen solution

For measurement of c-reactive protein (CRP).

Application

Use PAB<CRP>S IgG for turbidimetric / nephelometric assays. It can be coupled to latex surfaces and used in the respective assays.

Benefits

Detect aggregates of PAB<CRP>S IgG with CRP by turbidimetric measurements (λ 340/700 nm).

Product Description

Immunogen: Human C-reactive protein

Properties

The polyclonal antibody IgG directed to c-reactive protein is produced in sheep. It is prepared as solution containing protein (≥ 40 g/l); Tris buffer with NaN_3 , 0.09% (w/v).

Recommended working concentration: 15 mg/ml

Remark: When stored over longer periods at $+4^\circ\text{C}$, a slight turbidity may occur which can easily be removed by centrifugation. No alteration of antibody properties occurs thereby.

Specification

Appearance: Clear to slightly opalescent yellowish solution

pH value ($+25^\circ\text{C}$): 7.8-8.2

Protein (Biuret): 50-60 mg/ml

Purity (HPLC): ≥ 90 area%

Bioburden: ≤ 1000 CFU/ml

Function: Calibration curve characterization defined by turbidimetric measurement A_{340}/A_{700}

δA (Standard 1mg/dl): ≥ 0.025

δA (Standard 10mg/dl): ≥ 0.230

δA (Standard 25mg/dl): ≥ 0.480

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Background Information

C-reactive protein (CRP) is part of the β -globulin family found in human plasma. Increased levels of CRP are involved in a variety of inflammatory diseases. Furthermore elevated CRP serum levels indicate tissue injury, transplant rejection, carcinogenesis and acute myocardial infarction.

Cat. No. **Pack Size**

11 888 714 103 custom fill

Will be supplied as "PAK<CRP>S-IgG *SQ". Unit of Measure is "l".



For further processing only.

PAB<T3>S IgG

lyophilizate

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

PAB<T3>S IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on the optimal design for competitive assay formats.

Product Description

PAB<T3>S IgG is a polyclonal antibody directed to triiodothyronine produced in sheep.

Cat. No. **Pack Size**

10 907 332 103 10 mg (samples),
 ≥ 50 mg (custom fill)

Will be supplied as "PAK<T3>S-IGG(DE),(ES 3G)". Unit of Measure is "g active ingredient".

For further processing only.

Immunogen: Trijodothyronine derivative.

Specification

Appearance: White lyophilizate

Protein (A_{280}): ≥ 0.7 mg/mg lyophilizate

HPLC (HPLC / TSK 3000): ≥ 90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

**PAB<T4>S IgG
 lyophilizate**

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

PAB<T4>S IgG is highly qualified for heterogeneous immunoassays.

Benefits

Rely on the optimal design for competitive assay formats.

Product Description

PAB<T4>S IgG is a polyclonal antibody directed to trijodothyronine produced in sheep.

Immunogen: Trijodothyronine derivative.

Specification

Appearance: White lyophilizate

Protein (A_{280}): ≥ 0.7 mg/mg lyophilizate

HPLC (HPLC / TSK 3000): ≥ 90 area%

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
10 767 794 103	10 mg (samples), ≥ 50 mg (custom fill)

Will be supplied as "PAK<T4>S-IGG(DE),(ES)". Unit of Measure is "mg active ingredient".
 For further processing only.

Streptavidin, recombinant from *Streptomyces avidinii*, expressed in *E. coli*, lyophilizate

Biotin/Streptavidin Portfolio

Application

Use Streptavidin, recombinant as a tool for solid phase technology and universal detection systems in immunology and molecular diagnostics.

Benefits

- Rely on a complete solution using highly qualified reagents in combination with activated Biotin Esters.

Specification

Appearance: White lyophilizate

Protein (A_{282} ; factor 3.1): 0.6-0.8 mg/mg lyophilizate

Specific activity/Biotin binding capacity: ≥ 17 U/mg protein

Proteasen (incubation with Azocoll for up to 24 hours at +25°C): ≤ 0.001 U/mg

Absorption (A_{405} against repurified water): ≤ 0.01

Water (K. Fischer): $\leq 12\%$

IEF (pH 6-9): Two main bands between 6.8 and 7.5

SDS-PAGE: Chromatographically homogeneous

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Background Information

Streptavidin consists of four subunits with a molecular weight of 13 kD, each containing a single biotinbinding site. Each subunit has six thyrrosine residues. The protein is carbohydrate free.

Streptavidin + 4 biotin \rightarrow streptavidin \cdot (biotin)₄

The formation of the complex is measured at 233 nm.

Cat. No. **Pack Size**

11 520 679 103 custom fill

Will be supplied as "Streptavidin Special Quality". Unit of Measure is "g active ingredient".

For further processing only.

D-Biotin-N-hydroxysuccinimide ester crystalline powder

Biotin/Streptavidin Portfolio

Application

Use D-Biotin-N-hydroxysuccinimide ester as biotinylating reagent for proteins and aminolabeled oligonucleotides.

Benefits

- Rely on a complete solution using highly qualified reagents in combination with Streptavidin.

CAS: 35013-72-0

Properties

Nomenclature: D-Biotinyl-N-hydroxy-succinimide ester

Formula: C₁₄H₁₉N₃O₅S

Molecular weight: 341.4 D

Remark: Under mild conditions the activated ester reacts with amino groups. The aminocaproic acid spacer is useful if biotinylated macromolecules are coupled because steric hindrance is minimized.

Specification

Appearance: White crystallize

Biotin ester (from N): 98-103%

N (elementary analysis): 12.0-12.73%

Purity (TLC: silica gel, 1-butanol/glacial acetic acid/H₂O= 2/1/1; a) in UV; b) with KMnO₄): Chromatographically homogeneous

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No.	Pack Size
10 734 250 103	custom fill

Will be supplied as "D-Biotin-N-hydroxy-succinimide Ester". Unit of Measure is "g".

For further processing only.

D-Biotinoyl- ϵ -aminocaproic acid-N-hydroxysuccinimide ester powder

Biotin/Streptavidin Portfolio

Application

Use D-Biotinoyl- ϵ -aminocaproic acid-N-hydroxysuccinimide ester as biotiny-lating reagent for proteins and aminolabeled oligonucleotides.

Benefits

- Rely on a complete solution using highly qualified reagents in combination with Streptavidin.

CAS: 72040-63-2

Properties

Nomenclature: D-Biotinyl- ϵ -amido caproic acid N-hydroxysuccinimid ester

Formula: $C_{20}H_{30}N_4O_6S$

Molecular weight: 454.5 D

Remark: Under mild conditions the activated ester reacts with amino groups. The aminocaproic acid spacer is useful if biotinylated macromolecules are coupled because steric hindrance is minimized.

Specification

Appearance: White to beige powder

Biotin ester (from N): $\geq 97\%$

N (elementary analysis): 11.9-12.8%

C (elementary analysis): 51.0-54.7%

H (elementary analysis): 6.4 to 6.9%

Purity (TLC: silica gel, 1-butanol/glacial acetic acid/ H_2O = 50/15/25, iodide stream/UV): Chromatographically homogeneous

Hydrolysis product (NMR): $\leq 20\%$

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No.

Pack Size

11 003 933 103 custom fill

Will be supplied as "D-Biotinyl- ϵ -aminocap. AcidN-Hydroxy Succ.". Unit of Measure is "g".

For further processing only.

Streptavidin R-Phycoerythrin LumiGrade Reagent

Ready to use solution

Standard for highly sensitive fluorescent detection.

Application

Conjugated reporter dyes such as Streptavidin R-Phycoerythrin (SA-PE) are well established for Luminex's xMAP Assay Kits or array-based applications due to their excellent spectral characteristics. Their dedicated and reproducible design ensure high performance multiplex assays using antibodies, receptors, peptides and oligonucleotides.

Benefits

- Improve detection and quantification due to high signal-to-noise combined with very low background .
- Rely on premium production quality that results in:
 - excellent lot-to-lot consistency
 - high purity
 - ready to-use solution
 - reagent stability up to 24 months at +2 to +8°C

Product Description

The Streptavidin R-Phycoerythrin LumiGrade conjugate contains Phycoerythrin from red algae (RPE) and Streptavidin (SA), recombinant.

Properties

High molecular weight conjugate size distribution (700-1200 kD).

Specification

Appearance: Reddish solution in potassium phosphate, 50 mmol/l, pH 6.8, sodium azide, 0.05% (c=1 mg/ml)

Streptavidin (A_{280}): 16-21 weight%

Absorption ratio $A_{566/280}$: >3.3

Protein: 1.00 ± 0.10 mg/ml

Content of color relating to SA: 0.95-1.40 (molar ratio RPE : SA)

Purity (HPLC/ TSK 3000 XL): free SA <1%

Fluorescence emission: Maximum/intensity

488 nm excitation: 576 nm ± 5 / ≥200

545 nm excitation : 576 nm ± 5 / ≥250

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.

05 065 925 103

Pack Size

1, 5, 100 ml, custom fill

Will be supplied as "Streptavidin R-Phycoerythrin Lumi Grade Reagent". Unit of Measure is "mg active ingredient".

For further processing only.

Streptavidin R-Phycoerythrin LumiGrade Ultrasensitive Reagent

Ready to use solution

Standard for ultrasensitive fluorescent detection.

Application

Conjugated reporter dyes such as Streptavidin R-Phycoerythrin (SA-PE) are well established for Luminex's xMAP[®] Assay Kits and array-based applications due to their excellent spectral characteristics. Their dedicated and reproducible design ensure high performance multiplex assays using antibodies, receptors, peptides and oligonucleotides.

Benefits

- Improve detection and quantification due to high signal-to-noise combined with very low background .
- Rely on the superior production quality that results in:
 - excellent lot-to-lot consistency
 - high purity
 - ready to-use solution
 - reagent stability 12 months at +2 to +8°C

Product Description

The Streptavidin R-Phycoerythrin LumiGrade conjugate contains Phycoerythrin (PE) from red algae and Streptavidin (SA), recombinant.

Properties

High molecular weight conjugate size distribution (1,500-50,000 kD).

Specification

Appearance: Reddish solution in potassium phosphate, 50 mmol/l, pH 6.8, sodium azide, 0.05% (c=1 mg/ml)

A₅₆₆: 0.73-0.81

SA-R-PE from R-Phycoerythrin (A₅₆₆/7.7): 0.95-1.05 mg/ml


Purity (HPLC / TSK 6000): ≥ 99.7 area%

Contamination (HPLC / TSK 6000): ≤ 0.3 area%

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
05 351 693 103	1, 100 ml, custom fill

Will be supplied as "Streptavidin R-Phycoerythrin Lumi Grade ultrasensitive". Unit of Measure is "mg active ingredient".

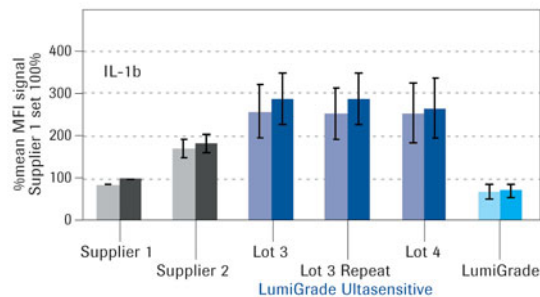

Excellent lot-to-lot consistency: Streptavidin R-Phycoerythrin LumiGrade Ultrasensitive conjugates of lot 3 and 4 show very consistent readouts with high signal-to-noise in the Interleukin 1b xMAP[®] Assay. Streptavidin R-PE of Supplier 1 serves as reference, normalized to 100% for each concentration. Mean values of ± Standard Deviation excluding blank values are shown.

For further processing only.

Multi Analyte Stripe universal device

Universal test stripe

Application

Use the Multi Analyte Stripes for a range of applications, such as classical sandwich immunoassays to NA/Oligo-detection.

Benefits

- Receive qualitative results in less than 10 minutes.

Product Description

Test stripes are delivered in boxes, each containing 50 stripes.

Specification

Components of the test stripe (4.6 nm):

MAB<Dig>IgG on gold conjugate: 0.2 µg per test stripe

Poly-Streptavidin (result line): 0.8 µg per test stripe

PAB<MouseFc>IgG (control line): 0.1 µg per test stripe

Sensitivity (analytical):

A biotin/-digoxigenin-peptide in 700µl buffer solution at a concentration of 50 pg/ml is visually detected as positive on the basis of the result line in the read out zone after chromatography.

Specificity (analytical):

The test stripes don't show a visible result line after chromatography of 700 µl buffer without adding biotin/-digoxigenin-peptide (negative control).

Stability: At +2 to +8°C from date of manufacturing for 18 months.

Remark: The Multi Analyte Stripe is an immunoassay test stripe employing anti-biotin/anti-digoxigenin.

Cat. No.

05 354 358 103

Pack Size

50 stripes in 1 box

Will be supplied as "Multi Analyte Stripe". Unit of Measure is "piece".

Minimum order size: 85 000 single strips.

For further processing only.

Streptavidin Magnetic Particles suspension

Streptavidin-coated magnetic particles

Application

Use Streptavidin Magnetbeads for the fast and simple separation of a variety of biotin-labeled molecules in solution. This includes single mRNA, oligo-nucleotides, DNA, DNA fragments, glycoconjugates and protein isolation. Streptavidin Magnetbeads can also be used as a solid phase within liquid ELISA systems.

Benefits

- Use the highly suitable streptavidin-coated magnetic particles for sensitive isolation and detection systems.

Properties

The beads show very high lot-to-lot consistency and are already included in standard applications. They are offered in a protease- and DNase/RNase-free format.

Specification

Appearance: Brown suspension

Content: 9-11 mg/ml (solid binding)

Performance test in mRNA-HS-kit: Function corresponds

Specific activity/Biotin binding capacity: ≥ 1800 pmol/mg

Proteases (incubation for up to 30 minutes at +37°C, casein-resorufin-marked): Not detectable

Cat. No.

11 636 502 103

Pack Size

custom fill

Will be supplied as "Streptavidin Magnetic-Particles". Unit of Measure is "g active ingredient".

For further processing only.

RNases (incubation for up to 4 hours at +37°C): Not detectable
DNases (incubation for up to 4 hours at +37°C): Not detectable
Unspecific binding of protein: Not detectable
Stability: At +2 to +8°C within specification range for 24 months.

StreptaWell, 384 plate transparent, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8 –10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C1 (standard bind)
SA-coated area: ≥90 µl/well
Blocked volume: >90 µl/well
Coating variance (CV): ≤8%
Total biotin binding capacity (competition assays): ≥1.5 ng/well
Homogeneity [VK] of series: ≤15%
Bleeding: <5 ngSA/well
Stability: At +2 to +8°C within specification range for 36 months.

Cat. No.	Pack Size
11 974 998 103	custom package

Will be supplied as "TRSA-SA MTP 384-well, clear". Unit of Measure is "piece".
 For further processing only.

StreptaWell, 384 plate white, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8 –10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C1 (standard bind)
SA-coated area: ≥90 µl/well
Blocked volume: >90 µl/well

Cat. No.	Pack Size
11 974 980 103	custom package

Will be supplied as "TRSA-SA MTP 384-well, white". Unit of Measure is "piece".
 For further processing only.

Coating variance (CV): ≤8%

Total biotin binding capacity (competition assays): ≥1.5 ng/well

Homogeneity [VK] of series: ≤15%

Bleeding: <5 ngSA/well

Stability: At +2 to +8°C within specification range for 36 months.

StreptaWell, C1, breakapart transparent, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8–10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C1 (standard bind)

SA-coated area: ≥250 μl/well

Blocked volume: >250 μl/well

Coating variance (CV): <5%

Total biotin binding capacity (competition assays): ≥5 ng/well

Homogeneity [VK] of series: ≤10%

Bleeding: >5 ngSA/well

Stability: At +2 to +8°C within specification range for 48 months.

Cat. No.	Pack Size
03 246 507 103	custom package

Will be supplied as "SA-MTP (N-breakap. transp./C1)". Unit of Measure is "piece".

For further processing only.

StreptaWell, C8 module, high binding capacity white, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8–10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C2 (high bind)

Cat. No.	Pack Size
11 975 021 103	custom package

Will be supplied as "TRSA-Bi/SACP) C8 plate, white". Unit of Measure is "piece".

For further processing only.

SA-coated area: $\geq 300 \mu\text{l}/\text{well}$
Blocked volume: $> 300 \mu\text{l}/\text{well}$
Coating variance (CV): $< 5\%$
Total biotin binding capacity (competition assays): $\geq 25 \text{ ng}/\text{well}$
Homogeneity [VK] of series: $\leq 10\%$
Stability: At +2 to +8°C within specification range for 36 months.

StreptaWell, C8, breakapart, high binding capacity

transparent, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8–10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C2 (high bind)
SA-coated area: $\geq 250 \mu\text{l}/\text{well}$
Blocked volume: $> 250 \mu\text{l}/\text{well}$
Coating variance (CV): $< 5\%$
Total biotin binding capacity (competition assays): $\geq 20 \text{ ng}/\text{well}$
Homogeneity [VK] of series: $\leq 10\%$
Bleeding: $> 5 \text{ ngSA}/\text{well}$
Stability: At +2 to +8°C within specification range for 36 months.

Cat. No.	Pack Size
11 986 694 103	custom package

Will be supplied as "SA-MTP (N-breakap.C8/C2 plus)". Unit of Measure is "piece".
 For further processing only.

StreptaWell, C8, lockwell

transparent, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8 –10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Cat. No.	Pack Size
04 869 532 103	custom package

Will be supplied as "SA-MTP (NUNC LOCKWELL C8 TRANSP/C1)". Unit of Measure is "piece".
 For further processing only.

Specification**Type of coating:** C1 (standard bind)**SA-coated area:** ≥ 250 $\mu\text{l}/\text{well}$ **Blocked volume:** > 250 $\mu\text{l}/\text{well}$ **Coating variance (CV):** $< 5\%$ **Total biotin binding capacity** (competition assays): ≥ 5 ng/well**Homogeneity [VK] of series:** $\leq 10\%$ **Bleeding:** > 2 ngSA/well**Stability:** At $+2$ to $+8^\circ\text{C}$ or at $+15$ to $+25^\circ\text{C}$ within specification range for 36 months.**StreptaWell, C96 plate
white, coated with recombinant streptavidin**

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid assay amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8 –10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification**Type of coating:** C1 (standard bind)**SA-coated area:** ≥ 300 $\mu\text{l}/\text{well}$ **Blocked volume:** > 300 $\mu\text{l}/\text{well}$ **Coating variance (CV):** $< 5\%$ **Total biotin binding capacity** (competition assays): ≥ 5 ng/well**Homogeneity [VK] of series:** $\leq 10\%$ **Stability:** At $+2$ to $+8^\circ\text{C}$ within specification range for 36 months.**StreptaWell, C96 plate, high binding
capacity
transparent, coated with recombinant streptavidin**

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid assay amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8 –10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Cat. No.	Pack Size
11 975 005 103	custom package

Will be supplied as "TRSA-SA MTP C96 Plate, white". Unit of Measure is "piece".

For further processing only.

Cat. No.	Pack Size
11 975 030 103	custom package

Will be supplied as "TRSA-Bi/SACP) C96 plate, clear". Unit of Measure is "piece".

For further processing only.

Specification

Type of coating: C2 (high bind)
SA-coated area: ≥300 µl/well
Blocked volume: >300 µl/well
Coating variance (CV): <5%
Total biotin binding capacity (competition assays): ≥25 ng/well
Homogeneity [VK] of series: ≤10%
Stability: At +2 to +8°C within specification range for 36 months.

StreptaWell, C96 plate, high binding capacity white, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8 –10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C2 (high bind)
SA-coated area: ≥300 µl/well
Blocked volume: >300 µl/well
Coating variance (CV): <5%
Total biotin binding capacity (competition assays): ≥25 ng/well
Homogeneity [VK] of series: ≤10%
Bleeding: >5 ngSA/well
Stability: At +2 to +8°C within specification range for 36 months.

Cat. No.	Pack Size
11 975 013 103	custom package

Will be supplied as "TRSA-Bi/SACP) C96 Plate, white". Unit of Measure is "piece".
 For further processing only.

StreptaWell, F8 module

transparent, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid assay amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8–10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C1 (standard bind)

SA-coated area: $\geq 300 \mu\text{l/well}$

Blocked volume: $> 300 \mu\text{l/well}$

Coating variance (CV): $< 5\%$

Total biotin binding capacity (competition assays): $\geq 5 \text{ ng/well}$

Homogeneity [VK] of series: $\leq 10\%$

Bleeding: $> 2 \text{ ngSA/well}$

Stability: At $+2$ to $+8^\circ\text{C}$ within specification range for 36 months.

Cat. No.	Pack Size
11 940 279 103	custom package

Will be supplied as "SA coated MTP Nunc F8". Unit of Measure is "piece".

For further processing only.

StreptaWell, F8 module, high binding capacity

transparent, coated with recombinant streptavidin

Streptavidin-coated microwell plates

Application

Use Streptavidin-coated microplates in all kinds of immunoassays as well as in assay systems for detecting protein-nucleic acid interactions and nucleic acid assay amplification and hybridization. Various types of transparent and white streptavidin-coated plates (96- and 384-well) are available for colorimetric, chemiluminescent and fluorescent applications. The minimum order size is 500 plates/order with a delivery time of 8–10 weeks.

Benefits

- Rely on Roche StreptaWell plates combining excellent signal to noise ratio with increased sensitivity and extremely low bleeding properties.
- Use Roche StreptaWell plates to achieve superb intra- and inter-assay precision.

Specification

Type of coating: C2 (high bind)

SA-coated area: $\geq 300 \mu\text{l/well}$

Blocked volume: $> 300 \mu\text{l/well}$

Coating variance (CV): $> 5\%$

Total biotin binding capacity (competition assays): $\geq 25 \text{ ng/well}$

Homogeneity [VK] of series: $\leq 10\%$

Bleeding: $< 5 \text{ ngSA/well}$

Stability: At $+2$ to $+8^\circ\text{C}$ within specification range for 36 months.

Cat. No.	Pack Size
11 965 875 103	custom package

Will be supplied as "SA-MTP (Nunc F8 transp./C2+)". Unit of Measure is "piece".

For further processing only.

Colloidal Gold 20 nm suspension

Colloidal Gold is well known as an established labeling tool for a broad range of blotting and diagnostic applications and for electron- /light microscopy.

Application

Use Colloidal Gold, 20 nm as conjugation partner for all kind of antibodies, proteins and macromolecules. It is recommended for low to medium sensitive assays.

Benefits

- Rely on the highly lot-to-lot consistent quality: Colloidal Gold 20 nm is manufactured in unique production lots of up to 250 l.

Product Description

The 20 nm Goldsol quality is red and spherical. The ready to use "gold suspension" is pH adjusted.

Properties

"Citrate Gold" obtained from reduction of Tetrachloro-auric acid (HAuCl_4) with citric acid.

Specification

Appearance: Clear, light red liquid

Particle size: 19-23 nm

Particle concentration (A_{520} unit): 0.85-1.00

λ_{max} : 516.0-518.5 nm

Stability: At +2 to +8°C within specification range for 12 months.

Background Information

Due to its intense red color Colloidal Gold is one of the basic components for test strip development and manufacturing.

Cat. No.	Pack Size
05 418 291 103	1, 5, 25 l

Will be supplied as "Colloidal Gold 20 nm". Unit of Measure is "l".
For further processing only.

Colloidal Gold 40 nm suspension

Colloidal Gold is well known as an established labeling tool for a broad range of blotting and diagnostic applications and for electron- /light microscopy.

Application

Use Colloidal Gold, 40 nm as conjugation partner for all kind of antibodies, proteins and macromolecules. It is recommended for high sensitive assays due to its unique shape and color.

Benefits

- Rely on the highly lot-to-lot consistent quality: Colloidal Gold 40 nm is manufactured in unique production lots of up to 250 l.

Product Description

The 40 nm Goldsol quality is red/violet ("potato shape"). The ready to use "gold suspension" is pH adjusted.

Properties

"Citrate Gold" obtained from reduction of Tetrachloro-auric acid (HAuCl_4) with citric acid.

Specification

Appearance: Turbid, slightly opalescent raspberry red liquid

Particle size: 38.0-43.0 nm

190 Particle concentration (A_{520} unit): 1 ± 0.2

Cat. No.	Pack Size
05 416 744 103	1, 5, 25 l

Will be supplied as "Colloidal Gold 40 nm". Unit of Measure is "l".
For further processing only.

λ_{max} : 531±1 nm

Stability: At +2 to +8°C within specification range for 12 months.

Background Information

Due to its intense red/violet color Colloidal Gold is one of the basic components for test strip development and manufacturing.

Interference Eliminating Proteins

Interference elimination is an important differentiation factor and challenging task for state-of-the-art immunoassay development and kit manufacturing. Test applications with increased sensitivity, better precision, improved test accuracy and optimized handling convenience demand powerful tools for efficient interference elimination.

Roche's product line of interference eliminating proteins (non-specific and specific for antibody interferences) can be combined to meet the specific test requirements.

Framework IEP lyophilizate

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Application

Use the monoclonal MAB33 Framework-IEP to eliminate monomeric (immuno-response) and highly specific interferences directed against the framework regions of antibodies. Its interference eliminating strength is based on identical immunoreactive amino acid sequences in the framework regions of the test antibody.

Benefits

- Eliminate exclusively interference to produce improved assay performance.
- Increase the assay sensitivity and accuracy.
- Profit from higher precision, greater accuracy and optimized handling convenience.

Product Description

Immunogen: PTH

Spleen donor: Mouse Balb/c

Antibody class: IgG1, kappa

Preparation: MAB33 Framework IEP is lyophilized from a solution containing protein, potassium phosphate and NaCl. No further preservatives are added.

Properties

Molecular structure: IgG1, monomer

Remark: Cannot be used in test systems for determination of PTH.

Serum concentrations of CK-MM: ≤ 3 U/ml have shown neither influence on interference elimination properties nor on recovery of analyte.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless to slightly opalescent solution in NaCl, 0.9% (c=10 mg/ml)

Purity (HPLC): $\geq 90\%$ IgG of total protein

Functional activity (relative titer based on masterlot determined by MTP assay): $\geq 80\%$

Recommended working concentration: 20-5,000 $\mu\text{g/ml}$ incubation buffer

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Literature

See **MAB33-IgG1**.

Cat. No.

03 369 846 103

Pack Size

5, 50, 250 mg

Will be supplied as "Framework IEP *SQ". Unit of Measure is "mg active ingredient".

For further processing only.

HAMA Serum, Type 1 lyophilizate

Cat. No.	Pack Size
11 767 275 103	custom fill

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Will be supplied as "Hama-Serum 1- Qual. Standard *SQ". Unit of Measure is "piece".
For further processing only.

Application

HAMA Serum Type 1 is primarily intended for polyvalent and spontaneous HAMA interference occurring in healthy donors.

Benefits

- Rely on a reproducible source for positive HAMA interference samples in immunoassays.
- Use HAMA Serum Type I as positive control for test development and quality control purposes.

Product Description

HAMA Serum Type 1 is a lyophilized human serum with serum ingredients within normal range. No preservative are added. The product must be handled just as carefully as patient specimens.

Properties

pH value (+25°C): 7-8

Protein (Biuret): ≥65 mg/vial

Analyte concentrations determined by Roche Elecsys® -Tests:

Alpha fetoprotein (AFP): 2.0 IU/ml

Carcinoembryonic antigen (CEA): 1.3 ng/ml

Follicle stimulating hormone (FSH): 10.0 mIU/ml

Human chorionic gonadotropin (HCG): 2.8 mIU/ml

Luteinizing hormone (LH): 8.1 mIU/ml

Prolactin: 160.7 µIU/ml

Prostate specific antigen (PSA): 0.3 ng/ml

PSA free: 0.06 ng/ml

Thyroid stimulating hormone (TSH): 2.5 µIU/ml

Troponin T (TN-T): Not detectable

Specification

Appearance: Yellowish lyophilizate

Interference effect: Corresponds to specification

Infectious parameters (determined by FDA approved methods):

HbsAg: negative

Anti HIV 1+2: Negative

Anti HCV: Negative

HIV 1 Ag: Negative

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Background Information

HAMA serum interferences in immunoassays can vary within a broad range depending on the person's immune system.

Literature

Roche Applied Science, Interference-Eliminating Proteins for the Diagnostics Industry, 7 ed., July 2009

HAMA Serum, Type 2 lyophilizate

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Application

HAMA Serum Type 2 primarily represents mono-/bivalent and specific HAMA interference occurring after treatment with monoclonal antibodies.

Benefits

- Rely on a reproducible source for HAMA interference material in immunoassays.
- Use HAMA Serum Type 2 as positive control development and quality control purposes.

Product Description

HAMA Serum Type 2 is a lyophilized human serum with serum ingredients within normal range. No preservative are added. The product must be handled just as carefully as patient specimens.

Properties

pH value (+25°C): 7-8

Protein (Biuret): ≥65 mg/vial

Analyte concentrations determined by Roche Elecsys® -Tests:

Alpha fetoprotein (AFP): 3.5 IU/ml

Carcinoembryonic antigen (CEA): 1.2 ng/ml

Follicle stimulating hormone (FSH): 7.0 mIU/ml

Human chorionic gonadotropin (HCG): 3.0 mIU/ml

Luteinizing hormone (LH): 6.7 mIU/ml

Prolactin: 114 µIU/ml

Prostate specific antigen (PSA): 0.3 ng/ml

PSA free: 0.05 ng/ml

Thyroid stimulating hormone (TSH): 1.6 µIU/ml

Troponin T (TN-T): 0.01 ng/ml

Specification

Appearance: Yellowish lyophilizate

Interference effect: Corresponds to specification

Infectious parameters:

HbsAg: Negative

Anti HIV 1+2: Negative

Anti HCV: Negative

HIV 1 Ag: Negative

Stability: At -15 to -25°C within specification range for 24 months.

Background Information

HAMA serum interferences in immunoassays can vary within a broad range depending on the person's immune system.

Literature

Roche Applied Science, Interference-Eliminating Proteins for the Diagnostics Industry, 7 ed., July 2009

Cat. No. **Pack Size**

05 167 060 103 custom fill

Will be supplied as "HAMA Serum 2L *SQ". Unit of Measure is "piece".

For further processing only.

MAB IgG2b/Fab2a Poly lyophilizate

Cat. No.	Pack Size
11 355 830 103	5, 50, 250 mg

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Application

Use the polymer MAB IgG2b/Fab2a Poly to reduce polymeric interference against IgG2a and/or IgG2b antibodies. It also covers Fab neo-epitopes.

Benefits

- Eliminate exclusively interference to produce improved assay performance.
- Increase the assay sensitivity and accuracy.
- Profit from higher precision, greater accuracy and optimized handling convenience.

Product Description

Immunogen: DPH (2b), human S-AMY (2a)

Spleen donor: Mouse Balb/c

Antibody class: IgG2b, kappa/IgG2a, kappa

Preparation: Lyophilized from a solution containing potassium phosphate and NaCl and 6% sucrose. No further preservatives are added.

Properties

Molecular structure: IgG2a-Fab, polymerized with IgG2b; defined molecular range distribution

Remarks: Cannot be used in test systems for determination of DPH, human S-AMY.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless to slightly opalescent solution in NaCl, 0.9% (c=5 mg/ml)

Functional activity (relative titer based on master lot determined by MTP assay): ≥80%

Recommended working concentration: 0.5-500 µg/ml incubation buffer

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Literature

- 1) E. Mössner, H. Lenz, G. Bienhaus, Poster: AACCC/IFCC (1990)
- 2) E. Mössner, H. Lenz, G. Bienhaus, Clin. Chem. **36**, 1093 (1990)
- 3) R. Valdes, J. Clin. Imm. **15/2**, 87 (1992)
- 4) H. Vaidya, Clin. Chem. **38/9**, 1737 (1992)
- 5) St. Avrameas, Mol. Imm. **30/12**, 1133 (1993)
- 6) E. Wilkinson, Clin. Chem. **39/10**, 2166 (1993)
- 7) L.J. Kricka, Clin. Chem of Acta **215**, 153 (1993)
- 8) H.J. Hansen, J. Clin. Imm. **16/4**, 294 (1993)
- 9) M. Kuroki, J. Imm. Meth. **180**, 81 (1995)
- 10) R. Sapin, Clin. Chem. **41/1**, 117 (1995)
- 11) H. Schlebusch et al., Hybridoma **14**, 167, 74 (1995)
- 12) P. Mikrosch, Eur. J. Chem. Biochem. **35/11**, 881 (1997)
- 13) U. Hasholzner, Anticancer Research, **17**, 3055 (1997)
- 14) St. Levison, J. Clin. Ligand Assay **20/2**, 180 (1997)
- 15) L.J. Kricka, Clin. Chem. **45/7**, 942 (1999)
- 16) Roche Applied Science, Interference-Eliminating Proteins for the Diagnostics Industry, 7 ed., July 2009

Will be supplied as "MAB-IgG(2b)/Fab(2a) Polymer, PolyMAB2b/2a". Unit of Measure is "mg active ingredient". For further processing only.

MAB33 IgG1 lyophilizate

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Application

Use MAB33 IgG1 for test formulations employing intact IgG1. MAB33 IgG1 is especially suitable for the elimination of monomeric and specific interference.

Benefits

- Eliminate exclusively interference to produce improved assay performance.
- Increase the assay sensitivity and accuracy.
- Profit from higher precision, greater accuracy and optimized handling convenience.

Product Description

Immunogen: h CK-MM

Spleen donor: Mouse Balb/c

Antibody class: IgG1, kappa

Preparation: Lyophilized from a solution containing potassium phosphate and NaCl. No further preservatives are added.

Properties

Molecular structure: IgG1, monomer

Remarks: Cannot be used in test systems for determination of CK-MM and CK-MB.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless to slightly opalescent solution in NaCl, 0.9% (c=10 mg/ml)

Protein (Biuret): 0.7 mg protein/mg lyophilizate

Purity (HPLC / Mono Q): ≥95 area% IgG of total protein

Functional activity (relative titer based on master lot determined by MTP assay): ≥80%

Recommended working concentration: 50-5000 µg/ml incubation buffer

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Background Information

MAB33 IgG1 is a monoclonal antibody with defined specificity.

Literature

- 1) E. Mössner, H. Lenz, G. Bienhaus, Poster: AACC/IFCC (1990)
- 2) E. Mössner, H. Lenz, G. Bienhaus, Clin. Chem. **36**, 1093 (1990)
- 3) R. Valdes, J. Clin. Imm. **15/2**, 87 (1992)
- 4) H. Vaidya, Clin. Chem. **38/9**, 1737 (1992)
- 5) St. Avrameas, Mol. Imm. **30/12**, 1133 (1993)
- 6) E. Wilkinson, Clin. Chem. **39/10**, 2166 (1993)
- 7) L.J. Kricka, Clin. Chem of Acta **215**, 153 (1993)
- 8) H.J. Hansen, J. Clin. Imm. **16/4**, 294 (1993)
- 9) M. Kuroki, J. Imm. Meth. **180**, 81 (1995)
- 10) R. Sapin, Clin. Chem. **41/1**, 117 (1995)
- 11) H. Schlebusch et al., Hybridoma **14**, 167, 74 (1995)
- 12) P. Mikrosch, Eur. J. Chem. Biochem. **35/11**, 881 (1997)
- 13) U. Hasholzner, Anticancer Research, **17**, 3055 (1997)
- 14) St. Levison, J. Clin. Ligand Assay **20/2**, 180 (1997)
- 15) L.J. Kricka, Clin. Chem. **45/7**, 942 (1999)
- 16) Roche Applied Science, Interference-Eliminating Proteins for the Diagnostics Industry, 7 ed., July 2009

Cat. No.

Pack Size

11 200 941 103 custom fill

Will be supplied as "MABM-33-IgG(DE),SQ MAB 33". Unit of Measure is "g active ingredient".

For further processing only.

MAB33 IgG1/Fab1 Poly lyophilizate

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Application

Use the polymer MAB33 IgG1/Fab1 Poly for assays working with Fab-conjugates.

Benefits

- Eliminate exclusively interference to produce improved assay performance.
- Increase the assay sensitivity and accuracy.
- Profit from higher precision, greater accuracy and optimized handling convenience.

Product Description

Immunogen: h CK-MM

Spleen donor: Mouse Balb/c

Antibody class: IgG1

Preparation: Lyophilized from a solution containing potassium phosphate and NaCl and 6% sucrose. No further preservatives are added.

Properties

Molecular structure: IgG1-Fab, polymerized with IgG1; defined molecular range distribution

Remarks: Cannot be used in test systems for determination of CK-MM and CK-MB.

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless to slightly opalescent solution in NaCl, 0.9% (c=10 mg/ml)

Functional activity (relative titer based on master lot determined by MTP assay): ≥80% **Recommended working concentration:** 0.5-500 µg/ml incubation buffer

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Background Information

The polymer MAB33 IgG1/Fab1 Poly incorporates elements comparable with those of MAB33 IgG1/IgG1 Poly, and also covers interferences against Fab neo-epitopes.

Literature

- 1) E. Mössner, H. Lenz, G. Bienhaus, Poster: AACC/IFCC (1990)
- 2) E. Mössner, H. Lenz, G. Bienhaus, Clin. Chem. **36**, 1093 (1990)
- 3) R. Valdes, J. Clin. Imm. **15/2**, 87 (1992)
- 4) H. Vaidya, Clin. Chem. **38/9**, 1737 (1992)
- 5) St. Avrameas, Mol. Imm. **30/12**, 1133 (1993)
- 6) E. Wilkinson, Clin. Chem. **39/10**, 2166 (1993)
- 7) L.J. Kricka, Clin. Chem of Acta **215**, 153 (1993)
- 8) H.J. Hansen, J. Clin. Imm. **16/4**, 294 (1993)
- 9) M. Kuroki, J. Imm. Meth. **180**, 81 (1995)
- 10) R. Sapin, Clin. Chem. **41/1**, 117 (1995)
- 11) H. Schlebusch et al., Hybridoma **14**, 167, 74 (1995)
- 12) P. Mikrosch, Eur. J. Chem. Biochem. **35/11**, 881 (1997)
- 13) U. Hasholzner, Anticancer Research, **17**, 3055 (1997)
- 14) St. Levison, J. Clin. Ligand Assay **20/2**, 180 (1997)
- 15) L.J. Kricka, Clin. Chem. **45/7**, 942 (1999)
- 16) Roche Applied Science, Interference-Eliminating Proteins for the Diagnostics Industry, 7 ed., July 2009

Cat. No.	Pack Size
11 368 338 103	5, 50, 250 mg

Will be supplied as "MAB-IgG/Fab (Polymer), SQ Poly MAB 33".
Unit of Measure is "g active ingredient".
For further processing only.

MAB33 IgG1/IgG1 Poly frozen solution

Part of the modular MAB33 family toolbox for highly specific interference elimination.

Application

Use the polymer MAB33 IgG1/IgG1 Poly for formulations employing intact IgG1. MAB33 IgG1/IgG1 Poly is more efficient for polymeric and less specific types of interference (compared to monoclonal MAB 33 IgG1).

Benefits

- Eliminate exclusively interference to produce improved assay performance.
- Increase the assay sensitivity and accuracy.
- Profit from higher precision, greater accuracy and optimized handling convenience.

Product Description

Immunogen: h CK-MM

Spleen donor: Mouse Balb/c

Antibody class: IgG1

Preparation: Lyophilized from a solution containing potassium phosphate and NaCl and 6% sucrose. No further preservatives are added.

Properties

Molecular structure: Molecular structure: IgG1, polymerized with IgG1; defined molecular range distribution

Remarks: Cannot be used in test systems for determination of CK-MM and CK-MB.

Specification

Appearance: Frozen liquid

Solubility: Yellowish clear to slightly opalescent solution, containing K-phosphate, NaCl and 4% sucrose, pH 7,5

Protein (Biuret): ≥30 mg protein/mg lyophilizate

Purity (HPLC): ≥90 area% IgG of total protein

Functional activity (relative titer based on master lot determined by MTP assay): ≥80 %

Recommended working concentration: 0.5-500 µg/ml incubation buffer

pH 5.5 treatment (60 minutes): Corresponds to specification

Turbidity properties δA_{334} : Corresponds to specification

Bioburden: ≤250 CFU/ml

Stability: At -60 to -90°C within specification range for 24 months. Avoid repeated freezing and thawing.

Background Information

MAB33 IgG1/IgG1 Poly is the polymerized chemical version of MAB33 IgG1.

Literature

- 1) E. Mössner, H. Lenz, G. Bienhaus, Poster: AACC/IFCC (1990)
- 2) E. Mössner, H. Lenz, G. Bienhaus, Clin. Chem. **36**, 1093 (1990)
- 3) R. Valdes, J. Clin. Imm. **15/2**, 87 (1992)
- 4) H. Vaidya, Clin. Chem. **38/9**, 1737 (1992)
- 5) St. Avrameas, Mol. Imm. **30/12**, 1133 (1993)
- 6) E. Wilkinson, Clin. Chem. **39/10**, 2166 (1993)
- 7) L.J. Kricka, Clin. Chem of Acta **215**, 153 (1993)
- 8) H.J. Hansen, J. Clin. Imm. **16/4**, 294 (1993)
- 9) M. Kuroki, J. Imm. Meth. **180**, 81 (1995)
- 10) R. Sapin, Clin. Chem. **41/1**, 117 (1995)
- 11) H. Schlebusch et al., Hybridoma **14**, 167, 74 (1995)
- 12) P. Mikrosch, Eur. J. Chem. Biochem. **35/11**, 881 (1997)

Cat. No.

11 939 661 103

Pack Size

5, 50, 250, 1000 mg

Will be supplied as "MAB-33-IgG-Polymer *SQ". Unit of Measure is "g active ingredient".

For further processing only.

Bovine IgG (PAB<->R-IgG) lyophilizate

Part of the product portfolio for nonspecific interference elimination.

Application

Use Bovine IgG to reduce nonspecific adsorption of antibodies to the solid phase and other cross-reactive, nonspecific antibody interactions.

Benefits

- Block the nonspecific linkage of test components to the solid phase.
- Minimize background signals.
- Improve assay sensitivity and dynamics.

Product Description

The polyclonal antibody Bovine IgG is lyophilized from a solution containing potassium phosphate and sodium chloride. No preservatives are added. IgG fraction is produced in bovine, purified by anion-exchange chromatography.

Properties

Recommended working concentration: 0.5-2.5 mg/ml

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in NaCl, 0.9% (c=10 mg/ml)

Turbidimetric measurement (A_{546} , against water): ≤ 100 mE

Protein (Biuret): ≥ 0.8 mg/mg lyophilizate

Aggregated IgG (HPLC / TSK3000): $\leq 5\%$

Country of origin: USA or NZL

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No. **Pack Size**

11 293 621 103 custom fill

Will be supplied as "PAB<->R-IgG(DE) Bovine IgG". Unit of Measure is "g active ingredient".



For further processing only.

Bovine Serum Albumin I lyophilizate

Part of the product portfolio for nonspecific interference elimination.

Application

Use BPLA Type I to reduce non-specific adsorption to the solid phase or to saturate unoccupied binding sites. It is recommended for assays with higher demands on sensitivity (e.g., thyroid, tumor markers).

Benefits

- Rely on the premium quality for blocking the non-specific linkage of test components to the solid phase.
- Minimize background signals.
- Improve assay sensitivity and dynamics.
- Profit from a higher purity grade with low levels of steroids, haptens and lipids than BSA V.

Product Description

BPLA Type I is lyophilized from a solution containing potassium phosphate and sodium chloride. No preservatives are added.

Specification

Appearance: Yellowish lyophilizate

Solubility: Clear to slightly turbid yellowish solution in water (c=60 mg/ml)

A_{405} (against water): ≤ 0.250

pH value: 6.5-7.5

Protein (from N, according to Dumas, factor 6.25): $\geq 95\%$

Cat. No. **Pack Size**

11 726 536 103 custom fill

Will be supplied as "Albumin RPLA 1 Assay Quality". Unit of Measure is "kg".

Additional products: BPLA Type IV, Catalog No. 11726544103; BPLA Type IV, new, with minimized content of proteases, is available for sampling under Catalog No. 03535240103.

For further processing only.

Purity (HPLC / TSK 3000): ≥95% (monomer)

Water (K. Fischer): ≤5%

Bioburden: ≤50 CFU/g

Ca: ≤0.1%

Fe: ≤0.005%

Cu: ≤0.002%

Complexing agent:

Recovery of Fe: 100±20%

Recovery of Cu: 100 ± 20%

Heavy metals (as Pb): ≤0.002%

P_i: ≤0.005%

Octanoic acid (GC): ≤0.5%

Analysis of T3, T4, Estradiol, Testosterone, Progesterone (for information only): Values stated on certificate of analysis.

Country of origin: USA

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Bovine Serum Albumin IV lyophilizate

Part of the product portfolio for nonspecific interference elimination.

Application

Use BPLA Type IV to reduce non-specific adsorption to the solid phase or to saturate unoccupied binding sites. BPLA Type IV represents a very high purity grade (≥ 95%) and is recommended for assays with higher demands on sensitivity (*e.g.*, thyroid, tumor markers).

Benefits

- Rely on the premium quality for blocking the non-specific linkage of test components to the solid phase.
- Minimize background signals.
- Improve assay sensitivity and dynamics.

Product Description

BPLA Type IV is lyophilized from a solution containing potassium phosphate and sodium chloride. No preservatives are added.

Specification

Appearance: Yellowish lyophilizate

Solubility: Clear to slightly turbid yellowish solution in water (c=60 mg/ml)

A₄₀₅ (against water): ≤0.200

pH value: 6.5-7.5

Protein (from N, according to Dumas, factor 6.25): ≥95%

Purity (HPLC / TSK 3000): ≥95% (monomer)

Water (K. Fischer): ≤5%

Bioburden: ≤50 CFU/g

Ca: ≤0.1%

Fe: ≤0.005%

Cu: ≤0.002%

pH 4.5 treatment (up to 3 hours): Corresponds to specification

Country of origin: USA, NZL

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No. **Pack Size**

11 726 544 103 custom fill

Will be supplied as "Albumin RPLA 4 Assay Quality". Unit of Measure is "kg".

For further processing only.

Poly BSA Type I frozen solution

Part of the product portfolio for nonspecific interference elimination.

Application

Use Polymeric BSA Type I for the elimination of nonspecific interferences.

Benefits

- Rely on the premium quality for effective elimination of multivalent interferences (*e.g.*, polyvalent antibodies/conjugates or antigens) in infectious disease assays and tumor marker tests.

Product Description

Poly BSA Type I is produced from bovine serum albumine (BSA) by polymerization. Its surface is chemically modified by acetylation, and the negative charge of the Poly BSA Type I eliminates hydrophobic interactions. The Poly BSA Type I solution contains protein, potassium-phosphate buffer and preservatives (chloroacetamide, methylisothiazolone and sucrose).

Properties

Recommended working concentration: 0.1-20 mg/ml incubation buffer.

Specification

Appearance: Yellowish, clear to slightly opalescent solution (frozen)

pH value (+25°C): 6.8-7.2

Protein (Biuret): ≥40 mg/ml

Particle size (Photon correlations spectrometry): 15-45 nm

Country of origin: USA, NZL

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Cat. No.	Pack Size
11 866 737 103	1, 5, 20 g

Will be supplied as "Poly BSA Type I *SQ". Unit of Measure is "g active ingredient".



For further processing only.

Poly BSA Type II frozen solution

Part of the product portfolio for nonspecific interference elimination.

Application

Use Polymeric BSA Type II for the elimination of nonspecific interferences.

Benefits

- Rely on the premium quality for effective elimination of multivalent interferences (*e.g.*, polyvalent antibodies/conjugates or antigens) in infectious disease assays and tumor marker tests.

Product Description

Poly BSA Type II is produced from bovine serum albumine (BSA) by polymerization. Its surface is chemically modified by succinylation, and the negative charge of the Poly BSA Type II eliminates hydrophobic interactions. The Poly BSA Type II solution contains protein, potassium-phosphate buffer, and preservatives (chloroacetamide, methylisothiazolone and sucrose).

Properties

Recommended working concentration: 0.1-20 mg/ml incubation buffer

Specification

Appearance: Yellowish, clear to slightly opalescent solution (frozen)

pH value (+25°C): 6.8-7.2

Protein (Biuret): ≥40 mg/ml

Cat. No.	Pack Size
11 816 438 103	1, 5, 20 g

Will be supplied as "poly BSA Type II *SQ". Unit of Measure is "g active ingredient".



For further processing only.

Interference Eliminating Proteins (IEPs)

Unspecific Interference

Particle size (Photon correlations spectrometry): 25-55 nm

Country of origin: USA, NZL

pH 5.5 treatment (30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months. Avoid repeated freezing and thawing.

Rabbit IgG (PAB<->K-IgG) lyophilizate

Part of the product portfolio for nonspecific interference elimination.

Application

Use Rabbit IgG to reduce the nonspecific antibody interference in assays using rabbit antibodies.

Benefits

- Block the non specific linkage of test components to the solid phase.
- Minimize background signals.
- Improve assay sensitivity and dynamics.

Product Description

The polyclonal antibody Rabbit IgG is lyophilized from a solution containing potassium phosphate and sodium chloride. No preservative are added. IgG fraction is produced in rabbit, purified by anion-exchange chromatography.

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl 0.9% (c=10mg/ml)

Protein (Biuret) : ≥0,7 mg/mg lyophilizate

Purity (HPLC / TSK 3000): Corresponds to specification (in comparison to master lot)

Stability: At -15 to -25°C within specification range for 36 months. Avoid repeated freezing and thawing.

Cat. No. **Pack Size**

10 912 280 103 custom fill

Will be supplied as "PAB<->K-IgG(DE-FF)". Unit of Measure is "g active ingredient".

For further processing only.

Sheep IgG (PAB<->S-IgG) lyophilizate

Part of the product portfolio for nonspecific interference elimination.

Application

Use Sheep IgG to reduce the nonspecific antibody interference in assays using sheep antibodies.

Benefits

- Block the nonspecific linkage of test components to the solid phase.
- Minimize background signals.
- Improve assay sensitivity and dynamics.

Product Description

The polyclonal antibody Sheep IgG is lyophilized from a solution containing potassium phosphate and sodium chloride. No preservatives are added. IgG fraction is produced in sheep, purified by anion-exchange chromatography.

Properties

Recommended working concentration: 0.5-1.5 mg/ml

Specification

Appearance: White lyophilizate

Solubility: Clear to slightly opalescent solution in NaCl, 0.9% (c= 5 mg/ml)

Cat. No. **Pack Size**

10 717 606 103 custom fill

Will be supplied as "PAB<->S-IgG". Unit of Measure is "g".

For further processing only.

Turbidimetric measurement (A_{546} , against water): ≤ 100 mE Protein (Biuret):

0.8 mg/mg lyophilizate

Aggregated IgG (HPLC / TSK3000): $\leq 10\%$

Stability: At -15 to -25°C within specification range for 24 months.

Alkaline Phosphatase Mutein, recombinant from calf intestine, expressed in *Pichia pastoris*, lyophilizate

Part of marker enzyme portfolio

Application

Use Alkaline Phosphatase Mutein (AP Mutein) to eliminate human serum derived AP directed assay interferences.

Benefits

- Rely on the special design for reducing alkaline phosphatase-related assay interference interactions.

Product Description

In addition to the AP Mutein, recombinant the lyophilizate contains a mixture of proteins from the *Pichia* system supporting the interference elimination. AP Mutein is lyophilized from a solution containing in NaCl, 0.2 mol/l; ZnCl₂, 0.1 mmol/l; Tea, 30 mmol/l; MgCl₂, 1 mmol/l; raffinose, 50% (w/v); pH approximately 7.6. Production is done according to the procedures of the active enzyme.

EC 3.1.3.1

Specification

Appearance: White to yellowish lyophilizate

Solubility: Clear to light yellowish solution in 50% glycerol solution (c=10 mg/ml)

pH value: 7.0-8.0

Protein (A₂₈₀, 1 mg/ml=1, against water): ≥0.2 mg protein/mg lyophilizate

Specific activity (+37°C, pNPP): ≤10 U/mg protein

SDS-gel (qualitative comparison of the gel bands in reference to the bands of a standard): Corresponds to specification

Stability: At +2 to +8°C within specification range for 24 months.

Background Information

Alkaline Phosphatase Mutein (AP Mutein), recombinant, is the inactive form of recombinant highly active AP, expressed in *Pichia pastoris*. The inactivation of AP Mutein is based on one single point mutation located in the active site of the alkaline phosphatase.

Cat. No. **Pack Size**

04 781 007 103 custom fill

Will be supplied as "AP-Mutein, rec.". Unit of Measure is "g".
For further processing only.

Alkaline Phosphatase, EIA Grade from calf intestine, solution

Part of marker enzyme portfolio

Application

Marker enzyme for preparing antibody- / antigen-enzyme conjugates incorporated in immunoassay reagents for colorimetric, fluorimetric and luminometric detection.

Benefits

- Synthesize stable and reproducible Alkaline Phosphatase (AP) conjugates with AP, sourced from NZL intestines, that serves as a reliable origin.

EC 3.1.3.1

Properties

Nomenclature: Orthophosphoric-monoester phosphohydrolase (alkaline

Cat. No. **Pack Size**

10 556 602 103 custom fill

Will be supplied as "Phosphatase, Alkaline, Calf Intestine". Unit of Measure is "g".
For further processing only.

optimum)

Molecular weight: ≥57 kD

Inhibitors: P_i, metal chelating agents, divalent heavy metal ions (e.g., Be²⁺, Zn²⁺), many amino acids (e.g., L-phenylalanine, L-tryptophan, L-cysteine), iodo-sobenzoate, iodoacet-amide.

Activators: Mg²⁺, Co²⁺, Mn²⁺

pH optimum: 9.8

pH stability: 8.0

Thermal stability: Up to +40°C

Specification

Appearance: Clear, colourless solution in NaCl, 3 mol/l; MgCl₂, 1 mmol/l;

ZnCl₂, 0.1 mmol/l; Tea, 30 mmol/l

pH value: 7.0-8.0

Protein (A₂₈₀, 1 mg/ml=1, against water): ≥10 mg/ml

Specific activity (+37°C, pNPP): ≥3000 U/mg

Alkaline Phosphatase (HPLC): ≥95 area% (HPLC profile added to certificate)

Amino groups: 8-16 mol/mol

Carbohydrates, n=2: No limit

pH 5.5 treatment (30 minutes): Corresponds to specification

Origin of bovine intestine: NZL

Stability: At +2 to +8°C within specification range for 15 months.

Background Information

Alkaline phosphatase catalyzes the hydrolysis of numerous phosphate esters, such as esters of primary and secondary alcohols, sugar alcohols, cyclic alcohols, phenols and amines. Phosphodiesterases do not react with Alkaline Phosphatase, EIA Grade. The enzyme hydrolyzes PP_i. The kinetic properties of the enzyme depend on many factors, such as purity of enzyme, concentration of enzyme in the assay, buffer, pH etc.

Literature

- 1) A.P. Schaap, H. Akhavan, L.J. Romano, Clin. Chem. **35**, No. 9 1863-1864 (1989)
- 2) M.S. Urdea, J. Kolberg, J. Clyne, J.A. Running, D. Besemer, B. Warner, R. Sanchez-Pescader, Clin. Chem. **35**, No. 8, 1571-1575 (1989)
- 3) E. Jablonski, E.W. Moomaw, R.H. Tullis, J.L. Ruth, Nucleic Acids Res. **14**, 6115-28 (1986)
- 4) D.G. Williams, J. of Immunological Methods **72**, 261-268 (1984)
- 5) Roche Applied Science, Alkaline Phosphatase, 4. ed., 2007.

Alkaline Phosphatase, recombinant, highly active

from calf intestine, expressed in *Pichia pastoris*, solution

Part of marker enzyme portfolio

Application

Marker enzyme for preparation of antibody- / antigen-enzyme conjugates incorporated in highly sensitive immunoassay reagents for colorimetric, fluorimetric and luminometric detection. Alkaline Phosphatase is recommended for conjugation via carbohydrate groups (content approximately 30%).

Benefits

- Rely on excellent superior product quality of Alkaline Phosphatase (AP) recombinant, highly active.
- Synthesize stable, highly active and reproducible AP conjugates.
- Eliminate the risk of BSE contamination: No animal-derived components are used in the production process.

EC 3.1.3.1

Cat. No.	Pack Size
03 137 031 103	custom fill

Will be supplied as "AP, Yeast, high act., rec., EIA, NaCl". Unit of Measure is "g".

For further processing only.

Properties

IEP (IEF, CE): 3.6-4.7

MALDI-TOF MS:

Total molecular weight: 124±0 kD

Molecular weight protein: 104 kD (=84%)

Molecular weight carbohydrate: 20±10 kD (16±6%)

Accessible N-glycosylation sites: 2/subunit

Branching type: Higher branched type (hybrid type) GlcNAc, Mannose, no NeuAc detected

O-glycosylation sites: Not detected

Number of isoenzymes:

present: 3 (MS)

based on protein: 1 (MS)

Specification

Appearance: Clear, colorless solution in NaCl, 3 mol/l; MgCl₂, 5 mmol/l; ZnCl₂, 0.1 mmol/l; Tea, 30 mmol/l, pH approximately 7.6

pH value: 7.0-8.0

Protein (A₂₈₀, 1 mg/ml=1, against water): ≥20 mg/ml

Specific activity (+37°C, pNPP): ≥7000 U/mg

Alkaline Phosphatase (HPLC): ≥95% (HPLC profile added to certificate)

Amino groups: 5-13 mol/mol

Carbohydrates, n=2: No limit

Stability: At +2 to +8°C within specification range for 12 months.

Background Information

Alkaline Phosphatase recombinant, highly active catalyzes the hydrolysis of numerous phosphate esters, such as esters of primary and secondary alcohols, sugar alcohols, cyclic alcohols, phenols and amines. Phosphodiester do not react with Alkaline Phosphatase, recombinant. The enzyme hydrolyzes PP_i. The kinetic properties of the enzyme depend on many factors, such as purity of enzyme, concentration of enzyme in the assay, buffer, pH etc.

Literature

- 1) Th. Manes, M.F. Hoylaerts, R. Mueller, F. Lottspeich, W. Hoelke, J.L. Millán, Genetic Complexity, Structure and Characterisation of Highly Active Bovine Intestinal Alkaline Phosphatases, JBC, **273**, 36, 23353-23360 (1998)
- 2) R.K. Bretthauer, F.J. Castellino, Glycosylation of Pichia Pastoris derived Proteins, Biotechnol. Appl. Biochem. **30**, 193-200 (1999)
- 3) I. Ceveghino, I. Gregg, Heterologous Protein Expression in the Methylotropic Yeast Pichia Pastoris, FEMS Microbiology Reviews **24**, 45-66 (2000)
- 4) G. Gellissen, Heterologous Protein Production in Methylotropic Yeasts, Appl. Microbiol. Biotechnol. **54**, 741- 50 (2000)
- 5) Roche Applied Science, Alkaline Phosphatase, 4. ed., 2007.

Alkaline Phosphatase, recombinant, highly active, carbohydrate reduced from calf intestine, expressed in *Pichia pastoris*, solution

Part of marker enzyme portfolio

Application

Marker enzyme for preparation of antibody- / antigen- enzyme conjugates incorporated in highly sensitive immunoassay reagents for colorimetric, fluorimetric and luminometric detection. Recommended for conjugation via amino groups.

Benefits

- Rely on excellent product quality of Alkaline Phosphatase (AP) recombinant, highly active.

Cat. No. **Pack Size**

03 535 452 103 custom fill

Will be supplied as "AP, highly active, recombinant, CR". Unit of Measure is "g".

For further processing only.

- Synthesize stable, highly reactive, reproducible AP conjugates.
- Eliminate the risk of BSE contamination: No animal-derived components are used in the production process.

EC 3.1.3.1

Properties**IEP** (IEF, CE): 5.2-6.0**MALDI-TOF MS:**

Total molecular weight: 111 kD

Molecular weight protein: 104 kD (= 94%)

Molecular weight carbohydrates: 6.5 kD (= 6%)

Accessible N-glycosylation sites: 2/subunit**Branching type:** Reduced branched type GlcNAc, Mannose no NeuAc detected**O-glycosylation sites:** Not detected**Number of isoenzymes:**

present: 1 (MS)

based on protein: 1 (MS)

Specification**Appearance:** Clear, colourless solution in NaCl, 3 mol/l; MgCl₂, 5 mmol/l; ZnCl₂, 0.1 mmol/l; Tea, 30 mmol/l, pH approximately 7.6**pH value:** 7.0-8.0**Protein** (A₂₈₀; 1 mg/ml=1; against water): ≥20±1 mg/ml**Specific activity** (+37°C, pNPP): ≥ 7,000 U/mg**Alkaline Phosphatase** (HPLC): ≥95 area% (HPLC profile added to certificate)**Amino groups:** 5-13 mol/mol**Carbohydrates**, n=2: ≤7%**Stability:** At +2 to +8°C within specification range for 12 months.**Background Information**

AP recombinant, highly active CR catalyzes the hydrolysis of numerous phosphate esters, such as esters of primary and secondary alcohols, sugar alcohols, cyclic alcohols, phenols and amines. Phosphodiesterases do not react with Alkaline Phosphatase, recombinant. The enzyme hydrolyzes PP_i. The kinetic properties of the enzyme depend on many factors, such as purity of enzyme, concentration of enzyme in the assay, buffer, pH and others. The product contains a significantly reduced carbohydrate moiety.

β-Galactosidase from *E. coli*, lyophilizate

Part of marker enzyme portfolio

Application

Marker enzyme for the manufacturing of antibody- and antigen-enzyme conjugates incorporated in immunoassays for colorimetric and fluorimetric detection.

Benefits

- Use β-Galactosidase as a reliable source for hydrolysis of β-Gal conjugates.

EC 3.2.1.23

Properties**Nomenclature:** β-D-galactohydrolase**Molecular weight** (by sequence): 465 kD**Structure:** 4 identical subunits; β-galactosidase contains no carbohydrates

Cat. No.	Pack Size
11 291 963 103	custom fill

Will be supplied as "β-Galactosidase, Lyo.". Unit of Measure is "MU".

For further processing only.

Specification**Appearance:** White lyophilizate**pH value** (c=10 mg/ml, in water): 7.0-8.0**Protein** (Biuret): 0.25-0.5 mg/mg lyophilizate**Activity** (+37°C, 2-NP-β-D-galactoside): ≥120 U/mg lyophilizate**Specific activity** (+37°C, 2-nitrophenyl-β-D-galactopyranoside): ≥300 U/mg protein**Contaminants** (expressed as percentage of β-Galactosidase activity):

β-Fruktosidase: ≤0.001

α-Galactosidase: ≤0.001

Glucose-DH: <0.001

α-Glucosidase: ≤0.001

"NADH oxidase": ≤0.001

Na (flame photometric): ≤2500 ppm**Stability:** At +2 to +8°C: within specification range for 12 months. Store dry.**Background Information**

β-Galactosidase hydrolyzes β-D-galactosides.

β-Galactosidase Mutein
from *E. coli* overproducer, lyophilizate

Part of marker enzyme portfolio

Application

Use β-Galactosidase Mutein to eliminate β-galactosidase directed interferences in immunoassays derived from human sera.

Benefits

- **Eliminate the interference activity of β-galactosidase directed antibodies:** β-Galactose Mutein reacts with β-galactosidase directed antibodies in the same way as the active β-galactosidase.

EC 3.2.1.23

Properties

β-Galactosidase Mutein is identical to native β-galactosidase with respect to immuno-reactivity, conjugation, properties, surface charge, hydrophobicity, molecular weight, production procedure and down stream processing.

Specification**Appearance:** White lyophilizate**Solubility:** Clear, colorless to slightly opalescent solution in water (c=20 mg/ml)**Protein** (Biuret): 0.15-0.30 mg/mg lyophilizate**Specific activity** (10 mg/ml, +37°C, 2-nitrophenyl-β-D-galactopyranoside): ≤0.2 U/mg protein**Aggregated β-Galactose Mutein** (HPLC / TSK 4000): ≤10%**Immunoreactivity** (based on ML): 80-120%**Stability:** At -15 to -25°C within specification range for 24 months.**Background Information**β-Galactose Mutein from *E.coli* is constructed using site-directed mutagenesis of single amino acids in the active site.**Literature**

- 1) T. Kohno et al., J. Anal. **5**, 197 – 205 (1991)
- 2) E. Ishikawa et al., Scand. J. Immunol. **8** (Suppl. 7), 43 – 55 (1978)
- 3) M. Imagawa et al., J. Biochem. **960**, 1727 – 1735 (1983)
- 4) R. Armenta et al., Analytical. Biochemistry **146**, 211 – 219 (1985)

Cat. No. **Pack Size****11 184 024 103** custom fill

Will be supplied as "β-Galactosidase Mutein". Unit of Measure is "mg active ingredient".



For further processing only.

β-Galactosidase, recombinant, EIA Grade from *E. coli* overproducer, lyophilizate

Cat. No.	Pack Size
10 570 079 103	custom fill

Will be supplied as "β-Galactosidase, Escherichia coli". Unit of Measure is "g active ingredient".



For further processing only.

Part of marker enzyme portfolio

Application

Marker enzyme for the manufacturing of antibody- and antigen-enzyme conjugates incorporated in immunoassays for colorimetric and fluorimetric detection.

Benefits

- Synthesize stable, highly active and reproducible β-Gal antigen and antibody conjugates.
- Eliminate the risk of BSE contamination: No animal-derived components are used in the production process.

EC 3.2.1.23

Properties

Nomenclature: β-D-galactohydrolase

Molecular weight (by sequence): 465 kD

Structure: 4 identical subunits, β-galactosidase contains no carbohydrates

Isoelectric point: 4.61

Michaelis constants:

Tris buffer, pH 7.6, +20°C / relation rate:

2-nitrophenyl-β-galactoside: 9.50×10^{-4} mol/l / 1.00

phenyl-β-D-galactoside: 3.23×10^{-3} mol/l / 0.05

lactose: 3.85×10^{-2} mol/l / 0.06

4-nitrophenyl-β-galactoside: 4.45×10^{-4} mol/l / ~0.50

Activators: Mg²⁺ and Na⁺ (or other monovalent cations) are essential for activity.

pH optimum: 8.0

pH stability: 6.0

Thermal stability: Up to +37°C

Thiol groups: 64 SH groups, approximately 16 of these are accessible for SH-reactive reagents. 4 of these (Cys 76) take part in conjugation.

Specification

Appearance: White lyophilizate, stabilized with phosphate buffer and sucrose

Solubility: Clear, colorless solution in water (c=20 mg/ml)

Protein (Biuret): Approximately 0.25 mg/mg lyophilizate

Specific activity (+37°C, 2-nitrophenyl-β-D-galactopyranoside): ≥700 U/mg protein

SH-groups, free (after dialysis): ≥12 moles/mole enzyme (corresponds to 465 000 g)

Aggregated β-galactosidase (HPLC): ≤3% (dimer-part with a molar mass of 0.93×10^6 D)

Stability: At -15 to -25°C within specification range for 24 months. Store under nitrogen.

Background Information

β-Galactosidase hydrolyzes β-D-galactosides. Although the enzyme activity with 2-nitrophenyl-β-D-galactoside as substrate is higher than with the 4-isomer, the enzyme reaction with the 4-compound is more sensitive due to a higher absorption coefficient for 4-nitrophenol, ϵ_{405} : 18.5 [mmol⁻¹ x l x cm⁻¹].

Peroxidase (POD), EIA Grade from horseradish, lyophilizate

Part of marker enzyme portfolio

Application

Peroxidase (POD), EIA Grade is a marker enzyme enabling peroxidation of reduced dyes in the indicator reaction producing a color, fluorimetric or luminescent derivative of the labeled molecule for further detection and quantification.

Benefits

- Synthesize stable and reproducible antibody- and antigen-conjugates incorporated in enzyme immuno assays (EIA) representing outstanding homogeneity with respect to isoenzyme C distribution.

EC 1.11.1.7

Properties

Horseradish peroxidase is a 44,173.9 D glycoprotein with 4 lysine residue.

Specification

Appearance: Red-brown lyophilizate

Activity (+25°C, guaiacol, H₂O₂): ≥225 U/mg lyophilizate

Specific Activity (+25°C, ABTS, H₂O₂, pH 5.0): ≥900 U/mg lyophilizate

Purity number (A₄₀₃/A₂₇₅): 3.0-3.5

A₄₀₃ (0.2 mg/ml; against buffer): No limit

Contaminants (expressed as percentage of Peroxidase activity):

ATPase: ≤0.001

Catalase: ≤0.7

Phosphatase, acidic: ≤0.001

Isoenzyme distribution (HPLC): ≥90% (homogeneous with respect to isoenzyme C)

Amino groups: 2-3 moles/mole enzyme

Carbohydrates: 12.0-14.5% (w/w)

Stability: At -15 to -25°C within specification range for 24 months. Keep tightly sealed.

Literature

1) R. Presentini, B. Terrana, J. of Immunoassay **16** (3), 309 (1995)

2) P. Tijssen, E. Kurstak, Analytical Biochemistry **136**, 451 (1984)

Cat. No.	Pack Size
10 815 462 103	custom fill

Will be supplied as "Peroxidase (POD) from Horse-radish". Unit of Measure is "g".

For further processing only.

Peroxidase (POD), Grade I from horseradish, lyophilizate

Part of marker enzyme portfolio

Application

Peroxidase (POD), Grade I is a marker enzyme enabling peroxidation of reduced dyes in the indicator reaction producing a color, fluorimetric or luminescent derivative of the labeled molecule for further detection and quantification.

EC 1.11.1.7

Properties

Horseradish peroxidase is a 44,173.9 D glycoprotein with 4 lysine residue.

Specification

Appearance: Red-brown lyophilizate

Solubility: Clear, red-brown solution in water (c=10 mg/ml)

pH value (c=10 mg/ml): 6.0-7.0

Cat. No.	Pack Size
10 121 606 103	custom fill

Will be supplied as "Peroxidase (POD), Grade I, Horse-radish". Unit of Measure is "MU".

For further processing only.

Activity (+25°C, guaiacol, H₂O₂): ≥250 U/mg lyophilizate

Purity number (A₄₀₃/A₂₇₅): 3.0-3.5

Contaminants (expressed as percentage of Peroxidase activity):

ATPase: ≤0.001

Catalase: ≤0.7

Phosphatase, acidic: ≤0.001

Stability: At +2 to +8 within specification range for 24 months. Keep tightly sealed.

Poly Peroxidase (Poly POD), EIA Grade from horseradish, lyophilizate

Part of marker enzyme portfolio

Application

Poly Peroxidase (Poly POD) is a marker enzyme enabling peroxidation of reduced dyes in the indicator reaction producing a color, fluorimetric or luminescent derivative of the labeled molecule for further detection and quantification.

Benefits

- Enhance ELISA sensitivity by using Polymeric Peroxidase (Poly POD).

Product Description

Poly Peroxidase is lyophilized in 10 mmol/l potassium phosphate, 50 mmol/l NaCl, 1 mmol/l EDTA, pH 6.1 and saccharose as stabilizer.

EC 1.11.1.7

Properties

Molecular weight: 0.8 ± 0.2 x 10⁶ D (~ 20 POD-monomers)

Activation: Is accomplished by MHS (Maleimidohexanoyl-N-hydroxysuccinimide ester) ≥40 MH-groups per Poly POD (MH) are accessible for conjugation with sulfhydryl groups.

Specification

Appearance: Red-brown lyophilizate

Solubility: Clear, red-brown solution in water (c= 5 mg/ml)

Specific activity (+25°C, ABTS): ≥600 U/mg

MH-groups: ≥2 (mol MH/mol POD)

Stability: At -60 to -90°C within specification range for 48 months.

Cat. No.	Pack Size
11 578 545 103	custom fill

Will be supplied as "Peroxidase, Polymerized (MH)". Unit of Measure is "mg active ingredient".



For further processing only.

Chlorophenolred- β -D-galactopyranoside (CPRG) sodium salt, powder

Substrate for marker enzyme

Application

Use CPRG as a substrate for β -Galactosidase.

CAS: 99792-79-7

Properties

Formula: C₂₅H₂₁O₁₀Cl₂Na

Molecular weight: 607.4 D

Specification

Appearance: Orange-red powder

Solubility: Clear, red colored solution in water (c=20 mg/ml)

CPRG (A₄₀₅, ϵ =22.57 l x mmol⁻¹ x cm⁻¹): 80-110%

CPRG (HPLC): \geq 97.5%

Na (AA): 3-4%

Water (K. Fischer): \leq 15%

Chlorophenolred, free (from A₅₇₈): \leq 0.1%

Galactose, free (enzymatically) : \leq 5.0%

A₅₇₈ (c=5 mmol/l water) : \leq 0.200

A₆₅₀ (c=5 mmol/l water; turbidity): \leq 0.030

Thin layer chromatography: Corresponds to reference

Reaction rate (β -galactosidase) of sample/2-NP-galactoside: \geq 8.5%

Stability: At -15 to -25°C within specification range for 24 months. Store dry. Keep tightly closed. Ship in dry ice containers.

Literature

- 1) B. Porstmann, T. Porstmann in T.T. Ngo: Nonisotopic Immunoassay, Plenum Press, N.Y. (1988)
- 2) J. Backhaus, H. Buschek, R. Machat, M. Kuhr, W.F. Weckerle, Chromogenic Substrates for β -D-galactosidase, 4th European Carbohydrate Symposium (July 12 -17, 1987)

Cat. No. **Pack Size**

11 379 119 001 custom fill

Will be supplied as "Chlorophenolred-ss-D-galactopyranoside". Unit of Measure is "g".



For further processing only.

2,2'-azino-bis(3-ethylbenzthiazoline) 6-sulphonic acid (ABTS) solution

Substrate for marker enzyme

Application

Use ABTS as substrate solution for horse radish peroxidase (405nm).

CAS: 30931-67-0

Specification

Appearance: Slightly green solution

UV spectrum (280-450 nm): Maximum at 342 \pm 4 nm; 500 \pm 50 mE

Extinction (A₄₀₅, 1 cm, against water): \leq 60 mE

Content of ABTS: \geq 0.3 g/l

Content of hydrogen peroxide: \geq 0.01%

Function ELISA: Corresponds to specification

Stability: At +2 to +8°C within specification range for 30 months.

Cat. No. **Pack Size**

11 668 234 103 custom fill

Will be supplied as "ABTS Bulk solution". Unit of Measure is "l". For further processing only.

BM Blue POD Substrate, soluble 3,3'-5,5'-Tetramethylbenzidine (TMB), solution

Substrate for marker enzyme

Application

Use BM Blue POD Substrate as substrate solution for horse radish peroxidase.

CAS: 54827-17-7

Specification

Appearance of solution: Clear, slightly yellowish

Performance test according to TM (function test): Corresponds to specification

Stability: At +2°C to 8°C within specification range for 18 months.

Cat. No.	Pack Size
11 432 559 103	custom fill

Will be supplied as "BM blue POD Substrate, Soluble". Unit of Measure is "l".

For further processing only.

3,3',5,5'-Tetramethylbenzidine (TMB) crystalline powder

Substrate for marker enzyme

Application

Use TMB as substrate solution for horse radish peroxidase (450nm).

CAS: 54827-17-7

Properties

Formula: C₁₆H₂₀N₂

Molecular weight: 240.35 D

Specification

Appearance: Yellowish to light brown crystalline powder

Melting range: +168 to +171°C

TMB (GC): ≥99.5 area%

TMB (titrimetric, based on dry weight): ≥97.5%

Loss on drying (for 2 hours at +105°C): ≤1%

Stability: At +2 to +8°C within specification range for 24 months. Protect from light.

Cat. No.	Pack Size
10 203 700 103	custom fill

Will be supplied as "3,3',5,5'-Tetramethylbenzidine". Unit of Measure is "kg".

For further processing only.

4-Aminophenyl Phosphate (pAPP) disodium salt, powder

Substrate for marker enzyme

Application

Use 4-Aminophenyl Phosphate (pAPP), Disodium Salt as substrate for alkaline phosphatase.

CAS: 75966-16-4

Properties

Formula: C₆H₆NO₄PNa₂

Molecular weight: 233.07 D

Specification

Appearance: Off white to brownish powder

Solubility: Clear, fawn to brownish solution in water (c=100 mg/ml)

pH value: 7.0-8.0

ESI-MS: 188.011 ± 0.005 D

Cat. No.	Pack Size
05 642 965 103	custom fill

Will be supplied as "4-Aminophenyl Phosphate Disodium Salt". Unit of Measure is "kg".

For further processing only.

pAPP (HPLC): ≥90 area%

Water (K. Fischer): ≤15%

Na (flame photometric): 15-21%

p-Nitrophenylphosphate (HPLC): ≤0.3 area%

Cl: ≤0.1%

P_i (acid labile): ≤0.5%

P_i: ≤1%

Stability: At +2 to +8°C within specification range for 12 months. Protect from light.

4-Nitrophenyl Phosphate (pNPP) disodium salt, crystalline powder

Substrate for marker enzyme

Application

Use pNPP as a substrate for alkaline phosphatase.

CAS: 4264-83-9

Properties

Formula: C₆H₄NO₆PNa₂ x 6 H₂O

Molecular weight: 371.1 D (pNPP: 219.1 D)

Detection: at 405 nm

Specification

Appearance: White to slightly yellow crystalline powder

Solubility: Clear, colorless to slightly yellow solution in water (c=50 mg/ml)

pH value: 9±1

4-NPP-Na₂ x 6 H₂O (calculated from value found enzymatically): ≥95%

4-NPP (enzymatically): ≥56%

Na (flame photometric): 13 ± 1%

Water (K. Fischer): 28 ± 3%

4-Nitrophenol (free): ≤0.07%

P_i: ≤ 0.3%

Blank (with TC A_{Popl}, ΔA/30 min) : ≤0.015

Reactions rates (AP): 100 ± 5%

Stability: At +2 to +8°C within specification range for 24 months. Store dry. Protect from light.

Literature

B. Porstmann, T. Porstmann in T.T. Ngo: Nonisotopic Immunoassay, Plenum Press, N.Y. (1988)

Cat. No.	Pack Size
10 004 847 103	custom fill

Will be supplied as "4-Nitrophenyl Phosphate, Disodium Salt". Unit of Measure is "kg".

Additional formulation: Tablets are available on request.

For further processing only.

5-Bromo-4-chloro-3-indolyl-phosphate (BCIP)

disodium salt, crystalline powder

Substrate for marker enzyme

Application

Use BCIP as precipitating substrate for alkaline phosphatase.

CAS: 102185-33-1

Properties

Formula: C₈H₆NO₄BrClP x Na₂ x 1.5 H₂O

Molecular weight: 397.6 D (BCIP: 326.4 D)

Detection: Forms a blue precipitate

Cat. No.	Pack Size
10 997 846 103	custom fill

Will be supplied as "5-Br-4-chloro-3-indolyl-phosphate, Di-Na".

Unit of Measure is "kg".



For further processing only.

Specification**Appearance:** White to slightly bluish microcrystalline powder**BCIP** (from N): ≥79%**BCIP** (HPLC): ≥99 area%**N** (elementary analysis): ≥3.39%**Na** (flame photometric): 10.5-12.5%**Water** (K. Fischer) : 6.5-8.5%**Performance test** (incubation with alkaline phosphatase, aerial oxidation):

Blue precipitate

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Protect from light.

**5-Bromo-4-chloro-3-indolyl-phosphate
(BCIP)****toluidin, crystalline powder**

Substrate for marker enzyme

Application

Use BCIP as precipitating substrate for alkaline phosphatase.

CAS: 6578-06-09**Properties****Formula:** $C_8H_6NO_4BrCIP \times C_7H_9N$ **Molecular weight:** 433.6 D (BCIP: 326.4 D)**Detection:** Forms a blue precipitate**Specification****Appearance:** White to slightly yellowish microcrystalline powder**BCIP x toluidine** (from N): ≥99%**BCIP x toluidine** (HPLC): ≥99 area%**N** (elementary analysis): 6.39-6.60%**Performance test** (incubation with alkaline phosphatase, aerial oxidation):

Blue precipitate

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Protect from light.

Cat. No. **Pack Size****10 760 978 103** custom fill

Will be supplied as "5-Br-4-Cl-3-indolyl- phosphate". Unit of Measure is "g".

For further processing only.

Human Serum

frozen solution

Qualified for the Cobas® Core / Elecsys® Modular Platform.

Application

Basic matrix for manufacturing calibrators and controls for immunoassays and assays in clinical chemistry.

Benefits

Rely on the highly qualified and well characterized source of human serum.

Product Description

Frozen solution from pooled blood donors.

Specification

Appearance (optical check): Yellowish slight turbid liquid

Turbidity (A_{280} , against water): ≤ 0.600 E

pH value (+25°C): 7.0-7.5

Protein (Biuret): ≥ 61 mg/ml

Cholesterol (CHOD-PAP): ≥ 140 mg/dl

Triglyceride (GPO-PAP): 65-206 mg/dil

Ca (o-Cresolphthalein complexone): ≤ 2.2 mmol/l

Cholinesterase (Butyrylthiocholine Gen 2): ≥ 4700 U/l

Creatine kinase (IFCC): ≤ 250 U/l

Non-reactive in HBsAg: corresponds

Anti-HIV I+II: negative

Anti-HCV: negative

HCV NAT non-reactive: corresponds

HIV-1 NAT non-reactive: corresponds

Stability: At -15°C to -25°C within specification range for 12 months.

Background Information

All products derived from human blood donors are prepared exclusively from the blood donors tested individually and shown to be free from HBsAg and antibodies to HCV and HIV. The testing methods applied were FDA-approved or cleared in compliance with the European Directive 98/79/EC, Annex II, List A.

Cat. No.

Pack Size


11 758 225 103 9 L (samples 100 ml)

Will be supplied as "Serum <-> Human". Unit of Measure is "l".



For further processing only.



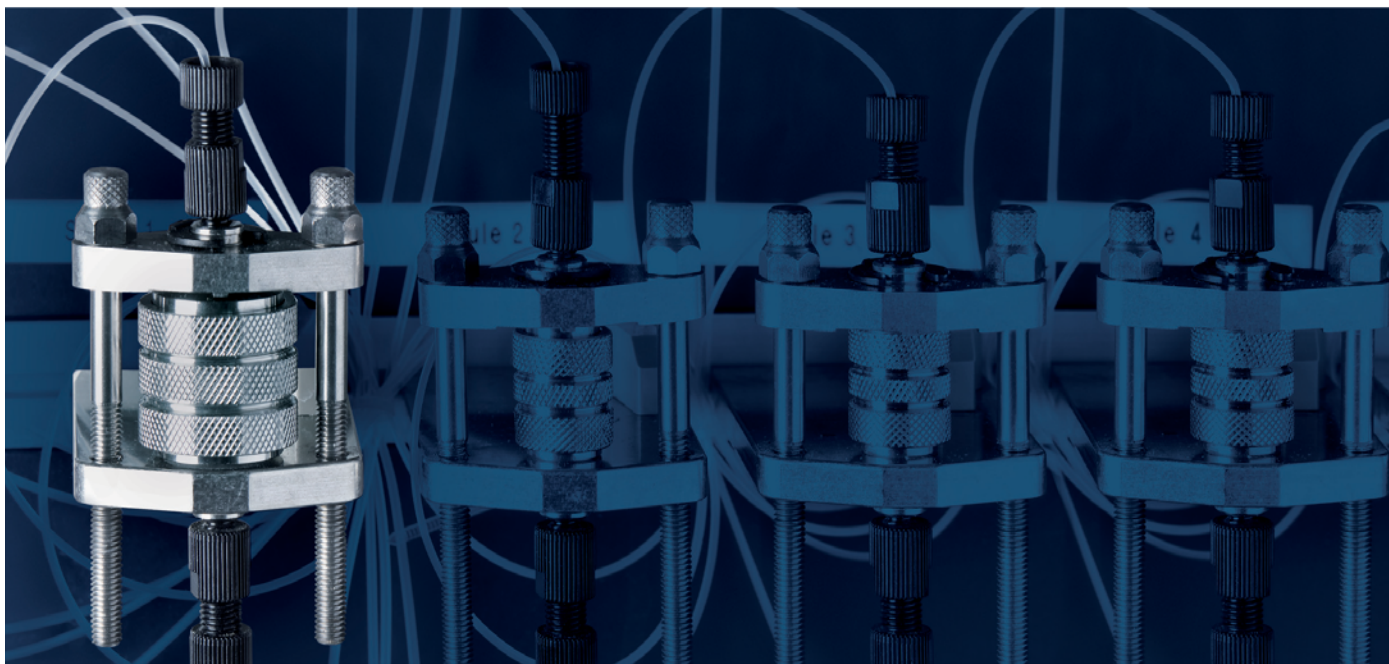

3 Molecular Diagnostics

Sample Preparation220	Labeling and Detection305
Chaotropic Salts220	Conjugates305
Enzymes220	Enzymes308
Amplification260	Labeled Nucleotides311
DNA Polymerases227	Carrier and Competitor Nucleic Acids323
Expand System	227	DNA323
T4 DNA Polymerase	232	RNA326
Taq DNA Polymerase	234	Glycogen328
Tth DNA Polymerase	238	Additional Reagents329
DNA Polymerases, Hot Start242	Enzymes229
ActiTaq Δexo DNA Polymerase	242	Proteins334
AptaTaq DNA Polymerase	243		
EagleTaq DNA Polymerase	252		
FastStart DNA Polymerase	253		
HawkZ05 DNA Polymerase	258		
HawkTaq DNA Polymerase	259		
DNA Master260		
Reverse Transcriptases267		
RNA Master271		
Nucleotides274		
deoxyNTPs	274		
dideoxyNTPs	286		
riboNTPs	290		
Additional Products296		

Guanidine Hydrochloride crystals

Application

Use Guanidine hydrochloride as a denaturing agent for proteins in a broad variety of nucleic acid purification applications.

Benefits

- **Perform high quality DNA and RNA purification.** Effectively denature all protein components in your nucleic acid.
- **Obtain intact mRNA and RNA.** Stabilize the RNA and mRNA molecules using highly concentrated Guanidinium Chloride solutions during the purification process.

CAS: 50-01-1

Properties

Formula: $\text{CH}_5\text{N}_3 \times \text{HCl}$

Molecular weight: 95.53 D

Specification

Appearance: White, crystalline powder

Chloride (qualitative): Positive

Melting range (Büchi): 183-188°C

Solubility: Clear, colorless to slightly yellow (c=764.4 mg/ml in water, 8 mol/l)

Guanidinium chloride (from N): $\geq 99\%$

Guanidinium chloride (from Cl): $\geq 99\%$

Nitrogen (elementary analysis): 43.5-44.5%

Chloride (argentometric titration): 36.7-37.5%

Heavy metals (as Pb): ≤ 10 ppm

Fe: ≤ 3 ppm

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No.

Pack Size

11 696 548 103 500 g

Will be supplied as "Guanidinium Chloride, Solid". Unit of Measure is "kg".

For further processing only.

DNase I, recombinant, Grade I from bovine pancreas, expressed in *Pichia pastoris*, lyophilizate

Recombinant produced DNase I in PCR grade, lyophilized quality, free of animal-derived materials, is an essential tool for all applications requiring DNase-free RNA templates.

Application

DNase I, recombinant, Grade I, is suitable for:

- Isolation of DNA-free RNA produced by *in vitro* transcription
- Producing DNA-free preparations of protein and RNA:
 - To ensure that RT-PCR templates are free of genomic DNA
 - To remove DNA templates after *in vitro* transcription of RNA
- Nick-translation labeling of DNA with added DNA polymerase I
- Determining the "footprint" of a DNA-binding protein
- Microarray analysis

Benefits

- **Be compliant with regulatory requirements.** Assure that your applications are free of animal-derived materials.
- **Achieve reliable results.** Experience excellent quality and greater lot-to-lot consistency due to advanced production processes in conjunction with rigorous analytical testing.

Product Description

DNase I, recombinant, Grade I, originally isolated from bovine pancreas, is a recombinant enzyme expressed in *Pichia pastoris*. It is a glycoprotein of a molecular weight of approximately 39 kD. DNase I, recombinant, Grade I, is a DNA-specific endonuclease that hydrolyzes phosphodiester linkages of double- and single-stranded DNA to a mixture of mono- and oligonucleotides. DNase I, recombinant, Grade I, is manufactured using state-of-the-art processes yielding animal component-free material.

EC 3.1.21.1

Properties

Nomenclature: DNase I

pH optimum: 7.0-8.0

Activators: DNase I requires bivalent cations for maximal activity.

Inhibitors: EDTA, EGTA, SDS

Specificity: Double-strand specific endonuclease that degrades DNA

Specification

Appearance: White to slightly yellowish lyophilizate

Activity (calf thymus DNA): ≥ 10 kU/vial lyophilizate

Activity (calf thymus DNA, modified buffer system): No limit

Unit definition: One unit according to Kunitz produces an increase in absorbance of 0.001/minute under assay conditions in 1 ml at 260 nm.

Proteases (resorufin-marked casein): Not detectable in up to 50 U after 17 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 2 U after 4 hours incubation at +37°C.

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
03 724 778 103	10 kU

Will be supplied as "DNase I rec RGI (10 KU)". Unit of Measure is "piece".

For further processing only.

DNase I, recombinant, RNase-free from bovine pancreas, expressed in *Pichia pastoris*, solution

Recombinant produced DNase I in PCR grade, solution quality, free of animal-derived materials, is an essential tool for all applications requiring DNase-free RNA templates.

Application

Use DNase I, recombinant, Grade I, for:

- Isolation of DNA-free RNA produced by *in vitro* transcription
- Producing DNA-free preparations of protein and RNA:
 - To ensure that RT-PCR templates are free of genomic DNA
 - To remove DNA templates after *in vitro* transcription of RNA
- Nick-translation labeling of DNA with added DNA polymerase I
- Determining the "footprint" of a DNA-binding protein
- Microarray analysis

Benefits

- **Be compliant with regulatory requirements.** Assure that your applications are free of animal-derived materials.
- **Achieve reliable results.** Experience excellent quality and greater lot-to-lot consistency due to advanced production processes in conjunction with rigorous analytical testing.
- **Obtain undegraded and stable RNA.** Rely on the highly purified and rigorously tested product that excludes RNase activity ensuring high sensitivity of your RT-PCR assay.

Product Description

DNase I, recombinant, RNase-free, originally isolated from bovine pancreas, is a recombinant enzyme expressed in *Pichia pastoris*. It is a glycoprotein of a molecular weight of approximately 39 kD. DNase I, recombinant, RNase-free, is a DNA-specific endonuclease that hydrolyzes the phosphodiester linkages of double- and single-stranded DNA to a mixture of mono- and oligonucleotides. DNase I, recombinant, is manufactured using state-of-the-art processes yielding animal component-free material.

The enzyme is highly purified and rigorously tested for contaminating RNase and protease activity for superb RT-PCR.

EC 3.1.21.1

Properties

Nomenclature: DNase I

pH optimum: 7.0-8.0

Activators: DNase I requires bivalent cations for maximal activity.

Inhibitors: EDTA, EGTA, SDS

Specificity: Double-strand specific endonuclease that degrades DNA

Specification

Appearance: Colorless to slightly yellow solution

Activity (calf thymus DNA): 9-14 kU/ml

Activity (calf thymus DNA, modified buffer system): No limit

Unit definition: One unit according to Kunitz produces an increase in absorbance of 0.001/minute under assay conditions in 1 ml at 260 nm.

Proteases (resorufin-marked casein): Not detectable in up to 50 U after 17 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 10 U after 4 hours incubation at +37°C.

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.

Pack Size

03 539 121 103 custom fill

Will be supplied as "DNase I rec RNase free in Glycerol". Unit of Measure is "kU".



For further processing only.

Proteinase K, recombinant, PCR Grade from *Tritirachium album*, expressed in *Pichia pastoris*, solution

Recombinant produced Proteinase K is free of animal-derived materials. In this PCR grade solution quality, this enzyme is a universal tool for a wide variety of template preparation applications.

Application

Proteinase K, recombinant, PCR Grade, digests native proteins very efficiently. This enzyme can be used to rapidly inactivate endogenous RNases and DNases during nucleic acid isolation. Proteinase K is particularly suited for the isolation of native RNA and DNA from tissues and cell lines.

The enzyme promotes cell lysis by activating a bacterial autolytic factor.

Proteinase K is also used for:

- Analysis of membrane structures by modifying proteins and glycoproteins on cell surfaces
- Removal of cellular debris during the preparation of colony lifts
- Treatment of tissue sections to ensure efficient probe infiltration during *in situ* hybridization

Benefits

- **Achieve reliable results.** Experience excellent quality and higher lot-to-lot consistency due to advanced production processes in conjunction with rigorous analytical testing.
- **Maximize the yield of target nucleic acids.** Proteinase K is rigorously tested for the absence of nucleases.
- **Effectively isolate low copy templates.** The DNA content of the enzyme preparation is strongly reduced, and the enzyme is tested for exogenous nucleic acids that may interfere in target amplification, potentially reducing sensitivity and test accuracy.

Product Description

Proteinase K, originally isolated from the mold *Tritirachium album*, is a recombinant enzyme expressed in *Pichia pastoris*. It is a highly active, subtilisin-related serine endopeptidases that does not exhibit any pronounced cleavage specificity. Thus, Proteinase K, recombinant, PCR Grade, is a universal tool for template preparation. Amino acid sequence and molecular structure of the recombinant enzyme and the native protease are identical. However, the production process of the recombinant Proteinase K guarantees an enzyme of outstanding reliability and purity meeting all the requirements of diagnostics' manufacturers.

Special emphasis has been placed on a low DNA-content of the enzyme preparation, making Proteinase K, recombinant, PCR Grade, ideally suited for isolating PCR and RT-PCR templates.

EC 3.4.23.1

Properties

Nomenclature: Proteinase K

Molecular weight: 28.8 kD

pH optimum: 7.5-10.5


Inhibitors: Proteinase K is inhibited by diisopropyl fluorophosphate and phenylmethylsulfonyl fluoride (PMSF) and is also totally inactivated by mercuric ions. Pefabloc® SC and Pefabloc® PLUS are specific, irreversible and nontoxic inhibitors.

Specificity: Proteinase K is one of the most active endopeptidases known and does not exhibit any pronounced cleavage specificity. Activity can be stimulated by addition of denaturing agents (SDS and urea).

Cat. No. **Pack Size**

03 654 672 103 850 ml, custom fill

Will be supplied as "Proteinase K, rec., PCR grade, solution". Unit of Measure is "IU".

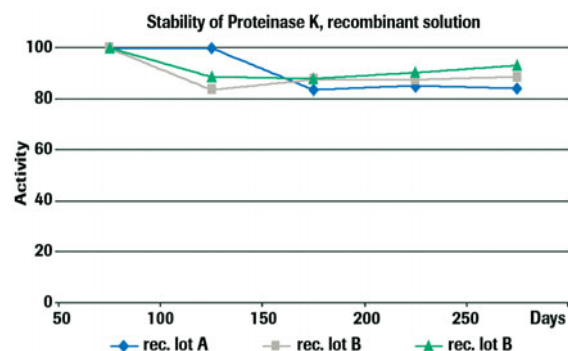

Temperature stability of recombinant Proteinase K. Accelerated stability tests at high temperature show the robustness of the recombinant enzyme. Three different lots of recombinant Proteinase K solution were tested for their temperature- stress stability at +35°C. Only minor activity loss is observed.

For further processing only.

Specification

Appearance: Clear, colorless solution

Volume activity (+37°C, hemoglobin): ≥600 U/ml

Specific activity (+37°C, hemoglobin): ≥30 U/mg protein

Unit definition (hemoglobin): One unit is the enzyme activity which releases folin positive amino acids and peptides equivalent to 1 μmol of tyrosine in 1 minute under the test conditions.

Volume activity (+25°C, Chromozym): ≥50 U/ml

Specific activity (+25°C, Chromozym): ≥2.5 U/mg protein

Unit definition (Chromozym): One unit is the enzyme activity which cleaves at +25°C in 1 minute 18 mmol Chromozym TRY.

Protein (Biuret): 14.0-22.0 mg/ml

Unspecific endonucleases (MWM III DNA): Not detectable in up to 200 μg after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 200 μg after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 40 μg after 16 hours incubation at +37°C.

DNA (Threshold): ≤10 pg/mg enzyme

Bioburden: ≤5 CFU/ml

Stability: At +2 to +8°C within specification range for 18 months.

Proteinase K, recombinant, PCR Grade from *Titrirachium album*, expressed in *Pichia pastoris*, lyophilizate

Recombinant produced Proteinase K is free of animal-derived materials. In PCR grade lyophilized quality, this enzyme is a universal tool for a wide variety of template preparation applications.

Application

Proteinase K, recombinant, PCR Grade, digests native proteins very efficiently. This enzyme can be used to rapidly inactivate endogenous RNases and DNases during nucleic acid isolation. Proteinase K is particularly suited for the isolation of native RNA and DNA from tissues and cell lines.

The enzyme promotes cell lysis by activating a bacterial autolytic factor.

Proteinase K is also used for:

- Analysis of membrane structures by modifying proteins and glycoproteins on cell surfaces
- Removal of cellular debris during the preparation of colony lifts
- Treatment of tissue sections to ensure efficient probe infiltration during *in situ* hybridization

Benefits

- **Achieve reliable results.** Experience excellent quality and higher lot-to-lot consistency due to advanced production processes in conjunction with rigorous analytical testing.
- **Maximize the yield of target nucleic acids.** Proteinase K is rigorously tested for the absence of nucleases.
- **Effectively isolate low copy templates.** The DNA content of the enzyme preparation is strongly reduced, and the enzyme is tested for exogenous nucleic acids that may interfere in target amplification, potentially reducing sensitivity and test accuracy.

Product Description

Proteinase K, originally isolated from the mold *Titrirachium album*, is a recombinant enzyme expressed in *Pichia pastoris*. It is a highly active, subtilisin-related serine endopeptidases that does not exhibit any pronounced cleavage specificity. Thus, Proteinase K, recombinant, PCR Grade, is a universal

Cat. No.	Pack Size
03 508 811 103	25 mg
03 450 376 103	50 mg
03 508 838 103	100 mg
03 450 384 103	250 mg
05 963 133 103	1 g
05 963 117 103	5 g

03450376103: Will be supplied as "Protein. K, rec PCR Grade Lyo. MPB 25mg". Unit of Measure is "piece".

03508811103: Will be supplied as "Protein. K, rec PCR Grade Lyo. MPB 50mg". Unit of Measure is "piece".

03508838103: Will be supplied as "Protein. K, rec PCR Grade Lyo. MPB 100mg". Unit of Measure is "piece".

03450384103: Will be supplied as "Protein. K, rec PCR Grade Lyo. MPB 250mg". Unit of Measure is "piece".

05963133103: Will be supplied as "Proteinase K, rec., PCR Grade, Lyo 1 g". Unit of Measure is "piece".

05963117103: Will be supplied as "Proteinase K, rec., PCR Grade, Lyo 5 g". Unit of Measure is "piece".

For further processing only.

tool for template preparation. Amino acid sequence and molecular structure of the recombinant enzyme and the native protease are identical. However, the production process of the recombinant Proteinase K guarantees an enzyme of outstanding reliability and purity meeting all requirements of diagnostics' manufacturers.

Special emphasis has been placed on a low DNA-content of the enzyme preparation, making Proteinase K, recombinant, PCR Grade, ideally suited for isolating PCR and RT-PCR templates.

EC 3.4.23.1

Properties

Nomenclature: Proteinase K

Molecular weight: 28.8 kD

pH optimum: 7.5-10.5

Inhibitors: Proteinase K is inhibited by diisopropyl fluorophosphate and phenylmethylsulfonyl fluoride (PMSF) and is also totally inactivated by mercuric ions. Pefabloc® SC and Pefabloc® PLUS are specific, irreversible and nontoxic inhibitors.

Specificity: Proteinase K is one of the most active endopeptidases known and does not exhibit any pronounced cleavage specificity. Activity can be stimulated by addition of denaturing agents (SDS and urea).

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=20 mg/ml)

Volume activity (+37°C, hemoglobin): ≥24 U/mg lyophilizate

Specific activity (+37°C, hemoglobin): ≥30 U/mg protein

Unit definition (hemoglobin): One unit is the enzyme activity which releases folin positive amino acids and peptides equivalent to 1 µmol of tyrosine in 1 minute under the test conditions.

Volume activity (+25°C, Chromozym): ≥2 U/mg lyophilizate

Specific activity (+25°C, Chromozym): ≥2.5 U/mg protein

Unit definition (Chromozym): One unit is the enzyme activity which cleaves at +25°C in 1 minute 18 mmol Chromozym TRY.

Unspecific endonucleases (MWM III DNA): Not detectable in up to 200 µg after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 200 µg after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 40 µg after 16 hours incubation at +37°C.

DNA (Threshold): ≤10 µg/mg enzyme

Bioburden: ≤125 CFU/g

Stability: At +2 to +8°C within specification range for 18 months.

RNase A**from bovine pancreas, lyophilizate, powder**

Standard RNase A in lyophilized quality is an essential tool for applications requiring RNA-free DNA templates.

Application

Use RNase A for isolation of genomic DNA. For this purpose, RNase A should be boiled before use.

Benefits

- **Rely on consistent quality.** Assure that your applications deliver consistent results due to the long-term stability of lyophilized RNase A.

EC 3.1.27.5

Properties

Molecular weight: 13.7 kD

Specificity: Pyrimidine-specific endoribonuclease that acts on single-stranded RNA

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=1 mg/ml)

Activity: ≥50 U/mg

Unit definition: One unit produces a decrease in absorbance from A_0 to A_1 in 1 minute under assay conditions (Kunitz). A_0 to A_1 corresponds to the total conversion, A_1 being the final absorbance.

Unspecific endonucleases (λ DNA): Not detectable in up to 1 μ g after 4 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 1 μ g after 4 hours incubation at +37°C.

Proteases ($\Delta A \leq 0.1$, 15 minutes, +37°C): Corresponds to reference

Turbidity, according to Maniatis ($\Delta A_{366} \leq 0.100$): Corresponds to reference

A_{280} (1 mg/ml water): 0.54-0.72

pH ≤5.0 treatment (≥3 hours): Corresponds to reference

Countries of origin: South Africa, Argentina, Australia, New Zealand, Uruguay or the United States

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Cat. No.

10 154 105 103

Pack Size

custom fill

Will be supplied as "Ribonuclease A from Bovine Pancreas". Unit of Measure is "g".

For further processing only.

Expand High Fidelity PCR System

For extremely accurate amplification of genomic DNA targets up to 5 kb using PCR.

Application

Use Expand High Fidelity PCR System for:

- Routine amplification of DNA fragments up to 5 kb from all DNA
- Amplification of DNA fragments up to 10 kb.
- Labeling of PCR products with modified nucleotides (e.g., DIG-dUTP, biotin-dUTP, fluorescein-dUTP)
- Combination with dUTP and Uracil-DNA Glycosylase for prevention of carryover contamination between PCR reactions
- Manufacture of amplification mixtures for regulated applications (e.g., *in vitro* diagnostics, quality control), including validation

Benefits

- **Improve fidelity of PCR.** Use this enzyme blend with its threefold greater accuracy than Taq Polymerase for more precise amplification of longer DNA templates.
- **Maximize target yield.** Minimize amplification of prematurely terminated products using an ideally formulated proofreading enzyme for increased full-length yields.
- **Save time using improved chemistry.** The balanced enzyme blend is optimized for high fidelity and yield.
- **Obtain consistent robust results.** The specially formulated buffer and magnesium concentration ensure high lot-to-lot consistency.

Product Description

Enzyme blend consisting of Taq DNA Polymerase and Tgo DNA Polymerase.

EC 2.7.7.7

Properties

Enzymes in Expand High Fidelity PCR System were originally isolated from the thermophilic eubacteria *Thermus aquaticus* (Taq) BM or *Thermococcus gorgonarius* (Tgo), both expressed in *E. coli*.

Enzyme activities:

Taq Polymerase: Highly processive 5'-3' DNA polymerase; double-strand-specific 5'-3' exonuclease; no 3'-5' exonuclease activity

Tgo Polymerase: Highly processive 5'-3' DNA polymerase; double-strand-specific 3'-5' exonuclease (also known as proofreading activity); no 5'-3' exonuclease activity.

pH optimum: Approximately 8.9 (+20°C)

Temperature optimum:

Fragment length <3 kb: Approximately +72°C

Fragment length >3 kb: Approximately +68°C

Substrates: Incorporates dNTP, dUPT, various labeled or modified nucleotides (200 µmol/l each is recommended of normal dNTP; increased concentrations of variants)

Divalent ion requirement: Mg²⁺ (1.5 mmol/l standard concentration)

Recommended usage per 50 µl reaction: 2.5 U (0.7 µl)

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; EDTA, 0.1 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50% (v/v); pH approximately 8.0 at +4°C

Volume activity: ≥3.5 U/µl

RNases (MSII-RNA): Not detectable in incubation with up to 30 U after 1 hour at +37°C.

Cat. No. **Pack Size**

03 310 256 103 custom fill

Will be supplied as "Expand High Fidelity". Unit of Measure is "kU". The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 48

- For the best fit reaction buffer, use Expand High Fidelity PCR Buffer, see page 228

Function test in PCR (200 ng human genomic DNA, 4.8 kb tPA fragment):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Expand High Fidelity PCR Buffer 10x conc., with MgCl₂

Standard reaction buffer for PCR using the Expand High Fidelity PCR System.

Application

Use this buffer together with Expand High Fidelity PCR System.

Benefits

- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer for fast and easy setup of highly reproducible PCR experiments.
- **Gain excellent performance.** Take full advantage of the Expand High Fidelity PCR System using a specially optimized reaction buffer.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 500 mmol/l; (NH₄)₂SO₄, 220 mmol/l; MgCl₂, 15 mmol/l; pH approximately 8.9 at +25°C

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Function test in PCR (human genomic DNA, 4.8 kb tPA fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Expand High Fidelity PCR Buffer 10x conc., without MgCl₂

Standard reaction buffer without MgCl₂ for optimization of MgCl₂ concentration in PCR using the Expand High Fidelity PCR System.

Application

Use this buffer together with Expand High Fidelity PCR System. whenever the amplification of difficult targets requires a MgCl₂ concentration that needs to be individually optimized.

Benefits

- **Amplify difficult targets.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer with an optimized MgCl₂ concentration for best results.
- **Gain excellent performance.** Take full advantage of the Expand High Fidelity PCR System using a specially optimized reaction buffer.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 500 mmol/l; (NH₄)₂SO₄, 220 mmol/l; pH approximately 8.9 at +25°C

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Cat. No.	Pack Size
05 917 131 103	1 ml 1 piece contains 1 ml.

Will be supplied as "Exp.HF Buffer 10x w MgCl₂ MPB". Unit of Measure is "piece".



DRY ICE

For further processing only.

Cat. No.	Pack Size
05 917 123 103	1 ml 1 piece contains 1 ml.

Will be supplied as "Exp.HF Buff. 10x w/o MgCl₂ MPB". Unit of Measure is "piece".



DRY ICE

For further processing only.

Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Function test in PCR (human genomic DNA, 4.8 kb tPA fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Expand Long Template PCR System

For extremely accurate amplification of genomic DNA targets up to 20 kb using polymerase chain reaction (PCR).

Application

Use Expand Long Template PCR System for:

- Routine amplification of DNA fragments up to 20 kb from all DNA
- Amplification of DNA fragments up to 40 kb from λDNA
- Labeling of PCR products with modified nucleotides (e.g., DIG-dUTP, biotin-dUTP, fluorescein-dUTP)
- Combination with dUTP and Uracil-DNA Glycosylase for prevention of carryover contamination between PCR reactions
- Manufacture of amplification mixtures for regulated applications (e.g., *in vitro* diagnostics, quality control), including validation

Benefits

- **Amplify longer templates than ever before.** Generate PCR products 5 to 20 kb in length from complex genomic DNA using this optimized enzyme blend.
- **Achieve higher yields and fidelity.** Three times higher fidelity with much higher yield compared to Taq DNA Polymerase.
- **Improve PCR efficiency.** More full-length product characterizing human gene loci, fingerprinting DNA, and isolating entire genes from cDNA or entire viral genomes.

Product Description

Enzyme blend consisting of Taq DNA Polymerase and Tgo DNA Polymerase.

EC 2.7.7.7

Properties

Enzymes in the Expand Long Template PCR System were originally isolated from the thermophilic eubacteria *Thermus aquaticus* (Taq) BM and *Thermococcus gorgonarius* (Tgo), both expressed in *E. coli*.

Enzyme activities:

Taq Polymerase: Highly processive 5'-3' DNA polymerase, double-strand specific 5'-3' exonuclease, no 3'-5' exonuclease activity

Tgo Polymerase: Highly processive 5'-3' DNA polymerase, double-strand specific 3'-5' exonuclease (also known as proofreading activity), no 5'-3' exonuclease activity

Temperature optimum:

Fragment length <3 kb: Approximately +72°C

Fragment length >3 kb: Approximately +68°C

Substrates: Incorporates dNTP, dUPT, various labeled or modified nucleotides

Divalent ion requirement: Mg²⁺ (1.75 mmol/l when using 350 µmol/l of each dNTP; 2.75 mmol/l when using 500 µmol/l of each dNTP)

Recommended usage per 50 µl reaction: 0.5-5.0 U (3.75 U standard concentration)

Specification

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; EDTA, 0.1 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50% (v/v); pH approximately 8.0 at +4°C

Cat. No. **Pack Size**

03 321 053 103 custom fill

Will be supplied as "Expand LT PCR Sys. Enzymmix, Bulk". Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 48

- For the best fit reaction buffer, use Expand Long Template PCR Buffer 1, see page 230
- For the best fit reaction buffer, use Expand Long Template PCR Buffer 2, see page 230
- For the best fit reaction buffer, use Expand Long Template PCR Buffer 3, see page 231

Volume activity: ≥ 5 U/ μ l

RNases (MSII-RNA): Not detectable in incubation with up to 30 U after 1 hour at +37°C.

Function test in PCR (200 ng human genomic DNA; 9.3, 12, and 15 kb fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Expand Long Template PCR Buffer 1 10x conc., with 17.5 mM MgCl₂

Reaction buffer 1 for PCR of fragments from 0.5 to 12 kb using the Expand Long Template PCR System.

Application

Use this buffer together with Expand Long Template PCR System.

Benefits

- **Gain excellent performance.** Take full advantage of the Expand Long Template PCR System using the specially optimized reaction buffer 1 for amplification of fragments from 0.5 to 12 kb.
- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled reaction buffer for fast and easy setup of highly reproducible PCR experiments.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 500 mmol/l; (NH₄)₂SO₄, 160 mmol/l; MgCl₂, 17.5 mmol/l; pH approximately 9.2 at +25°C

Unspecific endonucleases (λ DNA and MWM III DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Function test in PCR (200 ng human genomic DNA, 9.3 kb tPA fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
05 917 115 103	1 ml

Will be supplied as "Expand LT PCR Buffer 1 (10x), 1 ml". Unit of Measure is "piece".

For further processing only.

Expand Long Template PCR Buffer 2 10x conc., with 27.5 mM MgCl₂

Reaction buffer 2 for PCR of fragments from 12 to 15 kb using the Expand Long Template PCR System.

Application

Use this buffer together with Expand Long Template PCR System.

Benefits

- **Gain excellent performance.** Take full advantage of the Expand Long Template PCR System using the specially optimized reaction buffer 2 for amplification of fragments from 12 to 15 kb.
- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled reaction buffer for fast and easy setup of highly reproducible PCR experiments.

Cat. No.	Pack Size
05 420 075 103	1 ml

Will be supplied as "Expand LT PCR Buffer 2 (10x), 1 ml". Unit of Measure is "piece".

For further processing only.

Specification

Appearance: Clear, colorless solution
Contents: Tris/HCl, 500 mmol/l; (NH₄)₂SO₄, 160 mmol/l; MgCl₂, 27.5 mmol/l; DMSO, 20%; Tween 20, 1%; pH approximately 9.2 at +25°C
Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.
Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.
Function test in PCR (200 ng human genomic DNA, 12 kb tPA fragment): Corresponds to specification
Stability: At -15 to -25°C within specification range for 24 months.

Expand Long Template PCR Buffer 3
10x conc., with 27.5 mM MgCl₂

Reaction buffer 3 for PCR of fragments larger than 15 kb using the Expand Long Template PCR System.

Application

Use this buffer together with Expand Long Template PCR System.

Benefits

- **Gain excellent performance.** Take full advantage of the Expand Long Template PCR System using the specially optimized reaction buffer 3 for amplification of fragments larger than 15 kb.
- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled reaction buffer for fast and easy setup of highly reproducible PCR experiments.

Specification

Appearance: Clear, colorless solution
Contents: Tris/HCl, 500 mmol/l; (NH₄)₂SO₄, 220 mmol/l; MgCl₂, 27.5 mmol/l; DMSO, 20%; Tween 20, 1%; pH approximately 9.2 at +25°C
Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.
Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.
Function test in PCR (200 ng human genomic DNA, 15 kb tPA fragment): Corresponds to specification
Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
05 420 083 103	1 ml

Will be supplied as "Expand LT PCR Buffer 3 (10x), 1 ml". Unit of Measure is "piece".

For further processing only.

T4 DNA Polymerase

from T4 plasmid pTL43W infected *E. coli* 71-18, solution

Use T4 DNA Polymerase for 3' labeling of DNA.

Application

Use T4 DNA Polymerase for:

- Labeling of 3' termini of DNA. Extensive labeling is achieved using the replacement reaction in which the 3'-exonuclease activity of the enzyme first digests dsDNA to produce molecules with recessed 3' termini. After addition of high concentrations of labeled dNTP, the polymerase activity extends the 3' ends along the length of the template
- Gap-filling in site-directed mutagenesis experiments, in combination with T4 Gene 32 Protein

Benefits

- **Obtain improved performance.** Take advantage of this highly processive, contamination-controlled T4 DNA Polymerase.

EC 2.7.7.7

Properties

Enzyme activities: T4 DNA Polymerase is a DNA-dependent DNA polymerase which catalyzes the addition of dNTP to hydroxyl-termini of recessive ends of double-stranded DNA. With low dNTP concentrations, the enzyme has an extreme 3'-5' exonuclease activity. It has no 5'-3' exonuclease activity.

dNTP concentration: With <1 μmol/l dNTP, the exonuclease activity of T4 DNA Polymerase dominates and template DNA is degraded; 1 to 2 μmol/l dNTP are used for specific polymerase activity in labeling experiments. For maximal polymerase activity, exonuclease activity can be suppressed using high dNTP concentrations up to 100 μmol/l.

pH optimum: 8.0-9.0

Divalent ion requirement: Mg²⁺

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 10 mmol/l; MgCl₂, 2.5 mmol/l; NaCl, 100 mmol/l; DTE, 2 mmol/l; EDTA, 0.5 mmol/l; glycerol, 50% (v/v); pH approximately 8.0 at +4°C

Volume activity: ≥1 U/μl

Unit definition: One unit T4 DNA Polymerase is defined as the amount of enzyme which catalyzes the incorporation of 10 nmol [³H]dNTP into acid insoluble DNA in 30 minutes at +37°C.

Purity (SDS PAGE): ≥90%

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 40 U after 4 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 40 U after 4 hours incubation at +37°C.

Nuclease activity on single-stranded DNA (M13mp9 ssDNA): Not detectable in up to 40 U after 4 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 40 U after 4 hours incubation at +37°C.

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. **11 004 778 103** Pack Size **custom fill**

Will be supplied as "DNA Polymerase, T4". Unit of Measure is "kU". The enzyme is supplied without reaction buffer.



For further processing only.

- For the best fit reaction buffer, use T4 DNA Polymerase Incubation Buffer, see page 233

T4 DNA Polymerase Incubation Buffer 5x concentrated

Standard reaction buffer for 3' labeling of DNA using T4 DNA Polymerase.

Application

Use this T4 DNA Polymerase Incubation Buffer together with T4 DNA Polymerase. For a list of applications, refer to **T4 DNA Polymerase**.

Benefits

- **Gain excellent performance.** Take full advantage of T4 DNA Polymerase using a reaction buffer specially optimized for this enzyme.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 0.25 mmol/l; MgCl₂, 35 mmol/l; β-Mercaptoethanol, 50 mmol/l; EDTA, 0.5 mmol/l; (NH₄)₂SO₄, 75 mmol/l; BSA, 0.1 mg/ml; pH approximately 8.8 at +25°C **Unspecific endonucleases** (λDNA and pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +37°C.

Function test in combination with T4 DNA Polymerase: Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
05 187 168 103	1 ml

Will be supplied as "T4 DNA Polymerase Incubation Buffer 5x".
Unit of Measure is "piece".

For further processing only.

Taq DNA Polymerase, GMP Grade, 5 U/μl from *Thermus aquaticus* BM, expressed in *E. coli*, solution

Taq DNA Polymerase is the robust standard enzyme for the amplification of DNA fragments up to 3 kb in the polymerase chain reaction (PCR).

Application

Use Taq DNA Polymerase, GMP Grade, 5 U/μl, for:

- Routine PCR and RT-PCR applications
- Amplification of DNA fragments up to 3 kb from various sources of DNA
- Labeling of DNA with modified nucleotides (e.g., DIG-dUTP, biotin-dUTP, fluorescein-dUTP)
- Combination with dUTP and Uracil-DNA Glycosylase for prevention of carryover contamination between PCR reactions
- Manufacture of amplification mixtures for applications with regulatory requirements (e.g., *in vitro* diagnostics, quality control)

Benefits

- **Obtain consistent results.** Rely on the robust reaction performance and the lot-to-lot consistency of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to GMP (Good Manufacturing Practice) regulations.
- **Profit from cost efficiency.** Benefit from low cost per reaction.

EC 2.7.7.7

Properties

Taq DNA Polymerase is the recombinant full-length version of the thermostable enzyme from the eubacterium *Thermus aquaticus* BM, expressed in *E. coli*.

Enzyme activities: Highly processive 5'-3' DNA polymerase; double-strand specific 5'-3' exonuclease; no 3'-5' exonuclease activity

pH optimum: Approximately 9.0 (+20°C)

Temperature optimum: Approximately +75°C

Half life at +95°C: Approximately 40 minutes

Substrates: Incorporates dNTP, dUPT, dITP, various labeled or modified nucleotides (200 μmol/l each is recommended of normal dNTP; increased concentrations of variants)

Divalent ion requirement: Mg²⁺ (1.5 mmol/l standard concentration)

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50% (v/v); pH approximately 8.0 at +4°C.

Volume activity: ≥5 U/μl

Specific activity (Protein: A₂₈₀): ≥130,000 U/mg

Unit definition: One unit Taq DNA Polymerase is defined as the amount of enzyme that incorporates 10 nmol of total deoxyribonucleosidetriphosphates into acid precipitable DNA within 30 minutes at +75°C under standard assay conditions.

Purity (SDS PAGE): ≥98%

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 10 U after 1 hour incubation at +37°C.

Function test in PCR (10 pg λDNA, 0.5 kb fragment): Corresponds to

Cat. No. Pack Size

03 707 610 103 1 kU

03 707 628 103 5 kU

03 161 455 103 50 kU

03707610103: Will be supplied as "Taq DNA Polymerase Ind. GMP Grade, 1 kU". Unit of Measure is "piece".

03707628103: Will be supplied as "Taq DNA Polymerase Ind. GMP Grade, 5 kU". Unit of Measure is "piece".

03161455103: Will be supplied as "Taq DNA Polym GMP Grade 50ku". Unit of Measure is "piece".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 48

► For the best fit reaction buffer, use PCR Buffer ,see page 237

► For the best fit reaction buffer, use PCR Buffer Without MgCl₂, 10x concentrated, see page 237

reference

Function test in qPCR using LightCycler®

(human genomic DNA, β -globin gene): Corresponds to reference

(plasmid DNA, β -globin gene): Corresponds to reference

Bioburden: ≤ 50 CFU/ml

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Taq DNA Polymerase, 5 U/ μl from *Thermus aquaticus* BM, expressed in *E. coli*, solution

Taq DNA Polymerase is the robust standard enzyme for the amplification of DNA fragments up to 3 kb in the polymerase chain reaction (PCR).

Application

For applications see *Taq DNA Polymerase, GMP Grade, 5 U/ μl*

Benefits

- **Obtain consistent results.** Rely on the robust reaction performance and lot-to-lot consistency of this product.
- **Profit from cost efficiency.** Benefit from low cost per reaction.

EC 2.7.7.7

Properties

See *Taq DNA Polymerase, GMP Grade, 5 U/ μl*

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50% (v/v), pH approximately 8.0 at $+4^\circ\text{C}$

Volume activity: ≥ 5 U/ μl

Unit definition: One unit Taq DNA polymerase is defined as the amount of enzyme that

incorporates 10 nmol of total deoxyribonucleosidetriphosphates into acid precipitable DNA within 30 minutes at $+75^\circ\text{C}$ under standard assay conditions.

Unspecific endonucleases (λ DNA): Not detectable in up to 30 U after 16 hours incubation at $+37^\circ\text{C}$.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at $+37^\circ\text{C}$.

Exonucleases (^3H -DNA): Not detectable in up to 30 U after 4 hours incubation at $+65^\circ\text{C}$.

Function test in PCR using conventional blockcycler (10 μg λ DNA, 0.5 kb fragment): Corresponds to reference

Function test in qPCR using LightCycler® System

(human genomic DNA, β -globin gene): Corresponds to reference

(plasmid DNA, β -globin gene): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
11 147 633 103	custom fill

Will be supplied as "Taq DNA Polymerase". Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 48

- ▶ For the best fit reaction buffer, use PCR Buffer, see page 237
- ▶ For the best fit reaction buffer, use PCR Buffer Without MgCl_2 , 10x concentrated, see page 237

Taq DNA Polymerase, 50 U/μl from *Thermus aquaticus* BM, expressed in *E. coli*, glycerol-free solution

Taq DNA Polymerase is the robust standard enzyme for the amplification of DNA fragments up to 3 kb in the polymerase chain reaction (PCR).

Application

Use Taq DNA Polymerase, 50 U/μl, especially for:

- Setup of PCR master mixtures, when highly concentrated components are required
- Preparation of dried amplification mixtures for more convenience and increased stability at ambient temperature

For further applications see **Taq DNA Polymerase, GMP Grade, 5 U/μl**

Benefits

- **Prepare dried amplification mixtures.** Use this formulation for manufacture of dried-down reagents with high stability and convenience.
- **Obtain consistent results.** Rely on the robust reaction performance and the lot-to-lot consistency of this product.
- **Profit from cost efficiency.** Benefit from low cost per reaction.

Product Description

High concentrated, glycerol-free solution, ideal for preparation of dried-down amplification mixtures.

EC 2.7.7.7

Properties

See **Taq DNA Polymerase, GMP Grade, 5 U/μl**

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); pH approximately 8.0 at +4°C

Glycerol content: ≤0.1% (v/v)

Volume activity: 55±5 U/μl

Unit definition: One unit Taq DNA Polymerase is defined as the amount of enzyme that incorporates 10 nmol of total deoxyribonucleosidetriphosphates into acid precipitable DNA within 30 minutes at +75°C under standard assay conditions.

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 30 U after 4 hours incubation at +65°C.

Function test in qPCR using LightCycler® 480 System (human genomic DNA, tPA gene): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
04 827 007 103	custom fill

Will be supplied as "Taq DNA Pol., Glycerol-free". Unit of Measure is "KU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 48

PCR Buffer

10x conc., with 15 mM MgCl₂

Standard reaction buffer for PCR using Taq DNA Polymerase.

Application

Use this buffer together with Taq DNA Polymerase. For a list of applications, refer to *Taq DNA Polymerase, GMP Grade, 5 U/μl*.

Benefits

- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer for fast and easy setup of highly reproducible PCR experiments.
- **Achieve excellent PCR performance.** Take full advantage of Taq DNA Polymerase, GMP Grade, using this specially optimized reaction buffer.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 100 mmol/l; KCl, 500 mmol/l; MgCl₂, 15 mmol/l; pH approximately 8.3 at +20°C

Unspecific endonucleases (λDNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Function test in PCR (0.01 ng λDNA, 0.5 kb lambda fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
11 974 769 103	1 ml
11 271 326 103	custom fill

11271318103: Will be supplied as "PCR Buffer 10x 1 ml". Unit of Measure is "piece". 11271326103: Will be supplied as "PCR Buffer 10x". Unit of Measure is "ml".



For further processing only.

PCR Buffer

10x conc., without MgCl₂

Standard PCR reaction buffer without MgCl₂ using Taq DNA Polymerase for individual MgCl₂ optimization.

Application

Use this buffer without MgCl₂ together with Taq DNA Polymerase, GMP Grade, 5 U/μl for amplification of difficult targets requiring a MgCl₂ concentration that is optimized individually.

Benefits

- **Amplify difficult DNA targets using an improved buffer.** For best results use this premixed, pH-adjusted, contamination-controlled reaction buffer with individually optimized MgCl₂ concentrations.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 100 mmol/l; KCl, 500 mmol/l; pH approximately 8.3 at +20°C

Unspecific endonucleases (λDNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Function test in PCR (0.01 ng λDNA, 0.5 kb lambda fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
11 600 753 103	1 ml
11 600 761 103	custom fill

11600753103: Will be supplied as "PCR Buffer(10x) w/o MgCl₂ MPB". Unit of Measure is "piece".

11600761103: Will be supplied as "PCR Buffer (10x) w/o MgCl₂ Bulk". Unit of Measure is "ml".



For further processing only.

Tth DNA Polymerase

from *Thermus* species, expressed in *E. coli*, solution

Tth DNA Polymerase is a thermostable DNA Polymerase with intrinsic reverse transcriptase activity for RT-PCR amplification of RNA to a length of at least 1 kb.

Application

Use Tth DNA Polymerase for:

- One-step RT-PCR of single copy genes from eucaryotic genomes in the presence of Mn²⁺ ions
- One-step RT-PCR of a specific transcript or an entire population of transcripts
- PCR in the presence of Mg²⁺ ions
- Labeling of PCR products with modified nucleotides

Benefits

- **Amplify directly from RNA.** Benefit from the intrinsic reverse transcriptase activity of Tth Polymerase to directly amplify from RNA in one step.
- **Improve PCR yield.** Obtain more PCR product, because the Tth Polymerase is stable during prolonged repetitive high temperature incubations.
- **Enhance specificity of amplification.** Avoid loss of specificity due to RNA secondary structure using the higher annealing temperature Tth Polymerase allows compared to other reverse transcriptases.

EC 2.7.7.7

Properties

Tth Polymerase is the recombinant version of the thermostable enzyme from the thermophilic eubacterium *Thermus thermophilus* species, expressed in *E. coli*.

Enzyme activities: Highly processive 5'-3' DNA polymerase; no 3'-5' exonuclease activity; very efficient intrinsic reverse transcriptase (RT) activity in the presence of manganese ions; no RNase H activity

pH optimum: Approximately 9.0 (+25°C)

Temperature optimum for elongation: Approximately +72°C

Temperature optimum for reverse transcription: Approximately +60 to +70°C

Divalent ion requirement for PCR: Mg²⁺

Divalent ion requirement for RT activity and RT-PCR: Mn²⁺

Substrates: Incorporates dNTP, dUTP, dITP, various labeled or modified nucleotides (200 µmol/l each is recommended of normal dNTP; increased concentrations of variants)

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 10 mmol/l; KCl, 300 mmol/l; EDTA, 0.1 mmol/l; DTT, 1 mmol/l; Triton X-100, 0.1% (v/v); glycerol, 50% (v/v); pH approximately 7.5 at +25°C

Volume activity: ≥ 5 U/µl

Unit definition: One unit Tth DNA Polymerase is defined as the amount of enzyme that incorporates 10 nmol of total deoxyribonucleosidetriphosphates into acid precipitable DNA within 30 minutes at +75°C under standard assay conditions.

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 U after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 U after 16 hours incubation at +65°C.

Ribonucleases (MS2 RNA): Not detectable in up to 20 U after 4 hours

Cat. No. **11 485 954 103** Pack Size custom fill

Will be supplied as "Tth DNA Polymerase". Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 01

- For optimal reverse transcription activity of Tth DNA Polymerase, add Mn(OAc)₂ Stock Solution, see page 239
- For the best fit reaction buffer, use Tth DNA Polymerase Incubation Buffer, see page 239

incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 20 U after 4 hours incubation at +65°C.

Function test in PCR (10 ng human genomic DNA, 1.1 kb collagen fragment): Corresponds to specification

Function test in RT-PCR (10 ng human liver RNA, 630 bp MCAD fragment): Corresponds to specification

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 18 months.

Tth DNA Polymerase PCR Buffer 10x conc., with 15 mM MgCl₂

Standard reaction buffer for PCR using Tth DNA Polymerase.

Application

Use this Tth DNA Polymerase PCR Buffer in combination with Tth DNA Polymerase for PCR applications.

Benefits

- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer for fast and easy setup of highly reproducible PCR experiments.
- **Gain better performance.** Take full advantage of Tth DNA Polymerase using a reaction buffer specially optimized for this enzyme.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 100 mmol/l; KCl, 1 mol/l; MgCl₂, 15 mM; BSA, 500 µg/ml; Tween 20, 0.5%; pH approximately 8.9 at +25°C

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Ribonucleases (MS2 RNA): Not detectable in up to 20 µl after 4 hours incubation at +37°C.

Function test (10 ng human genomic DNA, 1.1 kb collagen fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
05 187 176 103	1 ml

Will be supplied as "Tth PCR Buffer 10x, 1ml". Unit of Measure is "piece".

For further processing only.

Mn(OAc)₂ Stock Solution 25 mM

RT-PCR grade Mn-acetate solution.

Application

Use this Mn(OAc)₂ Stock Solution in combination with the Tth DNA Polymerase RT-PCR Buffer to optimize the manganese concentration.

Benefits

- **Obtain reliable results.** Rely on the high lot-to-lot consistency of this product, thoroughly tested for constant quality.
- **Simplify your RT-PCR setup.** Save time producing a suitable, pure manganese solution by using this ready-to-use formulation.

Specification

Appearance: Clear, colorless to pinkish solution

Cat. No.	Pack Size
05 187 109 103	1 ml

Will be supplied as "Mn(OAc)₂ Stock Solution, 25 mM". Unit of Measure is "piece".

For further processing only.

DNA Polymerases

Tth DNA Polymerase

Contents: Mn-acetate, 25 mmol/l

Unspecific endonucleases (λ DNA and MWM III DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Function test (10 ng human liver RNA, 630 bp MCAD fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Tth DNA Polymerase RT-PCR Buffer 5x concentrated

Standard reaction buffer for RT-PCR using Tth DNA Polymerase.

Application

Use this Tth DNA Polymerase RT-PCR Buffer in combination with Tth DNA Polymerase for RT-PCR applications.

Benefits

- **Retranscribe and amplify difficult targets.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer with a customized manganese concentration for best results.
- **Gain better performance.** Take full advantage of Tth DNA Polymerase using a reaction buffer specially optimized for this enzyme.

Specification

Appearance: Clear, colorless solution

Contents: Bicine/KOH, 0.25 mol/l; potassium acetate, 575 mmol/l; glycerol, 40% (v/v); pH approximately 8.2

Unspecific endonucleases (λ DNA and MWM III DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Function test (10 ng human liver RNA, 630 bp MCAD fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
05 187 079 103	1 ml

Will be supplied as "Tth RT-PCR Buffer 5x, 1ml". Unit of Measure is "piece".

For further processing only.

RMS Z05 DNA Polymerase, 200 U/ μ L from *Thermus* species Z05, expressed in *E. coli*, solution

RMS Z05 DNA Polymerase is a one-step RT-PCR DNA polymerase.

Application

Use RMS Z05 DNA Polymerase for:

- Single buffer RT-PCR
- Reverse transcription of RNA targets up to 1kb
- Incorporation of modified nucleotides for labeling of PCR products
- Detection formats such as hydrolysis probes, hybridization probes and SYBR Green

Benefits

- Simplify PCR setup by using an RNA- and DNA-dependent DNA polymerase activity for PCR and RT-PCR.
- Rely on high lot-to-lot consistency and full traceability with GMP manufacturing.
- Obtain increased specificity, sensitivity, efficiency, and RT-PCR yield due to high temperature RT.

Cat. No.	Pack Size
05 206 979 190	20 kU
05 206 987 190	200 kU

05206979190: Will be supplied as "CMPNT RMS Z05 20KU, 200U/uL 0.1mL". Unit of Measure is "piece".

05206987190: Will be supplied as "CMPNT RMS Z05 200KU, 200U/uL 1.0mL". Unit of Measure is "piece".

For further processing only.

Patent and License Disclaimer(s): 66

- Control carryover contamination with a dUTP and Uracil-DNA Glycosylase treatment (UNG) compatible setup.
- Enjoy long product shelf life with 24 months stability.

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Tween 20, 0.2% (v/v); glycerol, 50.0% (v/v); pH approximately 8.0 at +4°C

Volume activity: 200-240 U/μl

Purity: One major band

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

DNA Polymerases, Hot Start

ActiQaq Δexo DNA Polymerase

ActiQaq Δexo DNA Polymerase

from *Thermus aquaticus* BM, expressed in *E. coli*, solution

ActiQaq Δexo DNA Polymerase is a chemically modified, N-terminal truncated Taq DNA polymerase without 5'-3' exonuclease activity. This modified DNA polymerase is ideal for specific, sensitive DNA amplification using hot start PCR with enzyme activation during initial heat denaturation to optimally detect mismatches.

Application

Use ActiQaq Δexo DNA Polymerase for:

- SNP analysis and genotyping
- Allele-specific PCR
- Multiplexing up to 350 bp
- Random primed PCR

Benefits

- **Optimize your SNP analysis.** Discriminate paired and unpaired primer ends using an enzyme optimized for allele-specific PCR.
- **Increase specificity of PCR.** The 5'-3' exonuclease activity without hot start will not amplify nonspecific, low-temperature primer-template hybrids.
- **Fine-tune your hybridization probe qPCR.** Avoid hybridization probe degradation using a 5'-3' exonuclease activity lacking polymerase.

EC 2.7.7.7

Properties

ActiQaq Δexo DNA Polymerase is designed for hot start PCR and must be heat-activated at the beginning of the PCR.

Enzyme activities: Highly processive 5'-3' DNA polymerase; no 5'-3' exonuclease activity; no 3'-5' exonuclease activity

Heat-activation: +95°C for 3-10 minutes (assay-dependent, recommendation is 10 minutes for full activation)

pH optimum: Approximately 8.3 (+25°C)

Temperature optimum: Approximately +72°C

Standard reaction buffer: Tris/HCl, 10 mmol/l; KCl, 6.25 mmol/l final concentration at pH 8.3

Substrates: Incorporates dNTP, dUPT, dITP, various labeled or modified nucleotides (200 μmol/l each is recommended of normal dNTP, increased concentrations of variants)

Divalent ion requirement: Mg²⁺ (2.0 mmol/l standard concentration)

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; EDTA, 0.1 mmol/l; DTT, 1 mmol/l; Tween 20; 0.5% (v/v); glycerol, 50% (v/v); pH approximately 8.5 at +25°C

Volume Activity: 4.4±0.4 U/μl

Residual activity prior reactivation: <5%

Performance test in qPCR using LightCycler® 480: Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. **03 788 075 103** Pack Size **custom fill**

Will be supplied as "ActiQaq delta exo DNA Polymerase". Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 49

AptaTaq DNA Polymerase, 5 U/μl from *Thermus aquaticus* BM, expressed in *E. coli*, solution

The novel AptaTaq hot start PCR technology combines the native speed and robustness of Taq DNA Polymerase with a fast hot start system ensuring sensitive, and specific amplification of the target DNA.

Application

Apply AptaTaq DNA Polymerase for:

- Fast PCR assays with no extra enzyme activation time and fast cycling protocols
- Single- and multiplex PCR and qPCR applications that require high specificity, sensitivity, and yield
- RT-PCR
- Difficult templates, such as complex secondary structures or GC-rich sequences
- Automated PCR workflows requiring high stability of the reaction mixtures during automated pipetting and prolonged handling at room temperature

Benefits

- **Reduce time to result.** Save up to 15 minutes per run by omitting the initial activation step required by chemically modified hot start polymerases, and reduce cycling time with fast protocols.
- **Maximize specificity, sensitivity, and yield.** Achieve reliable amplification of your target DNA from various sources (*e.g.*, genomic DNA, cDNA, plasmids).
- **Simplify PCR setup.** Store these highly stable polymerase for up to 1 month at +2° to +8°C and set up your hot start PCR reaction at room temperature.
- **Obtain consistent results.** Rely on Roche's standardized manufacturing processes, including extensive Quality Control release testing, resulting in high lot-to-lot consistency providing the perfect basis for (IVD) kit manufacturers and end users.

Product Description

AptaTaq DNA Polymerase is a blend of Taq DNA Polymerase and a specific oligonucleotide (aptamer) providing hot start features.

EC 2.7.7.7

Properties

AptaTaq DNA Polymerase is reversibly inhibited below +55°C and becomes active at temperatures over +60°C. This hot start feature eliminates the risk of nonspecific primer extension during PCR setup.

Enzyme activities: Highly processive 5'-3' DNA polymerase; double-strand-specific 5'-3' exonuclease; no 3'-5' exonuclease activity

pH optimum: Approximately 9.0 (+20°C)

Activation temperature: Active at ≥+60°C

Temperature optimum: Approximately +75°C

Half life at +95°C: Approximately 40 minutes

Substrates: Incorporates dNTP and various labeled or modified nucleotides (200 μmol/l each is recommended of normal dNTP, increased concentrations of variants).

Divalent ion requirement: Mg²⁺ (1.5 mmol/l standard concentration)

Specification


Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50.0% (v/v); pH approximately 8.0 at +4°C

Cat. No. **Pack Size**


05 457 882 103 custom fill

Will be supplied as "AptaTaq DNA Polymerase, 5 U/μl". Unit of Measure is "kU".

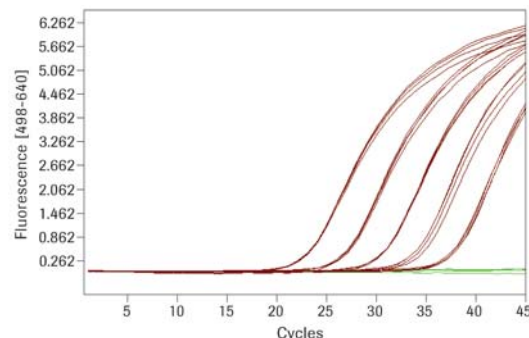

Sensitivity of AptaTaq DNA Polymerase, 5 U/μl, on a Real-Time PCR Instrument.

Various amounts of plasmid DNA (5000 fg to 0.5 fg) were used for the amplification of a Factor V wild-type fragment using HybProbe probe format. Even 5 fg can be detected without difficulties.

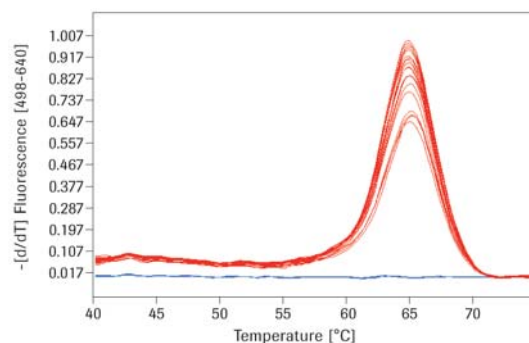

Specificity of AptaTaq DNA Polymerase, 5 U/μl, on a Real-Time PCR Instrument.

Melting curve analysis of Factor V wild-type fragments amplified from plasmid DNA (5000 fg to 0.5 fg) using HybProbe probe format results in a sharp T_m peak at about 65°C.

For further processing only.

Patent and License Disclaimer(s): 63

DNA Polymerases, Hot Start

AptaTaq DNA Polymerase

Volume activity: 5.5±0.5 U/μl

Aptamer concentration (HPLC): 3.6 μmol/l ±10%

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +36°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 30 U after 4 hours incubation at +65°C.

Performance test in qPCR using LightCycler® 480 (≥0.03 ng human genomic DNA, 339 bp tPA fragment): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

The aptamer/polymerase mixture is a hot start system with reversible inhibition of the polymerase activity at lower temperatures. Polymerase inactivation is achieved by a tight bond of the folded aptamer-oligonucleotide to the active site of the polymerase at lower temperatures. Upon heating above +60°C, the aptamer acts like a molecular switch, changing its temperature-dependent tertiary structure and releasing the active polymerase. Dropping the temperature below +55°C shuts off the polymerase activity again. Similar to antibody-based methods, the enzyme is much more quickly activated by heating, than chemically modified polymerases. In contrast to antibodies, the aptamer-oligonucleotide is much more stable, allowing longer storage at room temperature.

AptaTaq DNA Polymerase, 50 U/μl from *Thermus aquaticus* BM, expressed in *E. coli*, glycerol-free solution

Novel AptaTaq hot start PCR technology preserves the native speed and robustness of Taq DNA Polymerase in combination with a fast hot start system to ensuring sensitive and specific target DNA amplification. High concentration for dried-down amplification mixes.

Application

Apply AptaTaq DNA Polymerase for:

- Fast PCR assays with no extra enzyme activation time and fast cycling protocols
- Single- or multiplex PCR and qPCR applications requiring high specificity, sensitivity, and yield
- RT-PCR
- Difficult templates with secondary structures or GC-rich sequences
- Formulation of dried-down amplification reagents


Benefits

- **Reduce time to result.** Save up to 15 minutes per run by omitting the initial activation step required by chemically modified hot start polymerases, and reduce cycling time with fast protocols.
- **Maximize specificity, sensitivity, and yield.** Achieve reliable amplification of your target DNA from various sources (*e.g.*, genomic DNA, cDNA, plasmids).
- **Simplify PCR setup.** Store these highly stable polymerase for up to 1 month at +2° to +8°C and setup your hot start PCR reaction at room temperature.
- **Obtain consistent results.** Roche standardized manufacturing processes include extensive Quality Control release testing for high lot-to-lot consistency ideal for (IVD) kit manufacturers and end users.
- **Prepare stable amplification mixes in dry format.** Use this formulation for producing dried-down amplification mixes stable at room temperature.

Cat. No. **Pack Size**

05 187 605 103 custom fill

Will be supplied as "AptaTaq DNA Pol., Glycerol-free, 50 U/μl". Unit of Measure is "kU".

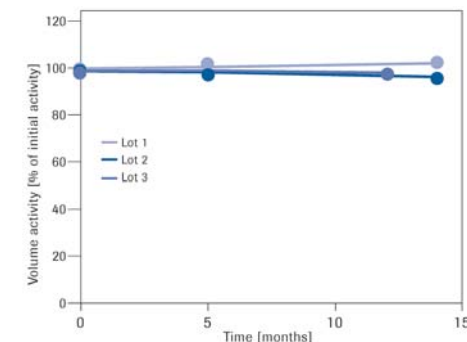

Real-time stability of AptaTaq DNA Polymerase, glycerol-free, 50 U/μl.

Volume activity was determined by a radioactive test after storage at -20°C for different time periods, starting at 100%. AptaTaq DNA Polymerase can be stored in the freezer for more than 15 months without any activity loss.

Product Description

AptaTaq DNA Polymerase is a blend of Taq DNA Polymerase and a specific oligonucleotide (aptamer) providing hot start features. The concentrated formulation does not contain glycerol and is suitable for the preparation of dry amplification mix preparations.

EC 2.7.7.7

Properties

AptaTaq DNA Polymerase is active at temperature above +60 to +65°C and inactive below +55°C. This hot start feature eliminates the risk of nonspecific primer extension. Taq DNA Polymerase is a highly processive 5'-3' DNA Polymerase that lacks 3'-5' exonuclease activity. Taq DNA Polymerase is stable during prolonged incubations at elevated temperatures (+95°C). The enzyme exhibits highest activity at a pH of approximately 9 (adjusted at +20°C) and temperatures approximately +75°C. The inherent stability of Taq DNA Polymerase is shown by the high storage stability in refrigerator and freezer (24 months at +2 to +8°C and -25 to -25°C). Taq DNA Polymerase accepts dNTP analogs as substrates.

pH optimum: Approximately 9.0 (+20°C)

Temperature optimum for elongation: Approximately +75°C

Half life at +95°C: Approximately 40 minutes

Divalent ion requirement: Mg²⁺ (standard concentration, 1.5 mmol/l)

dNTP requirement: Approximately 200 µmol/l for each dNTP

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); pH approximately 8.0 at +4°C

Volume activity: 55±5 U/µl

Glycerol content: ≤0.1% (v/v)

Aptamer concentration (HPLC): 35.75 µmol/l ±10%

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 30 U after 4 hours incubation at +65°C.

Performance test in qPCR using LightCycler® 480 (≥0.03 ng human genomic DNA, 339 bp tPA fragment): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

See *AptaTaq DNA Polymerase, 5 U/µl*


AptaTaq DNA Polymerase LDx, 5 U/µl from *Thermus aquaticus* BM, expressed in *E. coli*, solution

Taq DNA Polymerase with novel hot start system and extremely low DNA background for maximum sensitivity, speed and robustness in liquid assay formulations.

Application

Select AptaTaq DNA Polymerase LDx to perform microbial testing and other assays where the absence of contaminating bacterial, fungal, and/or human DNA is crucial. AptaTaq DNA LDx Polymerase is ideal for:

- Fast PCR assays with no extra enzyme activation time and fast cycling protocols

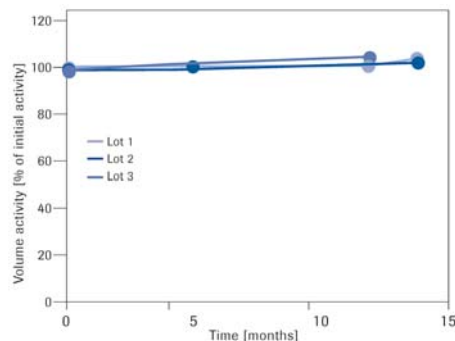

Real-time stability of AptaTaq DNA Polymerase, glycerol-free, 50 U/µl.

Volume activity was determined by a radioactive test after storage at +4°C for different time periods, starting at 100%. AptaTaq DNA Polymerase can be stored in the refrigerator for more than 15 months without any activity loss.

For further processing only.

Patent and License Disclaimer(s): 63

Cat. No.	Pack Size
05 884 314 103	custom fill

Will be supplied as "AptaTaq DNA Polymerase LDx, 5 U/µl". Unit of Measure is "kU".

For further processing only.

Patent and License Disclaimer(s): 63

- Single- and multiplex PCR and qPCR applications that require high specificity, sensitivity, and yield
- RT-PCR
- Difficult templates with secondary structures or GC-rich sequences
- Automated PCR workflows requiring high stability of the reaction mixtures during automated pipetting and prolonged handling at room temperature

Benefits

- **Minimize risks from contaminating nucleic acids.** AptaTaq DNA Polymerase LDx is extensively tested using ultra sensitive tests for contaminating nucleic acids from bacteria and fungi.
- **Rely on sophisticated manufacturing practices.** Roche has developed a nucleic acid-free workflow with clearly defined, highly consistent manufacturing processes to offer a product with very low nucleic acid background.
- **Enjoy the benefits of the advanced AptaTaq hot start system.** Use AptaTaq DNA Polymerase for additional benefits including speed, easy handling and consistent results.

Product Description

AptaTaq DNA Polymerase LDx is a blend of Taq DNA Polymerase and a specific oligonucleotide (aptamer) with hot start features, optimized for applications detecting lowest levels of DNA.

EC 2.7.7.7

Properties

AptaTaq DNA Polymerase LDx is active at temperature above +60 to +65°C and inactive below +55°C. This hotstart feature eliminates the risk of nonspecific primer extension. Taq DNA Polymerase is a highly processive 5'-3' DNA Polymerase that lacks 3'-5' exonuclease activity. Taq DNA Polymerase is stable during prolonged incubations at elevated temperatures (+95°C). The enzyme exhibits highest activity at a pH of approximately 9 (adjusted at +20°C) and temperatures approximately +75°C. Taq DNA Polymerase accepts dNTP analogs as substrates.

pH optimum: Approximately 9.0 (+20°C)

Temperature optimum for elongation: Approximately +75°C

Half life at +95°C: Approximately 40 minutes

Divalent ion requirement: Mg²⁺ (standard concentration, 1.5 mmol/l)

dNTP requirement: Approximately 200 µmol/l for each dNTP

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50% (v/v); pH approximately 8.0 at +4°C

Volume activity: 5.5±0.5 U/µl

Aptamer concentration (HPLC): 3.6 µmol/l ±10%

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 30 U after 4 hours incubation at +65°C.

Tests for the absence of contaminating nucleic acids

(human genomic DNA, β-Globin fragment): Corresponds to specification (LightCycler® UniTOOL ResoLight assay, detecting grampositive and gramnegative bacterial

DNA and fungal DNA, <1 genome equivalent/20 U enzyme): Corresponds to specification

246 Performance test in qPCR using LightCycler® 480 (≥0.03 ng human

genomic DNA, 339 bp tPA fragment): Corresponds to reference
Stability: At -15 to -25°C within specification range for 12 months.

Quality

AptaTaq DNA Polymerase LDx quality contains a very low DNA background, verified using an ultra-sensitive LightCycler® assay for the absence of gram(+), gram(-) bacteria, and fungal DNA. To pass this Quality Control test, the level of contaminating nucleic acid must be <1 genome equivalent per 20 units of DNA polymerase. Furthermore, AptaTaq DNA Polymerase LDx is analyzed for the absence of contaminating human DNA using a LightCycler® test specific for β-globin.

Background Information

Contaminating nucleic acids from various sources can affect PCR due to non-specific amplification, leading to reduced sensitivity and specificity, and false positive results. To minimize the risk of contamination and provide a product with very low nucleic acid background, Roche developed a nucleic acid-free workflow with defined, consistent manufacturing processes and ultra-sensitive quality control methods:

- Our raw materials have reduced DNA content.
- All equipment, buffers, and solutions are decontaminated.
- Highly trained staff and dedicated rooms ensure clean production.
- Remaining traces of DNA contamination are removed using chromatography.
- The final product is extensively characterized and tested for the absence of contaminating DNA.

For additional information on the AptaTaq hot start system, see **AptaTaq DNA Polymerase, 5 U/μl**

AptaTaq DNA Polymerase LDx, 50 U/μl from *Thermus aquaticus* BM, expressed in *E. coli*, glycerol-free solution

Taq DNA Polymerase with novel hot start system and extremely low DNA background for maximum sensitivity, speed and robustness. High concentration for dried-down amplification mixes and processing.

Application

Select AptaTaq DNA Polymerase LDx to perform microbial testing and other assays where the absence of contaminating bacterial, fungal, and/or human DNA is crucial. AptaTaq DNA LDx Polymerase is ideal for:


- Fast PCR assays with no extra enzyme activation time and fast cycling protocols
- Single- and multiplex PCR and qPCR applications that require high specificity, sensitivity, and yield
- RT-PCR
- Difficult templates with complex secondary structures or GC-rich sequences
- Formulation of dried-down amplification reagents

Benefits

- **Minimize risk of contaminating nucleic acids.** AptaTaq DNA Polymerase LDx is extensively evaluated using ultra sensitive tests for detecting contaminating nucleic acids from bacteria and fungi.
- **Rely on high quality manufacturing practices.** Roche has developed a nucleic acid-free workflow with clearly defined, highly consistent manufacturing processes resulting in a product with very low nucleic acid background.

Cat. No.	Pack Size
05 447 895 103	custom fill

Will be supplied as "AptaTaq DNA Pol. LDx, Glyc.-free, 50 U/μl". Unit of Measure is "kU".

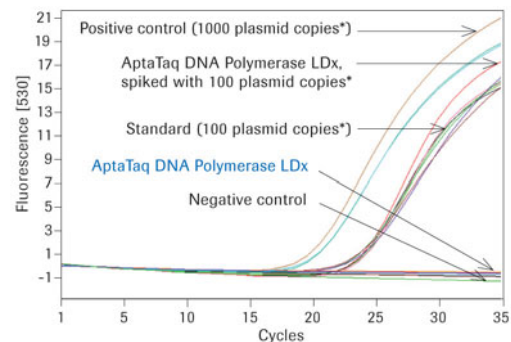

LightCycler® UniTool ResoLight quality control release assay for AptaTaq DNA Polymerase LDx on a Real-Time PCR Instrument.

Test of thirty units of AptaTaq DNA Polymerase LDx, glycerol-free, 50 U/μl shows no contaminating gram(+) or gram(-) bacterial DNA or fungal DNA. The Quality Control release value is defined as <1 genome equivalent/20 units DNA polymerase.

For further processing only.

Patent and License Disclaimer(s): 63

DNA Polymerases, Hot Start

AptaTaq DNA Polymerase

- **Prepare stable amplification mixes in dry format.** Use this formulation for producing dried-down amplification mixes stable at room temperature.
- **Enjoy the benefits of the advanced AptaTaq hot start system.** Refer to AptaTaq DNA Polymerase for additional benefits like speed, easy handling and consistent results.

Product Description

AptaTaq DNA Polymerase LDx is a blend of Taq DNA Polymerase and a specific oligonucleotide (aptamer) with hot start features, optimized for applications detecting the lowest levels of DNA.

EC 2.7.7.7

Properties

AptaTaq DNA Polymerase LDx is active at temperature above +60 to +65°C and inactive below +55°C. This hot start feature eliminates the risk of nonspecific primer extension. Taq DNA Polymerase is a highly processive 5'-3' DNA Polymerase lacking 3'-5' exonuclease activity. Taq DNA Polymerase is stable during prolonged incubations at elevated temperatures (+95°C). This enzyme exhibits highest activity at a pH of approximately 9 (adjusted at +20°C) and temperatures approximately +75°C. Taq DNA Polymerase accepts dNTP analogs as substrates.

pH optimum: Approximately 9.0 (+20°C)

Temperature optimum for elongation: Approximately +75°C

Half life at +95°C: Approximately 40 minutes

Divalent ion requirement: Mg²⁺ (standard concentration, 1.5 mmol/l)

dNTP requirement: Approximately 200 µmol/l for each dNTP

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); pH approximately 8.0 at +4°C

Volume activity: 55±5 U/µl

Glycerol content: ≤0.1% (v/v)

Aptamer concentration (HPLC): 35.75 µmol/l ±10%

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 30 U after 4 hours incubation at +37°C.

Tests for the presence of contaminating nucleic acids

(human genomic DNA, β-Globin fragment): Corresponds to specification (LC UniTool Resolight assay, specific for grampositive and gramnegative bacterial DNA and fungi DNA, <1 genome equivalent/20 U enzyme): Corresponds to specification

Performance test in qPCR using LightCycler® 480 (≥0.03 ng human genomic DNA, 339 bp tPA fragment): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Quality

LDx quality contains a very low DNA background, as verified using an ultra-sensitive LightCycler® assay for the absence of gram(+), gram(-) bacteria, and fungal DNA. To pass this Quality Control test, the level of contaminating nucleic acid must be <1 genome equivalent per 20 units of DNA polymerase. Furthermore, it is analyzed for the absence of contaminating human DNA with a LightCycler® test, specific for β-globin.

Background Information

For information on LDx refer to *AptaTaq DNA Polymerase LDx, 5 U/μl*
For additional information on the AptaTaq hot start system, see *AptaTaq DNA Polymerase, 5 U/μl*

AptaTaq Δexo DNA Polymerase, 5 U/μl from *Thermus aquaticus* BM, expressed in *E. coli*, solution

N-terminal truncated Taq DNA Polymerase with novel hot start system and no 5'-3' exonuclease activity for optimal detection of mismatches, speed and robustness in liquid assay formulations.

Application

Use AptaTaq Δexo DNA Polymerase for:

- SNP analysis and genotyping
- Allele-specific PCR
- Multiplexing
- Arbitrarily primed PCR
- Automated PCR requiring prolonged handling at room temperature

When time to result matters, this novel hot start technology is ideal as it does not require any activation time.

Benefits

- **Optimize your SNP analysis.** Discriminate between paired and unpaired primer ends using an enzyme optimized for allele-specific PCR.
- **Obtain reliable results fast.** Benefit from the general features of the AptaTaq DNA Polymerase System with the differentiating capabilities of a 5'-3' exonuclease activity-lacking Taq DNA Polymerase.

Product Description

This novel optimized mixture of high-quality N-terminal-deleted Taq DNA Polymerase and a specific oligonucleotide (aptamer) provides improved discrimination against misextension. As with the AptaTaq DNA Polymerase System, the AptaTaq Δexo DNA Polymerase-based assay shows high specificity and a broad dynamic range of products.

EC 2.7.7.7

Properties

AptaTaq Δexo DNA Polymerase is active at temperature above +60 to +65°C and inactive below +55°C. This hotstart feature eliminates the risk of nonspecific primer extension. Taq DNA Polymerase is a highly processive 5'-3' DNA Polymerase that lacks 5'-3' and 3'-5' exonuclease activity. Taq DNA Polymerase is stable during prolonged incubations at elevated temperatures (+95°C). The enzyme exhibits highest activity at a pH of approximately 9 (adjusted at +20°C) and temperatures approximately +75°C. Taq DNA Polymerase accepts dNTP analogs as substrates.

pH optimum: Approximately 8.3 (+20°C)

Temperature optimum for elongation: Approximately +72°C

Half life at +95°C: Approximately 40 minutes

Divalent ion requirement: Mg²⁺ (standard concentration, 2 mmol/l)

dNTP requirement: Approximately 200 μmol/l for each dNTP

Specification


Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50% (v/v); pH approximately 8.0 at +4°C

Volume activity: 5.5±0.5 U/μl


Cat. No.	Pack Size
05 458 030 103	custom fill

Will be supplied as "AptaTaq delta exo DNA Polymerase, 5 U/μl".
Unit of Measure is "kU".

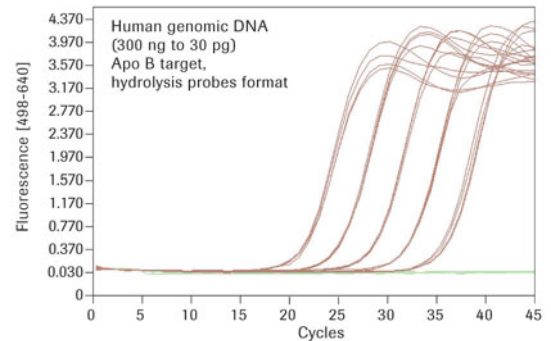

Amplification of Apo B: Sensitivity of AptaTaq Δexo DNA Polymerase, 5 U/μl, on a Real-Time PCR Instrument.

Thermal cycling conditions: Denaturation: 30 seconds at 95°C. Amplification: 5 seconds at 95°C, 15 seconds at 60°C, 10 seconds at 72°C, 45 cycles. Cooling: 60 seconds at 40°C.

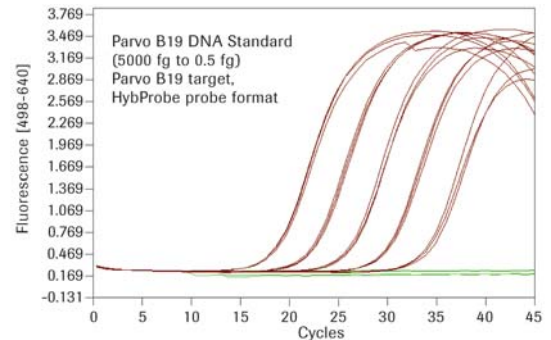

Amplification of Parvo B19: Sensitivity of AptaTaq Δexo DNA Polymerase, 5 U/μl, on a Real-Time PCR Instrument.

Thermal cycling conditions: Denaturation: 30 seconds at 95°C. Amplification: 5 seconds at 95°C, 15 seconds at 60°C, 10 seconds at 72°C, 45 cycles. Cooling: 60 seconds at 40°C.

For further processing only.

Patent and License Disclaimer(s): 63

DNA Polymerases, Hot Start

AptaTaq DNA Polymerase

Aptamer concentration (HPLC): 24.0 $\mu\text{mol/l} \pm 10\%$

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Exonucleases ($^3\text{H-DNA}$): Not detectable in up to 30 U after 4 hours incubation at +65°C.

Performance test in qPCR using LightCycler® 480 (≥ 0.03 ng human genomic DNA, 339 bp tPA fragment): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

See *AptaTaq DNA Polymerase, 5 U/ μl*

AptaTaq Δexo DNA Polymerase, 50 U/ μl from *Thermus aquaticus* BM, expressed in *E. coli*, glycerol-free solution

N-terminal truncated Taq DNA Polymerase with novel hot start system, and no 5'-3' exonuclease activity, for optimal detection of mismatches, speed and robustness. High concentration for dried-down amplification mixes and processing.

Application

Use AptaTaq Δexo DNA Polymerase for:

- SNP analysis and genotyping
- Allele-specific PCR
- Multiplexing
- Arbitrarily primed PCR
- Formulation of dried-down amplification reagents

When time to result matters, this novel hot start technology is ideal as it does not require any activation time.

Benefits

- **Optimize your SNP analysis.** Discriminate between paired and unpaired primer ends using an enzyme optimized for allele-specific PCR.
- **Obtain reliable results fast.** Benefit from the general features of the AptaTaq DNA Polymerase System with the discriminating capabilities of a 5'-3' exonuclease activity-lacking Taq DNA Polymerase.
- **Prepare stable amplification mixes in dried-down format.** Use this formulation for producing dried-down amplification mixes stable at room temperature.

Product Description

AptaTaq DNA Δexo Polymerase is a blend of N-terminal truncated Taq DNA Polymerase and a specific oligonucleotide (aptamer) with hot start features, optimized for excellent discrimination against misextension. This concentrated formulation contains no glycerol and is suitable for the preparation of dried-down amplification mix preparations.

EC 2.7.7.7

Properties

AptaTaq Δexo DNA Polymerase is active at temperatures above +60 to +65°C and inactive below +55°C. This hot start feature eliminates the risk of non-specific primer extension. Taq DNA Polymerase itself is a highly processive 5'-3' DNA Polymerase lacking 5'-3' and 3'-5' exonuclease activity. Taq DNA Polymerase is stable during prolonged incubations at elevated temperatures (+95°C). The enzyme exhibits highest activity at a pH of approximately

Cat. No.	Pack Size
05 364 086 103	custom fill

Will be supplied as "AptaTaq delta exo DNA Pol., Glyc.-free". Unit of Measure is "kU".



For further processing only.

Patent and License Disclaimer(s): 63

9 (adjusted at +20°C) and temperatures approximately +75°C. Taq DNA Polymerase accepts dNTP analogs as substrates.

pH optimum: Approximately 8.3 (+20°C)

Temperature optimum for elongation: Approximately +72°C

Half life at +95°C: Approximately 40 minutes

Divalent ion requirement: Mg²⁺ (standard concentration, 2 mmol/l)

dNTP requirement: Approximately 200 µmol/l for each dNTP

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Casein 1 g/l; glycerol-free; pH approximately 8.0 at +4°C

Volume activity: 55±5 U/µl

Aptamer concentration: 240 µmol/l ±10% (HPLC)

Unspecific endonucleases (λDNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 30 U after 16 hours incubation at +37°C.

Function test in qPCR using LightCycler® 480 (≥0.03 ng human genomic DNA, 339 bp tPA fragment): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

See *AptaTaq DNA Polymerase, 5 U/µl*

EagleTaq DNA Polymerase, 5 U/ μ L from *Thermus aquaticus*, expressed in *E. coli*, solution

EagleTaq DNA Polymerase is a chemically modified hot start DNA polymerase.

Application

Apply EagleTaq DNA Polymerase for:

- Hot start activated amplification
- Incorporation of modified nucleotides for generating labeled PCR products
- Detection formats such as hydrolysis probes, hybridization probes and SYBR Green

Benefits

- Obtain high specificity, sensitivity, and yield for genomic targets.
- Simplify PCR setup for automation by using hot start activation at +95°C with a chemically modified enzyme.
- Rely on lot-to-lot consistency and full traceability with GMP manufacturing.
- Control carryover contamination with a dUTP and Uracil-DNA Glycosylase treatment (UNG) compatible setup.

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Tween 20, 0.5% (v/v); glycerol, 50.0% (v/v); pH approximately 9.0 at +20°C

Volume activity: 5.4–5.9 U/ μ l

Function test: At least 1.5×10^5 fold amplification of λ DNA after 25 cycles. One band at 500 bp on an agarose gel.

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
05 206 944 190	1 kU
05 206 952 190	25 kU
05206944190: Will be supplied as "CMPNT EAGLETAQ 1 KU, 5U/ μ L 0.2mL". Unit of Measure is "piece".	
05206952190: Will be supplied as "CMPNT EAGLETAQ 25 KU, 5U/ μ L 5mL". Unit of Measure is "piece".	
For further processing only.	
Patent and License Disclaimer(s): 64	

FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl

from *Thermus aquaticus* BM, expressed in *E. coli*, solution

FastStart Taq DNA Polymerase enables specific and sensitive amplification of DNA fragments in a hot start polymerase chain reaction (PCR) with activation of the enzyme in the initial heat denaturation phase.

Application

Use FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl, for:

- Hot start PCR and RT-PCR with high specificity, sensitivity and yield
- Specific amplification of DNA fragments from various sources of DNA and for diverse down-stream applications
- Labeling of DNA with modified nucleotides (e.g., DIG-dUTP, biotin-dUTP, fluorescein-dUTP)
- The prevention of carryover contamination between PCR reactions in combination with dUTP and Uracil-DNA Glycosylase
- Manufacture of amplification mixtures for regulated applications (e.g., *in vitro* diagnostics, quality control) with requests for more stringent validation

Benefits

- **Maximize specificity and target yield.** Minimize the extension of non-specifically bound primers using this hot start reaction that maximizes the amplification of specific product.
- **Achieve highest sensitivity.** Prevent nonspecific priming and detect amplification product from as little as one copy of your target DNA.
- **Simplify PCR setup.** Pipette and handle the hot start reaction mix at ambient temperature.
- **Simplify assay design.** Create robust PCR assays with a minimum of optimization efforts also suited for multiplex PCR applications.
- **Obtain reliable results.** Rely on the robust reaction performance, and high lot-to-lot consistency of this product, thoroughly tested for a reproducible quality. Manufacturing and documentation are according to GMP (Good Manufacturing Practice) regulations.

EC 2.7.7.7

Properties

FastStart Taq DNA Polymerase is designed for hot start PCR and has to be heat-activated in the beginning of the reaction protocol.

Enzyme activities: Highly processive 5'-3' DNA polymerase; double-strand specific 5'-3' exonuclease; no 3'-5' exonuclease activity

Heat activation: +95°C for 3-10 minutes (assay-dependent; recommendation is 10 minutes for full activation)

pH optimum: Approximately 9.0 (+25°C)

Temperature optimum: Approximately +75°C

Half life at +95°C: Approximately 40 minutes

Substrates: Incorporates dNTP, dUPT, dITP, various labeled or modified nucleotides (200 μmol/l each is recommended of normal dNTP, increased concentrations of variants).

Divalent ion requirement: Mg²⁺ (1.5 mmol/l standard concentration)

Specification

Appearance: Clear, opalescent solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Tween 20, 0.2% (v/v); glycerol, 50% (v/v); pH 9.0±0.1 at +25°C

Volume activity: ≥5 U/μl

Unit definition: One unit Taq DNA Polymerase is defined as the amount of heat-activated enzyme that incorporates 10 nmol of total deoxyribonucleoside-triphosphates into acid precipitable DNA within 30 minutes at +75°C under

Cat. No. Pack Size

04 659 163 103 5 kU

Will be supplied as "FastStart Taq DNA Pol. Ind. GMP Grd, 5KU".
Unit of Measure is "piece".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 49

- ▶ For the best fit reaction buffer, use PCR Buffer, see page 256
- ▶ For the best fit reaction buffer, use FastStart PCR Buffer, see page 238

DNA Polymerases, Hot Start

FastStart DNA Polymerase

standard assay conditions.

Unspecific endonucleases (λ DN): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 25 U after 1 hour incubation at +37°C.

Function test in PCR

(human genomic DNA, tPA gene): Corresponds to reference

(human genomic DNA, ApoE gene): Corresponds to reference

Function test in qPCR using the LightCycler® System

(human genomic DNA, β -globin gene): Corresponds to reference

(plasmid DNA, β -globin gene): Corresponds to reference

(reverse transcribed cDNA, PBGD gene): Corresponds to reference

Bioburden: ≤ 50 CFU/ml

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Background Information

FastStart Taq DNA Polymerase is a chemically inactivated form of recombinant Taq DNA Polymerase. It remains inactive at temperatures up to +75°C. At higher temperatures, the modification is cleaved off and the polymerase acquires its enzymatic activity. Using FastStart Taq DNA Polymerase, PCR setup can be done conveniently at ambient temperature with no risk of nonspecific priming. The polymerase will not be activated until the initial denaturation step of the PCR protocol, at which point nonspecific hybridization can no longer occur.

FastStart Taq DNA Polymerase, 5 U/ μ l from *Thermus aquaticus* BM, expressed in *E. coli*, solution

FastStart Taq DNA Polymerase enables specific and sensitive amplification of DNA fragments in a hot start polymerase chain reaction (PCR) with activation of the enzyme in the initial heat denaturation phase.

Application

For applications see **FastStart Taq DNA Polymerase, GMP Grade, 5 U/ μ l**

Benefits

See benefits of **FastStart Taq DNA Polymerase, GMP Grade, 5 U/ μ l**

EC 2.7.7.7

Properties

See **FastStart Taq DNA Polymerase, GMP Grade, 5 U/ μ l**

Specification

Appearance: Clear, opalescent solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Tween 20; 0.2% (v/v); glycerol, 50% (v/v); pH approximately 9.0 at +25°C

Volume activity: ≥ 5 U/ μ l

Unit definition: One unit Taq DNA Polymerase is defined as the amount of heat-activated enzyme that incorporates 10 nmol of total deoxyribonucleoside-triphosphates into acid precipitable DNA within 30 minutes at +75°C under standard assay conditions.

Unspecific endonucleases (λ DN): Not detectable in up to 25 U after 16

Cat. No.	Pack Size
12 161 508 103	custom fill

Will be supplied as "Fast Start Taq DNA Polymerase". Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 49

- For the best fit reaction buffer, use PCR Buffer, see page 256
- For the best fit reaction buffer, use FastStart PCR Buffer, see page 238

hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 25 U after 1 hour incubation at +37°C.

Exonucleases (calf thymus DNA): Not detectable in up to 15 U after 4 hours incubation at +65°C.

Function test in PCR

(50 pg human genomic DNA, 365 bp tPA fragment): Corresponds to reference
(200 ng human genomic DNA, 284 bp ApoE fragment): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

See **FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl**

FastStart Taq DNA Polymerase, 100 U/μl
from *Thermus aquaticus* BM, expressed in *E. coli*,
solution

FastStart Taq DNA Polymerase enables specific and sensitive amplification of DNA fragments in a hot start polymerase chain reaction (PCR) with activation of the enzyme in the initial heat denaturation phase.

Application

Use FastStart Taq DNA Polymerase, 100 U/μl, especially for:

- Setup of PCR master mixtures, when highly concentrated components are requested
- Preparation of stabilized dried-down formulations of reaction mixtures

For further applications, see **FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl**

Benefits

See benefits of **FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl**

EC 2.7.7.7

Properties

See **FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl**

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; Tween 20, 0.2% (v/v); glycerol, 50% (v/v); pH 9.0 at +25°C ± 0.1

Volume activity: ≥ 100 U/μl

Unit definition: One unit Taq DNA Polymerase is defined as the amount of heat-activated enzyme that incorporates 10 nmol of total deoxyribonucleoside-triphosphates into acid precipitable DNA within 30 minutes at +75°C under standard assay conditions.

Unspecific endonucleases (λDNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 25 U after 1 hour incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 15 U after 4 hours incubation at +65°C.

Cat. No. **Pack Size**

04 433 785 103 custom fill

Will be supplied as "FastStart Taq DNA Pol. 100 U/μl". Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Patent and License Disclaimer(s): 49

DNA Polymerases, Hot Start

FastStart DNA Polymerase

Function test in PCR

(human genomic DNA, tPA gene): Corresponds to reference

(human genomic DNA, ApoE gene): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

See *FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl*

FastStart PCR Buffer 10x conc., with 20 mM MgCl₂

Standard reaction buffer for PCR using FastStart Taq DNA Polymerase.

Application

Use this buffer together with FastStart Taq DNA Polymerase. For applications refer to *FastStart Taq DNA Polymerase, 5 kU, GMP Grade*.

Benefits

- **Simplify PCR setup.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer for fast and easy setup of highly reproducible PCR experiments.
- **Gain better performance.** Take full advantage of FastStart Taq DNA Polymerase using a reaction buffer specially optimized for this enzyme.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 500 mmol/l; (NH₄)₂SO₄, 50 mmol/l; KCl, 100 mmol/l; MgCl₂, 20 mmol/l; pH approximately 8.3 at +25°C

Unspecific endonucleases (λDNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 20 μl after 1 hour incubation at +37°C.

Function test in PCR (50 pg human genomic DNA, 365 bp tPA fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
12 161 567 103	1 ml
12 161 516 103	custom fill
12161516103:	Will be supplied as "PCR buffer (10X) w MgCl ₂ ". Unit of Measure is "ml".
12161567103:	Will be supplied as "FastStart PCR Buffer (10x), 1 ml". Unit of Measure is "piece".
	
For further processing only.	



FastStart PCR Buffer 10x conc., without MgCl₂

Standard reaction buffer without MgCl₂ for optimization of the MgCl₂ concentration in PCR using FastStart Taq DNA Polymerase.

Application

Use this buffer together with FastStart Taq DNA Polymerase whenever the amplification of difficult target requires a specific MgCl₂ concentration.

Benefits

- **Amplify difficult targets.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer with an optimized MgCl₂ concentration for best results.
- **Gain excellent performance.** Take full advantage of FastStart Taq DNA Polymerase using a reaction buffer specially optimized for this enzyme.

Specification

256 Appearance: Clear, colorless solution

Cat. No.	Pack Size
05 917 166 103	1 ml
12 161 494 103	custom fill
12161559103:	Will be supplied as "PCR buffer (10xconc.) without MgCl ₂ MPB". Unit of Measure is "piece".
12161494103:	Will be supplied as "PCR buffer (10X) w/o MgCl ₂ ". Unit of Measure is "ml".
	
For further processing only.	



Contents: Tris/HCl, 500 mmol/l; $(\text{NH}_4)_2\text{SO}_4$, 50 mmol/l; KCl, 100 mmol/l; pH approximately 8.3 at +25°C

Unspecific endonucleases (λ DNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μl after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 20 μl after 1 hour incubation at +37°C.

Function test in PCR (50 pg human genomic DNA, 365 bp tPA fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

HawkZ05 DNA Polymerase, 40 U/μl from *Thermus* species Z05, expressed in *E. coli*, solution

The HawkZ05 DNA Polymerase is an one-step RT/PCR DNA polymerase.

Application

Apply HawkZ05 DNA Polymerase for:

- Efficient, high temperature cDNA synthesis and subsequent DNA amplification of RNA templates
- Use in multiplex PCR and qPCR applications that require high specificity, sensitivity, and yield
- Incorporation of modified nucleotides for labeling of PCR products
- Detection formats such as hydrolysis probes, hybridization probes and SYBR Green
- Fast-cycling diagnostic applications and other routine amplification of low-copy targets

Benefits

- Obtain increased detection sensitivity using this highly specific one- or two-step RT-PCR method.
- Enjoy fast PCR cycling without RNA degradation by using this fast activating, reversible aptamer-based hot start method.
- Control carryover contamination with a dUTP and Uracil-DNA Glycosylase treatment (UNG) compatible setup.
- Rely on lot-to-lot consistency and full traceability with GMP manufacturing.
- Simplify setup of PCR reactions using this stable aptamer hot start DNA polymerase at room temperature.

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Tween 20, 0.5% (v/v); glycerol, 50.0% (v/v); aptamer; pH approximately 8.0

Volume activity: 40-55 U/μl

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
05 230 322 190	5 kU
05 230 349 190	200 kU

05230322190: Will be supplied as "CMPNT HAWKZ05 DNA Polymerase 5KU 0.125mL". Unit of Measure is "piece".

05230349190: Will be supplied as "CMPNT HAWK Z05 DNA Polymerase 200KU 5mL". Unit of Measure is "piece".

For further processing only.

Patent and License Disclaimer(s): 63

HawkTaq DNA Polymerase, 5 U/ μ L from *Thermus aquaticus*, expressed in *E. coli*, solution

HawkTaq DNA Polymerase is a hot start DNA polymerase.

Application

Apply HawkTaq DNA Polymerase for:

- Fast-cycling diagnostic applications and other routine amplification of low-copy targets
- Use in multiplex PCR and qPCR applications that require high specificity, sensitivity, and yield
- Incorporation of modified nucleotides for labeling of PCR products
- Detection formats such as hydrolysis probes, hybridization probes and SYBR Green

Benefits

- Enjoy increased specificity, sensitivity, and yield compared to non-hot start methods.
- Utilize fast PCR cycling with the fast activating, reversible aptamer-based hot start.
- Control carryover contamination with a dUTP and Uracil-DNA Glycosylase treatment (UNG) compatible setup.
- Rely on high lot-to-lot consistency and full traceability with GMP manufacturing.
- Simplify setup of PCR reactions using this highly stable enzyme at room temperature.

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; EDTA, 0.1 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); glycerol, 50.0% (v/v); aptamer; pH approximately 8.0

Volume activity: 5.0-6.6 U/ μ l

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
05 230 357 190	5 kU
05 230 365 190	125 kU

05230357190: Will be supplied as "CMPNT HAWK TAQ DNA Polymerase 5KU 1mL". Unit of Measure is "piece".

05320365190: Will be supplied as "CMPNT HAWK TAQ DNA Polymerase 125KU 25mL". Unit of Measure is "piece".

For further processing only.

Patent and License Disclaimer(s): 63

AptaTaq DNA Master

5x concentrated

Ready-to-use, fast and robust 5x reaction mix with novel hot start technology for real-time and endpoint PCR, especially suited for high-throughput applications on a variety of system platforms and detection formats.

Application

Use the master mix to amplify targets efficiently with high specificity and sensitivity. The AptaTaq DNA Master is developed to match a broad range of applications. It is ideally suited for high-throughput applications with low reaction volume due to its 5x concentration and high stability at room temperature. In combination with appropriate dyes, it can be used on various instrument platforms for endpoint analysis and real-time PCR. Due to the use of dUTP, DNA carryover contamination can be prevented when adding Uracil-DNA Glycosylase (UNG).

Benefits

- **Use the PCR instrument of your choice.** Perform sensitive and efficient endpoint analysis and real-time PCR on a variety of platforms and formats.
- **Obtain results quickly.** Omit the enzyme activation and save up to 15 minutes per run.
- **Gain flexibility.** The 5x concentrated master mix allows you to vary reaction volumes and sample input for best results.
- **Benefit from high stability.** The master mix is stable during setup and in high-throughput instrument platform ready to be processed. The AptaTaq DNA Master can be kept in the refrigerator (+2 to +8°C) for at least 4 weeks without loss of activity and performance.

Product Description

AptaTaq DNA Master is a 5x concentrated, ready-to-use, one component hot start PCR mix. It contains AptaTaq DNA Polymerase, reaction buffer including an optimized Mg²⁺ concentration, and a dNTP mix with dUTP instead of dTTP.

EC 2.7.7.7

Properties

The master mix is very stable and can be stored in the refrigerator (+2 to +8°C) for at least 4 weeks without loss of activity and performance. It is stable at room temperature for at least 2 days.

Specification

Appearance: Clear, colorless solution

Performance test in qPCR using ABI 7500:

(human genomic DNA, CycA fragment): Corresponds to specification

(human genomic DNA, β-globin fragment): Corresponds to specification

(human genomic DNA, ApoE fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Background Information

For information on the AptaTaq hot start system, see **AptaTaq DNA**


Polymerase, 5 U/μl

Cat. No.

Pack Size


05 537 533 103 custom fill

Will be supplied as "AptaTaq DNA Master". Unit of Measure is "ml".

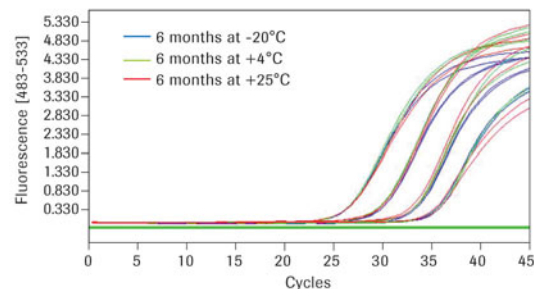

Storage stability of AptaTaq DNA Master.

Plasmid DNA from *Eurotium (Aspergillus) amstelodami* was detected in varying amounts (105 to 102 copies) using hydrolysis probes on a Real-Time PCR Instrument. AptaTaq DNA Master is stable for 6 months at +4°C and at least 4 weeks at +25°C.

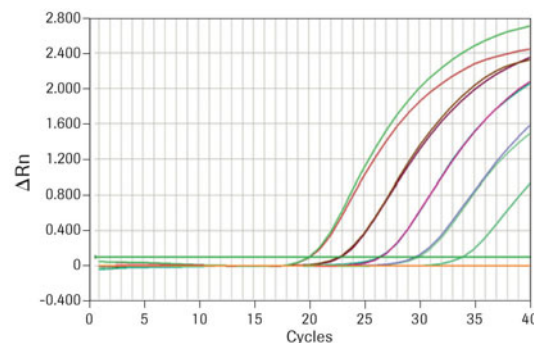

Performance of AptaTaq DNA Master with FRET-ROX on an ABI 7500 real-time PCR Instrument.

To demonstrate its sensitivity, the AptaTaq DNA Master has been tested in combination with hydrolysis probes detecting a specific human single-copy gene on an ABI 7500 real-time PCR Instrument. The result shows a broad dynamic range down to 5 pg (approximately 2 copies).

For further processing only.

Patent and License Disclaimer(s): 63

AptaTaq DNA Master without Mg²⁺ 5x concentrated

Ready-to-use, fast and robust 5x reaction mix without Mg²⁺, using novel hot start technology for real-time and endpoint PCR, especially suited for high-throughput applications on various system platforms and detection formats.

Application

Use this master mix to optimize Mg²⁺ concentration. Since it does not contain Mg²⁺, this master mix serves as basis for applications where a low Mg²⁺ concentration is required.

Benefits

- **Finetune AptaTaq DNA Master.** Use this master mix to optimize for best performance in your specific application.

Product Description

This master mix is identical to the AptaTaq DNA Master except that it does not contain Mg²⁺. The master mix is also part of the AptaTaq DNA Master Optimization Kit.

EC 2.7.7.7

Properties

The master mix is very stable and can be stored in the refrigerator (+2 to +8°C) for at least 4 weeks without loss of activity and performance. It is stable at room temperature for at least 2 days.

Specification

Appearance: Clear, colorless solution

Performance test in qPCR using ABI 7500:

(human genomic DNA, CycA fragment): Corresponds to specification

(human genomic DNA, β-globin fragment): Corresponds to specification

(human genomic DNA, ApoE fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Background Information

For information on the AptaTaq hot start system, see **AptaTaq DNA Polymerase, 5 U/μl**

AptaTaq DNA Master Optimization Kit

Kit for the adaptation of the ready-to-use, fast and robust 5x AptaTaq DNA Master for a variety of real-time PCR system platforms and detection formats.

Application

Use AptaTaq DNA Master Optimization Kit to design a master mix with best performance in your application, and to optimize reagents for your real-time PCR system platform and detection format. For high-throughput applications, a custom AptaTaq DNA Master can be provided. Please contact your sales representative.

Benefits

- **Optimize your master mix to your needs.**
- **Maximize versatility.** Use FRET-ROX dye on ABI instruments.
- **Generate higher fluorescence signals.** Detect dsDNA with high sensitivity using the next-generation LightCycler[®] ResoLight dye.

Product Description

The AptaTaq DNA Master Optimization Kit contains all reagents (except

Cat. No. **Pack Size**

05 548 802 103 custom fill

Will be supplied as "AptaTaq DNA Master w/o Mg²⁺". Unit of Measure is "ml".



For further processing only.

Patent and License Disclaimer(s): 63

Cat. No. **Pack Size**

05 537 568 001 1 Kit
for optimization of up to 750 reactions at
20 μl final reaction volume

Will be supplied as "AptaTaq DNA Master Optimization kit". Unit of Measure is "piece".



Contents

01. AptaTaq DNA Master
02. AptaTaq DNA Master without Mg²⁺, 5x concentrated
03. LightCycler[®] 480 ResoLight Dye
04. FRET-ROX Dye
05. MgCl₂ Stock Solution , 25 mmol/l
06. GC-RICH Solution, 5x concentrated
07. Water, PCR Grade

For life science research only. Not for use in diagnostic procedures.

Patent and License Disclaimer(s): 05

assay-specific components, such as primers and template) for optimizing PCR assays with AptaTaq DNA Polymerase. LightCycler® ResoLight dye allows the detection of dsDNA in real time using an appropriate PCR instrument (*e.g.*, the LightCycler® 480 Instrument), or in gel electrophoresis after the PCR has finished. A special Rox reference dye (FRET-ROX) enables you to design assays for all available real-time PCR instruments in which the Rox reference dye is required for quantitative analysis.

EC 2.7.7.7

Specification

Contents/Appearance:

Bottle 1: AptaTaq DNA Master, 5x concentrated; clear, colorless solution

Bottle 2: AptaTaq DNA Master without Mg²⁺, 5x concentrated; clear, colorless solution

Bottle 3: LightCycler® 480 ResoLight Dye; yellowish solution

Bottle 4: FRET-ROX Dye; slightly purple solution

Bottle 5: MgCl₂ Stock Solution, 25 mmol/l; clear, colorless solution

Bottle 6: GC-RICH Solution, 5x concentrated; clear, colorless solution

Bottle 7: Water, PCR Grade; clear, colorless solution

Performance test in qPCR using ABI 7500 (human genomic DNA, β-globin fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Background Information

For information on the AptaTaq hot start system, see **AptaTaq DNA Polymerase, 5 U/μl**

AptaTaq Genotyping Master 5x concentrated

Ready-to-use, fast and robust 5x reaction mix using novel hot start technology for high-throughput genotyping applications on real-time PCR instruments not requiring Rox normalization. Ideal PCR reagent mix for crude DNA extractions and robust against inhibitors.

Application

Use AptaTaq Genotyping Master in genotyping applications with all real-time PCR instruments that do not require Rox normalization. AptaTaq Genotyping Master is ideal for high-throughput applications using low reaction volumes. The master mix is very resistant to inhibitors and can be dried-down without loss of performance.


Benefits

- **Achieve results quickly.** Omit the enzyme activation and obtain results in less than half an hour.
- **Rely on a robust master mix.** Obtain reliable results from crude samples and even using low reaction volumes.
- **Ready for robotics.** Rely on the stability of the AptaTaq Genotyping Master mix for PCR automation. The viscosity of the master mix is optimized for accurate pipetting. The mix is stable during setup and on the stacker for more than 24 hours.
- **Gain flexibility.** The 5x concentrated master mix enables you to vary reaction volume and sample input for outstanding results. Use AptaTaq Genotyping Master mix for all real-time PCR instruments not requiring Rox normalization. For instruments requiring Rox normalization, use AptaTaq Genotyping Master (Rox).
- **Benefit from high stability.** Keep the master mix in the refrigerator for up to 4 weeks and profit from a quick setup without thawing first.

Cat. No.	Pack Size
05 955 807 103	10 ml for up to 5,000 reactions at 10 μl final reaction volume
05 890 152 103	custom fill

05955807103: Will be supplied as "AptaTaq Genotyping Master, 10ml". Unit of Measure is "piece".

05890152103: Will be supplied as "AptaTaq Genotyping Master". Unit of Measure is "ml".

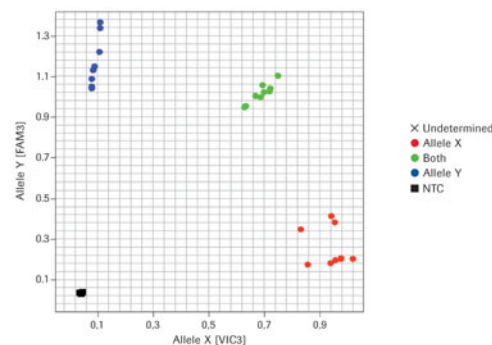

Clear allele separation by endpoint genotyping with hydrolysis probes using low reaction volumes on a Real-Time PCR Instrument.

Wild type and mutant TGFβ fragments were amplified using a 2-step protocol. Endpoint fluorescence was plotted in the VIC and FAM channels; reaction volume: 1 μl.

Product Description

AptaTaq DNA Master is a 5x concentrated, ready-to-use, one component hot start PCR mix, containing AptaTaq DNA Polymerase in an optimized concentration for the amplification of difficult sample types, reaction buffer, and a dNTP mix using dUTP instead of dTTP (for prevention of DNA contamination by PCR carryover by pretreatment with Uracil-DNA Glycosylase).

EC 2.7.7.7

Properties

The master mix is very stable and can be stored in the refrigerator (+2 to +8°C) for at least 4 weeks without loss of activity and performance. It is stable at room temperature for at least 2 days.

Specification

Appearance: Clear, colorless solution


Performance test in qPCR

(human genomic DNA, CypA fragment): Corresponds to specification
(human genomic DNA, β-globin fragment): Corresponds to specification
(human genomic DNA, ApoE fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.


Background Information

For information on the AptaTaq hot start system, see *AptaTaq DNA Polymerase, 5 U/μl*

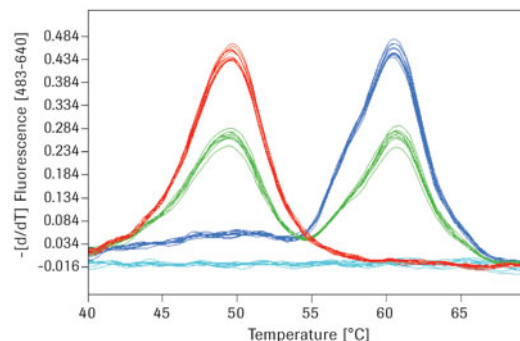

SNP detection with HybProbe probes using a Real-Time PCR Instrument.

Different SNPs of human Factor II were genotyped by melting curve analysis after amplification using a 3-step real-time PCR protocol.

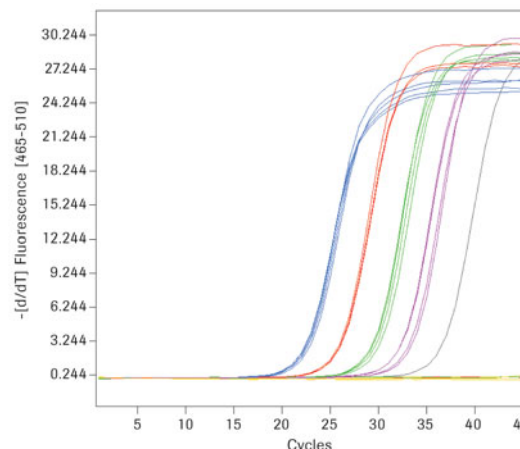

Robust amplification after DNA extraction that includes PCR inhibitors in crude plant material.

An *Arabidopsis* Gene E fragment was amplified by crude DNA extraction with chloroform and analyzed using a Real-Time PCR Instrument.

For further processing only.

Patent and License Disclaimer(s): 63

AptaTaq Genotyping Master (Rox)
5x concentrated

Ready-to-use, fast and robust 5x reaction mix with novel hot start technology for high-throughput genotyping applications on real-time PCR instruments requiring normalization with Rox. Ideal PCR reagent mix for crude DNA extractions and robust against inhibitors.

Application

Use AptaTaq Genotyping Master (Rox) in genotyping applications on instruments requiring normalization with Rox. AptaTaq Genotyping Master (Rox) is optimized for high-throughput applications using low reaction volumes. The master mix is very resistant to inhibitors and can be dried-down without loss of performance.

Benefits

- **Achieve results quickly.** Omit the enzyme activation and obtain results in less than half an hour.
- **Rely on a robust master mix.** Obtain reliable results from crude

Cat. No.	Pack Size
05 955 823 103	10 ml for up to 5,000 reactions at 10 μl final reaction volume
05 890 144 103	custom fill

05955823103: Will be supplied as "AptaTaq Genotyping Master (ROX), 10 ml". Unit of Measure is "piece".
05890144103: Will be supplied as "AptaTaq Genotyping Master (ROX)". Unit of Measure is "ml".
For further processing only.

Patent and License Disclaimer(s): 05, 63

samples even using low reaction volumes.

- **Ready for robotics.** Rely on the stability of the master mix for RT-PCR automation. The viscosity of AptaTaq Genotyping Master (Rox) is optimized for accurate pipetting. The mix is stable during setup and on the stacker for more than 24 hours.
- **Gain flexibility.** The 5x concentrated master mix allows you to vary reaction volume and sample input for best result. Use the same master mix for all types of real-time PCR instruments requiring Rox normalization.
- **Benefit from high stability.** Keep the master mix in the refrigerator for up to 4 weeks and profit from a quick setup without thawing first.

Product Description

AptaTaq DNA Master is a 5x concentrated, ready-to-use, one component hot start PCR mix, containing AptaTaq DNA Polymerase in an optimized concentration for the amplification of difficult sample types, reaction buffer, and a dNTP mix with dUTP (for prevention of DNA contamination by PCR carryover by pretreatment with Uracil-DNA Glycosylase). The special Rox reference dye (FRET-ROX) enables you to run assays for all real-time PCR instruments in which Rox reference dye is required for quantitative analysis.

EC 2.7.7.7

Properties

The PCR master mix is very stable and can be stored in the refrigerator (+2 to +8°C) for at least 4 weeks without loss of activity and performance. It is stable at room temperature for at least 2 days.

Specification

Appearance: Clear, slightly pink solution

Performance test in qPCR using ABI 7900

(human genomic DNA, CycA fragment): Corresponds to specification

(human genomic DNA, β-globin fragment): Corresponds to specification

(human genomic DNA, ApoE fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

For information on the AptaTaq hot start system, see *AptaTaq DNA Polymerase, 5 U/μl*

EagleTaq Master Mix

Eagle Taq Master Mix is a real-time PCR master mix.

Application

Apply EagleTaq Master Mix for:

- Gene expression analysis of cDNA
- Efficient amplification of rare cDNA and low copy DNA targets
- Fast thermal cycling for high-throughput real-time PCR applications

Benefits

- Enjoy efficient detection of rare targets.
- Rely on the stability of this master mix in automation and high-throughput plate stacking handling.
- Obtain high yield and a robust 5' nuclease based detection through use of the robust EagleTaq DNA Polymerase enzyme.
- Run the master mix on a variety of real-time PCR platforms that do not require Rox as reference dye.

Cat. No.	Pack Size
05 529 069 190	1 ml
05 529 085 190	50 ml

05529069190: Will be supplied as "KIT EAGLETAQ MMX 1mL RUO". Unit of Measure is "piece".

05529085190: Will be supplied as "KIT EAGLETAQ MMX 50mL". Unit of Measure is "piece".

For further processing only.

Patent and License Disclaimer(s): 64

- Minimize carryover PCR contamination with this 2x master mix containing dUTP, compatible with Uracil-DNA Glycosylase treatment (UNG).

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Function test: Average CT value of positive controls tested is between 18 and 28 cycles with starting template of 10 pg λDNA. Average CT value of test is within ±2 cycles of proven.

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

EagleTaq Master Mix (Rox)

Eagle Taq Master Mix (Rox) is a real-time PCR master mix for systems requiring a reference dye.

Application

Apply EagleTaq Master Mix (Rox) for:

- Gene expression analysis of cDNA
- Efficient amplification of rare cDNA and low copy DNA targets
- Fast thermal cycling for high-throughput real-time PCR applications

Benefits

- Enjoy efficient detection of rare targets.
- Rely on the stability of this master mix in automation and high-throughput plate stacking handling.
- Obtain high yield and a robust 5' nuclease based detection through use of the robust EagleTaq DNA Polymerase enzyme.
- Run the master mix on a variety of real-time PCR platforms that require Rox as reference dye.
- Minimize carryover PCR contamination with this 2x master mix containing dUTP, compatible with Uracil-DNA Glycosylase treatment (UNG).

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Function test: Average CT value of positive controls tested is between 18 and 28 cycles with starting template of 10 pg λDNA. Average CT value of test is within ±2 cycles of proven.

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
05 529 034 190	1 ml
05 529 026 190	50 ml

05529034190: Will be supplied as "KIT EAGLETAQ MMX WITH ROX 1mL RUO". Unit of Measure is "piece".
 05529026190: Will be supplied as "KIT EAGLETAQ MMX WITH ROX 50mL". Unit of Measure is "piece".
 For further processing only.
 Patent and License Disclaimer(s): 64

FastStart PCR Master

2x concentrated

Good value master (2x) for reliable results on block cyler systems and simple real-time monoplex applications.

Application

Use FastStart PCR Master for block cyler applications and for simple real-time monoplex applications. Containing the FastStart enzyme as hot start polymerase, mispriming by-products can be avoided in a broad range of applications.

Benefits

- **Profit from high cost-to-profit value.**
- **Avoid mispriming by-products using hot start enzymes.** Containing the FastStart enzyme as a hot start polymerase with a specially formulated buffer system, the master mix is ideally suited to avoid unwanted mispriming by-products in your PCR reaction.
- **Obtain fast results.** Profit from the speed of the master which allows fast cycling protocols.

Product Description

FastStart PCR Master is a ready-to-use, 2x concentrated master mix containing all reagents (except primers, probes, and template) required for block cyler assays and simple real-time RNA-detection assays, using either hydrolysis probe detection or HybProbe format. The FastStart PCR Master effectively suppresses mispriming by-products.

EC 2.7.7.7

Properties

Activation time (during initial PCR denaturation step): ≥2 minutes at +95°C

Temperature optimum for elongation:

Fragment length <3 kb: +72°C

Fragment length >3 kb: +68°C

The master mix can be stored in the refrigerator (+2 to +8°C) for at least 1 week without loss of activity and performance. It is stable at room temperature for at least one day.

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 100 mmol/l; KCl, 20 mmol/l; $(\text{NH}_4)_2\text{SO}_4$, 10 mmol/l; MgCl_2 , 4 mmol/l; dATP, dCTP, dGTP, dTTP (each 0.4 mmol/l)

pH value: 8.3±0.1 at +25°C

Volume activity: ≥50 U/ml

Unspecific endonucleases (λ DNA and MWM II DNA): Not detectable in up to 25 μl after 16 hours incubation at +37°C.

Exonucleases (^3H -DNA): Not detectable in up to 6 μl after 4 hours incubation at +37°C and +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 25 μl after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 25 μl after 1 hour incubation at +37°C

Performance test in PCR

(50 pg human genomic DNA, 365 bp tPA fragment): Corresponds to specification

(2 ng human genomic DNA, 1.1 kb collagen fragment): Corresponds to specification

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Background Information

For information on the FastStart hot start system, see **FastStart DNA Polymerase**

Cat. No.

04 659 155 103

Pack Size

custom fill

Will be supplied as "FastStart PCR Master IB". Unit of Measure is "ml".

For customized FastStart PCR Masters (e.g., with dUTP), please inquire.



For further processing only.

Patent and License Disclaimer(s): 49

AMV Reverse Transcriptase, recombinant, GMP Grade from Avian Myeloblastosis Virus, expressed in *E. coli*

AMV Reverse Transcriptase, recombinant, GMP Grade, transcribes RNA fragments up to 12 kb, providing high sensitivity in conjunction with high thermostability.

Application

Use AMV Reverse Transcriptase for synthesis of cDNA from total RNA or mRNA for:

- Two-step RT-PCR applications for amplification from RNA targets
- RT-PCR for detection of viral RNA
- TMA and NASBA nucleic acid amplification methods
- Synthesis of full-length cDNA for libraries or cloning
- Rapid amplification of cDNA end (RACE)
- Manufacture of amplification mixtures for applications with regulatory requirements (e.g., *in vitro* diagnostics, quality control)

Benefits

- **Achieve high sensitivity.** Use AMV Reverse Transcriptase for high sensitive RT-PCR applications in conventional thermal cyclers and real-time PCR instruments.
- **Obtain more full-length transcripts up to 12 kb.** Generate cDNA libraries with large inserts.
- **Reverse transcribe difficult templates.** Make use of the increased thermostability of this enzyme enabling reverse transcription at elevated temperatures to overcome RNA secondary structures (e.g., in GC-rich RNA templates) and achieve optimal reaction conditions for specific cDNA synthesis.
- **Efficiently label cDNA.** DIG-, biotin-, or dye-labeled nucleotides can be incorporated during cDNA synthesis.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to GMP (Good Manufacturing Practice) regulations.

EC 2.7.7.49

Properties

AMV Reverse Transcriptase, recombinant, GMP Grade, provides higher thermostability compared to the native forms of AMV or M-MLV reverse transcriptase, allowing higher temperatures for reverse transcription, thus achieving better performance with GC-rich RNA fragments and difficult secondary structures.

Enzyme activities: RNA-dependent DNA polymerase, DNA-dependent DNA polymerase, unwinding activity, RNase H (degrading RNA in RNA:DNA hybrids)

Recommended reaction temperature: +42 to +65°C

Substrates: Incorporates dNTP, ddNTP, dUPT, various labeled or modified nucleotides

Divalent ion requirement: Mg²⁺

Specification

Appearance: Clear, colorless solution

Storage buffer: Potassium phosphate, 200 mmol/l; DTT, 2 mmol/l; Triton X-100, 0.2% (v/v); glycerol, 50% (v/v); pH approximately 7.2

Volume activity: ≥20 U/μl

Specific activity: ≥50 kU/mg protein

Unit definition: One unit AMV Reverse Transcriptase, recombinant, GMP Grade, is defined as the amount of enzyme which incorporates 1 nmol of [³H] TMP into an acid insoluble product in 10 minutes at +37°C with poly(A)_x(dT)₁₅

Cat. No. Pack Size

03 203 166 103	20 kU 200 kU
-----------------------	-----------------

Will be supplied as "Rev.Transcriptase, AMV rec.". Unit of Measure is "kU".

The enzyme is filled as 20 or 200 KU per vial. Specify in the order, which filling.

The enzyme is supplied without reaction buffer.



For further processing only.

- For the best fit reaction buffer, use Transcriptor RT Buffer, see page 270

as substrate.

Purity (SDS PAGE): ≥90%

Unspecific endonucleases (MWM III DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 40 U after 4 hours incubation at +37°C.

Function test in RT-PCR (human skeletal muscle total RNA, 10 kb dystrophin gene fragment): Corresponds to reference

Bioburden: ≤50 CFU/ml

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

M-MLV Reverse Transcriptase, GMP Grade from Moloney Murine Leukemia Virus, expressed in *E. coli*

M-MLV Reverse Transcriptase generates full-length cDNA with high efficiency. It has a lower RNase H activity than AMV Reverse Transcriptase and lacks endonuclease activity.

Application

Use M-MLV Reverse Transcriptase for synthesis of cDNA from total RNA or mRNA for:

- Two-step RT-PCR applications for amplification from RNA targets
- RT-PCR for detection of viral RNA
- TMA and NASBA nucleic acid amplification methods
- Synthesis of full-length cDNA for libraries or cloning
- Rapid amplification of cDNA end (RACE)
- Manufacture of amplification mixtures for applications with regulatory requirements (*e.g.*, *in vitro* diagnostics, quality control)

Benefits

- **Achieve high sensitivity.** M-MLV Reverse Transcriptase lacks endonuclease activity and has much lower RNase H activity than AMV Reverse Transcriptase.
- **Obtain full-length transcripts up to 10 kb.** Generate cDNA libraries.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to GMP (Good Manufacturing Practice) regulations.

EC 2.7.7.49

Properties

M-MLV Reverse Transcriptase, GMP Grade, is highly processive and generates full length cDNA with high efficiency.

Enzyme activities: RNA-dependent DNA polymerase, DNA-dependent DNA polymerase, low RNase H activity, no endonuclease activity

Recommended reaction temperature: +37°C

Substrates: Incorporates dNTP, ddNTP, dUPT, various labeled or modified nucleotides

Divalent ion requirement: Mg²⁺

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 25 mmol/l; NaCl, 100 mmol/l; DTT, 10 mmol/l; EDTA,

Cat. No. **Pack Size**

04 707 486 103 200 kU

Will be supplied as "M-MLV RT Industrial GMP Grade, 200 KU".

Unit of Measure is "piece".

The enzyme is supplied without reaction buffer.



For further processing only.

0.1 mmol/l; Triton X-100, 0.01% (v/v); glycerol, 50% (v/v); pH approximately 8.4

Volume activity: ≥200-300 U/μl

Specific activity: ≥100 kU/mg protein

Unit definition: One unit M-MLV Reverse Transcriptase, GMP Grade, is defined as the amount of enzyme which incorporates 1 nmol of [³H]TMP into an acid insoluble product in 10 minutes at +37°C with poly(A)x(dT)₁₅ as substrate.

Purity (SDS PAGE): ≥90%

Unspecific endonucleases (MWM III DNA): Not detectable in up to 100 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 100 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 200 U after 1 hour incubation at +37°C.

Bioburden: ≤50 CFU/ml

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Transcriptor Reverse Transcriptase recombinant, expressed in *E. coli*

Transcriptor Reverse Transcriptase is the robust recombinant reverse transcriptase with thermostability up to +60°C, for transcription of RNA fragments up to 14 kb in two-step RT-PCR applications.

Application

Use Transcriptor Reverse Transcriptase for synthesis of cDNA from total RNA or mRNA for:

- Two-step RT-PCR applications using conventional thermal cyclers or real-time PCR instruments
- RT-PCR for detection of viral RNA
- TMA and NASBA nucleic acid amplification methods
- Synthesis of full-length cDNA up to 14 kb for libraries or cloning
- Rapid amplification of cDNA end (RACE)

Benefits

- **Achieve high sensitivity.** Use Transcriptor Reverse Transcriptase for high sensitive RT-PCR applications in conventional thermal cyclers and real-time PCR instruments.
- **Obtain more full-length transcripts up to 14 kb.** Generate cDNA libraries with large inserts.
- **Reverse transcribe difficult templates.** Make use of the increased thermostability of this enzyme enabling reverse transcription at elevated temperatures to overcome RNA secondary structures (e.g., in GC-rich RNA templates), and achieve optimal reaction conditions for specific cDNA synthesis.
- **Efficiently label cDNA.** DIG-, biotin-, or dye-labeled nucleotides can be incorporated during cDNA synthesis.

EC 2.7.7.49

Properties

Transcriptor Reverse Transcriptase offers higher thermostability compared to the native forms of AMV or M-MLV reverse transcriptase, allowing higher temperatures for reverse transcription, achieving high performance with GC-rich RNA fragments and difficult secondary structures.

Enzyme activities: RNA-dependent DNA polymerase, DNA-dependent DNA polymerase, unwinding activity, RNase H (degrading RNA in RNA:DNA

Cat. No.	Pack Size
03 531 252 103	custom fill

Will be supplied as "Transcriptor Bulk". Unit of Measure is "kU". The enzyme is supplied without reaction buffer.



For further processing only.

- For the best fit reaction buffer, use Transcriptor RT Buffer, see page 270

hybrids)

Recommended reaction temperature: +42 to +65°C

Substrates: Incorporates dNTP, ddNTP, dUPT, various labeled or modified nucleotides

Divalent ion requirement: Mg²⁺

Specification

Appearance: Clear, colorless solution

Storage buffer: Potassium phosphate, 200 mmol/l; DTT, 2 mmol/l; Triton X-100, 0.2% (v/v); glycerol, 50% (v/v), pH approximately 7.2

Volume activity: ≥20 U/μl

Specific activity: ≥50 kU/mg protein

Unit definition: One unit Transcriptor Reverse Transcriptase is defined as the amount of enzyme which incorporates 1 nmol of [³H]TMP into an acid insoluble product in 10 minutes at +37°C with poly(A)x(dT)15 as substrate.

Purity (SDS PAGE): ≥90%

Unspecific endonucleases (MWM III DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 25 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 40 U after 4 hours incubation at +37°C. **Function test in RT-PCR** (human skeletal muscle total RNA, 10 kb dystrophin gene fragment): Corresponds to reference

Function test in real-time RT-qPCR using the LightCycler® instrument (PBGD gene fragment from RNA standards): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Transcriptor RT Buffer 5x concentrated

Standard reaction buffer for Transcriptor Reverse Transcriptase.

Application

Use Transcriptor RT Buffer as an optimized reaction buffer for Transcriptor Reverse Transcriptase for synthesis of cDNA from total RNA or mRNA.

Benefits

- **Simplify reaction setup.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer for fast and easy setup of highly reproducible experiments.
- **Achieve excellent performance.** Take full advantage of Transcriptor Reverse Transcriptase using a reaction buffer specially optimized for this enzyme.

Specification

Appearance: Clear, colorless solution

Contents: Tris/HCl, 250 mmol/l; KCl, 150 mmol/l; MgCl₂, 40 mmol/l; pH approximately 8.5 at +25°C

Unspecific endonucleases (MWM III DNA): Not detectable after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable after 4 hours incubation at +37°C.

Function test in 2-step RT-PCR (human total RNA, 10 kb dystrophin gene fragment): Corresponds to reference

Function test in real-time RT-qPCR using the LightCycler® Instrument (RNA standards, PBGD gene fragment): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
03 531 325 103	1 ml

Will be supplied as "Transcriptor RT Buffer". Unit of Measure is "piece".



For further processing only.

AllStart RNA Master

High-end one-step RT-PCR master mix (2x) for high sensitivity and assay robustness.

Application

Use AllStart RNA Master for all one-step RNA PCR applications using instruments not requiring normalization with Rox. Due to its novel hot start system, AllStart RNA Master avoids mispriming by-products, making it ideal for multiplex assays with a large number of different primers. For high sensitivity, AllStart RNA Master is the one-step RT-PCR reagent mix of choice.

Benefits

- **Rely on a RT-PCR master mix suited for all target types.** The new hot start system avoids undesired by-products even for the most difficult targets.
- **Detect low abundant genes with ease.** Find that gene with less than 10 copies relying on AllStart Master's high sensitivity.
- **Obtain results quickly.** Profit from the fast RT-Step requiring only 3 to 10 minutes.
- **Ready for robotics.** Rely on the stability of the master mix ideal for RT-PCR automation. The viscosity of the master mix is optimized for accurate pipetting. AllStart Master mix is stable during setup and on the stacker for more than 12 hours.
- **Gain flexibility.** The AllStart Master mix concept allows you to use PCR instruments not requiring Rox normalization. AllStart RNA Master (Rox) is specifically for real-time PCR instruments requiring Rox normalization.

Product Description

AllStart RNA Master is a ready-to-use, 2x concentrated RT-PCR master mix containing all reagents (except primers, probes, and template) required for running quantitative, real-time RNA-detection assays in the hydrolysis probe detection format. AllStart RNA Master enables very sensitive detection and quantification of defined RNA sequences.

EC 2.7.7.7

Properties

This RT-PCR master mix can be stored in the refrigerator (+2 to +8°C) for at least 2 weeks without loss of activity and performance. It is stable at room temperature for at least one day.

Specification

Appearance: Clear, colorless solution

Performance test in qPCR

(human RNA, Rantes fragment, multiples reaction): Corresponds to specification

(human RNA, B2M fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
05 895 286 001	5 ml for up to 1,000 reactions at 10 µl final reaction volume
05 895 367 001	50 ml for up to 10,000 reactions at 10 µl final reaction volume

05895286001: Will be supplied as "AllStart RNA Master, 5 ml". Unit of Measure is "piece".

05895367001: Will be supplied as "AllStart RNA Master, 50 ml". Unit of Measure is "piece".

Contents

01. AllStart RNA Master Mix 2x concentrated
02. AllStart RNA Enzyme Mix 20x concentrated

05895286001: For life science research only. Not for use in diagnostic procedures.

05895367001: For further processing only.

Patent and License Disclaimer(s): 04

AllStart RNA Master (Rox)

High-end one-step RT-PCR master mix (2x) for high sensitivity and assay robustness.

Application

Use AllStart RNA Master (Rox) in all one-step RNA RT-PCR applications using instruments requiring normalization with Rox. Due to its novel hot start system, AllStart RNA Master (Rox) avoids mispriming by-products, making it ideal for multiplex assays with a large number of different primers. For high sensitivity, AllStart RNA Master (Rox) is the one-step RT-PCR master mix of choice.

Benefits

- **Rely on a RT-PCR master mix for all target types.** The new hot start system avoids undesired by-products even for the most difficult targets.
- **Detect low abundant genes with ease.** Find that gene with less than 10 copies due to AllStart RNA Master's high sensitivity.
- **Obtain results quickly.** Profit from the fast RT-Step requiring only 3 to 10 minutes.
- **Ready for robotics.** Rely on the stability of the master mix for RT-PCR automation. The viscosity of the master mix is optimized for accurate pipetting. AllStart RNA Master (Rox) mix is stable during setup and on the stacker for more than 12 hours.
- **Gain flexibility.** The AllStart RNA Master mix concept allows you to use PCR machines with or without Rox normalization. For real-time PCR not requiring Rox normalization, use AllStart RNA Master instead of AllStart RNA Master (Rox).

Product Description

AllStart RNA Master (Rox) is a ready-to-use, 2x concentrated RT-PCR master mix containing all reagents (except primers, probes, and template) required for running quantitative, real-time RNA-detection assays in the hydrolysis probe detection format. AllStart RNA Master (Rox) enables very sensitive detection and quantification of defined RNA sequences. A special Rox reference dye (FRET-ROX) enables you to run assays for all real-time PCR instruments in which Rox reference dye is needed for quantitative analysis.

EC 2.7.7.7

Properties

The RT-PCR master mix can be stored in the refrigerator (+2 to +8°C) for at least 2 weeks without loss of activity and performance. It is stable at room temperature for at least one day.

Specification

Appearance: Clear, colorless solution

Performance test in qPCR

(human RNA, Rantes fragment, multiples reaction): Corresponds to specification

(human RNA, B2M fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. Pack Size

05 999 375 001 5 ml
for up to 1,000 reactions at 10 µl final reaction volume


05 895 359 001 50 ml
for up to 10,000 reactions at 10 µl final reaction volume

05999375001: Will be supplied as "AllStart RNA Master (ROX), 5 ml". Unit of Measure is "piece".


05895359001: Will be supplied as "AllStart RNA Master (ROX), 50 ml". Unit of Measure is "piece".

Contents

01. AllStart RNA Master Mix (Rox), 2x concentrated
02. AllStart RNA Enzyme Mix, 20x concentrated

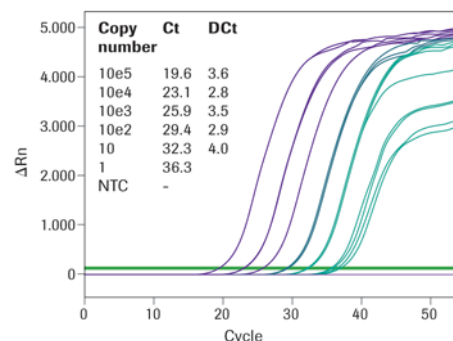

Highly sensitive detection of influenza A virus over a broad range of copy numbers. A serial dilution from 10^5 to 1 copy of viral RNA was reversely transcribed and amplified in one-step reactions of a final volume of 10 µl in 384-well plates on the Applied Biosystems 7900 Real-Time PCR System. Copy numbers of 10^5 and 10^4 were investigated in duplicate, copy numbers of 10^3 and lower were performed in quadruplicate. The following detection profile was used: 50°C 10 minutes; 95°C 2 minutes 55x (95°C 15 seconds; 60°C 35 seconds). Reactions without templates (NTC) generated no detectable fluorescence.

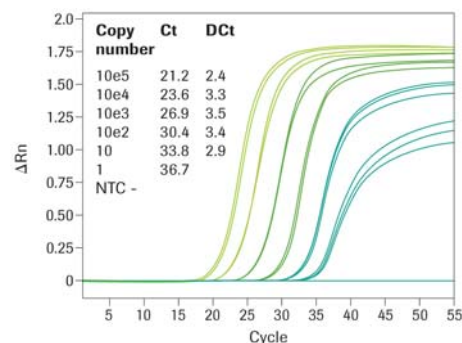

Efficient detection of hepatitis A virus. A serial dilution from 10^5 to 1 copy of hepatitis A viral RNA was reversely transcribed and amplified in one-step reactions of a final volume of 20 µl in 96-well plates on the Applied Biosystems® Real-Time PCR System. Copy numbers of 10^5 and 10^4 were investigated in duplicate, copy numbers of 10^3 and lower were performed in triplicate. The following detection profile was used: 50°C 10 minutes; 95°C 10 minutes; 55x (95°C 15 seconds; 60°C 35 seconds). Reactions without templates (NTC) generated no detectable fluorescence.

05999375001: For life science research only. Not for use in diagnostic procedures.

05895359001: For further processing only.

Patent and License Disclaimer(s): 04

HawkZ05 Fast One-Step RT-PCR Kit

HawkZ05 Fast One-Step RT-PCR Kit is a fast and robust, one tube, one enzyme RT-PCR reagent.

Application

Apply HawkZ05 Fast One-Step RT-PCR Kit for:

- High-throughput quantitative gene expression analysis
- Target detection and quantification
- Detection of rare transcripts
- Reverse transcription and amplification of RNA from limited samples

Benefits

- Enjoy increased specificity by employing hot start capability for both RT and PCR steps.
- Obtain reproducible detection of low copy numbers.
- Enjoy a broad dynamic range of eight logarithmic scales.
- Simplify the RT-PCR by using just one enzyme, buffer, and reaction step.
- Rely on this master mix for detection of target RNA with minor mismatches to primer sequences.

EC 2.7.7.7

Specification

Appearance: Clear, colorless solution

Function test: Average CT value of positive controls tested is between 20 and 30 cycles using starting template of 1×10^4 copy pAW 109 per reaction. Average CT value of real-time PCR test is within ± 2 cycles of the proven specification.

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
05 987 687 190	1 Kit 1x 5 ml master mix and 1x 1 ml Manganese acetate
05 987 695 190	1 Kit 2x 5 ml master mix and 2x 1 ml Manganese acetate

05987687190: Will be supplied as "KIT 1-STEP MMX W/ROX 1 x 5mL RUO". Unit of Measure is "piece".

05987695190: Will be supplied as "KIT 1-STEP MMX W/ROX 2 x 5mL LUO". Unit of Measure is "piece".

Contents

01. HawkZ05 Fast One-Step RT-PCR Master Mix
02. RMS Manganese Acetate (25 mM)

05987687190: For life science research only. Not for use in diagnostic procedures.

05987695190: For general laboratory use.

Patent and License Disclaimer(s): 65

dATP, PCR Grade

sodium salt, 100 mM

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTPs from the leading manufacturer of nucleotides.

Application

dATP, PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these nucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

dATP, PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

CAS: 1927-31-7

Properties

Nomenclature: 2'-Deoxy-adenosine-5'-triphosphate

Formula: C₁₀H₁₆N₅O₁₂P₃

Molecular weight: 491.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

dATP (1 μmol ± 15.0 A₂₆₀ units, pH 7.0): 100-110 mmol/l

dATP (high resolution HPLC method): ≥99 area%

dADP (HPLC): ≤0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 0.78±0.02

A₂₈₀/A₂₆₀: 0.15±0.01

A₂₉₀/A₂₆₀: ≤0.02

Function test in RT-PCR (RNA, human dystrophin, and mouse β-actin gene):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 42 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
04 631 056 103	20 ml (2,000 μmol)
11 889 516 103	100 ml (10,000 μmol)

04631056103: Will be supplied as "dATP PCR Grade, Sodium Solution, 20 ml". Unit of Measure is "μmol".

11889516103: Will be supplied as "dATPNa, Solution, (PCR Grade)". Unit of Measure is "μmol".



For further processing only.

dCTP, PCR Grade sodium salt, 100 mM

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTPs from the leading manufacturer of nucleotides.

Application

dCTP, PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these nucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

dCTP, PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

CAS: 2056-98-6

Properties

Nomenclature: 2'-Deoxy-cytidine-5'-triphosphate

Formula: C₉H₁₆N₃O₁₃P₃

Molecular weight: 467.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

dCTP (1 μmol ± 9.6 A₂₇₂ units, pH 7.0): 100-110 mmol/l

dCTP (high resolution HPLC method): ≥99 area%

dCDP (HPLC): ≤0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 0.82±0.02

A₂₈₀/A₂₆₀: 0.97±0.02

A₂₉₀/A₂₆₀: 0.30±0.02

Function test in RT-PCR (RNA, human dystrophin, and mouse β-actin gene):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 42 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

α-S-dCTP, Molecular Diagnostic Grade S-isomer, sodium salt, 100 mM

Modified dCTP for Strand Displacement Amplification and other special applications.

Cat. No. Pack Size

04 631 072 103 20 ml
(2,000 μmol)

11 889 508 103 100 ml
(10,000 μmol)

04631072103: Will be supplied as "dCTP PCR Grade, Sodium Solution, 20 ml". Unit of Measure is "μmol".

11889508103: Will be supplied as "dCTP, Na, Solution (PCR Grade)". Unit of Measure is "μmol".



For further processing only.

Cat. No. Pack Size

12 207 095 103 custom fill

Will be supplied as "Alpha-thio-d-CTP, Solution (Mol-DIA)". Unit of Measure is "μmol".



Application

Use α -S-dCTP in alternative amplification technologies such as Strand Displacement Amplification (SDA). In addition, these nucleotides can be used in applications, including amplification reactions when a higher resolution in capillary electrophoresis is desired. This product is the purified S-isomer, prepared using Roche's biocatalytical production process. Use of purified S-isomer eliminates the introduction of another isomer which does not participate in the amplification reaction as it is not a substrate for the polymerase.

Benefits

- Use only the reactive isomer in your amplification reaction.

Product Description

α -S-dCTP, Molecular Diagnostic Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice.

Properties

Nomenclature: α -Thio-2'-deoxy-cytidine-5'-triphosphate

Formula: C₉H₁₆N₃O₁₂SP₃ x H₂O

Molecular weight: 483.2 D

Specification

Appearance: Colorless solution

pH value: 8.1-8.5

α -S-d-CTP (1 μ mol \pm 9.6 A₂₇₂ units): 90-100 mmol/l

α -S-d-CTP, S-Isomer (HPLC): \geq 98.0 area%

α -S-d-CDP (HPLC): \leq 1.5 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 0.80-0.84

A₂₈₀/A₂₆₀: 0.95-1.00

A₂₉₀/A₂₆₀: 0.28-0.32

Stability: At -15 to -25°C within specification range for 12 months.

Quality

This α -S-dCTP is the purified S-isomer, not a diastereomeric mixture.

For further processing only.

**dGTP, PCR Grade
sodium salt, 100 mM**

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTPs from the leading manufacturer of nucleotides.

Application

dGTP, PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these nucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Cat. No. Pack Size

04 631 129 103 20 ml
(2,000 μ mol)

11 889 524 103 100 ml
(10,000 μ mol)

04631129103: Will be supplied as "dGTP PCR Grade, Sodium Solution, 20 ml". Unit of Measure is " μ mol".

11889524103: Will be supplied as "dGTP, Na, Solution (PCR Grade)". Unit of Measure is " μ mol".



For further processing only.

Product Description

dGTP, PGR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

CAS: 2564-35-4

Properties

Nomenclature: 2'-Deoxy-guanosine-5'-triphosphate

Formula: C₁₀H₁₆N₅O₁₃P₃

Molecular weight: 507.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

dGTP (1 μmol ± 13.7 A₂₅₂ units, pH 7.0): 100-110 mmol/l

dGTP (high resolution HPLC method): ≥99 area%

dGDP (HPLC): ≤0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 1.15±0.03

A₂₈₀/A₂₆₀: 0.67±0.02

A₂₉₀/A₂₆₀: 0.28±0.02

Function test in RT-PCR (RNA, human dystrophin, and mouse β-actin gene):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 42 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

**7-Deaza-2'-dGTP
lithium salt, 10 mM**

Use 7-Deaza-dGTP instead of dGTP in primer-extension reactions or PCR and receive a better resolution of GC-rich regions.

Application

Use 7-Deaza-dGTP as a substrate for most DNA polymerases, including Taq DNA polymerase.

7-Deaza-dGTP is used in the dideoxy-chain termination sequencing methods, in place of dGTP to overcome compression problems in gel electrophoresis when sequencing GC-rich stretches of DNA.

Benefits

- Achieve better sensitivity and resolution in primer-extension reactions.

Properties

Nomenclature: 7-Deaza-2'-deoxy-guanosine-5'-triphosphate

Formula: C₁₁H₁₇N₄O₁₃P₃

Molecular weight: 506.2 D

Specification

Appearance: Clear, colorless solution

pH value: 6.8-7.2

7-Deaza-2'-dGTP (1 μmol ± 13.4 A₂₅₀ units): 10.0-11.0 mmol/l

7-Deaza-2'-dGTP (HPLC): ≥95 area%

7-Deaza-2'-dGDP (HPLC): ≤4 area%

A₂₅₀/A₂₆₀: 0.84±0.04

A₂₈₀/A₂₆₀: 0.65±0.03

Cat. No. Pack Size

10 987 891 103 custom fill

Will be supplied as "7-Deaza-2'-deoxy-GTP, Di-Li". Unit of Measure is "μmol".



For further processing only.

Nucleotides

deoxyNTPs

A_{290}/A_{260} : 0.53±0.03

Stability: At -15 to -25°C within specification range for 30 months.

Background Information

Comparison of 7-Deaza-dGTP with dGTP and dITP showed that 7-Deaza-dGTP gives enhanced resolution compared with dGTP, resulting increased readability over long sequence regions compared with dITP. For sequencing reactions, dGTP is replaced by the same amount of 7-Deaza-dGTP in all four dideoxy-NTP solutions.

7-Deaza-dGTP performs equally well to all the other types of dideoxy sequencing and polymerization techniques. Partial substitution of 7-Deaza-dGTP for dGTP in PCR can improve the yield of reaction products for GC-rich templates containing strong secondary structures. Elimination of spurious GC-hydrogen bonding and relaxation of the secondary structure results in more efficient and specific PCR-product synthesis.

Incorporation of 7-Deaza-dGTP into DNA alters the fluorescent staining and electrophoretic mobility of the DNA.

dITP, PCR Grade sodium salt, 100 mM

High quality dITP from the leading manufacturer of nucleotides for the preparation of polynucleotides.

Application

dITP, PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics. Use dITP for the preparation of poly(dI) x poly(dC) and poly[d(I-C)] with DNA polymerase and dCTP.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these nucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

dITP, PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

CAS: 95648-77-4

Properties

Nomenclature: 2'-Deoxy-inosine-5'-triphosphate

Formula: C₁₀H₁₅N₄O₁₃P₃

Molecular weight: 492.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

dITP (1 μmol ± 12.3 A₂₄₉ units, pH 7.0): 100-110 mmol/l

Cat. No.	Pack Size
12 158 124 103	100 ml (10,000 μmol)

Will be supplied as "Desoxy-ITP, Na-Lsg., (PCR-Grade)". Unit of Measure is "μmol".



For further processing only.

dITP (high resolution HPLC method): ≥ 99 area%

dIDP (HPLC): ≤ 0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

A_{250}/A_{260} : 1.67 ± 0.03

A_{280}/A_{260} : 0.25 ± 0.03

A_{290}/A_{260} : 0.03 ± 0.02

Function test in RT-PCR (RNA, human dystrophin, and mouse β -actin gene):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 42 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

dTTP, PCR Grade sodium salt, 100 mM

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTPs from the leading manufacturer of nucleotides.

Application

dTTP, PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these nucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

dTTP, PCR Grade, is supplied in sealed and CO_2 -proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

CAS: 365-08-2

Properties

Nomenclature: 2'-Deoxy-thymidine-5'-triphosphate

Formula: $\text{C}_{10}\text{H}_{17}\text{N}_2\text{O}_{14}\text{P}_3$

Molecular weight: 482.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

dTTP ($1 \mu\text{mol} \pm 9.5 A_{267}$ units, pH 7.0): 100-110 mmol/l

dTTP (high resolution HPLC method): ≥ 99 area%

dTDP (HPLC): ≤ 0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

A_{250}/A_{260} : 0.64 ± 0.02

A_{280}/A_{260} : 0.74 ± 0.02

Cat. No. Pack Size

04 631 137 103 20 ml
(2,000 μmol)

11 889 559 103 100 ml
(10,000 μmol)

04631137103: Will be supplied as "dTTP PCR Grade, Sodium Solution, 20 ml". Unit of Measure is " μmol ".

11889559103: Will be supplied as "dTTP, Na, Solution (PCR Grade)". Unit of Measure is " μmol ".



For further processing only.

Nucleotides

deoxyNTPs

A_{290}/A_{260} : 0.24±0.02

Function test in RT-PCR (RNA, human dystrophin, and mouse β -actin gene):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 42 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

dUTP, PCR Grade sodium salt, 100 mM

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTPs from the leading manufacturer of nucleotides.

Application

dUTP, PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics. Avoid DNA carryover contamination between PCRs that can be a source of false positives.

To decontaminate PCR and RT-PCR mixes, dUTP is incorporated in place of dTTP into the PCR product. Reaction mixes are then pre-treated with Uracil-DNA Glycosylase (UNG) before amplification to degrade contaminating dUTP-containing PCR products.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Avoid false positive results.** Decontaminate your reaction mixes using Uracil-DNA Glycosylase before a new PCR amplification by incorporating dUTP instead of dTTP into PCR products.
- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these nucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

dUTP, PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

CAS: 1173-82-6

Properties

Nomenclature: 2'-Deoxy-uridine-5'-triphosphate

Formula: C₉H₁₅N₂O₁₄P₃

Molecular weight: 468.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

dUTP (1 μ mol \triangleq 9.9 A₂₆₀ units, pH 7.0): 100-110 mmol/l

dUTP (high resolution HPLC method): \geq 99 area%

dUDP (HPLC): \leq 0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

Cat. No.	Pack Size
04 631 145 103	20 ml (2,000 μ mol)
11 889 532 103	100 ml (10,000 μ mol)

04631145103: Will be supplied as "dUTP PCR Grade, Sodium Solution, 20 ml". Unit of Measure is " μ mol".

11889532103: Will be supplied as "dUTP, Na, Solution (PCR Grade)". Unit of Measure is " μ mol".



For further processing only.

A_{250}/A_{260} : 0.74±0.02

A_{280}/A_{260} : 0.38±0.02

A_{290}/A_{260} : 0.04±0.01

Function test in RT-PCR (RNA, human dystrophin, and mouse β -actin gene):

Corresponds to specification

Stability: At -15 to -25°C within specification range for 42 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

NucleoMix (dTTP), PCR Grade sodium salt, 40 mM (10 mmol/l each dNTP)

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTP NucleoMixes from the leading manufacturer of nucleotides.

Application

NucleoMix (dTTP), PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics. The ready-to-use mix eliminates process steps.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Save time and resources.** Eliminate process steps using a ready-to-use dNTP mix.
- **Rely on quality design and obtain reproducible results.** From the production process to formulation and final packaging, product release and long shelf life, lot after lot, NucleoMix meets all the requirements of the demanding kit manufacturer.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

NucleoMix (dTTP), PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

Concentration (sum, 1 μ mol \pm 10.7 A₂₆₀ units, pH 7.0): 40-44 mmol/l

Identity (HPLC diode array detector):

dATP: 9-12 mmol/l

dCTP: 9-12 mmol/l

dGTP: 9-12 mmol/l

dTTP: 9-12 mmol/l

dNTP (HPLC, sum of 4 peaks): \geq 99 area%

dNDP (HPLC, sum of 4 peaks): \leq 0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

Function test in RT-PCR: Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
03 186 083 103	100 ml

Will be supplied as "Nucleomix (10 mmol/l, with dTTP)". Unit of Measure is "ml".



For further processing only.

NucleoMix (dUTP), PCR Grade sodium salt, 40 mM (10 mmol/l each dNTP)

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTP NucleoMixes from the leading manufacturer of nucleotides.

Application

NucleoMix (dUTP), PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics. The ready-to-use mix eliminates process steps. Avoid carryover contamination between PCRs, a significant source of false positives. To decontaminate PCR and RT-PCR reagent mixes, dUTP is incorporated in place of dTTP. Subsequent reactions should be treated with Uracil-DNA Glycosylase (UNG) before amplification to degrade dUTP-containing PCR products.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Save time and resources.** Eliminate process steps using a ready-to-use dNTP mix.
- **Avoid false positive results.** Decontaminate before starting a new amplification reaction using Uracil-DNA Glycosylase to digest the dUTP-containing NucleoMix before PCR.
- **Rely on quality design and obtain reproducible results.** From the production process to formulation and final packaging, product release and long shelf life, lot after lot, NucleoMix meets all the requirements of the demanding kit manufacturer.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

NucleoMix (dUTP), PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

Concentration (sum, 1 μmol ± 11.0 A₂₆₀ units, pH 7.0): 40-44 mmol/l

Identity (HPLC diode array detector):

dATP: 9-12 mmol/l

dCTP: 9-12 mmol/l

dGTP: 9-12 mmol/l

dUTP: 9-12 mmol/l

dNTP (HPLC, sum of 4 peaks): ≥99 area%

dNDP (HPLC, sum of 4 peaks): ≤0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

Function test in PCR: Corresponds to specification

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No. **03 186 075 103** Pack Size **100 ml**

Will be supplied as "Nucleomix (10 mmol/l, with dUTP)". Unit of Measure is "ml".



For further processing only.

NucleoMix (dTTP), PCR Grade sodium salt, 100 mM (25 mmol/l each dNTP)

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTP NucleoMixes from the leading manufacturer of nucleotides.

Application

NucleoMix (dTTP), PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics. The ready-to-use mix eliminates process steps.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Save time and resources.** Eliminate process steps using a ready-to-use dNTP mix.
- **Rely on quality design and obtain reproducible results.** From the production process to formulation and final packaging, product release and long shelf life, lot after lot, NucleoMix meets all the requirements of the demanding kit manufacturer.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

NucleoMix (dTTP), PGR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

Concentration (sum, 1 $\mu\text{mol} \pm 10.7 A_{260}$ units, pH 7.0): 100-110 mmol/l

Identity (HPLC diode array detector):

dATP: 23-28 mmol/l

dCTP: 23-28 mmol/l

dGTP: 23-28 mmol/l

dTTP: 23-28 mmol/l

dNTP (HPLC, sum of 4 peaks): ≥ 99 area%

dNDP (HPLC, sum of 4 peaks): ≤ 0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

Function test in RT-PCR: Corresponds to specification

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.	Pack Size
04 920 171 103	20 ml
03 991 016 103	100 ml

03991016103: Will be supplied as "Nucleomix, 25mM, Na, St. (PCR-Gr.)". Unit of Measure is "ml".

04920171103: Will be supplied as "NucleoMix PCR Grade, 25 mmol/l, 20 ml". Unit of Measure is "ml".



For further processing only.

NucleoMix (dUTP), PCR Grade sodium salt, 100 mM (25 mmol/l each dNTP)

For outstanding, consistent performance in amplification reactions, use GMP-manufactured dNTP NucleoMixes from the leading manufacturer of nucleotides.

Application

NucleoMix (dUTP), PCR Grade, is designed for amplification reactions where high-quality reagents are required, such as for *in vitro* diagnostics. The ready-to-use mix eliminates process steps. Avoid carryover contamination between PCRs, a significant source of false positives. To decontaminate PCR and RT-PCR reagent mixes, dUTP is incorporated in place of dTTP. Subsequent reactions should be treated with Uracil-DNA Glycosylase (UNG) before amplification to degrade dUTP-containing PCR products.

Benefits

- **Insist on the highest purity for optimal results even with difficult targets.** PCR grade dNTPs are rigorously controlled for amplification inhibitors.
- **Save time and resources.** Eliminate process steps using a ready-to-use dNTP mix.
- **Avoid false positive results.** Decontaminate before starting a new amplification reaction using Uracil-DNA Glycosylase to digest the dUTP-containing NucleoMix before PCR.
- **Rely on quality design and obtain reproducible results.** From the production process to formulation and final packaging, product release and long shelf life, lot after lot, NucleoMix meets all the requirements of the demanding kit manufacturer.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

NucleoMix (dUTP), PCR Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for amplification.

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

Concentration (sum, 1 μmol ± 11.0 A₂₆₀ units, pH 7.0): 100-110 mmol/l

Identity (HPLC diode array detector):

dATP: 23-28 mmol/l

dCTP: 23-28 mmol/l

dGTP: 23-28 mmol/l

dUTP: 23-28 mmol/l

dNTP (HPLC, sum of 4 peaks): ≥99 area%

dNDP (HPLC, sum of 4 peaks): ≤0.9 area%

DNases/RNases: Negative

Nicking activity: Negative

Function test in PCR: Corresponds to specification

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No.


Pack Size

04 980 905 103 20 ml

Will be supplied as "NucleoMix PCR Grd (dU), 25 mmol/l, 20 ml".
Unit of Measure is "ml".



For further processing only.

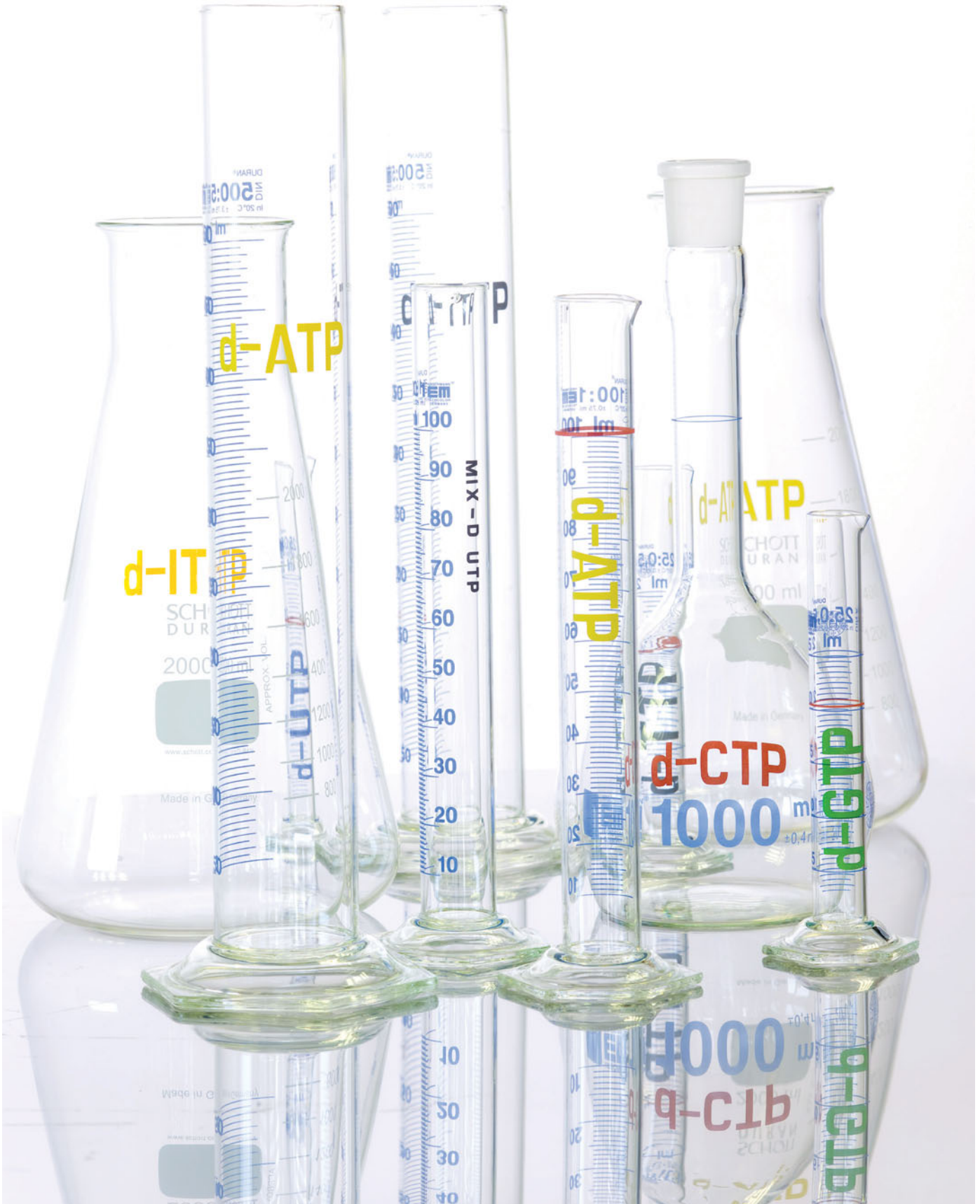

Amplification

Molecular Diagnostics

ddATP, Sequencing Grade

sodium salt, 10 mM

For consistent performance in sequencing applications and primer extension reactions, use dideoxyribonucleotides from the leading manufacturer of nucleotides.

Application

ddATP, Sequencing Grade, in dNTP mixes acts as chain terminator in sequencing and primer extension reactions. High purity and proven manufacturing processes produce a product ideal for the manufacturing of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these dideoxyribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the high lot-to-lot performance of this product.

Product Description

ddATP, Sequencing Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The pH is adjusted to match conditions for amplification.

Properties

Nomenclature: 2',3'-Dideoxy-adenosine-5'-triphosphate

Formula: C₁₀H₁₆N₅O₁₁P₃

Molecular weight: 475.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

ddATP (1 μmol ± 15.3 A₂₆₀ units, pH 7.0): 10.0-11.0 mmol/l

ddATP (HPLC): ≥98 area%

ddADP (HPLC): ≤1.5 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 0.76±0.02

A₂₈₀/A₂₆₀: 0.15±0.01

A₂₉₀/A₂₆₀: ≤0.02

Stability: At -15 to -25°C within specification range for 24 months.

Quality

All dideoxynucleotides can be GMP manufactured upon request.

Cat. No.	Pack Size
12 158 175 103	100 ml (1,000 μmol)

Will be supplied as "Dideoxy-ATP, Na, Sequencing Grade, Solution".
Unit of Measure is "μmol".



For further processing only.

ddCTP, Sequencing Grade sodium salt, 10 mM

For consistent performance in sequencing applications and primer extension reactions, use dideoxyribonucleotides from the leading manufacturer of nucleotides.

Application

ddCTP, Sequencing Grade, in dNTP mixes acts as chain terminator in sequencing and primer extension reactions. High purity and proven manufacturing processes produce a product ideal for the manufacturing of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these dideoxyribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the high lot-to-lot performance of this product.

Product Description

ddCTP, Sequencing Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The pH is adjusted to match conditions for amplification.

Properties

Nomenclature: 2',3'-Dideoxy-cytidine-5'-triphosphate

Formula: C₉H₁₆N₃O₁₂P₃

Molecular weight: 452.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

ddCTP (1 μmol \triangleq 9.6 A₂₇₂ units, pH 7.0): 10.0-11.0 mmol/l

ddCTP (HPLC): \geq 98 area%

ddCDP (HPLC): \leq 1.5 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 0.76 \pm 0.02

A₂₈₀/A₂₆₀: 1.06 \pm 0.02

A₂₉₀/A₂₆₀: 0.33 \pm 0.02

Stability: At -15 to -25°C within specification range for 24 months.

Quality

All dideoxynucleotides can be GMP manufactured upon request.

Cat. No. Pack Size

12 158 183 103 100 ml
(1,000 μmol)

Will be supplied as "Dideoxy-CTP, Na, Sequencing Grade, Solution". Unit of Measure is "μmol".



For further processing only.

ddGTP, Sequencing Grade

sodium salt, 10 mM

For consistent performance in sequencing applications and primer extension reactions, use dideoxyribonucleotides from the leading manufacturer of nucleotides.

Application

ddGTP, Sequencing Grade, in dNTP mixes acts as chain terminator in sequencing and primer extension reactions. High purity and proven manufacturing processes produce a product ideal for the manufacturing of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these dideoxyribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the high lot-to-lot performance of this product.

Product Description

ddGTP, Sequencing Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The pH is adjusted to match conditions for amplification.

Properties

Nomenclature: 2',3'-Dideoxy-guanosine-5'-triphosphate

Formula: C₁₀H₁₆N₅O₁₂P₃

Molecular weight: 491.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

ddGTP (1 μmol ± 13.7 A₂₅₂ units, pH 7.0): 10.0-11.0 mmol/l

ddGTP (HPLC): ≥98 area%

ddGDP (HPLC): ≤1.5 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 1.16±0.02

A₂₈₀/A₂₆₀: 0.67±0.02

A₂₉₀/A₂₆₀: 0.27±0.01

Stability: At -15 to -25°C within specification range for 24 months.

Quality

All dideoxynucleotides can be GMP manufactured upon request.

Cat. No.	Pack Size
12 158 191 103	100 ml (1,000 μmol)

Will be supplied as "DideoxyGTP, Na, Sequencing Grade, Solution".
Unit of Measure is "μmol".



For further processing only.

ddTTP, Sequencing Grade sodium salt, 10 mM

For consistent performance in sequencing applications and primer extension reactions, use dideoxyribonucleotides from the leading manufacturer of nucleotides.

Application

ddTTP, Sequencing Grade, in dNTP mixes acts as chain terminator in sequencing and primer extension reactions. The high purity and carefully designed manufacturing process are ideally suited for the manufacturing of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these dideoxyribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the high lot-to-lot performance of this product.

Product Description

ddTTP, Sequencing Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The pH is adjusted to match conditions for amplification.

Properties

Nomenclature: 2',3'-Dideoxy-thymidine-5'-triphosphate

Formula: C₁₀H₁₇N₂O₁₃P₃

Molecular weight: 466.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

ddTTP (1 μmol ≙ 8.4 A₂₆₀ units, pH 7.0): 10.0-11.0 mmol/l

ddTTP (HPLC): ≥98 area%

ddTDP (HPLC): ≤1.5 area%

DNases/RNases: Negative

Nicking activity: Negative

A₂₅₀/A₂₆₀: 0.60±0.02

A₂₈₀/A₂₆₀: 0.83±0.03

A₂₉₀/A₂₆₀: 0.31±0.02

Stability: At -15 to -25°C within specification range for 24 months.

Quality

All dideoxynucleotides can be GMP manufactured upon request.

Cat. No. Pack Size

12 158 205 103 100 ml
(1,000 μmol)

Will be supplied as "Dideoxy-TTP, Na, Sequencing Grade, Solution".
Unit of Measure is "μmol".



For further processing only.

ATP, Molecular Diagnostic Grade

sodium salt, 100 mM

For outstanding, consistent performance in transcription reactions, use GMP-manufactured ribonucleotides from the leading manufacturer of nucleotides.

Application

ATP, Molecular Diagnostics Grade, is designed for *in vitro* transcription reactions and alternative amplification reactions where high-quality reagents are required, such as for the manufacture of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these ribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

ATP, Molecular Diagnostic Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for reverse transcription and amplification.

CAS: 56-65-5

Properties

Nomenclature: Adenosine-5'-triphosphate

Formula: C₁₀H₁₆N₅O₁₃P₃

Molecular weight: 507.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

ATP (1 μmol ± 15.0 A₂₆₀ units, pH 7.0): 100-110 mmol/l

ATP (HPLC): ≥98 area%

ADP (HPLC): ≤1.5 area%

AMP (HPLC): ≤0.5 area%

Functional transcription assay: Corresponds to specification

DNases/RNases: Negative

Nicking activity: Negative

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No. _____ **Pack Size** _____

04 980 824 103 100 ml

Will be supplied as "ATP Mol Dia Grade, 100 mmol/l, 100 ml". Unit of Measure is "ml".



For further processing only.

CTP, Molecular Diagnostic Grade sodium salt, 100 mM

For outstanding, consistent performance in transcription reactions, use GMP-manufactured ribonucleotides from the leading manufacturer of nucleotides.

Application

CTP, Molecular Diagnostics Grade, is designed for *in vitro* transcription reactions and alternative amplification reactions where high-quality reagents are required, such as for the manufacture of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these ribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

CTP, Molecular Diagnostic Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for reverse transcription and amplification.

CAS: 65-47-4

Properties

Nomenclature: Cytidine-5'-triphosphate

Formula: C₉H₁₆N₃O₁₄P₃

Molecular weight: 483.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

CTP (1 μmol ± 7.4 A₂₆₀ units, pH 7.0): 100-110 mmol/l

CTP (HPLC): ≥98 area%

CDP (HPLC): ≤1.5 area%

CMP (HPLC): ≤0.5 area%

Functional transcription assay: Corresponds to specification

DNases/RNases: Negative

Nicking activity: Negative

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No. **Pack Size**

04 980 875 103 100 ml

Will be supplied as "CTP Mol Dia Grade, 100 mmol/l, 100 ml". Unit of Measure is "ml".



For further processing only.

GTP, Molecular Diagnostic Grade

sodium salt, 100 mM

For outstanding, consistent performance in transcription reactions, use GMP-manufactured ribonucleotides from the leading manufacturer of nucleotides.

Application

GTP, Molecular Diagnostics Grade, is designed for *in vitro* transcription reactions and alternative amplification reactions where high-quality reagents are required, such as for the manufacture of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these ribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

GTP, Molecular Diagnostic Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for reverse transcription and amplification.

CAS: 86-01-1

Properties

Nomenclature: Guanosine-5'-triphosphate

Formula: C₁₀H₁₆N₅O₁₄P₃

Molecular weight: 523.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

GTP (1 μmol ± 13.7 A₂₅₂ units, pH 7.0): 100-110 mmol/l

GTP (HPLC): ≥98 area%

GDP (HPLC): ≤1.5 area%

GMP (HPLC): ≤0.5 area%

Functional transcription assay: Corresponds to specification

DNases/RNases: Negative

Nicking activity: Negative

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No. _____ **Pack Size** _____

04 980 859 103 100 ml

Will be supplied as "GTP Mol Dia Grade, 100 mmol/l, 100 ml". Unit of Measure is "ml".



For further processing only.

UTP, Molecular Diagnostic Grade sodium salt, 100 mM

For outstanding, consistent performance in transcription reactions, use GMP-manufactured ribonucleotides from the leading manufacturer of nucleotides.

Application

UTP, Molecular Diagnostics Grade, is designed for *in vitro* transcription reactions and alternative amplification reactions where high-quality reagents are required, such as for the manufacture of *in vitro* diagnostics.

Benefits

- **Rely on quality by design.** From the production process to formulation and final packaging, product release and long shelf life, these ribonucleotides meet all the requirements of the demanding kit manufacturer.
- **Obtain consistent results.** Rely on the robust lot-to-lot performance of this product.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to Good Manufacturing Practice (GMP) regulations.

Product Description

UTP, Molecular Diagnostic Grade, is supplied in sealed and CO₂-proof bottles to ensure a stable pH during shipment on dry ice. The manufacturing process is fully validated including the final filling step. The pH is adjusted to match conditions for reverse transcription and amplification.

CAS: 63-39-8

Properties

Nomenclature: Uridine-5'-triphosphate

Formula: C₉H₁₅N₂O₁₅P₃

Molecular weight: 484.2 D

Specification

Appearance: Clear, colorless solution

pH value: 8.1-8.5

UTP (1 μmol ± 9.9 A₂₆₀ units, pH 7.0): 100-110 mmol/l

UTP (HPLC): ≥98 area%

UDP (HPLC): ≤1.5 area%

UMP (HPLC): ≤0.5 area%

Functional transcription assay: Corresponds to specification

DNases/RNases: Negative

Nicking activity: Negative

Stability: At -15 to -25°C within specification range for 30 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Cat. No. **Pack Size**


04 979 818 103 100 ml

Will be supplied as "UTP Mol Dia Grade, 100 mmol/l, 100 ml". Unit of Measure is "ml".



For further processing only.

Custom-Made Primers and Probes Produced According to cGMP

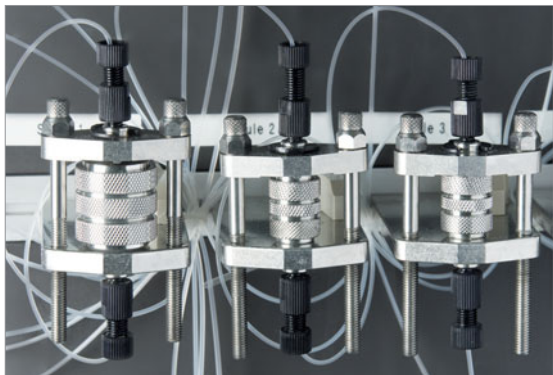

Columns for oligonucleotide synthesis



Final filling step takes place in a laminar flow box

DNA- and RNA-based diagnostics require unmodified or labeled oligonucleotides. Roche has the expertise, capability, and capacity to routinely develop and produce high-quality customized oligonucleotides in large quantities.

Roche has extensive experience in the development and production of custom primers and probes for diagnostics manufacturers, with a focus on filling and labeling according to customer needs:

- Oligonucleotides, modified or conjugated (*e.g.*, biotin, amino, or 3'-phosphate groups)
- Probes labeled with dyes, such as LightCycler® dyes
- Oligonucleotides coupled to customer-designed labels
- Modified amidites and CPG supports
- Additional analytical methods

The Roche oligonucleotide production facility allows synthesis scales from 15 μ moles to a few mmoles per batch. With a variety of proven purification protocols optimized for specific types of oligonucleotides, you are assured the high purity, quality, regulatory standards, and safety standards during synthesis.

- Choose from several production-scale synthesizers: 15 – 2000 μ mol.
- Select from a variety of preparative HPLC purification lines: 15 – 400 μ mol.
- Utilize ultrafiltration (UF) cross-contamination-free concentration/desalting capabilities: 15 – 2000 μ mol.
- Take advantage of proven oligo-protein-conjugate capabilities.
- Use comprehensive in-process control equipment: analytical HPLC systems, spectrometers, mass spectrometry.

Professional Support for Manufacturers of Molecular Diagnostics

Our personalized service features the detailed technical advice of Roche oligonucleotide experts right from the beginning.

We are happy to discuss project-specific details (*e.g.*, scale, specifications, sequence, regulatory requirements, time frame, patents) in-depth, and your project will be handled with complete confidentiality, with confidential disclosure agreements provided upon request.

- A pre-production test synthesis run may be required to ensure we are meeting your needs.
- Ensure quality with your own eyes. After your large scale synthesis run under regulatory demands, a representative sample is sent to you for function testing.
- Individual specifications will be set for each product.
- Production will proceed in complete accordance to your needs and regulatory requirements.



During purification by chromatography, each fraction is analyzed for purity.

GMP Grade Oligonucleotides

Roche is one of a few suppliers in the world that can produce customized oligonucleotides according to GMP. In answer to the rising importance of GMP manufacturing for diagnostic reagents, we built a state-of-the-art oligonucleotide production facility in early 2001.

In our processes, all raw materials must meet the highest quality and regulatory demands, and all vendors of building blocks and solvents must be qualified. Even the water used for the production process is monitored for microbial and RNase/DNase contamination. All production processes for oligonucleotides follow strict regulatory requirements and are designed to keep oligonucleotides in solution at all times, minimizing the risk of cross-contamination.

- All production areas are Class 100,000 with reference to non-viable particle counts ($> 0.5 \mu\text{cft}$) and less than 200 cfu/m^3 . Filling of bulk solutions takes place in a laminar flow box (Class 100).
- Documentation includes batch records with full traceability of raw material and equipment.
- A strict product segregation policy prevents all cross-contamination.
- The “dedication/single-use” concept ensures that materials having direct contact with the product (*e.g.*, graduated cylinders, fraction vials) are sterilized and used only once. Equipment in direct contact with the product (*e.g.*, prep. HPLC columns, UF membranes) are also dedicated for a single product or parameter.
- Processes and methods of production of oligonucleotides are validated and controlled by an independent quality measure.

Production process controls also apply to in-house personnel.

- A defined flow of personnel, material, and equipment maintains traceability and the highest quality standards.
- Access to dedicated laboratories is controlled for security and safety.

Amplification

Additional Products

GC-RICH Solution

5x concentrated

Additive for efficient amplification of difficult DNA templates as GC-rich targets up to 5 kb in the polymerase chain reaction (PCR).

Application

Use GC-RICH Solution in combination with PCR buffer to improve the amplification of GC-rich template DNA, GC-rich targets and repetitive sequences, as well as mixtures of nucleic acids with varying GC content.

Benefits

- **Enhance performance with difficult targets.** Use GC-RICH Solution to quickly optimize amplification of difficult templates.
- **Obtain reliable results.** Rely on the high lot-to-lot consistency of this product, thoroughly tested for a consistent quality.

Specification

Appearance: Clear, colorless solution.

Unspecific endonucleases (λ DNA and MWM III DNA): Not detectable in up to 20 μ l after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μ l after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 20 μ l after 1 hour incubation at +37°C.

Function test in PCR (200 ng human genomic DNA, 284 bp ApoE fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
----------	-----------

05 917 158 103	1 ml
----------------	------

12 161 583 103	custom fill
----------------	-------------

12161575103: Will be supplied as "GC-rich solution, MBF". Unit of Measure is "piece".

12161583103: Will be supplied as "GC-rich solution". Unit of Measure is "l".



For further processing only.

LightCycler® 480 Multiwell Plate 384

white, 4 barcodes

High-performance reaction device tailor-made for the LightCycler® 480 Instruments, 384-well version.

Application

Use LightCycler® 480 Multiwell Plates 384 for real-time PCR and melting curve analysis applications with the LightCycler® 480 Instrument, 384-well version.

Benefits

- **Gain excellent performance.** Benefit from plates, which fully exploit the special thermal and optical characteristics of the LightCycler® 480 Instrument, 384-well.
- **Simplify PCR setup.** Take advantage of the plate design with cut-away corner for error-free positioning.
- **Simplify PCR sample tracking.** Rely on identical barcodes fixed on each side of the plate for easy assay and sample tracking.
- **Profit from the cost efficiency.** The plate design is compatible with a variety of robotic systems. Save costs for sealing foils using proven automated sealing.

Product Description

LightCycler® 480 Multiwell Plates 384 are full-skirted, white polypropylene reaction devices supplied with a barcoded label on each side of the plate. Their dimensions (inner/outer well geometry) and material are designed to achieve optimal results using the LightCycler® 480 Instrument, 384-well, for reaction volumes from 5 to 20 μ l.

Cat. No.	Pack Size
----------	-----------

05 217 555 001	5 x 10 plates without sealing foils
----------------	-------------------------------------

Will be supplied as "LC480 Multiwell Plate 384, 4 bar codes". Unit of Measure is "piece". The plates are supplied without sealing foils.

For further processing only.

Patent and License Disclaimer(s): 51

Specification

Appearance (visual control): Corresponds to specification
Sample evaluation (visual control): Corresponds to specification
Stability: Store at +15 to +25°C.

LightCycler® 480 Multiwell Plate 96
white, 4 barcodes

High-performance reaction device tailor-made for the LightCycler® 480 Instruments, 96-well version.

Application

LightCycler® 480 Multiwell Plates 96 are used for real-time PCR and melting curve analysis applications with the LightCycler® 480 Instrument, 96-well version.

Benefits

- **Gain excellent performance.** Benefit from plates, which fully exploit the special thermal and optical characteristics of the LightCycler® 480 Instrument, 96-well.
- **Simplify PCR setup.** Take advantage of the plate design with cut-away corner for error-free positioning.
- **Simplify PCR sample tracking.** Rely on identical barcodes fixed on each side of the plate for easy assay and sample tracking.
- **Profit from the cost efficiency.** The plate design is compatible with a variety of robotic systems. Save costs for sealing foils using proven automated sealing.

Product Description

LightCycler® 480 Multiwell Plates 96 are half-skirted, white polypropylene reaction devices supplied with a barcode label on each side of the plate. Their dimensions (inner/outer well geometry) and material are designed to achieve optimal results using the LightCycler® 480 Instrument, 96-well, for reaction volumes from 10 to 100 µl.

Specification

Appearance (visual control): Corresponds to specification
Sample evaluation (visual control): Corresponds to specification
Stability: Store at +15 to +25°C.

MgCl₂ Stock Solution
25 mM

PCR Grade MgCl₂ Stock Solution.

Application

Use this MgCl₂ Stock Solution in combination with any PCR buffer without MgCl₂ to optimize the magnesium concentration.

Benefits

- **Obtain reliable results.** Rely on the high lot-to-lot consistency of this product, thoroughly tested for a constant quality.
- **Simplify your PCR setup.** Save time producing a suitable, pure MgCl₂ solution with this ready-to-use formulation.

Specification

Appearance: Clear, colorless solution
Contents: MgCl₂, 25 mmol/l; pH approximately 8.3 at +20°C

Cat. No.	Pack Size
05 220 319 001	5 x 10 plates without sealing foils

Will be supplied as "LC480 Multiwell Plate 96, 4 bar codes". Unit of Measure is "piece". The plates are supplied without sealing foils. For further processing only.

Patent and License Disclaimer(s): 51

Cat. No.	Pack Size
11 600 770 103	1 ml

11 600 788 103	custom fill
----------------	-------------

11600770103: Will be supplied as "MgCl₂-Slit. 25mM MPB". Unit of Measure is "piece".

11600788103: Will be supplied as "MgCl₂-Solution 25mM/Bulk". Unit of Measure is "ml".



For further processing only.

Unspecific endonucleases (λ DNA and MWM III DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μ l after 16 hours incubation at +65°C.

Function test in PCR (0.01 ng λ DNA, 0.5 kb lambda fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Protector RNase Inhibitor from rat lung, expressed in *E. coli*

Protector RNase Inhibitor is useful in any application where RNA degradation by RNases is critical.

Application

Use Protector RNase Inhibitor to:

- Protect RNA in cDNA synthesis reactions, *in vitro* transcription/translation systems, *in vitro* RNA synthesis, and *in vitro* virus replication
- Protect RNA during RNA isolation and purification
- Produce RNase activity-free antibodies

Benefits

- **Protect your RNA over a wide range of reaction conditions.** Use Protector RNase Inhibitor also with more thermostable reverse transcriptases for cDNA synthesis at elevated temperatures. It is stable at temperatures up to +55°C and at pH 5.0 to 9.0.
- **Don't worry about interference.** Even at high concentrations, Protector RNase Inhibitor does not interfere with other enzymes and reagents commonly used to analyse RNA.
- **Insist on a highly-purified preparation.** Each batch is function tested for the absence of nucleic acid modifying activities.

Properties

Protector RNase Inhibitor is a protein of 50 kD which inhibits enzymatic activity of RNases by noncovalently binding to the active site.

Activity: A minimum of 1 mmol/l DTT is required to maintain the inhibitor in its active state; a pH value between 5.0 and 9.0 is recommended (isoelectric point is at pH 4.7).

Inactivation: At temperatures > +65°C or under severe denaturing conditions the inhibitory activity disappears.

Specification

Appearance: Clear, colorless solution

Storage buffer: Hepes-KOH, 20 mmol/l; KCl, 50 mmol/l; DTT, 8 mmol/l; glycerol, 50% (v/v); pH approximately 7.6

Volume activity: ≥ 40 U/ μ l

Specific activity: ≥ 80 kU/mg protein

Unit definition: One unit Protector RNase Inhibitor is defined as the amount which inhibits 50% of the activity of 5 ng RNase A (activity measured according to Blackburn by inhibition of hydrolysis of cyclic cytidine-monophosphoric acid).

Purity (SDS PAGE): $\geq 95\%$

Unspecific endonucleases (MWM III DNA): Not detectable in up to 400 U after 1 hour incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 400 U after 1 hour incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 400 U after 1 hour incubation at +37°C.

Ribonuclease activity after thermal inactivation of inhibitor (10 minutes

Cat. No. **Pack Size**

03 335 429 103 custom fill

Will be supplied as "Protector RNase Inhibitor". Unit of Measure is "MU".



For further processing only.

at +65°C, MS2 RNA): Not detectable in up to 160 U after 1 hour incubation at +37°C.

Function test in RT-PCR (Titan One Tube RT-PCR Kit with control RNA): Corresponds to reference

Function test in real-time RT-qPCR using the LightCycler® instrument (PBGD gene fragment from RNA standards): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

RNase Inhibitor, recombinant, GMP Grade from rat lung, expressed in *E. coli*

RNase Inhibitor, recombinant, GMP Grade, is useful in any application where RNA degradation by RNases is critical.

Application

Use RNase Inhibitor, recombinant, GMP Grade, to:

- Protect RNA in cDNA synthesis reactions, *in vitro* transcription/translation systems, *in vitro* RNA synthesis, and *in vitro* virus replication
- Protect RNA during RNA isolation and purification
- Produce RNase activity-free antibodies
- Manufacture reaction mixtures for applications with regulatory requirements (*e.g.*, *in vitro* diagnostics, quality control)

Benefits

- **Protect your RNA over a wide range of reaction conditions.** Use RNase Inhibitor, recombinant, GMP Grade, also with more thermostable reverse transcriptases for cDNA synthesis at elevated temperatures. It is stable at temperatures up to +55°C and at pH 5.0 to 9.0.
- **Don't worry about interference.** Even at high concentrations, RNase Inhibitor, recombinant, GMP Grade, does not interfere with other enzymes and reagents commonly used to analyze RNA.
- **Insist on a highly-purified preparation.** Each batch is extensively function tested for the absence of nucleic acid modifying activities.
- **Stay ahead of regulatory requirements.** High quality manufacturing, quality control and documentation according to GMP (Good Manufacturing Practice) regulations.

Properties

RNase Inhibitor, recombinant, GMP Grade, is a protein of 50 kD which inhibits enzymatic activity of RNases by noncovalently binding to the active site.

Activity: A minimum of 1 mmol/l DTT is required to maintain the inhibitor in its active state; a pH value between 5.0 and 9.0 is recommended (isoelectric point is at pH 4.7).

Inactivation: At temperatures >+65°C or under severe denaturing conditions the inhibitory activity disappears.

Specification

Appearance: Clear, colorless solution

Storage buffer: Hepes-KOH, 20 mmol/l; KCl, 50 mmol/l; DTT, 8 mmol/l; glycerol, 50% (v/v); pH approximately 7.6

Volume activity: ≥40 U/μl

Specific activity: ≥80 kU/mg protein

Unit definition: One unit RNase Inhibitor, recombinant, GMP Grade, is defined as the amount which inhibits 50% of the activity of 5 ng RNase A (activity measured according to Blackburn by inhibition of hydrolysis of cyclic cytidine-monophosphoric acid).

Purity (SDS PAGE): ≥95%

DNA content: <100 pg/mg protein

Cat. No.	Pack Size
04 488 105 103	100 kU
03 247 058 103	1.4 MU

04488105103: Will be supplied as "RNase Inhibitor Ind. GMP Grade, 100 KU". Unit of Measure is "piece".

03247058103: Will be supplied as "RNase Inh. GMP, BIK 1.4M U/ bottle". Unit of Measure is "MU".



For further processing only.

Amplification

Additional Products

Unspecific endonucleases (MWM III DNA): Not detectable in up to 400 U after 1 hour incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 400 U after 1 hour incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 400 U after 1 hour incubation at +37°C.

Ribonucleases after thermal inactivation of inhibitor (10 minutes at +65°C, MS2 RNA): Not detectable in up to 160 U after 1 hour incubation at +37°C.

Function test in RT-PCR (Titan One Tube RT-PCR Kit with control RNA): Corresponds to reference

Function test in real-time RT-qPCR using the LightCycler® instrument (PBGD gene fragment from RNA standards): Corresponds to reference

Bioburden: <50 CFU/ml

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Manufactured under GMP (Good Manufacturing Practice) regulations.

Random Octamer Primer

2.5x concentrated solution

Use Random Octamer Primer for cDNA synthesis and to prime synthesis in oligonucleotide labeling.

Application

Use Random Octamer Primer for priming oligonucleotide labeling and cDNA synthesis.

Benefits

- Obtain consistent results in your RT-PCR reaction using the Random Octamer Primer.
- Produce high specific activity labeled DNA probes for *in situ* hybridization, northern and Southern blots.

Product Description

Random Octamer Primer is a mixture of 8-mer primers comprised of all possible sequences.

Specification

Appearance: Colorless solution

Contents: Tris/HCl, 120 mmol/l; MgCl₂, 12.5 mmol/l; octamer primer, 750 µg/ml; β-mercaptoethanol, 25 mmol/l; pH 7.0

Function test in octamer primed labeling (Cy5-dCTP):

Amount of DNA: ≥2.8 µg

Amount of fluorescently labeled dye: ≥100 pmol

Base to dye ratio: 40-80

Stability: At -15 to -25°C within specification range for 18 months.

Cat. No.	Pack Size
05 109 701 103	7 ml

Will be supplied as "Random Octamer (2.5X), MPB". Unit of Measure is "piece".



For further processing only.

T4 Gene 32 Protein, recombinant

recombinant from T4 phage, expressed in *E. coli*, solution

T4 Gene 32 Protein, recombinant, is a DNA-binding protein specific for single-stranded DNA, which can be used to improve the outcome of reverse transcription and PCR.

Cat. No.	Pack Size
05 793 262 103	custom fill

Will be supplied as "T4 Gene 32 Protein, rec.". Unit of Measure is "mg".

For further processing only.

Application

Use T4 Gene 32 Protein, recombinant, for:

- Optimization of reverse transcription and PCR (addition of T4 Gene 32 Protein to the reaction mixture can increase yield, specificity and efficiency of cDNA synthesis and DNA amplification)
- Stimulation of *in vitro* DNA synthesis
- Stabilization of single-stranded regions of DNA and RNA
- Sequencing of DNA with strong secondary structures
- Site-specific mutagenesis experiments using T4 DNA Ligase or T4 DNA Polymerase
- Complete digestion of DNA by restriction enzymes

Benefits

- Add T4 Gene 32 Protein to reaction mixes to improve performance of PCR and RT-PCR reactions.
- Rely on a contamination-controlled recombinant product with high lot-to-lot consistency.

Properties

Molecular weight: 35 kD

pH optimum: About 8.0

Isoelectric point: pH 5.5

Inactivation: After 20 minutes heat denaturation at +65°C the DNA binding activity is abolished.

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; NaCl, 100 mmol/l; EDTA, 1.0 mmol/l; DTT, 0.5 mmol/l; glycerol, 50% (v/v), pH approximately 8.0

ssDNA binding (shift in gelelectrophoresis): Corresponds to reference

Unspecific endonucleases (λ DNA): Not detectable in up to 50 μ g after 1 hour incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 50 μ g after 1 hour incubation at +37°C.

Exonucleases (3 H-DNA): Not detectable in up to 50 μ g after 1 hour incubation at +37°C.

Single-strand specific exonucleases (M13mp9 DNA): Not detectable in up to 50 μ g after 1 hour incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 50 μ g after 1 hour incubation at +37°C.

Absence of contaminating nucleic acids

(LightCycler® UniTOOL Resolight assay for detection of bacterial and fungal DNA): <1 genome equivalent/50 μ g

(LightCycler® PCR assay for detection of human β -Globin gene): <1 gene copy/50 μ g

Stability: At -15 to -25°C within specification range for 12 months.

Uracil-DNA Glycosylase, heat-labile from marine bacterium BMTU 3346, expressed in *E. coli*

Uracil-DNA Glycosylase (UNG) with an increased heat intolerance is the enzyme of choice for prevention of PCR carryover contamination.

Application

Use Uracil-DNA Glycosylase, heat-labile, (UNG) for prevention of carryover contamination with DNA amplification products in PCR. Always use dUTP containing PCR mixtures to enable decontamination by UNG treatment. In contrast to the UNG variant from *E.coli*, this heat-labile enzyme is completely

Cat. No.	Pack Size
11 780 565 103	custom fill

Will be supplied as "Uracil-DNA Glycosylase, heatlab". Unit of Measure is "kU".



For further processing only.

Amplification

Additional Products

inactivated in the initial heat denaturation step of a common PCR protocol and the formed PCR product will not be degraded.

Note: For high sensitive real-time PCR, specially optimized LightCycler® Uracil-DNA Glycosylase is recommended.

Benefits

- **Prevent carryover contamination.** Avoid the risk of false-positive PCR results by simple addition of this reagent to your PCR mixture.
- **Secure your PCR products.** No degradation of your amplified DNA due to the fast inactivation of this heat-labile enzyme in the beginning of the PCR protocol.

EC 3.2.2.15

Properties

Uracil-DNA Glycosylase hydrolyzes uracil-glycosidic bonds in DNA, creating abasic sites where the DNA is cleaved by heat, alkali, or endonuclease treatment. This heat-labile enzyme is easily inactivated by heat denaturation.

Specificity: Hydrolyzes uracil-glycosidic bonds in single- and double-stranded DNA; no activity on dU-free natural DNA and RNA.

Incubation: +15 to +25°C for 10 minutes are recommended for treatment of PCR mixtures; at higher temperatures the enzyme stability decreases.

Half life at +40°C: About 2 minutes

Heat inactivation: +95°C for 2 minutes are sufficient for inactivation

pH optimum: 8.3-8.9

Inhibition: Activity does not depend on metal ions; no inhibition in presence of EDTA or other chelating reagents.

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.1 mmol/l; glycerol, 50% (v/v); Nonidet P40, 0.5% (v/v); Tween 20, 0.5% (v/v); pH 8.0±0.1 at +4°C

Volume activity: ≥ 1 U/μl

Unit definition: One unit Uracil-DNA Glycosylase, heat-labile, is defined as the amount of enzyme required to completely degrade 1 μg purified single-stranded uracil-containing DNA (bacteriophage M13, grown in *E.coli* CJ236 dut-ung-) at +37°C within 60 minutes. For comparison, one Lindahl unit is comparable to 520,000 units based on our unit definition. One Lindahl unit is defined as the amount of enzyme required to release 1 mol uracil at +37°C in 1 minute.

Unspecific endonucleases (MWM III DNA and M13mp9 ssDNA): Not detectable in up to 20 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 U after 16 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 10 U after 4 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 20 U after 4 hours incubation at +65°C.

Function test, DNA decontamination (complete elimination of 10,000 copies of uracil-containing template DNA in a PCR assay): Corresponds to specification

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 18 months.

Water, PCR Grade

Specially purified, double-distilled, deionized, and autoclaved water.

Application

This product is specially tested and qualified for use in all PCR or RT-PCR. It can be used in conjunction with all endpoint or real-time PCR assays whenever highest quality water is required. All LightCycler® kits and many of the PCR amplification products, available from Roche Applied Science, contain Water, PCR Grade.

Benefits

- For high quality PCR, use pre-tested Water, PCR Grade.
- Trust this specially purified, double distilled, deionized, autoclaved water.
- There are no detectable endonucleases, ribonucleases and nicking activities.

CAS: 7732-18-5

Specification

Appearance: Clear, colorless solution

Content: Water

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +65°C.

Ribonucleases (MS2 RNA): Not detectable in up to 20 µl after 1 hour incubation at +37°C.

Performance test in PCR (0.01 ng λDNA, 0.5 kb λDNA fragment): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Function test: Each lot of Water, PCR Grade is tested using PCR amplification of a 0.5 kb target from 0.01 ng λDNA. When the PCR product is analyzed using gel electrophoresis, a clear, defined band is obtained.

Absence of nucleases: Water, PCR Grade, is tested using nucleic acids serving as nuclease substrates. No endonucleases, ribonucleases, or nicking activities are detected according to the strict Quality Control procedures.

Cat. No.	Pack Size
03 315 959 001	25 ml 1 x 25 ml
03 315 843 001	100 ml 4 x 25 ml

Will be supplied as "Water, PCR-Grade". Unit of Measure is "piece".



For life science research only. Not for use in diagnostic procedures.

Water, PCR Grade

Specially purified, double-distilled, deionized, and autoclaved water.

Application

Use Water, PCR Grade, for the preparation for solutions used in molecular diagnostics or molecular biology.

Benefits

- **Protect your expensive raw materials.** Use Water, PCR Grade, which is extensively tested to avoid the introduction of contaminants.
- **Obtain process security.** Use the same high quality water in all production lots.

CAS: 7732-18-5

Specification

Appearance: Clear solution

pH value: 5.0-7.0

Bioburden: <1 CFU/100 ml

TOC: <0.5 ppm

Conductivity: <1.3 µS/cm

DNases/RNases: Negative

Endotoxins (LAL assay): <0.25 EU/ml

Stability: At +15 to +25°C within specification range for 24 months.

Quality

Water, PCR Grade, is specially purified, double-distilled, deionized, and autoclaved. Nucleases that degrade DNA and RNA are not detectable. The product is tested for total organic content (TOC) and bioburden to avoid contamination with microorganisms and nucleic acids. In addition, Water, PCR Grade, is tested for endotoxins.

Cat. No.

Pack Size

03 036 430 103 0.5 Liter

Will be supplied as "Water PCR grade". Unit of Measure is "l".
For further processing only.

Anti-Digoxigenin-AP, Fab fragments polyclonal antibodies from sheep, lyophilizate

Cat. No.	Pack Size
11 082 736 103	custom fill

Fab fragments of a polyclonal anti-digoxigenin antibody conjugated to alkaline phosphatase for nonradioactive detection.

Will be supplied as "Anti-digoxigenin AP-conjugate". Unit of Measure is "kU".

For further processing only.

Application

Use Anti-Digoxigenin-AP, Fab fragments, for the detection of digoxigenin-labeled compounds using:

- cDNA array and nonradioactive DNA sequencing blot
- Colony/plaque hybridization
- Dot blot and ELISA
- Gel shift assay
- Immunohistochemistry
- *In situ* hybridization
- Northern, Southern, western blot
- RNase protection assay

Benefits

- Benefit from a hazard free labeling and detection system without the risks associated with radioactive assays.
- Achieve higher sensitivity and faster results compared to isotopic procedures.
- Profit from single-step detection using any substrate of alkaline phosphatase.
- Use Anti-Digoxigenin-AP, Fab fragments, for all your applications in combination with other Roche products.

Properties

Molecular weight: 50 kD (Fab fragment), 140 kD (alkaline phosphatase)

Specificity: Specific to digoxigenin and digoxin.

Cross-reactivity:

Digitoxin and digitoxigenin: < 1%

None to other human estrogen and androgen steroids, *e.g.*, estradiol or testosterone.

Digoxina: Not known

Binding:

Conjugate does not bind to itself.

Normally one molecule conjugate binds to one molecule digoxigenin, although there are two possible binding sites for digoxigenin.

Unspecific binding to RNA has not been reported.

Specification

Appearance: Clear, colorless solution

Volume activity (+37°C): 750-825 U/ml

Sensitivity: ≤0.1 pg

Stability: At +2 to +8°C within specification range for 24 months.

Anti-Digoxigenin-POD, Fab fragments polyclonal antibodies from sheep, lyophilizate

Fab-fragments of a polyclonal anti-digoxigenin antibody conjugated to horseradish peroxidase for nonradioactive detection.

Application

Use Anti-Digoxigenin-POD, Fab fragments, for the detection of digoxigenin-labeled compounds using:

- Dot blots and ELISA
- Immunohistochemistry
- *In situ* hybridization
- Western blot

Benefits

- Benefit from a hazard free labeling and detection system without the health risks associated with radioactive assays.
- Achieve higher sensitivity and faster results compared to isotopic procedures.
- Profit from the broad Roche product portfolio related to the DIG system.

Properties

Molecular weight: Molecular weight is within a broad range of 200 kD to several million kD depending on how large the complexes are.

Specificity: Specific to digoxigenin and digoxin.

Cross-reactivity: None to other human estrogen or androgen steroids, *e.g.*, estradiol or testosterone.

Specification

Appearance: Clear, brown solution

Volume activity (+25°C, ABTS): 3,400-4,600 U/ml

Storage buffer: Potassium phosphate 10 mmol/l; NaCl 100 mmol/l; succrose, 7.5%; RPLA4, 5 mg/ml; pH 7.5

pH 5.5 treatment (≥30 minutes): Corresponds to specification

Performance test (immune reaction signal height): Corresponds to specification

Stability: At -196°C within specification range for 48 months.

Cat. No. **Pack Size**

11 210 360 103 custom fill

Will be supplied as "Anti-digoxigenin POD-conjugate". Unit of Measure is "kU".



For further processing only.

Anti-Digoxigenin-Rhodamine, Fab fragments

polyclonal antibodies from sheep, lyophilizate

Fab-fragments of a polyclonal anti-digoxigenin antibody conjugated to rhodamine for nonradioactive detection.

Application

Use Anti-Digoxigenin-Rhodamine, Fab fragments, for the detection of digoxigenin-labeled compounds using:

- Digoxigenin-labeled sugars in glycoconjugate research
- Fluorescent *in situ* hybridization (FISH)
- Immunohistochemistry
- *In situ* hybridization

Benefits

- Simplify detection by direct visualization using a fluorescence microscope.
- Benefit from hazard free labeling and detection system without the risk associated with radioactive assays.
- Achieve higher sensitivity and faster results compared to isotopic procedures.
- Use Anti-Digoxigenin-Rhodamine, Fab fragments, for all your applications in combination with other Roche products.

Properties

Molecular weight: 50 kD (not conjugated) ≤ MW ≤ 54.22 kD (fully conjugated). Two to three rhodamine molecules are fused to one anti-digoxigenin, Fab fragment.

Specificity: Specific to digoxigenin and digoxin.

Cross-reactivity: None to other human estrogen and androgen steroids, *e.g.*, estradiol or testosterone.

Specification

Appearance: Clear, red, fluorescent solution

Performance test (*in situ* hybridization): Corresponds to specification

Stability: At -60 to -90°C within specification range for 36 months.

Cat. No.	Pack Size
11 210 378 103	custom fill

Will be supplied as "Anti-digoxigenin-rhodamine". Unit of Measure is "mg active ingredient".



For further processing only.

Klenow Enzyme, Labeling Grade from *E. coli* lysogenic NM 964, solution

Klenow enzyme is a DNA-dependent 5'-3' polymerase with 3'-5' exonuclease activity to synthesize DNA complementary to single-stranded DNA templates.

Application

Use Klenow Enzyme for:

- Random-primed DNA labeling using random oligonucleotides as primers for the incorporation of nonradioactively labeled and ³²P-labeled nucleotides
- Fill-in reaction for blunt-end formation of 3'-recessed (staggered) ends

Benefits

- Obtain consistent performance compared to the native enzyme due to high lot-to-lot consistency of the recombinant enzyme.
- Rely on the secure supply as calf thymus is no longer used.

Product Description

Klenow Enzyme is a DNA-dependent 5'-3' polymerase with 3'-5' exonuclease activity. It lacks the 5'-3' exonuclease activity of the native enzyme. Klenow Enzyme catalyzes the addition of mononucleotides to the 3'-OH terminus of DNA. This activity is used to synthesize DNA complementary to single-stranded DNA templates.

EC 2.7.7.7

Properties

Molecular weight: 68 kD

Specification

Appearance: Colorless solution

Storage buffer: Potassium phosphate, 50 mmol/l; DTE, 1 mmol/l; glycerol, 50% (v/v); pH approximately 7.0 at +4°C

Volume activity: ≥2,000 U/ml

Specific activity: ≥5,000 U/mg

Unit definition: One unit is the enzyme activity which incorporates 10 nmol of total nucleotides into an acid-precipitable fraction in 30 minutes under assay conditions at +37°C with poly [d(A-T)] as primer (Richardson, C.C. & Kornberg, A. (1994) J. Biol. Chem. **244**, 2996).

Purity (SDS gel electrophoresis, using up to 50 U enzyme): ≥90%

Protein: ≤0.4 mg/ml

Nicking activity (pBR322 DNA): Not detectable in up to 50 U after 16 hours incubation at +37°C.

Function test using Random Primed DNA Labeling Kit (>50% incorporation radioactive nucleotids after 30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No. _____ **Pack Size** _____

11 010 492 103 custom fill

Will be supplied as "Klenow Enzyme, Labeling Grade". Unit of Measure is "kU".



For further processing only.

Terminal Transferase, recombinant from calf thymus, expressed in *E. coli*, solution

Terminal Transferase catalyzes the addition of deoxy- and dideoxynucleoside triphosphates to the 3'-OH ends of double- and single-stranded DNA fragments and oligonucleotides.

Application

Use Terminal Transferase to add homopolymer tails to DNA fragments in cloning experiments, such as addition of overhanging ends onto cDNAs for easier cloning and labeling of 3' ends of double- and single-stranded DNA (e.g., oligonucleotides) with radioactive labeled nucleotides or nucleotides labeled with haptens, e.g., digoxigenin or biotin.

Benefits

- Profit from high lot-to-lot consistency due to the recombinant production process.
- Obtain outstanding performance compared to the native enzyme.

Product Description

Terminal Transferase catalyzes the template independent addition of deoxy- and dideoxynucleoside triphosphates to the 3'-OH ends of double- and single-stranded DNA fragments and oligonucleotides. Terminal Transferase incorporates digoxigenin-, biotin-, and fluorochrome-labeled deoxy- and dideoxynucleoside triphosphates, as well as radioactively labeled deoxy- and dideoxynucleoside triphosphates. The supplied 5x concentrated reaction buffer facilitates optimal tailing of all types of double-stranded DNA ends: blunt ended, with 3' overhang, or with 5' overhang.

EC 2.7.7.31

Properties

The enzyme catalyzes a template-independent addition of dNTPs or of a single ddNTP to the 3'-OH ends of double- or single-stranded DNA.

It accepts radioactively labeled nucleotides and nucleotides labeled with haptens, such as digoxigenin or biotin.

For activity, the enzyme requires the presence of divalent metal ions, preferably Co^{2+} .

Specification

Appearance: Clear, colorless solution

Storage buffer: Potassium phosphate, 60 mmol/l; KCl, 150 mmol/l; 2-mercaptoethanol, 1 mmol/l; Triton X-100, 0.5% (v/v); glycerol, 50% (v/v); pH approximately 7.2 at +4°C

Volume activity (Co^{2+}): $\geq 400 \times 10^3$ U/ml

Specific activity (Co^{2+}): $\geq 200 \times 10^3$ U/mg

Unit definition: One unit is the enzyme activity that leads to an incorporation of 1 nmol

dTMP into acid insoluble products within 30 minutes at +37°C under assay conditions (cacodylate, 200 mmol/l; Co^{2+} , 1 mmol/l) using $\text{d}(\text{pT})_6$ as primer.

Unspecific endonucleases (MWM II DNA): Not detectable in up to 400 U after 4 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 400 U after 4 hours incubation at +37°C.

Exonucleases (^3H -DNA): $\leq 0.1\%$

Function test (tailing reaction on a 30-mer oligonucleotide): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Background Information

The enzyme is shipped without reaction buffer. Please inquire to obtain an optimized buffer system.

Cat. No. **Pack Size**

03 289 869 103 custom fill

Will be supplied as "Terminale Transferase, recombi". Unit of Measure is "kU".



For further processing only.

- For the appropriate divalent metal ion, use Cobalt Chloride Solution, see page 310
- For the best fit reaction buffer, use Terminal Transferase Reaction Buffer, see page 310

Terminal Transferase Reaction Buffer 5x concentrated

Reaction buffer for Terminal Transferase, recombinant.

Application

Use Terminal Transferase Reaction Buffer for optimal results with Roche's Terminal Transferase, recombinant.

Specification

Appearance: Clear, colorless solution

Contents: Potassium cacodylate, 1 mol/l; Tris/HCl, 125 mmol/l; BSA, 1.25 mg/ml; pH 6.6±0.1 at +20°C

pH value (+20°C): 6.50-6.70

Unspecific endonucleases (MWM II DNA): Not detectable in up to 20 µl after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +37°C.

Performance test (tailing of a 30-mer oligonucleotide): Corresponds to specification

Stability: At -15 to -25°C within specification range for 18 months.

Cat. No.	Pack Size
06 280 927 103	1 ml

Will be supplied as "React buffer, term.transf., Sar. 1 ml". Unit of Measure is "piece".



For further processing only.

Cobalt Chloride Solution 25 mM solution

Cobalt chloride solution for the activation of terminal transferase and other enzymes.

Application

Use Cobalt Chloride Solution as a cofactor for many enzymes.

Specification

Appearance: Clear, reddish solution

Unspecific endonucleases (MWM II DNA): Not detectable in up to 20 µl after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 µl after 16 hours incubation at +37°C.

Performance test (tailing of a 30-mer oligonucleotide): Corresponds to specification

Stability: At -15 to -25°C within specification range for 18 months.

Cat. No.	Pack Size
05 895 391 103	1 ml

Will be supplied as "Cobalt Chloride Sol., Sarstedt 1 ml". Unit of Measure is "piece".



For further processing only.

5-Aminoallyl-UTP powder, lithium salt

Aminoallyl-modified nucleotides are used for efficient target labeling and a wide variety of dyes and haptens for subsequent post-labeling reactions.

Application

5-Aminoallyl-16-UTP is a substrate for SP6, T3, and T7 RNA polymerase. It can replace UTP in the *in vitro* transcription reaction for RNA labeling.

Linearized template DNA with T7, SP6, or T3 promoter is *in vitro* transcribed with the corresponding RNA polymerases using ATP, GTP, CTP, UTP, and aminoallyl-16-UTP, respectively.

Lacking bulky dye groups, aminoallyl-modified nucleotides are incorporated at extremely high and reproducible levels. Reaction of the amine-modified nucleic acid with an excess of amine-reactive dye or hapten results in high and consistent labeling efficiencies, regardless of the dye or hapten chosen. The two-step labeling method eliminates the need to optimize the labeling reaction.

Benefits

- Generate stronger signals using efficient incorporation of amino groups for post-labeling.
- Label targets with high uniformity and consistency.
- Increase the options for target labels by post-labeling.
- Save time with ready-to-use solutions.

Properties

Nomenclature: 5-Aminoallyl-uridine-5'-triphosphate

Formula: C₁₂H₁₆N₃O₁₅P₃Li₄

Molecular weight: 563 D

Specification

Appearance: White powder

5-Aminoallyl-UTP (1 μmol ≅ 9.73 A₂₄₀ units; phosphate buffer, 0.1 mol/l, pH 7.0): ≥90%

Purity (HPLC): ≥75 area%

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Cat. No. **Pack Size**

11 230 271 103 custom fill

Will be supplied as "5-Aminoallyl-UTP Tetra-Li". Unit of Measure is "g".



For further processing only.

5-Aminoallyl-dUTP powder, lithium salt

Aminoallyl-modified nucleotides are used for efficient target labeling and a wide variety of dyes and haptens for subsequent post-labeling reactions.

Application

5-Aminoallyl-16-dUTP replaces dTTP in labeling reactions using:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)
- T4 and T7 DNA polymerase
- Taq DNA polymerase and reverse transcriptases

Lacking bulky dye groups, aminoallyl-modified nucleotides can be incorporated at extremely high and reproducible levels. Reaction of the amine-modified nucleic acid with an excess of amine-reactive dye or hapten results in high and consistent labeling efficiencies, regardless of the dye or hapten chosen. The two-step labeling method eliminates the need to optimize the labeling reaction.

Cat. No. **Pack Size**

11 115 502 103 custom fill

Will be supplied as "5-Allylamino-2'-deoxy-UTP, Lithium Salt". Unit of Measure is "μmol".



For further processing only.

Benefits

- Generate stronger signals using efficient incorporation of amino groups for post-labeling.
- Label targets with high uniformity and consistency.
- Increase the options for target labels by post-labeling.
- Save time with ready-to-use solutions.

CAS: 1173-82-6

Properties

Nomenclature: 5-Aminoallyl-2'-deoxy-uridine-5'-triphosphate

Formula: C₁₂H₁₆N₃O₁₄P₃Li₄

Molecular weight: 547 D

Specification

Appearance: White powder

5-Aminoallyl-dUTP (1 μmol Δ 9.73 A₂₄₀ units; phosphate buffer, 0.1 mol/l, pH 7.0): ≥90%

Purity (HPLC): ≥75 area%

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

5-Aminoallyl-ddUTP powder, lithium salt

Aminoallyl-modified nucleotides are used for efficient target labeling and a wide variety of dyes and haptens for subsequent post-labeling reactions.

Application

5-Aminoallyl-11-ddUTP is used in 3'-end labeling reactions, serving as a substrate for:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)
- T4 and T7 DNA polymerase
- Taq DNA polymerase and reverse transcriptases

Lacking bulky dye groups, aminoallyl-modified nucleotides are incorporated at extremely high and reproducible levels. Reaction of the amine-modified nucleic acid with an excess of amine-reactive dye or hapten results in high and consistent labeling efficiencies, regardless of the dye or hapten chosen. The two-step labeling method eliminates the need to optimize the labeling reaction.

Benefits

- Generate stronger signals using efficient incorporation of amino groups for post-labeling.
- Label targets with high uniformity and consistency.
- Increase the options for target labels by post-labeling.
- Save time with ready-to-use solutions.

Cat. No.	Pack Size
11 858 114 103	custom fill

Will be supplied as "5'-Aminoallyl-2,3'-ddUTP, Li". Unit of Measure is "g".



For further processing only.

Properties

Nomenclature: 5-Aminoallyl-2',3'-dideoxy-uridine-5'-triphosphate

Formula: C₁₂H₁₆N₃O₁₃P₃Li₄

Molecular weight: 531 D

Specification

Appearance: White to yellowish, crystalline substance

5-Aminoallyl-ddUTP (1 μmol ± 9.73 A₂₄₀ units; phosphate buffer, 0.1 mol/l, pH 7.0): ≥85%

Purity (HPLC): ≥70 area%

Stability: At -15 to -25°C within specification range for 12 months. Store dry.

Biotin-11-CTP
10 mM solution

Biotin-labeled nucleotides are used for the efficient generation of biotinylated targets which can be subsequently captured using streptavidin coated solid phases or detected by streptavidin conjugates.

Application

Biotin-11-CTP is a substrate for SP6, T3, and T7 RNA Polymerase. It can replace CTP in the *in vitro* transcription reaction for RNA labeling. Double-stranded DNA (*i.e.*, linearized plasmid or cDNA) containing a T7, SP6, or T3 promoter serves as template for *in vitro* transcription with the corresponding RNA polymerase using ATP, GTP, CTP, UTP and Biotin-11-CTP. For a higher and/or more uniform label density Biotin-11-CTP can be used together with Biotin-16-UTP in the *in vitro* transcription reaction. The labeled RNA can be detected using a fluorescent streptavidin conjugate or by ELISA using a Streptavidin-AP conjugate.

Benefits

- Save time with ready-to-use solutions.

Properties

Nomenclature: Biotin-11-cytidine-5'-triphosphate

Formula: C₂₈H₄₄N₇O₁₇P₃SN₂

Molecular weight: 921.7 D

Specification

Appearance: Clear, colorless solution

Biotin-11-CTP (1 μmol ± 9.1 A₂₇₁ units, phosphate buffer, 0.1 mol/l, pH 7.0): 10.0-11.0 mmol/l

Purity (HPLC): 85-100 area%

Function test (incorporation using T7 Polymerase): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Biotin-16-UTP
10 mM solution

Biotin-labeled nucleotides are used for the efficient generation of biotinylated targets which can be subsequently captured using streptavidin coated solid phases or detected by streptavidin conjugates.

Application

Biotin-16-UTP is a substrate for SP6, T3, and T7 RNA Polymerase. It can replace UTP in the *in vitro* transcription reaction for RNA labeling. Linearized template DNA with T7, SP6 or T3 promoter is *in vitro* transcribed

Cat. No. **Pack Size**

04 762 924 103 custom fill

Will be supplied as "Biotin-11-CTP". Unit of Measure is "nmol".



For further processing only.

Cat. No. **Pack Size**

11 413 201 103 custom fill

Will be supplied as "Biotin-16-UTP". Unit of Measure is "nmol".



For further processing only.

with the corresponding RNA polymerases using ATP, GTP, CTP, UTP and Biotin-16-UTP, respectively.
The labeled RNA can be subsequently detected with a fluorescent streptavidin conjugate or by ELISA using a Streptavidin-AP conjugate.

Benefits

- Save time with ready-to-use solutions.

CAS: 86303-26-6

Properties

Nomenclature: Biotin-16-uridine-5'-triphosphate

Formula: C₃₂H₄₈N₇O₁₉P₃SLi₄

Molecular weight: 987.5 D

Specification

Appearance: Clear, colorless solution

Biotin-16-UTP (1 μmol ≙ 10.7 A₂₄₀ units, phosphate buffer, 0.1 mol/l, pH 7.0): 10.0-11.0 mmol/l

Purity (HPLC): ≥85 area%

Function test using DIG RNA Labeling Kit (SP6/T7): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Biotin-16-dUTP 20 mM solution

Biotin-labeled nucleotides are used for the efficient generation of biotinylated targets which can be subsequently captured using streptavidin coated solid phases or detected by streptavidin conjugates.

Application

Use Biotin-16-dUTP for nonradioactive DNA labeling such as random priming, PCR labeling or nick translation.

Biotin-16-dUTP is used as a substrate for:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)
- Taq DNA polymerase
- Reverse Transcriptase (*e.g.*, Transcriptor)

Biotin-16-dUTP replaces dTTP in the random-primed DNA labeling reaction or in nick translation in a ratio of 35% Biotin-16-dUTP and 65% dTTP, as well as in PCR.

The nucleotide also serves as a substrate for Terminal Transferase in 3'-end labeling.

Biotin labeled DNA can be detected with:

- Streptavidin-alkaline phosphatase conjugate and a chemiluminescent substrate (CSPD, CDP-*Star*) or a color substrate
- Biotin Luminescence Detection Kit

Benefits

- Save time with ready-to-use solutions.

Properties

Nomenclature: Biotin-16-2'-deoxy-uridine-5'-triphosphate

Formula: C₃₂H₄₈N₇O₁₈P₃SLi₄

Molecular weight: 971.5 D

Specification

Appearance: Clear, colorless solution

Cat. No. **Pack Size**

04 889 665 103 custom fill

Will be supplied as "Biotin-16-dUTP, 20 mM". Unit of Measure is "nmol".



For further processing only.

Biotin-16-dUTP (1 $\mu\text{mol} \triangleq 10.7 A_{240}$ units, phosphate buffer, 0.1 mol/l, pH 7.0):
20-25 mmol/l
Purity (HPLC): 85-100 area%
Function test using Biotin-High Prime (0.3 μg): Corresponds to reference
Stability: At -15 to -25°C within specification range for 12 months.

Biotin-16-dUTP 1 mM solution

Biotin-labeled nucleotides are used for the efficient generation of biotinylated targets which can be subsequently captured using streptavidin coated solid phases or detected by streptavidin conjugates.

Application

Use Biotin-16-dUTP for nonradioactive DNA labeling such as random priming, PCR labeling or nick translation.

Biotin-16-dUTP is used as a substrate for:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)
- Taq DNA polymerase
- Reverse Transcriptase (*e.g.*, Transcriptor)

Biotin-16-dUTP replaces dTTP in the random-primed DNA labeling reaction or in nick translation in a ratio of 35% Biotin-16-dUTP and 65% dTTP, as well as in PCR.

The nucleotide also serves as a substrate for Terminal Transferase in 3'-end labeling.

Biotin labeled DNA can be detected with:

- Streptavidin-alkaline phosphatase conjugate and a chemiluminescent substrate (CSPD, CDP-*Star*) or a color substrate
- Biotin Luminescence Detection Kit

Benefits

- Save time with ready-to-use solutions.

Properties

Nomenclature: Biotin-16-2'-deoxy-uridine-5'-triphosphate

Formula: $C_{32}H_{48}N_7O_{18}P_3SLi_4$

Molecular weight: 971.5 D

Specification

Appearance: Clear, colorless solution

Biotin-16-dUTP (1 $\mu\text{mol} \triangleq 10.7 A_{240}$ units, phosphate buffer, 0.1 mol/l, pH 7.0):
1.0-1.1 mmol/l

Purity (HPLC): 85-100 area%

Function test using Biotin-High Prime (0.3 μg): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Biotin-16-ddUTP 1 mM solution

Biotin-labeled nucleotides are used for the efficient generation of biotinylated targets which can be subsequently captured using streptavidin coated solid phases or detected by streptavidin conjugates.

Application

Biotin-16-ddUTP is used as a substrate for:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)

Cat. No.	Pack Size
11 093 711 103	custom fill

Will be supplied as "Biotin-16-dUTP Solution". Unit of Measure is "nmol".



For further processing only.

Cat. No.	Pack Size
11 431 692 103	custom fill

Will be supplied as "Biotin-16-ddUTP". Unit of Measure is "nmol".



For further processing only.

- T4 DNA polymerase
 - Taq DNA polymerase and reverse transcriptase (*e.g.*, Transcriptor)
- Biotin-16-ddUTP is preferentially used for 3'-end labeling of oligonucleotides with Terminal Transferase, recombinant.
- The biotin-labeled oligonucleotide can be used as a hybridization probe for:

- DNA and RNA transfers
- Colony and plaque screening
- *In situ* hybridization

In addition to common hybridization techniques, biotin-labeled oligonucleotides are specifically useful for screening expression libraries for sequence-specific DNA-binding proteins, for example, transcription factors.

The labeled oligomer can be subsequently detected by ELISA using the Streptavidin-AP conjugate for nucleic acid detection.

Oligonucleotides are enzymatically labeled at their 3' end with Terminal Transferase by incorporation of a single biotin-labeled dideoxyuridine-triphosphate (biotin-ddUTP).

Benefits

- Save time with ready-to-use solutions.

CAS: 422268-45-9

Properties

Nomenclature: Biotin-16-2',3'-dideoxy-uridine-5'-triphosphate

Formula: C₃₂H₄₈N₇O₁₇P₃SLi₄

Molecular weight: 955.5 D

Specification

Appearance: Clear, colorless solution

Biotin-16-ddUTP (1 μmol ≅ 10.7 A₂₄₀ units, phosphate buffer, 0.1 mol/l, pH 7.0): 1.0-1.1 mmol/l

Purity (HPLC): 85-100 area%

Function test using DIG Oligonucleotide 3'-End Labeling Kit, 2nd generation (DIG is replaced by Biotin-16-ddUTP): Corresponds to reference

Stability: At -15 to -25°C within specification range for 24 months.

Digoxigenin-11-UTP 10 mM solution

Digoxigenin-11-UTP is used to label RNA with digoxigenin. DIG-labeling permits fast and sensitive detection without the need of radioactive material.

Application

Use DIG-11-UTP as a substrate for SP6, T3, and T7 RNA polymerases. It replaces UTP in the *in vitro* transcription reaction for DIG labeling of RNA in a ratio of 35:65%. Linearized template DNA with T7, SP6, or T3 promoter is *in vitro* transcribed with the corresponding RNA polymerases using ATP, GTP, CTP, UTP and DIG-11-UTP, respectively. Labeled RNA can be detected using Anti-Digoxigenin-AP or Fab fragments.

Benefits

- **Eliminate radioactive material from your laboratory.** Only nonhazardous material is required for your assays.
- **Obtain results faster.** The required exposure time is much shorter compared to radioactive assays.
- **Profit from a stable probe.** DIG-labeled probes are stable for more

Cat. No.	Pack Size
11 230 409 103	custom fill

Will be supplied as "DIG-11-UTP". Unit of Measure is "nmol". For other than the listed concentrations please inquire.



For further processing only.

than a year and can be easily stripped.

- **Take advantage of the large number of published protocols.**

CAS: 186033-10-3

Properties

Nomenclature: Digoxigenin-11-uridine-5'-triphosphate

Formula: C₄₃H₆₁N₄O₂₂P₃Li₄

Molecular weight: 1106.7 D

pH stability: Digoxigenin is bound to the nucleotide portion via an alkali-resistant ether linkage. The preparation is stable to 0.1-0.5 M NaOH at +15 to +25°C.

Specification

Appearance: Clear, colorless solution

Digoxigenin-11-UTP (1 μmol ± 8.78 A₂₉₀ units, phosphate buffer, 0.1 mol/l, pH 7.0): 10.0-11.0 mmol/l

Purity (HPLC): 85 area%

Function test using Dig RNA Labeling Kit: Corresponds to reference

Stability: At -15 to -25°C within specification range for 18 months.

Quality

Function tested with the DIG RNA Labeling Kit.

Digoxigenin-11-UTP 3.5 mM solution

Digoxigenin-11-UTP is used to label RNA with digoxigenine. DIG-labeling permits fast and sensitive detection without the need of radioactive material.

Application

Use DIG-11-UTP as a substrate for SP6, T3, and T7 RNA polymerases. It replaces UTP in the *in vitro* transcription reaction for DIG labeling of RNA in a ratio of 35:65%. Linearized template DNA with T7, SP6, or T3 promoter is *in vitro* transcribed with the corresponding RNA polymerases using ATP, GTP, CTP, UTP, and DIG-11-UTP, respectively. Labeled RNA can be detected using Anti-Digoxigenin-AP or fab fragments.

Benefits

- **Eliminate radioactive material from your laboratory.** Only nonhazardous material is required for your assays.
- **Obtain results faster.** The required exposure time is much shorter compared to radioactive assays.
- **Profit from a stable probe.** DIG-labeled probes are stable for more than a year and can be easily stripped.
- **Take advantage from the large number of published protocols.**

Properties

Nomenclature: Digoxigenin-11-uridine-5'-triphosphate

Formula: C₄₃H₆₁N₄O₂₂P₃Li₄

Molecular weight: 1106.7 D

pH stability: Digoxigenin is bound to the nucleotide portion via an alkali-resistant ether linkage. The preparation is stable to 0.1-0.5 M NaOH at +15 to +25°C.

Specification

Appearance: Clear, colorless solution

Digoxigenin-11-UTP (1 μmol ± 8.78 A₂₉₀ units, phosphate buffer, 0.1 mol/l, pH 7.0): 3.3-3.7 mmol/l

Purity (HPLC): 85-100 area%

Cat. No.	Pack Size
03 359 239 103	custom fill

Will be supplied as "Digoxigenin-11-UTP". Unit of Measure is "nmol".

For other than the listed concentrations please inquire.



For further processing only.

Function test using DIG RNA Labeling Kit: Corresponds to reference

Stability: At -15 to -25°C within specification range for 18 months.

Quality

Function tested with the DIG RNA Labeling Kit.

Digoxigenin-11-dUTP, alkali-stable 1 mM solution

Digoxigenin-11-dUTP, alkali-stable, replaces radioactive assays for even higher sensitivity. The required exposure time is drastically reduced compared to radioactive assays.

Application

Use Digoxigenin-11-dUTP, alkali-stable, for nonradioactive DNA labeling (*e.g.*, random priming, PCR labeling, tailing, or nick translation). DIG-11-dUTP replaces dTTP in the random-primed DNA labeling reaction or in nick translation in a ratio of 35% DIG-11-dUTP and 65% dTTP. It is ideal as a substrate for DNA polymerase, Taq DNA polymerase, Terminal Transferase and reverse transcriptase.

Note: For labeling of probes which are preferentially used in hybridization experiments where stripping and reprobing of the membrane is intended, use DIG-11-dUTP, alkali-labile.

Benefits

- **Eliminate radioactive material from your laboratory.** Only nonhazardous material is required for your assays.
- **Obtain results faster.** The required exposure time is much shorter compared to radioactive assays.
- **Profit from a stable probe.** DIG-labeled probes are stable for more than a year and can be easily stripped.
- **Take advantage of the large number of published protocols.**

CAS: 1173-82-6

Properties

Nomenclature: Digoxigenin-11-2'-deoxy-uridine-5'-triphosphate

Formula: C₄₃H₆₁N₄O₂₁P₃Li₄

Molecular weight: 1090.7 D

Specification

Appearance: Clear, colorless solution

Digoxigenin-11-dUTP (1 μmol ± 8.78 A₂₉₀ units, phosphate buffer, 0.1 mol/l, pH 7.0): 0.9-1.1 mmol/l

Purity (HPLC): 85-100 area%

Function test using DIG DNA Labeling and Detection Kit: Corresponds to reference

Stability: At -15 to -25°C within specification range for 18 months.

Quality

Function tested by RPL.

Cat. No. **Pack Size**

11 093 681 103 custom fill

Will be supplied as "DIG-11-dUTP". Unit of Measure is "nmol".
For other than the listed concentrations please inquire.



For further processing only.

Digoxigenin-11-dUTP, alkali-labile 1 mM solution

Digoxigenin-11-dUTP, alkali-labile, is a replacement for radioactive assays for even higher sensitivity. The DIG-label can be cleaved under basic conditions resulting in an unlabeled probe for further processing.

Application

Use Digoxigenin-11-dUTP, alkali-labile, for nonradioactive DNA labeling (e.g., random-primed or nick translation). The alkali-labile compound should be used for labeling of probes which are preferentially used in hybridization experiments where stripping and reprobing of the membrane is intended. DIG-11-dUTP replaces dTTP in the random-primed DNA labeling reaction or in nick translation in a ratio of 35% DIG-11-dUTP and 65% dTTP. It is ideal as a substrate for DNA polymerase, Taq DNA polymerase, Terminal Transferase and reverse transcriptase.

Note: Do not use in experiments where alkaline treatment is required. For this application, use DIG-11-dUTP, alkali-stable.

Benefits

- **Eliminate radioactive material from your laboratory.** Only nonhazardous material is required for your assays.
- **Obtain results faster.** The required exposure time is much shorter compared to radioactive assays.
- **Profit from a stable probe.** DIG-labeled probes are stable for more than a year and can be easily stripped.
- **Obtain the unlabeled fragment after detection.** The DIG-label can be cleaved under basic conditions resulting in an unlabeled probe for further processing.
- **Take advantage of the large number of published protocols.**

CAS: 1173-82-6

Properties

Nomenclature: Digoxigenin-11-2'-deoxy-uridine-5'-triphosphate

Formula: C₄₅H₆₃N₄O₂₂P₃Li₄

Molecular weight: 1132.7 D

Specification

Appearance: Clear, colorless solution

Digoxigenin-11-dUTP (1 μmol ± 8.78 A₂₉₀ units, phosphate buffer, 0.1 mol/l, pH 7.0): 0.9-1.1 mmol/l

Purity (HPLC): 85-100 area%

Function test using DIG DNA Labeling and Detection Kit: Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Function tested by RPL.

Digoxigenin-11-ddUTP 1 mM solution

Digoxigenin-11-ddUTP permits replacement of radioactive assays for even higher sensitivity. The required exposure time is drastically reduced compared to a radioactive assay.

Application

Use Digoxigenin-11-ddUTP as a substrate for: Terminal Transferase, DNA polymerase I (holoenzyme and Klenow fragment), T4 and T7 DNA polymerase

Cat. No. **Pack Size**

11 579 541 103 custom fill

Will be supplied as "DIG-11-dUTP, Alkali-labile, Solution". Unit of Measure is "nmol".

For other than the listed concentrations please inquire.



For further processing only.

Cat. No. **Pack Size**

11 365 207 103 custom fill

Will be supplied as "DIG-11-ddUTP". Unit of Measure is "nmol".
For other than the listed concentrations please inquire.



For further processing only.

or Taq DNA polymerase and reverse transcriptase (e.g., Transcriptor). It is preferentially used for 3'-end labeling of oligonucleotides with Terminal Transferase, recombinant. The DIG-labeled oligonucleotide can be used as a hybridization probe for:

- DNA and RNA transfers
- Colony and plaque screening
- *In situ* hybridization

In addition to common hybridization techniques, DIG labeled oligonucleotides are specifically useful for screening expression libraries for sequence specific DNA-binding proteins such as transcription factors.

Benefits

- **Eliminate radioactive material from your laboratory.** Only nonhazardous material is required for your assays.
- **Obtain results faster.** The required exposure time is much shorter compared to radioactive assays.
- **Profit from a stable probe.** DIG-labeled probes are stable for more than a year and can be easily stripped.
- **Take advantage of the large number of published protocols.**

Properties

Nomenclature: Digoxigenin-11-2'-dideoxy-uridine-5'-triphosphate

Formula: $C_{43}H_{61}N_4O_{20}P_4Li$

Molecular weight: 1074.7 D

Specification

Appearance: Clear, colorless solution

Digoxigenin-11-ddUTP (1 μ mol \triangleq 8.78 A_{290} units, phosphate buffer, 0.1 mol/l, pH 7.0): 1.0-1.1 mmol/l

Purity (HPLC): 85-100 area%

Function test using DIG Oligonucleotide 3'-End Labeling Kit, 2nd

generation: Corresponds to reference

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Function tested by end labeling.

Fluorescein-12-UTP 10 mM solution

Enzymatic nonradioactive labeling reagent for *in vitro* transcription reactions.

Application

Use Fluorescein-12-UTP to add a nonradioactive label to RNA during *in vitro* transcription. The labeled RNA can easily and safely be detected either directly or with an enzyme-conjugated anti-fluorescein. Fluorescein-12-UTP can replace UTP as a substrate for the following enzymes:

- SP6 RNA polymerase
- T3 RNA polymerase
- T7 RNA polymerase

Labeled RNA can be subsequently detected by *in situ* hybridization and direct fluorescence detection and detection by ELISA using Anti-Fluorescein-AP, Fab fragments or Anti-Fluorescein-POD, Fab fragments.

Benefits

- **Save time** with a ready-to-use solution for *in situ* hybridization and ELISA.

Cat. No.	Pack Size
11 431 706 103	custom fill

Will be supplied as "Fluorescein-12-UTP". Unit of Measure is "nmol".

For other than the listed concentrations please inquire.



For further processing only.

CAS: 134367-01-4

Properties

Nomenclature: Fluorescein-12-uridine-5'-triphosphate

Formula: C₃₉H₃₇N₄O₂₂P₃Li₄

Molecular weight: 1034.4 D

Specification

Appearance: Clear, yellow solution

Fluorescein-12-UTP (1 μmol ± 63.3 A₄₉₅ units, phosphate buffer, 0.1 mol/l, pH 9.0): 10.0-11.0 mmol/l

Purity (HPLC, including isomere): >85-100 area%

Function test using DIG RNA Labeling Kit (SP6/T7) (DIG is replaced by Fluorescein-12-UTP): Corresponds to reference

Ribonucleases: Negative

Stability: At -15 to -25°C within specification range for 12 months. Store dry and protect from light.

Fluorescein-12-dUTP

1 mM solution

Enzymatic nonradioactive labeling reagent for cDNA synthesis, PCR, random primed labeling, nick-translation or primer extension.

Application

Use Fluorescein-12-dUTP as a substrate for:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)
- Taq DNA polymerase
- Reverse transcriptase (*e.g.*, from AMV and M-MuLV)

Fluorescein-12-dUTP replaces dTTP in the random-primed DNA labeling reaction or in nick translation reactions, as well as in PCR. The nucleotide also serves as a substrate for Terminal Transferase in 3'-end labeling.

Fluorescein-labeled probes can be used for *in situ* hybridization with direct fluorescence detection and detection by ELISA using Anti-Fluorescein-AP, Fab fragments.

Repeated fluorescence labeling using Tetramethylrhodamine-6-dUTP (red) and AMCA-6-dUTP (bright blue) is possible.

Benefits

- **Save time** with a ready-to-use solution for *in situ* hybridization and ELISA.
- **Be flexible.** The Fluorescein-12-dUTP labeled samples can be co-labeled using either Tetramethylrhodamine-dUTP or AMCA-6-dUTP.

CAS: 134344-32-4

Properties

Nomenclature: Fluorescein-12-deoxyuridine-5'-triphosphate

Formula: C₃₉H₃₇N₄O₂₁P₃Li₄

Molecular weight: 1018.4 D

Specification

Appearance: Clear, yellow solution

Fluorescein-12-dUTP (1 μmol ± 63.3 A₄₉₅ units, phosphate buffer, 0.1 mol/l, pH 9.0): 1.0-1.1 mmol/l

Purity (HPLC, including isomere): >85-100 area%

Function test using In Situ Cell Death Detection Kit, Fluorescein (*in situ* hybridization): Corresponds to reference

Stability: At -15 to -25°C within specification range for 12 months. Store dry and protect from light.

Cat. No. **Pack Size**

11 375 601 103 custom fill

Will be supplied as "Fluorescein-12-dUTP". Unit of Measure is "nmol".

For other than the listed concentrations please inquire.



For further processing only.

Tetramethylrhodamine-5-dUTP**1 mM solution, lithium salt**

Tetramethylrhodamine-5-dUTP is a substitute for dTTP in nick-translation and in the random-primed labeling reactions.

Application

Tetramethylrhodamine-5-dUTP is used for nonradioactive labeling of DNA. This modified nucleotide is a substrate for:

- Terminal Transferase
- DNA polymerase I (holoenzyme and Klenow fragment)
- Taq DNA polymerase
- Reverse transcriptase (e.g., Transcriptor Reverse Transcriptase and other reverse transcriptases)

Benefits

- Benefit from a hazard free labeling and detection system without the risks associated with radioactive assays.
- Profit from directly measuring the fluorescence for *in situ* hybridization.

Product Description

Tetramethylrhodamine-5-dUTP is a substitute for dTTP in nick-translation reactions and in the random-primed labeling technique for DNA labeling, as well as in PCR. The nucleotide also serves as a substrate for Terminal Transferase in 3'-end labeling. Tetramethylrhodamine-labeled probes show red fluorescence and are suitable for use in *in situ* hybridization for direct fluorescence detection. Multiple fluorescence labeling using Fluorescein-12-dUTP (yellow fluorescence) or other dye-labeled deoxynucleotides is possible.

Properties

Formula: C₃₇H₃₆N₅O₁₈P₃Li₄

Molecular weight: 959.4 D

Specification

Appearance: Clear, red solution

Tetramethylrhodamine-5-dUTP (1 μmol ± 70.0 A₅₅₁ units, phosphate buffer, 0.1 mol/l, pH 9.0): 1.0-1.1 mmol/l

Purity (HPLC): 85.0-100.0 area%

Function test (*in situ* hybridization): Corresponds to specification

Stability: At -15 to -25°C within specification range for 12 months. Protect from light.

Cat. No.**11 542 907 103****Pack Size**

custom fill

Will be supplied as "Tetramethylrhodamine-5-dUTP". Unit of Measure is "nmol".



For further processing only.

COT Human DNA, CGH Grade

from human male placenta DNA, enriched for repetitive sequences, solution

Obtain best results for *in situ* suppression (CISS) hybridization, DNA microarray application and many other hybridization applications. With the COT Human DNA, CGH Grade, reproducible and sensitive measurements of dsDNA concentration are possible.

Application

Use the COT Human DNA in the following applications:

- Nucleic acid labeling and detection
- DNA microarray applications such as comparative genome hybridization and sequence capture
- Complex hybridization of human nucleic acids like FISH

In microarray applications, COT Human DNA is used in hybridization solutions to block repetitive DNA sample sequences from nonspecific hybridizations. In filter and other hybridization techniques, COT Human DNA is also used in prehybridization solutions to inactivate nonspecific target binding sites.

Benefits

- Ensure specificity by blocking unspecific nucleic acid motifs using high concentrations of double-stranded, repetitive (COT) human DNA.
- Achieve reproducible, reliable results using fluorometrically quantified COT human DNA, when hybridizing to complex dsDNA.
- Use COT Human DNA that has been pretested using DNA microarrays and comparative genome hybridization (CGH).
- Suppress the background noise effectively in your experiments using high concentrations of repetitive dsDNA.

Product Description

The COT fraction of human genomic DNA consists largely of rapidly annealing repetitive elements. These interspersed repetitive sequences (IRS) such as SINEs (small interspersed repetitive elements, *e.g.*, Alu-elements) and LINEs (large interspersed repetitive elements, *e.g.*, L1-elements) are distributed ubiquitously throughout the genome. COT Human DNA is prepared from human placental DNA by shearing, denaturing, and reannealing under conditions that enrich these repetitive elements.

CAS: 99675-55-5

Specification

Appearance: Clear, colorless solution

COT Human DNA (A_{260} , water, 1 AB=50 µg/ml): ≥1.0 mg/ml

Fluorometrical determination of concentration: 1.0-1.5 mg/ml

Y-Chromosom (recovered exclusively from male human placenta): Corresponds to specification

A_{260}/A_{280} : 1.6-2.0

Absence of HIV 1/2 and HCV/HBV: Corresponds to specification

Performance test using gel electrophoretic separation (4% agarose gel without RE cleavage): middle chain length: 50-300 bp

Comparable intensity to previous lot: Corresponds to specification

Function test in the CGH array: Corresponds to specification

Stability: At -15 to -25°C within specification range for 18 months.

Quality

The product is HIV tested.

Background Information

Repetitive elements (IRS) present in a probe (*e.g.*, cosmids, YACs, chromosome painting probes) generate nonspecific hybridization signals that are distributed

Cat. No. **Pack Size**

05 111 854 103 custom fill

Will be supplied as "COT Human DNA, CGH Grade". Unit of Measure is "ml".



For further processing only.

over the whole chromosome or genome. To enable specific hybridization of the probe to the chromosomal target site (*e.g.*, single-copy sequences or low-copy repeats), the probe must be denatured in the presence of excess unlabeled COT Human DNA. This DNA serves as a competitor. In a subsequent preannealing step, the repetitive probe elements rapidly hybridize to excess repeats in the Cot Human DNA, while most of the specific probe sequences remain single stranded, and can thus hybridize to their chromosomal targets. This technique is known as chromosomal *in situ* suppression (CISS) hybridization.

COT Human DNA from human placenta DNA, enriched for repetitive sequences, solution

For suppression of cross-hybridization to human repetitive sequences in filter and *in situ* hybridizations.

Application

COT Human DNA is used in chromosome *in situ* suppression (CISS) hybridization. Cosmid or YAC probes contain repetitive elements that result in monospecific hybridization signals distributed over the entire chromosome. To enable specific hybridization to the chromosomal target site, the probe is denatured together with an excess of unlabeled COT Human DNA as a competitor. COT Human DNA can be used to suppress nonspecific hybridization to human repetitive sequences in microarray analysis, and in filter and *in situ* hybridization experiments.

Benefits

- Rely on the effective suppression due to the high concentration of repetitive dsDNA.

Product Description

The COT fraction of human genomic DNA consists largely of rapidly annealing repetitive elements. These interspersed repetitive sequences (IRS) such as SINEs (small interspersed repetitive elements, *e.g.*, Alu-elements) and LINEs (large interspersed repetitive elements, *e.g.*, L1-elements) are distributed ubiquitously throughout the genome. COT Human DNA is prepared from human placental DNA by shearing, denaturing, and reannealing under conditions that enrich these repetitive elements.

CAS: 99675-55-5

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 10 mmol/l; EDTA, 1 mmol/l; pH 7.4

COT Human DNA (A_{260} , water): 1.0-2.0 mg/ml

Y-Chromosom (recovered exclusively from male human placenta): Corresponds to specification

A_{260}/A_{280} : 1.6-2.0

Absence of HIV 1/2 and HCV/HBV: Corresponds to specification

Performance test using gel electrophoretic separation (4% agarose gel, without RE cleavage): middle chain length: 50-300 bp

Comparable intensity of Cot 1 DNA to masterlot: Corresponds to specification

Stability: At -15 to -25°C within specification range for 18 months.

Quality

The product is HIV tested.

Cat. No. **Pack Size**

11 582 011 103 custom fill

Will be supplied as "DNA,Cot-1,human". Unit of Measure is "mg".



For further processing only.

DNA from fish sperm, solution

This DNA preparation can serve for prevention of nonspecific binding in hybridization experiments.

Application

Use this preparation of single-stranded genomic DNA fragments to prevent nonspecific binding in membrane or *in situ* DNA hybridization experiments. It can be added directly to the hybridization mix with no need for sonification or denaturation.

CAS: 9007-49-2

Properties

The DNA is sonicated. The length of the DNA fragments is mostly in the range of 50 to 600 bp (not verified). UV absorption maximum is at 258 nm.

Specification

Appearance: White to slightly grey lyophilizate

DNA content (A_{260}): ≥ 15 AB/mg lyophilizate

DNA content (based on P_{organic}): $\geq 70\%$ of lyophilizate

P_{organic} ($P_{\text{total}} - P_{\text{p}}$): ≥ 60 $\mu\text{g}/\text{mg}$ lyophilizate

P: ≤ 3 $\mu\text{g}/\text{mg}$ lyophilizate

Na (flame photometric): $6 \pm 2\%$

Protein content (Lowry): ≤ 50 $\mu\text{g}/\text{mg}$ lyophilizate

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No.	Pack Size
10 223 638 103	custom fill

Will be supplied as "DNA, Sodium Salt from fish sperm". Unit of Measure is "g".

For further processing only.

Poly [d(A-T)] powder

Poly[d(A-T)] is a suitable template for RNA synthesis.

Application

Use Poly[d(A-T)] as a template for RNA polymerases.

Specification

Appearance: White lyophilizate

Melting range: 62.0-69.0°C

Content (A_{260} , content of one vial resolved in 1 ml water): 50-60 OD/ml

Mean strand length (gel electrophoresis): 500-100,000 bp

Ratio dA/dT: $1:1 \pm 10\%$

A_{250}/A_{260} : 0.70-0.76

A_{280}/A_{260} : 0.38-0.58

A_{290}/A_{260} : 0.09-0.16

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
11 336 312 103	50 OD[260] units

Will be supplied as "Poly[d(A-T)], Sodium Salt". Unit of Measure is "piece".



For further processing only.

RNA**from yeast, powder, free acid**

This product is a preparation of total RNA.

Application

Use this RNA preparation for studies which use natural RNA in a *in vivo* and *in vitro* protein-synthesizing system. It can also be used as carrier RNA in *in situ* hybridization experiments.

CAS: 63231-63-0

Properties

Total RNA from *Saccharomyces cerevisiae*.

Specification

Appearance: Yellowish to brown powder

Purity (A_{260}): $\geq 95.0\%$

P_i: $\leq 0.3\%$

A_{250}/A_{260} : 0.86-0.92

A_{280}/A_{260} : 0.44-0.48

A_{290}/A_{260} : 0.15-0.19

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No.**Pack Size**

10 153 320 103 custom fill

Will be supplied as "RNA from Yeast". Unit of Measure is "g".
For further processing only.

Poly(A)**potassium salt, solution**

Poly(A) supports precipitation of DNA and RNA.

Application

Use Poly(A) as a carrier for quantitative precipitation of DNA and RNA.

CAS: 26763-19-9

Properties

Poly(A) is a suitable carrier for quantitative precipitation of DNA and RNA, especially to improve recovery of low amounts of nucleic acid or of short fragments <200 bp.

Molecular weight: 100-500 kD

Specification

Appearance: Clear, colorless solution

pH value: 6.5±0.5

Content (A_{257}): 4.5-5.5 mg/ml

Mean strand length (gel electrophoresis): 3,000-10,000 nucleotides

Ribonucleases (fluorescence polarisation): Negative

A_{250}/A_{260} : 0.86-0.90

A_{280}/A_{260} : 0.28-0.32

A_{290}/A_{260} : 0.03-0.05

Stability: At -15 to -25°C within specification range for 36 months.

Cat. No.**Pack Size**

12 159 074 103 custom fill

Will be supplied as "Poly (A) Potassium Salt, Solution". Unit of Measure is "g".



For further processing only.

Poly(A) potassium salt, lyophilizate

Poly(A) supports precipitation of DNA and RNA.

Application

Use Poly(A) as a carrier for quantitative precipitation of DNA and RNA.

CAS: 26763-19-9

Properties

Poly(A) is a suitable carrier for quantitative precipitation of DNA and RNA, especially to improve recovery of low amounts of nucleic acid or of short fragments <200 bp.

Molecular weight: 100-500 kD

Specification

Appearance: White lyophilizate

Content (A_{257}): 2.1 $\mu\text{mol/mg}$

Mean strand length (gel electrophoresis): 2,100-10,000 nucleotides

Ribonucleases (fluorescence polarization): Negative

A_{250}/A_{260} : 0.86-0.90

A_{280}/A_{260} : 0.28-0.32

A_{290}/A_{260} : 0.03-0.05

Stability: At +2 to +8°C within specification range for 36 months.

Cat. No.	Pack Size
10 041 955 103	1 g

Will be supplied as "Polyadenylic Acid, Poly (A), K-Salt". Unit of Measure is "g".

For further processing only.

Glycogen, Molecular Biology Grade from mussels, solution

Glycogen, Molecular Biology Grade, supports precipitation of DNA and RNA.

Application

Use Glycogen, Molecular Biology Grade, as a carrier for quantitative precipitation of DNA and RNA.

Benefits

- **Gain excellent performance.** Take advantage of a contamination-controlled reagent with high lot-to-lot consistency.

CAS: 9005-79-2

Properties

Glycogen, Molecular Biology Grade, is a suitable carrier for DNA and RNA in ethanol precipitation and phenol/chloroform extraction, especially to increase sensitivity with low amounts of total nucleic acid. In contrast to carrier DNA and RNA, glycogen is inert in nucleic acid modifying processes. It has no influence on enzymatic treatment of nucleic acids or on gel electrophoresis. Glycogen does not bind to nucleic acids, and can be easily removed by gel electrophoresis or gel filtration.

Specification

Appearance: Clear, colorless solution

Concentration: ≥20 mg/ml

Unspecific endonucleases (λDNA and MWM III DNA): Not detectable in up to 200 µg after 4 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 200 µg after 4 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 200 µg after 4 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 200 µg after 4 hours incubation at +37°C.

Proteinases (colorimetric): Not detectable in up to 200 µg after 2 hours incubation at +37°C.

Nucleic acid (gel electrophoresis): Not detectable in up to 200 µg

Stability: At -15 to -25°C within specification range for 36 months.

Cat. No. _____ **Pack Size** _____

10 899 232 103 custom fill

Will be supplied as "Glycogen, from Mussels SQ for Mol.Biol.". Unit of Measure is "g".



For further processing only.

Alkaline Phosphatase, recombinant, 1 U/μl from bovine intestine, expressed in *Pichia pastoris*, solution

Alkaline Phosphatase for 5' dephosphorylation of DNA and RNA.

Application

Use Alkaline Phosphatase, recombinant, 1 U/μl for:

- Dephosphorylation of 5' phosphate from DNA and RNA
- Coupling to other proteins via its amino or carbohydrate groups, for detection in immunoassays and western blot analysis. The increased heat instability of the recombinant enzyme allows quick deactivation by heating after dephosphorylation

Benefits

- **Rely on fast end efficient performance.** Recombinant AP allows fast dephosphorylation and quick deactivation.
- **Gain excellent performance.** Take advantage of this highly processive Alkaline Phosphatase, recombinant, 1 U/μl, specifically tested for molecular biology applications.

EC 3.1.3.1

Properties

Enzyme activity: Alkaline Phosphatase, recombinant, 1 U/μl catalyzes the dephosphorylation of 5' phosphate from DNA and RNA.

pH activity optimum: 9.8

pH stability optimum: 8.0

Cofactor: Zn²⁺

Aktivators: Mg²⁺, Mn²⁺, Co²⁺

Inactivation: Complete inactivation after 5 minutes at +75°C.

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 25 mmol/l; MgCl₂, 1 mmol/l; ZnCl₂, 0.1 mmol/l; glycerol, 50% (v/v), pH approximately 7.6 at +4°C

Volume activity: ≥1 U/μl

Specific activity: ≥5 kU/mg protein

Unit determination: Photometric test with 4-nitrophenyl phosphate

Unspecific endonucleases (λDNA): Not detectable in up to 100 U after 4 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 100 U after 4 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 60 U after 4 hours incubation at +37°C. **Ribonucleases (MS2 RNA):** Not detectable in up to 100 U after 1 hour incubation at +37°C. **Function test, 5' labeling (MWM III DNA, incorporation rate ≥30%):** Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
04 571 550 103	custom fill

Will be supplied as "Alkaline Phosphatase rec., MB Grade, 1 U".
Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

- For the best fit reaction buffer, use Dephosphorylation Buffer for Alkaline Phosphatase, see page 331

Alkaline Phosphatase, recombinant, 20 U/μl

from bovine intestine, expressed in *Pichia pastoris*,
solution

Alkaline Phosphatase for 5' dephosphorylation of DNA and RNA.

Application

Use Alkaline Phosphatase, recombinant for:

- Dephosphorylation of 5' phosphate from DNA and RNA
- Coupling to other proteins via its amino or carbohydrate groups, for detection in immunoassays and western blot analysis

Benefits

- **Gain excellent performance.** Take advantage of this highly processive, contamination-controlled Alkaline Phosphatase, recombinant.

EC 3.1.3.1

Properties

Enzyme activity: Alkaline Phosphatase, recombinant, 1 U/μl catalyzes the dephosphorylation of 5' phosphate from DNA and RNA.

pH activity optimum: 9.8

pH stability optimum: 8.0

Cofactor: Zn²⁺

Aktivators: Mg²⁺, Mn²⁺, Co²⁺

Inactivation: Complete inactivation after 5 minutes at +75°C.

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 25 mmol/l; MgCl₂, 1 mmol/l; ZnCl₂, 0.1 mmol/l; glycerol, 50% (v/v); pH approximately 7.6 at +4°C

Volume activity: ≥20 U/μl

Specific activity: ≥5 kU/mg protein

Unit determination: Photometric test with 4-nitrophenyl phosphate.

Unspecific endonucleases (λDNA): Not detectable in up to 100 U after 4 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 100 U after 4 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 60 U after 4 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 100 U after 1 hour incubation at +37°C.

Function test, 5' labeling (MWM III DNA, incorporation rate ≥30%): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No. **Pack Size**

04 571 363 103 custom fill

Will be supplied as "Alkaline Phosphatase rec., MB Grade, 20 U".
Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

- For the best fit reaction buffer, use Dephosphorylation Buffer for Alkaline Phosphatase, see page 331

Dephosphorylation Buffer for Alkaline Phosphatase

10x concentrated

Standard reaction buffer for Alkaline Phosphatase.

Application

Use Dephosphorylation Buffer for Alkaline Phosphatase as optimized reaction buffer for Alkaline Phosphatase, recombinant, 1 U/ μ l.

Benefits

- **Simplify reaction setup.** Use this premixed, pH-adjusted and contamination-controlled standard reaction buffer for fast and easy setup of highly reproducible experiments.
- **Gain excellent performance.** Take full advantage of this Dephosphorylation Buffer for Alkaline Phosphatase, using it as reaction buffer specially optimized for Alkaline Phosphatase, recombinant, 1 U/ μ l.

Specification

Appearance: Clear, colorless solution

Contents: 10x concentrated solution, comprising Tris/HCl, 500 mmol/l; EDTA, 1 mmol/l; pH approximately 8.5 at +20°C

Unspecific endonucleases (λ DNA): Not detectable in up to 20 μ l after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 20 μ l after 16 hours incubation at +37°C.

Ribonucleases (MWM II RNA): Not detectable in up to 20 μ l after 1 hour incubation at +50°C.

Performance test (5' labeling on MWM III DNA, incorporation rate \geq 30%): Corresponds to reference

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No.	Pack Size
05 989 639 103	1 ml

Will be supplied as "Dephosphorylation Buffer, AP (10x), 1 ml".

Unit of Measure is "piece".

For further processing only.

T4 DNA Ligase**recombinant form of the enzyme from T4 phage, solution**

Use T4 DNA Ligase for ligation of DNA fragments.

Application

Use T4 DNA Ligase to ligate DNA fragments with blunt or overlapping ends.

Benefits

- **Obtain consistent performance.** Take advantage of this highly processive, contamination-controlled T4 DNA Ligase.

EC 6.5.1.1

Properties

Enzyme activities: T4 DNA Ligase catalyzes the formation of phosphodiester bonds between neighbouring 3'-hydroxyl and 5'-phosphate ends in double-stranded DNA. Sticky- and blunt-ended DNA fragments are ligated. Single-stranded nicks in double-stranded DNA are also closed.

Appropriate ligation buffer, 10x concentrated: Tris/HCl, 660 mmol/l; MgCl₂, 50 mmol/l; DTT, 50 mmol/l; ATP, 10 mmol/l, pH 7.5 at +20°C (Note: ATP is not stable).

pH optimum: 7.2-7.8

Divalent ion requirement: Mg²⁺

Inactivation: After 10 minutes heat denaturation at +65°C ligase activity is stopped.

Specification

Appearance: Clear, colorless solution

Storage buffer: Tris/HCl, 20 mmol/l; KCl, 60 mmol/l; DTE, 5 mmol/l; EDTA, 1 mmol/l; glycerol, 50% (v/v); pH approximately 7.5 at +4°C

Volume activity: ≥5 U/μl

Unit definition: One unit T4 DNA Ligase is defined as the amount of enzyme which converts 1 nmol of [³²P] from pyrophosphate into Norit-absorbable material in 20 minutes at +37°C.

Glycosylases (M13mp11(U) ssDNA): Not detectable in up to 10 U after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 10 U after 16 hours incubation at +37°C.

Exonucleases (³H-DNA): Not detectable in up to 15 U after 4 hours incubation at +37°C.

Animal-derived additives: None

Stability: At -15 to -25°C within specification range for 18 months.

Cat. No.

Pack Size

10 909 246 103 custom fill

Will be supplied as "DNA Ligase, T4 from T4-infected EcoNM989".
Unit of Measure is "kU".

The enzyme is supplied without reaction buffer.



For further processing only.

Pyrophosphatase, inorganic (PPase) from yeast

Pyrophosphatase, inorganic (PPase) catalyzes the hydrolysis of inorganic pyrophosphate to form orthophosphate.

Application

Pyrophosphatase, inorganic (PPase) enhances enzymatic reactions producing pyrophosphate, since these reactions typically have equilibrium constants and are inhibited by the generated pyrophosphate.

EC 3.6.1.1

Properties

Molecular weight: 63 kD

Cofactor: Mg²⁺ is required for enzymatic activity.

Isoelectric point: pH 4.8

Specification

Appearance: White suspension

Storage solution: Ammonium sulfate, 3.2 mol/l; pH approximately 6.0

Specific activity: ≥200 U/mg

Unit definition: One unit Pyrophosphatase, inorganic (PPase) generates 1.6 μmol orthophosphate per minute at +25°C by hydrolysis of inorganic pyrophosphate.

Standardized concentration: 5.0±0.5 mg/ml

Contaminating activities (as percentage of PPase activity):

Phosphatase: ≤0.01%

ATPase: ≤0.01%

Stability: At +2 to +8°C within specification range for 36 months.

Cat. No.	Pack Size
10 150 681 103	custom fill

Will be supplied as "Pyrophosphatase, Inorganic (PPase), Yeast".
Unit of Measure is "g".

For further processing only.

Bovine Serum Albumin, Molecular Biology Grade

2% solution

Bovine Serum Albumin, Molecular Biology Grade, supports enzyme stability.

Application

Use Bovine Serum Albumin, Molecular Biology Grade, for enzyme stabilization and for dilution of nucleic acid modifying enzymes.

Benefits

- **Gain excellent performance.** Take advantage of a contamination-controlled reagent with high lot-to-lot consistency.

CAS: 9048-46-8

Properties

Special quality for molecular biology

Molecular weight: 63 kD

Specification

Appearance: Clear, yellowish solution

Storage buffer: Tris/HCl, 50 mmol/l; NaCl, 100 mmol/l; 2-mercaptoethanol, 1 mmol/l; EDTA, 0.25 mmol/l; glycerol, 50% (v/v); pH approximately 7.5 at +25°C

Protein concentration (A_{280} , 1 mg/ml \triangleq 0.67 OD): ≥ 20 mg/ml

Unspecific endonucleases (λ DNA): Not detectable in up to 0.75 mg/ml after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 0.75 mg/ml after 16 hours incubation at +37°C.

Exonucleases (3 H-DNA): Not detectable in up to 0.75 mg/ml after 4 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 1 mg/ml after 4 hours incubation at +37°C.

Proteinases (colorimetric): Not detectable in up to 0.75 mg/ml after 2 hours incubation at +37°C.

pH ≤ 5.5 treatment (≥ 30 minutes): Corresponds to specification

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Prepared of bovine plasma from USA with official veterinary's certificate of health of the donor animals and of the deactivation of animal material at pH ≤ 5.5 for ≥ 30 minutes.

Cat. No.

Pack Size

10 715 859 103 custom fill

Will be supplied as "Albumin.(BSA) SQ for Molecular Biology".
Unit of Measure is "g".



For further processing only.

Bovine Serum Albumin, Molecular Biology Grade

10% solution

Bovine Serum Albumin, Molecular Biology Grade, supports enzyme stability.

Application

Use Bovine Serum Albumin, Molecular Biology Grade, for enzyme stabilization and for dilution of nucleic acid modifying enzymes.

Benefits

- **Gain excellent performance.** Take advantage of a contamination-controlled reagent with high lot-to-lot consistency.

CAS: 9048-46-8

Properties

Special quality for molecular biology

Molecular weight: 63 kD

Specification

Appearance: Clear, yellowish solution

Storage buffer: Tris/HCl, 50 mmol/l; NaCl, 100 mmol/l; DTT, 1 mmol/l; EDTA, 0.25 mmol/l; glycerol, 50% (v/v); pH approximately 7.5 at +25°C

Protein concentration (A_{280} , 1 mg/ml \pm 0.67 OD): 100 \pm 5 mg/ml

Unspecific endonucleases (λ DNA): Not detectable in up to 0.75 mg/ml after 16 hours incubation at +37°C.

Nicking activity (pBR322 DNA): Not detectable in up to 0.75 mg/ml after 16 hours incubation at +37°C.

Exonucleases (3 H-DNA): Not detectable in up to 0.75 mg/ml after 4 hours incubation at +37°C.

Ribonucleases (MS2 RNA): Not detectable in up to 1 mg/ml after 4 hours incubation at +37°C.

Proteinases (colorimetric): Not detectable in up to 0.75 mg/ml

pH \leq 5.5 treatment (\geq 30 minutes): Corresponds to specification

Country of origin: USA

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Prepared of bovine plasma from USA with official veterinary's certificate of health of the donor animals and of the deactivation of animal material at pH \leq 5.5 for \geq 30 minutes.

Cat. No.	Pack Size
05 931 665 103	1000 ml

Will be supplied as "Albumin Bovine Plasma, 10%, Mol Biol Grd".
Unit of Measure is "ml".

For further processing only.

Blocking Reagent powder

Blocking Reagent for nonradioactive filter hybridization and the detection of nucleic acid hybrids, especially for the Roche DIG system.

Application

Use Blocking Reagent to decrease the background in nonradioactive filter hybridization and the detection of nucleic acid hybrids.

Benefits

- **Obtain reliable results.** The Blocking Reagent is perfectly suited for the Roche DIG system.

Specification

Appearance: White, fine-grained powder

Solubility:

Clear, opalescent solution in buffer with 50% formamide (5% [w/v])

Clear, opalescent solution in maleic acid, 0.1 mol/l; NaCl, 0.15 mol/l; pH 7.5 at +20°C (10% [w/v])

Clear, opalescent solution in Tris/HCl, 0.1 mol/l; NaCl, 0.15 mol/l; pH 7.5 at +20°C (5% [w/v])

DNases/RNases: Negative

Nicking activity: Negative

Function test: Hybridization and detection of homologous RNA using DIG

Detection System: Corresponds to specification

Stability: At +15 to +25°C within specification range for 24 months.

Cat. No. **Pack Size**

10 057 177 103 custom fill

Will be supplied as "Blocking Reagent". Unit of Measure is "kg".
For further processing only.

Histone H3 from calf thymus, lyophilizate

Histone H3 forms complexes with DNA and can be of interest in studies of DNA-protein interactions.

Application

Use Histone H3 in studies of DNA-protein interactions, in investigations on structure and function of chromatin, and as substrate for protein kinases.

Properties

Histone H3 is electrophoretically homogeneous.

Molecular weight: 15.3 kD

Specification

Appearance: White lyophilizate

Solubility: Clear, colorless solution in water (c=1 mg/ml)

Protein content: ≥1 mg protein/mg lyophilizate

Purity (SDS gel electrophoresis): ≥95%

pH <5.5 treatment (≥30 minutes): Corresponds to reference

Stability: At +2 to +8°C within specification range for 24 months.

Quality

Prepared of calf thymus from *Belorussia* with official veterinary's certificate of health of the donor animals and of the deactivation of animal material at pH ≤5.5 for ≥30 minutes.

Cat. No. **Pack Size**

11 039 202 103 custom fill

Will be supplied as "Histone H3 from Calf Thymus". Unit of Measure is "mg".
For further processing only.





4 Pharma Biotech

Enzymes340
Enzymes for cell isolation340
Glycohydrolase344
Proteases348
Industrial Process Control359
Mycoplasma Testing359
Biochemicals361
Activated Sugars361
Cofactors365
Additional Reagents367

Liberase MNP-S**lyophilizate, sterile acc. to Ph. Eur.**

Liberase MNP-S is a highly purified enzyme blend for tissue dissociation where high purity, consistent quality, low levels of bacterial endotoxins and a sterile product (STERILE A) are required for reproducible high cell yield and viability.

Application

Use Liberase MNP-S for *in vitro* enzymatic dissociation of tissue, to isolate single cells from a broad range of tissue types. In particular intended for use in isolation of tissue based stem cells and chondrocytes from cartilage.

Benefits

- **Increase safety.** Liberase MNP-S is manufactured free of mammalian and avian tissue-derived raw materials.
- **Maximize viability and yield of isolated cells** due to reduced clostridial and trypsin activity, as well as reduced endotoxin content.
- **Save time and resources** with high experimental reproducibility due to higher lot-to-lot consistency.
- **Count on higher specific activity of these enzyme blends** due to higher Collagenase I + II purity as determined by HPLC analysis.

Product Description

Liberase MNP-S contains highly purified collagenase class I and class II from *Clostridium histolyticum*. The two collagenase isoforms are blended in a precise ratio with each other and with a medium concentration of highly purified thermolysin, a neutral protease isolated from *Bacillus thermoproteolyticus*. The blend is classified STERILE A according to European Pharmacopoeia. A product manufactured in accordance with current Good Manufacturing Practice (cGMP) is available upon request.

EC 3.4.24.-

Specification (for 5 mg pack size)**Appearance:** White lyophilized cake**Total protein content:** 41.5-62.3 mg/bottle (target: 51.9 mg protein)**Collagenase I b, c part of collagenase I:** ≤40 area%**Collagenase I content:** 16.8-25.2 mg/bottle (target: 21.0 mg protein)**Collagenase II content:** 11.2-16.8 mg/bottle (target: 14.0 mg protein)**Thermolysin/neutral protease content:** 13.5-20.3 mg/bottle (target: 16.9 mg protein)**Activity** (Wünsch, calculated): 142-237 U/bottle**Activity** (thermolysin/neutral protease, calculated): 147,030-275,000 U/bottle**Ratio collagenase II/total collagenase:** 0.3-0.5**Endotoxin:** ≤50 EU/mg protein**Sterility** (according to European Pharmacopoeia, current version): Corresponds to reference**Stability:** At -15 to -25°C within specification range for 12 months.**Cat. No.** **Pack Size****05 578 582 001** 5 mg**05 578 566 001** 35 mg

Will be supplied as "Liberase MNP-S". Unit of Measure is "piece".



Blended Proteolytic Enzyme for Tissue Dissociation. For evaluation purposes only.

Liberase MTF C/T, GMP Grade

0.2µm filtered, lyophilizate

Liberase MTF C/T GMP Grade kit is a highly purified enzyme blend for tissue dissociation where high purity, consistent quality and low levels of bacterial endotoxins are required in achieving reproducible high cell yield and viability.

Application

Use the Liberase MTF C/T, GMP Grade kit for *in vitro* enzymatic dissociation of tissue, to isolate single cells. In particular intended for use in isolation procedures of pancreatic islets from human and pig.

Benefits

- **Increase safety.** Liberase MTF C/T, GMP Grade, is manufactured free of mammalian and avian tissue-derived raw materials.
- **Maximize viability and yield of isolated cells** due to reduced clostridial and trypsin activity, as well as reduced endotoxin content.
- **Save time and resources** with high experimental reproducibility due to higher lot-to-lot consistency.
- **Count on higher specific activity of these enzyme blends** due to higher Collagenase I + II purity as determined by HPLC analysis.

Product Description

Liberase MTF C/T, GMP Grade kit contains highly purified collagenase class I and class II from *Clostridium histolyticum*. The two collagenase isoforms are blended in a precise ratio with each other. Highly purified thermolysin from *Bacillus thermoproteolyticus* is provided in a separate vial to allow individual adjustments of enzyme ratios for optimization of your tissue dissociation protocol.

EC 3.4.24.-

Specification

I. Collagenase I/II MTF Blend (≥2,000 U)

Appearance: White lyophilized cake

Protein content: 430-644 mg (target: 537 mg protein)

Collagenase I content: 258-386 mg (target: 322 mg protein)

Collagenase II content: 172-258 mg (target: 215 mg protein)

Collagenase I b, c part of Collagenase I: ≤10 area%

Purity (HPLC, total peak areas of Col I + Col II): ≥85%

Activity (Wünsch, calculated): 2,172-3,617 U/vial

Ratio collagenase II/total collagenase (HPLC): 0.3-0.5

Bioburden: ≤10 CFU/bottle

Endotoxins: ≤10 EU/mg protein

II. Thermolysin MTF (15 mg)

Appearance: White lyophilized cake

Protein content: 12.0-18.0 mg (target: 15.0 mg protein)

Activity neutral protease (calculated): 130,500-234,000 U/vial

Purity (HPLC, total peak area thermolysin): ≥85%

Bioburden: ≤10 CFU/bottle

Endotoxins: ≤50 EU/mg protein

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No. **Pack Size**

05 339 880 001 2x500 mg


Will be supplied as "Liberase MTF C/T (2:3) GMP Grade". Unit of Measure is "piece".



Contents

01. Collagenase I/II MTF Blend (≥2,000 U) GMP Grade, 2x

02. Thermolysin MTF (15mg) GMP Grade, 3x

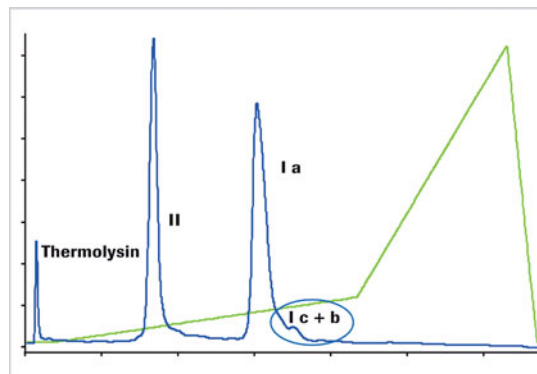

Unmatched purity: HPLC Analysis of Liberase MTF C/T GMP Grade shows only trace levels of Collagenase I b + I c, common impurities in other products.



Human Islets isolated by using Liberase MTF C/T GMP Grade.

Blended Proteolytic Enzyme for Tissue Dissociation.
For further processing only.




CUSTOMBIOTECH

Liberase Enzyme Blends

Leading Performance in Tissue Dissociation

- Mammalian Tissue Free (MTF)
- GMP Grade and/or sterile according to the European Pharmacopoeia
- Custom sizes and blends available, according to your specifications

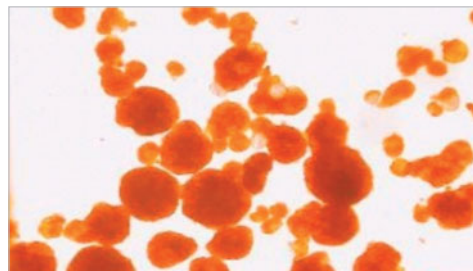

Human Islets isolated by using Liberase MTF

Blended proteolytic enzyme for tissue dissociation.

custombiotech.roche.com/liberasegmp

Your Roche Custom Biotech Customer Service

Europe, Middle East, Africa, and Latin America +49 621 759 8580
mannheim.custombiotech@roche.com

Japan +81 3 5443 5285
japan.custombiotech@roche.com

Asia Pacific +65 6371 6638
apac.custombiotech@roche.com

Canada +1 450 686 7050
custombiotech.can@roche.com

United States
+1 800 428 5433, ext. 14649 (toll-free)
custombiotech.us-sales@roche.com

Roche Diagnostics GmbH
Sandhofer Straße 116
68305 Mannheim, Germany

LIBERASE is a trademark of Roche.

© 2011 Roche Diagnostics.
All rights reserved.



Liberase T-Flex, Research Grade

0.2µm filtered, lyophilizate

The Liberase T-Flex, Research Grade kit enables individual blends of collagenase and thermolysin in low level ratios. The enzymes are manufactured for applications in research where high purity, consistent quality and low levels of bacterial endotoxins are required for reproducible high cell yield and viability.

Application

Use Liberase T-Flex, Research Grade, for *in vitro* enzymatic dissociation of a broad range of tissues, enabling isolation of single cells. In particular intended for use in isolation procedures of pancreatic islets from several animal types.

Benefits

- **Increase safety.** Liberase T-Flex, Research Grade, is manufactured free of mammalian and avian tissue-derived raw materials.
- **Maximize viability and yield of isolated cells** due to reduced clostridial and trypsin activity, as well as reduced endotoxin content.
- **Save time and resources** with high experimental reproducibility due to higher lot-to-lot consistency.
- **Count on higher specific activity of these enzyme blends** as a result of higher Collagenase I + II purity as determined by HPLC analysis.

Product Description

The Liberase T-Flex Research Grade kit contains highly purified collagenase class I and class II from *Clostridium histolyticum*. The two collagenase isoforms are blended in a precise ratio with each other. Highly purified thermolysin from *Bacillus thermoproteolyticus* is provided in a separate vial to allow individual adjustments of enzyme ratios for optimization of your tissue dissociation protocol.

EC 3.4.24.-

Specification

I. Collagenase I/II (500 mg)

Appearance: White lyophilized cake

Protein content: 430-644 mg protein (target: 537 mg protein)

Activity target (Wünsch, calculated): 2,895 U/bottle

Ratio collagenase II/total collagenase: 0.3-0.5

II. Thermolysin (15 mg)

Appearance: White lyophilized cake

Protein content: 12.0-18.0 mg protein (target 15.0 mg protein)

Activity target neutral protease (calculated): 180,000 U/bottle

Stability: At -15 to -25°C within specification range for 12 months.

Cat. No.	Pack Size
05 989 132 001	500 mg

Will be supplied as "Liberase T-Flex Research Grade". Unit of Measure is "piece".

Contents

01. Collagenase I/II Blend (500mg), 1x

02. Thermolysin (15mg), 2x

For life science research only. Not for use in diagnostic procedures.

Neuraminidase (exo- α -Sialidase) from *Vibrio cholerae*, solution

Neuraminidase is a glycohydrolase.

Application

Neuraminidase hydrolyzes terminal N- or O-acetylneuraminic acids which are α 2,3-, α 2,6-, α 2,8- linked to oligosaccharides, polysaccharides, mucopolysaccharides, glycoproteins, and glycolipids. For the hydrolysis of glycolipids, the presence of a detergent is necessary. Roche's Neuraminidase is very well suited for structural research studies on glycoconjugates and for hydrolytic cleavage of sialic acid from biological material.

Benefits

- Profit from broad substrate specificity and absence of proteases.

Product Description

Neuraminidase hydrolyzes terminal N- or O-acyl-neuraminic acids that are α 2,3-, α 2,6-, or α 2,8-linked to galactose, Hex, NAc, or N- or O-acylated neuraminyl residues in oligosaccharides/glycoconjugates or sialic acid. Relative rate of cleavage is α 2 \rightarrow 3 > α 2 \rightarrow 6 > α 2 \rightarrow 8, determined on bonds in tri- and tetrasaccharides.

EC 3.2.1.18

Properties

Molecular weight: Approximately 95 kD

pH optimum: 5.5-6.2

Specific activity: Approximately 20 U/mg total protein (approximately 40 U/mg enzyme protein at +37°C) and pH 5.5 with N-acetyl-neuraminosyl-D-lactose as the substrate.

Specification

Appearance: Clear, colorless solution

Contents (after blending with Micr-O-protect® and EDTA): Sodium acetate, 50 mmol/l; NaCl, 154 mmol/l; CaCl₂, 9 mmol/l; EDTA, 10 mmol/l; Micr-O-protect® 0.1% (w/v)

pH value: 5.0-6.0

Activity (+37°C, with N-acetylneuraminyl-lactose): \geq 1.0 U/ml

Stability: At +2 to +8°C within specification range for 18 months.

Cat. No. _____ **Pack Size** _____

11 087 096 103 custom fill

Will be supplied as "Neuraminidase (Sialidase)". Unit of Measure is "U".

For further processing only.

N-Glycosidase A (PNGase A) solution

Enzyme that cleaves N-glycans.

Application

Use N-Glycosidase A to cleave all types of asparagine bound N-glycans including high mannose-, hybrid-, biantennary-, triantennary- and tetraantennary complex types, provided that the amino group as well as the carboxyl group are present in peptide linkage. Use N-Glycosidase A also to cleave a single N-acetylglucosamine residue from the peptide.

Benefits

- Obtain a clean reaction by using an enzyme without contaminating glycohydrolase activities.

EC 3.5.1.52

Specification

Appearance: Clear, colorless solution

Activity: ≥ 50 mU/ml

Specific activity (Bradford): ≥ 0.5 U/mg protein

Contaminants (expressed as percentage of N-Glycosidase-A activity):

α -Galactosidase: ≤ 0.1

β -Galactosidase: ≤ 0.1

β -Glucosidase: ≤ 0.1

α -Mannosidase: ≤ 0.1

β -Mannosidase: ≤ 0.1

β -N-Acetyl-glucosaminidase: ≤ 0.1

β -Xylosidase: ≤ 0.1

α -Fucosidase: ≤ 0.1

Sialidase (after 17 hours at +37°C): ≤ 0.1

Proteases (casein-resorufin, for 17 hours at +37°C): $\Delta A/17h \leq 0.025$

Stability: At -15 to -20°C within specification range for 24 months.

Cat. No. **Pack Size**

11 646 583 103 custom fill

Will be supplied as "N-Glycosidases A, Solution". Unit of Measure is "milliU".



For further processing only.

Catalase

from *Corynebacterium glutamicum*, lyophilizate

Catalase can be used for removal of H_2O_2 .

Application

Catalase is used in eye care industry for neutralization of H_2O_2 for cleaning contact lenses.

Benefits

Increase cost efficiency by reducing enzyme input. In comparison to catalase from bovine liver, Roche Catalase from *Corynebacterium glutamicum* has:

- Higher thermal stability
- Higher pH stability
- More stability in the presence of H_2O_2

Product Description

Biological activity: 150 IU of Catalase are sufficient to completely remove 3% hydrogen peroxide in 10 minutes.

Source: Microbial from *Corynebacterium glutamicum*

EC 1.11.1.6

Specification

Appearance: Brown lyophilizate

Solubility: Clear solution in phosphate buffer, 0.05 mol/l, pH 7.5 (c=40 mg/ml)

pH value (c=10 mg/ml in water): 7.0-8.0

Activity (+25°C, H_2O_2): 11,000-17,000 U/mg lyophilizate

Sucrose (TC food, Ident No. 10 139 041): ≥60%

Purity (HPLC, A_{280}): ≥90 area%

Water (K. Fischer): ≤5%

Bioburden:

Total amount: ≤50 CFU/g lyophilizate

Germ differentiation according to DAB 10 (*E. coli*, *Salmonellae*, *Staphylococcus aureus*, *Pseudomonas aeruginosa*): Below the limit of detection


Stability: At +2 to +8°C within specification range for 18 months.

Cat. No.


Pack Size

11 650 645 103 custom fill


Will be supplied as "Catalase, Microbial Lyophilizate". Unit of Measure is "GU".

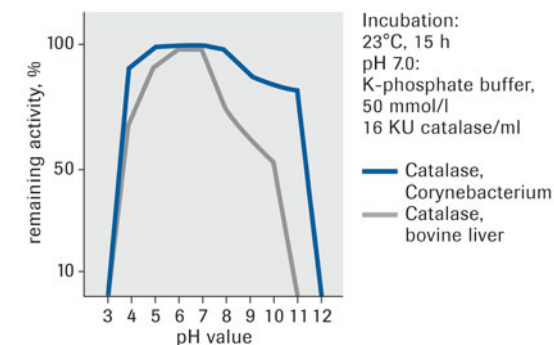


pH optimum

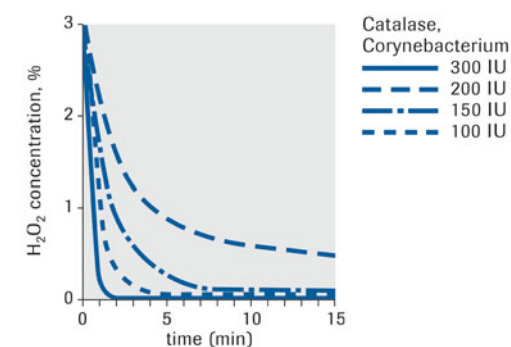

Reduction potential

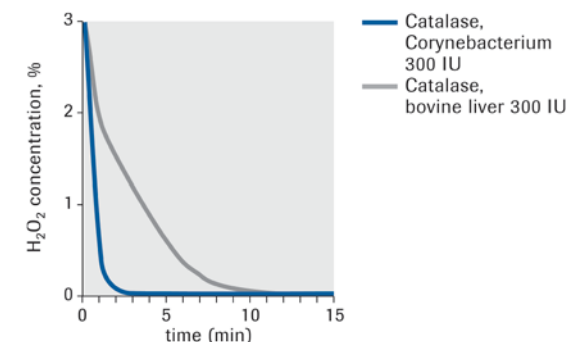

Removal efficiency of H_2O_2

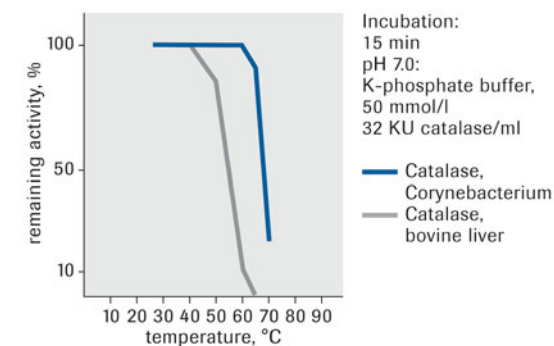

Thermal stability

For further processing only.

Catalase

from *Corynebacterium glutamicum*, solution

Catalase can be used for removal of H₂O₂.

Application

Catalase is used in eye care industry for neutralization of H₂O₂ for cleaning contact lenses.

Benefits

Increase cost efficiency by reducing enzyme input. In comparison to catalase from bovine liver, Roche Catalase from *Corynebacterium glutamicum* has:

- Higher thermal stability
- Higher pH stability
- More stability in the presence of H₂O₂

Product Description

Biological activity: 150 IU of Catalase are sufficient to completely remove 3% hydrogen peroxide in 10 minutes.

Source: Microbial from *Corynebacterium glutamicum*

EC 1.11.1.6

Specification

Appearance: Brownish-green, clear solution

pH value: 7.0-8.0

Density (according to DAB 10): No limit

Activity (+25°C, H₂O₂): ≥500,000 U/ml

Glycerol (enzymatically): 310-430 mg/ml (± 25-35%, v/v)

Ethanol (enzymatically): 71-87 mg/ml (± 9-11%, v/v)

Water: No limit

Purity (HPLC, A₂₈₀): ≥90 area%

Bioburden:

Total amount: ≤15 CFU/ml

Germ differentiation according to DAB 10 (*E. coli*, *Salmonellae*, *Staphylococcus aureus*, *Pseudomonas aeruginosa*): Below the limit of detection

Stability: At +2 to +8°C within specification range for 18 months.

Cat. No. **Pack Size**

11 668 153 103 custom fill

Will be supplied as "Catalase Microbial in Glycerol". Unit of Measure is "GU".

For further processing only.

Carboxypeptidase B, recombinant from rat pancreas, expressed in *Pichia pastoris*, solution

Carboxypeptidase B, recombinant, is intended to use in highly regulated production processes at pharmaceutical companies.

Application

Use the animal component-free and GMP-manufactured Carboxypeptidase B, recombinant, as critical raw material for the production of active pharmaceutical ingredients (API), *i.e.*, insulin.

Benefits

- **Obtain a cGMP-manufactured product** for use in highly regulated processes.
- **Eliminate the risk of virus contamination** and the risk of animal-related cross-contamination.
- **Rely on high purity.** Minimize the risk of host cell protein contamination in your final product.
- **Increase the safety of your production processes** with robust and reproducible performance and high lot-to-lot consistency.

Product Description

Carboxypeptidase B is a widely used metalloprotease, typically isolated from pancreas of different animals, that specifically releases arginine and lysine from the C-terminus of peptides and proteins. Roche has chemically synthesized a gene encoding for the amino acid sequence of the rat Carboxypeptidase B and has transformed the gene into the expression host *Pichia pastoris*, which expresses the recombinant Carboxypeptidase B as active protease with identical properties compared to the native Carboxypeptidase B. The product is manufactured under current good manufacturing practice (cGMP). No animal-derived products are used in the fermentation, purification and final formulation. The production process is validated resulting in a very high lot-to-lot consistency.

EC 3.4.17.2

Properties

Molecular weight: 34.6 kD

Specification

Appearance: Clear, colorless to slightly yellowish solution

Storage buffer: Tris/HCl, 33 mmol/l; ZnCl₂, 0.1 mmol/l; pH 7.5-8.5 at +25°C

Activity (hippurylarginine): ≥400 U/ml

Total activity: 30 kU ±15%

Specific activity: ≥210 U/mg

Protein: 2.8±0.5 mg/ml

Purity (RP-HPLC, according to master lot): ≥85%

Bioburden: ≤50 CFU/ml

Trypsin (Chromozym TRY): ≤0.005%

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Carboxypeptidase B, recombinant, is produced completely animal component-free and according to cGMP guidelines.

Background Information

For several years Roche has successfully pursued the strategy of replacing animal-derived enzymes, frequently used in pharmaceutical production processes with recombinant, animal component-free enzymes.

Related products are recombinant Trypsin, recombinant Proteinase K, recombinant DNase I and others.

Cat. No. _____ **Pack Size** _____

03 358 682 103 30 kU

Will be supplied as "CpB rec.". Unit of Measure is "MU".



For further processing only.

Endoproteinase Asp-N, Sequencing Grade from a mutant of *Pseudomonas fragi*, lyophilizate

Cat. No. **Pack Size**

11 058 541 103 2 µg

Will be supplied as "Endoproteinase Asp-N Sequencing Grade".
Unit of Measure is "piece".
For further processing only.

Endoproteinase Asp-N can be used for specific cleavage of peptides.

Application

Use Endoproteinase Asp-N, Sequencing Grade, for protein structure analysis and sequence analysis.

Benefits

- **Obtain consistent and clear peptide sequencing results.** Minimize the risk of unknown peptide impurities by using this highly purified quality.

Product Description

Endoproteinase Asp-N is a widely used metalloprotease, that specifically hydrolyzes peptide bonds at the amino side of asparagine and cysteine. If cysteine is reduced or alkylated only -↓-Asp-X is cleaved. Endoproteinase Asp-N is isolated from *Pseudomonas fragi* mutant. The protease is supplied as lyophilizate.

EC 3.4.24.33

Properties

Molecular weight: 27 kD

pH optimum: 7.0-8.0

Inhibitors: EDTA and α-phenanthroline

Specification

Appearance: White lyophilizate

Activity (+37°C, azocoll): ≥40 U/bottle

Specific activity: ≥20,000 U/mg protein

Purity (SDS PAGE, homogeneity; PhastSystem®, 8-25% reducing terms): ≥90%

Specificity (HPLC, glucagon, after 1 hour incubation): ≥90%

Unspecific cleavage peptides (HPLC, melittin, after 4 hours incubation): ≤10%

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Quality

Endoproteinase Asp-N, Sequencing Grade, is free of impurities, according to the current Quality Control procedures, that may interfere with separation of peptides in RP HPLC.

Background Information

Roche offers a broad portfolio of highly purified endoproteases suitable for peptide sequencing or other specific application. Related products are Endoproteinase Lys-C, Sequencing Grade, Endoproteinase Glu-C, Sequencing Grade, Endoproteinase Arg-C, Sequencing Grade, and others.

Endoproteinase Glu-C

from *Staphylococcus aureus* V8, lyophilizate, salt-free

Endoproteinase Glu-C can be used for specific cleavage of peptides.

Application

Use Endoproteinase Glu-C for protein structure analysis and sequence analysis.

Product Description

Endoproteinase Glu-C is a widely used serine protease, that specifically hydrolyzes peptide bonds at the carboxylic side of glutamic acid (in ammonium bicarbonate, pH 7.8, or ammonium acetate buffer, pH 4.0), or glutamic acid and aspartic acid (in phosphate buffer, pH 7.8). Endoproteinase Glu-C is isolated from *Staphylococcus aureus* and is supplied as a salt-free lyophilizate.

EC 3.4.21.19

Properties

Molecular weight: 30 kD

pH optimum: 4.0-7.8

Inhibitors: DFP, α_2 -macroglobulin, and TLCK

Specification

Appearance: White lyophilizate

Activity: ≥ 20 U/mg lyophilizate

Stability: At +2 to +8°C within specification range for 24 months.

Background Information

Roche offers a broad portfolio of highly purified endoproteinases suitable for hydrolyzing of peptides or proteins. Related products are Endoproteinase Lys-C, Endoproteinase Asp-N, Endoproteinase Arg-C, and others.

Cat. No.	Pack Size
10 787 906 103	50 μ g

Will be supplied as "Endoproteinase Glu-C (Staph.aureus V 8)".
Unit of Measure is "mg".

For further processing only.

Endoproteinase Glu-C, Sequencing Grade from *Staphylococcus aureus* V8, lyophilizate

Cat. No.	Pack Size
11 058 525 001	50 µg

Endoproteinase Glu-C, Sequencing Grade, can be used for specific cleavage of peptides.

Application

Use Endoproteinase Glu-C, Sequencing Grade, for protein structure analysis and sequence analysis.

Benefits

- **Obtain consistent and clear peptide sequencing results.** Minimize the risk of unknown peptide impurities by using this highly purified quality.

Product Description

Endoproteinase Glu-C is a widely used serine protease, that specifically hydrolyzes peptide bonds at the carboxylic side of glutamic acid (in ammonium bicarbonate, pH 7.8, or ammonium acetate buffer, pH 4.0), or glutamic acid and aspartic acid (in phosphate buffer, pH 7.8). Endoproteinase Glu-C is isolated from *Staphylococcus aureus* and is supplied as a salt-free lyophilizate.

EC 3.4.21.19

Properties

Molecular weight: 30 kD

pH optimum: 4.0-7.8

Inhibitors: DFP, α_2 -macroglobulin, and TLCK

Specification

Appearance: White lyophilizate

Specific activity: ≥ 15 U/mg protein

Protein (BCA): ≥ 50 µg/bottle

Contaminants (HPLC): $\leq 10\%$

Specificity (HPLC, insulin B_{0x}^{*} after 1 hour incubation): $\geq 90\%$

Unspecific cleavage peptides (after 18 hours incubation): $\leq 5\%$

Stability: At +2 to +8°C within specification range for 24 months. Store dry.

Quality

Endoproteinase Glu-C is free of impurities, according to the current Quality Control procedures, that may interfere with separation of peptides in RP HPLC.

Background Information

Roche offers a broad portfolio of highly purified endoproteases suitable for peptide sequencing or other specific application. Related products are Endoproteinase Lys-C, Sequencing Grade, Endoproteinase Asp-N, Sequencing Grade, Endoproteinase Arg-C, Sequencing Grade, and others.

Will be supplied as "Endoproteinase Glu-C Sequencing Grade".

Unit of Measure is "piece".

For further processing only.

Endoproteinase Lys-C

from *Lysobacter enzymogenes*, lyophilizate

Endoproteinase Lys-C can be used for specific cleavage of peptides.

Application

Use Endoproteinase Lys-C for protein structure analysis and for sequence analysis. It is suitable to digest proteins in solution, in polyacrylamide gels or on blotting membranes.

Product Description

Endoproteinase Lys-C is a widely used serine protease, that specifically hydrolyzes amide, ester, and peptide bonds at the carboxylic side of lysin. Endoproteinase Lys-C is isolated from *Lysobacter enzymogenes* and is supplied as a lyophilizate.

EC 3.4.21.50

Properties

Molecular weight: 33 kD (reduced), 30 kD (nonreduced)

pH optimum: 8.5-8.8

Inhibitors: DFP, TLCK, aprotinin, and leupeptin

Specification

Appearance: White lyophilizate

Activity: ≥3 U/bottle

Specific activity: ≥150 U/mg protein

Unspecific cleavage peptides (HPLC, after 6 hours incubation): ≤5 area%

Stability: At +2 to +8°C within specification range for 12 months.

Background Information

Roche offers a broad portfolio of endoproteases suitable for hydrolyzing of peptides or proteins. Related products are Endoproteinase Glu-C, Endoproteinase Asp-N, Endoproteinase Arg-C and others.

Cat. No.	Pack Size
10 476 978 103	3 U

Will be supplied as "Endoproteinase Lys-C, Lysobacter enzymog.". Unit of Measure is "U".

For further processing only.

Endoproteinase Lys-C, Sequencing Grade from *Lysobacter enzymogenes*, lyophilizate

Cat. No. **Pack Size**

11 058 533 103 5 µg

Endoproteinase Lys-C, Sequencing Grade, can be used for specific cleavage of peptides.

Application

Use Endoproteinase Lys-C, Sequencing Grade, for protein structure analysis and for sequence analysis. It is suitable to digest proteins in solution, in polyacrylamide gels or on blotting membranes.

Benefits

- **Obtain consistent and clear peptide sequencing results.** Minimize the risk of unknown peptide impurities by using this highly purified quality.

Product Description

Endoproteinase Lys-C, Sequencing Grade, is a widely used serine protease, that specifically hydrolyzes amide, ester, and peptide bonds at the carboxylic side of lysin. Endoproteinase Lys-C is isolated from *Lysobacter enzymogenes* and is supplied as a lyophilizate.

EC 3.4.21.50

Properties

Molecular weight: 33 kD (reduced), 30 kD (nonreduced)**pH optimum:** 8.5-8.8**Inhibitors:** DFP, TLCK, aprotinin, and leupeptin

Specification

Appearance: White lyophilizate**Activity:** ≥3 U/bottle**Specific activity:** ≥150 U/mg protein**Unspecific cleavage peptides** (HPLC, after 6 hours incubation; based on area sum): ≤5 area%**Stability:** At +2 to +8°C within specification range for 12 months.

Quality

Endoproteinase Lys-C, Sequencing Grade, is free of impurities, according to the current Quality Control procedures, that may interfere with separation of peptides in RP HPLC. Purity is checked by HPLC and SDS PAGE using silver staining.

Background Information

Roche offers a broad portfolio of highly purified endoproteases suitable for peptide sequencing or other specific application. Related products are Endoproteinase Glu-C, Sequencing Grade, Endoproteinase Asp-N, Sequencing Grade, Endoproteinase Arg-C, Sequencing Grade, and others.

Will be supplied as "Endoprot. Lys-C, Sequ., MPB 5 UG". Unit of Measure is "piece".

For further processing only.

Endoproteinase Lys-C, Sequencing Grade from *Lysobacter enzymogenes*, solution

Endoproteinase Lys-C, Sequencing Grade, can be used for specific cleavage of peptides.

Application

Use Endoproteinase Lys-C, Sequencing Grade, for protein structure analysis and for sequence analysis. It is suitable to digest proteins in solution, in polyacrylamide gels or on blotting membranes.

Benefits

- **Obtain consistent and clear peptide sequencing results.** Minimize the risk of unknown peptide impurities by using this highly purified quality.

Product Description

Endoproteinase Lys-C is a widely used serine protease, that specifically hydrolyzes amide, ester, and peptide bonds at the carboxylic side of lysin. Endoproteinase Lys-C is isolated from *Lysobacter enzymogenes* and is supplied as a solution.

EC 3.4.21.50

Properties

Molecular weight: 33 kD (reduced), 30 kD (nonreduced)

pH optimum: 8.5-8.8

Inhibitors: DFP, TLCK, aprotinin, and leupeptin

Specification

Appearance: Clear, colorless solution

Activity (+25°C, Chromozym PL): No limit (U/ml)

Specific activity: ≥200 U/mg protein

Protein (A_{280}): ≥0.1 mg/ml

Purity (SDS PAGE): ≥90%

Specificity (HPLC, melittin, after 1 hour incubation): ≥90%

Unspecific cleavage peptides (HPLC, Insulin Box, after 18 hours incubation): ≤5 area%

Stability: At -15 to -25°C within specification range for 12 months.

Quality

Endoproteinase Lys-C, Sequencing Grade, is free of impurities, according to the current Quality Control procedures, that may interfere with separation of peptides in RP HPLC. Purity is checked by HPLC and SDS PAGE using silver staining.

Background Information

Roche offers a broad portfolio of highly purified endoproteases suitable for peptide sequencing or other specific application. Related products are Endoproteinase Glu-C, Sequencing Grade, Endoproteinase Asp-N, Sequencing Grade, Endoproteinase Arg-C, Sequencing Grade, and others.

Cat. No.

Pack Size

11 051 199 103 custom fill

Will be supplied as "Endoproteinase Lys-C Sequencing Grade".
Unit of Measure is "mg".



For further processing only.

Papain

from *Carica papaya*, suspension

Cat. No.	Pack Size
10 154 644 103	custom fill

Papain can be used for complete proteolytic cleavage of peptides and proteins.

Will be supplied as "Papain from *Carica papaya*". Unit of Measure is "g".

For further processing only.

Application

Use Papain for:

- Complete proteolytic cleavage of peptides and proteins
- Limited hydrolysis of native immunoglobulins
- Tissue dissociation (together with collagenase, esterase, trypsin)
- Solubilization of integral membrane proteins
- Production of glycopeptides from purified proteoglycans

Benefits

- **Increase the safety of your processes** with reproducible performance and high lot-to-lot consistency.

Product Description

Papain is a cysteine endopeptidase, hydrolyzing proteins, peptides, amides, and esters of amino acids and peptides, especially at bonds involving Arg, Lys, Glu, His, Gly and Tyr. Upon prolonged incubation, further bonds are cleaved; shows also esterase and transamidase activity, and is used for peptide synthesis.

Papain is isolated from *Carica papaya* and is supplied as a lyophilizate.

EC 3.4.22.2

Properties

Molecular weight: 23 kD ±2

Specification

Appearance: White suspension

pH value (c=10 mg/ml): 4.0-5.0

Specific activity (+25°C, BAEE): ≥30 U/mg protein

Protein (ΔA_{280} , 1%₂₄): ≥10 mg/ml

Stability: At + 2 to +8°C within specification range for 12 months.

Papain**from *Carica papaya*, 0.2µm filtered, solution**

Papain can be used for complete proteolytic cleavage of peptides and proteins.

Application

Use Papain for:

- Complete proteolytic cleavage of peptides and proteins
- Limited hydrolysis of native immunoglobulins
- Tissue dissociation (together with collagenase, esterase, trypsin)
- Solubilization of integral membrane proteins
- Production of glycopeptides from purified proteoglycans

Benefits

- **Increase the safety of your processes** with reproducible performance and high lot-to-lot consistency.
- **Minimize the risk of microbial cross-contamination** in your product by very low bioburden.

Product Description

Papain is a cysteine endopeptidase, hydrolyzing proteins, peptides, amides, and esters of amino acids and peptides, especially at bonds involving Arg, Lys, Glu, His, Gly and Tyr. Upon prolonged incubation, further bonds are cleaved. Papain shows also esterase and transamidase activity, and is used for peptide synthesis.

Papain is isolated from *Carica papaya* and is supplied as a lyophilizate.

EC 3.4.22.2

Properties

Molecular weight: 23 kD ±2

Specification

Appearance: White suspension

pH value: 4.3-4.8

Solubility: Clear, colorless solution in water (c=12 mg/ml)

Specific activity (+25°C, BAEE): ≥16.0 U/mg protein

Protein (ΔA_{280} , 1% = 24, standardized, measured against sodium acetat, 50 mmol/l): 12 ± 1 mg/ml

Purity (SDS PAGE): ≥80%

Band profile based on reference lot: Corresponds to specification

Bioburden (according to DAB 10): ≤100 CFU/ml

Endotoxins: ≤100 EU/ml

Stability: At +2 to +8°C within specification range for 6 months.

Quality

0.2 µm filtered describes the state-of-the art process to reduce the risk of contamination by passing solutions through a filter with 0.22 µm or less in pore size. This reagent is not specified as being sterile according to pharmacopoeia monographs.

Cat. No.

Pack Size

11 720 422 103 custom fill

Will be supplied as "Papain filtered". Unit of Measure is "g".
For further processing only.

Pronase from *Streptomyces griseus*, lyophilizate

Pronase is used for total degradation of proteins.

Application

Use Pronase for:

- Total degradation of proteins during isolation of DNA and RNA
- Total hydrolysis of proteins for technical purposes
- Tissue dissociation in conjunction with collagenases and trypsin
- Production of glycopeptides from purified glycoproteins

Benefits

- Rely on a high proteolytic activity and a low nuclease contamination activity.
- Increase the safety of your application processes with reproducible performance and high lot-to-lot consistency.

Product Description

Pronase is a mixture of several nonspecific endo- and exoproteases that digest proteins down to single amino acids. Resolves carboxylic acids and alcohols.

EC 3.4.24.4

Properties

pH optimum: 6.0-7.5

Specification

Appearance: Near white lyophilizate

Activity (+40°C, casein): ≥6,000 U/g lyophilizate

Contaminants:

DNases/RNases: ≤0.002 U/mg lyophilizate

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Quality

Contaminants are ≤0.002 U DNase/RNase per mg lyophilizate.

Background Information

The preparation contains approximately 20% calcium acetate.

Cat. No. **Pack Size**

10 165 913 103 custom fill

Will be supplied as "Pronase from *Streptomyces griseus*". Unit of Measure is "g".

For further processing only.

Trypsin, recombinant from porcine pancreas, expressed in *Pichia pastoris*

Trypsin, recombinant, is intended to use in highly regulated production processes at pharmaceutical companies.

Application

Use the animal component-free and GMP-manufactured Trypsin, recombinant, as critical raw material for the production of active pharmaceutical ingredients (API), *i.e.*, insulin, vaccines and for cell dissociation.

Benefits

- **Obtain a cGMP-manufactured product** for use in highly regulated processes.
- **Eliminate the risk of virus contamination** and the risk of animal-related cross-contamination.
- **Rely on high purity.** Minimize the risk of host cell protein contamination in your final product.
- **Increase the safety of your production processes** with reproducible performance and high lot-to-lot consistency.

Product Description

Trypsin is a widely used serine protease, typically isolated from pancreas of different animals, that specifically cleaves at the C-terminus of arginine and lysine within a peptide chain. Roche has chemically synthesized a gene encoding for the amino acid sequence of the porcine Trypsin and has transformed the gene into the expression host *Pichia pastoris*, which expresses the recombinant Trypsin as active protease with identical properties compared to the native Trypsin.

The product is manufactured under current good manufacturing practice (cGMP). No animal-derived products are used in the fermentation, purification and final formulation. The production process is validated resulting in a very high lot-to-lot consistency.

EC 3.4.4.4

Properties

Molecular weight: 23.5 kD

pH optimum: 8.0

Inhibitors: TLCK, DFP, PMSF, leupeptin, soybean trypsin inhibitor, trypsin inhibitor from hen egg, aprotinin, α_2 -macroglobulin, α_1 -antitrypsin, APMSE, and antipain

Specification

Appearance: Clear, colorless to slightly yellowish solution

Storage buffer: HCl, 10 mmol/l; CaCl₂, 20 mmol/l

pH value: 2±0.5

Activity (Chromozym TRY): ≥10,800 U/ml

Total activity (Chromozym TRY): 3.5 MU ±10%

Specific activity: ≥180 U/mg

Protein: 70±10 mg/ml

Purity (RP-HPLC, according to master lot): ≤20% α -trypsin, ≥70% β -trypsin

Bioburden: ≤100 CFU/ml

Stability: At -15 to -25°C within specification range for 24 months.

Quality

Trypsin, recombinant, is produced completely animal component-free and according to cGMP guidelines.

Background Information

For several years Roche has successfully pursued the strategy of replacing animal-derived enzymes, frequently used in pharmaceutical production processes with recombinant, animal component-free enzymes.

Related products are recombinant Carboxypeptidase B, recombinant DNase I, and others.

Cat. No. **Pack Size**

03 358 658 103 3.5 MU

06 369 880 103 1 g

Will be supplied as "Trypsin rec. Bulk". Unit of Measure is "MU".



For further processing only.

MycoTOOL PCR Mycoplasma Detection Kit

Highly sensitive PCR test for the detection of Mycoplasmas, accepted for validation according to the E.P. 2.6.7 directive.

This test is suitable for release testing and in-process control. It can replace culture and indicator cell tests.

Application

Use the MycoTOOL test in the field of biopharmaceutical and vaccine production, cell therapy, transplantation, and veterinary testing.

Description of the workflow: 1 ml of sample is divided and purified in 2 x 200 µl eluate. From each eluate two Mycoplasma PCR reactions and two GAPDH PCR reactions are performed (20 µl reaction volume). A GAPDH PCR is performed using undiluted eluate to control for inhibition effects of the matrix. This control eliminates the risk of undetected intracellular Mycoplasma. A PCR with a 1:100 dilution of the eluate controls test sensitivity. The desired dilution can be adapted to the tested matrix. Analysis of the amplicons is performed using polyacrylamide gel electrophoresis. No poststaining is required since the detection dye is included in the PCR mix.

Benefits

- **Save time.** Replace traditional 28 days culture testing with the rapid MycoTOOL PCR Test.
- **Speed up your validation.** Rely on published Roche validation data (see Literature).
- **Detect traces of Mycoplasmas.** The sensitivity is 1 CFU/ml.
- **Minimize the risk of false negatives and false positives.** The matrix lysis control eliminates risk of undetected intracellular Mycoplasma. The diluted GAPDH and positive control identify potential PCR inhibition. False positives are prevented by nucleic acid free reagents, tested according to Roche's current ultrasensitive Quality Control procedures. Uracil-DNA Glycosylase pretreatment is performed to prevent contamination by PCR carryover.
- **Detect the broad panel of Mollicute species.** A universal primer design covers a whole range of Mycoplasmas.
- **Achieve process security.** Rely on test results using a consistent kit quality and established change control procedures.

Product Description

The MycoTOOL Mycoplasma Detection Prep Kit contains all necessary reagents for purifying the DNA, while the MycoTOOL Mycoplasma Detection Amplification Kit supplies all reagents for the amplifications step.

Specification

MycoTOOL, Mycoplasma Detection Amp Kit, Cat. No. 05 184 240

Bottle 1a: Clear, viscous solution

Bottle 1b, 2, 3, 4, 6, 7, 8, 9, 10: Clear, colorless solution

Bottle 5: Clear, slightly yellow solution

Bottle 8 (negative control): No amplicon is detectable using Roche reference material.

Bottle 9 (positive control of the PCR): 2 of 4 replicates are recognized at 360 bp using Roche reference material.

Stability: At -15 to -25°C within specification range for 12 months. Once the kit is opened, store all kit components at +15 to +25°C. In addition, store vial 5, the detection dye, protected from light.

MycoTOOL, Mycoplasma Detection Prep Kit, Cat. No. 05 184 592

Bottle 1, 2, 3, 4, 5, 6: Colorless solution

Packaging: 45 tubes with colorless cap in printed Zip-Lock bag

Stability: At -15 to -25°C within specification range for 12 months. Once the kit is opened, store all kit components at room temperature. Do not freeze or store the kit in a refrigerator.

Cat. No. Pack Size

05 200 709 001 1 Kit
for testing of 10 samples (1 ml each)

The MycoTOOL PCR Mycoplasma Detection Kit is shipped in 2 separate subkits, according to the components' shipping temperature requirements.

Subkit 1: MycoTOOL Mycoplasma Detection Prep Kit, Cat. No. 05 184 592 001; shipping at room temperature.

Subkit 2: MycoTOOL Mycoplasma Detection Amplification Kit, Cat. No. 05 184 240 001; shipping on dry ice. Unit of Measure is "piece".



Contents

Subkit 1: MycoTOOL Mycoplasma Detection Prep Kit, Cat. No. 05 184 592 001

01. Proteinase K
02. Lysis Buffer
03. Precipitation Reagent
04. Washing Buffer
05. Dissolution Buffer
06. Reaction vials

Subkit 2: MycoTOOL Mycoplasma Detection Amplification Kit, Cat. No. 05 184 240 001

- 01a. RM1a Master Mix
- 01b. RM1b Master Mix
02. MgCl₂-Solution
03. Primer Mix, Mycoplasma
04. Primer Mix, GAPDH
05. Detection Dye
06. Water, PCR Grade
07. Dilution Buffer
08. Negative Control
09. Positive Control
10. DNA Molecular Weight Marker

For use in quality control/manufacturing processes only.

Quality

False positives are prevented by nucleic acid free reagents, tested according to Roche's current ultrasensitive Quality Control procedures. Uracil-DNA Glycosylase pretreatment prevents contamination by PCR carryover. Change control procedures are used for all reagents and kits.

Background Information

Mycoplasmas are frequent causes of contamination in biopharmaceutical production, cell therapy and tissue engineering. Traditional detection methods, required by Pharmacopoeias and drug regulating agencies worldwide, use growth on culture media to detect contaminating organisms. These culture-based methods are time-consuming, requiring as much as 28 days to complete, laborious and difficult to interpret. Rapid methods, like the new MycoTOOL Test described here, for detecting Mycoplasma contamination, can help improve efficiency, quality and safety in the manufacturing of pharmaceutical and biological products.

Literature

Sven M. Deutschmann, Holger Kavermann, Yvonne Knack, Validation of a NAT-based Mycoplasma assay according European Pharmacopoeia, Biologicals, Volume **38**, Issue 2, Special Section: Mycoplasma (pp. 181-248), March 2010, Pages 238-248, ISSN 1045-1056

CMP-N-Acetylneuraminic Acid sodium salt

CMP-N-Acetylneuraminic Acid is an activated sugar.

Application

Use CMP-N-Acetylneuraminic Acid for glycosylation of target substances with the suitable glycotransferase.

CAS: 3063-71-6

Properties

Molecular weight: 659.39 D

Use the activated sugar to transfer sialic acid on an existing sugar chain together with the appropriate sialyltransferase.

Specification

Appearance: White, lyophilizate

Solubility: Clear, colorless solution in water

CMP-NANA (A_{260}): $\geq 87\%$ ($\epsilon=7,4 \text{ l x mmol}^{-1} \times \text{cm}^{-1}$)

CMP-NANA (HPLC): $\geq 90 \text{ area}\%$

CMP-NANA + CMP (HPLC): $\geq 95 \text{ area}\%$

Na (flame photometric): 6-8%

Heavy metals (as Pb): $\leq 25 \text{ ppm}$ (in validation)

CMP (HPLC): $\leq 10 \text{ area}\%$

CDP (HPLC): $\leq 1 \text{ area}\%$

CTP (HPLC): $\leq 0.5 \text{ area}\%$

Sum of unknown impurities (HPLC): $\leq 1\%$

Bioburden: $\leq 100 \text{ CFU/g}$

Stability: At -15 to -25°C within specification range for 12 months.

Quality

CMP-N-Acetylneuraminic Acid is a highly purified product suitable for pharmaceutical applications. It can be produced according to cGMP regulations or as a completely animal component-free product on request. Due to the special designed synthesis pathway, Roche can produce kilogram amounts of CMP-N-Acetylneuraminic Acid in a cost-effective manner.

Background Information

Several years ago, Roche developed a new cost-effective synthesis pathway for the production of activated sugars in large scale. For the production of all activated sugars a valid TSE/BSE certificate is available. Production can be completely animal component-free and according to cGMP regulations on request.

Related products are UDP-Gal (04589173103), UDP-Glc (10154938103), and UDP-GlcNAc (11787900103).

Cat. No. **Pack Size**

05 974 003 103 custom fill

Will be supplied as "CMP-NANA, Sodium salt". Unit of Measure is "kg".



For further processing only.

UDP-N-Acetylglucosamine disodium salt

UDP-N-Acetylglucosamin is an activated sugar.

Application

Use UDP-GlcNAc for glycosylation of target molecules with a suitable enzyme.

CAS: 7277-98-7

Properties

Formula: $C_{17}H_{25}N_3O_{17}P_2Na_2$

Molecular weight:

UDPGlcNAc: 607.4 D

UDPGlcNAc- Na_2 : 651.4

Specification

Appearance: White to slightly yellowish, amorphous powder

Solubility: Clear, colorless to slightly yellowish solution in water (c=100 mg/ml)

pH value (c=100 mg/ml, in water): 5.0-7.0

UDPGlcNAc- Na_2 (A_{260} , $\epsilon=9.9 \text{ l} \times \text{mmol}^{-1} \times \text{cm}^{-1}$): $\geq 80.0\%$

UDPGlcNAc (HPLC): ≥ 96.0 area%

Na (flame photometric): 6.0-8.0%

Water (K. Fischer): $\leq 6.0\%$

P_i: $\leq 0.4\%$

Heavy metals (as Pb): ≤ 10 ppm

Fe (AAS): ≤ 10 ppm

UDP (HPLC): ≤ 1.0 area%

UMP (HPLC): ≤ 1.0 area%

UDP-galactose (HPLC): ≤ 0.5 area%

UDPG (HPLC): ≤ 1.5 area%

Uridine (HPLC): ≤ 0.5 area%

AMP (HPLC): ≤ 1.0 area%

ATP (HPLC): ≤ 0.2 area%

UTP (HPLC): ≤ 0.2 area%

Sum of unknown impurities (HPLC): ≤ 1.0 area%

Ethanol (GC): $\leq 10.0\%$

Isopropanol (GC): ≤ 20 ppm

Aceton (GC): ≤ 20 ppm

Methanol (GC): ≤ 20 ppm

Bioburden: ≤ 100 CFU/g

Mould: ≤ 50 CFU/g

Yeast: ≤ 50 CFU/g

Stability: At +2 to +8°C within specification range for 30 months.

Cat. No.

Pack Size

11 787 900 103 custom fill

Will be supplied as "UDP-Glc-Nac". Unit of Measure is "kg".
For further processing only.

UDP-Galactose disodium salt, powder

UDP-Galactose is an activated sugar.

Application

Use UDP-Galactose for glycosylation of target substances using suitable transferase enzyme.

CAS: 137868-52-1

Properties

Formula: C₁₅H₂₂N₂O₁₇P₂Na₂

Molecular weight:

UDP-galactose: 566.3 D

UDP-galactose-Na₂: 610.3 D

Specification

Appearance: White powder

Solubility: Clear, colorless solution in water (c=100 mg/ml)

pH value (solution, c=100 mg/ml): 5.0-7.0

UDP-galactose, disodium salt (A₂₆₀, ε=9.9 l x mmol⁻¹ x cm⁻¹): ≥80.0%

UDP-galactose (HPLC): ≥95.0 area%

Na (flame photometric): 6.0-8.0%

Water (K. Fischer): ≤6.0%

UDP (HPLC): ≤1.5 area%

UMP (HPLC): ≤0.5 area%

UDP-GlcNAc (HPLC): ≤0.5 area%

UDPG (HPLC): ≤3.0 area%

Uridine (HPLC): ≤0.5 area%

AMP (HPLC): ≤0.5 area%

ATP (HPLC): ≤0.2 area%

UTP (HPLC): ≤0.2 area%

Total unknown impurities (HPLC): ≤0.5 area%

Heavy metals (as Pb): ≤10 ppm

P_i: ≤0.4%

Fe (AAS): ≤10 ppm

Ethanol (GC): ≤5.5%

2-Propanol (GC): ≤100 ppm

Methanol (GC): ≤50 ppm

Acetone (GC): ≤100 ppm

Bioburden:

Total amount: ≤100 CFU/g

Mould: ≤50 CFU/g

Yeasts: ≤50 CFU/g

Stability: At +2 to +8°C within specification range for 12 months.

Cat. No. **Pack Size**

04 589 173 103 custom fill

Will be supplied as "UDP-Galactose, Di-Na". Unit of Measure is "kg".

For further processing only.

UDP-Glucose

disodium salt, powder

UDP-Glucose is an activated sugar.

Application

Use UDP-Glucose for glycosylation of target molecules with suitable enzyme.

CAS: 28053-08-9

Properties

Formula: $C_{15}H_{22}N_2O_{17}P_2Na_2$

Molecular weight:

UDP-glucose: 566.3 D

UDP-glucose- Na_2 : 610.3 D

Specification

Appearance: White to yellowish powder

Solubility: Clear, colorless solution in water (c=50 mg/ml)

pH value (solution, c=50 mg/ml): 5.8-7.2

Identity (NIR spectrum): Corresponds to reference

UDP-Glucose (enzymatically): $\geq 83.5\%$

UDP-Glucose (A_{260} , $\epsilon=9.9 \text{ l x mmol}^{-1} \text{ x cm}^{-1}$): $\geq 83.5\%$

UDP-Glucose, disodium salt (HPLC, based on dry weight): $\geq 95.0 \text{ area}\%$

Na (flame photometric): 7.0-8.0%

Water (K. Fischer): $\leq 6.0\%$

UDP (HPLC): $\leq 0.5 \text{ area}\%$

UMP (HPLC): $\leq 0.5 \text{ area}\%$

UTP (HPLC): $\leq 0.2\%$

UDP-Galactose (HPLC): $\leq 1.5 \text{ area}\%$

Uridine (HPLC): $\leq 0.2 \text{ area}\%$

Singleunknown impurity (HPLC): $\leq 0.5 \text{ area}\%$

Total unknown impurities (HPLC): $\leq 3.0 \text{ area}\%$

Sum impurities (HPLC): $\leq 4.5 \text{ area}\%$

Heavy metals (as Pb): $\leq 10 \text{ ppm}$

P_i: $\leq 0.3 \%$

Fe (bathophenanthroline): $\leq 25 \text{ ppm}$

Ethanol (GC): $\leq 1,000 \text{ ppm}$

2-Propanol (GC): $\leq 100 \text{ ppm}$

Methanol (GC): $\leq 50 \text{ ppm}$

Acetone (GC): $\leq 100 \text{ ppm}$

A_{250}/A_{260} : 0.71-0.75

A_{250}/A_{280} : 0.38-0.42

A_{290}/A_{260} : 0.03-0.05

Stability: At +2 to +8°C within specification range for 36 months.

Cat. No. **Pack Size**

10 154 938 103 custom fill

Will be supplied as "UDP-glucose, Di-Na". Unit of Measure is "kg".
For further processing only.

β -NAD free acid, crystalline powder

β -NAD is a cofactor for redox enzymes.

Application

NAD is involved in redox reactions, carrying electrons from one reaction to another and can also act as a cofactor for enzymatic reactions such as dehydrogenases.

Benefits

- Use this animal component- and gluten-free product in the food and pharmaceutical industry.
- Rely on large lot sizes and high lot-to-lot consistency for use in technical applications.

CAS: 53-84-9

Properties

Nomenclature: Nicotinamide adenine dinucleotide

Formula: $C_{21}H_{27}N_7O_{14}P_2 \times 3 H_2O$

Molecular weight:

NAD: 663.4 D

NAD x 3 H₂O: 717.4 D

Specification

Appearance: White, crystalline powder

β -NAD x 3 H₂O (calculated from value, enzymatically): $\geq 99\%$

β -NAD (enzymatically, A_{340} , $\epsilon=6.3 \text{ l x mmol}^{-1} \times \text{cm}^{-1}$): $\geq 92\%$

NAD (CN complex): $\geq 92\%$

NAD (A_{260} , $\epsilon=17.6 \text{ l x mmol}^{-1} \times \text{cm}^{-1}$): $\geq 92\%$

NAD (HPLC): ≥ 98 area%

Water (K. Fischer): 4-8%

AMP (enzymatically): $\leq 0.1\%$

Acetone (GC): $\leq 0.05\%$

Ethanol (GC): $\leq 0.05\%$

Isopropanol (GC): $\leq 0.05\%$

Methanol (GC): $\leq 0.05\%$

A_{250}/A_{260} : 0.81-0.85

A_{280}/A_{260} : 0.20-0.24

Stability: At +2 to +8°C within specification range for 24 months.

Cat. No. **Pack Size**

10 768 197 103 custom fill

Will be supplied as "beta-NAD, free acid". Unit of Measure is "kg".
For further processing only.

NADH, Food Grade

disodium salt, lyophilizate

NADH, Food Grade, is a cofactor for redox enzymes.

Application

Use NADH, Food Grade, as a cofactor for redox enzymes for industrial biotechnology and also for the determination of pyruvate, LDH, NH_3 , GIDH, MDH, and aldehyde.

Benefits

- Use this animal component-free product in the food and pharmaceutical industry.
- Rely on large lot sizes and high lot-to-lot consistency for use in technical applications.

CAS: 606-68-8

Properties

Nomenclature: β -Nicotinamide-adenine-dinucleotide, reduced

Formula:

β -NADH- Na_2 : $\text{C}_{21}\text{H}_{27}\text{N}_7\text{O}_{14}\text{P}_2\text{Na}_2$

β -NADH: $\text{C}_{21}\text{H}_{27}\text{N}_7\text{O}_{14}\text{P}_2$

Molecular weight:

β -NADH- Na_2 : 709.4 D

β -NADH: 665.4 D

Specification

Appearance: White to slightly yellowish lyophilizate

Solubility: Clear, colorless to slightly yellowish solution in water (c=50 mg/ml)

pH value (c=100 mg/ml): 8-9

NADH- Na_2 (enzymatically, based on dry weight): $\geq 98\%$

NADH (enzymatically): $\geq 85\%$

NADH (A_{260} , $\epsilon = 14.3 \text{ l x mmol}^{-1} \text{ x cm}^{-1}$): $\geq 85\%$

NADH (HPLC): ≥ 95 area%

Na (flame photometric): 6-7%

Water (K.Fischer): $\leq 6\%$

Ethanol (GC): $\leq 0.2\%$

NAD (enzymatically): $\leq 2\%$

AMP (enzymatically): $\leq 0.2\%$

Heavy metals (as Pb): ≤ 10 ppm

Bioburden:

Total amount: $\leq 10^4$ CFU/g

Coliforme: $\leq 10^2$ CFU/g

Yeasts and moulds: $\leq 3 \times 10^2$ CFU/g

Germ differentiation (*Staphylococcus aureus*, *Salmonellae*, *E. coli*): Negative

Stability: At -15 to -25°C within specification range for 24 months.

Cat. No. **Pack Size**

03 277 372 103 custom fill

Will be supplied as "NADH, Di-Na, Food grade, Lyo.". Unit of Measure is "kg".



For further processing only.

BM Condimed H1 0.2µm filtered, solution

BM Condimed H1 is a culture medium supplement for proliferation of hybridoma cells.

Application

Use BM Condimed H1 as a supplement to the normal culture medium to support the growth of B-cell hybridomas after fusion and during cloning. BM Condimed H1 replaces feeder cells. It is also used to optimize the growth of hybridomas after the thawing of cells stored in liquid nitrogen.

Working concentration: 10% (v/v) in culture medium that contains 10-20% FBS.

Benefits

- **Obtain increased proliferation rate** of your freshly fused hybridoma cells.
- **Obtain increased size and numbers** of colonies.
- **Rely on the unique composition** that makes feeder cells in your application unnecessary.

Product Description

BM Condimed H1 is prepared from the supernatant of a mouse lymphoma cell line stimulated with phorbol myristate acetate (PMA). It contains a complex mixture of growth factors and cytokines that have a marked stimulatory effect on the growth of hybridoma cells after fusion and during cloning.

Specification

Appearance: Clear, reddish solution

Contents: BM Condimed H1 is supplied as 0.2 µm filtered solution in RPMI 1640. The solution also contains FCS (fetal calf serum), 15% (v/v); oxalacetate, 1mmol/l; sodium pyruvate, 1mmol/l; insulin, 0.2 µg/ml; hIL-6, 1 ng/ml; phorbol myristate acetate, 10 ng/ml; phenol red.

Function (cell culture): Corresponds to reference

0.2 µm filtration (tested according to approved microbiological methods): Corresponds to reference

Stability: At -15 to -25°C within specification range for 24 months.

Quality

BM Condimed H1 is assayed for its ability to promote the growth of freshly fused hybridoma cells.

Background Information

BM Condimed H1 is added as a supplement (10%, v/v) to normal culture medium (basal medium, *e.g.* RPMI 1640, DMEM, IMDM) that also contains 10-20% FCS. Such a medium can support the growth of B-cell hybridomas, both after fusion and during cloning.

BM Condimed H1 should not be used at higher concentration, as basal medium or as a replacement for serum.

Cat. No. **Pack Size**

04 155 645 001 custom fill

Will be supplied as "BM Condimed H1, Bulk". Unit of Measure is "Liter".



For further processing only.

Chromozym TRY powder

Chromozym TRY is used for the determination of trypsin activity.

Application

Substrate for the reliable photometric determination of activity of proteases which hydrolyze peptides at the carboxylic side of arginine (trypsin, endoproteinase Arg-C and others).

Benefits

- Perform precise trypsin activity determination.
- Decrease variance between different measurements.
- Enjoy the easy handling, especially in routine testing.

Product Description

Synthetic substrate showing low variance for different activity measurements. The increase in absorption at 405 nm is easily measured.

CAS: 52299-14-6

Properties

Nomenclature: Carbobenzoxy-L-valyl-L-glycyl-L-arginine-4-nitranilide acetate

Formula: $C_{27}H_{36}N_8O_7$

Molecular weight: 584.6 D

Formula: $C_{27}H_{36}N_8O_7 \times AcOH$

Molecular weight: 644.7 D

Specification

Appearance: White to yellowish powder

Chromozym TRY, acetate: $\geq 89\%$

Chromozym TRY (enzymatically): $\geq 80\%$

Water (K. Fischer): $\leq 5\%$

4-Nitraniline, free: $\leq 0.5\%$

Stability: At +15 to +25°C within specification range for 24 months.

Background Information

A detailed test procedure is available upon request.

Cat. No. _____ **Pack Size** _____

10 378 496 103 custom fill

Will be supplied as "Chromozym TRY". Unit of Measure is "g".
For further processing only.

Glutathione, oxidized form (GSSG) lyophilizate

Cat. No.	Pack Size
10 151 327 103	custom fill

Application

Use GSSG for *in vitro* renaturation of proteins in inclusion bodies, and also as a reaction partner of GSH.

Benefits

- Use an animal component-free product to simplify approval processes.
- Rely on high lot-to-lot consistency and large scale production.

CAS: 27025-41-8

Properties

Formula: C₂₀H₃₂N₆O₁₂S₂

Molecular weight: 612.6 D

Specification

Appearance: White lyophilizate

Purity (enzymatically): ≥90%

Water (K. Fischer): ≤5.0%

Glutathione, reduced form (enzymatically): ≤0.5%

Fe (AAS): ≤20 ppm

Heavy metals (as Pb, AAS): ≤10 ppm

Stability: At +2 to +8°C within specification range for 36 months.

Will be supplied as "Glutathione, Oxidized Form (GSSG)". Unit of Measure is "g".

For further processing only.

Glutathione, reduced form (GSH) crystalline powder

Cat. No.	Pack Size
10 002 801 103	custom fill

GSH is used as a food supplement or fermentation media compound.

Application

Use GSH as an antioxidative to significantly reduce free radicals.

Benefits

- **Rely on high lot to lot consistency** and large scale production.

CAS: 70-18-8

Properties

Formula: C₁₀H₁₇N₃O₆S

Molecular weight: 307.3 D

Specific rotation [α]_D²⁵ (water): -17° to -21°

Specification

Appearance: White, crystalline powder

Solubility: Clear, colorless in water (c=50 mg/ml)

A₄₀₅ (hydrous solution, alkaline; 4.0 ml + 0.2 ml NaOH, 2 mol/l): ≤0.050

Melting range: +189 to +200°C

Purity (enzymatically and iodometric): ≥98%

Loss on drying (+105°C, 2 hours): ≤1%

Glutathione, oxidized form (enzymatically): ≤1.5%

Fe (bathophenanthroline): ≤5 ppm

Heavy metals (as Pb, AAS): ≤5 ppm

Methanol (GC): ≤50 ppm

Aceton (GC): ≤50 ppm

Ethanol (GC): ≤0.1%

Stability: At +2 to +8°C within specification range for 24 months.

Will be supplied as "Glutathione, Reduced Form (GSH)". Unit of Measure is "kg".

For further processing only.

GTP**disodium salt, powder**

For best and consistent performance rely on long-term experience from the leading manufacturer of nucleotides.

Application

Use GTP as a coenzyme for protein biosynthesis in a cell-free system.

Benefits

- **Obtain consistent results.** Rely on the excellent lot-to-lot performance of this product.

CAS: 56001-37-7

Properties

Formula: $C_{10}H_{14}N_5O_{14}P_3Na_2$

Molecular weight:

GTP: 523.2 D

GTP- Na_2 : 567.1 D

Specification

Appearance: White powder

GTP (enzymatically): $\geq 74\%$

GTP (A_{252}^1 , $\epsilon = 14.9 \text{ l x mmol}^{-1} \text{ x cm}^{-1}$): $\geq 74\%$

GTP (HPLC): ≥ 90 area%

Na (flame photometric): $8 \pm 1\%$

Water (K. Fischer): $\leq 10\%$

P_i: $\leq 0.9\%$

A₂₅₀¹/A₂₆₀¹: 1.15 ± 0.03

A₂₈₀¹/A₂₆₀¹: 0.66 ± 0.02

A₂₉₀¹/A₂₆₀¹: 0.28 ± 0.01

Stability: At +2 to +8°C within specification range for 12 months. Store dry.

Cat. No.**10 150 398 103****Pack Size**

custom fill

Will be supplied as "Guanosine-5'-triphosphate (GTP), Di-Na". Unit of Measure is "g".



For further processing only.

Pepstatin

from *Streptomyces* species, lyophilizate

Pepstatin is a protease inhibitor.

Application

Use Pepstatin to inhibit potently the HIV protease and other aspartic proteases such as pepsin, renin, cathepsin D, chymosin, and many microbial acid proteases.

Benefits

- **Rely on an outstanding purity** of ≥98% according to HPLC analysis.
- **Increase the safety of your application processes** with reproducible performance and high lot-to-lot consistency.

Product Description

Pepstatin is a highly efficient inhibitor for many aspartic and acid proteases. It is a hexa-peptide containing the unusual amino acid statine (Sta, (3S,4S)-4-amino-3-hydroxy-6-methylheptanoic acid), having the sequence Iva-Val-Val-Sta-Ala-Sta. Pepstatin is known for its ability to inhibit pepsin in a picomolar concentration. The product is soluble in methanol and ethanol, but insoluble in water.

EC 3.4.14.9

Properties

Formula: C₃₄H₆₃N₅O₉

Molecular weight: 685.9 D

Pepstatin is soluble in methanol (1 mg/ml and 20 mg/ml as well) and ethanol (1 mg/ml) if you allowed to sit overnight, and to 300 µg/ml in 6 N acetic acid. Stable at least 1 week at +2 to +8°C or 1 month if stored in aliquots at -15 to -25°C.

Specification

Appearance: White powder

Melting range: 227-231°C

Solubility:

Clear to slightly opalescent, colorless solution in methanol (c=1 mg/ml)

Clear, colorless solution in methanol (c=20 mg/ml)

Clear to slightly opalescent, colorless solution in ethanol (overnight)

Purity (HPLC): ≥98 area%

Performance of inhibition of pepsin:

With 0.01 ml sample: ≥85%

With 0.1 ml sample: ≥100%

Values are taken from supplier certificate.

Stability: At +2 to +8°C within specification range for 24 months.

Quality

Pepstatin is highly purified according to HPLC analysis.

Cat. No. **Pack Size**

10 253 294 103 custom fill

Will be supplied as "Pepstatin from *Streptomyces* species". Unit of Measure is "g".

For further processing only.

Patent and License Disclaimers

01. Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under such patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patent claims require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
04. A license to perform the patented 5' Nuclease Process for research is obtained by the purchase of (i) both Authorized 5' Nuclease Core Kit and Licensed Probe, (ii) a Licensed 5' Nuclease Kit, or (iii) license rights from Applied Biosystems. This product is an Authorized 5' Nuclease Core Kit. Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the applicable claims of US patents and corresponding claims outside the United States. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patent claims require a separate license from Roche. Further information on purchasing licenses may be obtained from the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
05. Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
48. Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patent claims require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
49. Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, nontransferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patent claims require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
51. No license is conveyed with the purchase of this product under any US patent claims and corresponding patent claims outside the United States, or any other patents or patent applications, relating to the 5' Nuclease and dsDNA-Binding Dye Processes. For further information contact the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
62. NOTICE TO PURCHASER: DISCLAIMER OF LICENSE This product is optimized for use in the polymerase chain reaction (PCR) covered by patents owned by F. Hoffmann-La Roche Ltd ("Roche"). No license under these patents to use the PCR Process is conveyed expressly or by implication to the purchaser by the purchase of this product. A license to use the PCR Process for certain research and development activities accompanies the purchase of certain Roche, Applied Biosystems or other licensed suppliers' reagents when used in conjunction with an authorized thermal cycler, or is available from Applied Biosystems. Diagnostic purposes require a license from Roche. Further information on purchasing licenses may be obtained by contacting the Director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.
63. NOTICE TO PURCHASER: LIMITED LICENSE Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patent claims require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Licensing Department of Roche Molecular Systems, Inc., 4300 Hacienda Drive, Pleasanton, California 94588 or Roche Diagnostics GmbH, Sandhofer Strasse 116, 68305 Mannheim, Germany.

64. NOTICE TO PURCHASER: LIMITED LICENSE Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patent claims require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Licensing Department of Roche Molecular Systems, Inc., 4300 Hacienda Drive, Pleasanton, California 94588 or Roche Diagnostics GmbH, Sandhofer Strasse 116, 68305 Mannheim, Germany.
65. NOTICE TO PURCHASER: LIMITED LICENSE Use of this product is covered by US patent claims and corresponding patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses require a separate license from Roche. Further information on purchasing licenses for diagnostic applications may be obtained by contacting the Licensing Department of Roche Molecular Systems, Inc., 4300 Hacienda Drive, Pleasanton, California 94588, USA.
66. NOTICE TO PURCHASER: LIMITED LICENSE The purchase price of this product includes a limited, non-transferable license under US patent claims or non-U.S. counterparts owned by Roche Molecular Systems, Inc. or F. Hoffmann-La Roche Ltd to use only this amount of product solely for the research and development activities of the purchaser. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. Further information on purchasing licenses to practice these patented processes may be obtained by contacting the Licensing Department of Roche Molecular Systems, Inc., 4300 Hacienda Drive, Pleasanton, California 94588 or Roche Diagnostics GmbH, Sandhofer Strasse 116, 68305 Mannheim, Germany.

The Following trademarks are trademarks of Roche

ABTS	EAGLETAQ	HIGH PURE	MAGNA PURE LC	STREPTAWELL
ACTITAQ	ENZYMUN-TEST	HYBPROBE	MGRADE	TAQMAN
AMPLICOR	ELECSYS	LC	MYCOTOOL	UNITOOL
AMPLITAQ	EXPAND	LIBERASE	NUCLEOMIX	X-TREMEGENE
APTATAQ	FASTSTART	LIGHTCYCLER	QUICKTAG	
BM CONDIMED	GENOPURE	LUMIGRADE	REALTIME READY	
COBAS	HAWKTAQ	MAGNA LYSER	RESOLIGHT	
COMPLETE and design	HAWKZO5	MAGNA PURE	SIMPLEPROBE	

Other Trademarks

ABI PRISM, CDP-Star, CSPD, StepOnePlus are trademarks of Life Technologies Corporation or its subsidiaries in the US and/or certain other countries.

SYBR is a registered trademark of Molecular Probes, Inc.

The ATCC trademark and trade name and any and all ATCC catalog numbers are trademarks of the American Type Culture Collection.

Tween is a registered trademark of ICI Americans, Inc., USA

Chromozym is a registered trademark of Pentapharm AG, Basle, Switzerland

CY and Sephadex are trademarks of GE Healthcare.

Eppendorf is a registered trademark of Eppendorf-NetherlerHinz GmbH.

FPLC is a registered trademark of Pharmacia Biotech AB, Uppsala, Sweden.

Genapol is a registered trademark of Hoechst AG, Frankfurt/M., Germany.

Triton and Kathon are trademarks of Rohm & Haas Company, Philadelphia, PA, USA.

LNA, ProbeFinder, and ProbeLibrary are registered trademarks of Exiqon A/S, Vedbaek, Denmark.

Nonidet is a registered trademark of Shell International Petroleum Company Limited, U.K.

Pefabloc is a trademark of DSM IP Assets B.V., Heerlen, Netherlands.

Thesit a registered trademark of Desitin-Werk, Carl Klinke GmbH, Hamburg, Germany.

Threshold is a registered trademark of Molecular Devices, UK.

xMAP is a registered trademark of Luminex Corporation, Austin, USA

Other brands or product names are trademarks of their respective holders

All business related to the sale of the Products listed in this catalog shall be subject to the following general **TERMS OF DELIVERY AND PAYMENT**. These terms may be superseded by terms set by the sales organization in your country. It is the responsibility of the Purchaser to inquire about the respective terms.

Orders

All orders and agreements are binding. When ordering, please be certain to indicate both catalog number and description of product. Minimum order value is U.S.\$ 1,000.- or the equivalent in foreign currency. The Supplier retains the right to change the minimum order value in individual cases.

Prices

The prices in effect on the day the goods are ordered shall be applicable. Should the Purchaser object to this, the Supplier shall be entitled to unilaterally cancel the contract either in whole or in part, as far as quantities pending delivery are concerned.

Delivery

Any delivery time shall not be binding on the Supplier, but the Supplier shall use its best efforts to adhere to it. The Supplier shall be entitled to postpone or cancel his delivery obligations in whole or in part, in the event of interruptions of delivery caused by Force Majeure such as, but not limited to Acts of God, disturbances of company operation, strikes, or unavailability of transportation facilities. The Purchaser shall not be entitled to any claim for damages caused by delayed delivery.

Terms of Payment

Payment must be made net 30 days from the date of invoice. The Supplier retains the right to change these Terms of Payment upon mutual agreement with the Purchaser. All duties will be paid by the Purchaser.

Claims for Damage in Transit or Shortages

The Supplier must be informed about shipping discrepancies or damages in transit within 3 days of receipt of goods. In addition, the shipping company must be informed about damages in transit immediately upon receipt of goods, if such damages are concealed at the latest within 3 days upon receipt of goods.

In case of shipping discrepancies or damages in transit the Purchaser shall not be entitled to make deductions from the sales price or to hold back payment of same. Claims that are acknowledged by Supplier to fall under its responsibility shall be restricted to delivery of a satisfactory substitute.

Returning Material

No product may be returned without authorization from the Supplier. This authorization is necessary to ensure correct return of material and issuance of credit or exchange of material. Requests for permission to return products must be made within 10 days of receipt. Returned products should be clearly marked with sender's name and address, carefully packed and shipped prepaid. All re-

turned products must be directed to the Supplier in writing. Products must be returned in the same condition as received and within 15 days of original shipment by Supplier.

WARRANTY

The Supplier warrants its products to the original purchaser against defects in materials and workmanship under normal use and application. The Supplier's sole obligation under this warranty shall be to replace defective products.

ALL PRODUCTS ARE DEVELOPED, DESIGNED AND SOLD FOR RESEARCH PURPOSES AND/OR USE IN MANUFACTURING. This warranty applies only to products in the original container that have been handled and stored in accordance with the Supplier's recommendations or instructions, and does not apply to a product which has been the subject to alteration, misuse or abuse. All claims under this warranty must be directed to the Supplier in writing and must be accompanied by a copy of the Purchaser's invoice.

THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING THE WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, IN NO CASE SHALL THE SUPPLIER BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES. Any modification of the above terms will be valid only if confirmed in writing by the Supplier at its corporate offices.

Retention of Title

All merchandise shall remain the property of the Supplier until all debts outstanding, including all current accounts receivable from the Purchaser have been paid in full.

IMPORTANT INFORMATION

Our preparations listed in this catalog are intended exclusively for further use in manufacturing, unless indicated otherwise. They are not for use on humans and animals.

The absence of a toxicity warning with one of our products does not, however, preclude a possible health hazard. With all of our products, due care should be exercised to prevent human contact and ingestion. All preparations should be handled only by trained personnel. We reserve the right to change product specifications (e.g. buffer systems) as we see fit to improve product quality.

The sale of the Products to Purchaser does not grant Purchaser, unless indicated otherwise, any rights to or license on any patent application, patent or other intellectual proprietary right of Roche or its Affiliates. Thus, it is Purchaser's responsibility to identify and acquire such a right or license.

Further, it is Purchaser's obligation to check whether Purchaser's use or other commercialization of the Products requires any third party licenses and to acquire such a license. Purchaser shall indemnify Roche and its Affiliates from any third party claims resulting from a possible infringement of any third party intellectual property rights.

1

1,4-Dithiothreitol	
⊙ 1,4-Dithiothreitol (DTT) crystallizate	
1,4-Dithiothreitol (DTT) crystallizate	31

2

2,2'-azino-bis(3-ethylbenzthiazoline-6-sulphonic acid (ABTS) solution	212
2,2'-Azino-di-[3-ethylbenzthiazoline sulfonate (6)]	
⊙ ABTS 2,2'-Azino-di-[3-ethylbenzthiazoline sulfonate (6)] diammonium salt	
2-Chloroacetamide (CAA) crystalline powder	7
2-Hydroxypyridine-N-oxide (Oxy-PYRION) crystalline powder	7
2-Hydroxypyridine-N-oxide (Oxy-PYRION), reduced sodium crystalline powder	8

3

3,3',5,5'-Tetramethylbenzidine (TMB) crystalline powder	213
3,5-Dichlorophenolsulfonic Acid disodium salt	139
3-[(3-Cholamidopropyl)dimethyl-ammonio]-1-propane sulfonate, zwitterionic	
⊙ CHAPS zwitterionic detergent, crystalline powder	
3-Hydroxy-1,2,3,4-tetrahydrobenzo[h]quinoline crystalline powder	31

4

4-(2-Aminoethyl)-benzene-sulfonyl fluoride	
⊙ Pefabloc SC (AEBSF) 4-(2-Aminoethyl)-benzenesulfonyl fluoride hydrochloride, powder	
4-Aminoantipyrine (4-APP) crystalline powder	139
4-Aminophenyl Phosphate (pAPP) disodium salt, powder	213
4-Nitrophenyl Phosphate (4-NPP) di-Tris salt	141
4-Nitrophenyl Phosphate (pNPP) disodium salt, crystalline powder	214
4-Nitrophenyl-α-D-maltohexaoside powder	140
4-Nitrophenyl-α-D-maltopentaoside powder	140

5

5-Aminoallyl-ddUTP powder, lithium salt	312
5-Aminoallyl-dUTP powder, lithium salt	311
5-Aminoallyl-UTP powder, lithium salt	311
5-Bromo-4-chloro-3-indolyl-phosphate (BCIP) disodium salt, crystalline powder	214
5-Bromo-4-chloro-3-indolyl-phosphate (BCIP) toluidin, crystalline powder	215
5-Bromo-5-Nitro-1,3-Dioxane (BND), Grade I crystalline powder	9
5-Bromo-5-Nitro-1,3-Dioxane (BND), Grade II crystalline powder	10

6

6-Phosphogluconolactonase from <i>Leuconostoc mesenteroides</i>, lyophilizate	130
--	-----

7

7-Deaza-2'-dGTP lithium salt, 10 mM	277
--	-----

G

α-Glucosidase from yeast overproducer, multifunctional, lyophilizate	96
α-Ketoglutarate (2-Oxoglutarate) disodium salt	148
α-Ketoglutarate (2-Oxoglutarate) disodium salt, dihydrate	148
α-Ketoglutarate (2-Oxoglutarate) for potassium test free acid	149
α-Ketoglutarate (2-Oxoglutarate) free acid	147
α-S-dCTP, Molecular Diagnostic Grade S-isomer, sodium salt, 100 mM	275

A

ABTS 2,2'-Azino-di-[3-ethylbenzthiazoline sulfonate (6)] diammonium salt	141
Acetate-CoA Ligase (Acetyl-CoA Synthetase) from microorganism, lyophilizate	58
Acetate-CoA Ligase (Acetyl-CoA Synthetase) from yeast, lyophilizate	58
Acetyl-CoA	
⊙ Acetyl-Coenzyme A trilithium salt	
Acetyl-Coenzyme A trilithium salt	35
Acetyl-Coenzyme A trisodium salt	35
Acid Phosphatase from potato, lyophilizate	59
ActiQaq Dexo DNA Polymerase from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i>, solution	242
Acyl-CoA Oxidase from microorganisms, lyophilizate	59
Adenosin-5'-diphosphate	
⊙ ADP disodium salt	
Adenosine Deaminase from calf intestine, suspension	59
Adenosine-5'-O-(2-thiodiphosphate) trilithium salt	45
Adenosine-5'-O-(3-thiotriphosphate) tetralithium salt	45
ADP disodium salt	46
ADP for potassium test free acid	47
ADP potassium salt	46
AEBSF	
⊙ Pefabloc SC (AEBSF) 4-(2-Aminoethyl)-benzenesulfonyl fluoride hydrochloride, powder	
Alanine Aminotransferase	
⊙ Alanine Aminotransferase (ALT) (GPT) from pig heart, suspension	
Alanine Aminotransferase (ALT) (GPT) from pig heart, lyophilizate	60
Alanine Aminotransferase (ALT) (GPT) from pig heart, suspension	60
α-Albumins	
↑ Bovine Serum Albumin (BSA), Fraction V fatty acids ≤ 0.2 mg/g, lyophilizate	
↑ Bovine Serum Albumin (BSA), Fraction V lyophilizate	
↑ Bovine Serum Albumin, Molecular Biology Grade 2% solution	
↑ Bovine Serum Albumin, Molecular Biology Grade 10% solution	
Alcohol Dehydrogenase from yeast, lyophilizate	61
Alcohol Dehydrogenase from yeast, suspension	63

Alcohol Dehydrogenase, chemically modified from yeast, lyophilizate	62
Aldehyde Dehydrogenase from yeast, lyophilizate	64
Aldose 1-Epimerase (Mutarotase) from pig kidney, suspension	65
Alkaline Phosphatase Mutein, recombinant from calf intestine, expressed in <i>Pichia pastoris</i> , lyophilizate	204
Alkaline Phosphatase, Cloning of DNA Fragments ◉ Alkaline Phosphatase, recombinant, 1 U/μl from bovine intestine, expressed in <i>Pichia pastoris</i> , solution	
Alkaline Phosphatase, EIA Grade from calf intestine, solution	204
Alkaline Phosphatase, recombinant, 1 U/μl from bovine intestine, expressed in <i>Pichia pastoris</i> , solution	329
Alkaline Phosphatase, recombinant, 20 U/μl from bovine intestine, expressed in <i>Pichia pastoris</i> , solution	330
Alkaline Phosphatase, recombinant, highly active from calf intestine, expressed in <i>Pichia pastoris</i> , solution	205
Alkaline Phosphatase, recombinant, highly active, carbohydrate reduced from calf intestine, expressed in <i>Pichia pastoris</i> , solution	206
↔ Alkaline Phosphatase, shrimp, alternative ↑ Alkaline Phosphatase, recombinant, 1 U/μl from bovine intestine, expressed in <i>Pichia pastoris</i> , solution	
AllStart RNA Master	271
AllStart RNA Master (Rox)	272
Amino-2-(hydroxymethyl)-1,3-propanediol, 2- ◉ Tris crystallizate	
Aminoethyl-benzene-sulfonyl fluoride, 4-2- ◉ Pefabloc SC (AEBSF) 4-(2-Aminoethyl)-benzenesulfonyl fluoride hydrochloride, powder	
AMP disodium salt	49
AMP free acid	48
AMV Reverse Transcriptase, recombinant, GMP Grade from Avian Myeloblastosis Virus, expressed in <i>E. coli</i>	267
Anti-Digoxigenin-AP, Fab fragments polyclonal antibodies from sheep, lyophilizate	305
Anti-Digoxigenin-POD, Fab fragments polyclonal antibodies from sheep, lyophilizate	306
Anti-Digoxigenin-Rhodamine, Fab fragments polyclonal antibodies from sheep, lyophilizate	307
Aprotinin from bovine lung, lyophilizate	30
AptaTaq DNA Master 5x concentrated	260
AptaTaq DNA Master Optimization Kit	261
AptaTaq DNA Master without Mg ²⁺ 5x concentrated	261
AptaTaq DNA Polymerase LDx, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , solution	245
AptaTaq DNA Polymerase LDx, 50 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , glycerol-free solution	247
AptaTaq DNA Polymerase, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , solution	243
AptaTaq DNA Polymerase, 50 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , glycerol-free solution	244
AptaTaq Genotyping Master (Rox) 5x concentrated	263
AptaTaq Genotyping Master 5x concentrated	262
AptaTaq Dexo DNA Polymerase, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , solution	249

AptaTaq Dexo DNA Polymerase, 50 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , glycerol-free solution	250
Ascorbate Oxidase from <i>Cucurbita</i> species, lyophilizate	66
Ascorbate Oxidase from <i>Cucurbita</i> species, poor of Aspartate aminotransferase (AST/GOT), lyophilizate	67
Ascorbate Oxidase, chemically modified from <i>Cucurbita</i> species, lyophilizate	68
Aspartate Aminotransferase (AST) (GOT) from pig heart, lyophilizate	69
Aspartate Aminotransferase (AST) (GOT) from pig heart, suspension	69
Aspartate Aminotransferase - GOT ◉ Aspartate Aminotransferase (AST) (GOT) from pig heart, suspension	
ATP, Grade I disodium salt	49
ATP, Grade II disodium salt	50
ATP, Molecular Diagnostic Grade sodium salt, 100 mM	290
ATP: D-hexose 6-phosphotransferase <i>EC 2.7.1.1</i> ◉ Hexokinase (HK) from yeast overproducer, lyophilizate	
Azino-di-[3-ethylbenzthiazoline sulfonate (6)], 2,2'- ◉ ABTS 2,2'-Azino-di-[3-ethylbenzthiazoline sulfonate (6)] diammonium salt	

β

β-D-Glucose: oxygen-1-oxidoreductase <i>EC 1.1.3.4</i> ◉ Glucose Oxidase (GOD), Grade II from <i>Aspergillus niger</i> overproducer, lyophilizate	
β-Galactosidase from <i>E. coli</i> , lyophilizate	207
β-Galactosidase Mutein from <i>E. coli</i> overproducer, lyophilizate	208
β-Galactosidase, recombinant, EIA Grade from <i>E. coli</i> overproducer, lyophilizate	209
β-Glucuronidase from <i>E. coli</i> , solution	97
β-NAD free acid, crystalline powder	365

B

Benzylidene-4-NP-G7 4,6-Benzylidene-4-nitrophenyl-α-D-maltoheptaoside, lyophilizate	142
Biotin-11-CTP 10 mM solution	313
Biotin-16-ddUTP 1 mM solution	315
Biotin-16-dUTP 1 mM solution	315
Biotin-16-dUTP 20 mM solution	314
Biotin-16-UTP 10 mM solution	313
Blocking Reagent powder	336
BM Blue POD Substrate, soluble 3,3',5,5'-Tetramethylbenzidine (TMB), solution	213
BM Condimed H1, 0.2 μm filtered, solution	367
Bovine IgG (PAB<->R-IgG) lyophilizate	199
Bovine Serum Albumin (BSA), Fraction V fatty acids ≤0.2 mg/g, lyophilizate	3
Bovine Serum Albumin (BSA), Fraction V lyophilizate	2
Bovine Serum Albumin (BSA), reduced sodium and potassium lyophilizate	4

Bovine Serum Albumin I lyophilizate	199
Bovine Serum Albumin IV lyophilizate	200
Bovine Serum Albumin, Molecular Biology Grade 10% solution	335
Bovine Serum Albumin, Molecular Biology Grade 2% solution	334
~ Bromo-4-chloro-3-indolyl phosphate (BCIP) ↑ 5-Bromo-4-chloro-3-indolyl-phosphate (BCIP) toluidin, crystalline powder	
C	
CaA, grade II, trilitium salt	
○ Coenzyme A, Grade II trilitium salt	
Carboxypeptidase B, recombinant from rat pancreas, expressed in <i>Pichia pastoris</i> , solution	348
Catalase from <i>Corynebacterium glutamicum</i> , lyophilizate	346
Catalase from <i>Corynebacterium glutamicum</i> , solution	347
CHAPS zwitterionic detergent, crystalline powder	28
CHAPSO zwitterionic detergent, crystalline powder	28
Chlorophenolred-beta-D-galactopyranoside (CPRG) sodium salt, powder	212
Cholamidopropyl-dimethyl-ammonio]-1-propane sulfonate, 3-3-, zwitterionic	
○ CHAPS zwitterionic detergent, crystalline powder	
Cholate ionic detergent, sodium salt	22
Cholesterol Esterase from <i>Candida cylindracea</i> , lyophilizate	73
Cholesterol Esterase from <i>Candida cylindracea</i> , solution	74
Cholesterol Esterase from <i>Pseudomonas species</i> , lyophilizate	71
Cholesterol Esterase from <i>Pseudomonas species</i> , solution	72
Cholesterol Esterase, chemically modified from <i>Pseudomonas species</i> , lyophilizate	71
Cholesterol Oxidase from <i>Brevibacterium species</i> , expressed in <i>E.coli</i> , lyophilizate	75
Cholesterol Oxidase from <i>Nocardia erythropolis</i> , lyophilizate	76
Cholesterol Oxidase from <i>Nocardia erythropolis</i> , solution	78
Cholesterol Oxidase from <i>Streptomyces species</i> , lyophilizate	79
Chromogenic Substrate for Lipase	142
Chromozym TRY powder	368
CIP (Calf Intestinal Phosphatase)	
○ Alkaline Phosphatase, recombinant, 20 U/μl from bovine intestine, expressed in <i>Pichia pastoris</i> , solution	
Citrate Lyase from <i>Klebsiella pneumoniae</i> , lyophilizate	80
Citrate Synthase from pig heart, suspension	80
Cleland's Reagent	
○ 1,4-Dithiothreitol (DTT) crystallizate	
CMP-N-Acetylneuraminic Acid sodium salt	361
CoA, grade I, free acid	
○ Coenzyme A, Grade I free acid	
CoA, grade I, trilitium salt	
○ Coenzyme A, Grade I trilitium salt	
Cobalt Chloride Solution 25 mM solution	310
Coenzyme A, Grade I free acid	36

Coenzyme A, Grade I trilitium salt	36
Coenzyme A, Grade II trilitium salt	37
Colipase from porcine pancreas, lyophilizate	81
Colloidal Gold 20 nm suspension	190
Colloidal Gold 40 nm suspension	190
COT Human DNA from human placenta DNA, enriched for repetitive sequences, solution	324
COT Human DNA, CGH Grade from human male placenta DNA, enriched for repetitive sequences, solution	323
Creatinase from microorganism, lyophilizate	82
Creatine Phosphate disodium salt	150
Creatininase from <i>Pseudomonas species</i> , expressed in <i>E.coli</i> , lyophilizate	82
Creatinine Deaminase from <i>Corynebacterium lilium</i> , lyophilizate	83
CTP, Molecular Diagnostic Grade sodium salt, 100 mM	291
D	
D(-)-Lactate monolithium salt	151
D-Biotin-N-hydroxysuccinimide ester crystalline powder	179
D-Biotinoyl-ε-aminocaproic acid-N-hydroxysuccinimide ester powder	180
D-Glucose-6-phosphate:NADP ⁺ 1-oxidoreductase <i>EC 1.1.1.49</i>	
○ Glucose-6-phosphate Dehydrogenase (G6P-DH) from yeast, lyophilizate	
D-Lactate Dehydrogenase (D-LDH) from microorganism, lyophilizate	110
D-Lactate Dehydrogenase (D-LDH), Grade I from <i>Lactobacillus delbrückii</i> , lyophilizate	111
D-Lactate Dehydrogenase (D-LDH), Grade II from <i>Lactobacillus delbrückii</i> , lyophilizate	112
D-Mannitol reduced sodium	32
dATP, PCR Grade sodium salt, 100 mM	274
dCTP, PCR Grade sodium salt, 100 mM	275
ddATP, Sequencing Grade sodium salt, 10 mM	286
ddCTP, Sequencing Grade sodium salt, 10 mM	287
ddGTP, Sequencing Grade sodium salt, 10 mM	288
ddTTP, Sequencing Grade sodium salt, 10 mM	289
Deoxycholate ionic detergent, sodium salt	23
Dephosphorylation Buffer for Alkaline Phosphatase 10x concentrated	331
dGTP, PCR Grade sodium salt, 100 mM	276
Di(adenosine-5'-)penta-phosphate trilitium salt	151
Diethanolamine 85% solution	16
Digoxigenin-11-ddUTP 1 mM solution	319
Digoxigenin-11-dUTP, alkali-labile 1 mM solution	319
Digoxigenin-11-dUTP, alkali-stable 1 mM solution	318
Digoxigenin-11-UTP 10 mM solution	316
Digoxigenin-11-UTP 3.5 mM solution	317
Dilaurylglycerosulfate powder	22

Dithiothreitol, 1,4- ○ 1,4-Dithiothreitol (DTT) crystallizate	
dITP, PCR Grade sodium salt, 100 mM	278
DNA from fish sperm, solution	325
DNase I, recombinant, Grade I from bovine pancreas, expressed in <i>Pichia pastoris</i> , lyophilizate	221
DNase I, recombinant, RNase-free from bovine pancreas, expressed in <i>Pichia pastoris</i> , solution	222
dTTP, PCR Grade sodium salt, 100 mM	279
dUTP, PCR Grade sodium salt, 100 mM	280

E

EagleTaq DNA Polymerase, 5 U/μl from <i>Thermus aquaticus</i> , expressed in <i>E. coli</i> , solution	252
EagleTaq Master Mix	264
EagleTaq Master Mix (Rox)	265
Endoproteinase Asp-N, Sequencing Grade from a mutant of <i>Pseudomonas fragi</i> , lyophilizate	349
Endoproteinase Glu-C from <i>Staphylococcus aureus</i> V8, lyophilizate, salt-free	350
Endoproteinase Glu-C, Sequencing Grade from <i>Staphylococcus aureus</i> V8, lyophilizate	351
Endoproteinase Lys-C from <i>Lysobacter enzymogenes</i> , lyophilizate	352
Endoproteinase Lys-C, Sequencing Grade from <i>Lysobacter enzymogenes</i> , lyophilizate	353
Endoproteinase Lys-C, Sequencing Grade from <i>Lysobacter enzymogenes</i> , solution	354
Esterase ○ Cholesterol Esterase from <i>Candida cylindracea</i> , solution	
Ethylidene-4-NP-G7 Ethyliden-4-nitrophenyl-α-D-maltoheptaosid (EPS), powder	144
Expand High Fidelity PCR Buffer 10x conc., with MgCl ₂	228
Expand High Fidelity PCR Buffer 10x conc., without MgCl ₂	228
Expand High Fidelity PCR System	227
Expand Long Template PCR Buffer 1 10x conc., with 17.5 mM MgCl ₂	230
Expand Long Template PCR Buffer 2 10x conc., with 27.5 mM MgCl ₂	230
Expand Long Template PCR Buffer 3 10x conc., with 27.5 mM MgCl ₂	231
Expand Long Template PCR System	229

F

FAD disodium salt	37
FastStart PCR Buffer 10x conc., with 20 mM MgCl ₂	256
FastStart PCR Buffer 10x conc., without MgCl ₂	256
FastStart PCR Master 2x concentrated	266
FastStart Taq DNA Polymerase, 100 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , solution	255
FastStart Taq DNA Polymerase, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , solution	254
FastStart Taq DNA Polymerase, GMP Grade, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i> , solution	253
Fluorescein-12-dUTP 1 mM solution	321

Fluorescein-12-UTP 10 mM solution	320
Formate Dehydrogenase from yeast, lyophilizate	84
Framework IEP lyophilizate	192
Fructose-1,6-diphosphate trisodium salt	152

G

Galactose 1-Dehydrogenase from <i>E.coli</i> overproducer, lyophilizate	84
Galactose 1-Dehydrogenase from <i>E.coli</i> overproducer, suspension	85
Galactose 1-Dehydrogenase from <i>Pseudomonas fluorescens</i> , suspension	85
GC-RICH Solution 5x concentrated	296
Germall 115 crystalline powder	11
Germall 115, reduced sodium crystalline powder	11
Glucose Oxidase (GOD), chemically modified from <i>Aspergillus niger</i> overproducer, lyophilizate	88
Glucose Oxidase (GOD), Grade I from <i>Aspergillus niger</i> overproducer, lyophilizate	86
Glucose Oxidase (GOD), Grade II from <i>Aspergillus niger</i> overproducer, lyophilizate	87
Glucose-1,6-diphosphate tetra(cyclohexylammonium) salt	152
Glucose-6-phosphate Dehydrogenase (G6P-DH) from <i>Leuconostoc mesenteroides</i> , expressed in <i>E. coli</i> , reduced phosphate, lyophilizate	93
Glucose-6-phosphate Dehydrogenase (G6P-DH) from <i>Leuconostoc mesenteroides</i> , expressed in <i>E. coli</i> , solution	91
Glucose-6-phosphate Dehydrogenase (G6P-DH) from <i>Leuconostoc mesenteroides</i> , lyophilizate	92
Glucose-6-phosphate Dehydrogenase (G6P-DH) from <i>Leuconostoc mesenteroides</i> , suspension	94
Glucose-6-phosphate Dehydrogenase (G6P-DH) from <i>Leuconostoc mesenteroides</i> , expressed in <i>E. coli</i> , lyophilizate	89
Glucose-6-phosphate Dehydrogenase (G6P-DH) from yeast, lyophilizate	94
Glucose-6-phosphate Dehydrogenase (G6P-DH), chemically modified from <i>Leuconostoc mesenteroides</i> , expressed in <i>E. coli</i> , lyophilizate	90
Glucose-6-phosphate disodium salt	153
Glucose-6-phosphate Isomerase ○ Glucose-6-phosphate Isomerase from yeast, suspension	
Glucose-6-phosphate Isomerase from yeast, suspension	95
Glupa-carboxylate monoammonium salt	144
Glutamate Dehydrogenase (NAD(P)) from beef liver, chemically modified, lyophilizate	99
Glutamate Dehydrogenase (NAD(P)) from beef liver, lyophilizate	98
Glutamate Dehydrogenase (NAD(P)) from beef liver, solution	100
Glutamate Dehydrogenase (NAD(P)) from <i>E.coli</i> overproducer, lyophilizate	97
Glutathione, oxidized form (GSSG) lyophilizate	369
Glutathione, reduced form (GSH) crystalline powder	369
Glycerol Kinase (GK) from <i>Bacillus stearothermophilus</i> , lyophilizate	102

Product Index


Glycerol Kinase (GK) from <i>Bacillus stearothermophilus</i> , solution	103	L-Lactate Dehydrogenase (L-LDH) from pig muscle, for use of AST/GOT-Determination according to IFCC recommendations, lyophilizate	113
Glycerol Kinase (GK), concentrated from <i>Bacillus stearothermophilus</i> , solution	103	L-Lactate Dehydrogenase (L-LDH) from pig muscle, for use of AST/GOT-Determination according to IFCC recommendations, solution	114
Glycerol-3-phosphate Dehydrogenase from rabbit muscle, suspension	104	L-Lactate Dehydrogenase (L-LDH) from pig muscle, suspension	115
Glycerol-3-phosphate Oxidase from <i>E.coli</i> overproducer, lyophilizate	105	L-Lactate Dehydrogenase (L-LDH) from rabbit muscle, suspension	117
Glycerol-3-phosphate Oxidase, chemically modified from <i>E.coli</i> overproducer, lyophilizate	106	L-Lactate Dehydrogenase (L-LDH), chemically modified from pig heart, lyophilizate	116
Glycogen, Molecular Biology Grade from mussels, solution	328	Lactate 2-Monooxygenase (Lactate oxidase) from <i>Pedio-coccus</i> species, lyophilizate	119
Glycylglycine crystalline powder	16	Lactate 2-Monooxygenase (Lactate oxidase), Grade I from <i>Aerococcus viridans</i> , expressed in <i>E. coli</i> , lyophilizate	118
GTP disodium salt, powder	370	Lactate 2-Monooxygenase (Lactate oxidase), Grade II from <i>Aerococcus viridans</i> , expressed in <i>E. coli</i> , lyophilizate	119
GTP, Molecular Diagnostic Grade sodium salt, 100 mM	292	Liberase MNP-S lyophilizate, sterile acc. to Ph. Eur.	340
Guanidine Hydrochloride crystals	220	Liberase MTF C/T, GMP Grade, 0.2 µm filtered, lyophilizate	341
Guanosine-5'-O-(2-thiodiphosphate) trilithium salt	51	Liberase T-Flex, Research Grade, 0.2 µm filtered, lyophilizate	343
Guanosine-5'-O-(3-thiodiphosphate) tetralithium salt	51	LightCycler® 480 Multiwell Plate 384 white, 4 barcodes	296
H		LightCycler® 480 Multiwell Plate 96 white, 4 barcodes	297
HAMA Serum, Type 1 lyophilizate	193	Lipase from porcine pancreas, lyophilizate	120
HAMA Serum, Type 2 lyophilizate	194	Lipoprotein Lipase from <i>Pseudomonas</i> species, lyophilizate	120
HawkTaq DNA Polymerase, 5 U/µL from <i>Thermus aquaticus</i> , expressed in <i>E. coli</i> , solution	259	Lipoprotein Lipase, chemically modified from <i>Pseudomonas</i> species, lyophilizate	121
HawkZ05 DNA Polymerase, 40 U/µl from <i>Thermus</i> species Z05, expressed in <i>E. coli</i> , solution	258	Long Range PCR ○ Expand Long Template PCR System	
HawkZ05 Fast One-Step RT-PCR Kit	273	Long Template PCR ○ Expand Long Template PCR System	
Hepes crystalline powder	17	Lysozyme from hen egg white, crystalline powder	122
Hexokinase (HK) from yeast overproducer, lyophilizate	107	M	
Hexokinase (HK) from yeast overproducer, solution	109	M-MLV Reverse Transcriptase, GMP Grade from Moloney Murine Leukemia Virus, expressed in <i>E. coli</i>	268
Hexokinase (HK) from yeast, lyophilizate	110	MAB IgG2b/Fab2a Poly lyophilizate	195
Hexokinase (HK), chemically modified from yeast overproducer, lyophilizate	108	MAB<AFP>M-LJ738 IgG lyophilizate	161
High Fidelity PCR System, Expand ○ Expand High Fidelity PCR System		MAB<AFP>M-TU11 IgG lyophilizate	161
Histone H3 from calf thymus, lyophilizate	336	MAB<CEA>M-TU2 IgG lyophilizate	162
Human Serum frozen solution	216	MAB<CEA>M-TU3 IgG lyophilizate	162
Hybridoma Cloning Supplement ○ BM Condimed H1, 0.2 µm filtered, solution		MAB<CK-MB>M-6.12.47 IgG lyophilizate	163
I		MAB<CK-MB>M-7.4.5 IgG lyophilizate	160
Imidazole crystalline powder	18	MAB<CK-MM>Mix frozen solution	160
Inorganic Pyrophosphatase ○ Pyrophosphatase, inorganic (PPase) from yeast		MAB<DD>M-1.2.57 IgG lyophilizate	163
K		MAB<DD>M-2.1.16 IgG lyophilizate	164
Klenow Enzyme, Labeling Grade from <i>E. coli</i> lysogenic NM 964, solution	308	MAB<Ferr>M-3.170 IgG lyophilizate	165
Kryptofix 221 solution	33	MAB<Ferr>M-4.184 IgG lyophilizate	165
L		MAB<FSH>M-1.303 IgG lyophilizate	166
L(+)-Alanine crystalline powder	147	MAB<FSH>M-W3g IgG lyophilizate	166
L-Lactate Dehydrogenase (L-LDH) from pig heart, suspension	116	MAB<H-S-Amy>M-Tu66C7 IgG lyophilizate	168
		MAB<H-S-Amy>M-Tu88E8 IgG lyophilizate	169


MAB<HCG>M-INN2 IgG lyophilizate	167
MAB<HCG>M-INN22 IgG lyophilizate	167
MAB<IGE>M-323 IgG lyophilizate	169
MAB<IGE>M-7H8 IgG lyophilizate	170
MAB<INSULIN>M-BM1 IgG lyophilizate	170
MAB<INSULIN>M-ST3 IgG lyophilizate	171
MAB<LH>M-11412 IgG lyophilizate	171
MAB<LH>M-2.406 IgG lyophilizate	172
MAB<PRL>M-C4E4 IGG lyophilizate	172
MAB<PRL>M-H12G10 IgG lyophilizate	173
MAB<TSH>M-A8 IgG lyophilizate	174
MAB<TSH>M-TU1.20 IgG lyophilizate	174
MAB33-IgG1 lyophilizate	196
MAB33-IgG1/Fab1 Poly lyophilizate	197
MAB33-IgG1/IgG1 Poly frozen solution	198
Magnesium Chloride Solution	
○ MgCl ₂ Stock Solution 25 mM	
Magnetic Particles, Streptavidin	
○ Streptavidin Magnetic Particles suspension	
Malate Dehydrogenase from pig heart, lyophilizate	123
Malate Dehydrogenase, chemically modified from pig heart, lyophilizate	123
Malate Dehydrogenase, IFCC Quality from pig heart, lyophilizate	124
Malate Dehydrogenase, IFCC Quality from pig heart, solution	125
MES crystallizate	18
MgCl₂ Stock Solution 25 mM	297
Micr-O-protect solution	12
Microprotect	
○ Micr-O-protect solution	
Mn(OAc)₂ Stock Solution 25 mM	239
Multi Analyte Stripe universal device	183
MycotoOL PCR Mycoplasma Detection Kit	359
N	
N-Acetyl-L-Cysteine crystallizate	153
N-Carbamoylsarcosine Amidase from <i>E.coli</i> overproducer, lyophilizate	70
n-Dodecyl-β-D-maltoside nonionic detergent, powder	25
N-Glycosidase A (PNGase A) solution	345
N-Methylhydantoinase (ATP-hydrolyzing) from <i>Arthrobacter</i> species, expressed in <i>E. coli</i>, lyophilizate	126
N-Methylisothiazolone (MIT) crystalline powder	13
N-Methylisothiazolone (MIT) for potassium test crystalline powder	14
n-Octyl β-D-glucoside nonionic detergent, powder	25
NAD(P)H Dehydrogenase (quinone) (Diaphorase) from pig heart, suspension	126
NAD, Grade I free acid	38
NAD, Grade II free acid	38


NADH, Food Grade disodium salt, lyophilizate	366
NADH, Grade I disodium salt	39
NADH, Grade II disodium salt	40
NADH, Grade II for potassium test, disodium salt	41
NADP disodium salt	42
NADP monopotassium salt	42
NADPH tetrasodium salt	43
Neuraminidase (exo-α-sialidase) from <i>Vibrio cholerae</i>, solution	344
Nitrate Reductase from <i>Aspergillus</i> species, lyophilizate	127
Non-specific Protease	
○ Pronase from <i>Streptomyces griseus</i> , lyophilizate	
Nonidet P40 nonionic detergent, aqueous solution	26
↔ Nonradioactive Labels	
↑ Biotin-16-ddUTP 1 mM solution	
↑ Biotin-16-dUTP 1 mM solution	
↑ Digoxigenin-11-ddUTP 1 mM solution	
↑ Digoxigenin-11-dUTP, alkali-labile 1 mM solution	
↑ Digoxigenin-11-dUTP, alkali-stable 1 mM solution	
↑ Fluorescein-12-dUTP 1 mM solution	
↑ Tetramethylrhodamine-5-dUTP 1 mM solution, lithium salt	
↔ Nonradioactive RNA Labels	
↑ Biotin-16-UTP 10 mM solution	
↑ Digoxigenin-11-UTP 10 mM solution	
↑ Fluorescein-12-UTP 10 mM solution	
NucleoMix (dTTP), PCR Grade sodium salt, 100 mM (25 mmol/l each dNTP)	283
NucleoMix (dTTP), PCR Grade sodium salt, 40 mM (10 mmol/l each dNTP)	281
NucleoMix (dUTP), PCR Grade sodium salt, 100 mM (25 mmol/l each dNTP)	284
NucleoMix (dUTP), PCR Grade sodium salt, 40 mM (10 mmol/l each dNTP)	282
O	
Oxalate Oxidase from barley seedings, lyophilizate	128
P	
PAB<CRP>S IgG frozen solution	176
PAB<T3>S IgG lyophilizate	176
PAB<T4>S IgG lyophilizate	177
pancreatic trypsin inhibitor	
○ Aprotinin from bovine lung, lyophilizate	
Papain from <i>Carica papaya</i>, 0.2 μm filtered, solution	356
Papain from <i>Carica papaya</i>, suspension	355
PCR Buffer 10x conc., with 15 mM MgCl₂	237
PCR Buffer 10x conc., without MgCl₂	237
Pefabloc SC (AEBSF) 4-(2-Aminoethyl)-benzenesulfonyl fluoride hydrochloride, powder	30
Pepstatin from <i>Streptomyces</i> species, lyophilizate	371
Peroxidase (POD), EIA Grade from horseradish, lyophilizate	210
Peroxidase (POD), Grade I from horseradish, lyophilizate	210
Peroxidase (POD), Grade II from horse radish, lyophilizate	128


~ Peroxidase Substrates	
↑ 2,2'-azino-bis(3-ethylbenzthiazoline-6-sulphonic acid (ABTS) solution	
↑ ABTS 2,2'-Azino-di-[3-ethylbenzthiazoline sulfonate (6)] diammonium salt	
↑ BM Blue POD Substrate, soluble 3,3'-5,5'-Tetramethylbenzidine (TMB), solution	
Phosphatase, acid	
○ Acid Phosphatase from potato, lyophilizate	
Phosphoenolpyruvate (PEP) monosodium salt	155
Phosphoenolpyruvate (PEP) tri(cyclohexylammonium) salt	154
Phosphoenolpyruvate (PEP), for potassium test tri(cyclohexylammonium) salt	154
Phosphogluconate Dehydrogenase (decarboxylating) from yeast, lyophilizate	129
Pipes disodium salt	19
Pipes free acid	19
Polidocanol (Thesit)	26
Poly [d(A-T)] powder	325
Poly BSA Type I frozen solution	201
Poly BSA Type II frozen solution	201
Poly Peroxidase (Poly POD), EIA Grade from horseradish, lyophilizate	211
Poly(A) potassium salt, lyophilizate	327
Poly(A) potassium salt, solution	326
Pronase from <i>Streptomyces griseus</i>, lyophilizate	357
Protease V8	
○ Endoproteinase Glu-C from <i>Staphylococcus aureus</i> V8, lyophilizate, salt-free	
Protector RNase Inhibitor from rat lung, expressed in <i>E. coli</i>	298
Proteinase K Inhibitor	
○ Pefabloc SC (AEBSF) 4-(2-Aminoethyl)-benzenesulfonyl fluoride hydrochloride, powder	
Proteinase K, recombinant, PCR Grade from <i>Titrirachium album</i>, expressed in <i>Pichia pastoris</i>, lyophilizate	224
Proteinase K, recombinant, PCR Grade from <i>Titrirachium album</i>, expressed in <i>Pichia pastoris</i>, solution	223
Pyrophosphatase, inorganic (PPase) from yeast	333
Pyruvate Kinase from <i>Bacillus stearothermophilus</i>, lyophilizate	132
Pyruvate Kinase from rabbit muscle, suspension	133
Pyruvate monosodium salt	155
Pyruvate Oxidase from <i>E.coli</i> overproducer, lyophilizate	134
R	
Rabbit IgG (PAB<->K-IgG) lyophilizate	202
Random Octamer Primer 2.5x concentrated solution	300
Reverse Transcriptase, Transcriptor	
○ Transcriptor Reverse Transcriptase recombinant, expressed in <i>E. coli</i>	
RMS Z05 DNA Polymerase, 200 U/μL from <i>Thermus species</i> Z05, expressed in <i>E. coli</i>, solution	240
RNA from yeast, powder, free acid	326
RNase A from bovine pancreas, lyophilizate, powder	225
RNase Inhibitor, Protector	
○ Protector RNase Inhibitor from rat lung, expressed in <i>E. coli</i>	
RNase Inhibitor, recombinant, GMP Grade from rat lung, expressed in <i>E. coli</i>	299
S	
S-Butyrylthiocholine Iodide crystallizate	156
Sarcosine Oxidase from <i>E.coli</i> overproducer, lyophilizate	135
Sheep IgG (PAB<->S-IgG) lyophilizate	202
Streptavidin Magnetic Particles suspension	183
Streptavidin R-Phycoerythrin LumiGrade Reagent Ready to use solution	181
Streptavidin R-Phycoerythrin LumiGrade Ultrasensitive Reagent solution	182
Streptavidin, recombinant from <i>Streptomyces avidinii</i>, expressed in <i>E. coli</i>, lyophilizate	178
StreptaWell, 384 plate transparent, coated with recombinant streptavidin	184
StreptaWell, 384 plate white, coated with recombinant streptavidin	184
StreptaWell, C1, breakapart transparent, coated with recombinant streptavidin	185
StreptaWell, C8 module, high binding capacity white, coated with recombinant streptavidin	185
StreptaWell, C8, breakapart, high binding capacity transparent, coated with recombinant streptavidin	186
StreptaWell, C8, lockwell transparent, coated with recombinant streptavidin	186
StreptaWell, C96 plate white, coated with recombinant streptavidin	187
StreptaWell, C96 plate, high binding capacity transparent, coated with recombinant streptavidin	187
StreptaWell, C96 plate, high binding capacity white, coated with recombinant streptavidin	188
StreptaWell, F8 module transparent, coated with recombinant streptavidin	189
StreptaWell, F8 module, high binding capacity transparent, coated with recombinant streptavidin	189
T	
T4 DNA Ligase recombinant form of the enzyme from T4 phage, solution	332
T4 DNA Polymerase from T4 plasmid pTL43W infected <i>E. coli</i> 71-18, solution	232
T4 DNA Polymerase Incubation Buffer 5x concentrated	233
T4 Gene 32 Protein, recombinant recombinant from T4 phage, expressed in <i>E. coli</i>, solution	300
Taq DNA Polymerase, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i>, solution	235
Taq DNA Polymerase, 50 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i>, glycerol-free solution	236
Taq DNA Polymerase, GMP Grade, 5 U/μl from <i>Thermus aquaticus</i> BM, expressed in <i>E. coli</i>, solution	234
Taurodesoxycholol sodium salt	23
Terminal Transferase Reaction Buffer 5x concentrated	310
Terminal Transferase, recombinant from calf thymus, expressed in <i>E. coli</i>, solution	309


Tetramethylrhodamine-5-dUTP 1 mM solution, lithium salt	322
Thio-NAD free acid	43
Thrombin from human plasma, lyophilizate	136
TMB ○ 3,3',5,5'-Tetramethylbenzidine (TMB) crystalline powder	
TOOS (N-Ethyl-N-(2-hydroxy-3-sulfopropyl)-m-toluidine, monosodium salt, dihydrate)	145
Transcriptor Reverse Transcriptase recombinant, expressed in <i>E. coli</i>	269
Transcriptor RT Buffer 5x concentrated	270
Tribromo-hydroxybenzoic acid crystallizate	146
Triose-phosphate Isomerase from rabbit muscle, suspension	136
Tris crystallizate	20
Triton X-100 nonionic detergent, viscous liquid	27
Trypsin Inhibitor, Pancreas Type (BPTI) ○ Aprotinin from bovine lung, lyophilizate	
Trypsin Substrate ○ Chromozym TRY powder	
Trypsin, recombinant from porcine pancreas, expressed in <i>Pichia pastoris</i>	358
trypsin-kallikrein inhibitor ○ Aprotinin from bovine lung, lyophilizate	
Tth DNA Polymerase from <i>Thermus</i> species, expressed in <i>E. coli</i> , solution	238
Tth DNA Polymerase PCR Buffer 10x conc., with 15 mM MgCl ₂	239
Tth DNA Polymerase RT-PCR Buffer 5x concentrated	240
Tween 20 purified, solution	27
U	
UDP-Galactose disodium salt, powder	363
UDP-Glucose disodium salt, powder	364
UDP-N-Acetylglucosamine disodium salt	362
Uracil-DNA Glycosylase, heat-labile from marine bacterium BMTU 3346, expressed in <i>E. coli</i>	301
Urease from jack bean, lyophilizate	137
Uricase from <i>Arthrobacter protophormiae</i> , lyophilizate	138
UTP, Molecular Diagnostic Grade sodium salt, 100 mM	293
V	
V8 Protease ○ Endoproteinase Glu-C from <i>Staphylococcus aureus</i> V8, lyophilizate, salt-free	
Valinomycin crystallizate	33
W	
Water, PCR Grade	303, 304
Z	
Zwittergent 3-14 zwitterionic detergent, powder	29
Y	
γ-Glutamyltransferase from hog kidney, lyophilizate	101















Published by

Roche Diagnostics GmbH
Sandhofer Straße 116
68305 Mannheim
Germany

© 2011 Roche Diagnostics.
All rights reserved.

custombiotech.roche.com

03064395990  0311

